

Sensation OPEN

CT

Software

System

Installation of SOMARIS/5 VB39

© Siemens 2007

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

Print No.: CT02-023.816.13.01.02
Replaces: n.a.

English
Doc. Gen. Date: 11.07

Document revision level

The document corresponds to the version/revision level effective at the time of system delivery. Revisions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order current revision levels.

Disclaimer

The installation and service of equipment described herein is to be performed by qualified personnel who are employed by Siemens or one of its affiliates or who are otherwise authorized by Siemens or one of its affiliates to provide such services.

Assemblers and other persons who are not employed by or otherwise directly affiliated with or authorized by Siemens or one of its affiliates are directed to contact one of the local offices of Siemens or one of its affiliates before attempting installation or service procedures.

1	General	5
	Definitions and Abbreviations	5
	Software not approved by Siemens	7
	Overview	9
	Software kit VB39A for ICS, IRS and IES.	10
	Images and Raw Data	11
	Required system data.	12
	Service key	12
	Host name.	12
	Regional settings	12
	Report.htm and license.dat.	13
	Backup information.	14
	Creating backup files	14
2	SW Installation Navigator	15
	Preparation	15
	Required data media	16
	Starting the Installation.	17
	System type and serial numbers.	17
	Starting the Installation.	18
	Installation is running	20
	Configure the keyboard language and enable Windows XP Access	20
	Restoring Steps	22
	Restore 1	22
	Restore 2	22
	Completion of configuration	25
	Import Reference Images	25
	Regional settings	25
	Setting the time.	25
	Flexible Cupping Correction	25
	Set Auto Filming, Auto Viewing, or Auto 3D in the Wizard	25
	Siemens Remote Services	27
3	SW Installation of IRS	28
	Installation of IRS	28
	Additional steps after installing the TCR.	29
4	SW Installation Wizard	30
	Preparation	30
	Required data media	31
	Starting the Installation.	32
	Starting the Installation.	33

Installation running	35
Configure the keyboard language and enable Windows XP Access.	35
Restore Steps	37
Restore 1	37
Restore 2	37
Completion of configuration	39
Regional settings	39
Setting the time	39
Restart	39
Siemens Remote Services	40
5 Additional Worksteps	41
Software self test.	41
Map drive C: of ICS (Navigator) at Wizard	43
Backup at the Navigator	45
Disable Windows XP Access.	46
Backup at the Wizard	47
Disable Windows XP Access.	48
6 Changes to previous version	49
7 Index	50

Definitions and Abbreviations

Abbreviations	Definitions
ACR-NEMA	American College of Radiology - National Electrical Manufacturers Association
BIOS	Basic Input Output System, lowest PC configuration settings
AET	Application Entity Title - Name of DICOM service participating in the communication either as a user (SCU) or provider (SCP)
CD	Compact Disk
CDR	Compact Disc Recordable
CT	Computed Tomography
DICOM	Digital Image Communication in Medicine (ACR-NEMA standard)
DHCP	Dynamic Host Configuration Protocol - procedure for assigning numeric IP address to nodes in which blocks of addresses are leased to clients on a first-come, first-serve basis. Typically, DHCP is combined with WINS. This means that an IP address is both dynamically assigned and automatically associated with a domain name. This can be combined with DNS, resulting in a dynamically assigned numeric IP address that is automatically associated with both NetBIOS and Internet names.
DNS	Domain Name Services - address resolution procedure that associates numeric IP addresses with Internet-style names
FTP	File Transfer Protocol - a client-server application that transfers files from one computer to another
Gateway	A display that controls traffic between two networks using factors above the network layer
ICS = Navigator	Image Control System - control computer for the entire CT scanner, with interface to the acquisition system
IES = Wizard	Imaging Evaluation System
IMS	Imaging System, consisting of ICS, IRS, and console
IRS	Image Reconstruction System - reconstruction computer
LAN	Local Area Network (here, always Ethernet)
LMHOSTS	Address resolution procedure which statically associates numeric IP addresses with NetBIOS names. The LMHOSTS file is checked to see if the name is being addressed in this file.
MEP	Market Entrance Phase
MOD	Magneto-Optical Disk

Abbreviations	Definitions
NetBIOS	Network Basic Input Output System - a collection of network services and programming interfaces. NetBIOS services run on network protocols such as TCP/IP.
NIC	Network Interface Card - an expansion card that provides a physical connection to the network
NT	"New Technology" Operating System - trademark of Microsoft, Inc.
OI	Operating Instructions
PB	Patient Browser
PPP	Point-to-Point Protocol - a data link and network protocol designed for use over serial connections (such as dial-up modems)
RAS	Remote Access Services - a Windows NT service that supports dial-in and dial-out connections using modems and ISDN.
Router	A multiport computer or piece of dedicated hardware that examines network layer addresses of packets on a network and retransmits them on a different network, if necessary
SCU	Service Class User
SCP	Service Class Provider
TCP/IP	Transport Control Protocol/Internet Protocol - the most widely used standard for networking
UI	User Interface
USC	Uptime Service Center
WINS	Windows Internet Name Service - address resolution procedure that dynamically associates numeric IP addresses with NetBIOS names. It is possible to combine DNS and WINS.
XP	"eXPerience" Operating System - trademark of Microsoft, Inc.

Software not approved by Siemens

⚠ WARNING

Unauthorized or incorrect operation or changes to software can cause malfunctions in your system. This can result in injuries or system defects.

This can result in injuries or system defects.

- ⇒ For this reason, the use or installation of software not approved by Siemens is strictly prohibited.

Hazards and Warnings

⚠ CAUTION

[hz_serdoc_F13G14U01M01 (p30)]

Connecting non-validated equipment!

General incorrect treatment and/or physical injury.

- ⇒ Persons who connect additional equipment to the medical device are considered to be system configurators and are therefore responsible for ensuring compliance with the valid version of relevant standards (e.g. system standard IEC/EN 60601-1-1 and/or other applicable standards). If you have any questions, please contact Siemens Customer Service.

⚠ CAUTION

[hz_serdoc_F13G14U02M01 (p30)]

Unauthorized system modifications!

Incorrect patient treatment

- ⇒ Unauthorized modifications on the system hardware, the software, or the software configuration can lead to serious malfunctions of the system, to data corruption, or to incorrect patient treatment. Unauthorized installation and execution of software other than that specified for use on the scanner with VSIM may jeopardize patient safety. All software installations have to be performed by authorized personnel.

⚠ CAUTION

[hz_serdoc_F09G07U01M01 (p30)]

Switching off the computer without shutting down the system!

Data loss or system corruption is possible.

- ⇒ To turn off or shut down the computer properly, use Options > End Session > Shutdown System in the syngo user interface.

 CAUTION

[hz_serdoc_F09G08U01M01 (p30)]

Mismatch in DICOM conformance!

Incorrect patient treatment

⇒ **Confirm customer receipt of the validation summary as part of the ATP.**

Overview

ICS + IRS = syngo CT Acquisition Workplace

In the standard configuration, the ICS is the main console for interchangeable scanning and imaging functions. The IRS is equipped with four computers as well as a switchable keyboard, mouse, and monitor.

IES = syngo CT Workplace

The IES console may be installed in parallel to the main operator ICS console via a shared database link. The IES is used as a second console for post-processing applications.

NOTE

The IES can run only when the ICS is switched on; if the ICS is switched off, the shared database cannot be accessed.

Software kit VB39A for ICS, IRS and IES

The software installation kit for SOMATOM Sensation Open systems is included in the delivery volume and comprises the following items required for a complete SW installation:

- OEMsetup_1 english VB39A CD 1 of 5
- OEMsetup_2 english VB39A CD 2 of 5
- syngo CT 2008O Som/5_ICS/IES_SW/FW VB39A CD 3 of 5
- DICOM_SR_Viewer CD 4 of 5
- syngo CT 2008O Som/5_ICS/IES_SW/FW VB39A CD 5 of 5
- SOM/5_IRS_VB39A_V1.0_32bit

Images and Raw Data

NOTE

For safety reasons, the customer needs to archive all images and raw data prior to installing the SW, if possible.

Required system data

For installing the software, the following data should always be available:

Service key

- From the syngo user interface, select **Options -> Local Service**.
- Write down the service key

ICS	Service key (level=7)		
IES 1	Service key (level=7)		
IES 2	Service key (level=7)		

----- Service key = 20 characters ----- 14 characters ----- last 6 characters

Host name

- In the main menu of the service software, select **Configuration**.
- Click Next.
- Select **Local Host -> TCP/IP LAN**.

Write down the host name.

ICS	Computer name	
IES 1	Computer name	
IES 2	Computer name	

Regional settings

- From the syngo user interface, select -> **Options -> Configuration**.
- In the configuration panel, double-click the **Regional and language options** icon.

Write down the following data:

Regional Options	Language		
Customized	Time		Time separator:
	Date		Date separator:
Languages ->Details	Input language		

Report.htm and license.dat

In case a backup file becomes corrupted, it may be useful to have the system configuration and the license file available. For this reason, you have to save these files.

- Open **Windows Explorer**
- Copy the file **report.htm** in folder **C:\Somaris\service\html\report** to, e.g., a floppy
- Copy the file **license.dat** in folder **C:\Somaris\Config\Licensing** to, e.g., a floppy.

Backup information

Make sure that up-to-date backups of the system are always available.

The CD-writer (drive [-R-] CD-R) is the standard backup device for the Navigator and the Wizard. Alternatively, the MOD may be used, if available.

Make sure a separate backup medium (e.g. "CD Recordable") is on hand for both the Navigator and the Wizard (if available).

Creating backup files

If possible, you should generate new backup files before installing the software.

- Start the service software and select the **Backup & Restore** function in the main menu.
- Insert the backup CD-R.
- From Command, select **Backup**.
- From Drives, select the CD drive, e.g., [-R-] CD-R.
- From Packages, select **the corresponding name** and click **Go**. Perform the backup of...
 - SW-settings
 - CT-settings
 - VSIM- settings
 - Inspace settings (at the IES)
 - Security settings

⇒ **BackupRestore Running ... Wait until Ready is displayed.**
- When "Ready" is displayed, **remove the backup medium**.

Preparation

NOTE

If you don't have a current backup, create one now, if possible. Refer to [\(Creating backup files / p. 14\)](#).

Required data media

To perform the SW installation of the Navigator (ICS), the following data media are required:

OEMsetup_1 english VB39A CD 1 of 5	
OEMsetup_2 english VB39A CD 2 of 5	
syngo CT 2008O Som/5_ICS/IES_SW/FW VB39A CD 3 of 5	
DICOM_SR_Viewer CD 4 of 5	
syngo CT 2008O Som/5_ICS/IES_SW/FW VB39A CD 5 of 5	
Backup CD-R	

Starting the Installation

NOTE	During the software installation of the ICS, the IES has to be switched off.
NOTE	Unplug the network cable between the IES and the Ethernet switch.
NOTE	Switch off the gantry at the control box.

System type and serial numbers

Name	Identifier for backup files	Material No.	Serial no.	
SOMATOM Sensation Open	P30L	8872017	49xxx	

Starting the Installation

- From the **Wizard**, shut down the application and switch off the system.
- Insert the **OEMsetup_1 english VB39A CD 1 of 5** in the Navigator.
- From the **Navigator**, select **System/End** and **Shut down the system**.
- Switch off the system at the **control box**.

Disconnect the network cable of the ICS from the Ethernet switch.

- **Switch on the system at the control box** (computer ON).

⇒ The computer boots from the CD-ROM. After a short while, the following window appears:

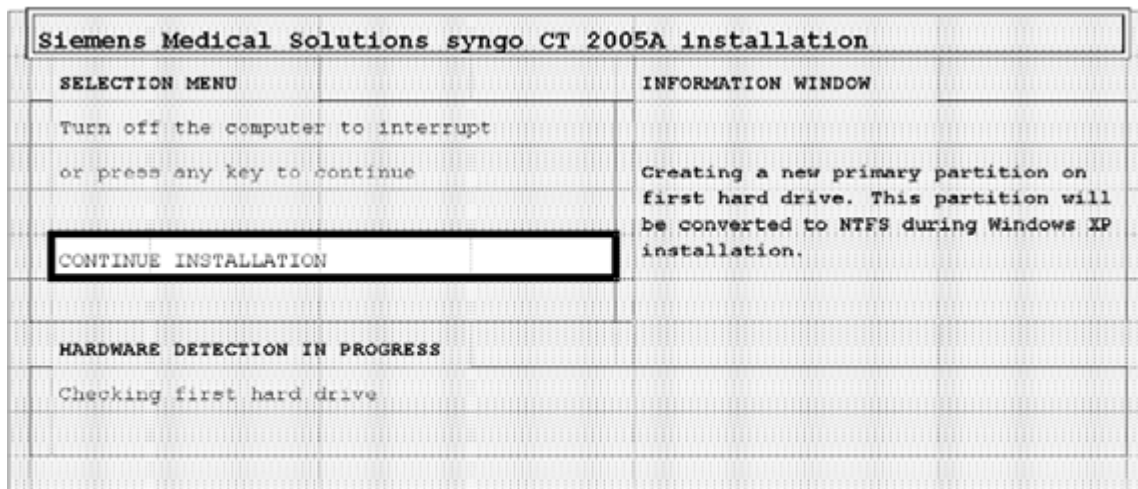


Fig. 1: Starting the installation screen

- Press **<Enter>** to continue installation (the ICS is beeping). This reboots the ICS.
- Select the distribution type: **CUSTOMER** (default).
- Enter the **machine name** (host name) as noted by replacing the default entry "MedPC" and pressing Enter.

NOTE

Please note that the computer name must not contain blanks or underscores. Use the following naming convention for the ICS computer name: **ct<serial number>** (e.g. **ct123456**).

- Select **Keep data on second harddrive** or select **Clear data on second harddrive** (deletes patient database).
Use ENTER to change the selection.

NOTE

Do not use 'Clear all harddrives 1-2(3)'. This option will format all hard drives only (removes all software). Reinstallation of the SW will not be started.

- Use **CONTINUE INSTALLATION** to confirm the start of installation (ICS is beeping if "Clear data" is selected).
- Press any key.
 - ⇒ All necessary files are copied from OEMsetup_1 english VB39A CD 1 of 5 to the hard drive.
- When the computer is beeping, take out OEMsetup_1 english VB39A CD 1 of 5 and insert **OEMsetup_2 english VB39A CD 2 of 5**, then **press any key**
- Insert the **syngo CT 2008O Som/5_ICS/IES_SW/FW VB39A CD 3 of 5**.
 - ⇒ The installation continues automatically.
 - ⇒ After approx. 40 minutes, the CD is ejected.
- Insert the **DICOM_SR_Viewer CD 4 of 5**.
 - ⇒ The installation continues automatically.
 - ⇒ After approx. 3 minutes, the CD is ejected.
- Insert the **syngo CT 2008O Som/5_ICS/IES_SW/FW VB39A CD 5 of 5**.
 - ⇒ The installation continues automatically.
 - ⇒ After approx. 30 minutes, the CD is ejected

NOTE

If you miss inserting the requested CD in time, a warning is displayed:

No package has been found. Please insert the CD "....." in the CD drive to continue the installation.

Insert the requested CD and press any key to continue.

Installation is running

- ⇒ The system reboots automatically after finishing the copy job.
- The following screen appears:

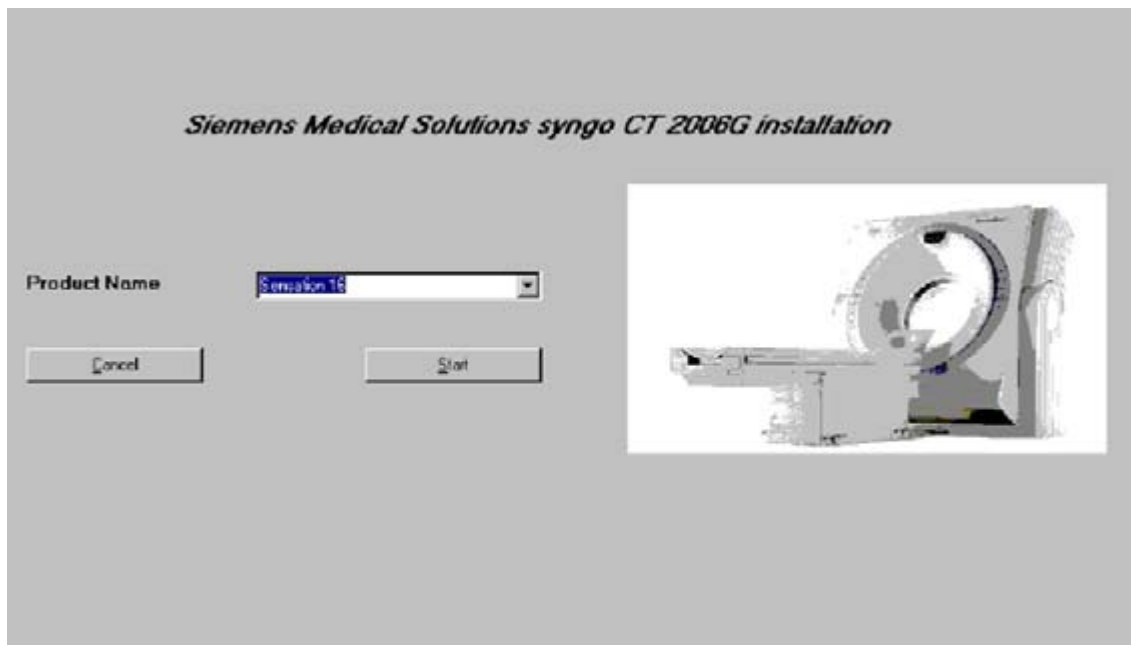


Fig. 2: Starting the installation

- Select the **system to be installed**.

Use the drop-down menu to change the product name.

- Click Start

Completion of the ICS installation (approx. 15 minutes) is indicated by a pop-up window with the message **Please enter license**.

- Click '**Click here to start the Local Service Configuration**' to start the service software (complete service license required).

Configure the keyboard language and enable Windows XP Access

Perform the following steps to configure keyboard languages **other than English**. Select **Start -> Settings -> Control Panel**

- Search for and click **Date, Time, Language and Regional and Options**
- Click **Regional and Language Options**.
- Select the **Languages -> Details** tab card
- Click **Add** and select the correct **Input Language** and confirm with **OK**.
- Select the English keyboard language, click **Remove** and confirm with **OK**.

⇒ A message window appears informing you that removal is only possible after reboot. Confirm with **OK**

- Click **OK** in the remaining windows or close the remaining windows.

Perform the following steps to **enable Windows XP** access at the time of the update.

- Enter the Service Key and open the Local Service main menu.
- From the service software home menu, select **Utilities**.
- Under Source, select **Escape to OS** and at Command **NT Command Interpreter**.
- Under Parameters, type **c:\Somaris\Tools\bin\KeyFilter.bat disable**.
- Activate with **Return** or **GO**.
 - ⇒ A reboot is required to make the change effective. This reboot will be performed after Restore 1.
- Click **Home**.

Restoring Steps

Restore 1

- In the main menu of the service software, select the **Backup & Restore** function.
- Insert the **Navigator** backup CD.
- From Command, select **Restore**.
- From Drives, select the **corresponding drive**.
- From Archive, select **<hostname>_SW-Settings02-PXXX_date_time_ok.ar**.

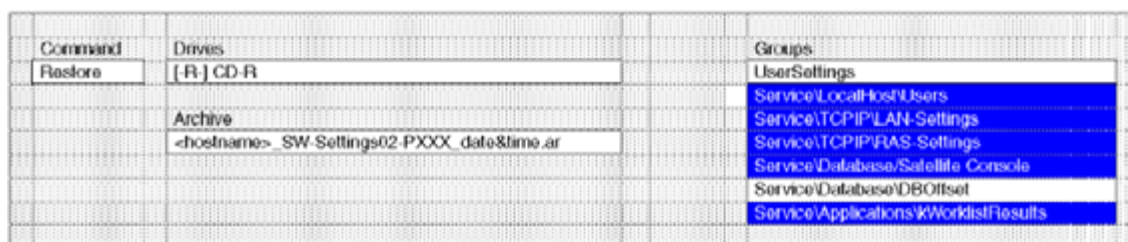


Fig. 3: Restoring 1 SW settings

- From Groups, select **all** except **User Settings** and **Service\Database\DBOffset** and click **Go**.
 - ⇒ **BackupRestore Running ...**; wait until **Ready** is displayed.
- From Archive, select **<hostname>_CT-Settings-PXXX_date_time_ok.ar**.



Fig. 4: Restoring 1 CT settings

- From Groups, select **all** except **IRS** and click **Go**.
 - ⇒ **BackupRestore Running ...**; wait until **Ready** is displayed.
- From Archive, select **<hostname>_SecuritySettings-PXXX_date_time_ok.ar**.
- From Groups, select **all** and click **Go**.
- Wait until "Ready" is displayed.
- Click **Home**.
- Connect the network cable again.

Restore 2

- In the main menu of the service software, select the **Backup & Restore** function.
- From Command, select **Restore**.
- From Drives, select the **corresponding drive**.

- From Archive, select **<hostname>_SW-Settings02-PXXX_date_time_ok.ar**.

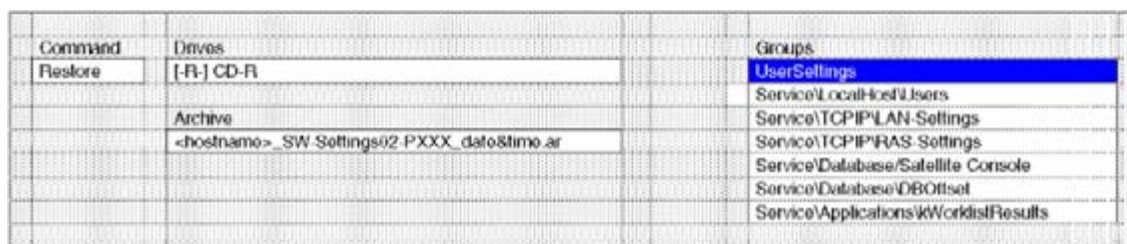


Fig. 5: Restoring 2 SW settings

NOTE

During the restore of SW settings, the following red error message can be ignored

```
cmd /C C:\Somaris\Config\BaReUp\CapPb\update.bat
```

- From Groups, select **User Settings** and click **Go**.
⇒ **BackupRestore Running ...**; wait until **Ready** is displayed.
- From Archive, select **<hostname>_SecuritySettings-PXXX_date_time_ok.ar**.
- From Groups, select **All** and click **Go**.
⇒ **BackupRestore Running ...**; wait until **Ready** is displayed.
- From Archive, select **<hostname>_CT-SettingsPXXX_date&time**.



Fig. 6: Restoring 2 CT settings

- From Groups, select **IRS** and click **Go**.

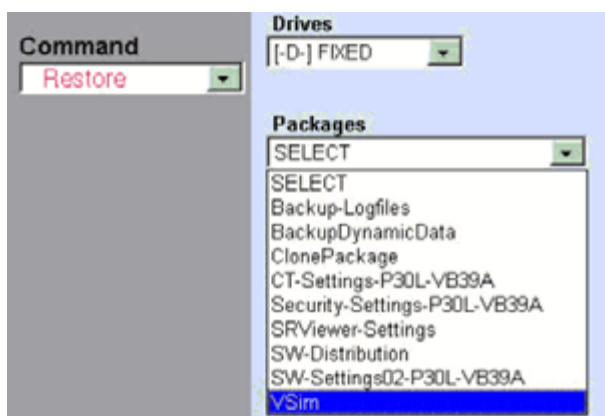


Fig. 7: Backup VSim settings

- From Groups, select **VSIM** and click **Go**.
⇒ **BackupRestore Running ...**; wait until **Ready** is displayed.

- Take out the backup CD.
- Click **Home**.

Completion of configuration

Import Reference Images

This is mandatory if "Clear data on second hard drive" was selected in the installation screen.

- From the patient browser, select **Transfer -> Import from Offline**.
- From the drop-down list, select the path **H:\SiteData\ReferenceImages**.
- Select all images and click **OK**.
 - ⇒ After the import, the images are available in the patient browser.

Regional settings

- From the syngo user interface, select -> **Options -> Configuration**.
- In the configuration panel, double-click the **Regional and language options** icon.
- In the Regional Options tab card, select a preferred country and click **Customize**.
- Change the date, time, and numbers as noted in ([Starting the Installation / p. 17](#)).
- In the **Languages** tab card, click **Details** and select the default input language type as noted in ([Starting the Installation / p. 17](#)).

Setting the time

- In the main menu of the service software, select the **Configuration** function.
- Select **Site Info** and click >. Check/correct the **Time zone and the time**
Note: do not activate the NTP time server or select automatic Daylight Savings Time
- If changes were made, click **Save**. Confirm the message displayed with OK.

Flexible Cupping Correction

- From the service software home menu, select **Utilities**.
- Under Source, select **Escape to OS** and at the Command select NT **Command Interpreter**.
- Under Parameters, type **c:\Somaris\tools\bin\configure_CupCor.exe**
- Activate with **Return** or **GO**.
- In the window that is displayed, keep the selected value and click **Start**

Set Auto Filming, Auto Viewing, or Auto 3D in the Wizard

If the customer requires Auto Filming in the first or second Wizard, he or she can select this under Options/Config/Examination. This checkbox, however, is only available after performing the following work step (determinewizard.pl) at the Navigator.

- Select **Programs / Accessories ->CommandPrompt**.
- In the open command window, type in **perl c:\Somaris\utils\determinewizard.pl**
- Activate with **Return**.
 - ⇒ This script will detect all Wizards that are active and connected to the shared database and make them available to the customer in the Configuration menu.
- Press **Start**. Select **Shut down... -> Restart**.

Siemens Remote Services

System Management installation
• Open Local Service Home menu and go to Configuration. • At System Mgmt-> Master check and correct the IP addresses. • Click Save and then Next • At System Mgmt-> Agent Control, see whether the Installed checkbox is set. • Click Go • At System Mgmt-> Agent Control set the Activated checkbox. ⇒ This may take a few minutes • Click Go • Click Finish
Set Remote Service Access Level
• Go to Options -> Service -> Remote Service • Set the Remote Service Access Control Mode to Limited access permanent

Installation of IRS

After replacing one of the 4 computers in the IRS, the software has to be installed at the computer that was replaced.

To install the software:

- Insert the **CD SOM/5_IRS_VB39A** into the CD drive of the replaced computer
- In the command window, select IRS -> Reboot
 - ⇒ The computer boots from CD and starts the SW installation.
- When the installation is finished, take out the installation CD and reboot all 4 computers

Additional steps after installing the TCR

NOTICE

During the TCR software installation, the IRS "default tables" are deleted.

The IRS "default tables" have to be manually transferred to the IRS.

⇒ As a result, the batch file `lrsTablesInstall.bat` has to be executed. Afterward, the "site-specific tune-up tables" have to be restored.

- From the service software home menu, select **Utilities**.
- Under Source, select **Escape to OS** and at Command, **NT Command Interpreter**.
- Under Parameters, type `c:\Somaris\utils\lrsTablesInstall.bat`
- Activate with **Return** or **GO**.
 - ⇒ The "**default tables**" are transferred to the IRS.
- In the main menu of the service software, select the **Backup & Restore** function.
- Insert the **Navigator** backup CD.
- From Command, select **Restore**.
- From Drives, select the **corresponding drive**.
- From Archive, select **CT-Settings-P30XXX-<hostname>**.

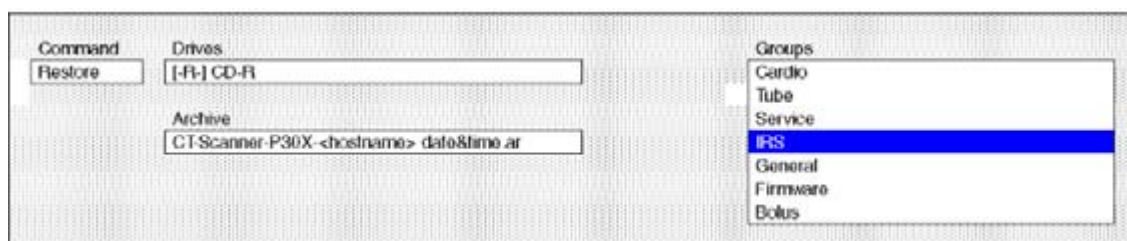


Fig. 8: Restoring CT settings after TCR exchange
(This is only an example!)

- From Groups, select **IRS** and click **Go**.
 - ⇒ The "**site-specific tune-up tables**" are transferred to the IRS.
- When "Ready" is displayed, remove the backup CD **Navigator** and restart the complete system.

Preparation

NOTE

If you don't have a current backup, create one now, if possible. Refer to [\(Creating backup files / p. 14\)](#).

Required data media

To perform the SW installation of the Navigator (ICS), the following data media are required:

OEMsetup_1 english VB39A CD 1 of 5	
OEMsetup_2 english VB39A CD 2 of 5	
syngo CT 2008O Som/5_ICS/IES_SW/FW VB39A CD 3 of 5	
DICOM_SR_Viewer CD 4 of 5	
syngo CT 2008O Som/5_ICS/IES_SW/FW VB39A CD 5 of 5	
Backup CD-R	

Starting the Installation

NOTE

During the software installation of the IES, the ICS has to be switched on.

NOTE

Unplug the network cable between the IES and the Ethernet switch.

Starting the Installation

- Click **Home** and terminate the service software by closing the browser window.
- Select **System/End** and **Shut down the system**.
- Insert the **OEMsetup_1 english VB39A CD 1 of 5**
- **Power off the IES.**
- **Disconnect the network** cable of the IES from the Ethernet switch.
- Power on the **IES**.
 - ⇒ The computer boots from CD-ROM. After a short while, the following window appears:

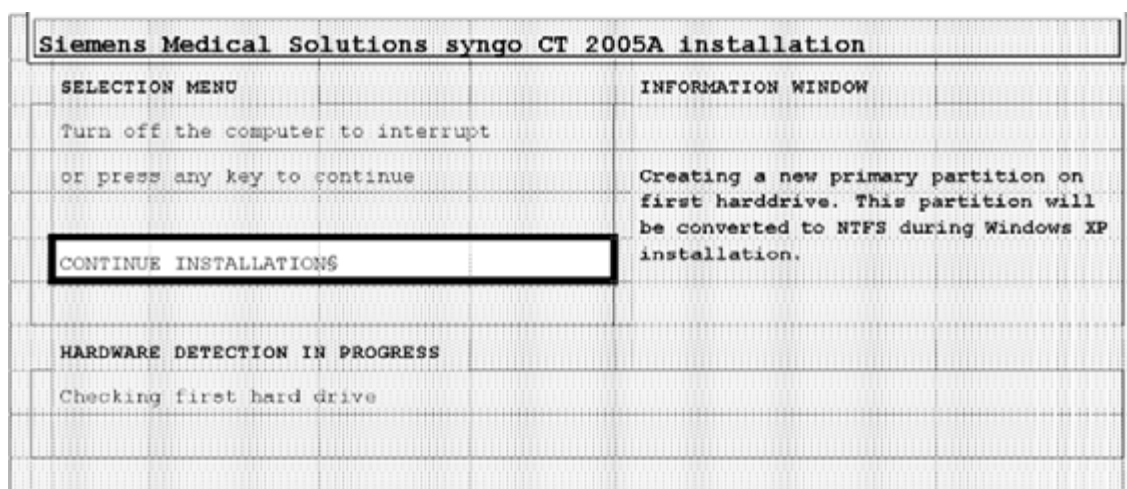


Fig. 9: SW installation Wizard

- Press **<Enter>** to continue installation (the IES is beeping). The IES restarts.
- Select distribution type **CUSTOMER** (default).
- Enter the **machine name** (host name) as noted by replacing the default entry "MedPC" and press Enter.

NOTE

Please note that the computer name must not contain blanks or underscores. Use the following naming convention for the IES computer name: ctw<serial number> (e.g. ctw12345, wizard connected to ct12345).

- Keep data on other hard drive(s) is the default and may not be modified (only one hard drive available)

NOTE

Do not use 'Clear all hard drive(s)'. This option will only format all hard drives (removes all software). No reinstallation of SW will be started.

- Use **CONTINUE INSTALLATION** to confirm the start of installation.
- Press any key.

- ⇒ All necessary files are copied from OEMsetup_1 english VB39A CD 1 of 5 to the hard drive.
- When the computer is beeping, take out OEMsetup_1 english VB39A CD 1 of 5 and insert **OEMsetup_2 english VB39A CD 2 of 5** and **press any key**
- Insert the **syngo CT 2008O Som/5_ICS/IES_SW/FW VB39A CD 3 of 5**.
 - ⇒ The installation continues automatically.
 - ⇒ After approx. 40 minutes, the CD is ejected.
- Insert the **DICOM_SR_Viewer CD 4 of 5**.
 - ⇒ The installation continues automatically.
 - ⇒ After approx. 3 minutes, the CD is ejected.
- Insert the **syngo CT 2008O Som/5_ICS/IES_SW/FW VB39A CD 5 of 5**.
 - ⇒ The installation continues automatically.
 - ⇒ After approx. 30 minutes, the CD is ejected

NOTE

If you miss inserting the requested CD in time, a warning is displayed:

No package has been found. Please insert the CD "....." in the CD drive to continue the installation.

Insert the requested CD and press any key to continue.

Installation running

- ⇒ The system reboots automatically after finishing the copy job.
- The following screen appears:

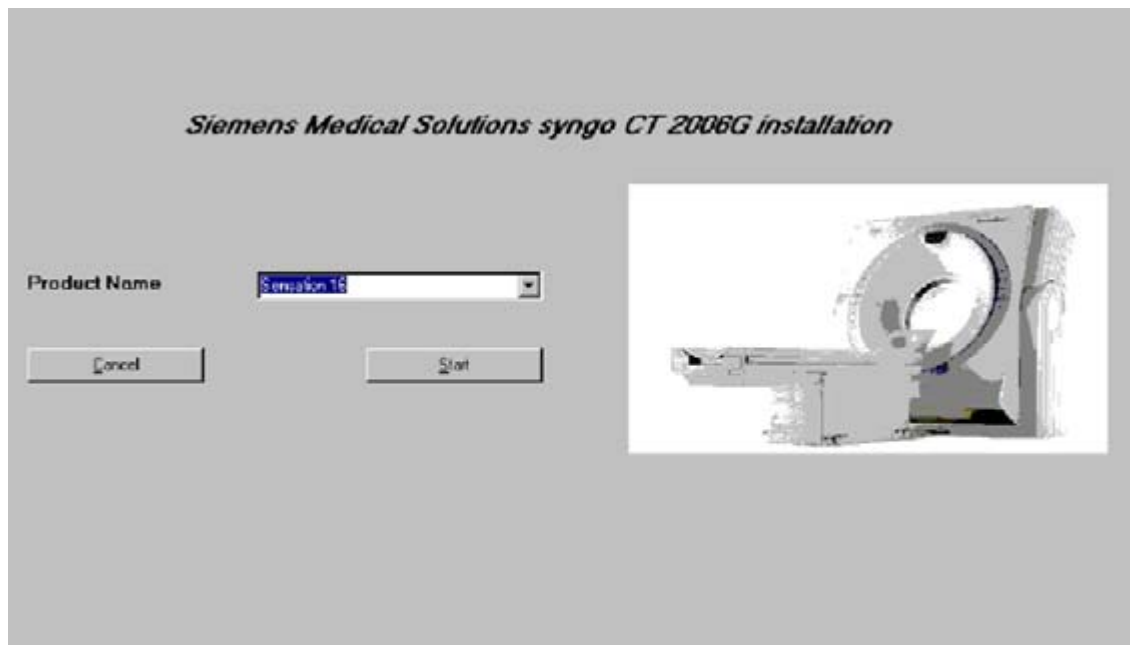


Fig. 10: Starting the installation

- Select **syngo CT Workplace** in the field "Product Name".
Use the drop down menu to change to syngo CT Workplace.
- Click Start

Completion of IES installation (approx. 15 minutes) is indicated by a pop-up window with the message **Please enter license**.

NOTE

Continue only when the Navigator is fully installed and started up.

- Click '**Click here to start the Local Service Configuration**' to start the service software (complete service license required).

Configure the keyboard language and enable Windows XP Access

Perform the following steps to configure keyboard languages **other than English**.

- Select **Start -> Settings -> Control Panel**
- Search for and click **Date, Time, Language and Regional and Options**
- Click **Regional and Language Options**.
- Select the **Languages -> Details** tab card
- Click **Add** and select the correct **Input Language** and confirm with **OK**.

- Select the English keyboard language, click **Remove** and confirm with **OK**.
 - ⇒ A message window appears informing you that removal is only possible after reboot. Confirm with **OK**
- Click **OK** in the remaining windows or close remaining windows.

Perform the following steps to **enable Windows XP** access at the time of the update.

- Enter the Service Key and open the Local Service main menu.
- From the service software home menu, select **Utilities**.
- Under Source, select **Escape to OS** and at Command **NT Command Interpreter**.
- Under Parameters, type **c:\Somaris\Tools\bin\KeyFilter.bat disable**.
- Activate with **Return** or **GO**.
 - ⇒ A reboot is required to make the change effective. This reboot will be performed after Restore 1.
- Click **Home**.

Restore Steps

NOTE

To perform the following steps, the Navigator has to be installed and running. Switch on the Navigator.

Restore 1

- In the main menu of the service software, select the **Backup & Restore** function.
- Insert the backup media for **Wizard**.
- From Command, select **Restore**.
- From Drives, select the **corresponding drive**.
- From Archive, select **SW-Settings02-Wiz_date_time.ar**.
- From Groups, select **all** except **User Settings** and click **Go**.
 - ⇒ **BackupRestore Running ...**; wait until **Ready** is displayed.
- From Archive, select **<hostname>_SecuritySettings-PXXX_date_time_ok.ar**.
- From Groups, select **all** and click **Go**.
 - ⇒ **BackupRestore Running ...**; wait until **Ready** is displayed.
- From Archive, select **CT-Settings-date_time_ok.ar**.
- From Groups, select **all** and click **Go**.
 - ⇒ **BackupRestore Running ...**; wait until **Ready** is displayed.
- Click **Home**.
- **Connect the network cable** of the IES to the Ethernet switch.

Restore 2

NOTE

The reboot takes a while, because the database is set up.

- Once started up, access the service software.
- In the main menu of the service software, select the **Backup & Restore** function.
- From Command, select **Restore**.
- From Drives, select the **corresponding drive**.
- From Archive, select **SW-Settings02-Wiz-<hostname>**.
- From Groups, select **User Settings** and click **Go**.
 - ⇒ **BackupRestore Running ...**; wait until **Ready** is displayed.
- From Archive, select **InSpace** (optional) and click **Go**.

- From Groups, select **%MEDHOME%\rt3d\preserve** and click **Go**.
 - ⇒ **BackupRestore Running ...**; wait until **Ready** is displayed.
- Click **Home** and terminate the service software by closing the browser window.

Completion of configuration

Regional settings

- From the syngo user interface, select -> **Options -> Configuration**.
- In the configuration panel, double-click the **Regional and Language Options** icon.
- In the **Regional Options** tab card, select a preferred country and click **Customize**.
- Change the date and time as noted in ([SW Installation Navigator / p. 15](#))

NOTE

The ARGUS option allows users to enter the parameters for "patient height" and "patient weight". For this reason, the correct measurement units have to be selected.

Number > Measurement System.

For cm/kg, the setting has to be set to "Metrics".

- In the Languages tab card, select the default input language and the keyboard type as noted in ([Starting the Installation / p. 17](#))

Setting the time

- In the main menu of the service software, select the **Configuration** function.
- Select **Site Info** and click >. Check/correct the **Time zone and the time**
Note: do not activate the NTP time server or select automatic Daylight Savings Time
- If changes were made, click **Save**. Confirm the message displayed with OK.

Restart

If no automatic restart was executed,

- Perform **System/End -> Restart the System** (restart is necessary to enable options and make the tab cards option available in the customer UI).

Siemens Remote Services

System Management installation
• Open the Local Service Home menu and go to Configuration. • At System Mgmt-> Master check and correct the IP addresses. • Click Save and then Next • At System Mgmt-> Agent Control, see whether the Installed checkbox is set. • Click Go • At System Mgmt-> Agent Control, set the Activated checkbox. ⇒ This may take as few minutes • Click Go • Click Finish
Set Remote Service Access Level
• Go to Options -> Service -> Remote Service • Set the Remote Service Access Control Mode to Limited access permanent

Software self test

Navigator

- Open the **Local Service Home** menu and go to **Test Tools**
- Click **Software SelfTest**



- Select the **Communication** group
- Press **Go**.



- Select the **Create SystemInfoFile**
- Press **Go**.

⇒ The file **SomSTSysInfo.info** will be created under
C:\Somaris\Tools\SelfTest\Communication\Local

Wizard

- Open the **Local Service Home** menu and go to **Test Tools**
- Click **Software SelfTest**



- Select the **Communication** group
- Press **Go**.




- Select the **Create SystemInfoFile**
- Press **Go**.

⇒ The file **SomSTSysInfo.info** will be created under
C:\Somaris\Tools\SelfTest\Communication\Local

Map drive C: of ICS (Navigator) at Wizard

Wizard

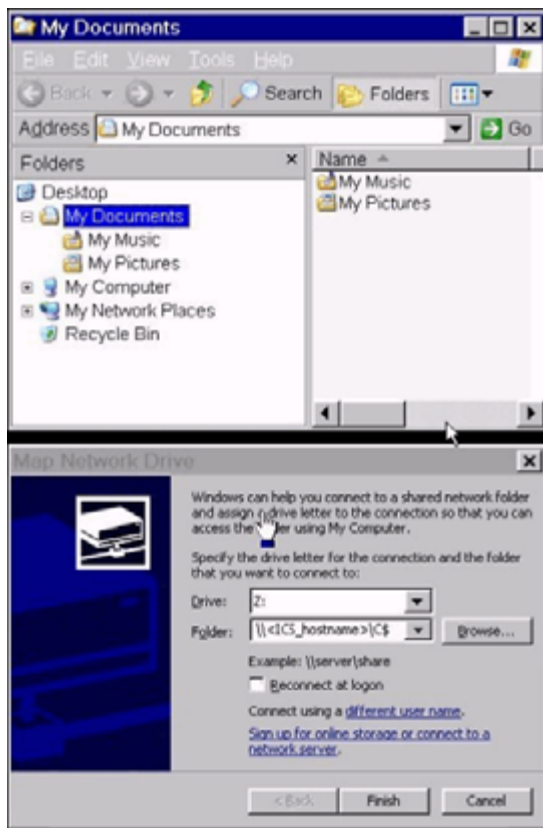

- Open Windows Explorer.
- Select **Tools->Map Network Drive**

- Select **Drive: Z:**
- Type

Folder: \\<ICS-Hostname>\C\$

⇒ In the window that is displayed, log on as the Administrator using the administrator password of the IES

Fig. 11:

- Copy the file **SomSTSysInfo.info**

from **C:\Somaris\Tools\SelfTest\Communication\Local**
to **Z:\Somaris\Tools\SelfTest\Communication\Remote**

and

from **Z:\Somaris\Tools\SelfTest\Communication\Local**
to **C:\Somaris\Tools\SelfTest\Communication\Remote**

Navigator and Wizard

- Open the **Local Service Home** menu and go to **Test Tools**
- Click **Software SelfTest** and execute all subtests.
 - Environment
 - Configuration
 - Security
 - Database
 - Misc
 - Communication (**only if the Wizard (IES) is connected**)
 - Performance

⇒ In case of errors, refer to the hints for troubleshooting

NOTE : If you have to perform troubleshooting, make sure that you recreate the **SomST-SysInfo.info** and copy and paste it between ICS and IES as described above. Otherwise, you will get the same error messages even though you have already solved the problem.

Navigator

- From the **Navigator**, select **System/End** and **Shut down the system**.
- Switch off the system at the **control box**.
- Wait a few seconds.
- Switch the system back on at the control box.

Note: Switching the system off and on is necessary to avoid network problems after the software self test.

Backup at the Navigator

- Start the service software and select the **Backup & Restore** function in the main menu.
- Insert the backup CD-R for the **Navigator**.
- From Command, select **Backup**.
- From Drives, select the CD drive, e.g., [-R-] CD-R.
- From Packages, select **CT-Settings-PXXX** and click **Go**.

⇒ **BackupRestore Running ...**

Command	Drives
Backup\$	[-R-] CD-R\$
	Packages
	CT-Scanner-PXXX\$

Fig. 12: Backup CT settings

⇒ Wait until **Ready** is displayed.

- Close the CD tray.
- From Packages, select **SW-Settings02-PXXX** and click **Go**.

⇒ **BackupRestore Running ...**

Command	Drives
Backup	[-R-] CD-R
	Packages
	SW-Settings02-PXXX

Fig. 13: Backup SW settings

⇒ Wait until **Ready** is displayed.

- Close the CD tray.
- From Packages, select **Security-Settings-PXXX** and click **Go**.

⇒ **BackupRestore Running ...**

Command	Drives
Backup	[-R-] CD-R
	Packages
	Security-Settings-PXXX

Fig. 14: Backup security settings

- From Packages, select **VSIM-Settings-PXXX** and click **Go**.

⇒ **BackupRestore Running ...**

⇒

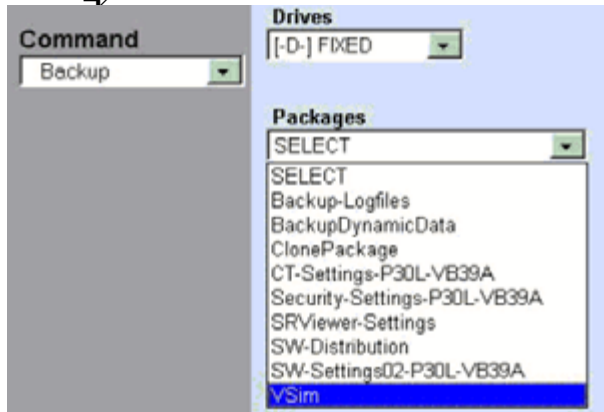


Fig. 15: Backup VSIM settings

⇒ Wait until **Ready** is displayed.

- When "Ready" is displayed, **remove the backup medium**.
- Click **Home** and terminate the service software by closing the browser window.

NOTE

The backup files have been saved correctly when **ok** is added at the end of the backup filename.

e.g. <hostname>-SW-Settings02-PXXX-date_time_ok.ar

Disable Windows XP Access

Perform the following steps to disable Windows XP access for the customer. Administrator privileges are required; you must use the command shell of the service software.

- From the service software home menu, select **Utilities**.
- Under Source, select **Escape to OS** and at Command **NT Command Interpreter**.
- Under Parameters, type **c:\Somaris\Tools\bin\KeyFilter.bat enable**.
- Activate with **Return** or **GO**.
- Reboot** to make the change effective.
 - To enable Windows XP access, type in
 - c:\Somaris\Tools\bin\KeyFilter.bat disable**

Backup at the Wizard

Wizard

- Start the service software and select the **Backup & Restore** function in the main menu.
- Insert the backup CD-R for the **Wizard**.
- From Command, select **Backup**.
- From Drives, select the corresponding drive.
- From Packages, select **CT-Settings-IES** and click **Go**.

⇒ **BackupRestore Running ...**

Command	Drives		
Backup	[-R-] CD-R		
	Packages		
	CT-Settings-IES		

Fig. 16: Backup CT settings

- ⇒ Wait until **Ready** is displayed.
- ⇒ Close the CD tray.
- From Packages, select **SW-Settings02** and click **Go**.
- ⇒ **BackupRestore Running ...**

Command	Drives		
Backup	[-R-] CD-R		
	Packages		
	SW-Settings02		

Fig. 17: Backup SW settings

- ⇒ Wait until **Ready** is displayed.
- ⇒ Close the CD tray.
- From Packages, select **Security-Settings-IES** and click **Go**.
- ⇒ **BackupRestore Running ..**

Command	Drives		
Backup	[-R-] CD-R		
	Packages		
	Security-Settings-IES		

Fig. 18: Backup security settings

- ⇒ Wait until **Ready** is displayed.
- ⇒ Close the CD tray.

- From Packages, select **InSpace** (optional) and click **Go**.
 ⇨ BackupRestore Running ..

Command	Drives			
Backup	[-R-] CD-R			
	Packages			
	InSpace			

Fig. 19: Backup InSpace

- ⇨ Wait until Ready is displayed.
- ⇨ Close the CD tray.
- When "Ready" is displayed, **remove the backup medium**.
- Click **Home** and terminate the service software by closing the browser window.

Disable Windows XP Access

Perform the following steps to disable Windows XP access for the customer. Administrator privileges are required; as a result, deactivation is performed using the command shell of the service software.

- From the service software home menu, select **Utilities**.
- Under Source, select **Escape to OS** and at Command **NT Command Interpreter**.
- Under Parameters, type **c:\Somaris\Tools\bin\KeyFilter.bat enable**.
- Activate with **Return** or **GO**.
- **Reboot** to make the change effective.
 - To enable Windows XP access, type in
 - **c:\Somaris\Tools\bin\KeyFilter.bat disable**

Chapter	Change	Reason
All	First version.	Change no. 511526

H	hz_serdoc_F09G07U01M01(p30).	7
	hz_serdoc_F09G08U01M01(p30).	8
	hz_serdoc_F13G14U01M01(p30).	7
	hz_serdoc_F13G14U02M01(p30).	7