


6-2-5 DEINSTALL PAC MODULE



The existing PAC Module will be replaced by a new PAC 2B Module with independent vector gating.

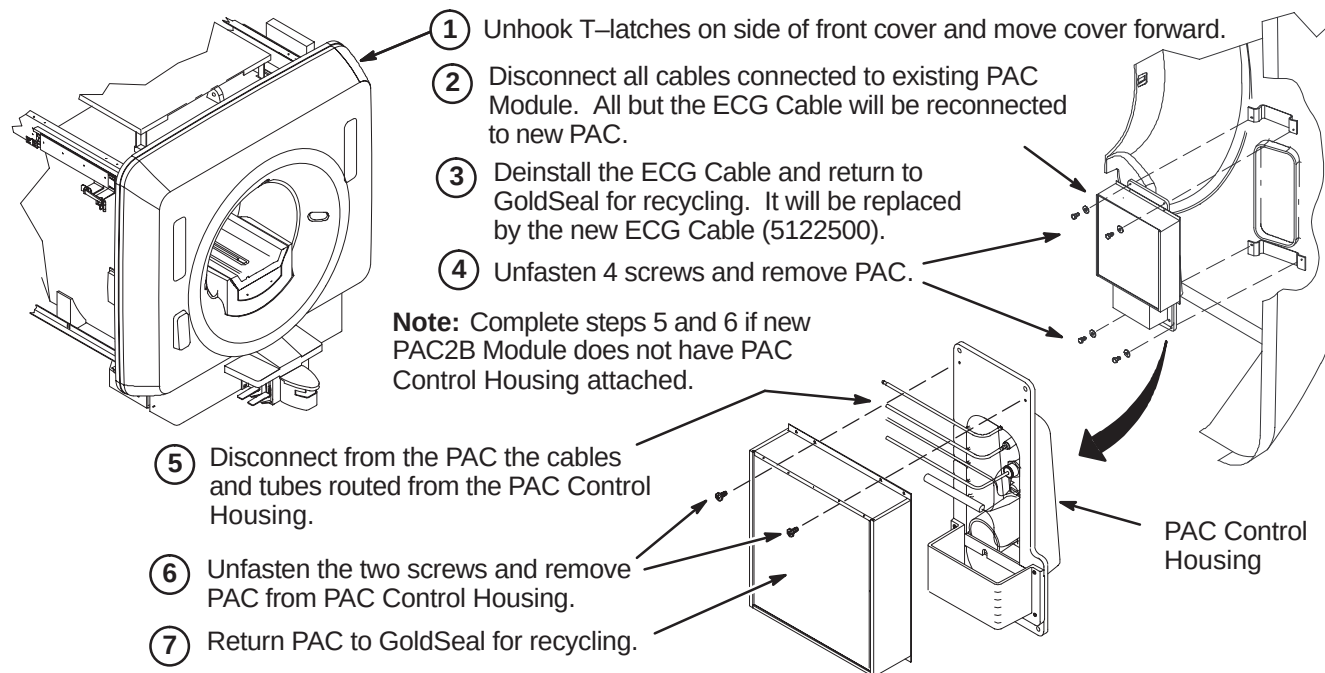

6-3 INSTALL NEW DUAL TEMPERATURE SENSOR CABLE ON FRONT END BELL


The existing Temperature Sensor Cable, Run 749, needs to be removed and must be replaced with new Dual Thermal Sensor cable, Run 1260.

6-3-1 UNFASTEN FLAT PANEL BORE LIGHTS

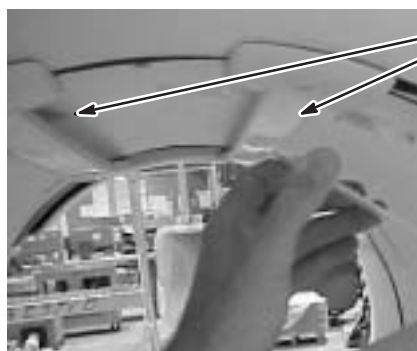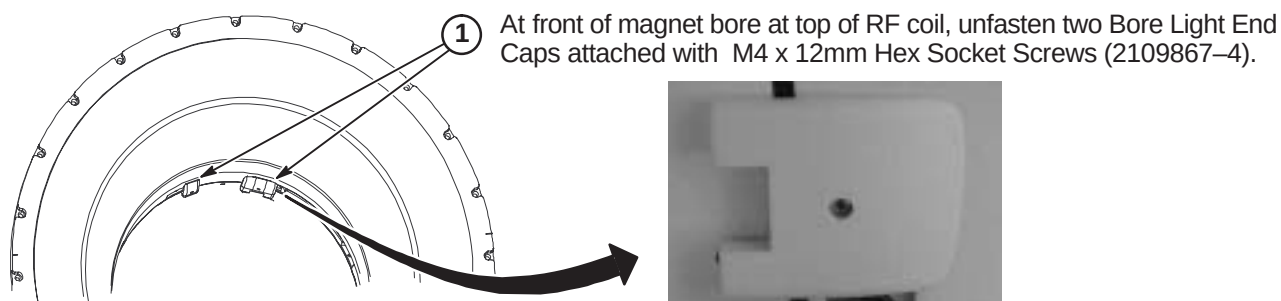


- ② Carefully lower the two Flat Panel lights as not to damage the fiber optic bundles they are attached to at the rear of the bore.

6-3-2 REMOVE FRONT END BELL

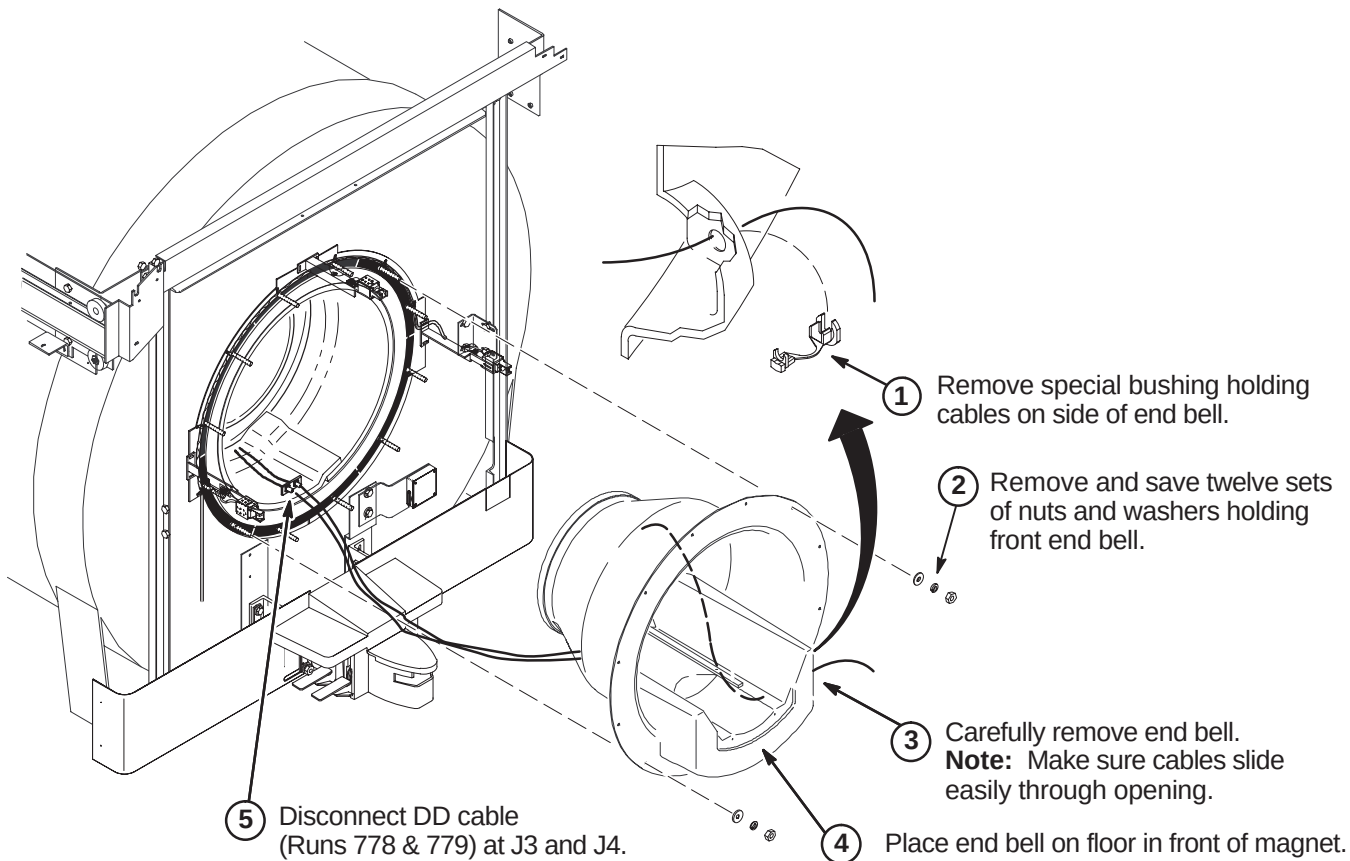

6-3-3 DE-INSTALL EXISTING RUN 749


The site may have upgraded to a Dual Thermal Sensor Run 749 recently, but this cable has to be removed and replaced with new Dual Thermal Sensor cable, Run 1260 (5123812). Signa HDx requires Run 1260 to be installed.

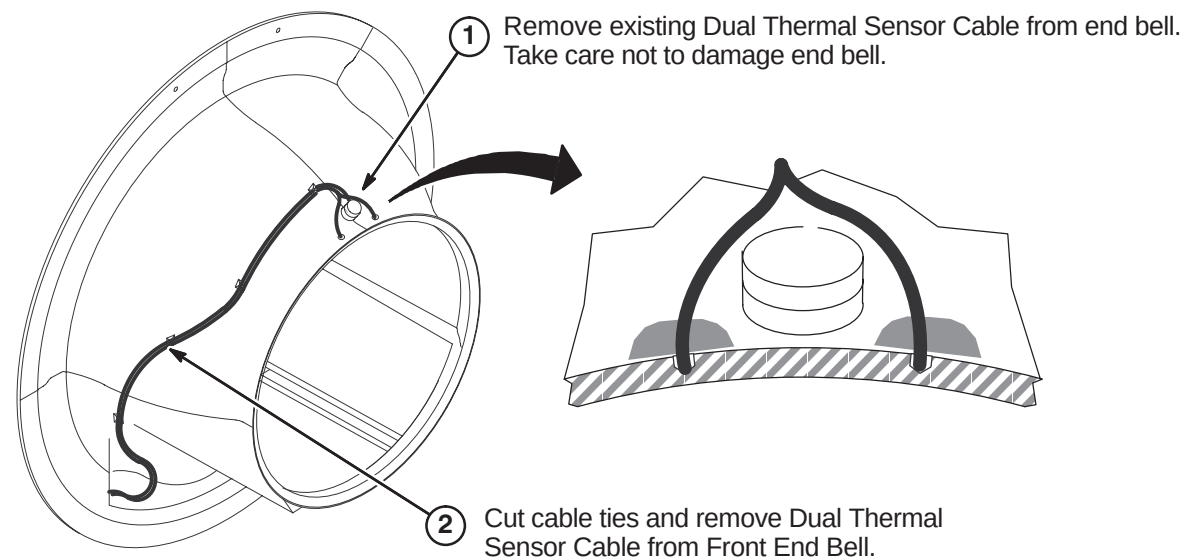


6-3-4 PREPARE END BELL FOR INSTALLATION OF DUAL THERMAL SENSOR CABLE

WARNING!

FERROUS MATERIAL HAZARD! DO NOT DRILL, HAMMER, CHISEL, ETC. OR USE ANY FERROUS TOOLS IN MAGNET ROOM UNLESS MAGNET IS RAMPED DOWN!

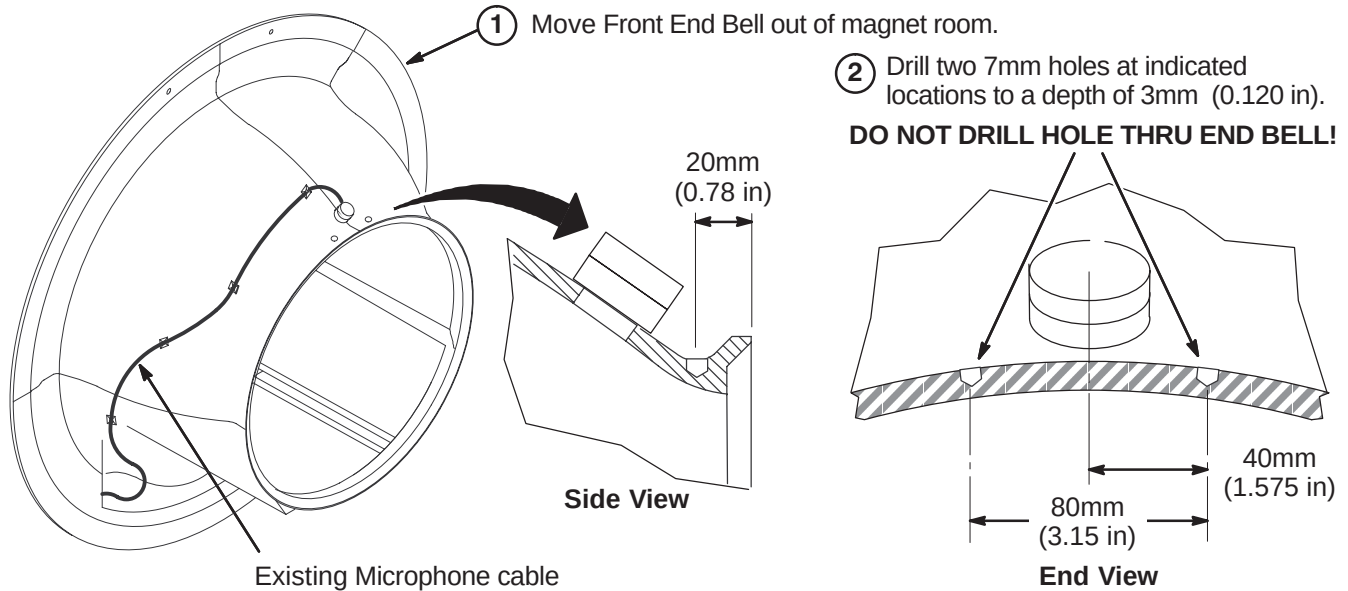


6-3-5 ATTACH DUAL THERMAL SENSOR CABLE

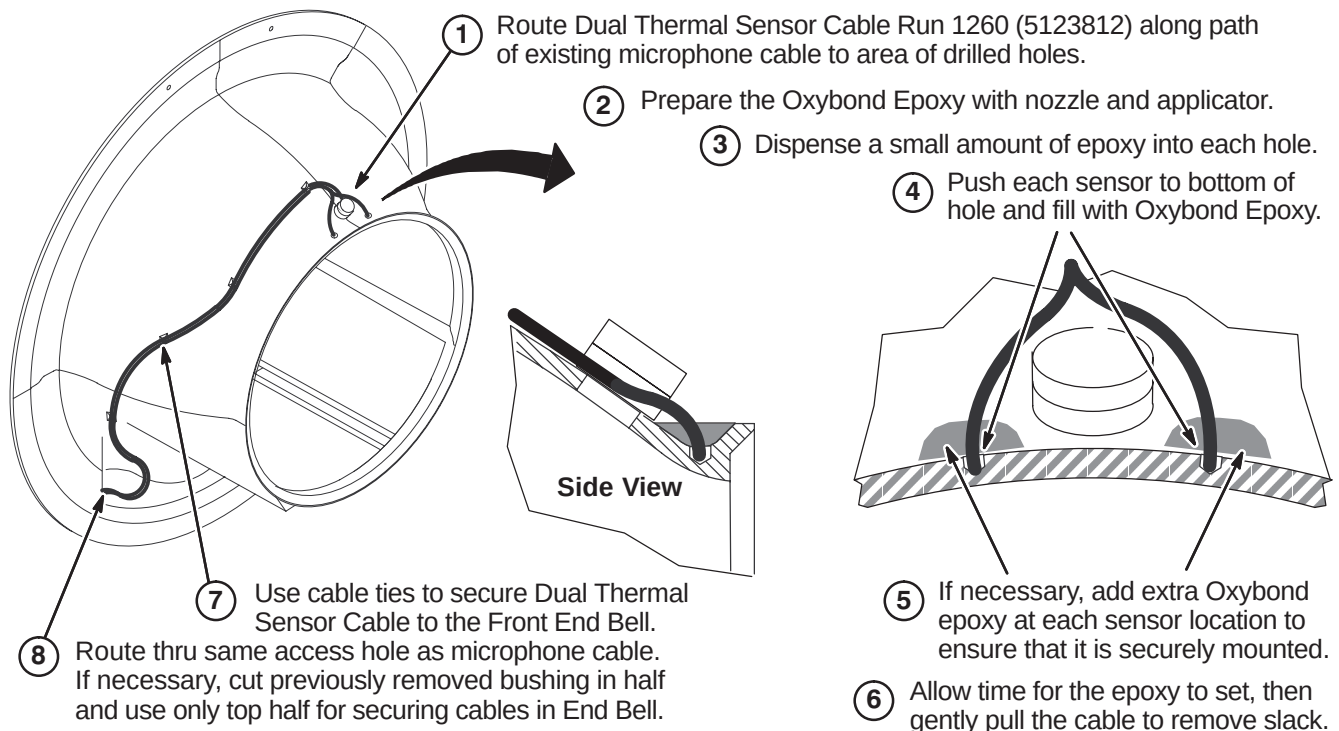

6-3-6 ATTACH FRONT END BELL TO MAGNET

Note: A good seal is required to prevent air leaks. Any air leaks between RF Tube (Gradient Coil) and End Bell joint can cause a loud humming noise during patient scanning. If possible, conduct an air leakage test prior to Bridge installation. This would require connecting power to the Blower Box and connecting 6" blower hose from Blower Box to Bulkhead.

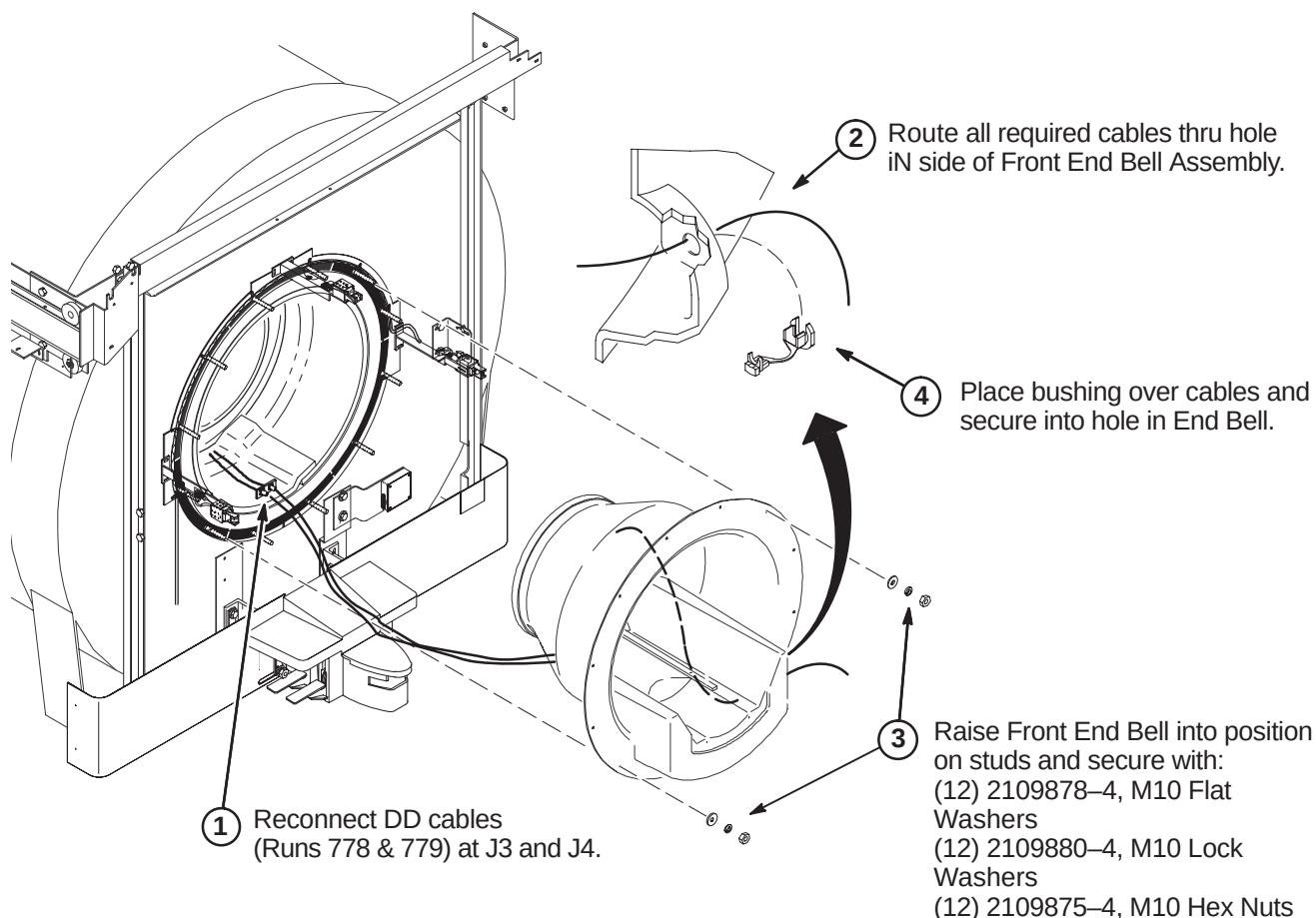

6-3-7 ATTACH BORE LIGHTS

