

eVO HP Tip for MR CSEs: Install or Update Software for VLM Gradient Supervision

By Dean Wilson, CSE
02-FEB-2018

This procedure describes step-by-step how to download and install the VLM software (hereon, “**the software**”) on eVO Windows 10 devices. A prerequisite is the installation of the US232-100 RS232-to-USB adapter (see SSN “MRI-CSE” stream) on Windows 10 eVO systems.



Figure 1 VLM Gradient Supervision Unit



Figure 2 VLM Control Software Platform

The tablet MUST be connected to the **Siemens Intranet**.

Browse to the Siemens MR Product Information page:

<https://productinformation.healthcare.siemens.com/ProductInformation.aspx?uuid=7abf6bb6-5daf-440c-9920-fab2d4640084>

If the link above fails, copy and paste the URL into an intranet browser.

Click on **Downloads**

Click on **Gradient Supervision**

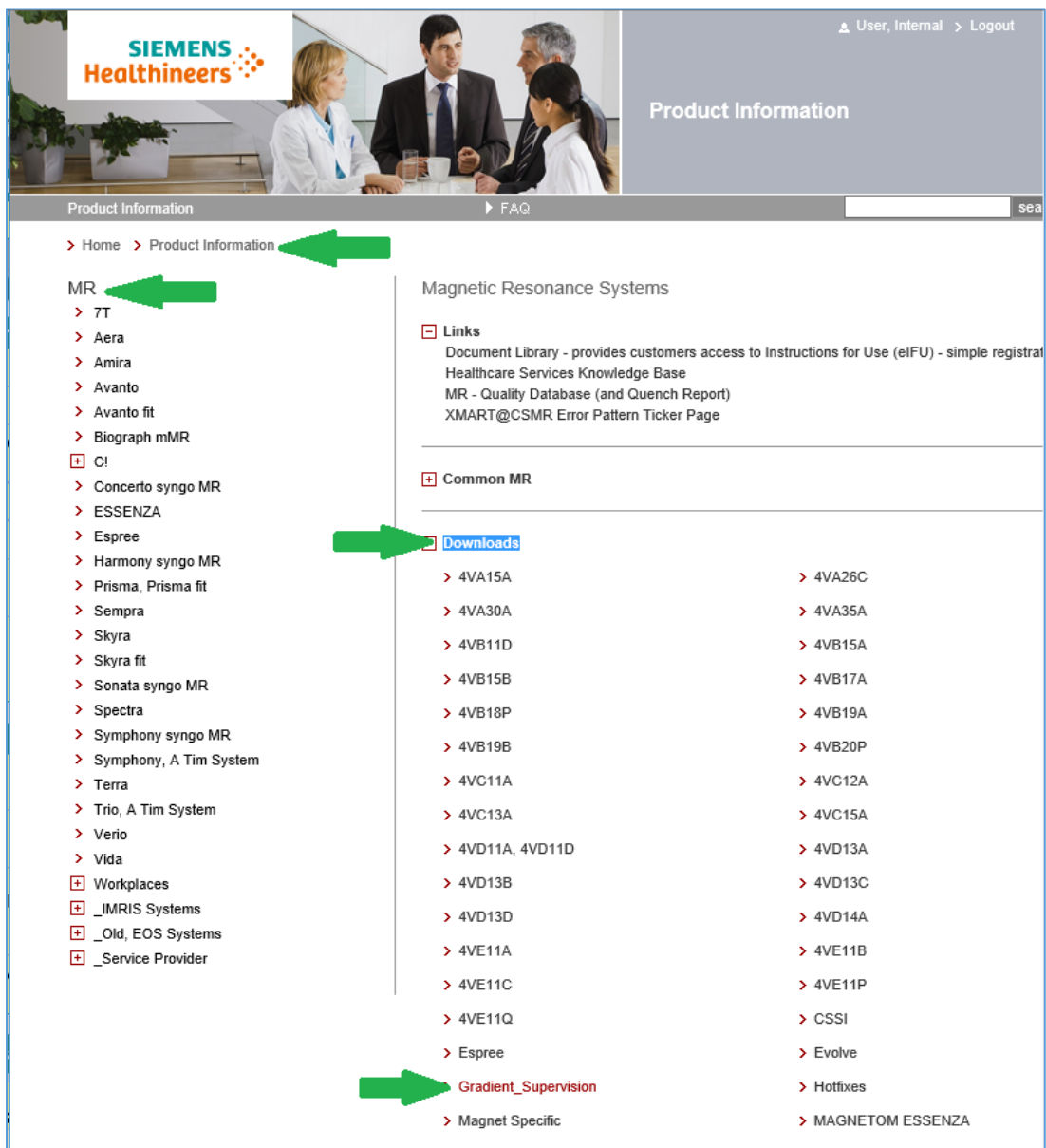


Figure 3 MR Product Information Page

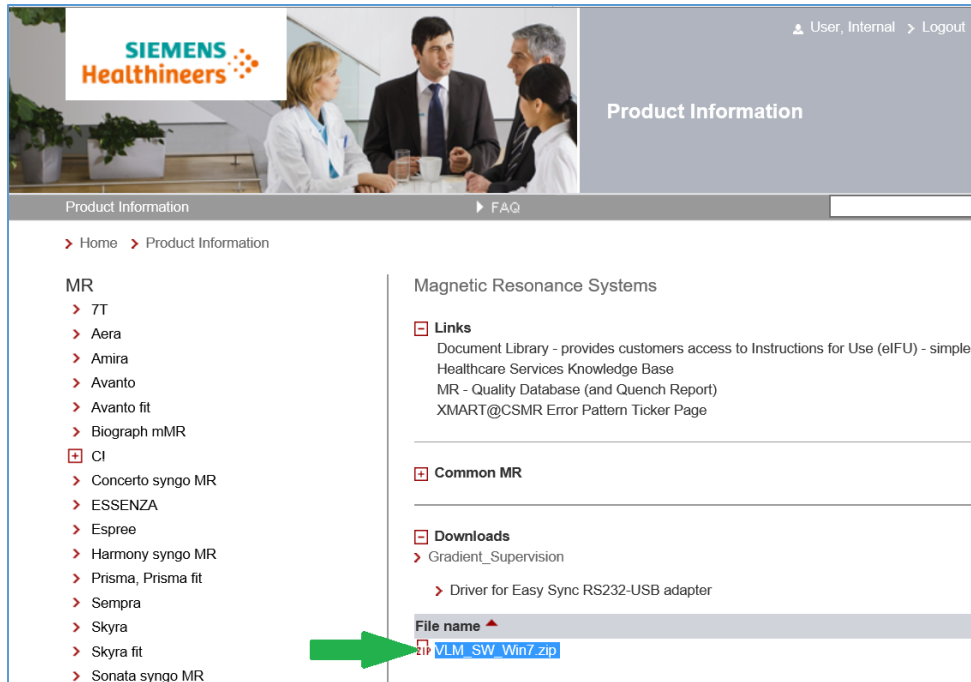


Figure 4 Downloads Menu - Gradient Supervision

Click on the VLM software ZIP file link

NOTE: The Windows7 (Win7) version functions in Windows 10.

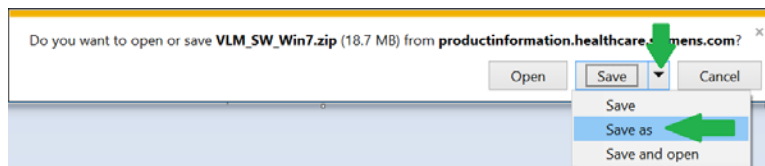


Figure 5 Download Popup

In the download popup bar, click **Save**, then **Save as**

Browse to a selected folder on your tablet to save the file, and click **Save**.

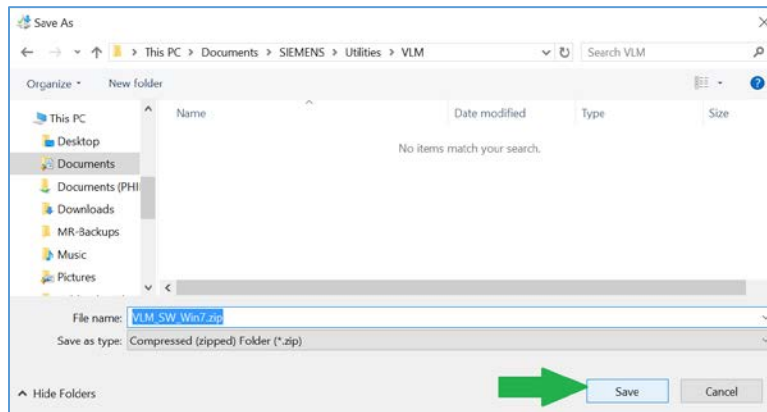


Figure 6 Save As to Folder

When the file download is complete, click **Open folder**

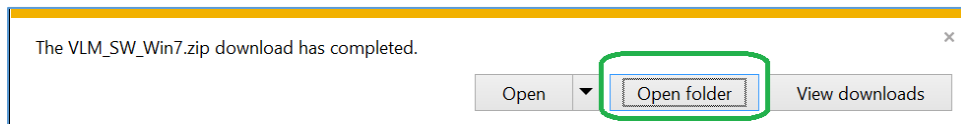


Figure 7 Download Popup Open Folder

Right-click on the file name and select 7-Zip, then Extract to "VLM_xx_xxx"

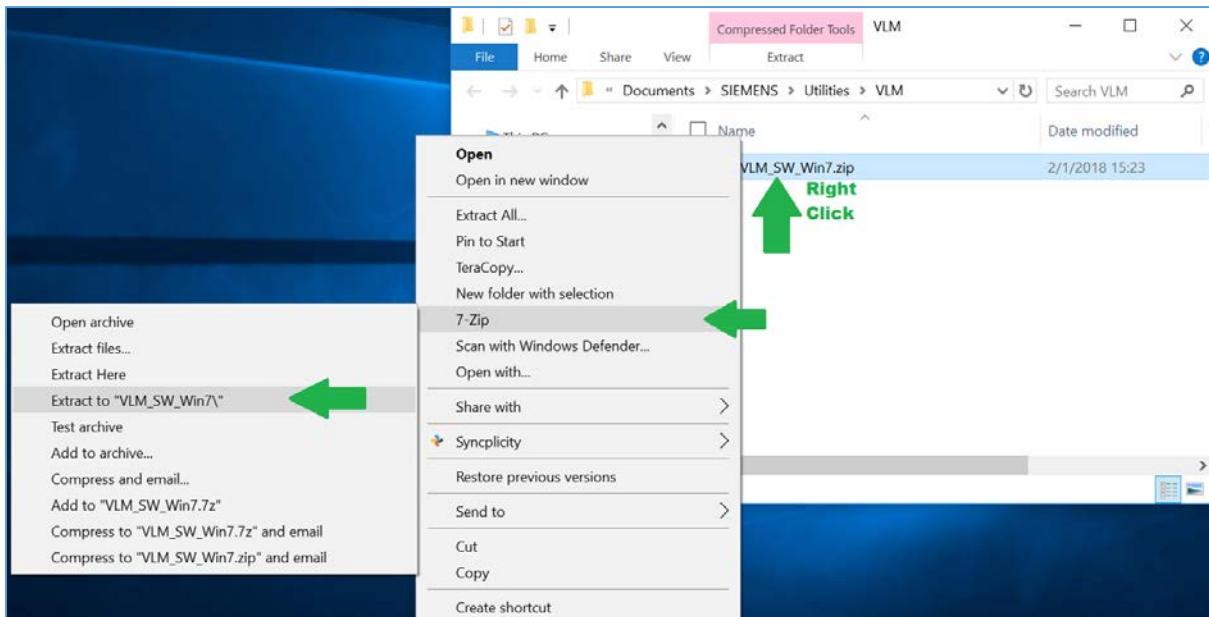


Figure 8 Extract Zip File

The folder and files are extracted

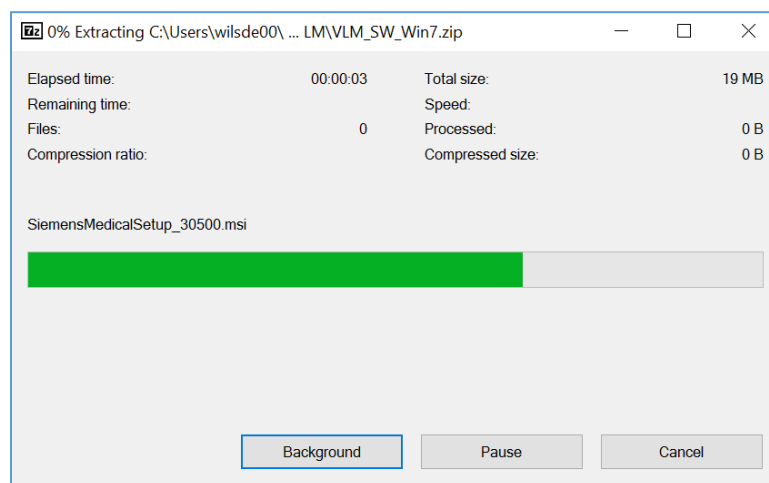


Figure 9 Zip File Extracting

Double-click on the folder to open it.

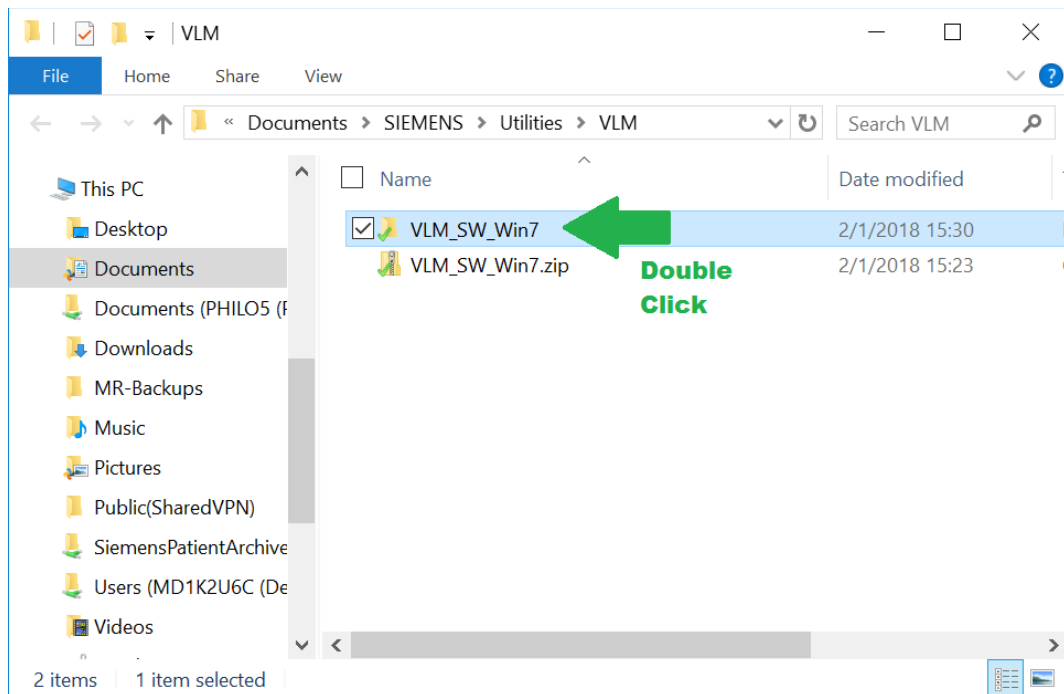


Figure 10 Open Extracted Folder

Double-click on the installation file to begin software installation.

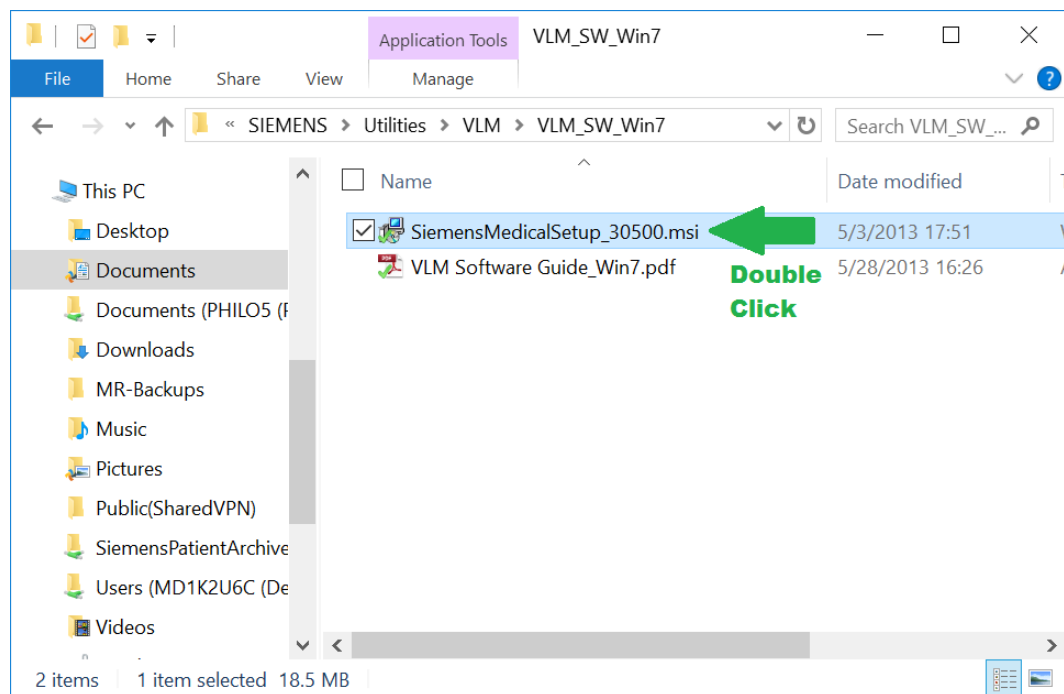


Figure 11 Begin Installation

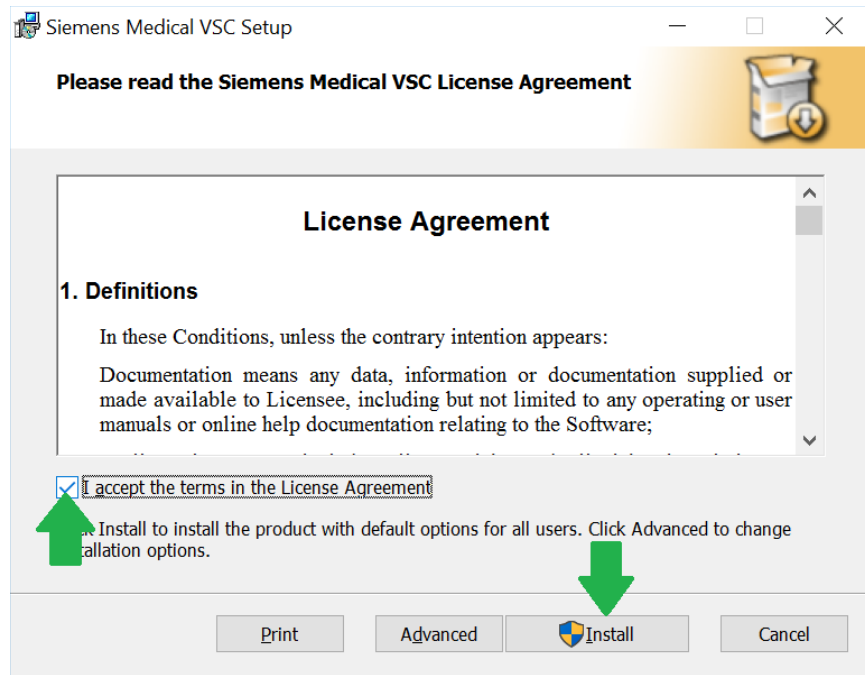


Figure 12 Accepting the EULA

The installation files extract

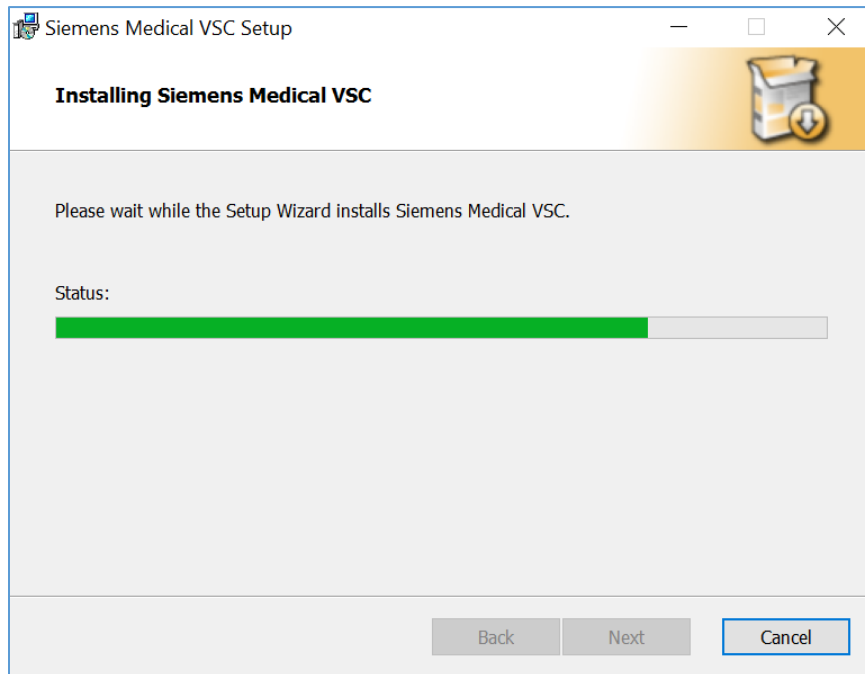


Figure 13 Installation Extracting

When complete, click **Finish** to close the Installation Wizard

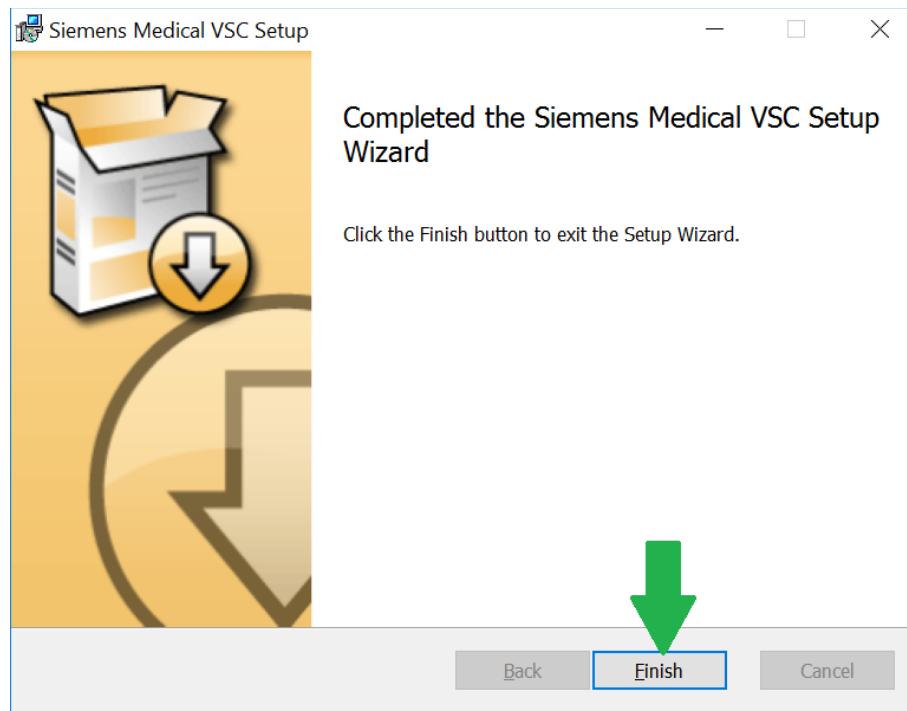


Figure 14 Closing the Wizard

An icon for the software appears on the desktop

Click on the START button and locate the icon for the software also.

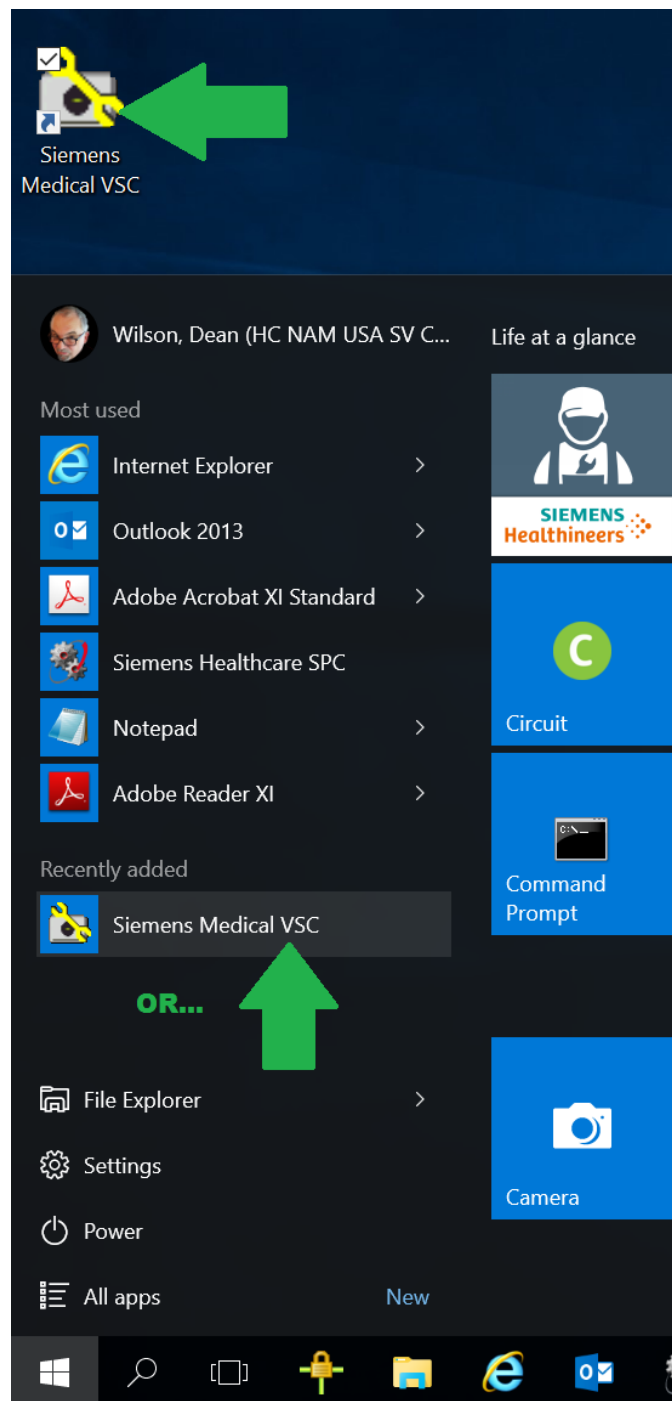


Figure 15 Software Icon on Desktop and START Menu

Double-click on the icon to launch the software.

In the software menu, click **C**onnection then **C**onnect and View



Figure 16 Software Connection Menu

Click on the **D**etails >>> button

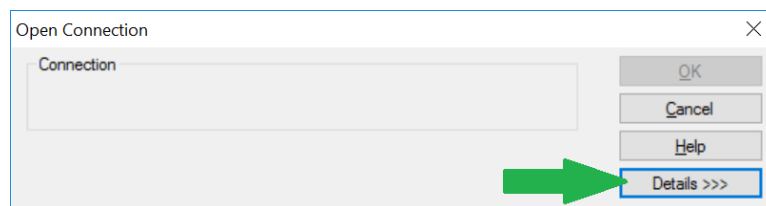


Figure 17 Open Connection Details Button

Connect the USB-to-RS232 adapter to the USB port and to the VLM monitor.
(this device driver must have already been installed previously)

In the Open Connection platform, click **Add...**

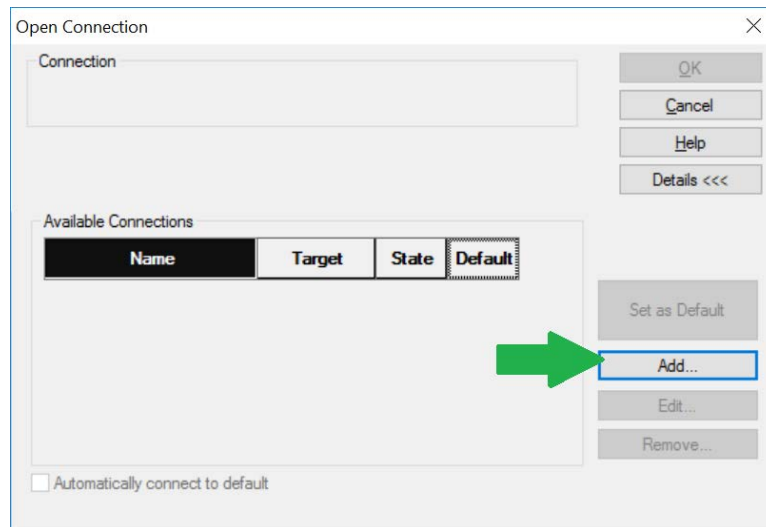


Figure 18 Add Open Connection

Select the pulldown menu to see a list of available COM ports. Select the port matching the USB device seen in Device Manager (not shown). In this example, COM4 is selected. This COMx port may vary.

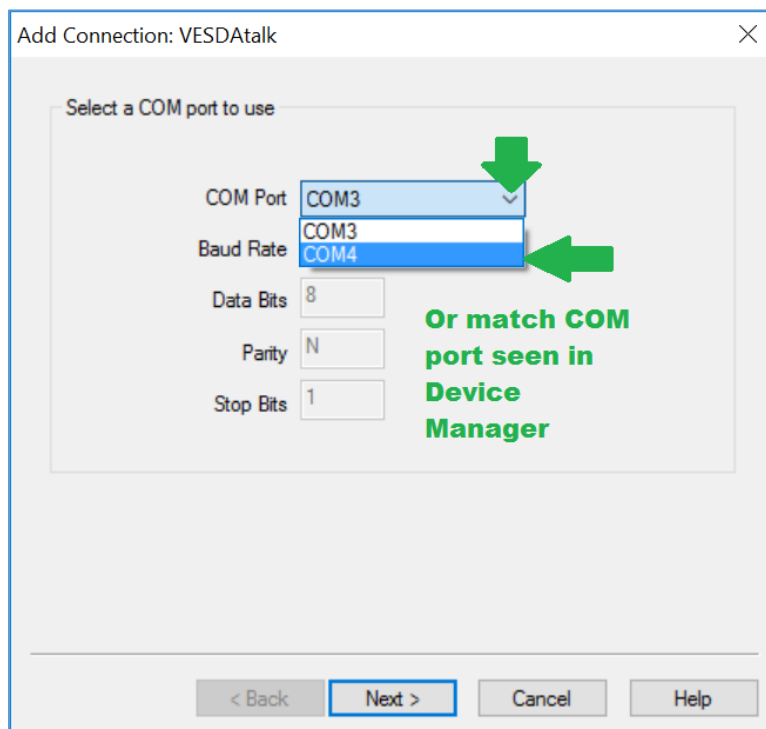


Figure 19 Select COM port

The remaining parameters are set by default. Click **Next >**

Click **Finish**.


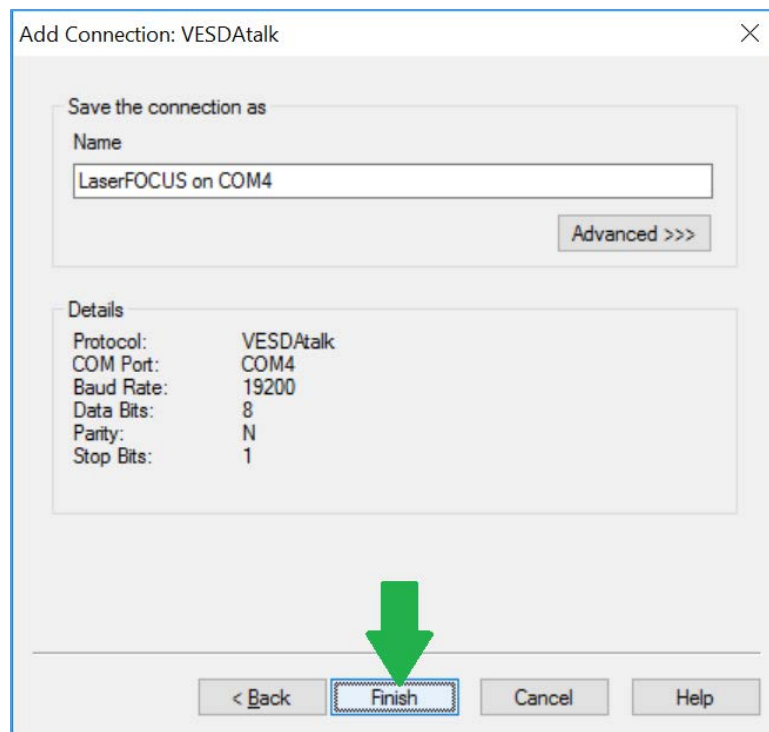


Figure 20 Finish Software Configuration

Click **OK**

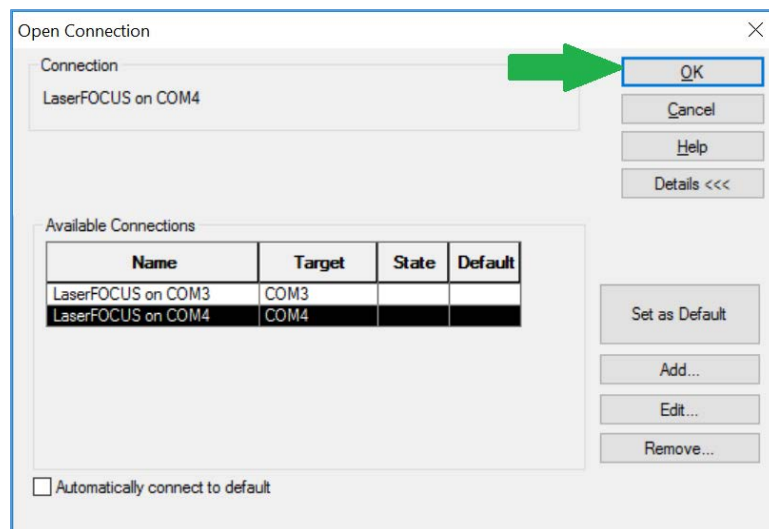

The software connects to the VLM monitor.

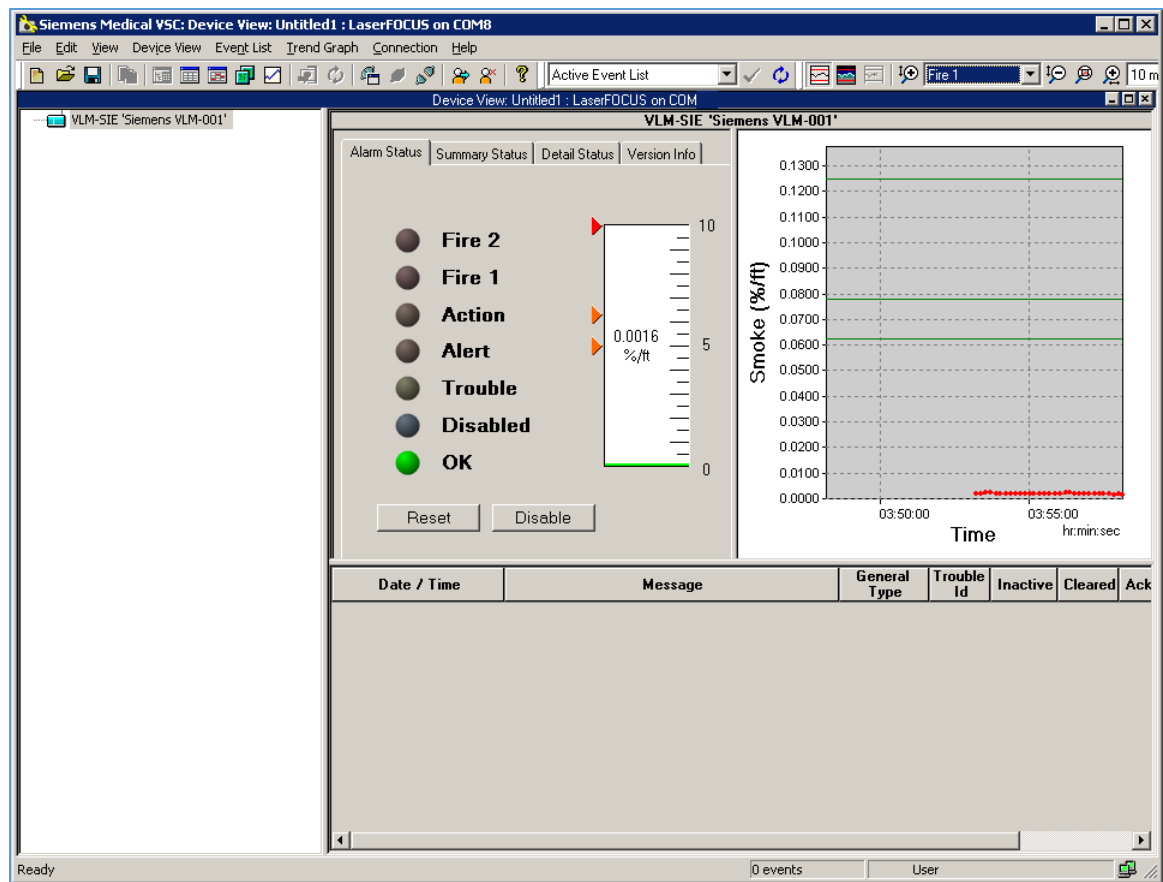


Figure 21 VLM Monitor Software Platform

TROUBLESHOOTING

If an error popup appears,

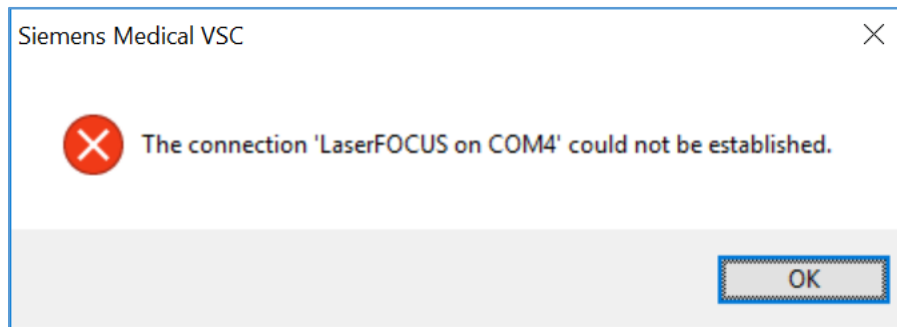


Figure 22 Connection Error

- Recheck the connections
- Verify COM port settings using Device Manager
 - Right-click on **START** button
 - Select **Device Manager**
 - Scroll down to **Ports (COM & LPT)**
 - Locate the port specific to the RS232 adapter and its port designation (i.e. **COM4**)
 - Verify the designated COM port is matched in the VLM configuration.
- If only one COM port appears, the USB-RS232 device is not connected properly, or the drivers are not installed properly.

Instructions on update/installation of the US232-100 device are available on the Siemens Social Network stream "MRI CSE" at [this link](#).