


TABLE OF CONTENTS

TABLE OF CONTENTS	1
1 - COMPONENT LOCATIONS	2
1-1 System Support Module (SSM).....	3
1-1-1 Front SSM Detail	4
1-1-2 Rear SSM.....	4
1-1-3 Rear SSM Connector Descriptions	5
1-2 RFI Module.....	6
1-2-1 RFI Module Front Connector Descriptions	6
1-2-2 RFI Module Rear Connector Descriptions	7
1-3 1.0T TELI Solid State RF Amps.....	8
1-3-1 1.0T TELI Solid State RF Amp Capacitors	9
1-3-2 1.0T TELI Solid State RF Amp Fans.....	10
1-3-3 1.0T TELI Solid State RF Amp Power Supplies.....	11
1-3-4 1.0T TELI Solid State RF Amp Fuse Location	13
1-3-5 1.0T TELI Solid State RF Amp Connector Descriptions	13
2 - INTERCONNECTS.....	14
2-1 1.0T SRF Cabinet Interconnects.....	14
2-2 System Support Module Interconnects	15
REVISION HISTORY	16

Description - This section applies to Signa systems installed with the 1.0T SRF cabinet.

1 - COMPONENT LOCATIONS


See Illustration 1-1.

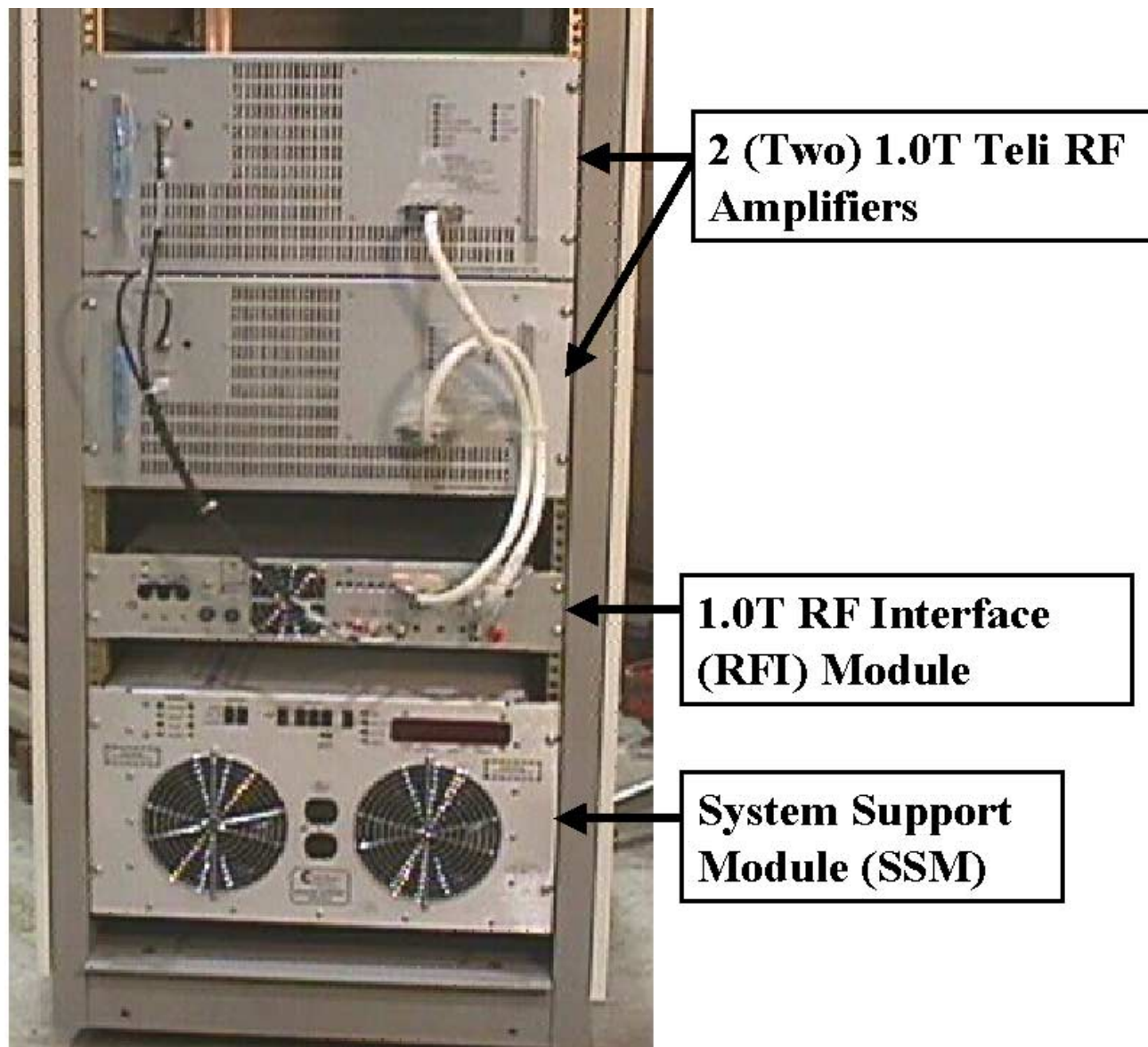

SRF CABINET LOCATIONS
ILLUSTRATION 1-1

1-1 System Support Module (SSM)


See Illustration 1-2.

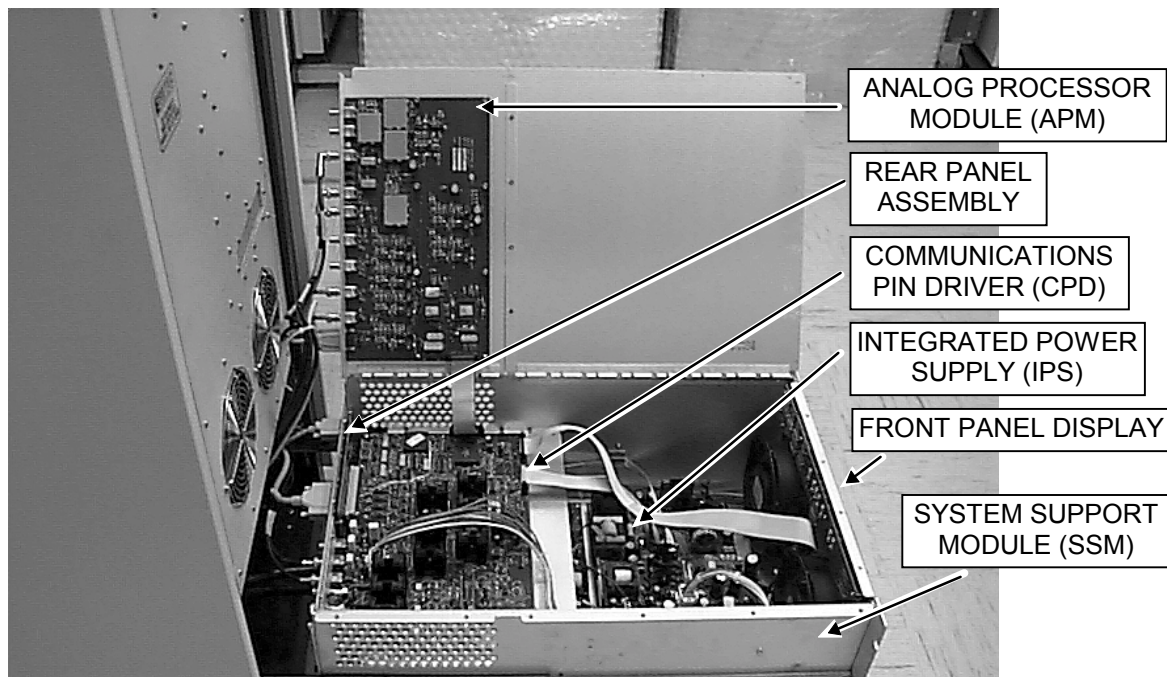

SSM COMPONENT LOCATIONS
ILLUSTRATION 1-2

1-1-1 Front SSM Detail


See Illustration 1-3.

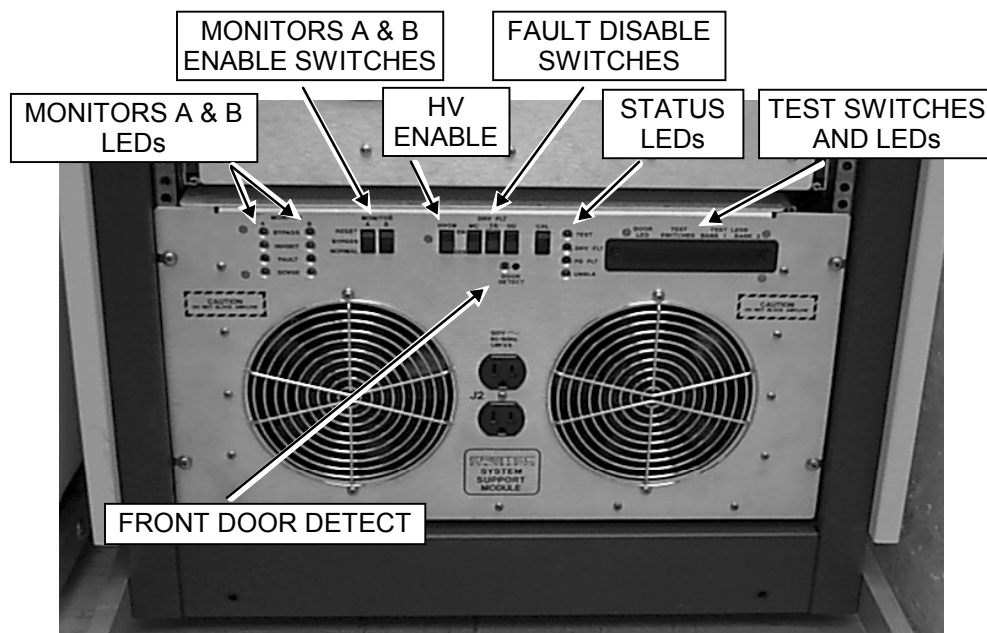

FRONT SSM LOCATIONS DETAIL
ILLUSTRATION 1-3

1-1-2 Rear SSM

See Illustration 1-4.

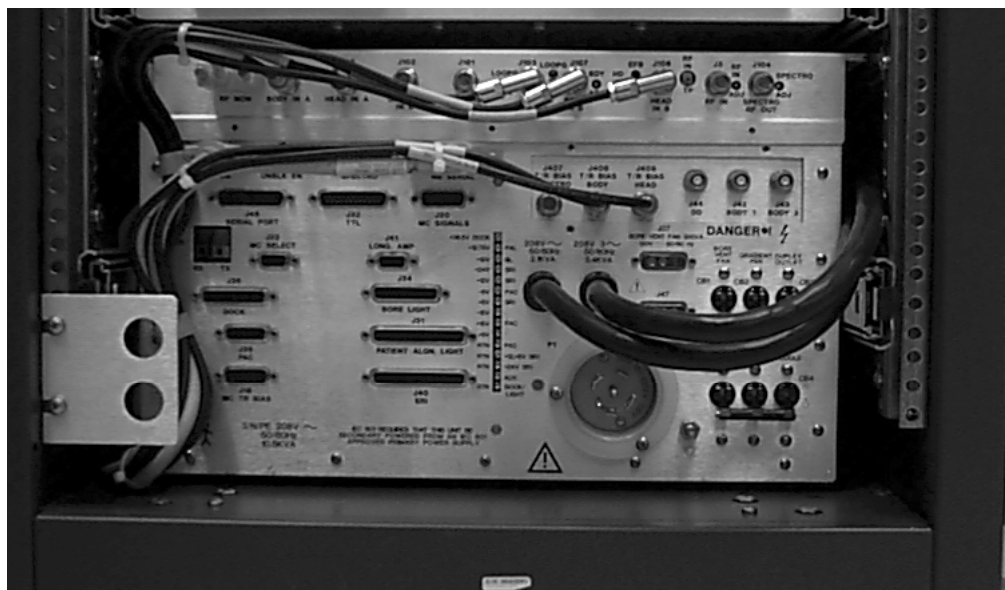

1-1-3 Rear SSM Connector Descriptions


See Table 1-1.

TABLE 1-1
REAR SSM CONNECTOR DESCRIPTIONS

CONNECTOR	DESCRIPTION	CONNECTOR	DESCRIPTION
J110	BODY IN A	CB4	SYSTEM SUPPORT MODULE MAIN
J108	HEAD IN A	P1	AC IN
J107	BODY IN B	J35	PAC
J106	HEAD IN A	J40	SRI
J104	SPECTRO RF OUT	J32	TTL
J103	EFB RF OUT	J41	LONGITUDINAL AMP
J102	SPECTRO IN B	J36	DOCK
J101	SPECTRO IN A	J34	BORE LIGHT
J43	BODY 2	J31	PATIENT ALIGNMENT LIGHT
J44	DD	J37	BORE VENT FAN
J42	BODY 1	J47	GRADIENT FAN
J407	T/R BIAS SPECTRO	J3	RF IN
J408	T/R BIAS BODY	J46	SERIAL PORT
J409	T/R BIAS HEAD	J5	UNBLANK ENABLE
J501	NB SERIAL PORT	J20	MULTICOIL SIGNALS
J503	RF MON	J22	MULTICOIL SELECT
J504	PCM	J18	MC TR BIAS
J507	SPECTRO	J2	HEAD/BODY DIRECTIONAL COUPLER OUTPUT
CB1	BORE VENT	RX	MDS RECEIVE
CB2	GRADIENT FAN	TX	MDS TRANSMIT
CB3	DUPLEX OUT		

1-2 RFI Module

See Illustration 1-5.

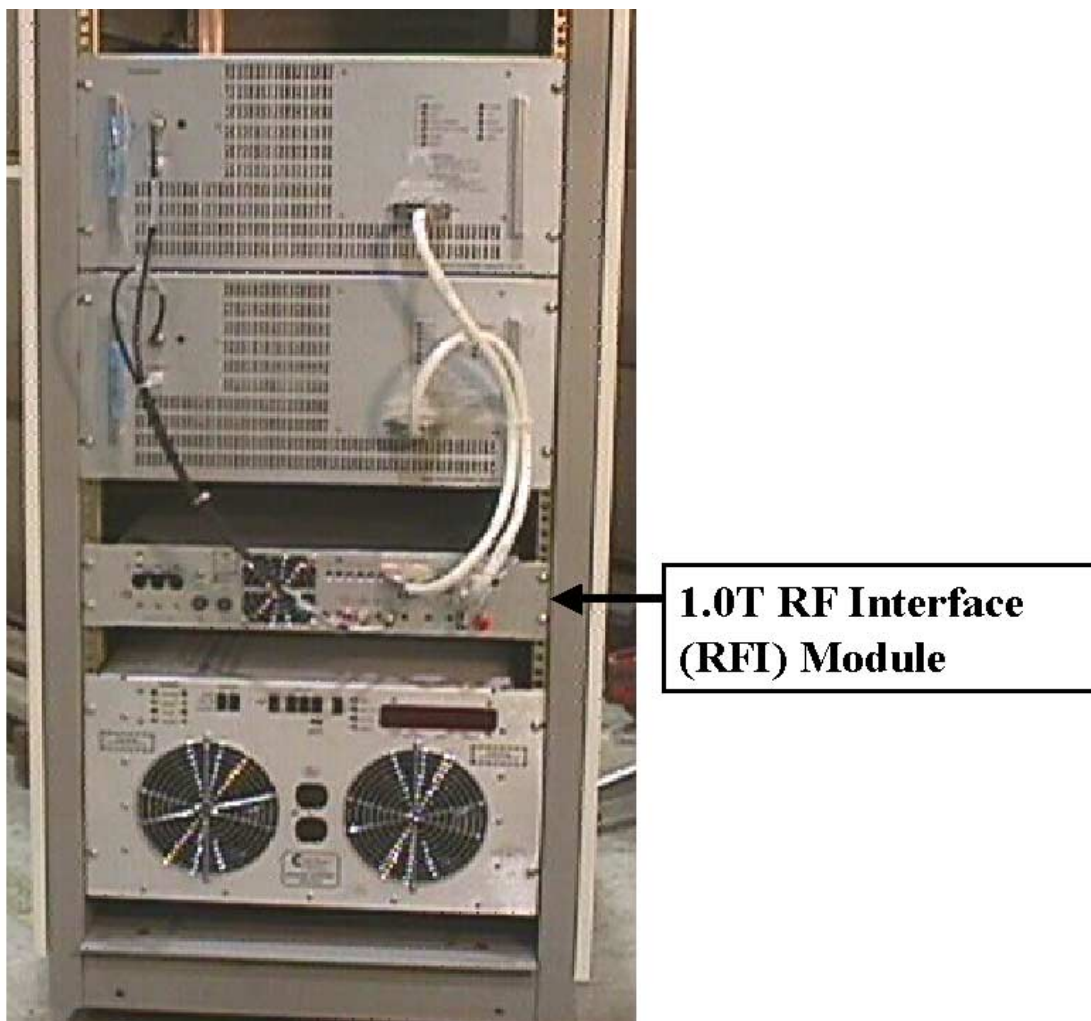

RFI MODULE LOCATION
ILLUSTRATION 1-5

1-2-1 RFI Module Front Connector Descriptions


See Table 1-2.

TABLE 1-2
FRONT RFI MODULE CONNECTOR DESCRIPTIONS

Connector	Description
J15	Balance Out
J11	Amp 1 Drive
J12	Amp 2 Drive
J14	RF Input
J13	RF Test Port
J16	Amp 1 Comm
J17	Amp 2 Comm

1-2-2 RFI Module Rear Connector Descriptions

See Illustration 1-6 for locations and Table 1-3 for descriptions.

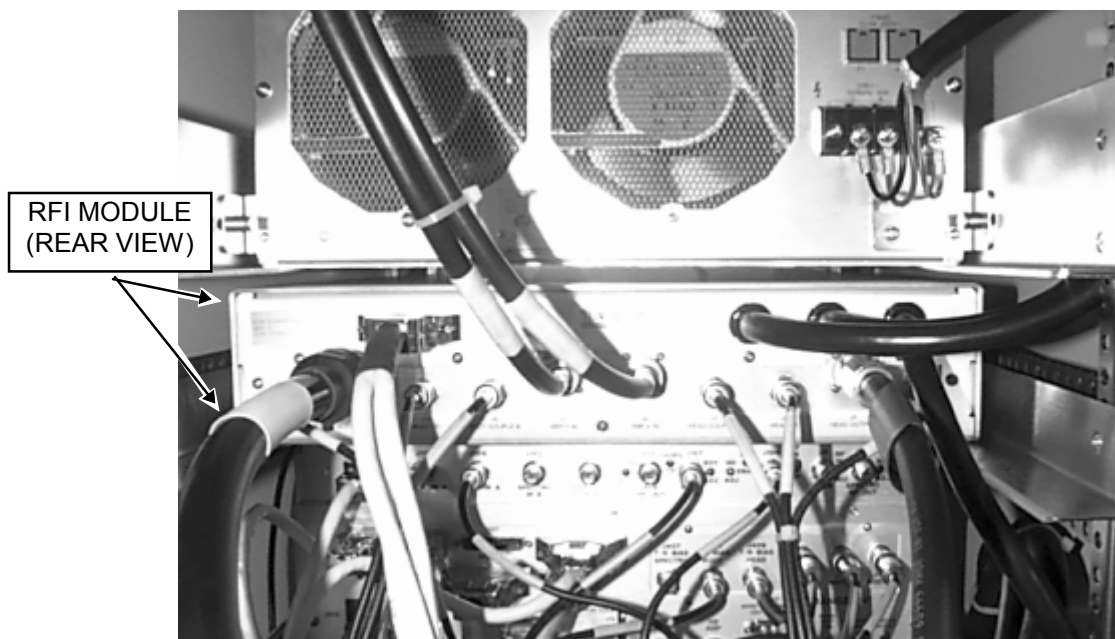


RFI MODULE REAR CONNECTIONS
ILLUSTRATION 1-6

TABLE 1-3
RFI MODULE REAR CONNECTOR DESCRIPTIONS

Connector	Description	Connector	Description
J8	Body Couple A	J3	Head Output
J4	Body Output	J5	Head Couple B
J10	Body TR	J18	General Comm
J1	Amp 1 In	J19	
J2	Amp 2 In	J20	
J7	Head Couple A	J21	Power from PDU
J9	Head TR		

1-3 1.0T TELI Solid State RF Amps

See Illustration 1-7.

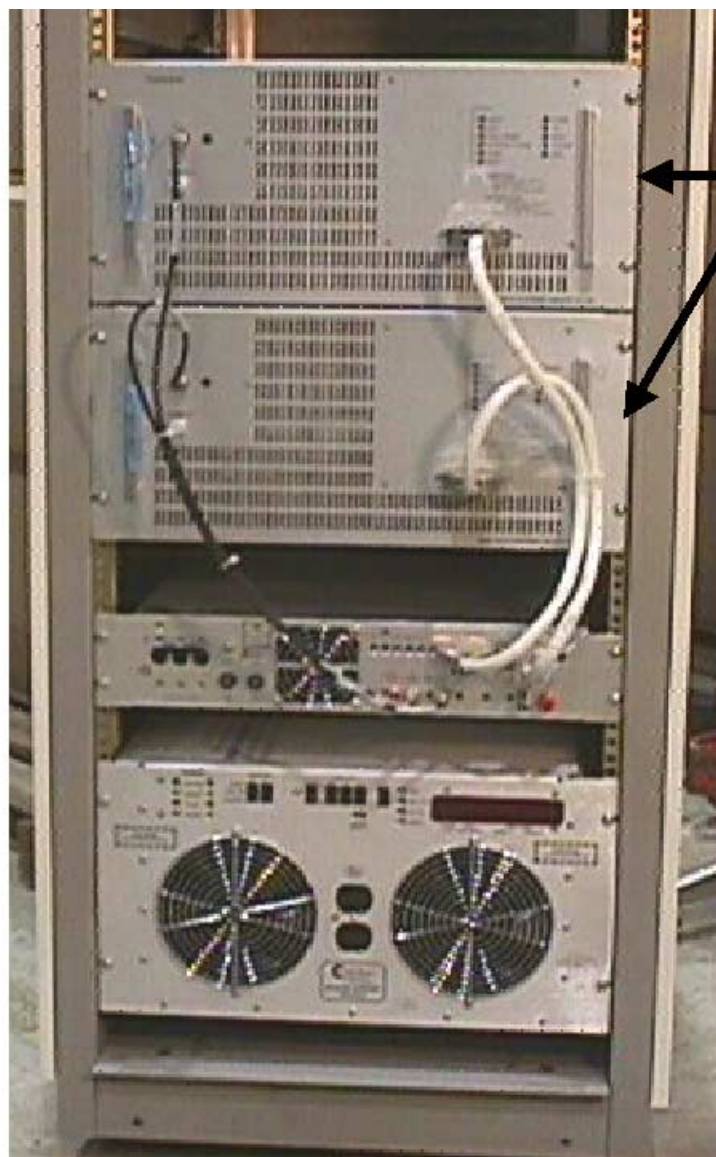


**2 (Two) 1.0T Teli RF
Amplifiers**

**RF AMP LOCATIONS
ILLUSTRATION 1-7**

1-3-1 1.0T TELI Solid State RF Amp Capacitors


See Illustration 1-8.

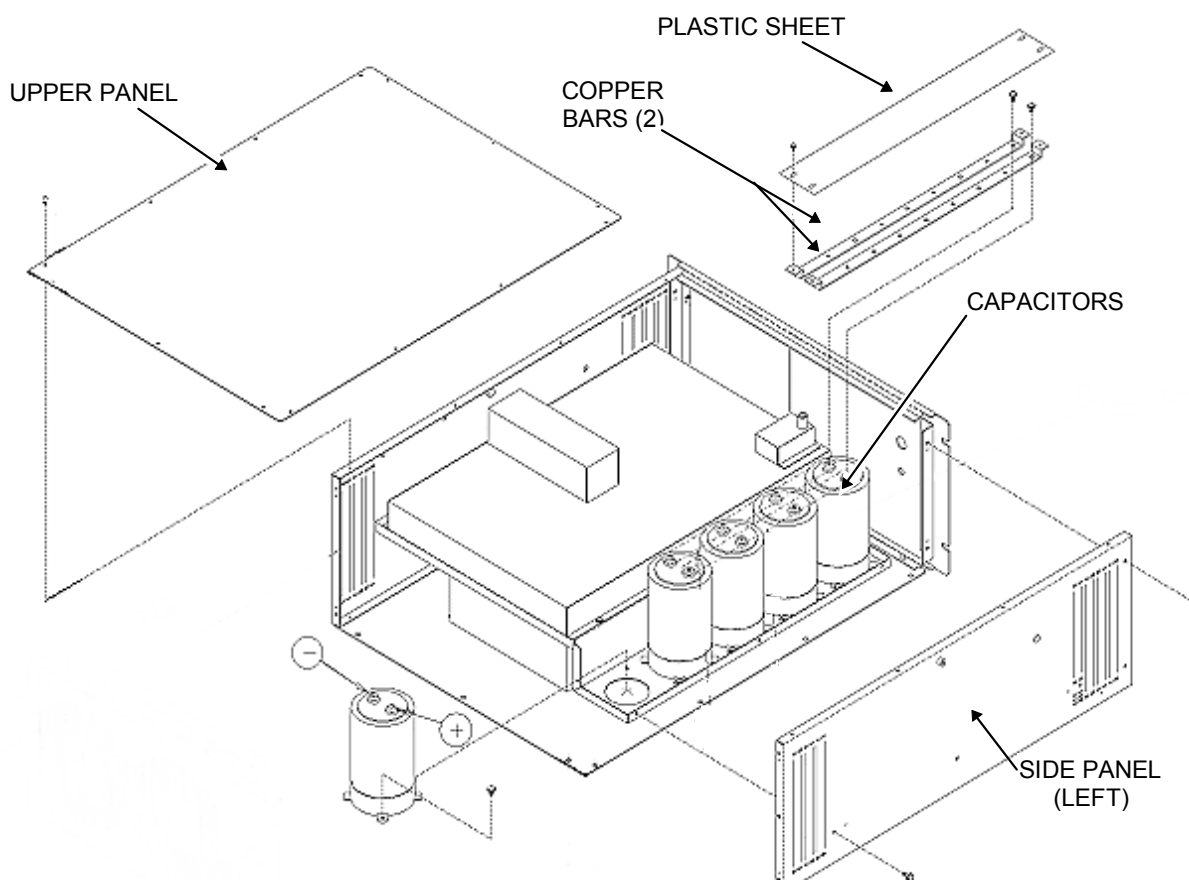

1.0T RF AMP CAPACITOR LOCATION
ILLUSTRATION 1-8

1-3-2 1.0T TELI Solid State RF Amp Fans

See Illustration 1-9.

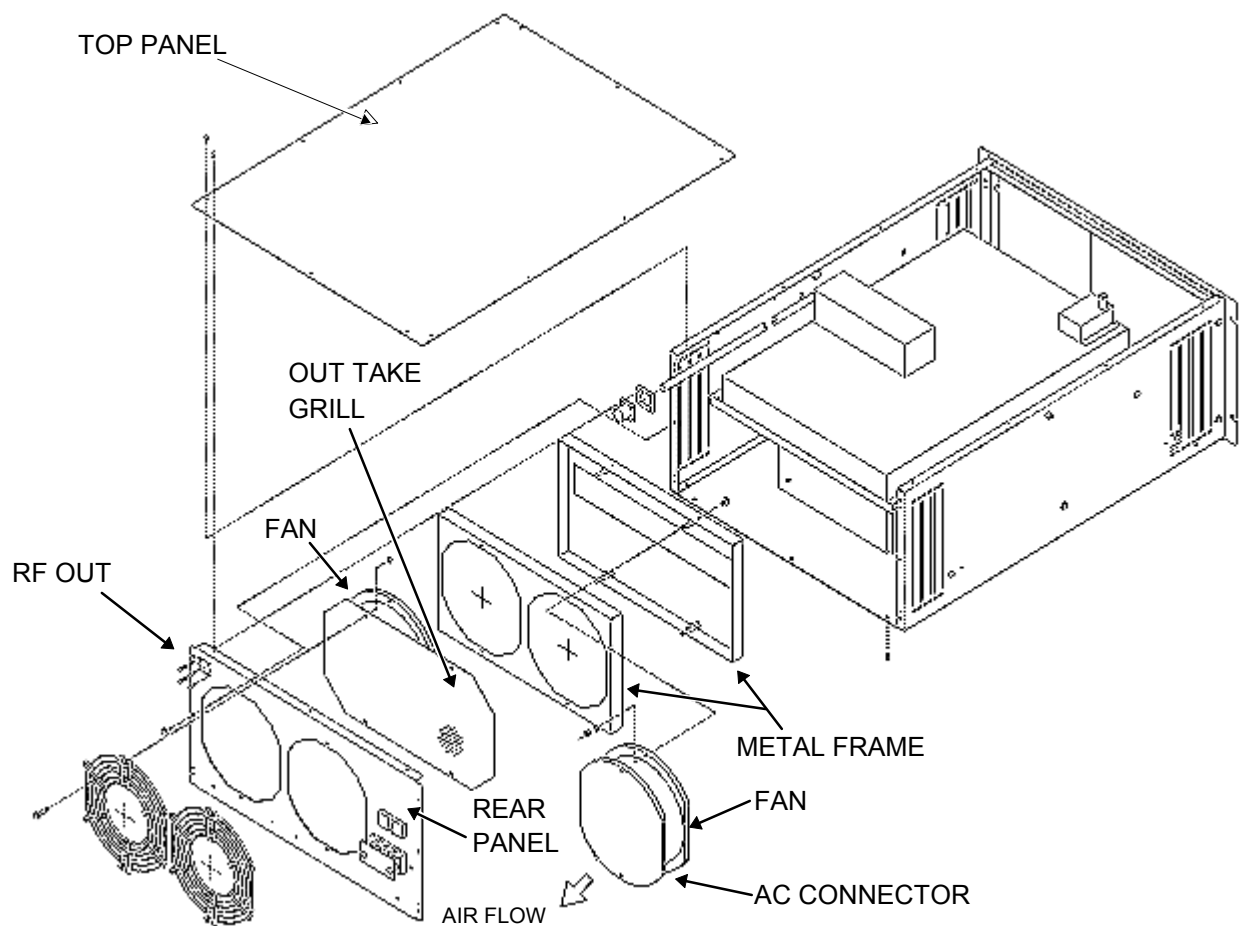


1.0T RF AMP FAN LOCATION
ILLUSTRATION 1-9

1-3-3 1.0T TELI Solid State RF Amp Power Supplies

48 Volt Power Supply

See Illustration 1-10.

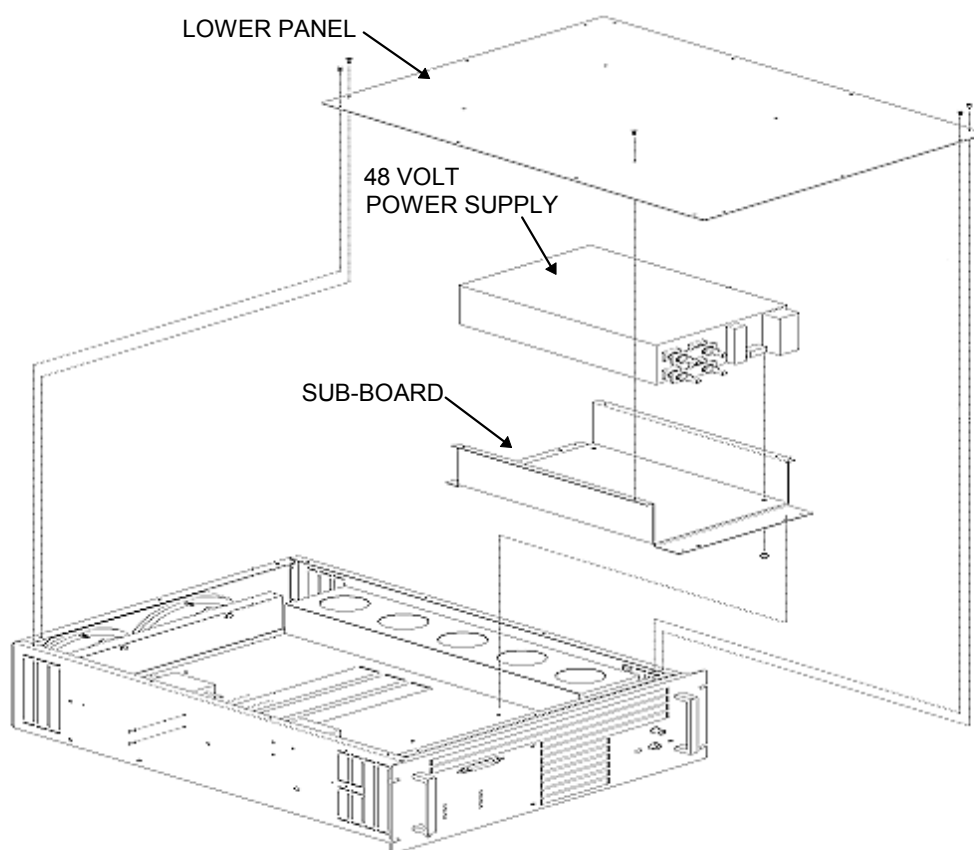


48 VOLT POWER SUPPLY LOCATION
ILLUSTRATION 1-10

1-3-3 1.0T TELI Solid State RF Amp Power Supplies (continued)

24 Volt Power Supply


See Illustration 1-11.

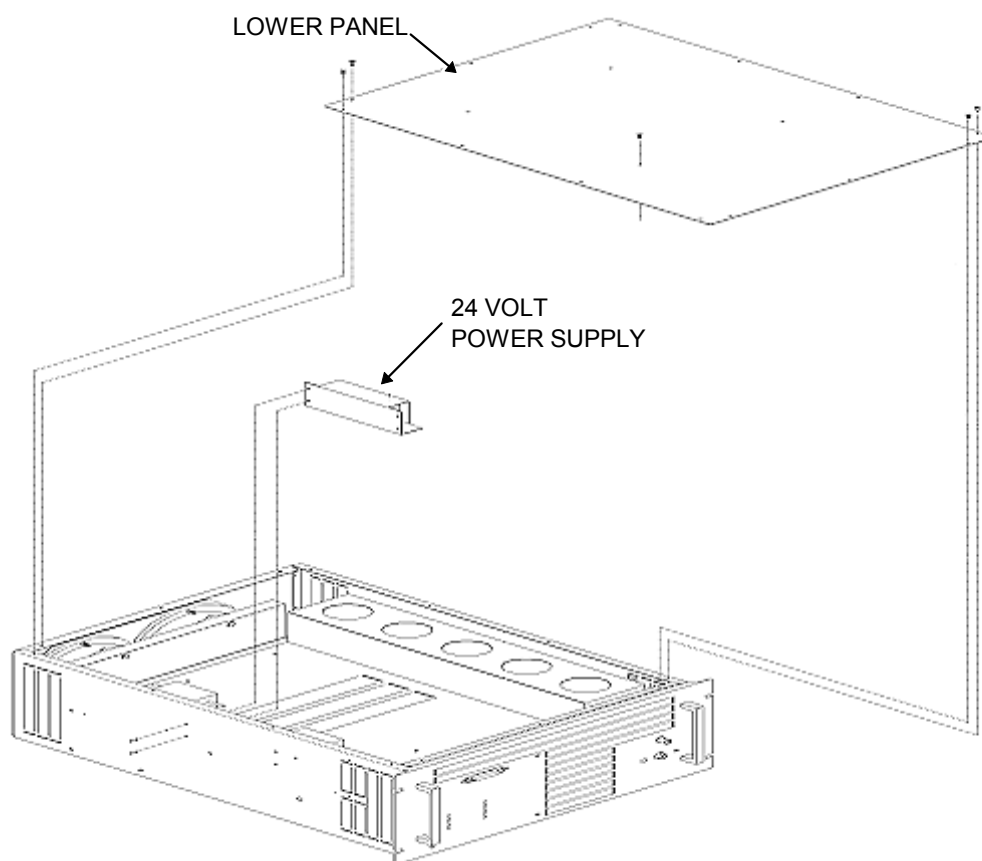

24 VOLT POWER SUPPLY LOCATION
ILLUSTRATION 1-11

1-3-4 1.0T TELI Solid State RF Amp Fuse Location

See Illustration 1-12.

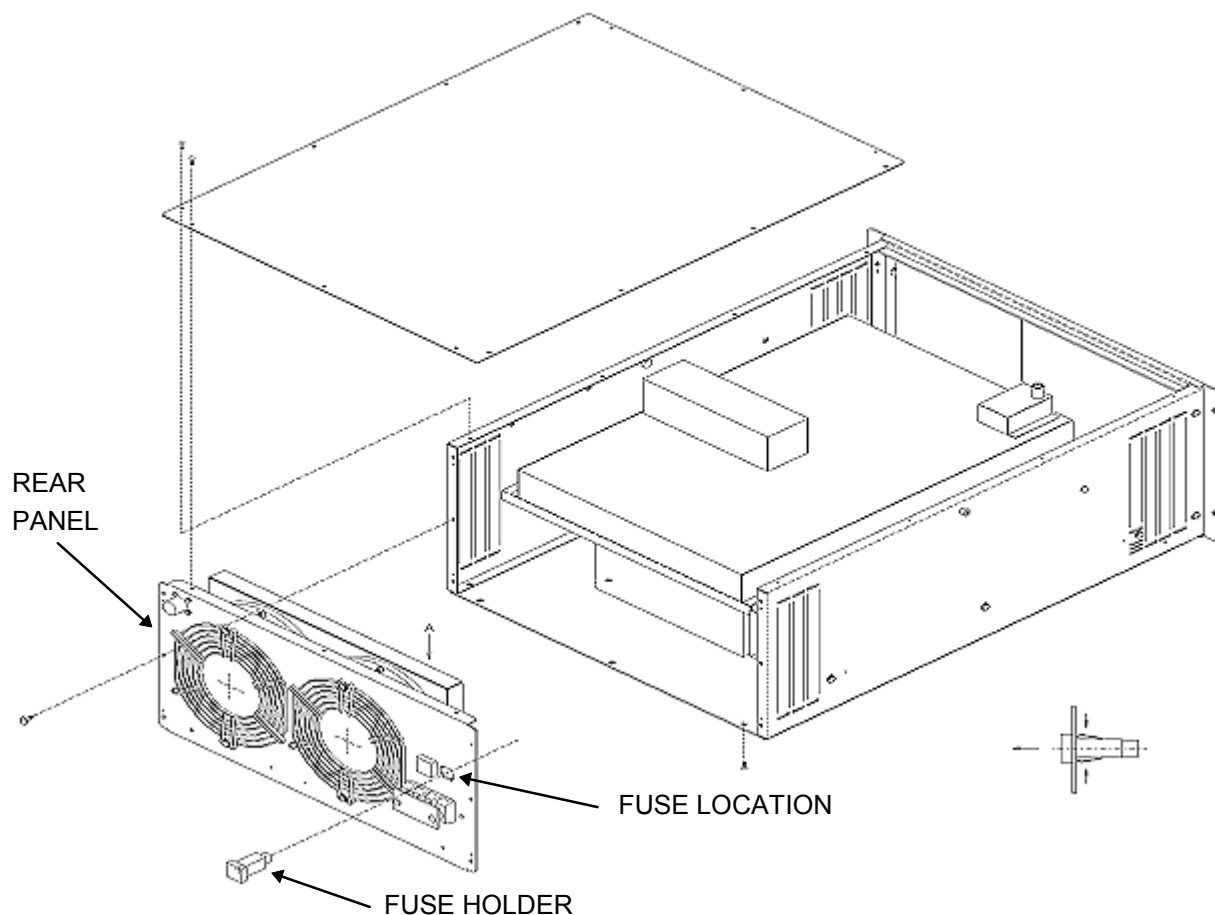

RF AMP FUSE LOCATION
ILLUSTRATION 1-12

1-3-5 1.0T TELI Solid State RF Amp Connector Descriptions

See Table 1-4.


TABLE 1-4
RF AMP CONNECTOR DESCRIPTIONS

Connector	Description
Front	
J3	RF In
J2	RF Sample
J1	Ext. Cont.
Rear	
J4	RF Out
T1	Phase L
T2	Phase N

2 - INTERCONNECTS

2-1 1.0T SRF Cabinet Interconnects


See Illustration 2-1 for rear cabinet interconnects.

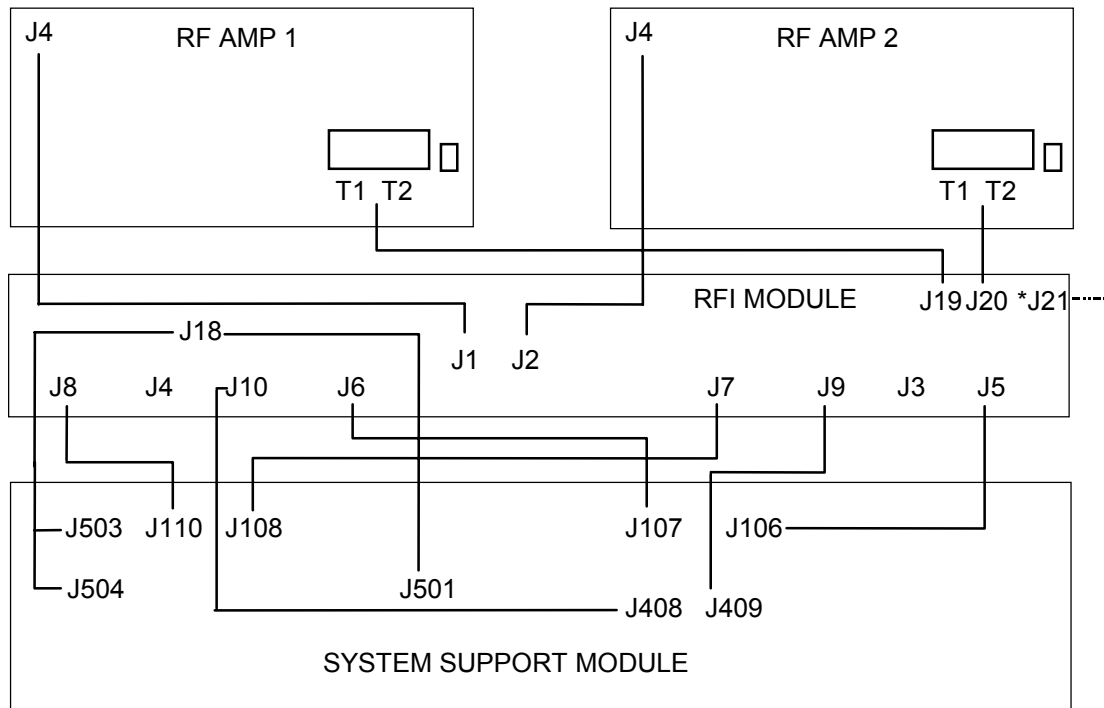

*J21 CONNECTS TO CIRCUIT BREAKER CB1 IN THE REAR OF THE SRF CABINET.

REAR CABINET RF INTERCONNECTS
ILLUSTRATION 2-1

See Illustration 2-2 for front cabinet interconnects.

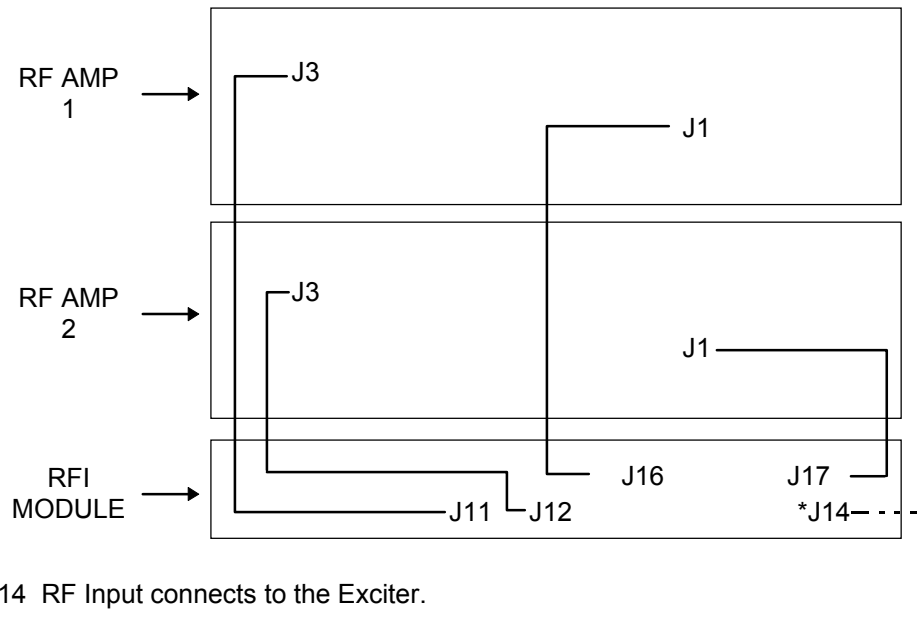


*J14 RF Input connects to the Exciter.

FRONT CABINET RF INTERCONNECTS
ILLUSTRATION 2-2

2-2 System Support Module Interconnects

See Illustration 2-3.

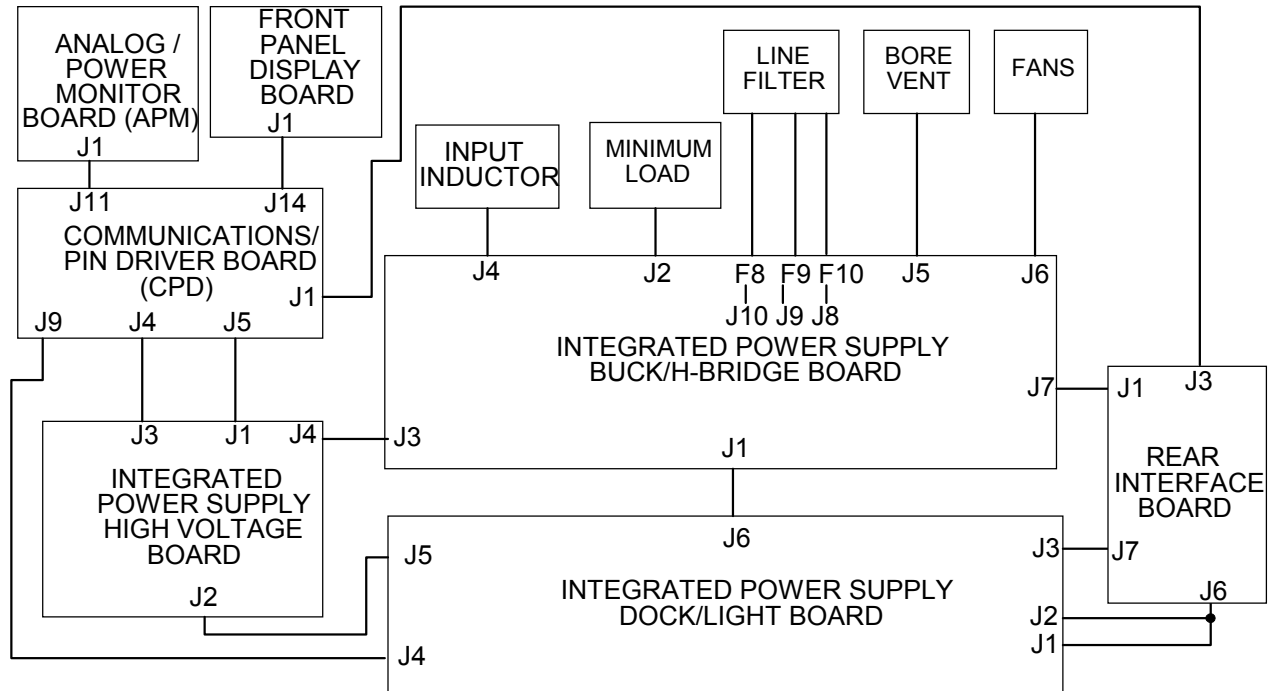

SYSTEM SUPPORT MODULE INTERCONNECTS
ILLUSTRATION 2-3

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	DEC. 20, 2000	D. Thome	Initial Release