



How To Burn Offline files

Correct Configuration VA25A, VB11D :

- Check that in Windows explorer folder CDR Offline exists under My Computer ->MED_SYSTEM (C:) ->MedCom ->temp -> CDR_OFFLINE
- CDR_OFFLINE is allowed to contain up to 2 subfolders (i.E. CDR_OFFLINE\ Sub1\ Sub2)

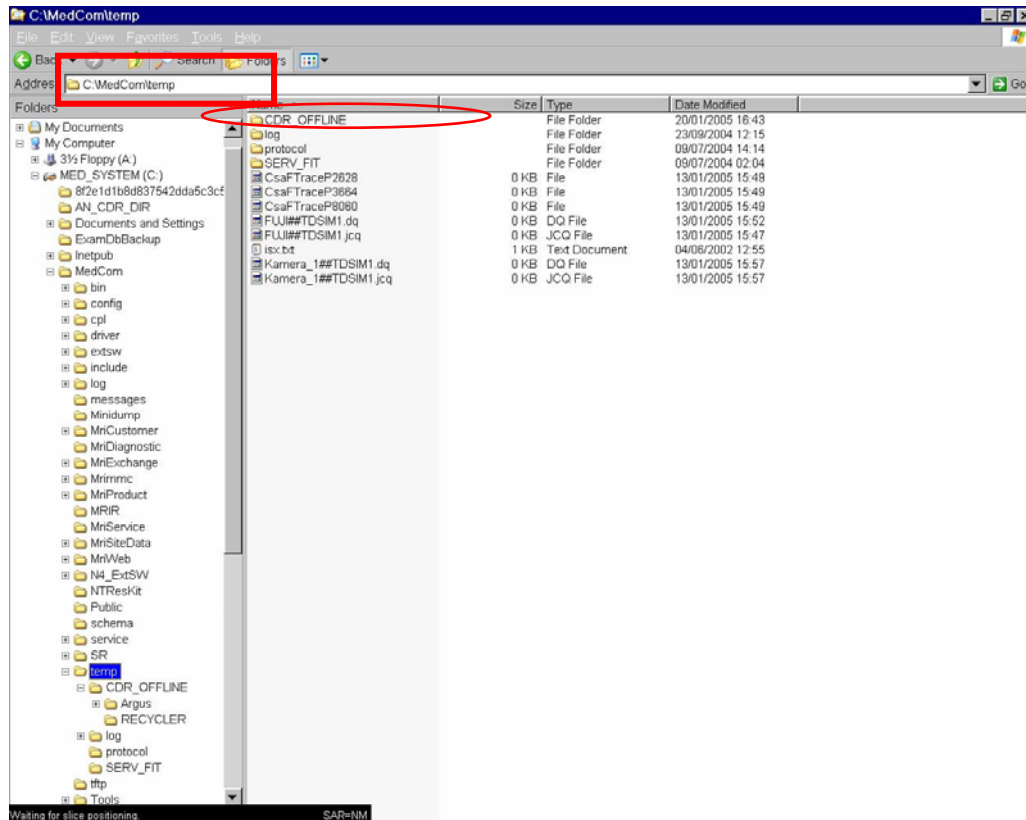


- Check that in Local Service path for Windows folder is correctly installed
- Open Local Service -> Configuration -> Image Import/Export has to be checked
-> Next -> Import/Export Directories -> Next -> Select logical HD Device
CDR_OFFLINE
- The Path shown here has to be identical with Windows Explorer Path

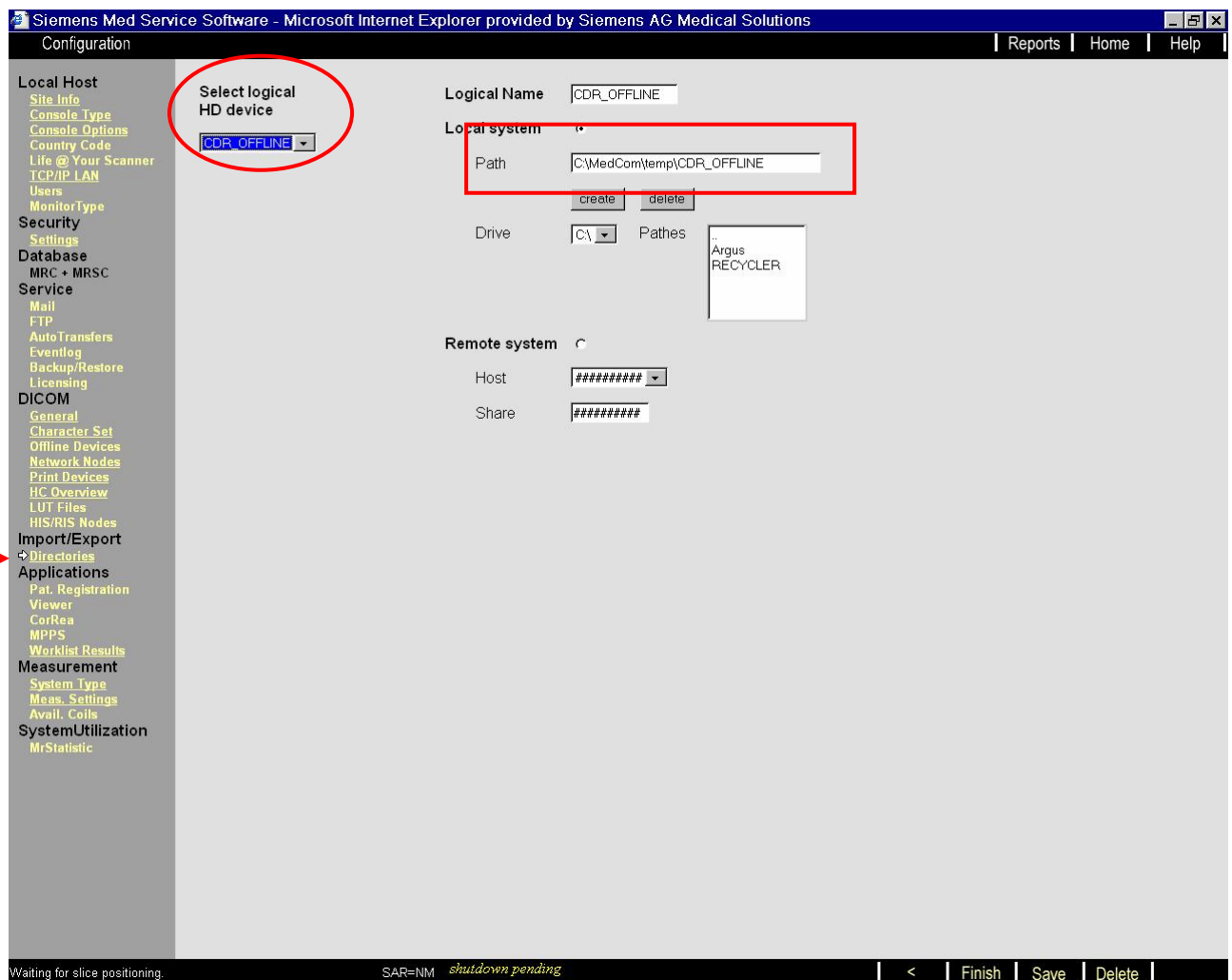

The configuration is delivered from the factory this way and it should not be changed by service or anybody else.

Correct Configuration Leonardo:

- Check that in Windows explorer folder SiteData exists under:
 - ➔ My Computer
 - ➔ ->CDR_OFFLINE (H:)
 - ➔ -> SiteData
- SiteData is allowed to contain up to 2 subfolders (i.E. SiteData\ MR\ Sub1)

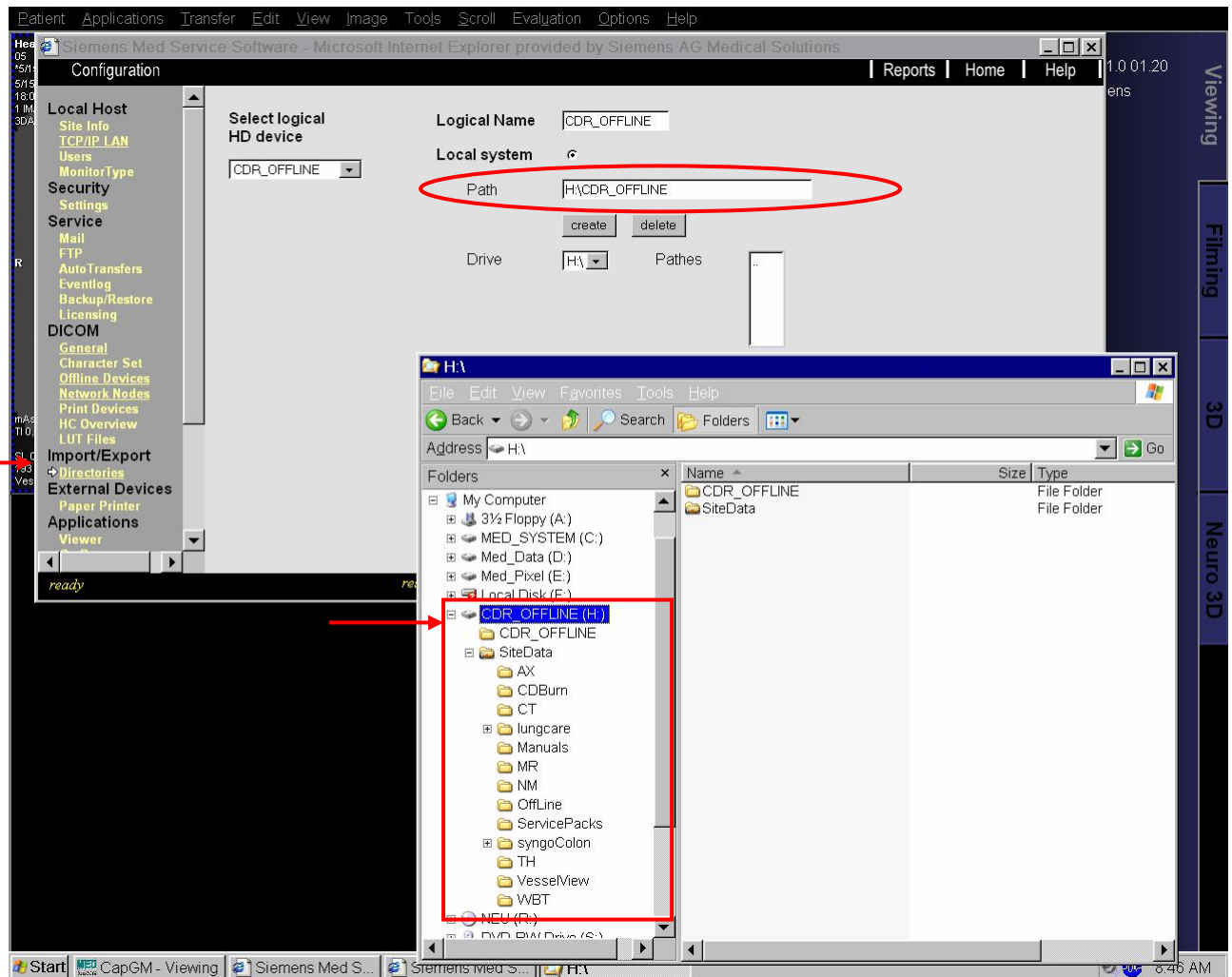


Saving data to CDR OFFLINE

- From Viewing/ Patient Browser -> select data -> Transfer -> Export to Offline
- From Argus-> Argus Tools -> Movie Control -> Save Movie As Avi
- From VesselView -> Movie card -> Camera icon


In all saving dialogs select the correct path C:\Medcom\temp\CDR_OFFLINE

In Leonardo H:\SiteData


! Please don't create any subfolders in this dialog !

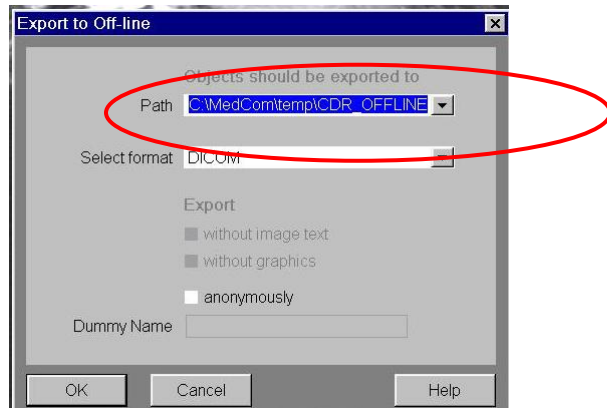

Argus:

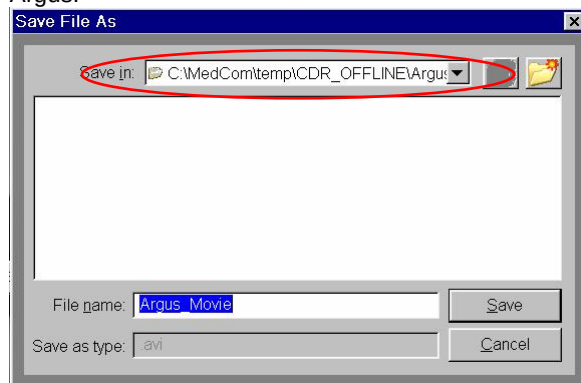


Vessel View:

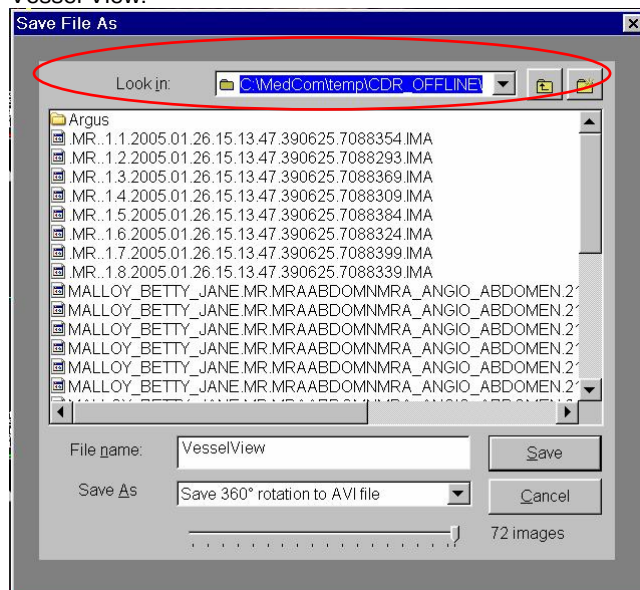

When you see text in the screen footer: "Export to off-line ... successfully completed"



you can check that the data is in Windows Explorer under
C:\Medcom\temp\CDR_OFFLINE
In Leonardo H:\SiteData

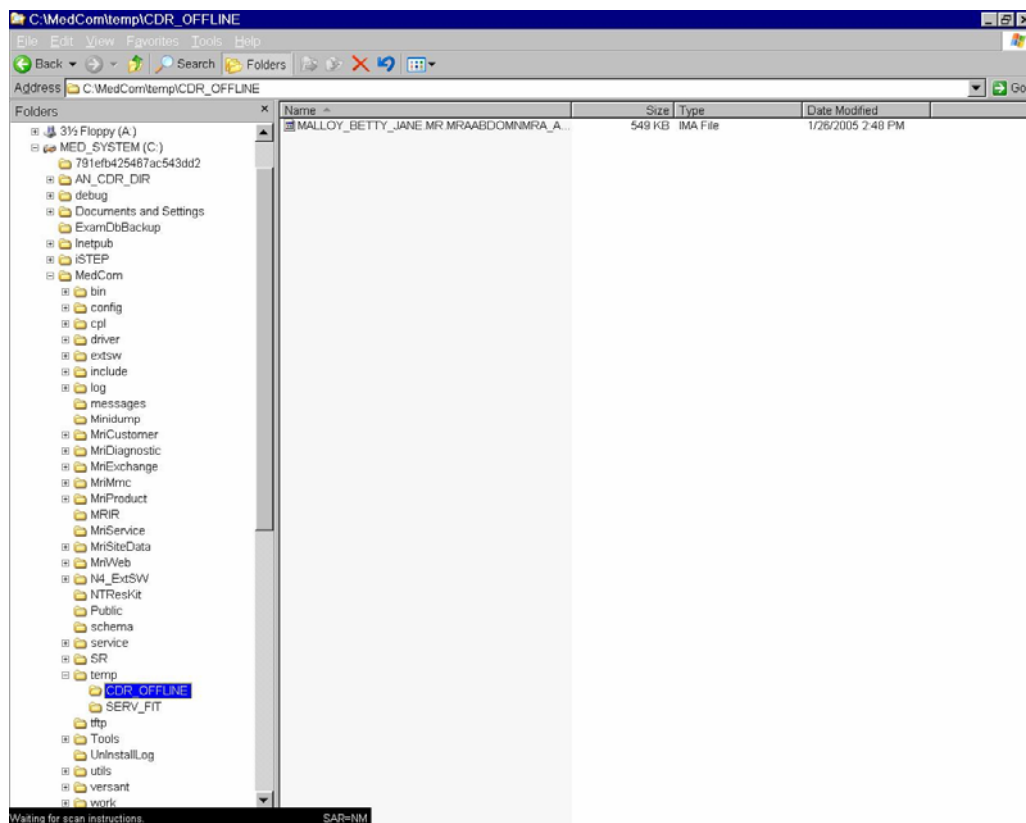

Clean Local Job Status

Local Job status has to be ! *completely empty* ! before burning Offline data is possible

- Open Transfer -> Local Job Status
- Delete all jobs
- If some jobs active, wait until they are completed and delete the completed jobs



Insert CD into system *before* selecting record command


Insert a new, empty CD. Multisession can not be used with Offline files, therefore the inserted CD has to be empty and cannot be used again after one burning session has been closed

Wait until the system has recognized the new CD, following text appears in screen footer:

Medium in CD-R inserted and contents will be read.

Select Record Offline files

Transfer -> Record Offline files

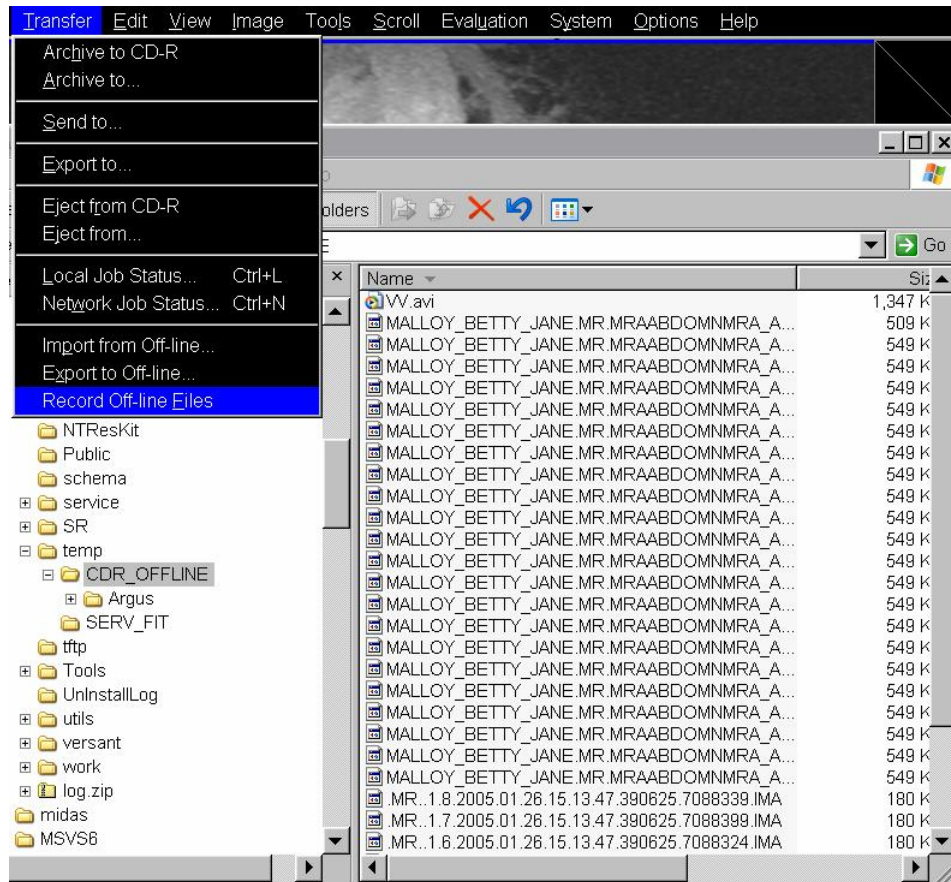

! Note the system doesn't give any information, that it has started burning !
In PC by the CD burner a yellow light starts blinking, this means that the burning has started.

When burning has finished the CD is automatically ejected and following text appears to screen footer:

Recording successfull.

To be successful with burning Offline files, please follow the procedure order described in this manual.

If you should have any problems succeeding, though following the description, please send us an e-mail to ApplikationMR@Siemens.com with complete savelog and problem description.