

Chapter 7

GOC2 (w/Linux PC)

Note: LightSpeed 3.X CT scanners may use any one of five Operator's Consoles. These include:

- **H3 Console** - See Chapter 5.
- **GOC1 (w/ Octane 2)** - See Chapter 6.
- **GOC2 (w/ Linux PC)** - See *This* Chapter.
- **GOC3 or GOC4** - See the *GOC3 IPL*, 2394950-100, or the *GOC4 IPL*, 5134893-100, respectively.


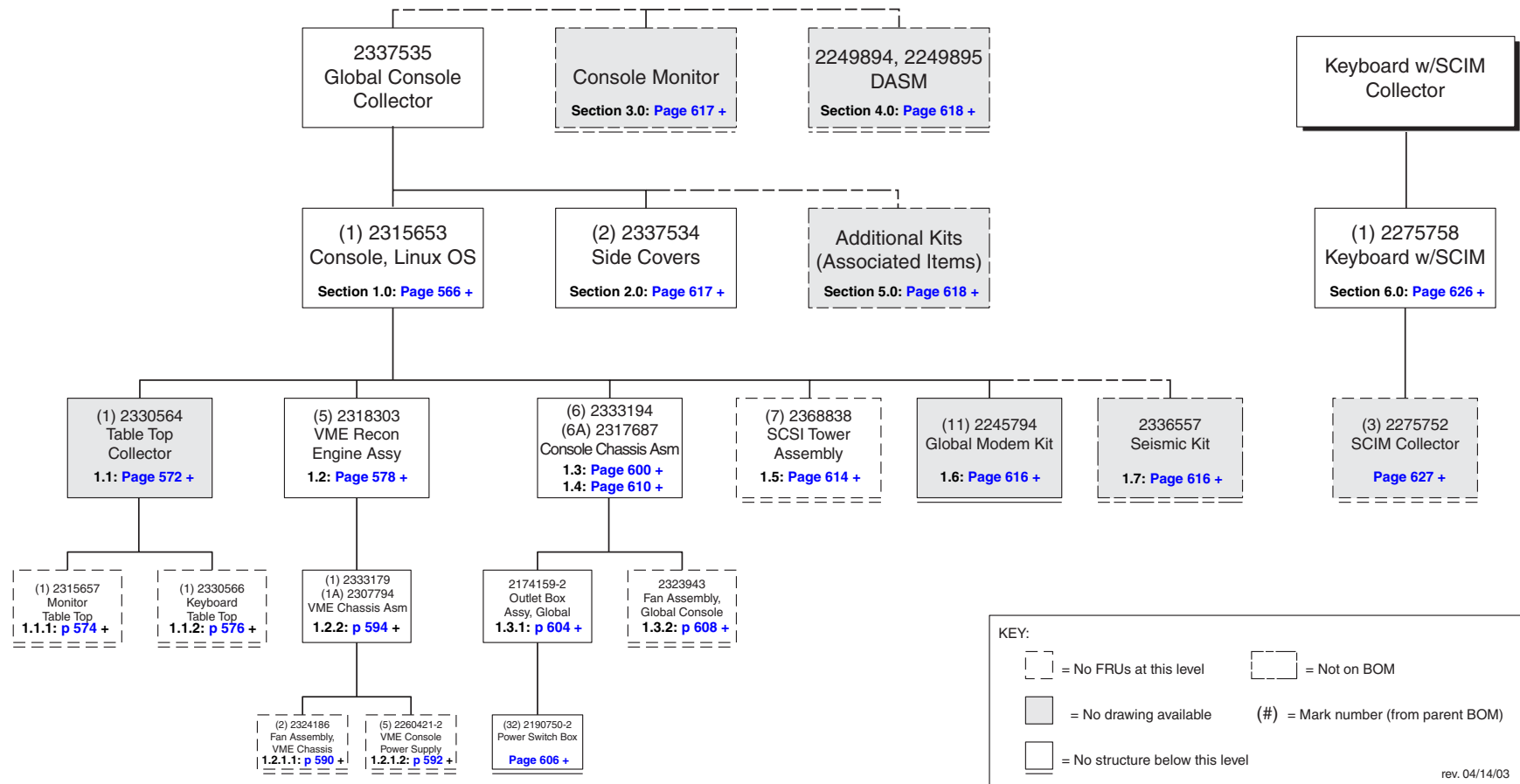


Figure 7-1 Global Console - Linux (GOC2) BOM Tree

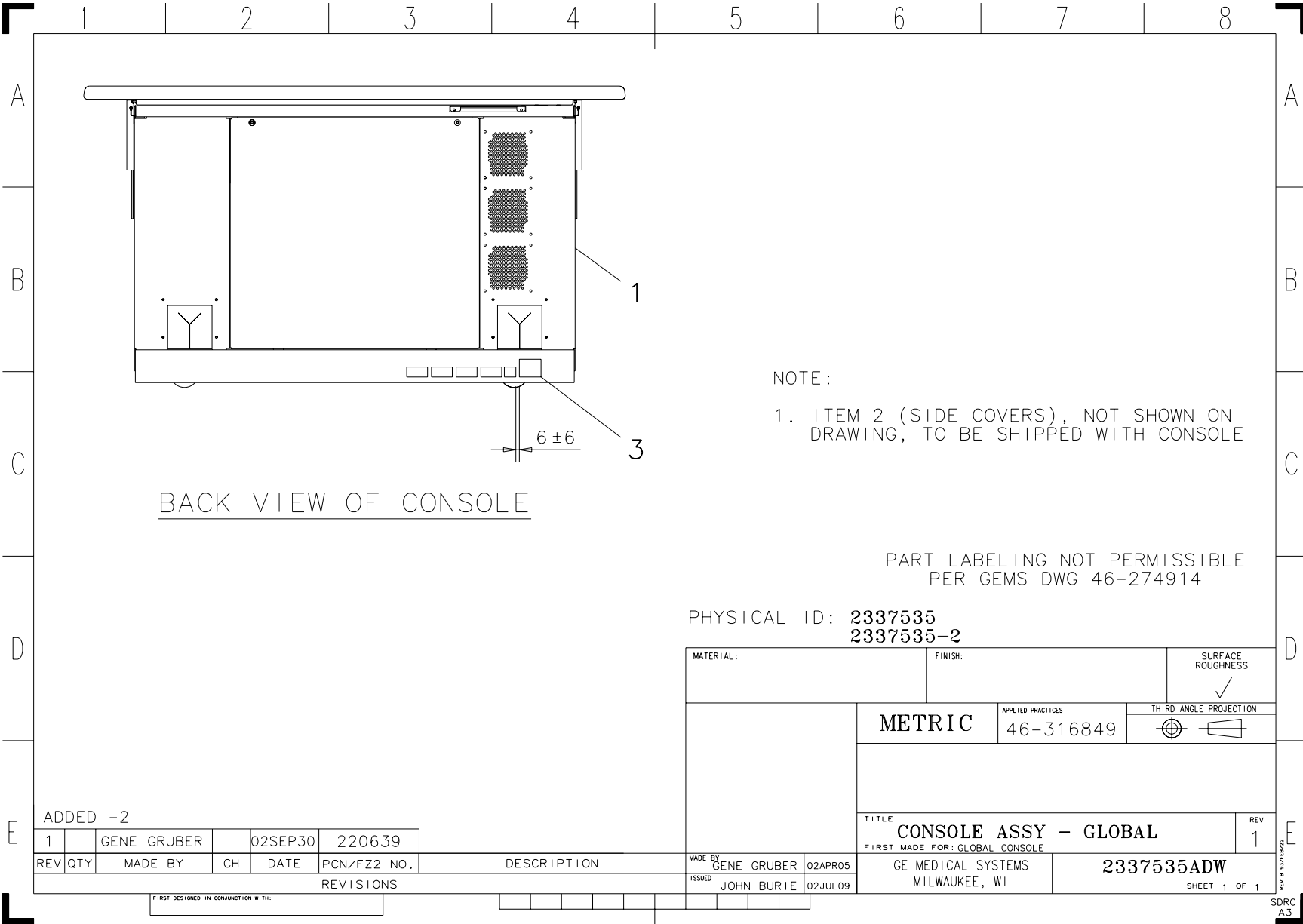

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2315653	1.000	Yes	CONSOLE-GOC, LINUX OPER. {For further detail, see Section 1.0 , beginning on page 566 .}
2	2337534	1.000	Yes	SIDE COVERS GOC {For further detail, see Section 2.0 , beginning on page 617 .}
3	46-302200P1	1.000	No	RATING PLATE W/ SN, MN_&

Table 7-1 Global Console Collector (2337535 BOM, rev 1)

Section 1.0 — Console Assembly (2315653)

FRU = YES

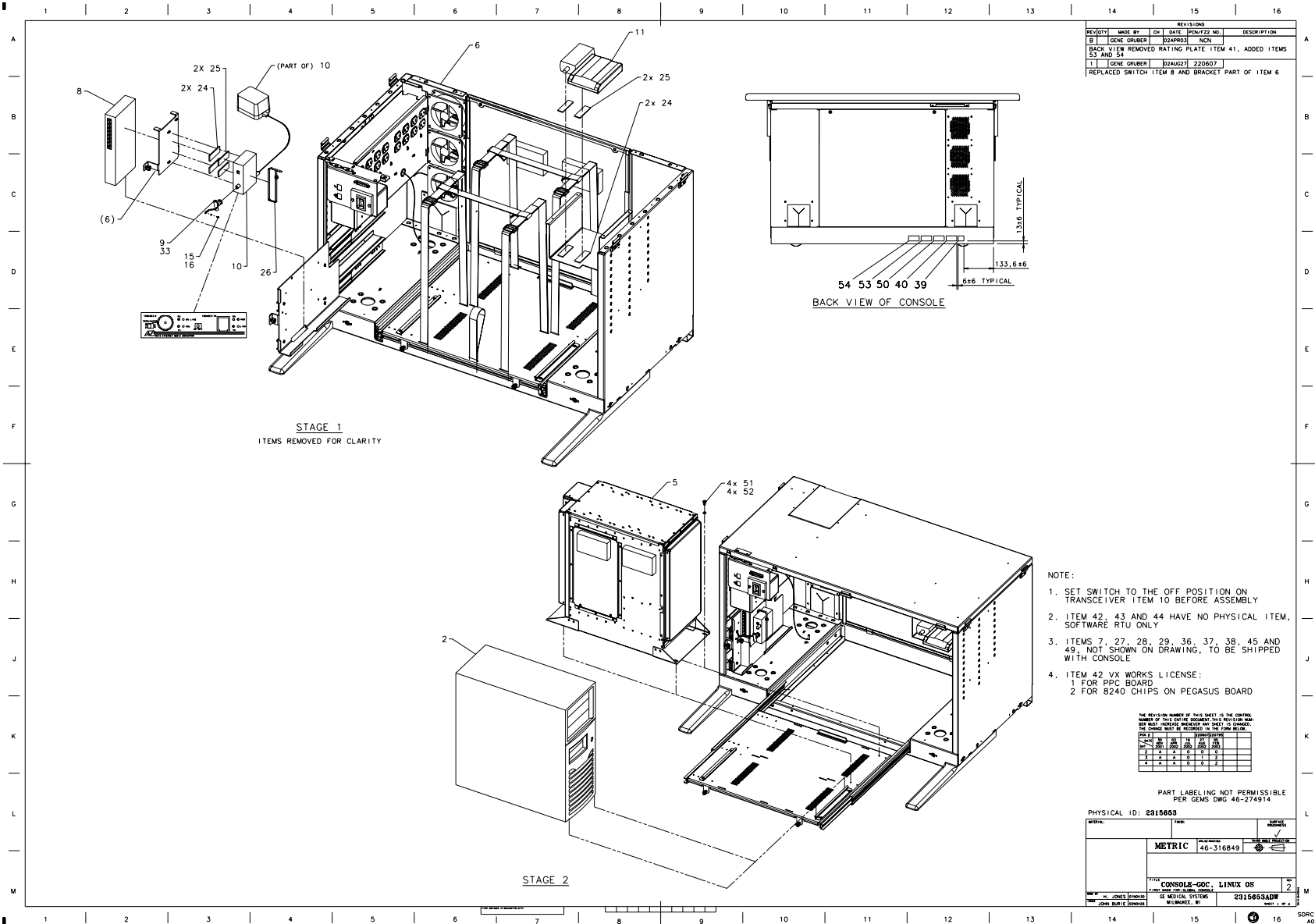


Figure 7-3 Console Assembly (2315653, sheet 1, rev 2)

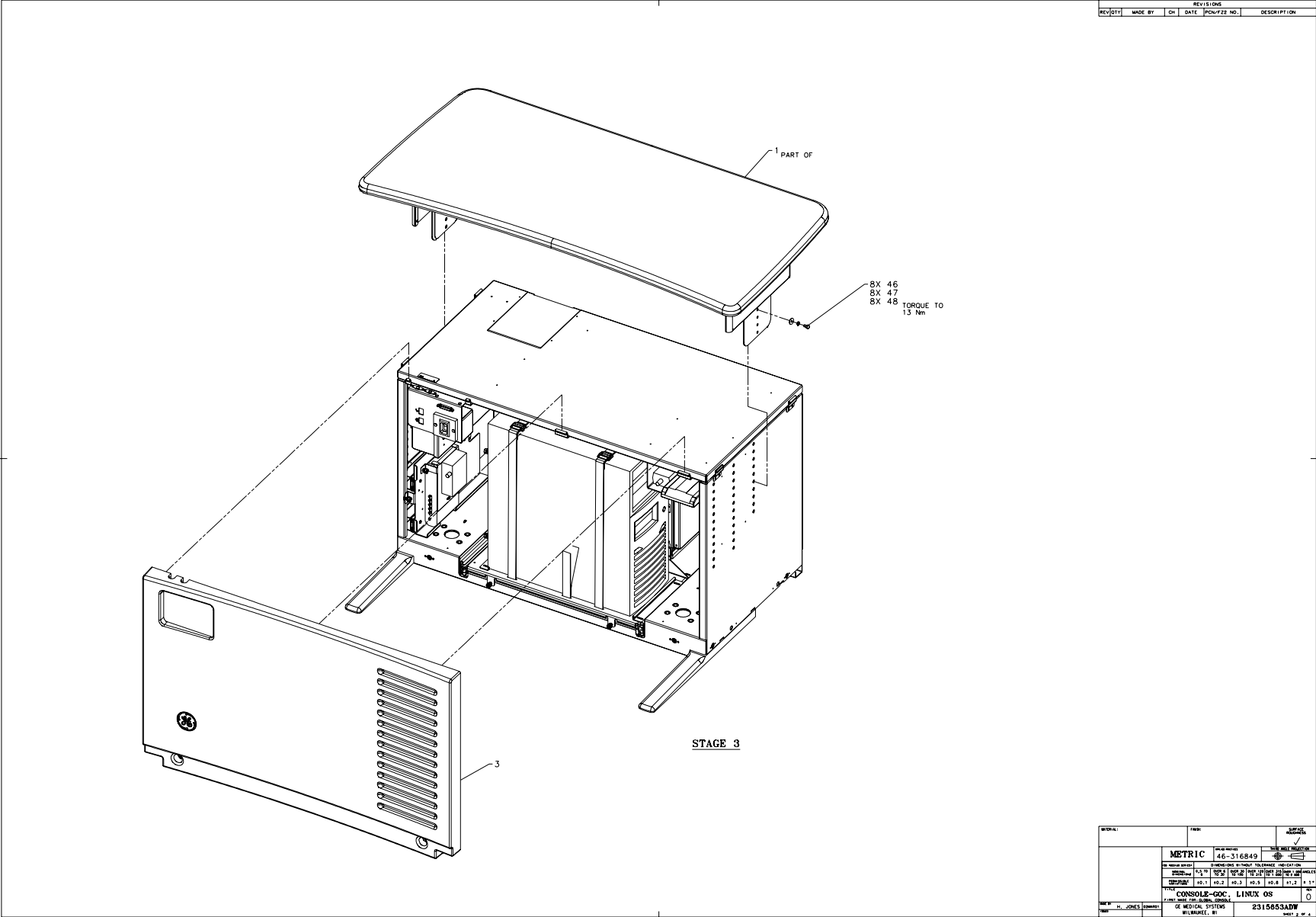


Figure 7-4 Console Assembly (2315653, sheet 2, rev 0)



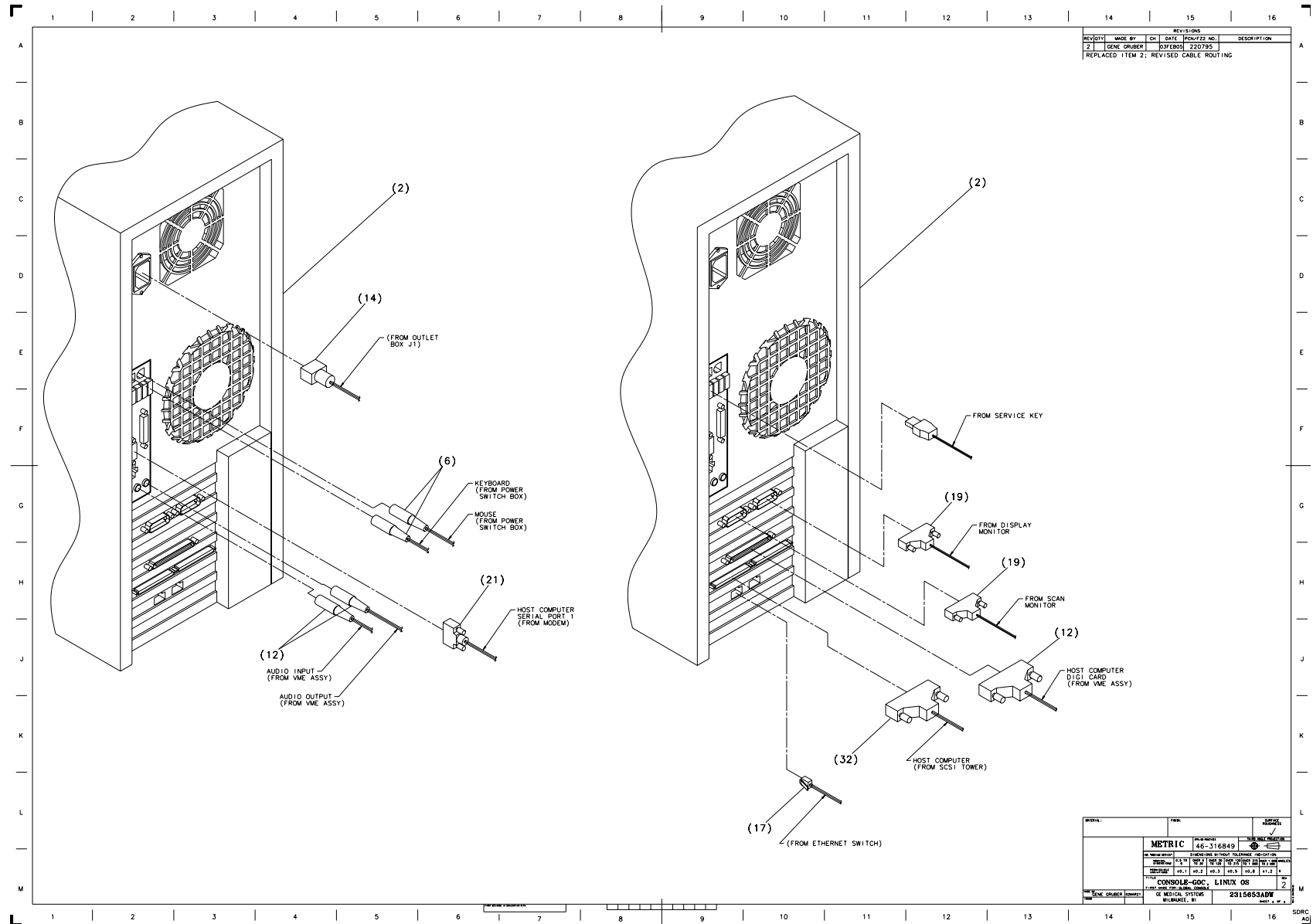


Figure 7-6 Console Assembly (2315653, sheet 4, rev 2)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
2	2368697	1.000	Yes	LINUX DDR WORKSTATION - G
3	2315655	1.000	Yes	FRONT COVER, CONSOLE CHAS
5	2318303	1.000	Yes	VME RECON ENGINE ASSEMBLY {For further detail, see Section 1.2 , beginning on page 578 .}
7	2368838	1.000	Yes	SCSI TOWER, DVD-RAM, MOD {For further detail, see Section 1.5 , beginning on page 614 .}
8	2197229-5	1.000	Yes	ETHERNET SWITCH, 8 PORT
9	2154019	1.000	Yes	50 OHM BNC GROUNDED TERMI
10	2148559	1.000	Yes	10BASE-2 TO UTP CONVERTER
11	2245794	1.000	Yes	GLOBAL MODEM KIT {For further detail, see Section 1.6 , beginning on page 616 .}
27	2336906	1.000	Yes	TRACKBALL FOR LINUX CONSO
28	2368697-2	2.000	Yes	LINUX DDR WORKST'N 3BUTTO
28A	2330676-2	ALT	Yes	COMMON WORKST'N 3-BUTTON
1	2330564	1.000	No	CONSOLE TABLE TOP COLLECTOR {For further detail, see Section 1.1 , beginning on page 572 .}
6	2333194	1.000	No	ASSY-CONSOLE CHASSIS {For further detail, see Section 1.3 , beginning on page 600 .}
6A	2317687	ALT	No	CONSOLE CHASSIS ASSEMBLY {For further detail, see Section 1.4 , beginning on page 610 .}
12	2317488	1.000	No	OCTOPUS CABLE, LINUX COMP
13	2142225	3.000	No	MONITOR POWER CORD
14	2142226	3.000	No	POWER CORD FOR EMC AREA
15	46-328421P4	1.000	No	SCREW PAN HEAD 2.5 MM 6 M
16	46-328432P1	1.000	No	WASHER LOCK - SPRING 2.6
17	2318311	2.000	No	ETHERNET CABLE, RJ45 TO R
18	2323942	1.000	No	CABLE, HD50 TO HD50
19	2304737-2	2.000	No	CABLE, DVI VIDEO CABLE -
21	2323945	1.000	No	CABLE - MODEM TO HOST COM

Table 7-2 Console Assembly (2315653 BOM, rev 4)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
24	46-136246P14	1.000 ft	No	2" WIDE, WHT, ADH BACK, V
25	46-136246P21	1.000 ft	No	2" WIDE X 25 YD WHITE VEL
26	46-208758P9	3.000	No	11.08"X .140" SELF-LOCKIN
29	2227726	1.000	No	MOUSE PAD FOR LIGHTSPEED
30	2318712	1.000	No	VME ETHERNET CABLE, LINUX
33	46-297332P1	1.000	No	50 OHM, BNC TEE ADAP
34	46-208781P2	7.000	No	WRITE-ON LABEL FOR MARKIN
35	2324697-2	1.000	No	SCSI CABLE, VHD68 TO CENT
36	2335376	1.000	No	CABLE - PDU TO CONSOLE FO
39	2132370	1.000	No	CSA LABEL
40	2117390	1.000	No	LABEL, CE MARKING COMPLIA
42	2284884	1.000	No	VXWORKS 5.4 LICENCE
43	46-269568P1	1.000	No	INFORMIX SOFTWARE-DATA BA
44	2240540	1.000	No	NETSCAPE LICENSE
45	2214164	1.000	No	5.25" 2.3GB 512 BYTES/SEC
46	2337083-2	8.000	No	SCREW, M6 X 20 SOC HEAD,
47	46-328430P3	8.000	No	WASHER PLAIN - NORMAL 6.4
48	46-328432P3	8.000	No	WASHER LOCK - SPRING 6.1
49	2337414	1.000	No	DVD RAM MEDIA, 9.4 GB
50	2132620	1.000	No	CE MARKING LABEL CLASS A
51	2103585-2	4.000	No	SCREW PAN HEAD 4 MM 10 MM
52	46-328431P1	4.000	No	WASHER PLAIN - LARGE 4.3
53	46-279769P1	1.000	No	IEC TYPE B SYMBOL
54	46-269337P6	1.000	No	LABEL OLD
55	46-208758P3	16.000	No	7.31" X .184" SELF-LOCKIN
56	46-208758P2	3.000	No	5.5" X .140" SELF-LOCKING

Table 7-2 Console Assembly (2315653 BOM, rev 4) (Continued)

1.1 Console Table Top Collector (2330564)

FRU = NO

USED ON:
2315653

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2315657	1.000	Yes	TABLE TOP ASSEMBLY, CONSOLE CHASSIS {For further detail, see Section 1.1.1 , beginning on page 574 .}
2	2330566	1.000	Yes	KEYBOARD TABLE TOP {For further detail, see Section 1.1.2 , beginning on page 576 .}

Table 7-3 Console Table Top Collector (2330564 BOM, rev 0)

This Page Intentionally Blank.

1.1.1 Table Top Assembly, Monitor (2315657)

FRU = YES


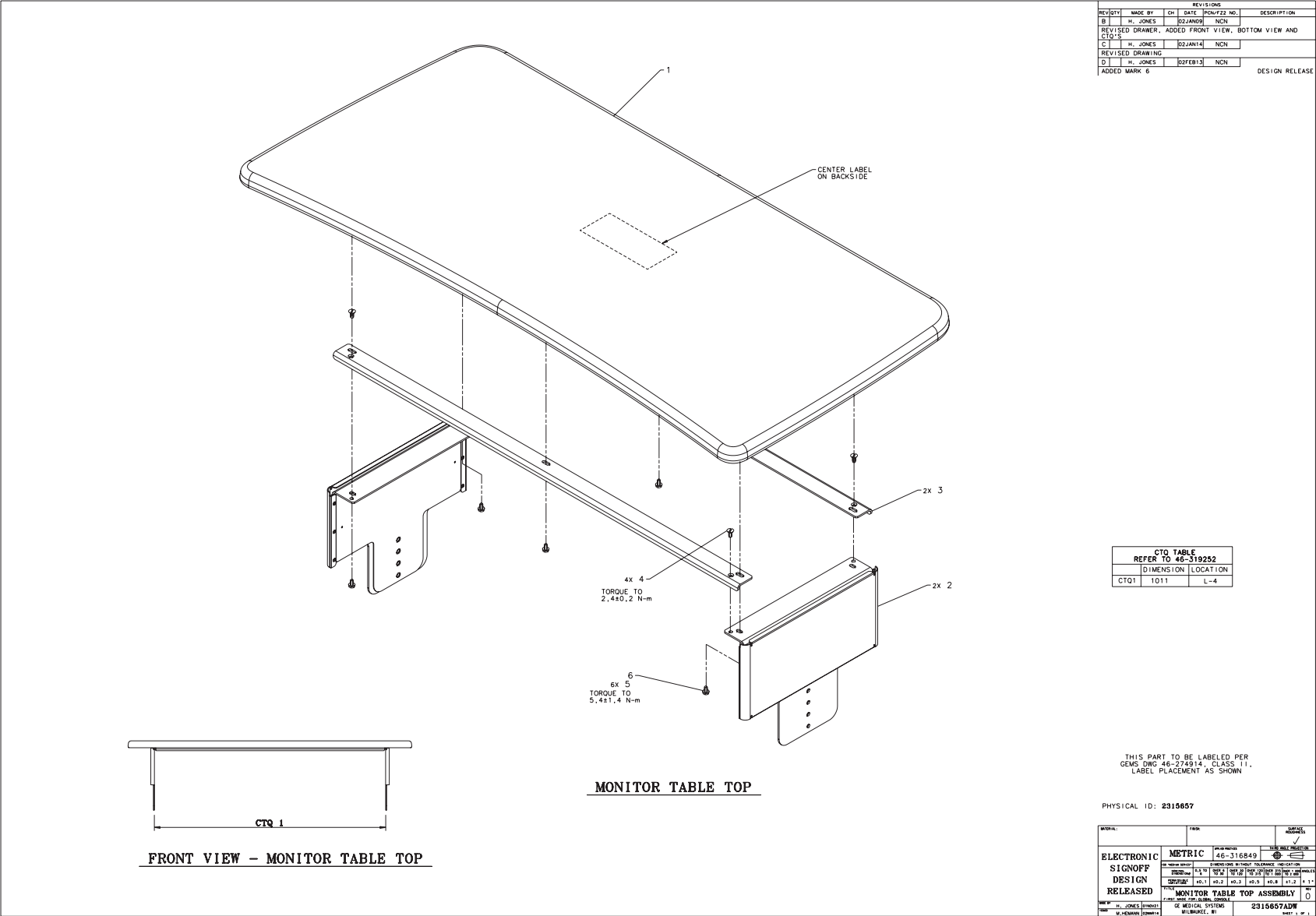


Figure 7-7 Table Top Assembly, Monitor (2315657, sheet 1, rev 0)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2317496	1.000	No	MONITOR TABLE TOP
2	2317588	2.000	No	ADJUSTMENT BRKT, MONITOR
3	2317503	2.000	No	SUPPORT BRACKET, MONITOR
4	46-328422P501	4.000	No	SCREW COUNTERSUNK HEA 5 M
5	46-170015P23	6.000	No	SCREW
6	46-170686P3	AR	No	LOCTITE 242, 50CC

Table 7-4 Table Top Assembly, Monitor (2315657 BOM, rev 0)

1.1.2 Table Top Assembly, Keyboard (2330566)
FRU = YES

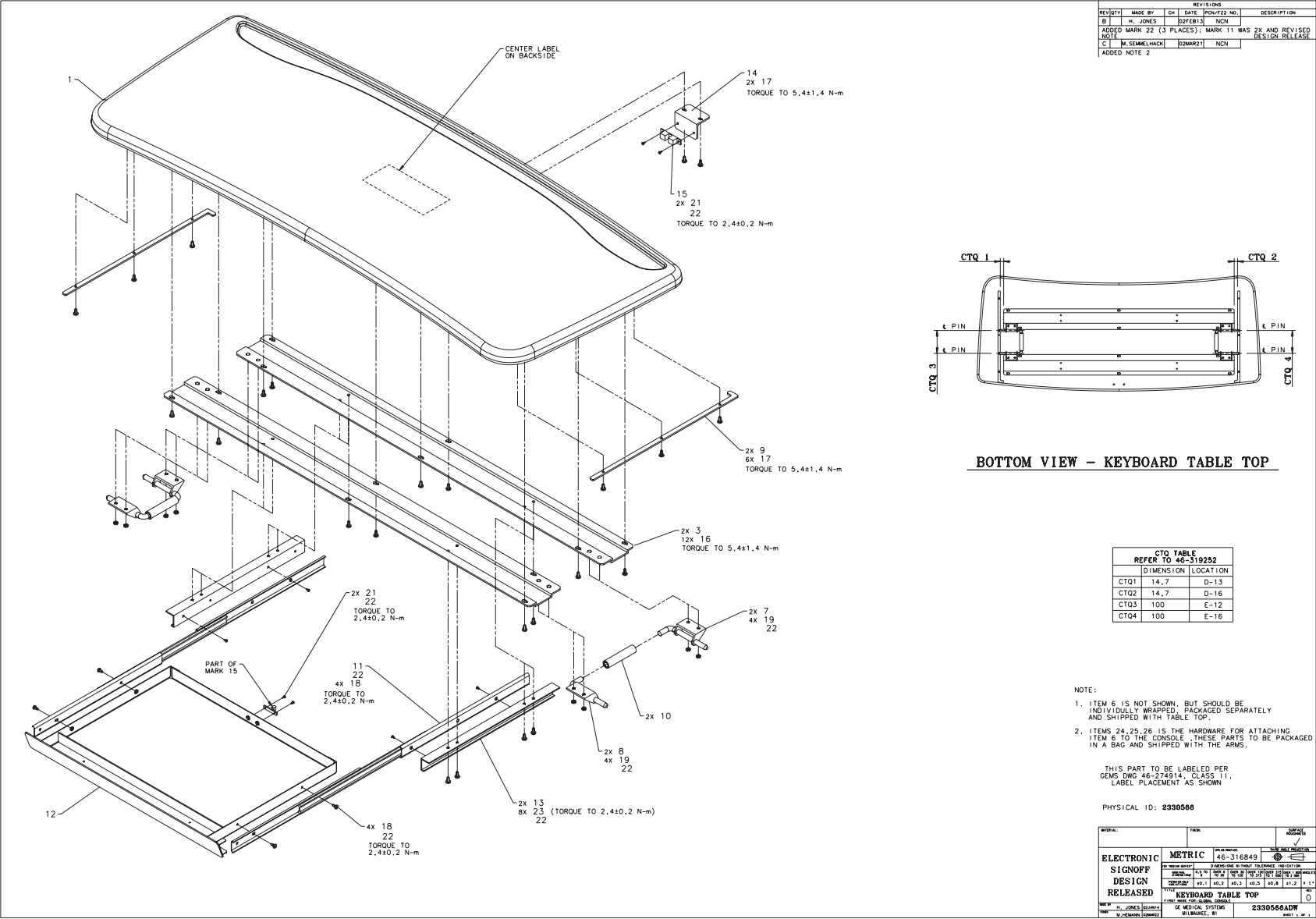

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2317494	1.000	No	KEYBOARD TABLE TOP
3	2317512	2.000	No	SUPPORT BRACKET, KEYBD TA
6	2317498	2.000	No	KEYBOARD ADJUSTMENT BRACK
7	2317593	2.000	No	LATCH HINGE, LEFT HANDED
8	2317594	2.000	No	LATCH HINGE, RIGHT HANDED
9	2317497	2.000	No	ALIGNMENT BRACKET
10	2317500	2.000	No	TUBE, KEYBD ALIGNMENT DEV
11	2317595	1.000	No	SLIDE RAIL FOR DRAWER ASS
12	2317499	1.000	No	DRAWER
13	2317513	2.000	No	DRAWER BRACKET
14	2317501	1.000	No	BRACKET, DRAWER BALL LATC
15	2317502	1.000	No	LATCH, DRAWER ASSEMBLY
16	46-170015P23	12.000	No	SCREW
17	46-170015P14	8.000	No	PAN HD SCREW 8-32 X 3/8
18	46-328421P11	8.000	No	SCREW PAN HEAD 4 MM 5 MM
19	46-312346P2	8.000	No	NUT, M4 X 0.7 HEX WITH LO
21	46-328421P3	4.000	No	SCREW PAN HEAD 2.5 MM 5 M
22	46-170686P3	AR	No	LOCTITE 242, 50CC
23	2103585	8.000	No	SCREW PAN HEAD 4 MM 8 MM
24	2337083	13.000	No	SCREW, M6 X 35 SOC HEAD
25	46-328430P3	13.000	No	WASHER PLAIN - NORMAL 6.4
26	46-328432P3	13.000	No	WASHER LOCK - SPRING 6.1

Table 7-5 Table Top Assembly, Keyboard (2330566 BOM, Rev 0)

1.2 VME Recon Engine Assembly (2318303)

FRU = YES

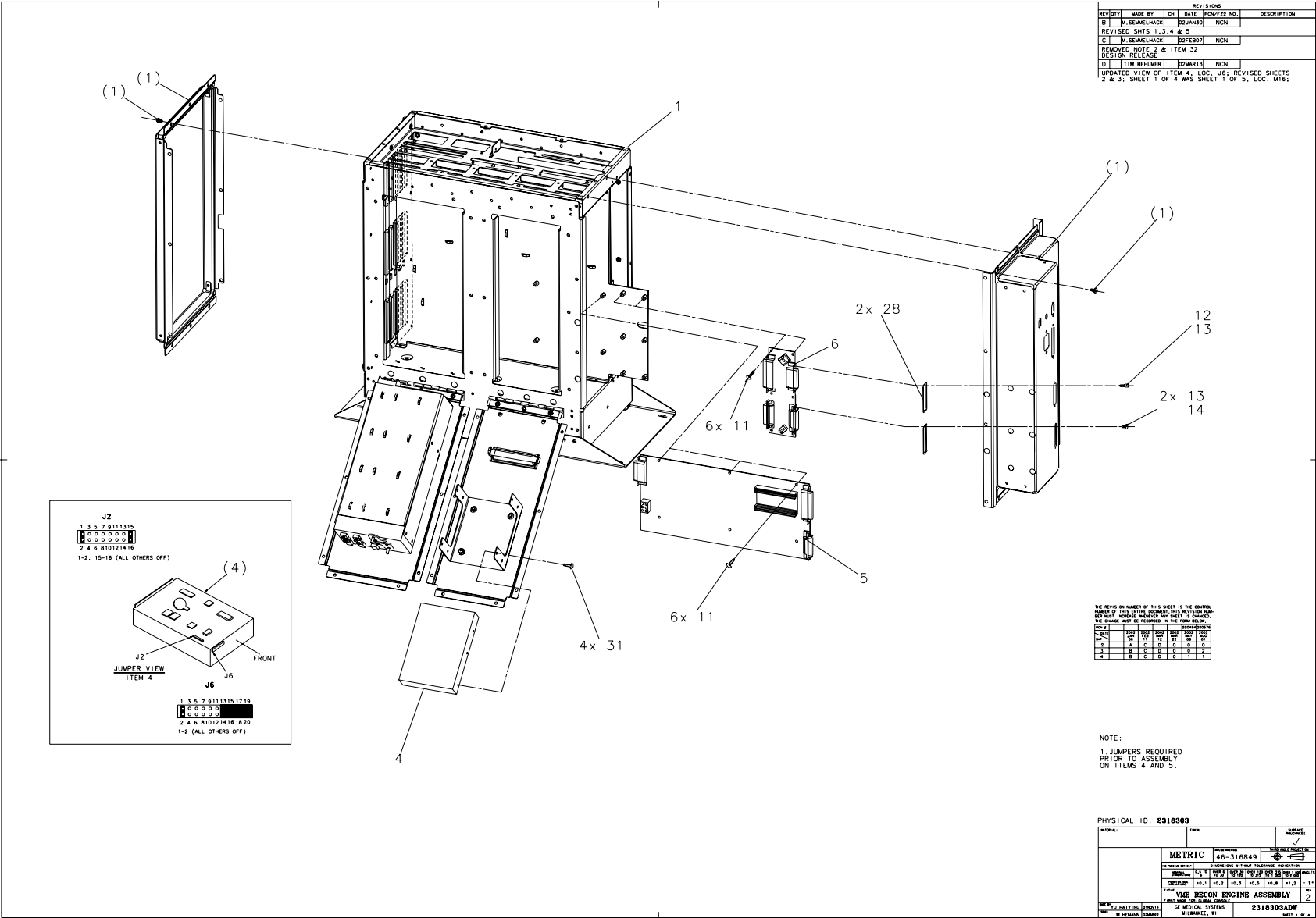


Figure 7-9 VME Recon Engine Assembly (2318303, sheet 1, rev 2)

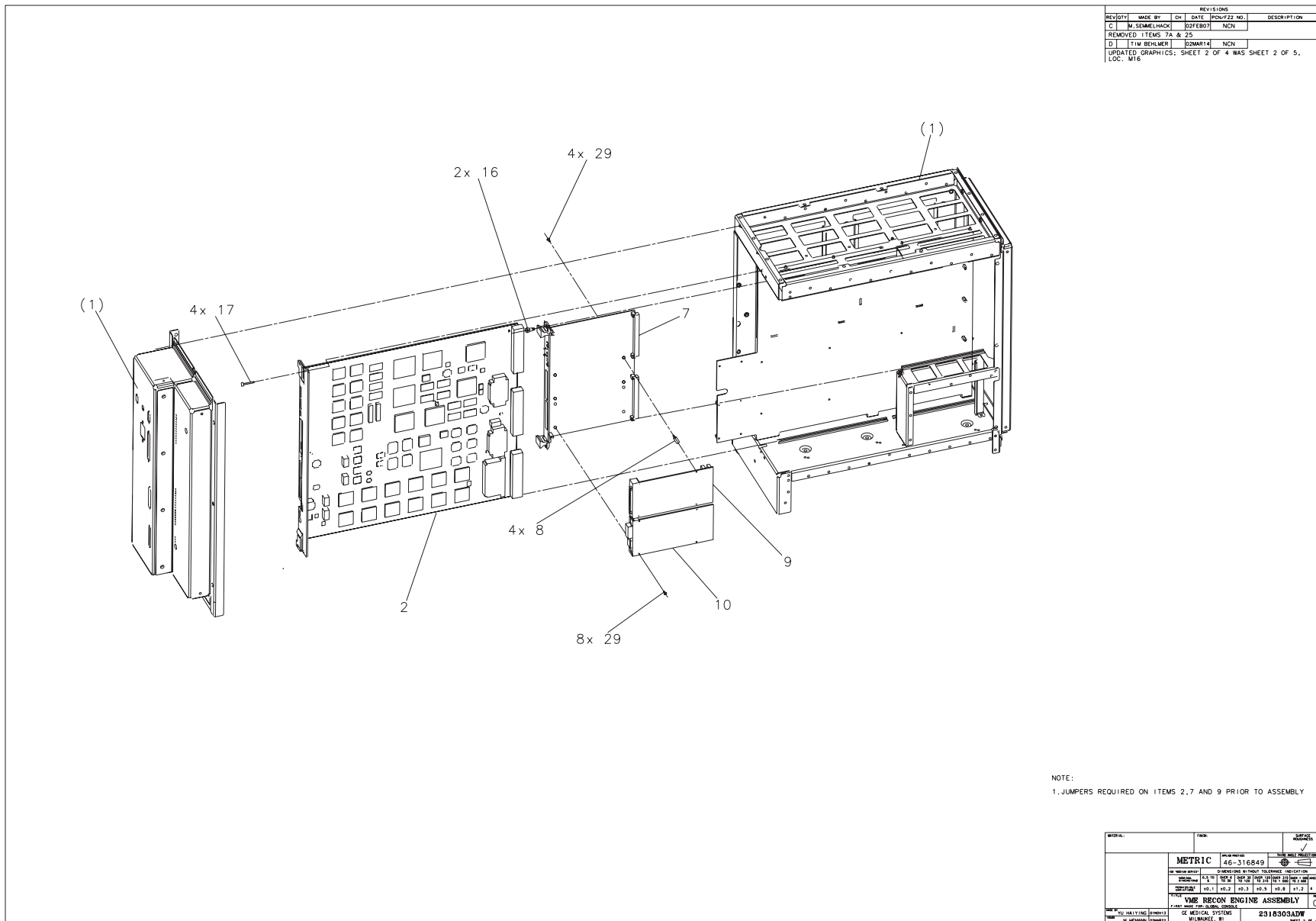


Figure 7-10 VME Recon Engine Assembly (2318303, sheet 2, rev 0)

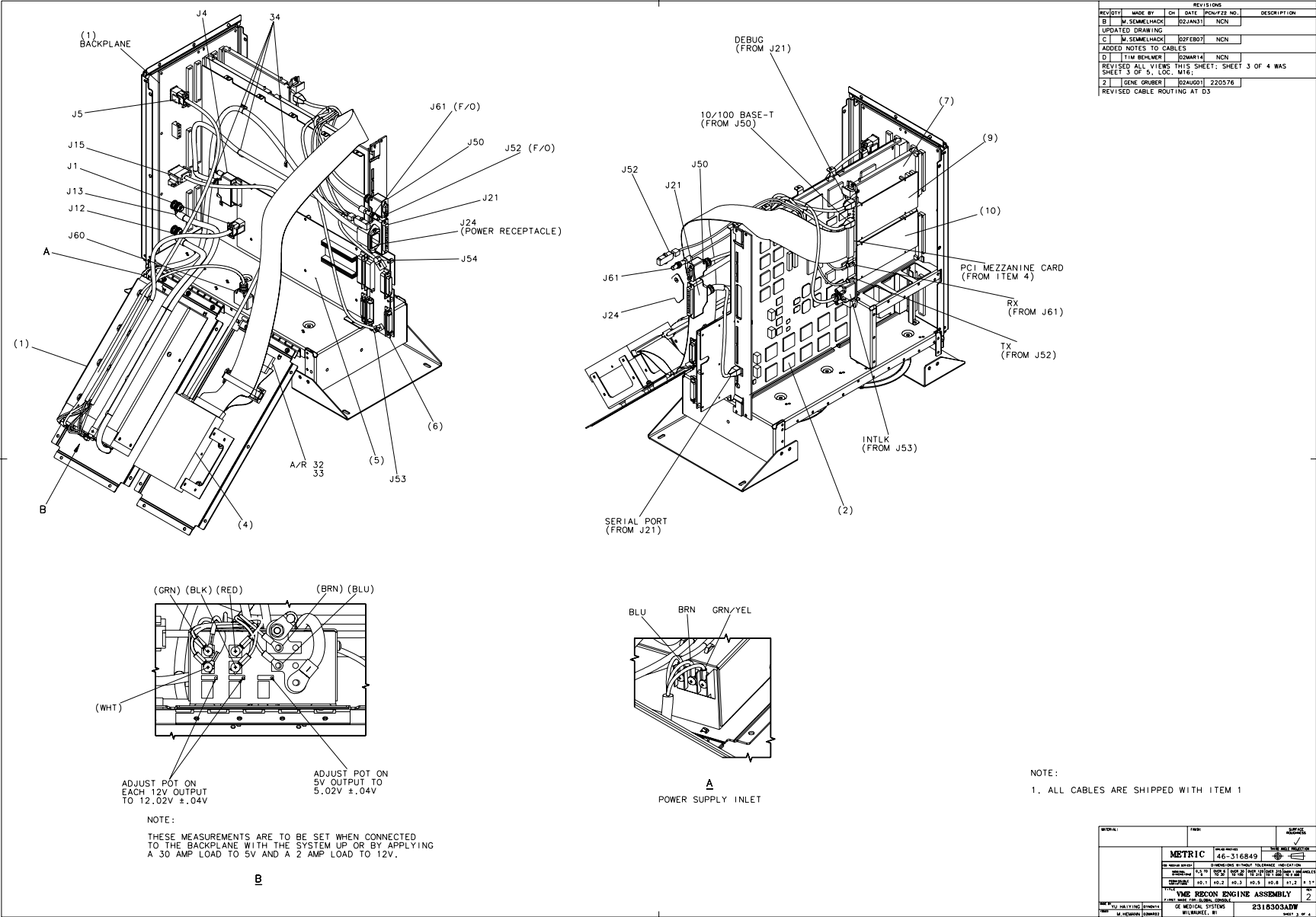


Figure 7-11 VME Recon Engine Assembly (2318303, sheet 3, rev 2)

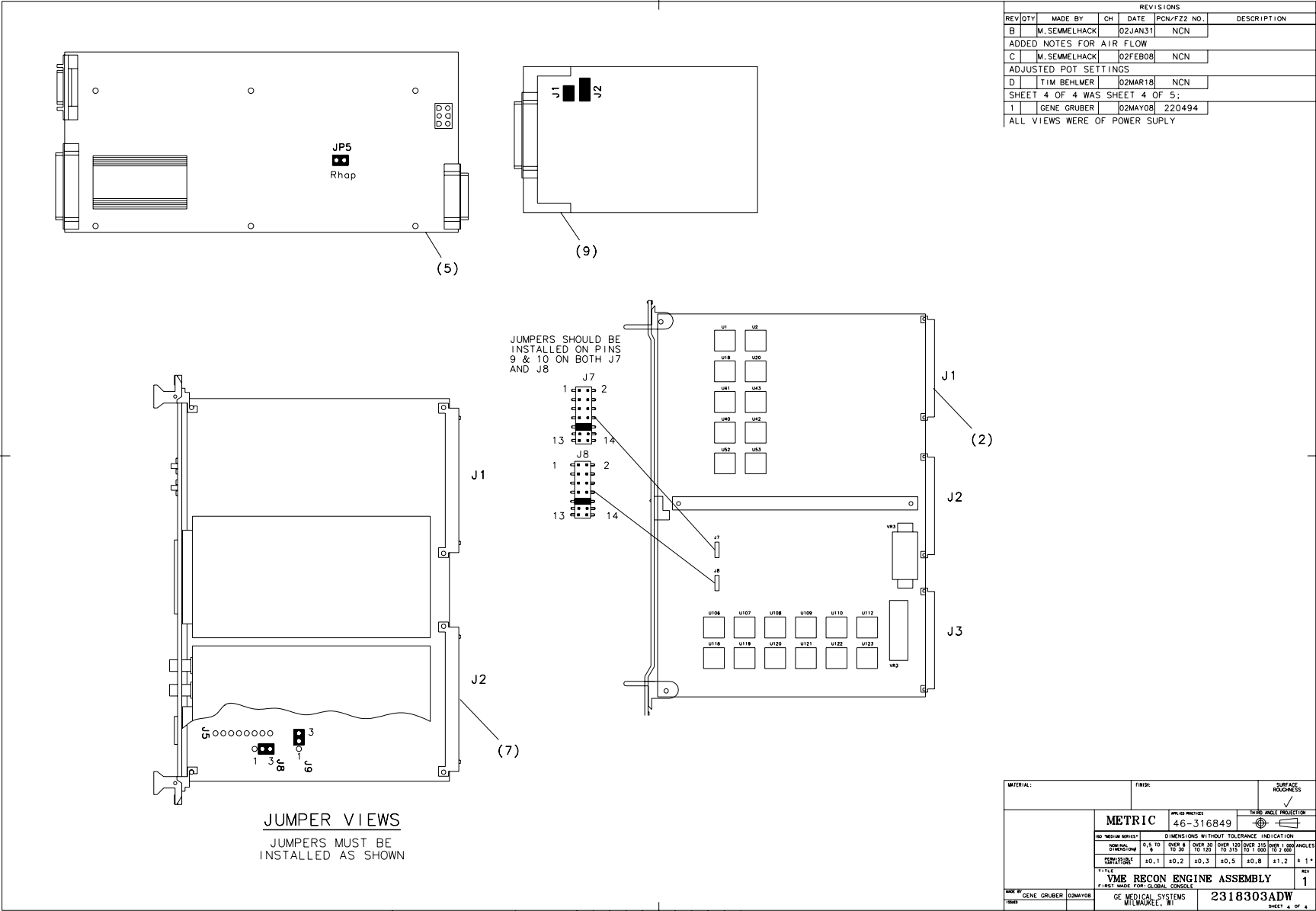


Figure 7-12 VME Recon Engine Assembly (2318303, sheet 4, rev 1)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2333179	1.000	Yes	CHASSIS,VME RECON BOX {For further detail, see Section 1.2.1 , beginning on page 584 .}
1A	2307794	ALT	Yes	VME CHASSIS ASSEMBLY {For further detail, see Section 1.2.2 , beginning on page 594 .}
2	2216467	1.000	Yes	PEGASUS IMAGE GENERATOR P
4	2225004-4	1.000	Yes	HARD DRIVE 18GB 15KRPM SC
5	2167014	1.000	Yes	RHAPSODE INTERCOM
6	2269601	1.000	Yes	PRESCRIBED TILT BOARD
7	2339284	1.000	Yes	PPC MVME2400 ASSY, HAWK3
9	2265396	1.000	Yes	CIRCUIT BOARD-SCSI DAUGHT
10	2259530	1.000	Yes	H3 PMC DIP
8	2209277	2.000	No	STANDOFF/SPACER FEMALE/FE
11	46-316540P2	12.000	No	SCREW PAN HEAD 3 MM 8 MM
12	46-221417P5	1.000	No	FEMALE SCREWLOCKS FOR HD
13	46-170686P3	AR	No	LOCTITE 242, 50CC
14	46-287187P4	2.000	No	
16	2280022	2.000	No	HANDSCREW, MOTOROLA BOARD
17	2142006	4.000	No	HEX SOC SCREW M2,5 X 16
28	2119591-2	2.000	No	25 PIN SUB-D EMI FILTER G
29	46-328421P4	10.000	No	SCREW PAN HEAD 2.5 MM 6 MM
31	46-170015P10	8.000	No	SEMS,6-32 X .25
32	46-136246P14	0.250 ft	No	2" WIDE, WHT, ADH BACK, VELCRO
33	46-136246P21	0.250 ft	No	2" WIDE X 25 YD WHITE VELCRO
34	46-208758P2	4.000	No	5.5" X .140" SELF-LOCKING CABLE TIE

Table 7-6 VME Recon Engine Assembly (2318303 BOM, Rev 4)

This Page Intentionally Blank.

1.2.1 Chassis, VME Recon Box (2333179)

FRU = YES

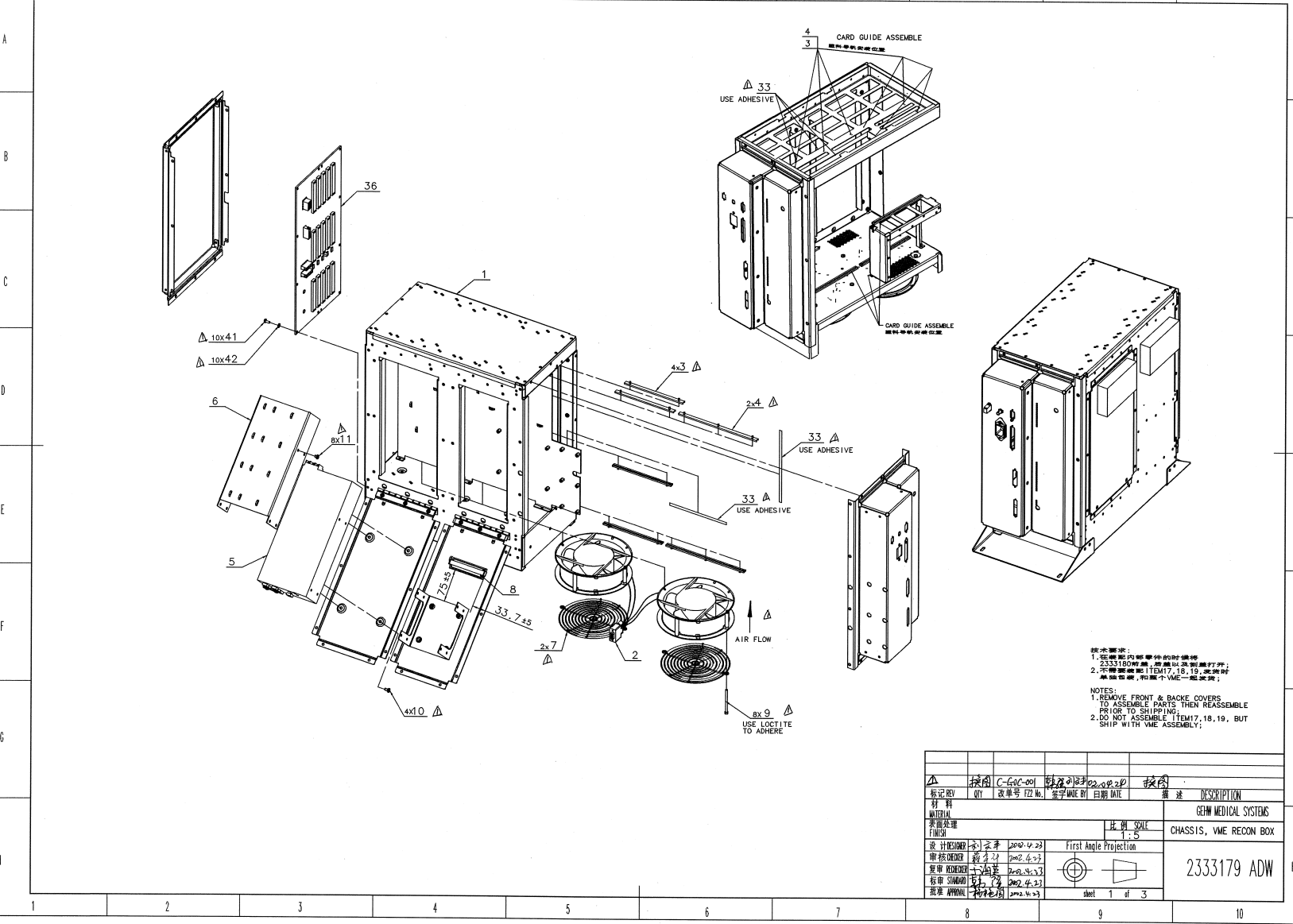


Figure 7-13 Chassis, VME Recon Box (2333179 sheet 1, rev 1)



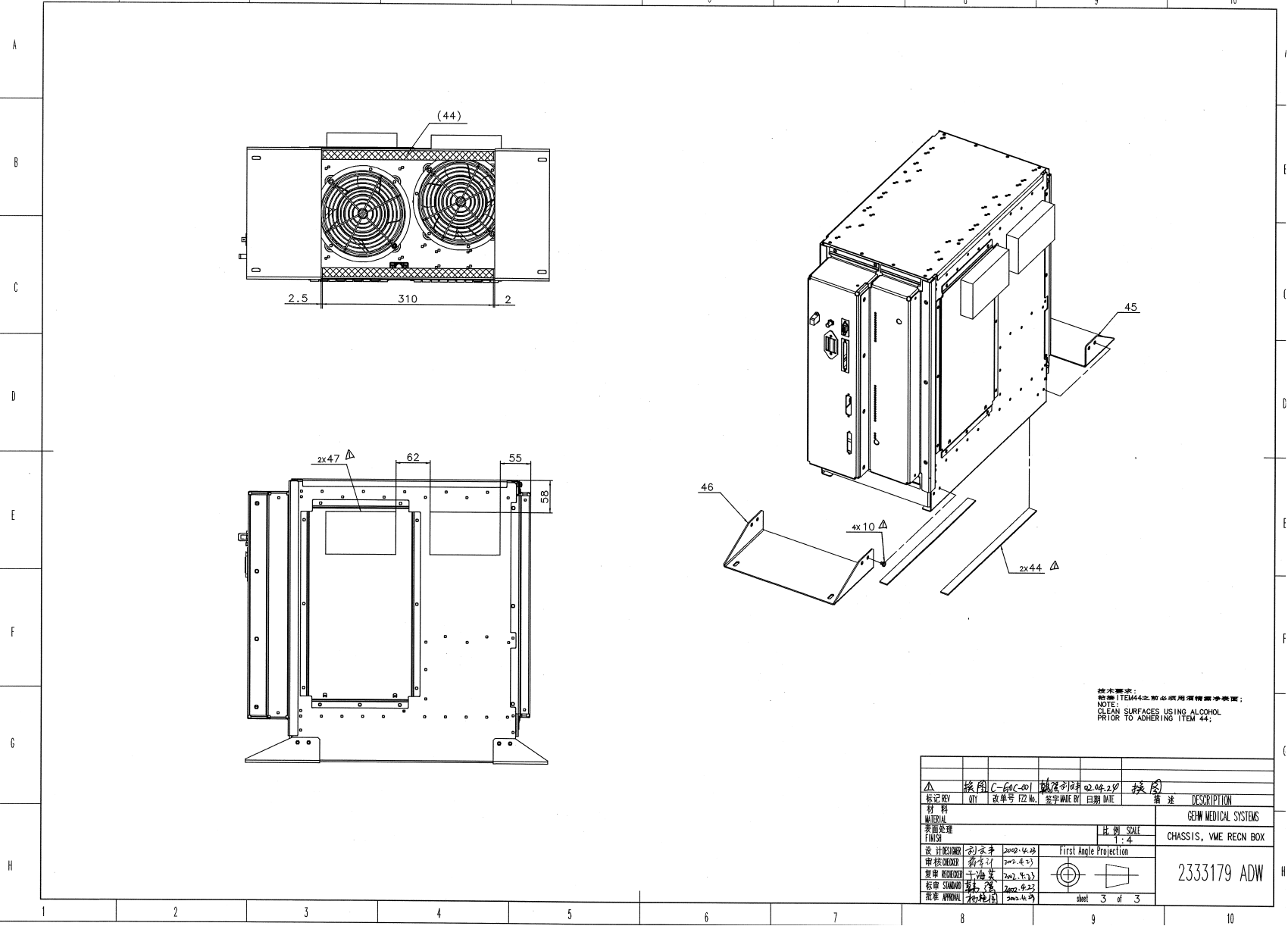


Figure 7-15 Chassis, VME Recon Box (2333179 sheet 3, rev 1)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2333180	1.000	No	VME RECON BOX
2	2324186	1.000	No	FAN ASSEMBLY, VME CHASSIS {For further detail, see Section 1.2.1.1 , beginning on page 590 .}
3	46-221513P3	4.000	No	CLAMP OR CLIP
4	46-221513P7	2.000	No	CLAMP OR CLIP
5	2260421-2	1.000	No	VME CONSOLE POWER SUPPLY {For further detail, see Section 1.2.1.2 , beginning on page 592 .}
6	2333181	1.000	No	POWER SUPPLY TIE DOWN BRA
7	2321164	2.000	No	FAN GUARD, 172 MM FANS
8	46-220191P5	1.000	No	PANDUIT FCM3-A-14 GRAY
9	2324209	8.000	No	SCREW SOCKET HEAD M4 X 60
10	2103585	16.000	No	SCREW PAN HEAD 4 MM 8 MM
11	2109873-11	8.000	No	SCREW PAN HEAD 4 MM X 5 M
12	2319051	1.000	No	5V POWER CABLE (+)
13	2319050	1.000	No	5V POWER CABLE, (-)
14	2318710	1.000	No	VME POWER CABLE, 110 VOLT
15	2307792	1.000	No	POWER CABLE, VME FANS
16	2337509	1.000	No	CABLE, ETHERNET
17	2322199	1.000	No	HD POWER CABLE, VME RECON
18	2315116	1.000	No	SCSI CABLE, 68 PIN RIBBON
19	2318714	1.000	No	CABLE, INTERNAL ABORT
20	2265406	1.000	No	FIBER-OPTIC COUPLING RECE
21	2142224	1.000	No	VME POWER SUPPLY TO RECEP
23	2318713	1.000	No	VME FIBER OPTIC CABLE
24	2318711	1.000	No	VME OCTOPUS CABLE
25	2319648	1.000	No	FIBER OPTIC CABLE, SC TO
26	2119591-3	1.000	No	37 PIN SUB-D EMI FILTER G
27	2119591-4	1.000	No	9 PIN SUB-D EMI FILTER GA

Table 7-7 VME Recon Box Chassis (2333179 BOM, rev 4)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
28	46-328421P3	2.000	No	SCREW PAN HEAD 2.5 MM 5 M
29	2319650	1.000	No	FIBER OPTIC BULKHEAD RECE
30	46-221417P4	5.000	No	FEMALE SCREWLOCKS FOR HD
31	N9406JS	2.000	No	SCREW
32	46-208758P2	8.000	No	5.5" X .140" SELF-LOCKING
33	U0507AC	1.000	No	PROTECTOR
36	2141161	1.000	No	RHAPSODE BACKPLANE
37	46-208935P8	4.000	No	1/4-20, 7/16" HEX X 3/16
38	46-220184P54	2.000	No	SCREW
39	1000904P475	8.000	No	ID 0.281 OD 0.625 THK .
40	46-221449P6	4.000	No	SCREW SIZE 1/4 (0.250) PH
41	2103585-2	10.000	No	SCREW PAN HEAD 4 MM 10 MM
42	46-328431P1	10.000	No	WASHER PLAIN - LARGE 4.3
43	46-221852P2	5.000	No	SWITCH
44	2327361	1.000	No	SINGLE SIDED SILICON RUBB
45	2323985	1.000	No	BRACKET, VME CHASSIS-BACK
46	2323987	1.000	No	BRACKET, VME CHASSIS, FRO
47	2332913	2.000	No	FOAM, POLYETHLENE, SELF A
48	46-312679P2	4.000	No	INTERNALLY RIBBED GROMMET
49	12-276	4.000	No	ID 0.2040 OD 0.5630 THK .

Table 7-7 VME Recon Box Chassis (2333179 BOM, rev 4) (Continued)

This Page Intentionally Blank.

1.2.1.1 Fan Assembly, VME Chassis (2324186)

FRU = NO

USED ON:
2333179,
2307794

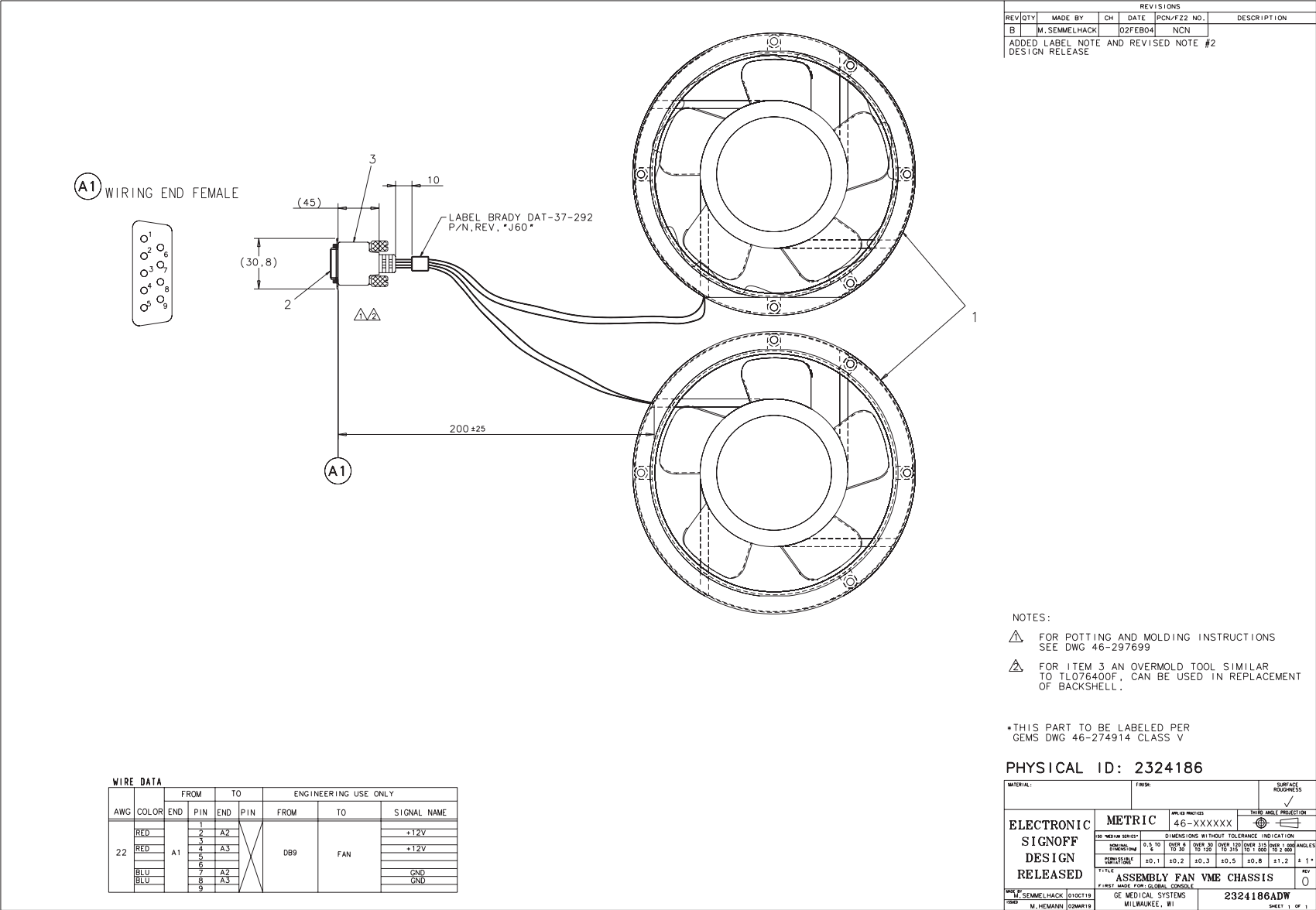


Figure 7-16 Fan Assembly, VME Chassis (2324186, sheet 1, rev 0)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2307791	2.000	No	FAN,172MMX150X51MM, 120 C
2	2324197	1.000	No	DB9 CONNECTOR, SOLDER CUP
3	2324196	1.000	No	DB9 BACKSHELL KIT WITH SC

Table 7-8 Fan Assembly, VME Chassis (2324186 BOM, Rev 0)

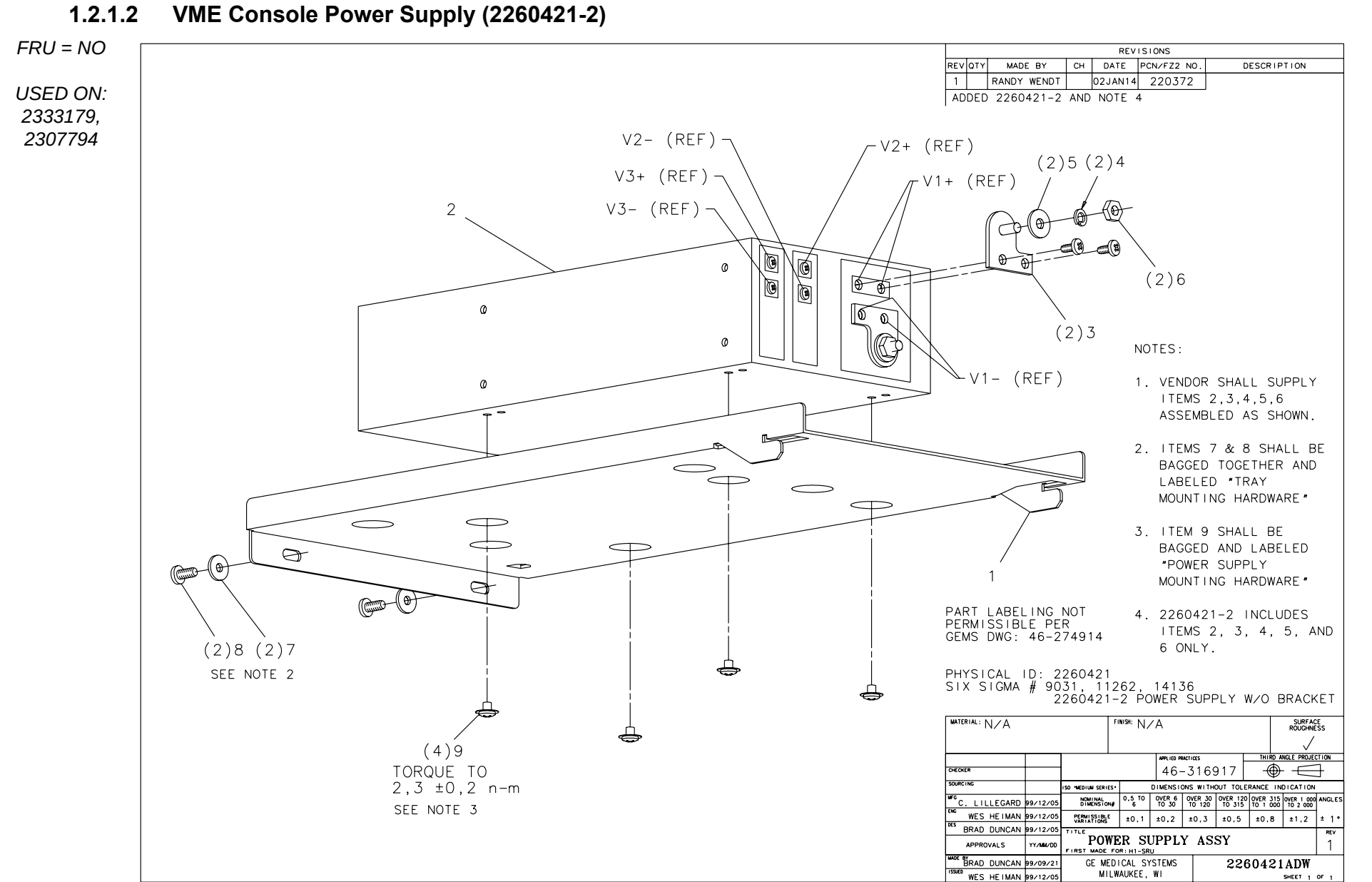


Figure 7-17 VME Console Power Supply (2260421, sheet 1, rev 1)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
2	2255986	1.000	No	POWER SUPPLY, MODULAR,600
3	2255660	2.000	No	5V POWER SUPPLY ADAPTER
4	46-220181P5	2.000	No	WASHER #10 (.190) SCREW S
5	12-276	2.000	No	ID 0.2040 OD 0.5630 THK .
6	46-208941P9	2.000	No	00010-32 .375 HEX X .124

Table 7-9 VME Console Power Supply (2260421-2 BOM, rev 1)

1.2.2 VME Chassis Assembly (2307794)
FRU = YES

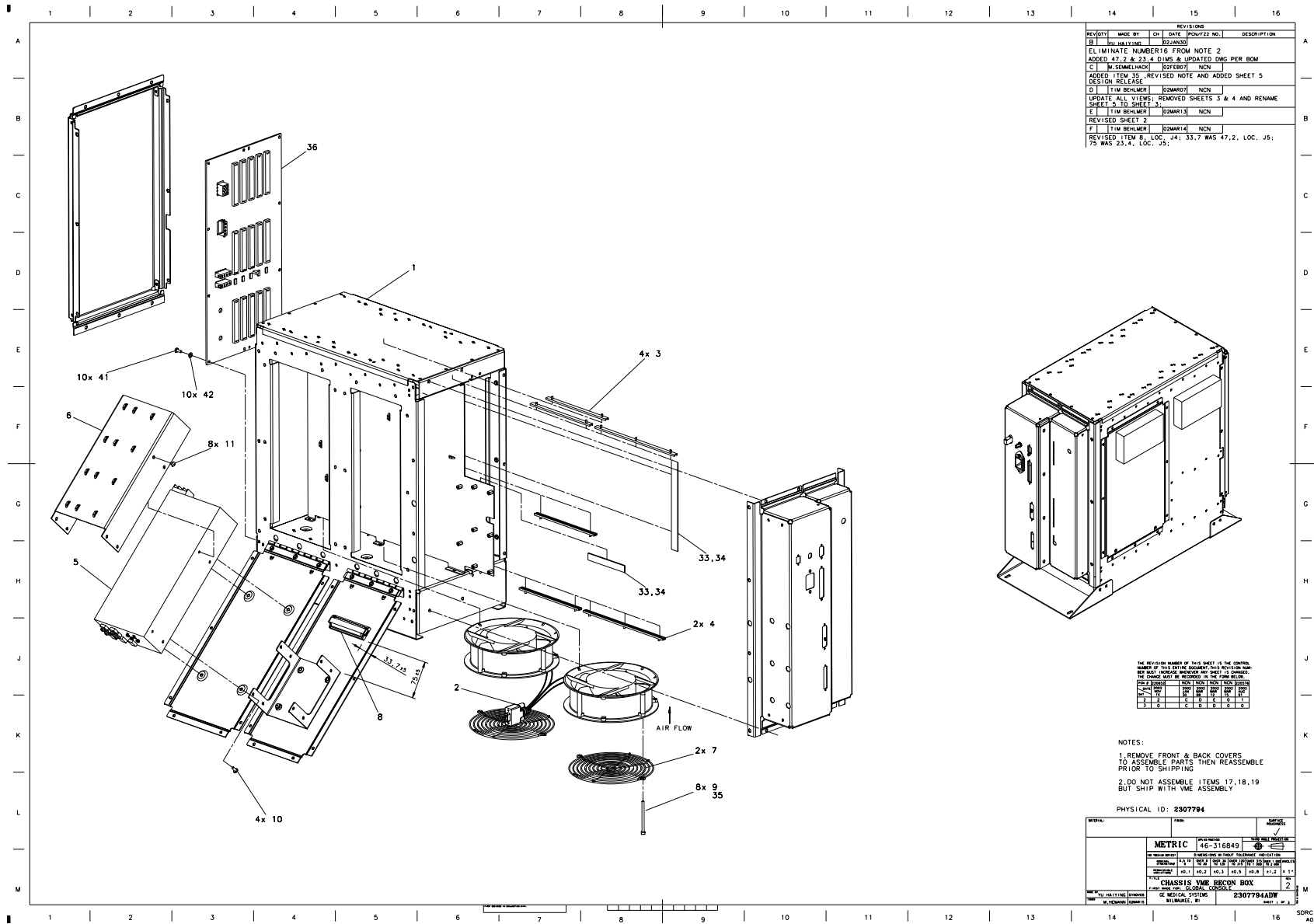


Figure 7-18 VME Chassis Assembly (2307794, sheet 1, rev 2)



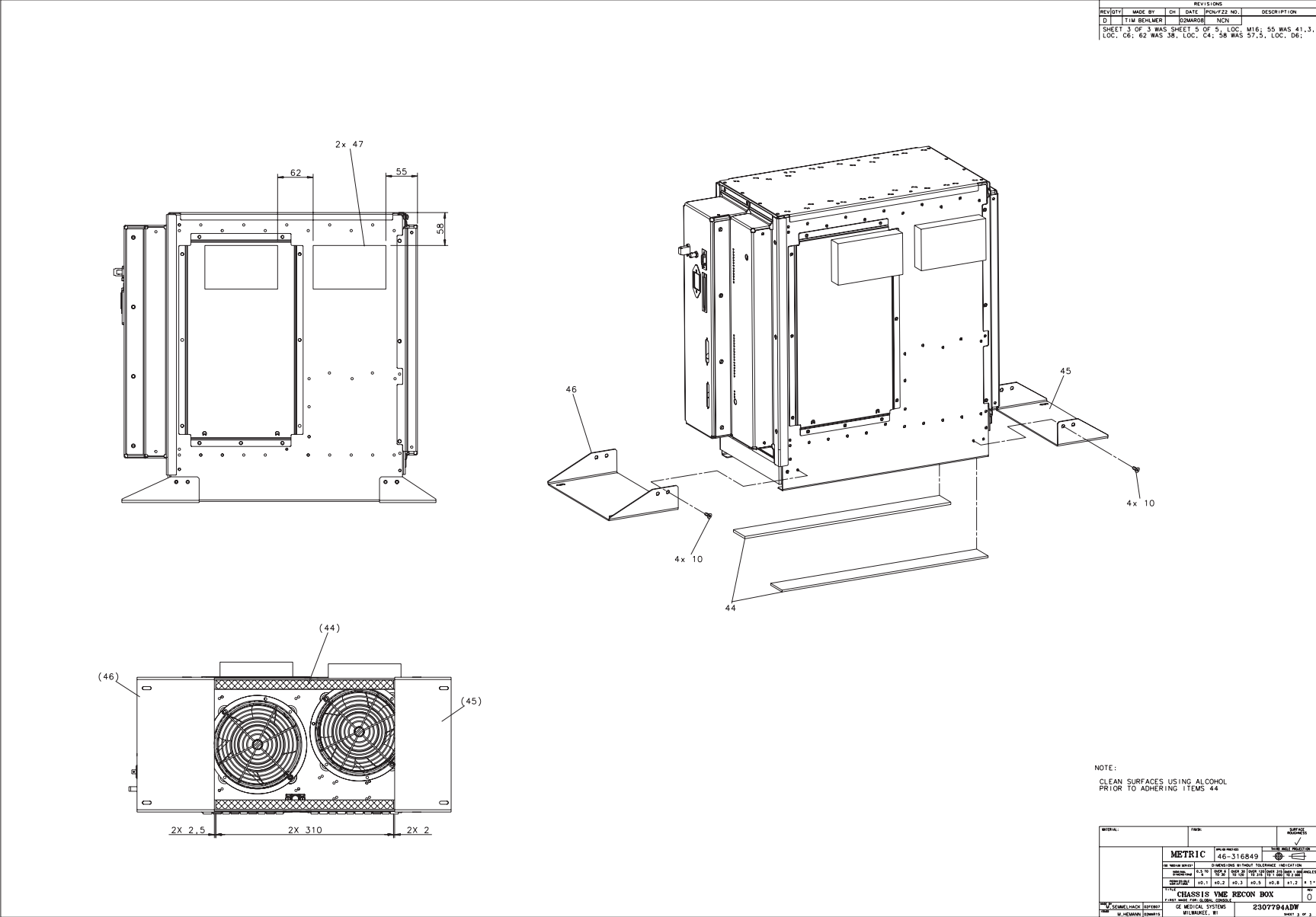


Figure 7-20 VME Chassis Assembly (2307794, sheet 3, rev 0)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
2	2324186	1.000	Yes	FAN ASSEMBLY, VME CHASSIS {For further detail, see Section 1.2.1.1 , beginning on page 590 .}
5	2260421-2	1.000	Yes	VME CONSOLE POWER SUPPLY {For further detail, see Section 1.2.1.2 , beginning on page 592 .}
1	2304975	1.000	No	CHASSIS, VME RECON BOX
3	46-221513P3	4.000	No	CLAMP OR CLIP
4	46-221513P7	2.000	No	CLAMP OR CLIP
6	2320003	1.000	No	POWER SUPPLY TIE DOWN BRA
7	2321164	2.000	No	FAN GUARD, 172 MM FANS
8	46-220191P5	1.000	No	PANDUIT FCM3-A-14 GRAY
9	2324209	8.000	No	SCREW SOCKET HEAD M4 X 60
10	2103585	16.000	No	SCREW PAN HEAD 4 MM 8 MM
11	2109873-11	8.000	No	SCREW PAN HEAD 4 MM X 5 M
12	2319051	1.000	No	5V POWER CABLE (+)
13	2319050	1.000	No	5V POWER CABLE, (-)
14	2318710	1.000	No	VME POWER CABLE, 110 VOLT
15	2307792	1.000	No	POWER CABLE, VME FANS
16	2337509	1.000	No	CABLE, ETHERNET
17	2322199	1.000	No	HD POWER CABLE, VME RECON
18	2315116	1.000	No	SCSI CABLE, 68 PIN RIBBON
19	2318714	1.000	No	CABLE, INTERNAL ABORT
20	2265406	1.000	No	FIBER-OPTIC COUPLING RECE
21	2142224	1.000	No	VME POWER SUPPLY TO RECEP
23	2318713	1.000	No	VME FIBER OPTIC CABLE
24	2318711	1.000	No	VME OCTOPUS CABLE
25	2319648	1.000	No	FIBER OPTIC CABLE, SC TO
26	2119591-3	1.000	No	37 PIN SUB-D EMI FILTER G
27	2119591-4	1.000	No	9 PIN SUB-D EMI FILTER GA

Table 7-10 VME Chassis Assembly (2307794 BOM, rev 0)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
28	46-328421P3	2.000	No	SCREW PAN HEAD 2.5 MM 5 M
29	2319650	1.000	No	FIBER OPTIC BULKHEAD RECE
30	46-221417P4	3.000	No	FEMALE SCREWLOCKS FOR HD
31	46-328422P12	2.000	No	SCREW COUNTERSUNK HEAD 4 MM
32	46-208758P2	8.000	No	5.5" X .140" SELF-LOCKING
33	46-220178P10	AR	No	GROMMET
34	46-170071P1	AR	No	MULTI-PURPOSE ADHESIVE US
35	46-170686P3	AR	No	LOCTITE 242, 50CC
36	2141161	1.000	No	RHAPSODE BACKPLANE
37	46-208935P8	4.000	No	1/4-20, 7/16" HEX X 3/16
38	46-220184P54	2.000	No	SCREW
39	1000904P475	8.000	No	ID 0.281 OD 0.625 THK .
40	46-221449P6	4.000	No	SCREW SIZE 1/4 (0.250) PH
41	2103585-2	10.000	No	SCREW PAN HEAD 4 MM 10 MM
42	46-328431P1	10.000	No	WASHER PLAIN - LARGE 4.3
43	46-221852P2	5.000	No	SWITCH
44	2327361	1.000	No	SINGLE SIDED SILICON RUBB
45	2323985	1.000	No	BRACKET, VME CHASSIS-BACK
46	2323987	1.000	No	BRACKET, VME CHASSIS, FRO
47	2332913	2.000	No	FOAM, POLYETHLENE, SELF A

Table 7-10 VME Chassis Assembly (2307794 BOM, rev 0)

This Page Intentionally Blank.

USED ON:
2315653

Figure 7-21 Console Chassis Assembly (2333194 sheet 1, rev 1)

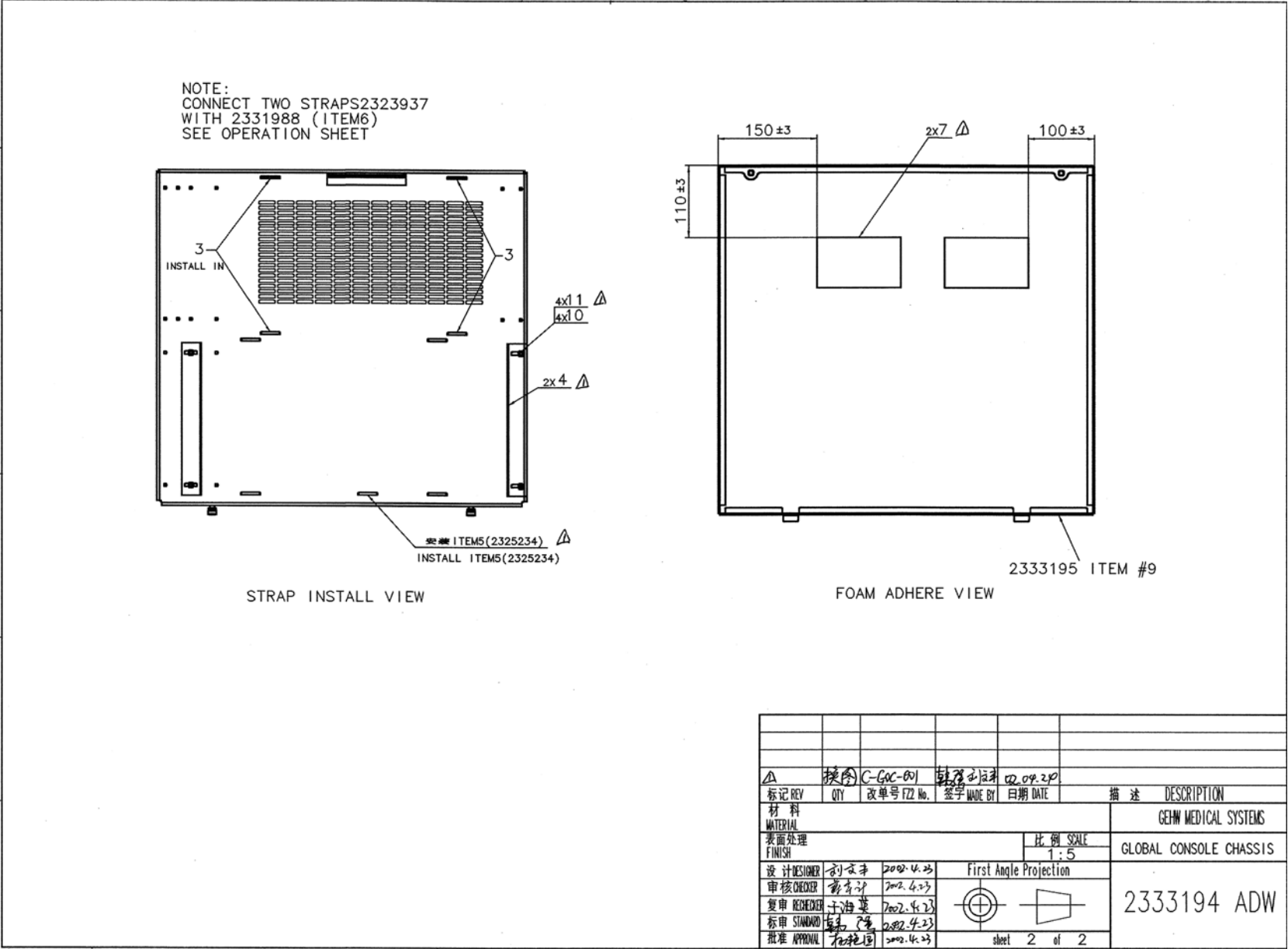


Figure 7-22 Console Chassis Assembly (2333194 sheet 2, rev 1)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2335086	1.000	No	CONSOLE CHASSIS FROM SUPP
2	2174159-2	1.000	No	OUTLET BOX ASSEMBLY GLOBA {For further detail, see Section 1.3.1 , beginning on page 604 .}
3	2323937	4.000	No	STRAP, HOLD DOWN
4	2323986	2.000	No	BRACKET, HOST COMPUTER
5	2325234	1.000	No	PULL STRAP, HOST COMPUTER
6	2331988	2.000	No	STRAP, HOLDING STRAP
7	2332913	2.000	No	FOAM, POLYETHLENE, SELF A
8	2323943	1.000	No	FAN ASSEMBLY, GLOBAL CONS {For further detail, see Section 1.3.2 , beginning on page 608 .}
9	46-208897P16	12.000	No	SCREW
10	46-328431P1	9.000	No	WASHER PLAIN - LARGE 4.3
11	2103585	9.000	No	SCREW PAN HEAD 4 MM 8 MM
12	2342018	1.000	No	LABEL, KEYBOARD/MOUSE I/O

Table 7-11 Chassis Assembly (2333194 BOM, rev 2)

This Page Intentionally Blank.

1.3.1 Outlet Box Assembly (2174159-2)

FRU = NO

USED ON:
2333194,
2317687

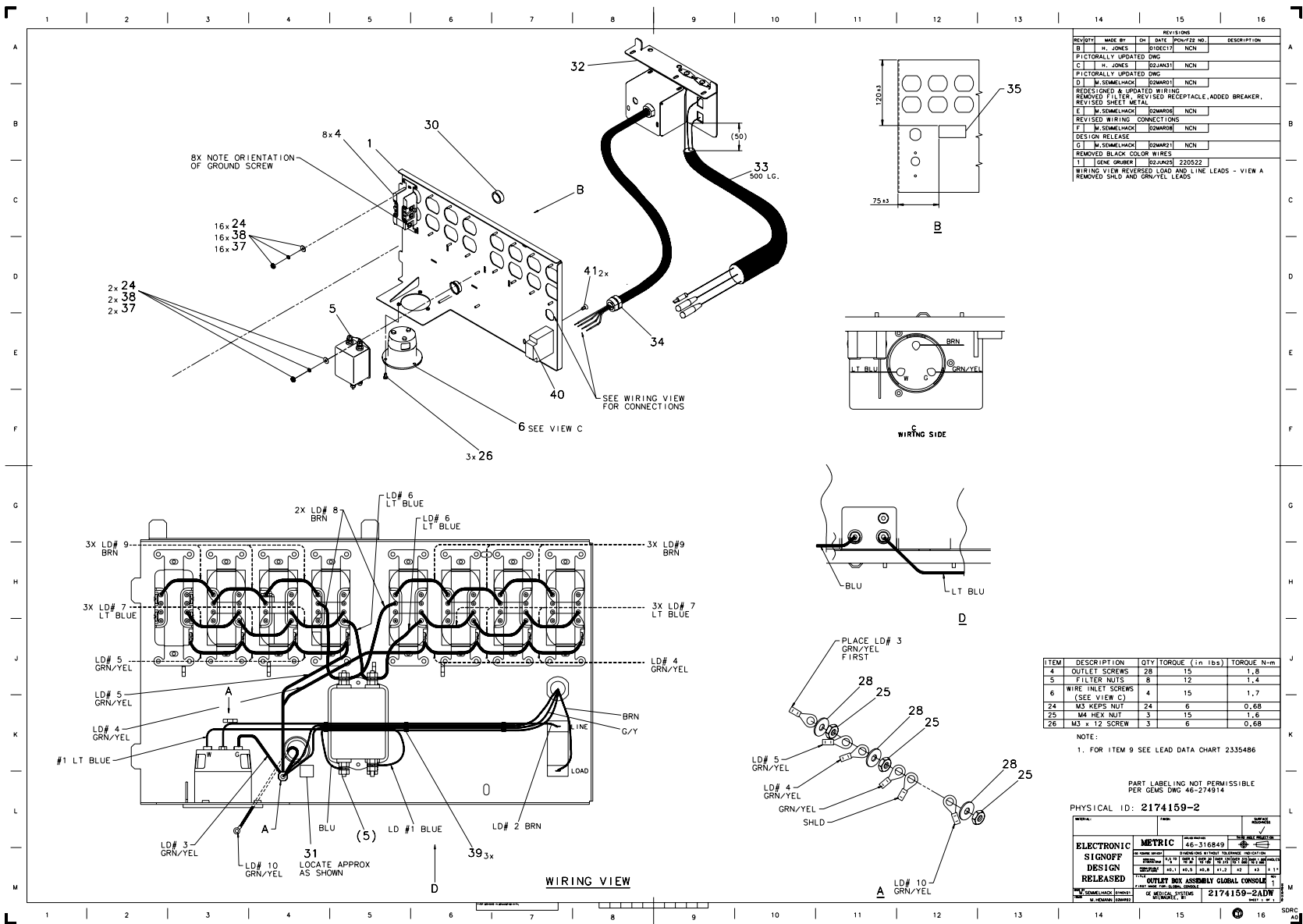


Figure 7-23 Outlet Box Assembly (2174159-2, sheet 1, rev 1)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2334322	1.000	No	OUTLET BOX COVER, GOC
4	46-278446P1	8.000	No	15A, 125VAC, 2 POLE, 3 WI
5	46-297335P3	1.000	No	16A, 250V 2-STAGE PWR LIN
6	46-243925P1	1.000	No	20A, 125V, 2-POLE, 3 WIRE
9	2335486	1.000	No	LEAD DATA FOR OUTLET BOX
24	46-328425P10	18.000	No	NUT HEXAGON 3 MM
25	46-328425P1	3.000	No	NUT HEXAGON 4 MM
26	46-312358P30	3.000	No	M3 X 0.5 X 12MM METRIC FR
28	46-328430P2	3.000	No	WASHER PLAIN - NORMAL 4.3
30	46-208695P61	1.000	No	HEYCO PART NO.2126 SB-87
31	46-271110P1	1.000	No	LABEL, PROTECTIVE EARTH (
32	2190750-2	1.000	No	POWER SWITCH BOX-GLOB CON {For further detail, see "Power Switch Box (2190750-2)," on page 606.}
33	46-313024P5	1.650	No	
34	46-208697P57	1.000	No	HEYCO PART NO.1247 SR-7W
35	2150839	1.000	No	CAUTION RECEPTACLES NOT F
37	46-311804P2	18.000	No	3.2X9X.8 ZP ST DIN
38	46-311805P10	18.000	No	HELICAL LOCK WASHER
39	46-208758P2	3.000	No	5.5" X .140" SELF-LOCKING
40	2336074	1.000	No	CIRCUIT BREAKER 15 AMPS 2
41	46-328421P501	2.000	No	SCREW PAN HEAD 3 MM 5 MM

Table 7-12 Outlet Box Assembly (2174159-2 BOM, rev 0)

Power Switch Box (2190750-2)

FRU = NO

USED ON:
2174159-2

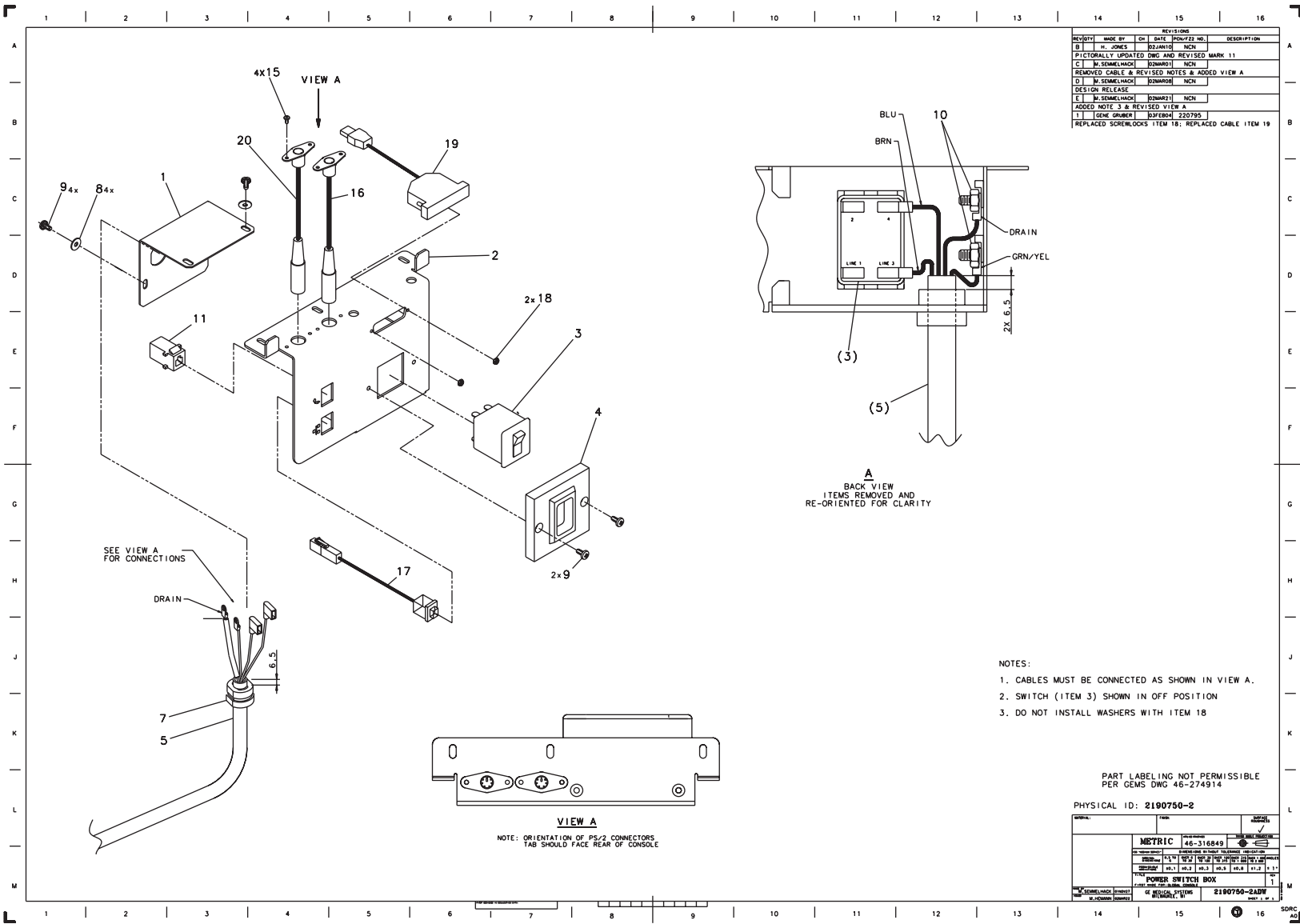


Figure 7-24 Power Switch Box (2190750-2, sheet 1, rev 1)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2317186	1.000	No	COVER-BACK COVER FOR SERV
2	2317185	1.000	No	BRACKET-FRONT SERVICE PAN
3	2336081	1.000	No	ROCKER SWITCH, 2PST
4	2150830	1.000	No	POWER BREAKER BEZEL
5	2200242-2	1.000	No	PWR CBL OUTLET BOX TO SER
7	46-208697P57	2.000	No	HEYCO PART NO.1247 SR-7W
8	46-328430P2	4.000	No	WASHER PLAIN - NORMAL 4.3
9	2103585	6.000	No	SCREW PAN HEAD 4 MM 8 MM
10	46-328425P1	2.000	No	NUT HEXAGON 4 MM
11	2108900	1.000	No	CONNECTOR 8 CONTACTS 150
15	46-328421P3	4.000	No	SCREW PAN HEAD 2.5 MM 5 M
16	2315117	1.000	No	PS/2 BULKHEAD CABLE
17	2318308	1.000	No	ETHERNET CABLE, CHASSIS S
18	46-170012P34	2.000	No	4-40, 1/4" HEX, 3/32 THIC
19	2368504	1.000	No	CONVERTER, DB25 TO USB
20	2315117-2	1.000	No	CABLE, PS/2 KEYBOARD CABL

Table 7-13 Power Switch Box (2190750-2 BOM, rev 2)

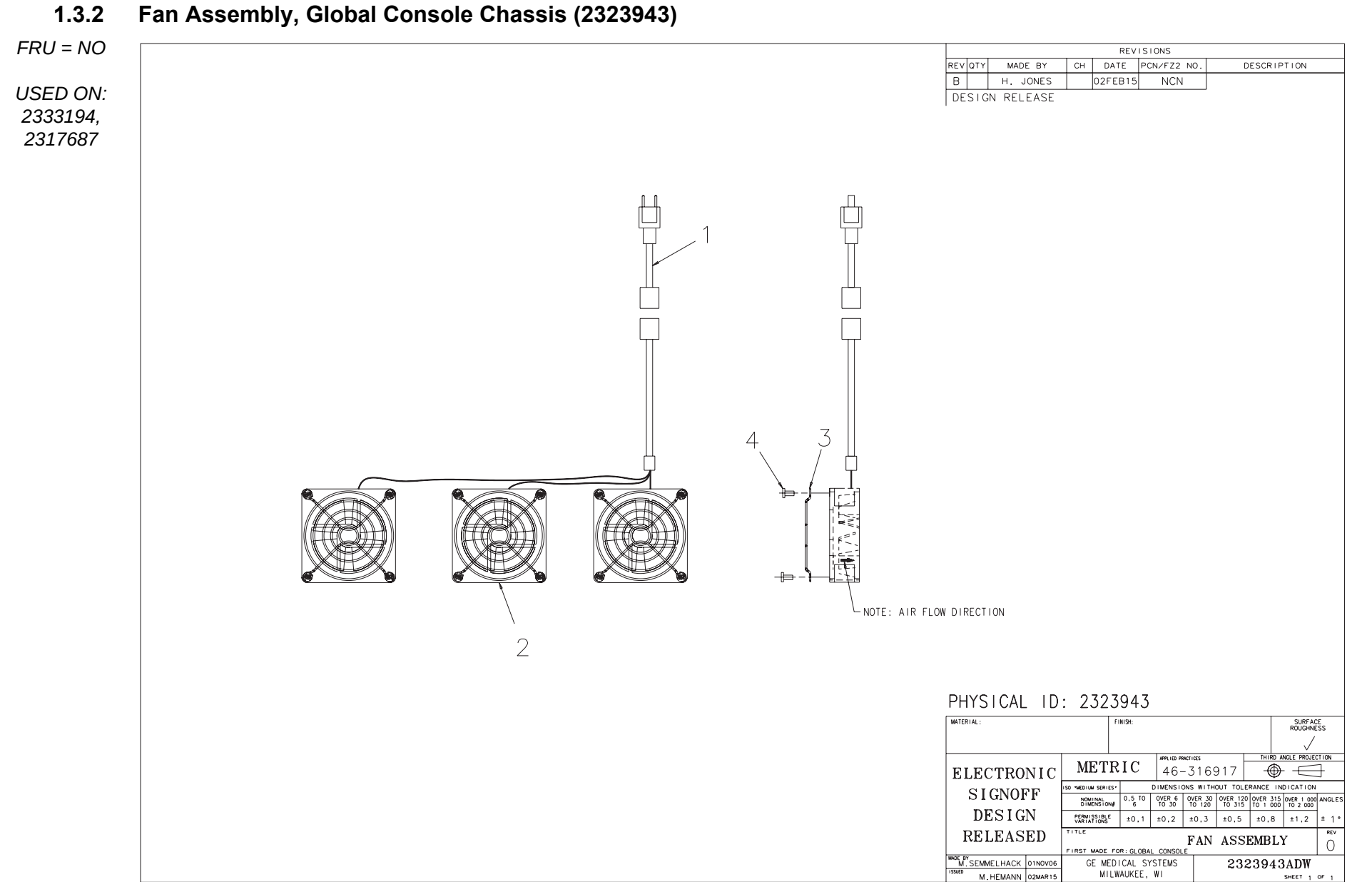


Figure 7-25 Fan Assembly, Console Chassis (2323943, sheet 1, rev 0)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2326685	1.000	No	POWER CABLE FOR 110V FAN
2	2326686	1.000	No	FAN ASSEMBLY, 110V, 120X120
3	2293755	3.000	No	TWIN COOLING FAN GRILL
4	46-208897P16	12.000	No	SCREW

Table 7-14 Fan Assembly, Global Console Chassis (2323943 BOM, rev 0)

USED ON:
2315653

Figure 7-26 Console Chassis Assembly (2317687, sheet 1, rev 1)

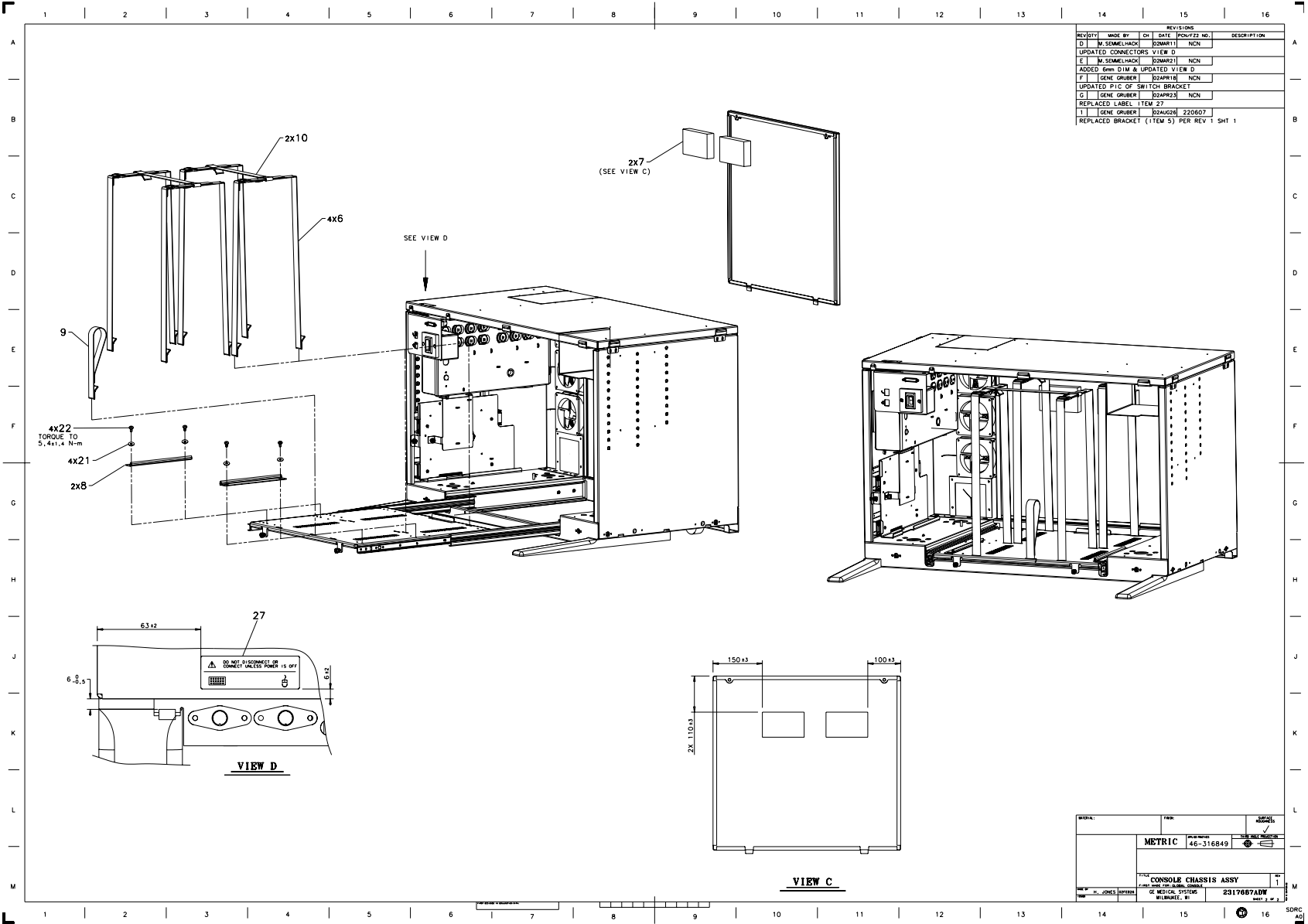


Figure 7-27 Console Chassis Assembly (2317687, sheet 2, rev 1)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
2	2174159-2	1.000	Yes	OUTLET BOX ASSEMBLY GLOBAL CONSOLE {For further detail, see Section 1.3.1 , beginning on page 604 .}
3	2318681	2.000	Yes	ANTI-TIP LEG, CONSOLE CHA
4	2319043	4.000	Yes	CASTER, 2 1/2" SWIVEL
11	2323943	1.000	Yes	FAN ASSEMBLY, GLOBAL CONSOLE {For further detail, see Section 1.3.2 , beginning on page 608 .}
1	2310119	1.000	No	GLOBAL CONSOLE CHASSIS
5	2317787	1.000	No	BRACKET, ETHERNET SWITCH
6	2323937	4.000	No	COMPUTER HOLD DOWN STRAP
7	2332913	2.000	No	FOAM, POLYETHYLENE, SELF A
8	2323986	2.000	No	BRACKET, HOST COMPUTER
9	2325234	1.000	No	PULL STRAP, HOST COMPUTER
10	2331988	2.000	No	STRAP, HOLDING STRAP
12	2109877-2	18.000	No	NUT HEXAGON 6 MM
13	46-208897P16	12.000	No	SCREW
14	2283156	2.000	No	RUBBER BOOT, REAR ACCESS
15	2331392	2.000	No	BRACKET FOR RUBBER BOOT
16	2331393	8.000	No	SCREW, POLYTOP SCREW, M4
17	2109867-25	4.000	No	SCREW HEXAGON SOCKET
18	2109878-4	4.000	No	WASHER PLAIN - NORMAL 10.
19	2109881-4	4.000	No	WASHER LOCK - INTERNAL 10
20	2325236	2.000	No	LEVELER, CONSOLE ANTI-TIP
21	46-328431P1	9.000	No	WASHER PLAIN - LARGE 4.3
22	2103585	9.000	No	SCREW PAN HEAD 4 MM 8 MM
24	46-328431P2	19.000	No	WASHER PLAIN - LARGE 6.4
25	2109881-3	18.000	No	WASHER LOCK - INTERNAL 6.
26	46-170686P3	AR	No	LOCTITE 242, 50CC
27	2342018	1.000	No	LABEL, KEYBOARD/MOUSE I/O

Table 7-15 Console Chassis Assembly (2317687 BOM, rev 1)

This Page Intentionally Blank.

FRU = YES

Figure 7-28 SCSI Tower Assembly (2368838, sheet 1, rev 0)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2319046	1.000	No	DVDRAM PLAYER
2	2318598	1.000	No	SONY MOD SMO-F551-SD
4	2315118	1.000	No	2 BAY EXTERNAL SCSI TOWER
5	46-328421P501	8.000	No	SCREW PAN HEAD 3 MM 5 MM
8	2350404	1.000	No	TERMINATOR, CENT50, ACTIV

Table 7-16 SCSI Tower, DVD-RAM, MOD (2368838 BOM, rev 0)

1.6 Global Modem Kit (2245794)

FRU = YES

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
3	2216697	1.000	Yes	CABLE (RJ11 25 FEET)
5	2249450	1.000	Yes	GLOBAL MODEM ADAPTER ASSY
1	2245773	1.000	No	MULTI-TECH GLOBAL MODEM
2	2253088	1.000	No	LABEL,CUSTOMER PROPERTY F
4	2255874	1.000	No	CABLE, RJ11 MODULAR, 6-CON

Table 7-17 Global Modem Kit (2245794 BOM, rev 3)

1.7 Seismic Kit (2336557)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2330594	2.000	No	SEISMIC BRACKET, GOC
2	2323984	1.000	No	BRACKET, SCSI TOWER MOUNT
3	2323937-2	1.000	No	STRAP, HOLD DOWN, 1" WIDE
4	46-136246P14	2.000 ft	No	2" WIDE, WHT, ADH BACK, V
5	46-136246P21	2.000 ft	No	2" WIDE X 25 YD WHITE VEL
6	46-208759P1	6.000	No	7.81" X .184" SELF-LOCKIN
7	46-328431P1	4.000	No	WASHER PLAIN - LARGE 4.3
8	46-328421P21	4.000	No	SCREW PAN HEAD 4 MM 40 MM
9	46-312346P2	4.000	No	NUT, M4 X 0.7 HEX WITH LO

Table 7-18 Seismic Kit (2336557 BOM, rev 0)

Section 2.0 — Side Covers (2337534)

FRU = YES

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2315656	1.000	No	RIGHT SIDE COVER, GLOBAL
2	2315656-2	1.000	No	LEFT SIDE COVER, GLOBAL C

Table 7-19 Side Covers (2337534 BOM, rev 0)

Section 3.0 — Console Display Monitors

3.1 CRT Monitors

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1A	2189192	2	Yes	PHILIPS 20 INCH MONITOR FOR LIGHTSPEED
1B	2244724	2	Yes	SONY GDM-500PS 21" COLOR MONITOR (RGB INPUTS)
1C	2254212	2	Yes	SONY GDM-5402T9 21" COLOR MONITOR (DB15 INPUT)

Table 7-20 Console Display Monitors

3.2 LCD Monitor & Cable (2309592)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2349792	1.000	Yes	18.1 IN LCD COLOR MONITOR
2	2295019	1.000	Yes	13W3-BD15 RGB LCD STRAIGHT CONN CBL

Table 7-21 LCD Monitor & Cable (2309592 BOM, rev 1)

Section 4.0 — DASMs

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1A	2249894	1	Yes	Analog DASM
1B	2249895	1	Yes	Digital DASM

Table 7-22 DASMs

Section 5.0 — Kits (Associated Items)

NUMBER	FRU	DESCRIPTION/NAME
2333829	Yes	GLOBAL CONSOLE, SCSI CABLE KIT {For further detail, see Section 5.1 , beginning on page 619.}
2333831	Yes	GLOBAL CONSOLE, ETHERNET / {For further detail, see Section 5.2 , beginning on page 619.}
2333832	Yes	GLOBAL CONSOLE, VIDEO/POWER CABLES {For further detail, see Section 5.3 , beginning on page 620.}
2333833	Yes	GLOBAL CONSOLE, ELECTRICAL COMPONENT {For further detail, see Section 5.4 , beginning on page 620.}
2333834	Yes	GLOBAL CONSOLE, BRACKETS & STRAPS {For further detail, see Section 5.5 , beginning on page 620.}
2333835	Yes	GLOBAL CONSOLE, MISC PM HARDWARE {For further detail, see Section 5.6 , beginning on page 621.}
2333836	Yes	GLOBAL CONSOLE, HARDWARE KIT {For further detail, see Section 5.7 , beginning on page 621.}
2333837	Yes	GLOBAL CONSOLE, PWR SWITCH & CABLES {For further detail, see Section 5.8 , beginning on page 623.}
2336912	Yes	LINUX CONSOLE SLIDE SERVICE KIT {For further detail, see Section 5.9 , beginning on page 623.}
2333829-2	Yes	SCSI CABLE KIT, GOC-LINUX {For further detail, see Section 5.10 , beginning on page 624.}

Table 7-23 Associated Items (2315653 AILs, rev. 2003-04-14)

5.1 SCSI Cable Kit (2333829)

FRU = YES

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2324697-2	1.000	No	SCSI CABLE, VHD68 TO CENT
3	2351053	1.000	No	CABLE, SCSI HD68 TO CENT5
4	2191429	1.000	No	NEXGEN CABLE FROM COMPUTE
5	2323942	1.000	No	CABLE, HD50 TO HD50
6	2317488	1.000	No	OCTOPUS CABLE, LINUX COMP
7	2315116	1.000	No	SCSI CABLE, 68 PIN RIBBON
8	2341433	1.000	No	CABLE, OCTOPUS - OCTANE 2
9	2341109	1.000	No	CABLE, SCSI, VHD68 TO HD6
10	2349539	1.000	No	SCSI CABLE, SCAN DATA, H1

Table 7-24 SCSI Cable Kit (2333829 BOM, rev 2)

5.2 Ethernet / Communication Cable Kit (2333831)

FRU = YES

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2318712	1.000	No	VME ETHERNET CABLE, LINUX
2	2318308	1.000	No	ETHERNET CABLE, CHASSIS S
3	2323945	1.000	No	CABLE - MODEM TO HOST COM
4	2318311	1.000	No	ETHERNET CABLE, RJ45 TO R
5	2318713	1.000	No	VME FIBER OPTIC CABLE
7	2319648	1.000	No	FIBER OPTIC CABLE, SC TO
8	2341436	1.000	No	CABLE, VME ETHERNET - OCT

Table 7-25 Communication Cable Kit (2333831 BOM, rev 0)

5.3 Video / Power Cables (2333832)

FRU = YES

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2142225	2.000	No	MONITOR POWER CORD
2	2142226	2.000	No	POWER CORD FOR EMC AREA
3	2304737-2	2.000	No	CABLE, DVI VIDEO CABLE -

Table 7-26 Video/Power Cables (2333832 BOM, rev 0)

5.4 Electrical Components Kit (2333833)

FRU = YES

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
4	46-297335P3	1.000	Yes	16A, 250V 2-STAGE PWR LIN
5	2350404	1.000	Yes	TERMINATOR, CENT50, ACTIV
1	46-269601P6	2.000	No	HIGH DENSITY 50-PIN SLICK
2	46-297332P1	2.000	No	50 OHM, BNC TEE ADAP
3	2154019	2.000	No	50 OHM BNC GROUNDED TERMI

Table 7-27 Electrical Components (2333833 BOM, rev 1)

5.5 Brackets and Straps (2333834)

FRU = YES

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2323984	1.000	No	BRACKET, SCSI TOWER MOUNT
2	2323937-2	1.000	No	STRAP, HOLD DOWN, 1" WIDE
3	2317787	1.000	No	BRACKET, ETHERNET SWITCH
4	2323937	4.000	No	STRAP, HOLD DOWN
5	2323986	2.000	No	BRACKET, HOST COMPUTER
6	2325234	1.000	No	PULL STRAP, HOST COMPUTER
7	2331988	2.000	No	STRAP, HOLDING STRAP

Table 7-28 Brackets & Straps (2333834 BOM, rev 0)

5.6 Miscellaneous PM Hardware (2333835)

FRU = YES

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2227726	1.000	No	MOUSE PAD FOR LIGHTSPEED
2	2325236	2.000	No	LEVELER, CONSOLE ANTI-TIP
3	2283156	2.000	No	RUBBER BOOT, REAR ACCESS
4	46-221513P3	4.000	No	CLAMP OR CLIP
5	46-221513P7	2.000	No	CLAMP OR CLIP
6	2119591-2	1.000	No	25 PIN SUB-D EMI FILTER G
7	2332913	1.000	No	FOAM, POLYETHYLENE, SELF A
8	2119670-2	4.000	No	EMC CONDUCTIVE FOAM TAPE
9	2336907	1.000	No	CONSOLE, HARD RUBBER EDGE

Table 7-29 Miscellaneous PM Hardware (2333835 BOM, rev 0)

5.7 Console Hardware Kit (2333836)

FRU = YES

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2103585	8.000	No	SCREW PAN HEAD 4 MM 8 MM
2	2104490	8.000	No	NYLON WASHER KEYSTONE #33
3	2142006	8.000	No	HEX SOC SCREW M2,5 X 16
4	2260836	8.000	No	SCREW PAN HEAD 3 MM 6 MM
5	2280022	8.000	No	HANDSCREW, MOTOROLA BOARD
6	2324209	8.000	No	SCREW SOCKET HEAD M4 X 60
7	2325236	2.000	No	LEVELER, CONSOLE ANTI-TIP
8	2331393	8.000	No	SCREW, POLYTOP SCREW, M4
9	1000904P475	8.000	No	ID 0.281 OD 0.625 THK .
10	2103585-2	8.000	No	SCREW PAN HEAD 4 MM 10 MM
11	2109867-25	8.000	No	SCREW HEXAGON SOCKET
12	2109873-11	8.000	No	SCREW PAN HEAD 4 MM X 5 M

Table 7-30 Hardware Kit (2333836 BOM, rev 0)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
13	2109877-2	8.000	No	NUT HEXAGON 6 MM
14	2109878-4	4.000	No	WASHER FLAT ISO7089 M10 A
15	2109881-3	4.000	No	WASHER LOCK - INTERNAL 6.
16	2109881-4	4.000	No	WASHER LOCK - INTERNAL 10
17	2119591-2	2.000	No	25 PIN SUB-D EMI FILTER G
18	46-170015P14	8.000	No	PAN HD SCREW 8-32 X 3/8
19	46-170015P23	4.000	No	SCREW, 1/4"-20 X 1/2"
20	46-208935P8	4.000	No	1/4-20, 7/16" HEX X 3/16
21	46-220178P10	2.000	No	GROMMET
22	46-220184P54	4.000	No	SCREW
23	46-221417P4	8.000	No	FEMALE SCREWLOCKS FOR HD
24	46-221449P6	4.000	No	SCREW SIZE 1/4 (0.250) PH
25	46-221852P2	8.000	No	SWITCH
26	46-287187P4	8.000	No	SCREW, CSNK 2.5MM X 6LG
27	46-311804P2	8.000	No	3.2X9X.8 ZP ST DIN
28	46-311805P10	8.000	No	HELICAL LOCK WASHER
29	46-312346P2	8.000	No	NUT, M4 X 0.7 HEX WITH LO
30	46-328421P11	8.000	No	SCREW PAN HEAD 4 MM 5 MM
31	46-328421P3	8.000	No	SCREW PAN HEAD 2.5 MM 5 M
32	46-328421P4	8.000	No	SCREW PAN HEAD 2.5 MM 6 M
33	46-328421P501	8.000	No	SCREW PAN HEAD 3 MM 5 MM
34	46-328422P501	8.000	No	SCREW COUNTERSUNK HEA 5 M
35	46-328431P1	8.000	No	WASHER PLAIN - LARGE 4.3
36	46-328431P2	8.000	No	WASHER PLAIN - LARGE 6.4
37	46-328432P1	8.000	No	WASHER LOCK - SPRING 2.6
38	46-312679P2	4.000	No	INTERNALLY RIBBED GROMMET
39	46-208695P61	2.000	No	HEYCO PART NO.2126 SB-87
40	46-170686P3	1.000	No	LOCTITE 242, 50CC

Table 7-30 Hardware Kit (2333836 BOM, rev 0) (Continued)

5.8 Power Switch & Cables (2333837)

FRU = YES

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2150833	1.000	No	BREAKER SWITCH
2	2108900	1.000	No	CONNECTOR 8 CONTACTS 150
3	2315117	1.000	No	CABLE, PS/2 MOUSE CABLE,
4	2315117-2	1.000	No	CABLE, PS/2 KEYBOARD CABL
5	2318308	1.000	No	ETHERNET CABLE, CHASSIS S
6	2191429	1.000	No	NEXGEN CABLE FROM COMPUTE
7	2341435	1.000	No	CABLE, KEYBOARD - OCTANE
8	2341434	1.000	No	CABLE, MOUSE - OCTANE 2

Table 7-31 Power Switch & Cables (2333837 BOM, rev 0)

5.9 Slide Service Kit (2336912)

FRU = YES

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2336909	1.000	No	COMMUNICATION TRAY SLIDE
2	2336911	1.000	No	RECON / COMPUTER TRAY SLI
3	2317595	1.000	No	SLIDE RAIL FOR DRAWER ASS
4	2317501	1.000	No	BRACKET, DRAWER BALL LATC
5	2317502	1.000	No	LATCH, DRAWER ASSEMBLY

Table 7-32 Slide Service Kit (2336912 BOM, rev 0)

5.10 SCSI Cable Kit, GOC-Linux (2333829-2)

FRU = YES

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2324697-2	1.000	No	SCSI CABLE, VHD68 TO CENT
5	2323942	1.000	No	CABLE, HD50 TO HD50
6	2317488	1.000	No	OCTOPUS CABLE, LINUX COMP
7	2315116	1.000	No	SCSI CABLE, 68 PIN RIBBON
8	2341433	1.000	No	CABLE, OCTOPUS - OCTANE 2
9	2341109	1.000	No	CABLE, SCSI, VHD68 TO HD6
10	2349539	1.000	No	SCSI CABLE, SCAN DATA, H1
11	2368504	1.000	No	CONVERTER, DB25 TO USB

Table 7-33 SCSI Cable Kit, GOC-Linux (2333829-2 BOM, rev 0)

This Page Intentionally Blank.

Section 6.0 — Keyboard w/SCIM Collector (2275758)



Figure 7-29 Keyboard & SCIM Assembly (Rev 11/15/00)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2275756	1.000	Yes	KEYBOARD, NORTH AMERICAN
3	2275752	1.000	Yes	SCIM COLLECTOR - ENGLISH {For further detail, see SCIM Collector – English (2275752) , below.}
4	2323942	1.000	Yes	CABLE, HD50 TO HD50
2	2275757	1.000	No	FUNCTION KEY TEMPLATE - E

Table 7-34 Keyboard w/SCIM Collector - English (2275758 BOM, rev 1)

SCIM Collector – English (2275752)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2262158	1.000	No	SCAN CONTROL INTERCOM MOD
2	2262158-11	1.000	No	OVERLAY, SCIM CT, TILT -
3	2262158-21	1.000	No	OVERLAY, SCIM CT, NO TILT

Table 7-35 SCIM Collector - English (2275752 BOM, rev 0)

End of Chapter