


REVISION	ECR/ECO	DATE	SIGN OFF
3	N/A	3/19/2008	M. ROSE
4	N/A	4/20/2008	M. ROSE
5	2056102	5/10/2008	M. ROSE

GOC6 CONSOLE OUTLET BOX SCHEMATIC

PART NUMBER 5272833

120/208Y VAC OR
120/240 VAC

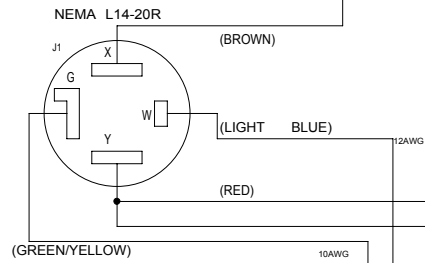


WIRE SIDE SHOWN
ALL WIRES 14AWG UNLESS NOTED

ANY EXPOSED CONDUCTIVE ELEMENT MUST
BE TIED TO CHASSIS GROUND, AND MUST NOT
BE DISCONNECTED WITHOUT THE USE OF A TOOL
"OUT" REFERS TO THE PINS ON THE
OUTSIDE OF THE OUTLET BOX

CIRCUIT BREAKER
SINGLE POLE, SERIES
15 AMPS 50/60HZ

CIRCUIT BREAKER
SINGLE POLE, SERIES
15 AMPS 50/60HZ

FRONT-SWITCH
PANEL

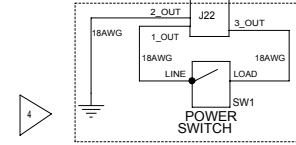

GROUND BAR
GB1

AREA NEAR GROUND BAR MUST HAVE
CENTRAL REFERENCE GROUND SYMBOL

CIRCUIT BREAKER
SINGLE POLE, SERIES
10 AMPS 50/60HZ

CIRCUIT BREAKER MUST BE
ACCESSIBLE TO THE CUSTOMER
MUST BE LABELED "O"

EXTERNAL OUTLETS MUST BE LABELED.
120VAC, 50/60HZ, 8.0 AMPS.
EQUIPOTENTIAL GROUND LABEL FOR STUDS

ACC 1 GROUND

ACC 2 GROUND

PERFORMANCE EMC FILTER
250VAC, 16 AMPS 50/60HZ

PERFORMANCE EMC FILTER
250VAC, 16 AMPS 50/60HZ

IEC320
SHEET
C13
F

GE Healthcare
Company Restricted

Tue Jun 03 10:13:02 2008

Block Path
goc6_console_outlet_box

Block Sheet
1 of 1

Dwg No
5272833SCH

Rev Sheet
5 1 of 1