

Effectivity

The following CT Equipment is affected:

2276776 LightSpeed Plus Table Assembly

Purpose

The purpose of this FMI is to replace table cradle P/N 2115992 with the new table cradle P/N 2115992-2. The new table cradle allows for the installation of a table handle (P/N 2287357).

Related FMIs

FMI 25280: Table Handle Installation

FMI 25283: IV Pole Installation

Furnished Materials

FMI Kit Part Number 2299384 includes the items listed below.

Qty	Description	Part #
1	FMI Document	2299187-100
1	Table Cradle	2115992-2

Required Tools

Qty	Description	Part #
1	3/16" Allen Wrench	

Disposal

U.S. ONLY

- 1.) Return all obsolete parts to:

GEMS Recycling Center
c/o ACE Warehouse, Bldg. 11
2200 College Avenue
Cudahy, WI 53110

- 2.) Destroy and/or dispose of all obsolete software and documentation according to GE policy.

INTERNATIONAL SITES (OUTSIDE THE U.S.)

Destroy and/or dispose of all obsolete hardware, software and documentation according to GE policy.

People Power

Allow 1.0 hour of labor for the installation of the new table cradle.

Procedure

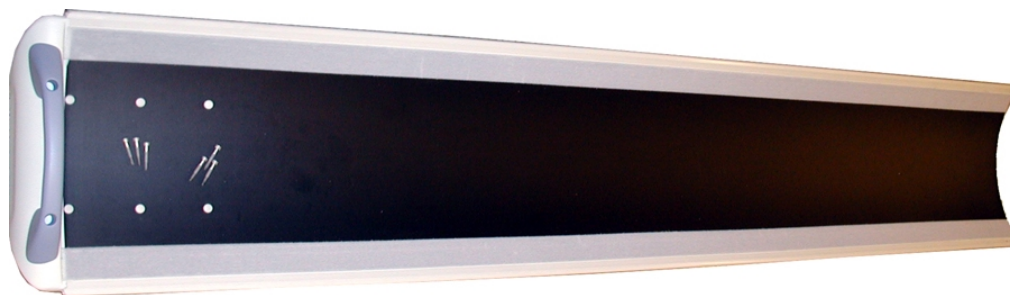


Figure 1 CT Table Cradle (PN 2115992-2) and six mounting screws

- 1.) Raise the table to a comfortable working height.
- 2.) Drive the cradle to the latched, home position, before attempting to remove the assembly.



NOTICE

An unlatched Cradle/carriage assembly could quickly move toward the gantry, and damage the longitudinal encoder assembly.

- 3.) Remove, and keep, the six plug buttons that cover the cradle bolt holes.
- 4.) Loosen and remove the six screws located beneath the plug buttons.
- 5.) Lift the Cradle upward, to remove the assembly from the table.
- 6.) Attach Table Handle to the new Cradle, per Procedure 2 of FMI25280 (Table Handle Installation).
- 7.) If Keeper is missing from new cradle, transfer Keeper from old cradle to new cradle, or install a new one (P/N 2222889-2).
- 8.) Install the new Cradle (with handle attached):
 - a.) Position the rear end of the Cradle over the carriage.
 - b.) Align the holes in the cradle to the threads in the carriage.



Figure 2 Threaded mounting holes in table carriage

- c.) Start, but do not tighten, all six screws.
- d.) Look through the gantry bore at the end of table.
- e.) Laterally position the Cradle on the drive rollers, until the same gap exists between the sloped sides of the Cradle and the guide rollers.
- f.) Tighten the six screws to 13 ft-lbs.
- g.) Replace the plug buttons.

End of Document