

[new...](#) | [new revision](#) | [page status](#) | [copy](#) | [move](#) | [delete](#)

welcome!
ct tech pubs web page
ct water cooler
ct tech pubs updates
surveys
9800/hla
highspeed family (us)
ct/i
lightspeed family
yms family
all ct products
index
old ct emails
communiques pending
approval
customize
members

search:

advanced search
news: daily | weekly

go offline | chat | notify |
print | tutorial | help

home> lightspeed family

Power Pan Circuit Breaker

christopher hardiman (christopher c hardiman), 07/11/2002 - 09:37 am



PROPRIETARY WARNING

This information is Property of GE - For GE Service Personnel Only.
Disclosure to third parties prohibited.

Communique Number

MQUE-5BXK3D

Equipment

LightSpeed Plus, Ultra & 16

Hints/Tips

A circuit breaker, CB1, exists in the gantry power pan that will shut off 120 VAC to the table. Up to now, documentation of this circuit breaker has been lacking. This will be corrected in the next release of the schematics manuals for these systems. Please refer to the attached diagrams for details on the location of this device and the updated schematic for the power pan.

You can download any of the files below by dragging them to your Windows desktop.
To open a file for reading, double-click it.




Author

Christopher C Hardiman

Revision

07/11/2002 09:37 AM


[new...](#) | [new revision](#) | [page status](#) | [copy](#) | [move](#) | [delete](#)

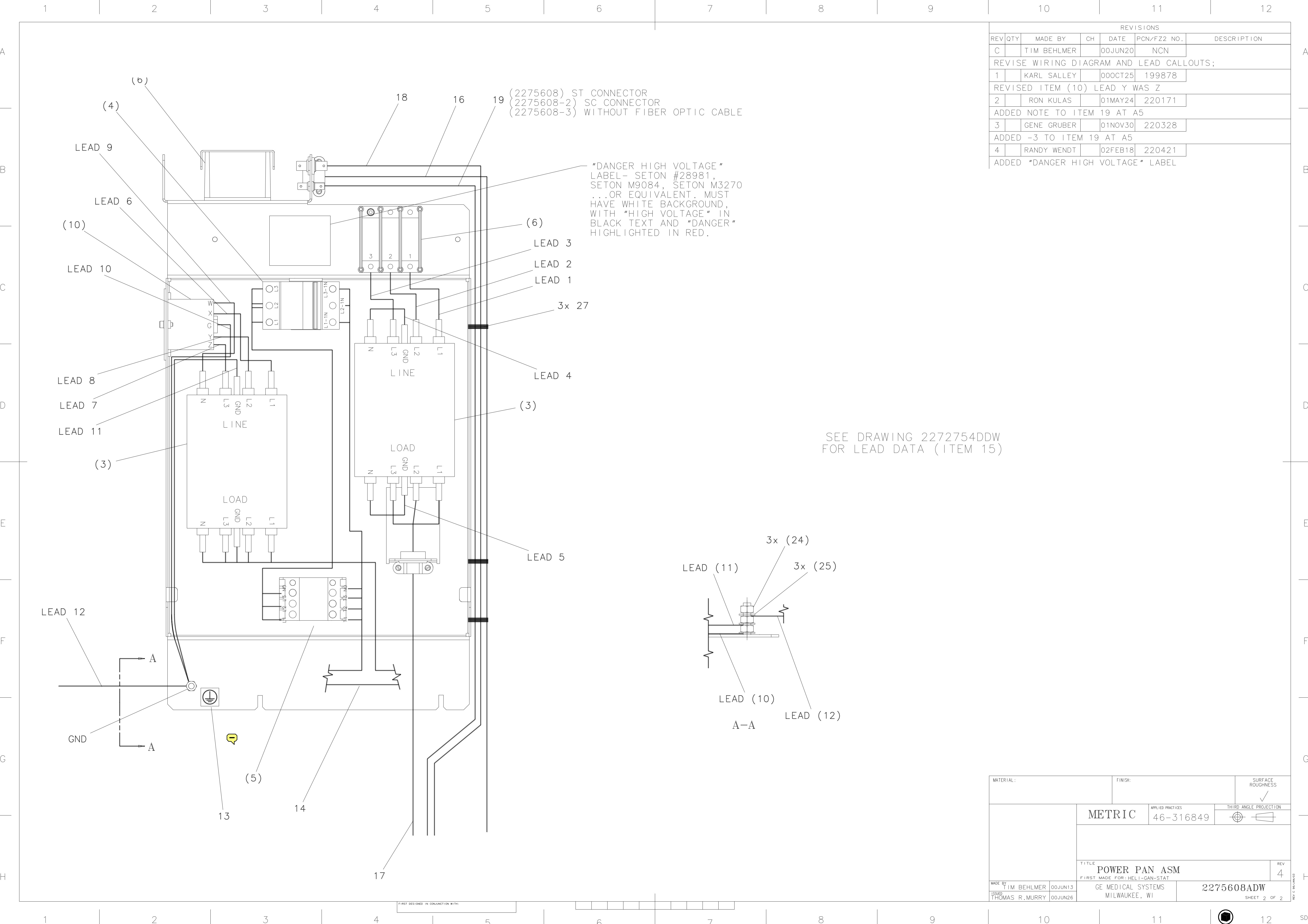

REVISIONS						
REV	QTY	MADE BY	CH	DATE	PCN/FZ2 NO.	DESCRIPTION
C		TIM BEHLMER		00JUN20	NCN	REVISE WIRING DIAGRAM AND LEAD CALLOUTS;
1		KARL SALLEY		00OCT25	199878	REVISED ITEM (10) LEAD Y WAS Z
2		RON KULAS		01MAY24	220171	ADDED NOTE TO ITEM 19 AT A5
3		GENE GRUBER		01NOV30	220328	ADDED -3 TO ITEM 19 AT A5
4		RANDY WENDT		02FEB18	220421	ADDED "DANGER HIGH VOLTAGE" LABEL

"DANGER HIGH VOLTAGE"
LABEL- SETON #28981,
SETON M9084, SETON M3270
...OR EQUIVALENT. MUST
HAVE WHITE BACKGROUND,
WITH "HIGH VOLTAGE" IN
BLACK TEXT AND "DANGER"
HIGHLIGHTED IN RED.

SEE DRAWING 2272754DDW
FOR LEAD DATA (ITEM 15)

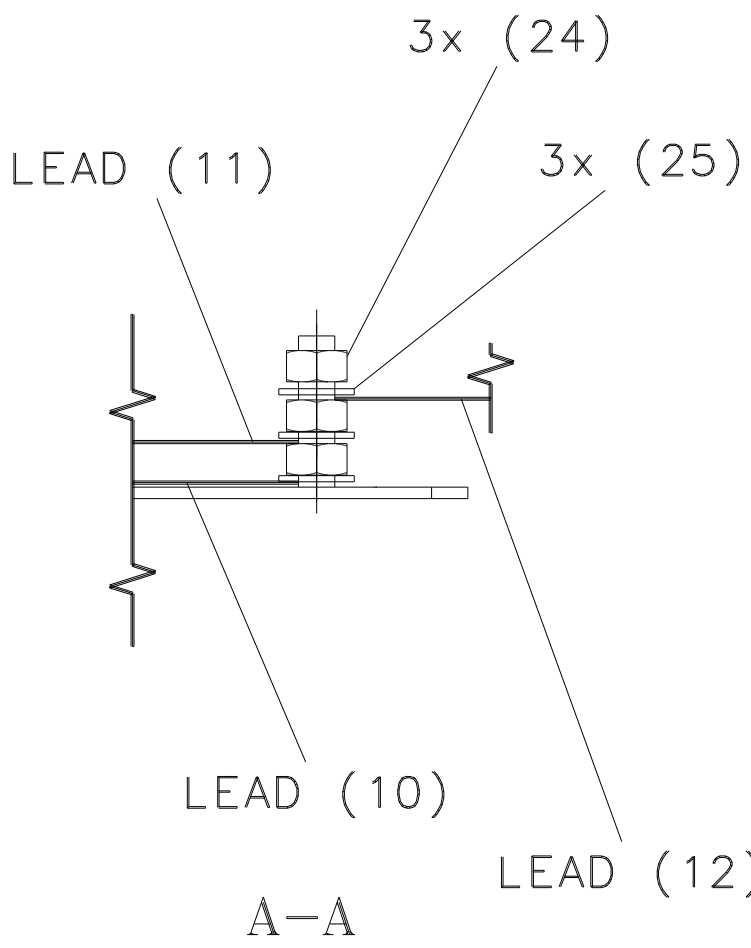


MATERIAL :			FINISH:		SURFACE ROUGHNESS ✓		
			METRIC		APPLIED PRACTICES 46-316849	THIRD ANGLE PROJECTION 	
			TITLE POWER PAN ASM FIRST MADE FOR: HELI-GAN-STAT				
MADE BY T. IM BEHLMER		00JUN13	GE MEDICAL SYSTEMS MILWAUKEE, WI			2275608ADW	
ISSUED THOMAS R. MURRY		00JUN26				SHEET 2 OF 2	


FIRST DESIGNED IN CONJUNCTION WITH:

REV C 06/20/2023

SDRC
A1

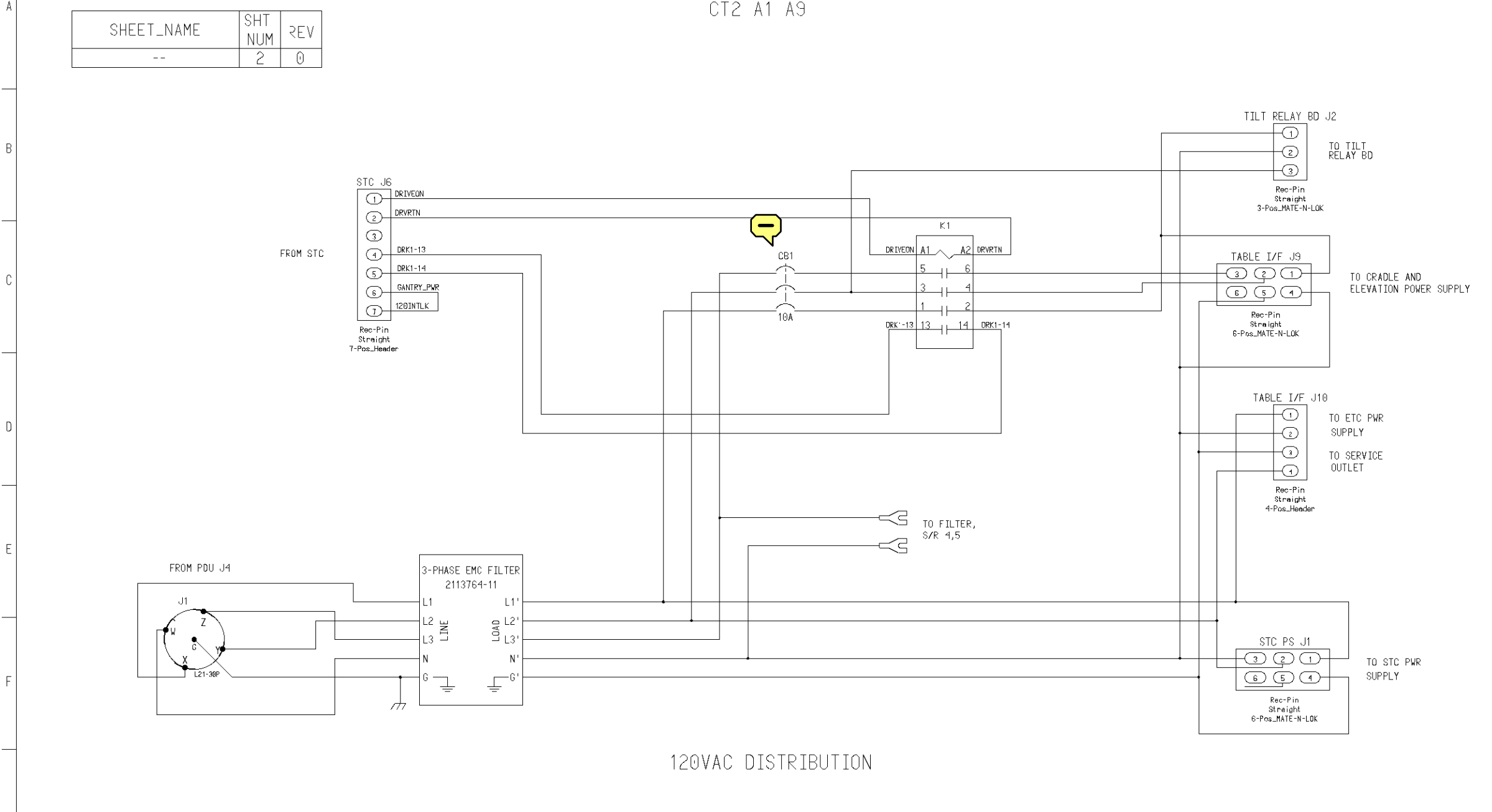

BLOCK PATHNAME \$GEMS_CT_XRAY_GEN/ASM_DOCS/2275608_ASM/H2_GAN_PWR_PAN					BLOCK SHT 1 OF 2		DESIGN TITLE H2 GANTRY POWER PAN				
REV	PCN	DESCRIPTION	DATE	APPR BY	DATE	MADE	FIRST MADE FOR HELI-GAN-STAT				
MADE BY						CHECK					
REV	PCN	DESCRIPTION	DATE	APPR BY	DATE	APPR	DRAWING NO. 2275608SCH				
MADE BY						OTHER	REV 0				
							SIZE 8		SHT/CONT. ON 1 / 2		

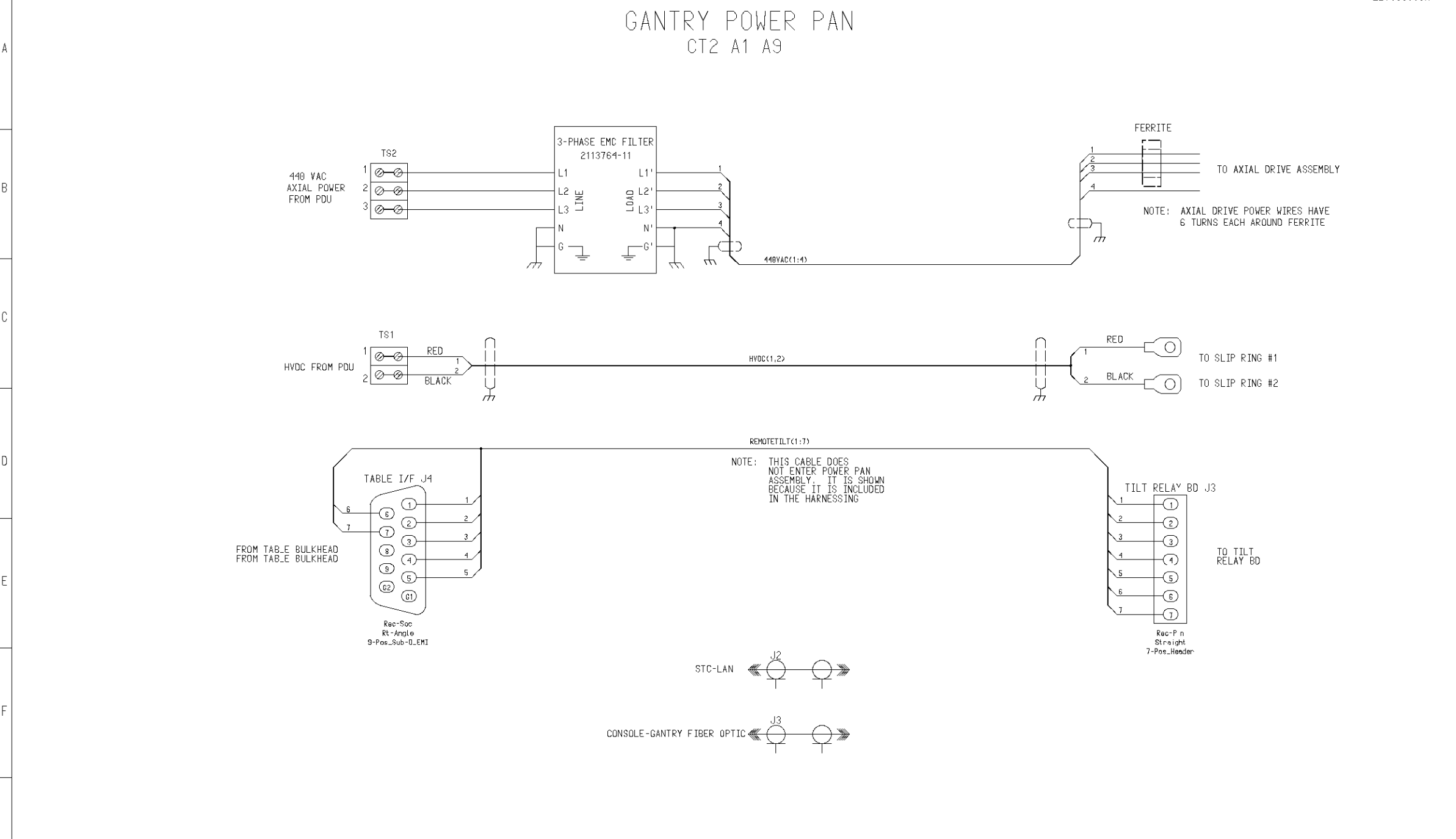

BLOCK PATHNAME \$GEMS_CT_XRAY_GEN/ASM_DOCS/2275608_ASM/H2_GAN_PWR_PAN					BLOCK SHT 2 OF 2		DESIGN TITLE H2 GANTRY POWER PAN			
REV PCN	DESCRIPTION	APPR BY	DATE	MADE R. VANDEN BOOGARD	DATE 04/JAN/2001	GE MEDICAL SYSTEMS MILWAUKEE, WI USA	FIRST MADE FOR HELI-GAN-STAT			
MADE BY	DATE			CHECK J. R. SCHMIDT	DATE 07/MAR/2001					
REV PCN	DESCRIPTION	APPR BY	DATE	APPR J. R. SCHMIDT	DATE 07/MAR/2001	DRAWING NO. 2275608SCH	REV 0	SIZE 8	SHT/CONT. ON 2 / -	
MADE BY	DATE			OTHER	DATE					