

MAGNETOM

System

Software-Handling Instructions

Software Installation syngo MR XA20A

M11-030.816.04.01.02 | 11.19

1 Software installation syngo MR XA20

1.1 Prerequisite

i

During the system shutdown an "Export package" with site specific data is created and stored at the service partition of the MaRS. To avoid loss of data before re-installing software **it's strongly recommended to get the site specific data from the MaRS service partition** if the MaRS is still accessible.

Mars is connected but powered off

Harddisks are empty at the MRAWP

Use the Blancco erasure tool [↗ Blancco link](#)

i

If the BIOS password was changed by local IT or CSE and the password is unknown, then what will happen is:

It is not possible to change the boot order or selecting a boot device for software installation or disaster recovery.

Workaround: Open the workstation and resetting the BIOS, afterwards the default BIOS settings will apply.

1.2 Installation procedure MRAWP

Shut down the MRAWP if it is not already turned off.

Plug in the SD card with the to-be-installed Numaris/X package

Turn on the MRAWP machine. Hold Esc on boot up.

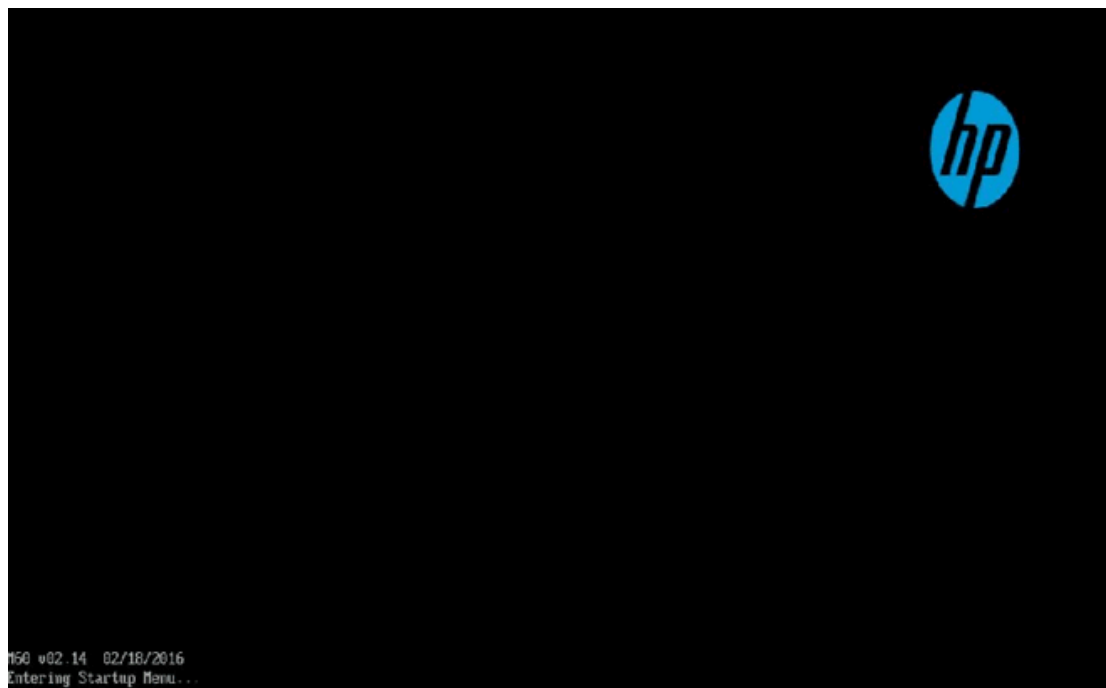


Fig. 1: Boot screen MRAWP

In the startup menu select Boot Menu.

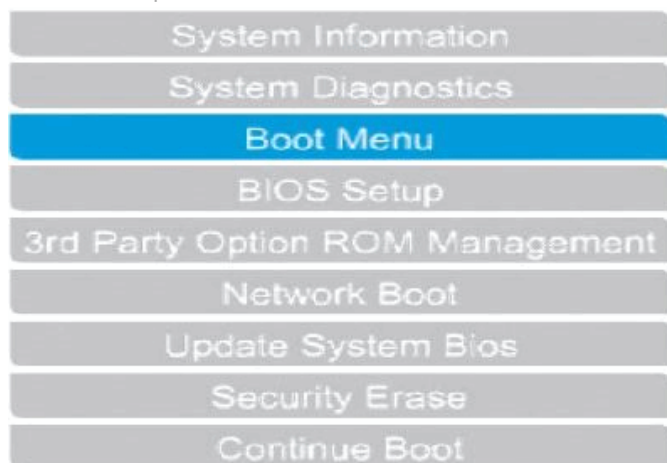


Fig. 2: Startup menu MRAWP

You will be prompted for a password.

In the boot menu select the SD card with the Numaris/X package. In this particular example the package is on a USB memory stick (SanDisk Cruzer Switch). Adapt to your scenario!

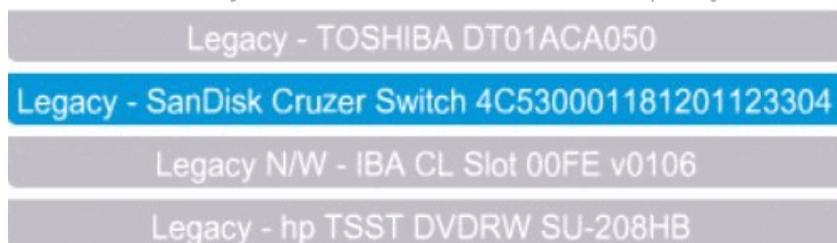


Fig. 3: Boot menu MRAWP

After selecting the correct storage device (SD-Card) the MR Deployment OS will boot up.

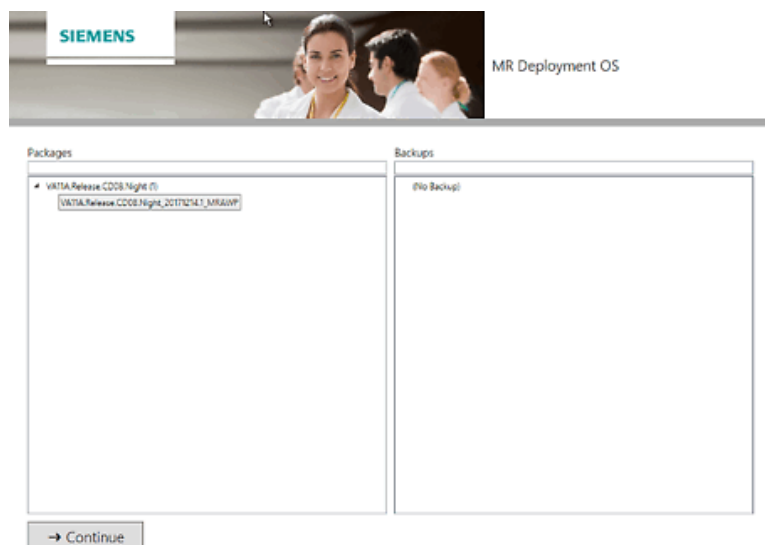


Fig. 4: Package selection screen

There will be one available package that can be installed. Select that package and press ⇒ Continue. Next you'll have to configure the MRAWP. Insert the correct hostname, serialnumber, system type and Network configuration.

Fig. 5: Deployment Portal example for a Vida configuration

Press ⇒ Continue . The installation will now begin.

i

The elapsed time values for total and between two installation images are only a hint. The installation process runs unattended and will take somewhere around 2 hours.



Fig. 6: Preparing system for the new package

Elapsed time since previous image: 30 min

Total elapsed time: 30 min

Do nothing if you see the Windows lock screen. The system will reboot automatically. After that the creation of native images will begin .

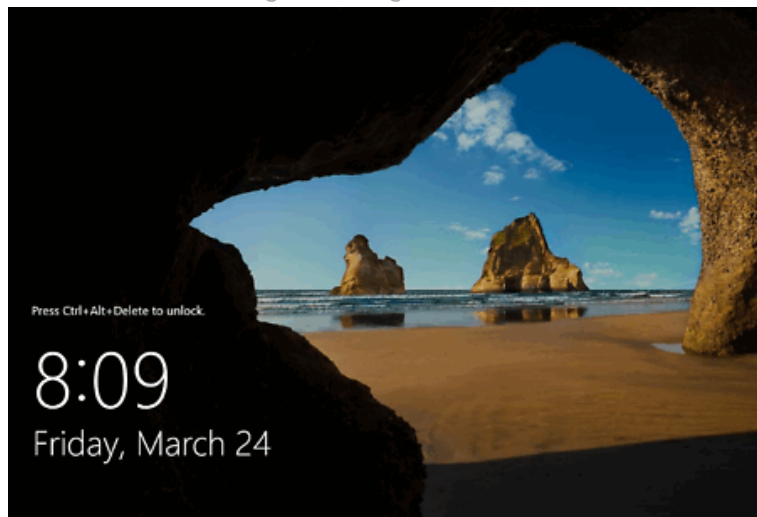


Fig. 7: Windows lock screen

Elapsed time since previous image: 45 min

Total elapsed time: 75 min

After that the Device Guard configuration and scan will start

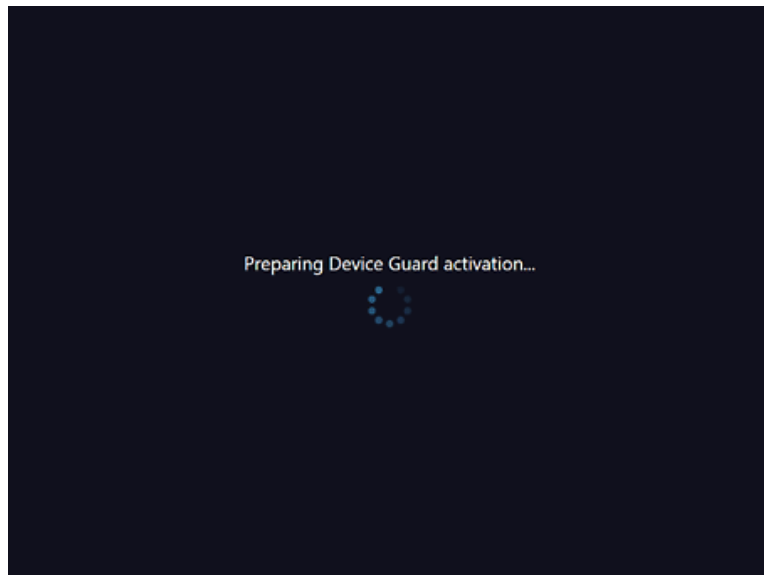


Fig. 8: Device Guard setup

i

At this step remove now the lan cable for the hospital network on the MRAWP.

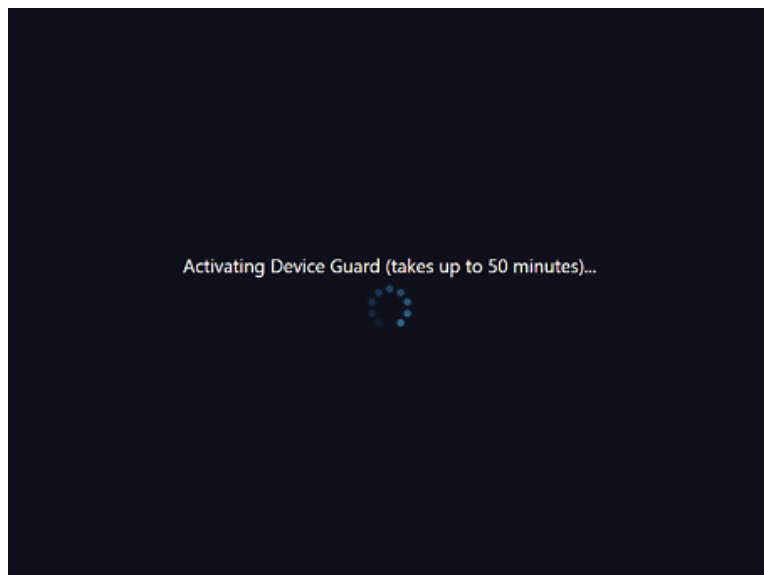


Fig. 9: Activating Device Guard

Elapsed time since previous image: 50 min

Total elapsed time: 125 min

The system will reboot twice after the Device Guard activation. Wait until the Bringup window appears than the installation is almost complete.

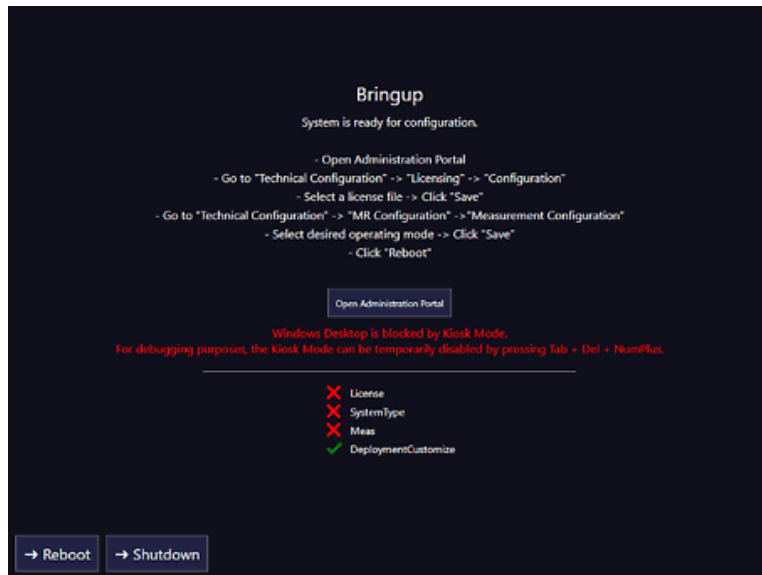


Fig. 10: Bringup screen

Reconnect the lan cable of the hospital network to the MRAWP
Windows 10 activation see [▶ Activating Windows 10 licensing](#)

1.3 Installation procedure MRWP

Shut down the MRWP if it is not already turned off.
Plug in the SD card with the to-be-installed Numaris/X package
Turn on the MRWP machine. Hold Esc on boot up.

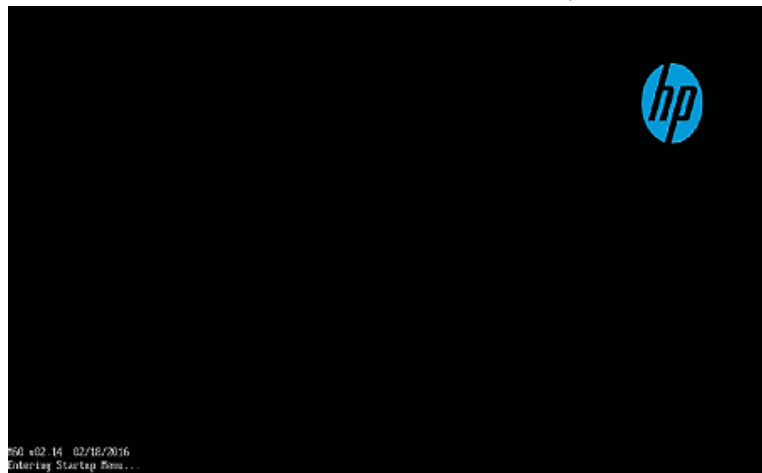


Fig. 11: Boot screen MRWP

In the startup menu select Boot Menu.