

GE Healthcare

LightSpeed 7.X Illustrated Parts List



OPERATING DOCUMENTATION



5116416-100
Rev 12

Chapter 2

GT Table

Section 1.0 — GT 2000 Table (5121647)

FRU = NO

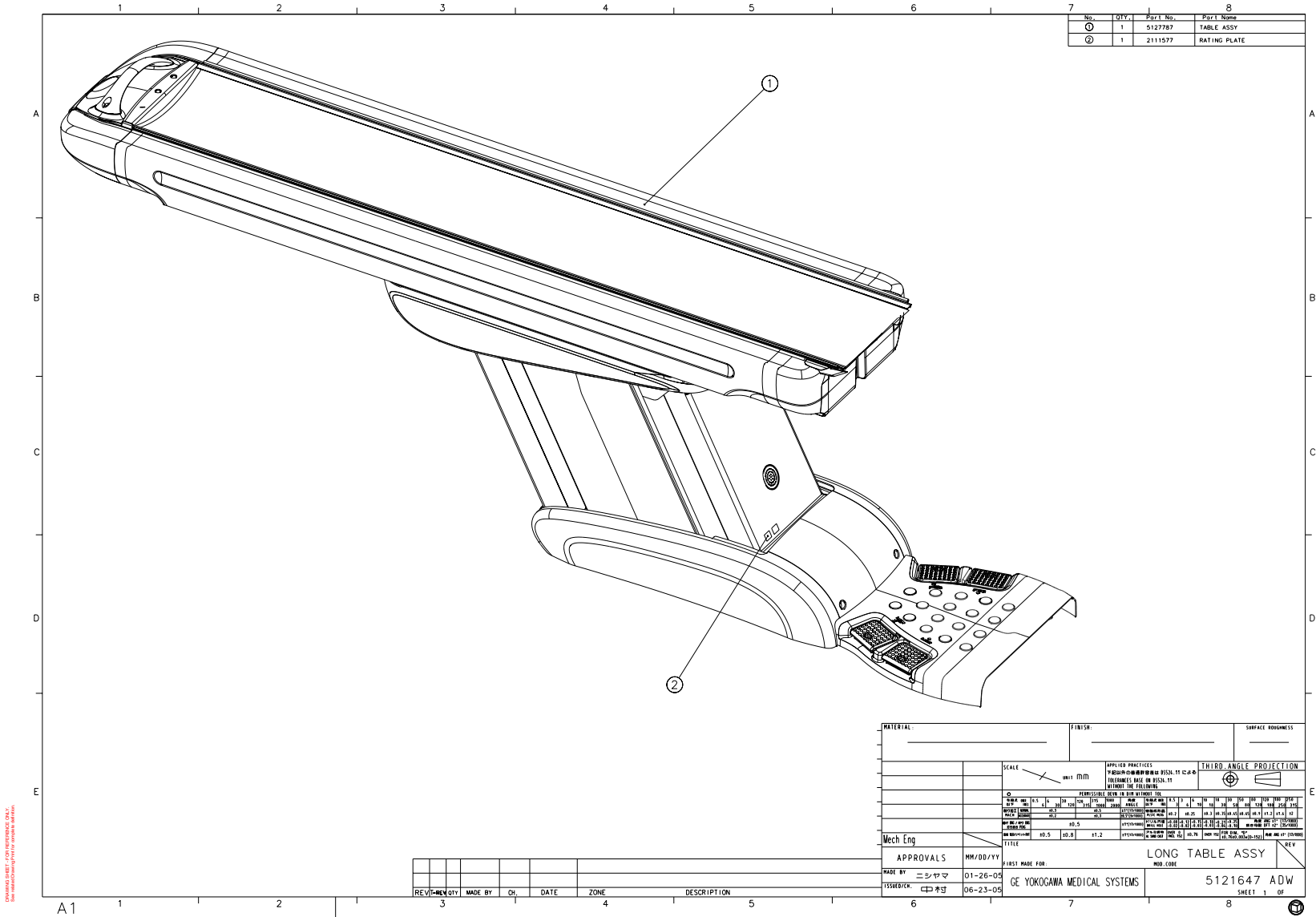


Figure 2-1 Positioning GT LONG TABLE ASSY (5121647ADW, sheet 1, rev 4)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
	5143974	1.0	Yes	GTCN SET
	5168472-2	1.0	Yes	IMS DRIVE REPLACEMENT SET Positioning GT
	5192989	1.0	Yes	IMS MOTOR ASSY Collector
	5237567	1.0	Yes	Lubricant U0303BD-2
2	46-303268P1	1.0	Yes	PRODUCT LOCATOR CARD FORM FOR RATING PLATES, PRINTED BY DUPLITECH
1	5127787	1.0	No	TABLE ASSY {For further detail, see Section 1.1, beginning on page 32. }
3	5129558	1.0	No	SHIP COLLECTOR LONG TABLE Positioning GT {For further detail, see Section 1.2, beginning on page 110. }

Table 2-1 Positioning GT LONG TABLE ASSY (5121647 BOM, Rev 18)

$FRU = NO$

Figure 2-2 Table Assembly (5127787ADW, sheet 2, rev 10)

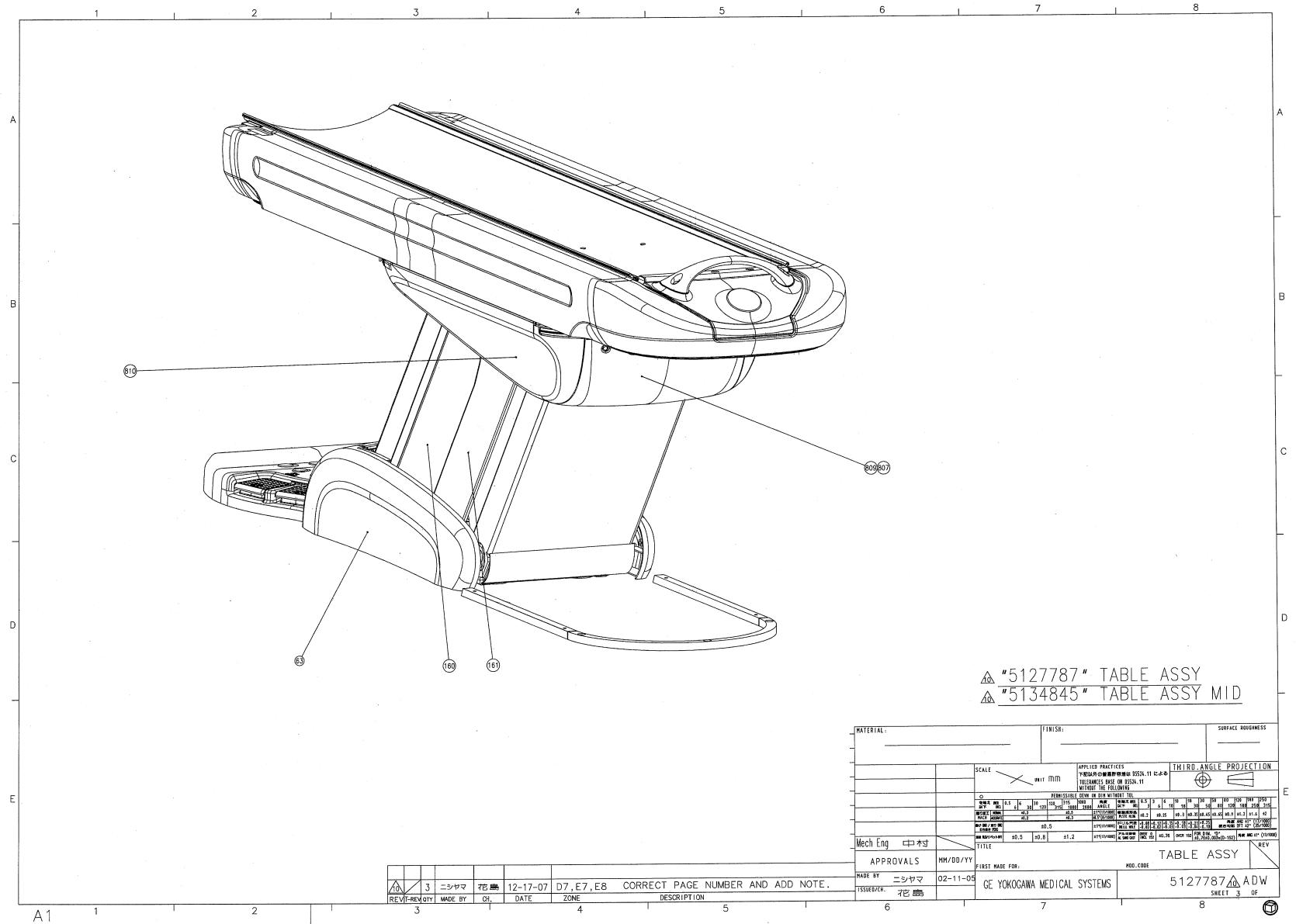
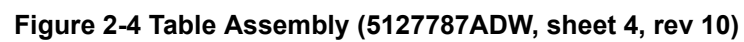
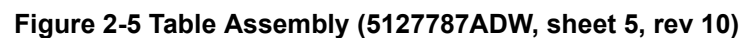
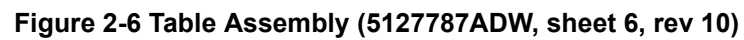


Figure 2-3 Table Assembly (5127787ADW, sheet 3, rev 10)







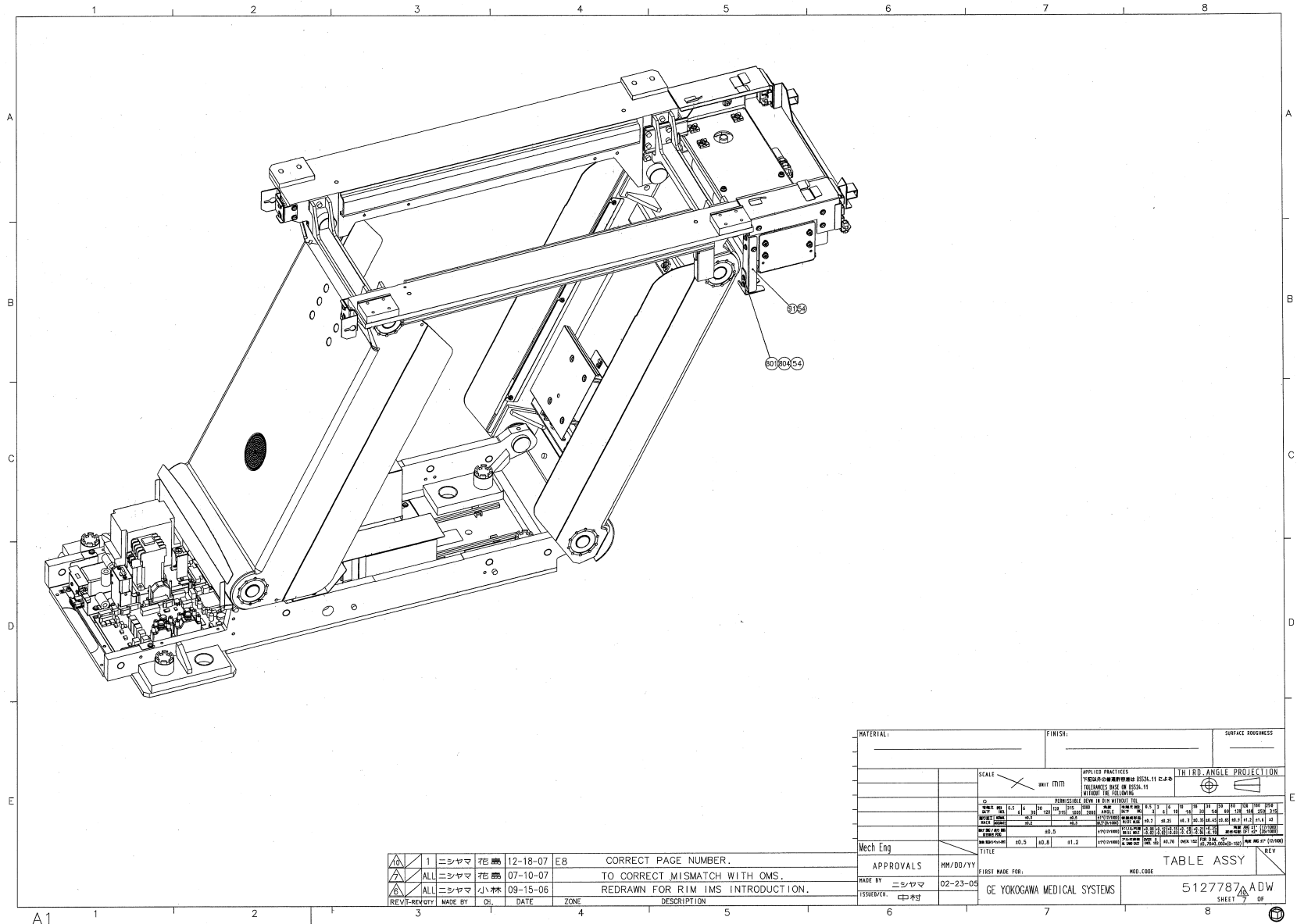


Figure 2-7 Table Assembly (5127787ADW, sheet 7, rev 10)

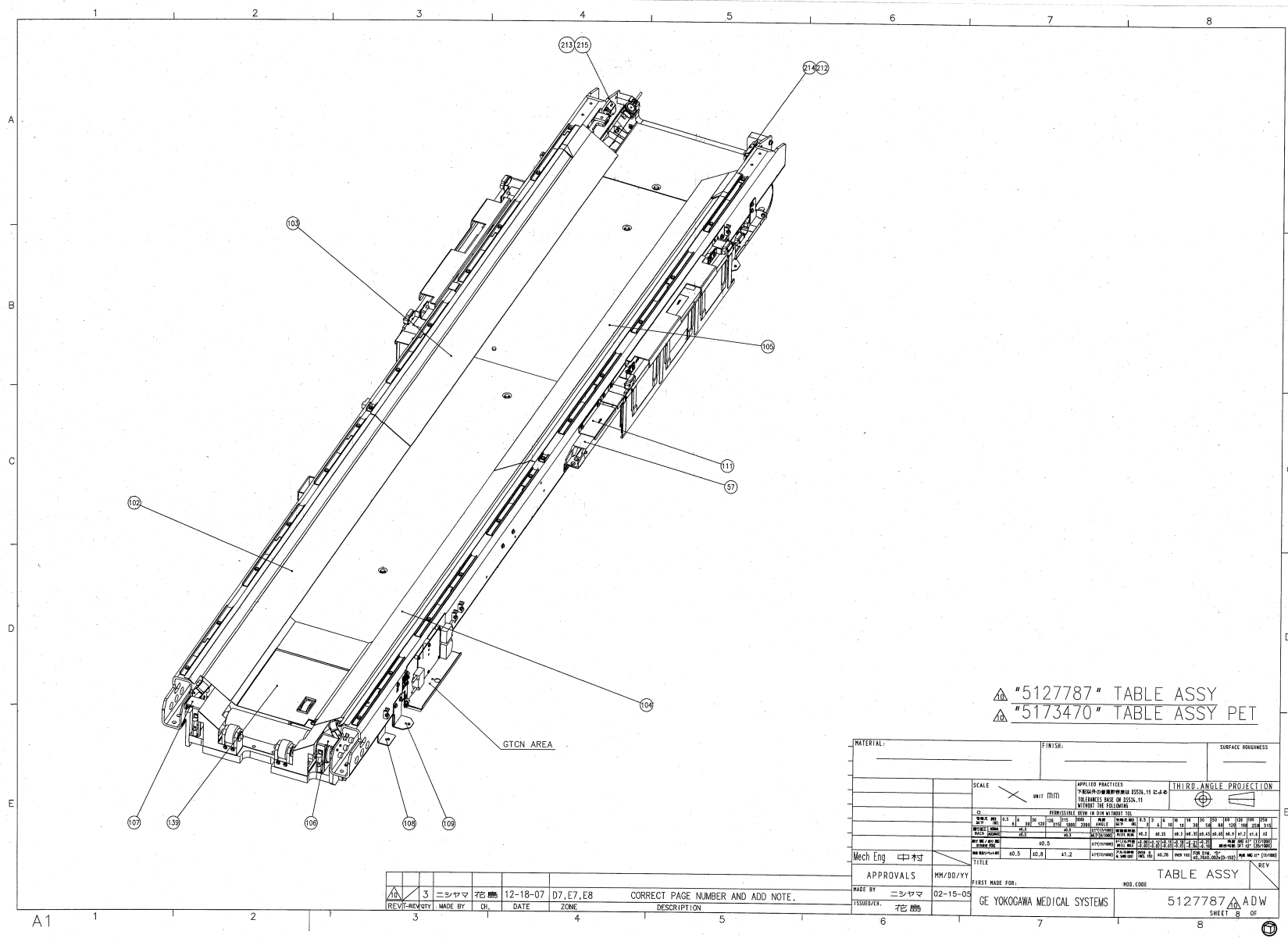


Figure 2-8 Table Assembly (5127787ADW, sheet 8, rev 10)

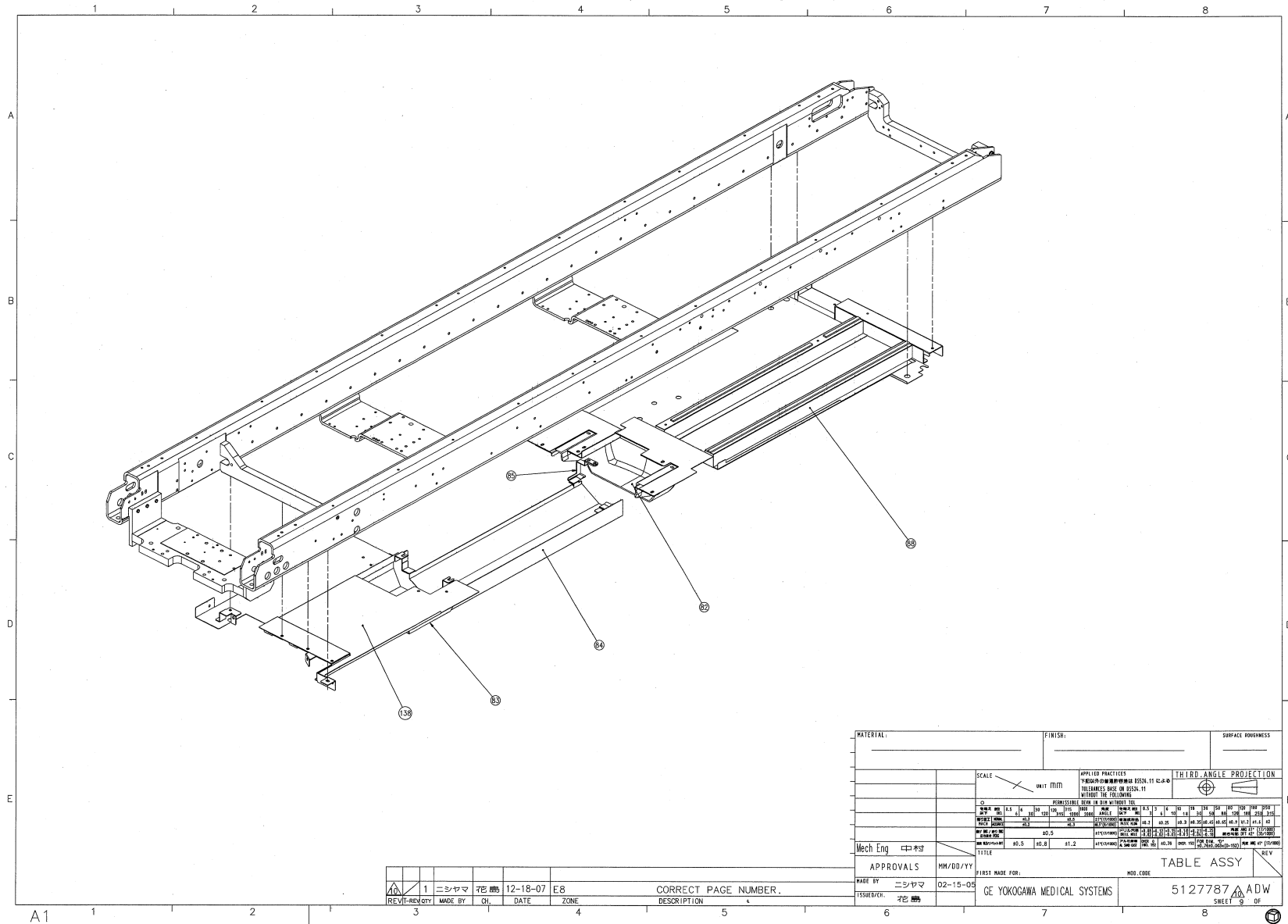


Figure 2-9 Table Assembly (5127787ADW, sheet 9, rev 10)

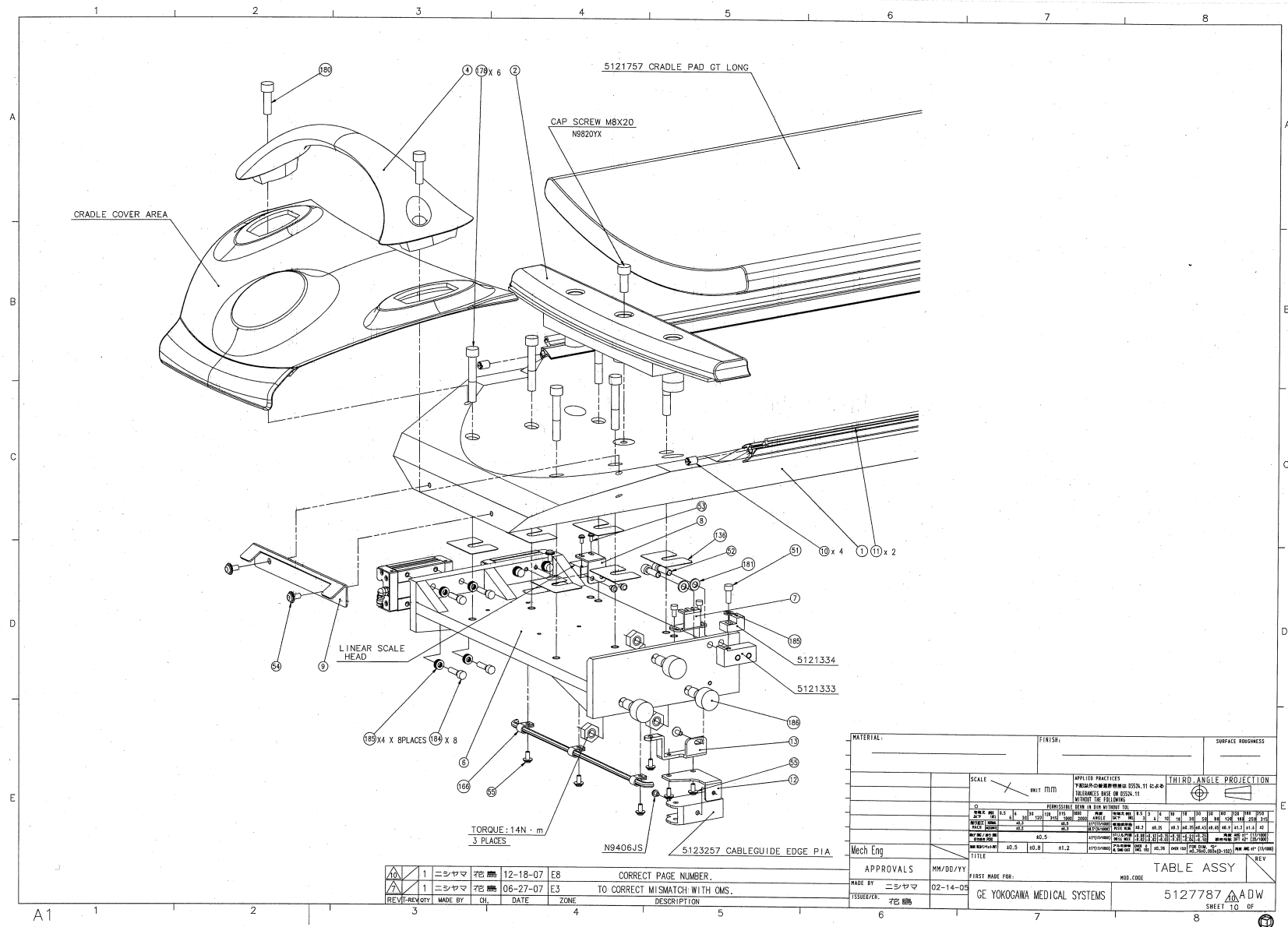


Figure 2-10 Table Assembly (5127787ADW, sheet 10, rev 10)

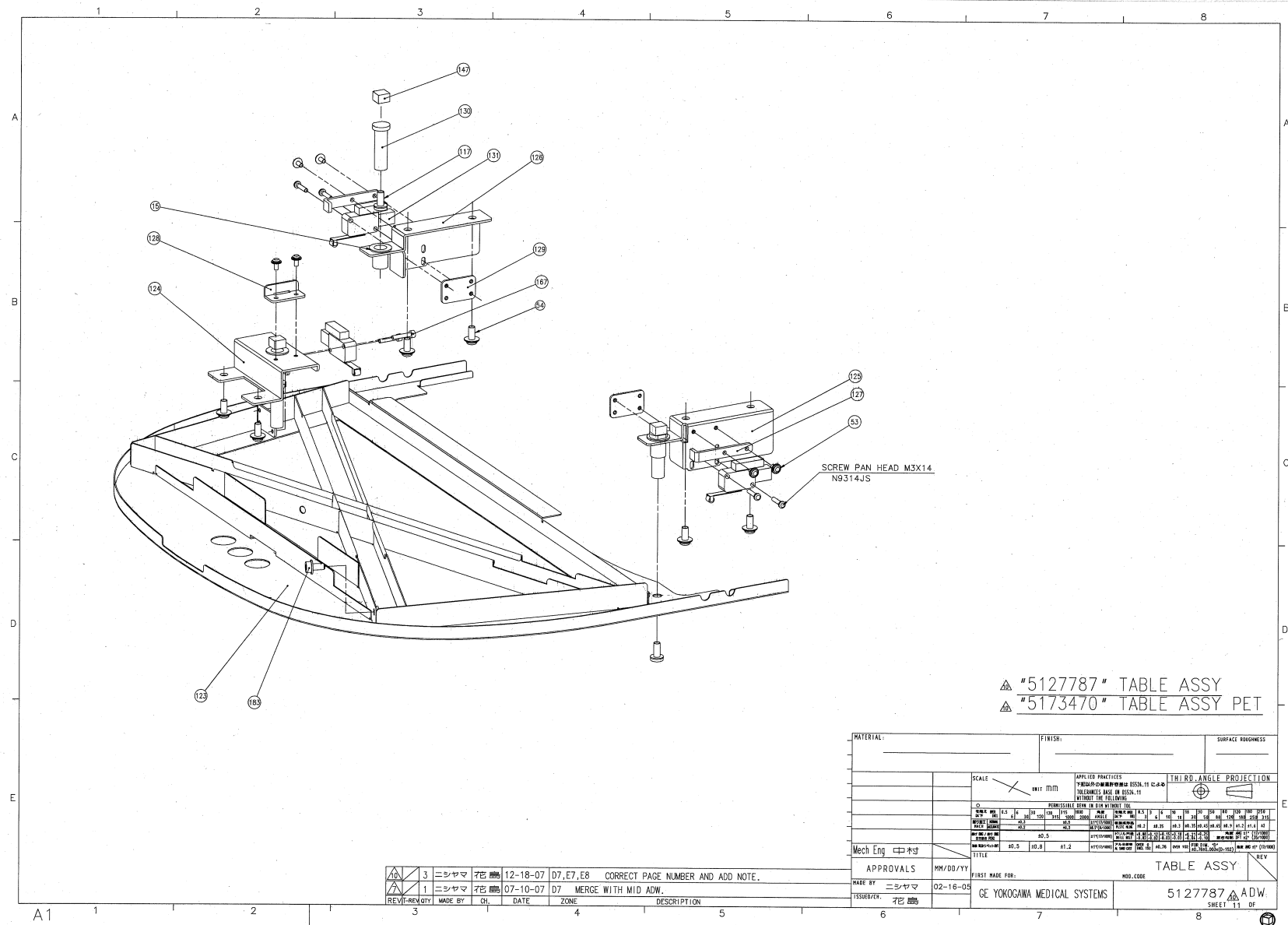


Figure 2-11 Table Assembly (5127787ADW, sheet 11, rev 10)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
	5127282	1.0	Yes	HYDRAULIC CYLINDER ASSY Positoning GT
	5127283	1.0	Yes	HYDRAULIC RETURN HOSE ASSY Positoning GT
	5127284	1.0	Yes	PROPORTIOAL VALVE ASSY Positoning GT
	5127285	1.0	Yes	SHUT OFF VALVE ASSY Positoning GT
	5127286	1.0	Yes	HYDRAULIC HIGH PRESUURE HOSE Positoning GT
	5127287	1.0	Yes	ELEVATION HYDRAULIC PUMP Positoning GT
	5138908	1.0	Yes	TABLE SUPPORT JIG ASSY
	5164938	1.0	Yes	TAPE SWITH 4 ASSY Positioning GT
	5166082-2	1.0	Yes	RIM BASE COVER LH REPLACEMENT SET positioning GT
	5166242-2	1.0	Yes	RIM BASE COVER RH REPLACEMENT SET positioning GT
	5181038	1.0	Yes	TAPE SW-5 Positioning GT
	5191991	1.0	Yes	RIM IMS COVER SET positioning GT
	5220568	1.0	Yes	INJECTOR GUARD KIT Positioning GT
1	5173098	1.0	Yes	CRADLE LONG ASSY Positioning GT
10	5123613	4.0	Yes	PLASTIC SET SCREW Positioning GT
11	5120311	2.0	Yes	ATTACHMENT-RAIL Positioning GT
25	5117286	1.0	Yes	GTCN J7 TO J221 Positioning GT
26	5114080	1.0	Yes	PE Positioning GT
27	5114118	1.0	Yes	GTCBJ3 TO J104 Positioning GT
28	5114120	1.0	Yes	GTCBJ10 TO SPEAKER Positioning GT
29	5114122	1.0	Yes	GTCBJ11 TO J200 Positioning GT
30	5114124	1.0	Yes	GTCBJ12 TO J201 Positioning GT
31	5114125	1.0	Yes	GTCBJ6 TO IMSDRIVER Positioning GT
32	5114126	1.0	Yes	GTCBJ5 TO CRADLEDRIVER Positioning GT
33	5114128	1.0	Yes	CRADLEDRIVER TO J204
34	5114129	1.0	Yes	IMSDRIVER TO J203 Positioning GT
35	5114135	1.0	Yes	GTCN TO J207 Positioning GT

Table 2-2 Table Assembly (5127787 BOM, Rev 35)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
36	5114146	1.0	Yes	GTCN TO J218 Positioning GT
37	5114148	1.0	Yes	GTCN TO LINERSCALE Positioning GT
38	5118961	1.0	Yes	LIMIT SW CABLE Positioning GT
39	5119130	1.0	Yes	J102 TO J112 Positioning GT
40	5120569	2.0	Yes	4PIN470DUM Positioning GT
41	5120570	1.0	Yes	6PIN470DUM Positioning GT
42	5120572	1.0	Yes	9PIN470DUM Positioning GT
43	5120574	1.0	Yes	TS_TOP_FRONT Positioning GT
44	5120577	1.0	Yes	TS_IMS_RIGHT Positioning GT
45	5120579	1.0	Yes	TS_IMS_LEFT Positioning GT
46	5120581	1.0	Yes	TS_IMS_REAR Positioning GT
47	5120773-2	1.0	Yes	GTCB TO GTTS Cable For EMC Positioning GT
48	5114146-2	1.0	Yes	LIMIT SW Positioning GT
50	N9408ZX	4.0	Yes	SCREW HEXAGON SOCKET 4 MM 8 MM CAP SCREW 4X8 FE
52	N9620YX	2.0	Yes	SCREW HEXAGON SOCKET 6 MM 20 MM CAP SCREW 6X20 FE
53	N9306HR	2.0	Yes	SCREW TP 3X6 FE ZN
54	N9510HR	120.0	Yes	SCREW TP
55	N9408HR	3.0	Yes	SCREW TP
57	5127780-2	1.0	Yes	IMS CABLE BEAR ASSY
62	5127458-2	1.0	Yes	BASE COVER CTR ASSY positioning GT
69	5127465-3	1.0	Yes	RIM TOP COVER FRONT ASSY positioning GT {For further detail, see Section 1.1.4, beginning on page 90. }
70	5154872	1.0	Yes	RIM TOP COVER REAR ASSY positioning GT
71	5127467	1.0	Yes	TOP COVER SIDE L ASSY {For further detail, see Section 1.1.6, beginning on page 94. }
72	5127468	1.0	Yes	TOP COVER SIDE R ASSY {For further detail, see Section 1.1.7, beginning on page 96. }
86	5113088	1.0	Yes	COVER-PLATE-L positioning GT
87	5113091	1.0	Yes	COVER-PLATE-R positioning GT

Table 2-2 Table Assembly (5127787 BOM, Rev 35) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
102	5140899	1.0	Yes	RAIL COVER 2000 RF Positioning GT
103	5140900	1.0	Yes	RAIL COVER 2000 RR Positioning GT
104	5140901	1.0	Yes	RAIL COVER 2000 LF Positioning GT
105	5140902	1.0	Yes	RAIL COVER 2000 LR Positioning GT
109	5113617	4.0	Yes	TOP COVER SIDE BRKT 11 positioning GT
119	N9612YX	16.0	Yes	SCREW HEXAGON SOCKET 6 MM 12 MM CAP SCREW 6X12 FE
121	5127469	1.0	Yes	TOUCH SENSOR L 2000 ASSY Positioning GT
122	5127470	1.0	Yes	TOUCH SENSOR R 2000 ASSY Positioning GT
123	5119698	1.0	Yes	TOUCH PLATE-REAR positioning GT
131	5125112	1.0	Yes	GT REAR TOUCH SWITCH CABLE Positioning GT
146	5219979	1.0	Yes	CRADLE COVER ASSY WITH LOGO LABEL Positioning GT of FRU
148	U0442BA	10.0	Yes	TUBE 3X0.25 BLACK
150	U0031AC	9.0	Yes	CLAMP
151	U0032AC	2.0	Yes	CLAMP OLD CLAMP
160	5130814	1.0	Yes	SIDE-PANEL-FRONT-LEFT-ASSY, Positioning, GT {For further detail, see Section 1.1.8, beginning on page 98. }
161	5130815	1.0	Yes	SIDE-PANEL-REAR-LEFT-ASSY, Positioning, GT {For further detail, see Section 1.1.9, beginning on page 100. }
162	5130816	1.0	Yes	SIDE-PANEL-FRONT-RIGHT-ASSY, Positioning, GT {For further detail, see Section 1.1.10, beginning on page 102. }
163	5130817	1.0	Yes	SIDE-PANEL-REAR-RIGHT-ASSY, Positioning, GT {For further detail, see Section 1.1.11, beginning on page 104. }
164	U0034BA	150.0	Yes	TAPE #5421 16Mug
165	U0047BA	80.0	Yes	TAPE
166	U0001AC	4.0	Yes	CLAMP NK-4N CLAMP
173	N9320JS	4.0	Yes	SCREW PAN HEAD 3 MM 20 MM SCREW osX 3X20 FE ZN
174	N9405JS	1.0	Yes	SCREW NABE 4X5 Fe Zn
176	N9508JS	19.0	Yes	SCREW NABE 5X8 FE ZN
180	N9825YX	2.0	Yes	SCREW CAP SCREW 8X25 Fe

Table 2-2 Table Assembly (5127787 BOM, Rev 35) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
181	N9601WS	4.0	Yes	WASHER PLAIN - LARGE 6.4 MM 12.5 MM WASHER 6.4X12.5X1.6FEZN
182	5126552	1.0	Yes	GND WIRE COVER Positioning GT
184	N9520ZX	8.0	Yes	SCREW HEXAGON SOCKET 5 MM 20 MM CAP SCREW 5X20 FE
185	N9501WS	32.0	Yes	WASHER PLAIN - LARGE 5.3 MM 10 MM WASHER M5
302	U0405BD	0.313	Yes	ADHESIVE YhjbN 271(af)
304	2221887	0.00625	Yes	U0427BD ADHESIVE
307	N9306ES	21.0	Yes	SCREW COUNTERSUNK HEA 3 MM 6 MM
308	N9306ZX	3.0	Yes	SCREW HEXAGON SOCKET 3 MM 6 MM CAP SCREW 3X6 FE
309	N9408JS	23.0	Yes	SCREW NABE 4X8 Fe Zn
310	N9414JS	2.0	Yes	SCREW PAN HEAD 4 MM 14 MM SCREW osX 4X14 FE ZN
311	N9510JS	8.0	Yes	SCREW NABE 5X10 Fe Zn
313	U0201AC	26.0	Yes	CLAMP CV_100
314	U0202AC	20.0	Yes	CLAMP CV_150
315	U0203AC	30.0	Yes	CLAMP CV_200
316	U0208AC	5.0	Yes	CLAMP TC-142 TAIRATSUPU KOTEI
317	U0212AC	24.0	Yes	CLAMP TMS_20
324	U0506AC	45.0	Yes	PROTECTOR EDGING CE-012
325	U0509AC	0.35	Yes	PROTECTOR CE_032
809	5142244	1.0	Yes	IMS-COVER-CTR-RIM positioning GT
810	5191722	1.0	Yes	RIM IMS COVER SIDE LH ASSY positioning GT {For further detail, see Section 1.1.12, beginning on page 106. }
811	5191736	1.0	Yes	RIM IMS COVER SIDE RH ASSY positioning GT {For further detail, see Section 1.1.13, beginning on page 108. }
2	5112083	1.0	No	CRADLE ATT BASE Positioning GT
4	5120588	1.0	No	cradle handle positioning GT
6	5112214	1.0	No	CRADLE-CARRIAGE Positioning GT
7	5112420	1.0	No	CARRIAGE WIRE HOLDER Positioning GT
8	5112424	1.0	No	CARRIAGE SENSOR SUP Positioning GT

Table 2-2 Table Assembly (5127787 BOM, Rev 35) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
9	5121378	1.0	No	CRADLE COVER BRKT positioning GT
12	5124107	1.0	No	CARRIAGE-CABLE-SUP-A Positioning GT
13	5124109	1.0	No	CARRIAGE-CABLE-SUP-B Positioning GT
15	5123640	3.0	No	OILLESS BUSH Bearing Special Parts U0934AH
51	N9514ZX	6.0	No	SCREW HEXAGON SOCKET 5 MM 14 MM
56	N9510ES	4.0	No	SCREW COUNTERSUNK HEA 5 MM 10 MM SCREW s[5X10 FE ZN
58	5127773	1.0	No	ELEVATION ASSY Positioning GT {For further detail, see Section 1.1.1, beginning on page 50. }
59	5117455	1.0	No	LINK BAR ASSEMBLY Positioning GT
60	5127782	1.0	No	TOP FRAME ASSY {For further detail, see Section 1.1.2, beginning on page 60. }
61	5126796-2	1.0	No	FOOT SW ASSY Positioning GT {For further detail, see Section 1.1.3, beginning on page 86. }
63	5154867-2	1.0	No	RIM BASE COVER LH ASSY positioning GT
64	5155337-2	1.0	No	RIM BASE COVER RH ASSY positioning GT
76	5113566	1.0	No	BASE BRKT 2R positioning GT
77	5113567	1.0	No	BASE BRKT 2L positioning GT
78	5113573	1.0	No	BRKT BASE FRAME L positioning GT
79	5113574	1.0	No	BRKT BASE FRAME R positioning GT
80	5159250	1.0	No	SIDE PANEL LWR BRKT L COMMON positioning GT
81	5159249	1.0	No	SIDE PANEL LWR BRKT R COMMON positioning GT
82	5112108-3	1.0	No	IMS Driving Cover Positioning GT
83	5112217	1.0	No	IMS D COVER FRONT Positioning GT
84	5112218	1.0	No	IMS D COVER MID Positioning GT
85	5118042	1.0	No	IMS D CV M BRKT 1 positioning GT
88	5126184	1.0	No	IMS-D-COVER-RR positioning GT
90	5136785	2.0	No	IMS COVER BRKT 18 positioning GT
91	5113591	1.0	No	IMS COVER BRKT 13L positioning GT
92	5113592	1.0	No	IMS COVER BRKT 13R positioning GT

Table 2-2 Table Assembly (5127787 BOM, Rev 35) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
95	5113596	2.0	No	IMS COVER BRKT 3 positioning GT
98	5136786	1.0	No	IMS COVER BRKT 19L positioning GT
99	5136787	1.0	No	IMS COVER BRKT 19R positioning GT
106	5113610	1.0	No	TOP FRONT BRKT 2L positioning GT
107	5113611	1.0	No	TOP FRONT BRKT 2R positioning GT
108	5113612	1.0	No	TOP FRONT BRKT 3 positioning GT
111	5119714	1.0	No	CABLE BEAR GUIDE positioning GT
114	N9408ES	4.0	No	SCREW COUNTERSUNK HEA 4 MM 8 MM
117	N9512LK	12.0	No	SCREW BIND UTIBATUKI 5X12 BsNi
118	N9514ZX	4.0	No	SCREW HEXAGON SOCKET 5 MM 14 MM
124	5119692	1.0	No	L-SW-BRKT-1 positioning GT
125	5119693	1.0	No	L-SW-BRKT-2 positioning GT
126	5119694	1.0	No	L-SW-BRKT-3 positioning GT
127	5119695	2.0	No	L-SW-BRKT-4 positioning GT
128	5119696	1.0	No	L-SW-BRKT-5 positioning GT
129	5119697	3.0	No	NUT-PLATE LS positioning GT
130	5119711	3.0	No	SHAFT TSP positioning GT
136	5127957	12.0	No	CARRIAGE SHIM02 Positioning GT
138	5123743	1.0	No	IMS-D-COVER-FF positioning GT
139	5127759	1.0	No	TRAY-F positioning GT
141	5127761	1.0	No	AMP-COVER positioning GT
147	5129565	3.0	No	SHAFT CUSHION Potisioning GT
167	N9316ZX	2.0	No	SCREW CAP SCREW 3X16 Fe * CAP SCREW 3X16 FE
171	2262896	10.0	No	U0465AA FLANGE SOCKET
175	N9430JS	2.0	No	SCREW PAN HEAD 4 MM 30 MM SCREW osX 4X30 FE ZN
177	N9801BS	4.0	No	NUT HEXAGON 8 MM FE ZN
178	N9850YX	6.0	No	SCREW HEXAGON SOCKET 8 MM 50 MM CAP SCREW 8X50 FE
179	U0241AA	14.0	No	NUT HEXAGON 4 MM NUT

Table 2-2 Table Assembly (5127787 BOM, Rev 35) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
183	2299437	12.0	No	U0472AA_FLANGE_SOCKET
186	5121904	3.0	No	CAM FOLLOWER BEARING U0443AH
210	5139862	1.0	No	IMS-COVER-BRKT-24L positioning GT
211	5139863	1.0	No	IMS-COVER-BRKT-24R positioning GT
212	5151800	1.0	No	TOP REAR BRKT COMMON 6L Positioning GT
213	5151801	1.0	No	TOP REAR BRKT COMMON 6R Positioning GT
214	5151802	1.0	No	TOP REAR BRKT COMMON 7L Positioning GT
215	5151803	1.0	No	TOP REAR BRKT COMMON 7R Positioning GT
300	5160028	0.0025	No	LUBRI PLATE,U0304BD
301	U0303BD	6.25	No	MORICORT M-77 PASTE
303	U0409BD	0.313	No	ADHESIVE Yhjb #601
305	U0046BA	180.0	No	TAPE
306	5160029	0.00625	No	KURE CRC,U0305BD
500	5128983	1.0	No	ETL LABEL Positioning GT
502	5129157	1.0	No	PLACE LABEL Potisioning GT
503	5129016	2.0	No	SHIPPING LABEL Positioning GT
506	U0036BN	12.0	No	NAME PLATE
507	5139485	1.0	No	WEEE LABEL WITH BIN ONLY
801	5142259	1.0	No	IMS COVER BRKT 31L positioning GT
802	5142261	1.0	No	IMS COVER BRKT 31R positioning GT
803	5142269	1.0	No	IMS COVER BRKT 33 positioning GT
804	5175971	2.0	No	PLATE NUT Positioning GT
805	5175384	1.0	No	IMS COVER BRKT 33R Positioning GT
806	5175137	1.0	No	IMS COVER BRKT 33L Positioning GT
807	5120580-2	1.0	No	TS_IMS_REAR Positioning GT
808	5189420	1.0	No	IMS DUMMY PLUG LABEL positioning GT

Table 2-2 Table Assembly (5127787 BOM, Rev 35) (Continued)

1.1.1 Elevation Assembly (5127773)

FRU = NO

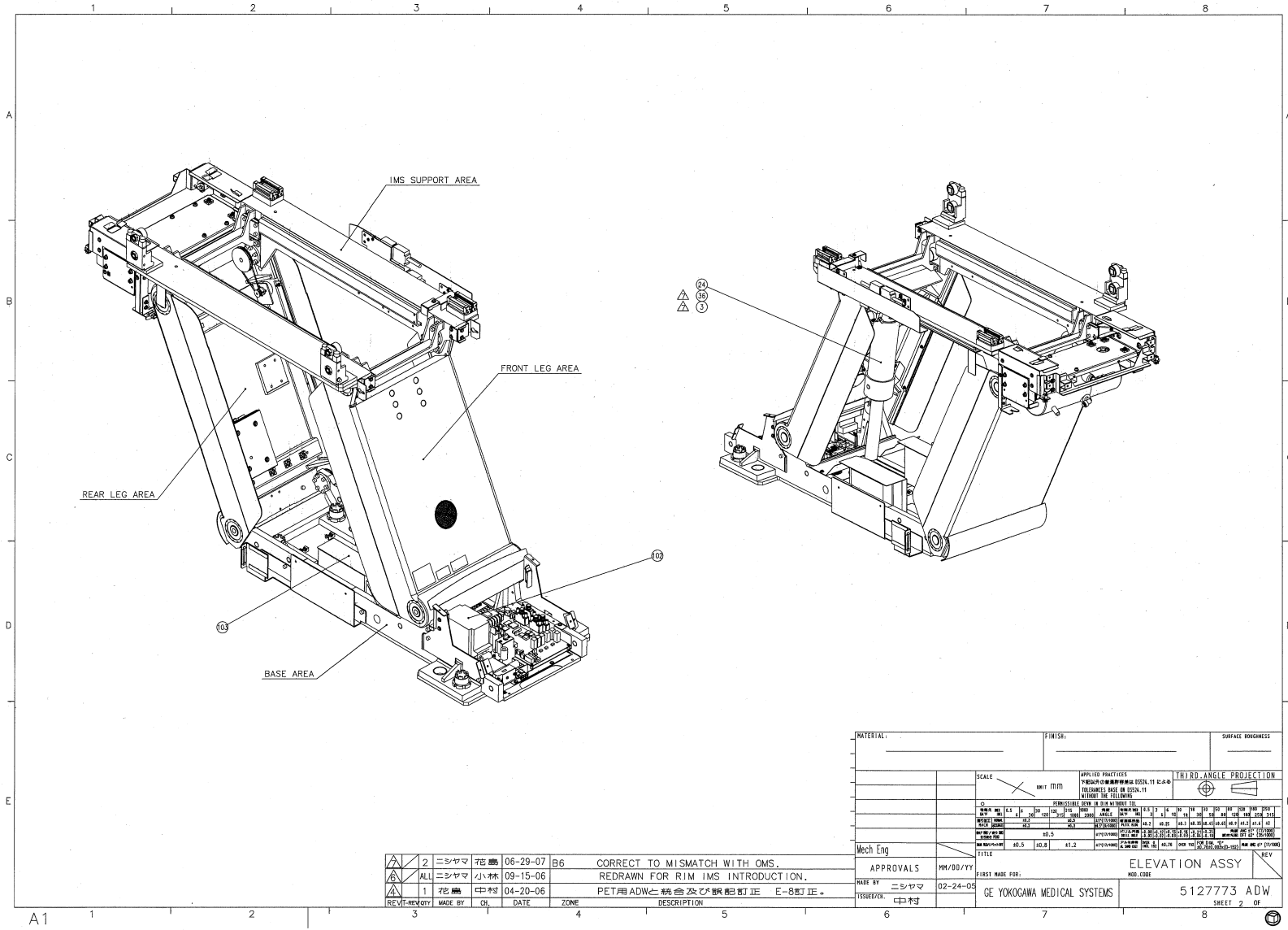


Figure 2-12 Elevation Assembly (5127773ADW, sheet 2, rev 8)

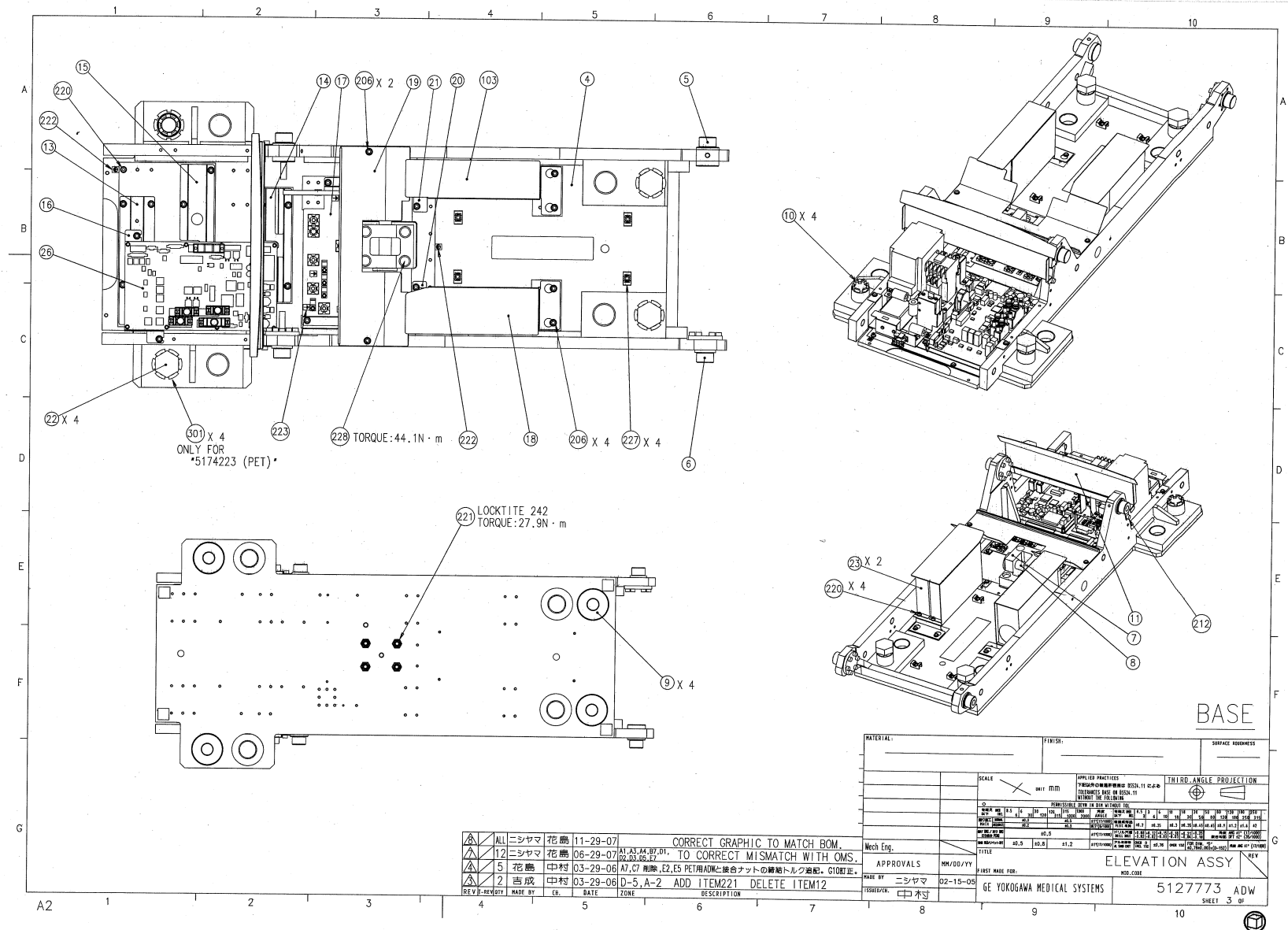


Figure 2-13 Elevation Assembly (5127773ADW, sheet 3, rev 8)

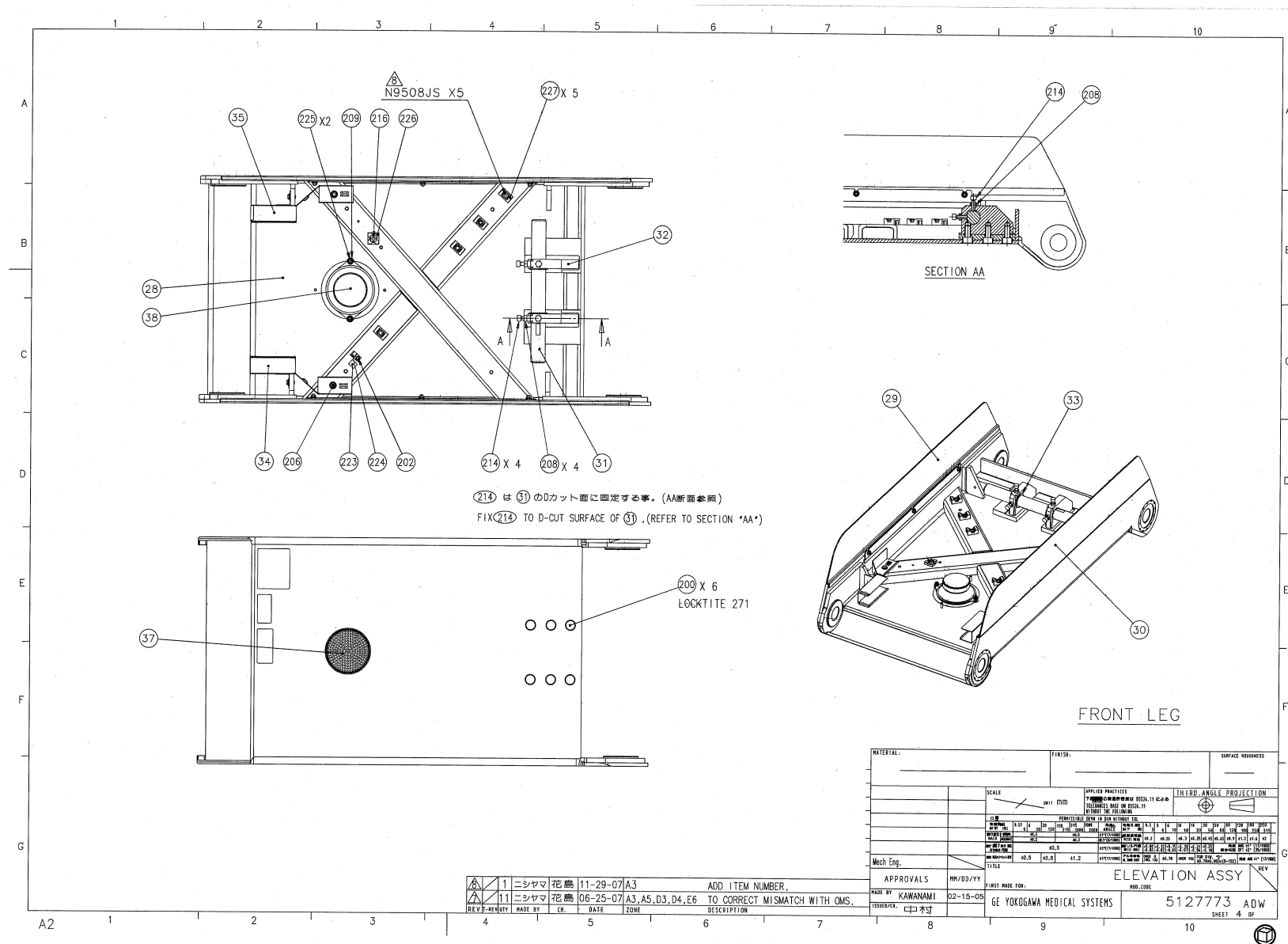


Figure 2-14 Elevation Assembly (5127773ADW, sheet 4, rev 8)



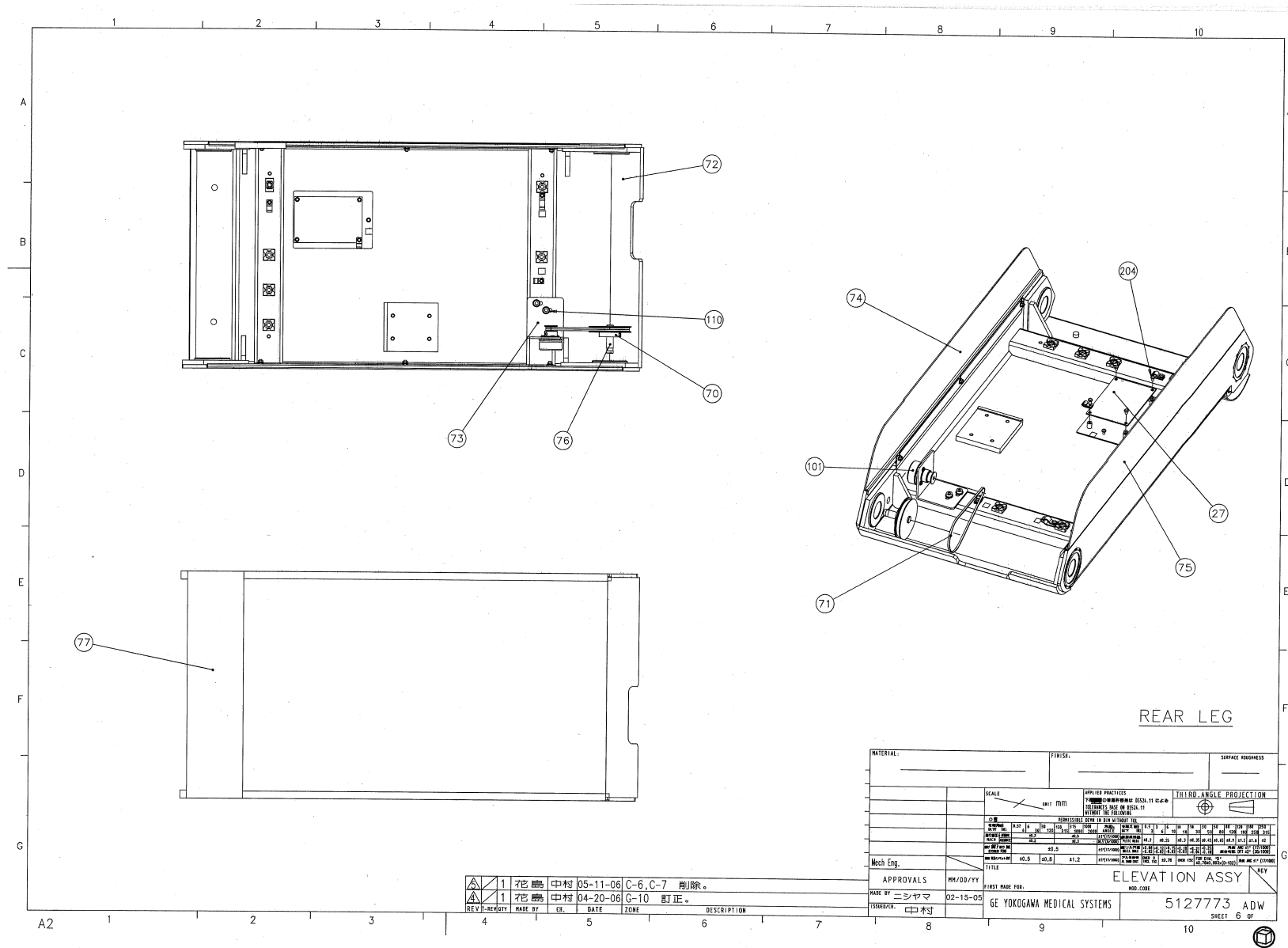


Figure 2-16 Elevation Assembly (5127773ADW, sheet 6, rev 8)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
2	U0097AG	8.0	Yes	STANDOFF SPACER MALE_FEMALE 4 MM 10
23	5132669-3	2.0	Yes	STEP MOTOR DRIVER Positioning GT
26	5142177-3	1.0	Yes	GTCB2 BOARD ASSY Positioning GT
27	5116570	1.0	Yes	GTTS BOARD ASSY Positioning GT
38	46-297687P1	1.0	Yes	8 OHM, 4W, 3-1/2~ SPEAKER. 2.9 OZ CERAMIC MAGNET, 90 DB. (0.187~ Q.C. TERMINALS) (DWG REV 3/92, JFK - OK).
71	5125804	1.0	Yes	BELT 350-2GT-3 Positioning GT
101	5127774	1.0	Yes	ELEVATION POTENTIOMETER ASSY
102	5127381-2	1.0	Yes	INVERTER ASSY Positioning GT
103	5127384	1.0	Yes	POWER SUPPLY ASSY Positioning GT
200	N9020YX	1.0	Yes	SCREW HEXAGON SOCKET 10 MM 20 MM CAP SCREW 10X20 FE
202	N9408JS	1.0	Yes	SCREW NABE 4X8 Fe Zn
203	N9401WS	1.0	Yes	WASHER PLAIN _ LARGE 4.3 MM 9 MM WASHER 4.3x9x0.8 FE ZN
204	N9408HR	8.0	Yes	SCREW TP
205	N9508SJ	3.0	Yes	SCREW SET SCREW 5 MM 8 MM SET SCREW M5X0.8X8htXvsg
206	N9510HR	38.0	Yes	SCREW TP
207	N9510ZX	5.0	Yes	SCREW HEXAGON SOCKET 5 MM 10 MM CAP SCREW 5X10 FE
208	N9601BS	2.0	Yes	NUT HEXAGON 6 MM NUT M6 2 FE ZN
209	N9601WS	16.0	Yes	WASHER PLAIN - LARGE 6.4 MM 12.5 MM WASHER 6.4X12.5X1.6FEZN
212	N9612YX	23.0	Yes	SCREW HEXAGON SOCKET 6 MM 12 MM CAP SCREW 6X12 FE
214	N9620YX	14.0	Yes	SCREW HEXAGON SOCKET 6 MM 20 MM CAP SCREW 6X20 FE
216	N9306ES	1.0	Yes	SCREW COUNTERSUNK HEA 3 MM 6 MM
221	46-328426P3	4.0	Yes	NUT HEXAGON_THIN 10 MM
224	U0030AC	1.0	Yes	CLAMP OR CLIP CLAMP
226	U0212AC	1.0	Yes	CLAMP TMS_20
227	U0208AC	9.0	Yes	CLAMP TC-142 TAIRATSUPU KOTEI
229	N9801WS	2.0	Yes	WASHER PLAIN - LARGE 8.4 MM 17 MM

Table 2-3 ELAVATION ASSY (5127773 BOM, Rev 25)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
3	5123754	1.0	No	FINE THREAD NUT U0276AA
4	5112074	1.0	No	TABLE BASE Positioning GT
5	5112928	2.0	No	BASE SHAFT R positioning GT
6	5112929	2.0	No	BASE SHAFT L positioning GT
7	5116290	1.0	No	CYLINDER BASE MOUNT 22 Positioning GT
8	5116291	1.0	No	CYLINDER BASE SHAFT 22
9	5112943	4.0	No	LEVERING PAD positioning GT
10	5112944	4.0	No	LEVERING ADJUSTER positioning GT
11	5112945	1.0	No	LIQUID RECIEVER positioning GT
13	5117581	1.0	No	TRAY-A Positioning GT
14	5117582	1.0	No	TRAY-B Positioning GT
15	5117584	1.0	No	TRAY-C Positioning GT
16	5117585	1.0	No	GTCB-BASE Positining GT
17	5121246	1.0	No	CABLE BASE Positioning GT
18	5121247	1.0	No	DRIVER BASE COVER Positioning GT
19	5121248	1.0	No	LIQUID RECEIVER COVER Positioning GT
20	5121249	1.0	No	TRAY D Positioning GT
21	5122886	1.0	No	TRAY E Positioning GT
22	5123854	4.0	No	SHIPPING BOLT M27 P2 Positioning GT
24	5112219	1.0	No	CYLINDER ASSY Positioning GT
28	5112107	1.0	No	FRONT LEG Positioning GT
29	5112949	1.0	No	FRONT LEG PLATE L positioning GT
30	5112950	1.0	No	FRONT LEG PLATE R positioning GT
31	5117114	1.0	No	CYLINDER ROD D25 Positioning GT
32	5117124	2.0	No	CYLINDER BASE HIGH D25 Positioning GT
33	5117129	2.0	No	WASHER 253215 Positioning GT
34	5125241	1.0	No	CABLE DUCT FLEG R Positioning GT
35	5125242	1.0	No	CABLE DUCT FLEG L Positioning GT

Table 2-3 ELAVATION ASSY (5127773 BOM, Rev 25) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
36	5121918	1.0	No	ROD END SCREW TYPE U0910AP
37	2261674	1.0	No	SPEAKER GRILL, (REF 46-297684P1)
39	5123254	4.0	No	RUBBER DUMPER CBU Positioning GT
40	5121924	4.0	No	FINETHREAD NUT U0275AA
41	5121913	4.0	No	ROLLER FOLLOWER BEARING U0444AH
46	5112109	1.0	No	LEFT CHANNEL SUPPORT Positioning GT
47	5112110	1.0	No	RIGHT CHANNEL SUPPORT Positioning GT
48	5113098	2.0	No	SUP REIN positioning GT
49	5113084	2.0	No	PUMP BASE BRKT R positioning GT
50	5113085	1.0	No	PUMP SPRT BRKT UPPER positioning GT
51	5113086	1.0	No	CHANNEL SUPPORT SHAFT positioning GT
52	5113087	1.0	No	CHANNEL SPRT SHAFT POT positioning GT
53	5113092	4.0	No	BEARING SHAFT positioning GT
54	5113093	2.0	No	LOWER SUP IMS positioning GT
55	5113095	2.0	No	UPPER SUP IMS positioning GT
56	5112928	2.0	No	BASE SHAFT R positioning GT
57	5113580-2	1.0	No	PUMP-SPRT-BRKT-L Positioning GT
58	5113581-2	1.0	No	PUMP-SPRT-BRKT-R Positioning GT
59	5113582	1.0	No	CABLE BEAR SUP BASE positioning GT
61	5116961	1.0	No	CABLE GUIDE BRKT1 Positioning GT
62	5116962	1.0	No	CABLE GUIDE BRKT2 Positioning GT
63	5117028	1.0	No	CONNECTOR BRKT Positioning GT
64	5117316	2.0	No	IMS SHIM01 Positioning GT
65	5117317	2.0	No	IMS SHIM02 Positioning GT
66	5123881	1.0	No	IMS SHIPPING LOCK positioning GT
67	5122745	4.0	No	CABLE-MOUNT-13 positioning GT
68	5125877	2.0	No	PUMP NUT PLATE Positioning GT
69	5180482	1.0	No	FINGER GUARD Positioning GT

Table 2-3 ELAVATION ASSY (5127773 BOM, Rev 25) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
70	5125802	1.0	No	POT PULLY BIG Positioning GT
72	5112202	1.0	No	REAR LEG Positioning GT
73	5113104	1.0	No	HEIGHT POT SPRT positioning GT
74	5113107	1.0	No	REAR LEG PLATE L positioning GT
75	5113108	1.0	No	REAR LEG PLATE R positioning GT
76	5113109	1.0	No	SHAFT POT 1 positioning GT
77	5115669	1.0	No	REAR LEG COVER positioning GT
99	5127955	2.0	No	IMS SUPPORT SHIM01 Positioning GT
100	5127956	2.0	No	IMS SUPPORT SHIM02 Positioning GT
110	2263506	2.0	No	U0466AA SET SCREW
111	5130886	2.0	No	IMS-SUP-L-SHIM005S positioning GT
112	5130887	8.0	No	IMS-SUP-L-SHIM01S positioning GT
113	5130881	4.0	No	IMS-SUPPORT-SHIM05 positioning GT
114	5131440	2.0	No	IMS-SUPPORT-SHIM01S positioning GT
201	N9030YX	4.0	No	SCREW HEXAGON SOCKET 10 MM 30 MM CAP 10X30
217	N9820YX	6.0	No	SCREW HEXAGON SOCKET 8 MM 20 MM CAP SCREW 8X20 FE
218	Z9817DL	1.4	No	SUPER HIGHLAND
219	2262896	8.0	No	U0465AA FLANGE SOCKET
220	2299437	5.0	No	U0472AA_FLANGE_SOCKET
222	U0022BN	6.0	No	NAME PLATE NAME PLATE-EARTH MARK
223	U0036BN	2.0	No	NAME PLATE
225	U0241AA	2.0	No	NUT HEXAGON 4 MM NUT
228	N9025YX	4.0	No	SCREW CAP SCREW 10X25FE * CAP SCREW 10X25FE

Table 2-3 ELAVATION ASSY (5127773 BOM, Rev 25) (Continued)

1.1.2 TOP FRAME ASSY (5127782)

FRU = NO

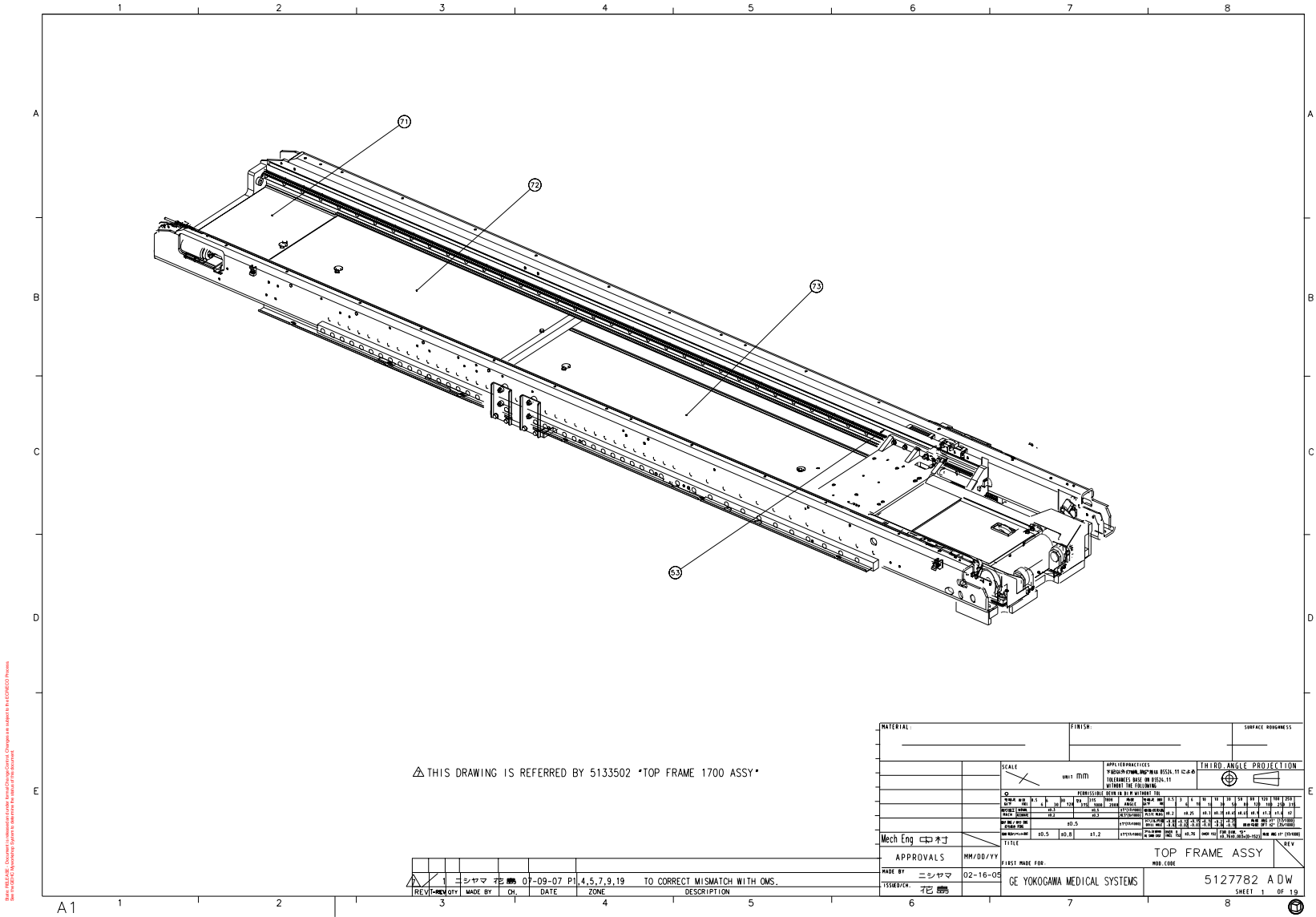


Figure 2-17 Top Frame Assembly (5127782ADW, sheet 1, rev 7)

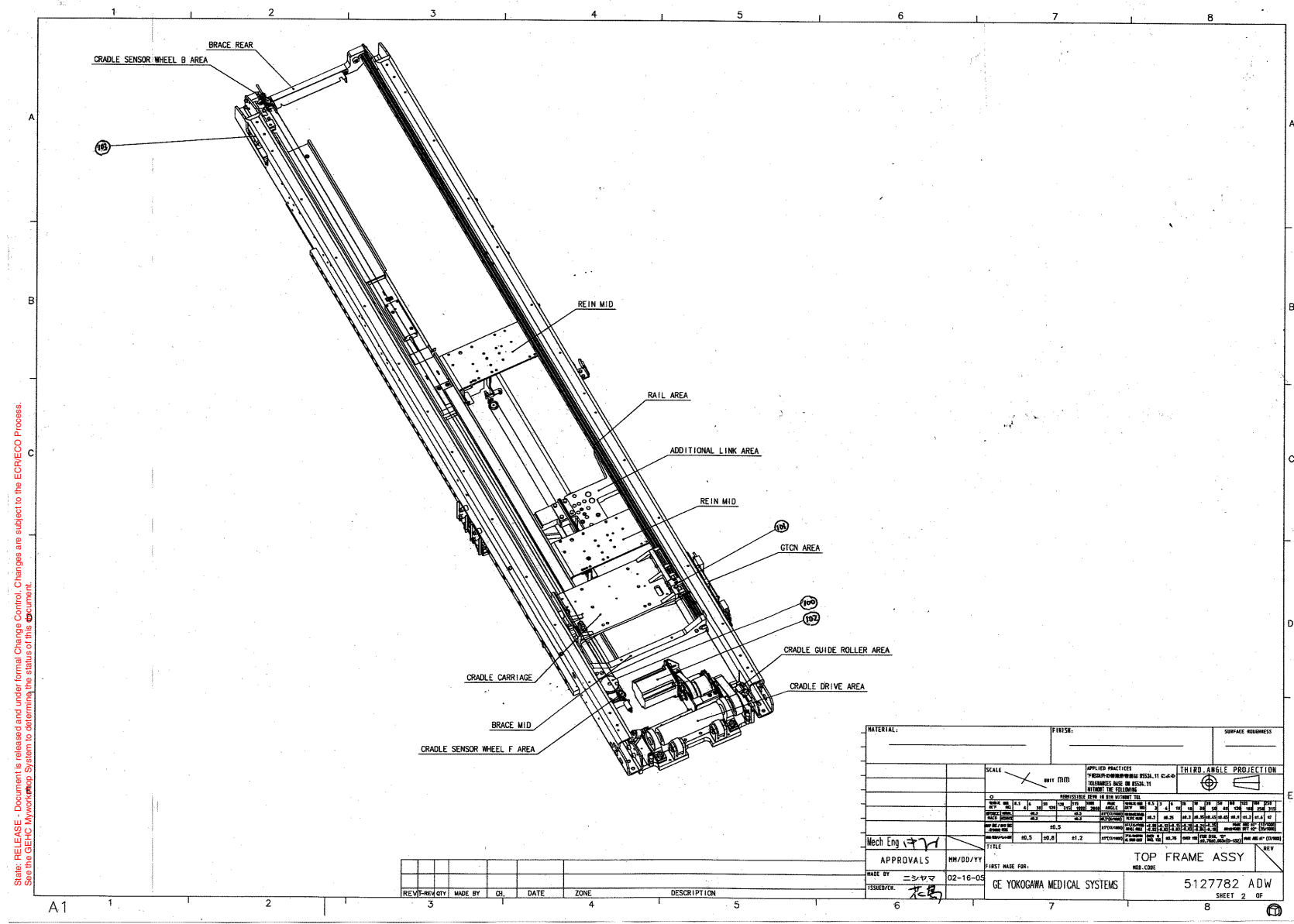


Figure 2-18 Top Frame Assembly (512778ADW2, sheet 2, rev 7)



Section 1.0 — GT 2000 Table (5121647)

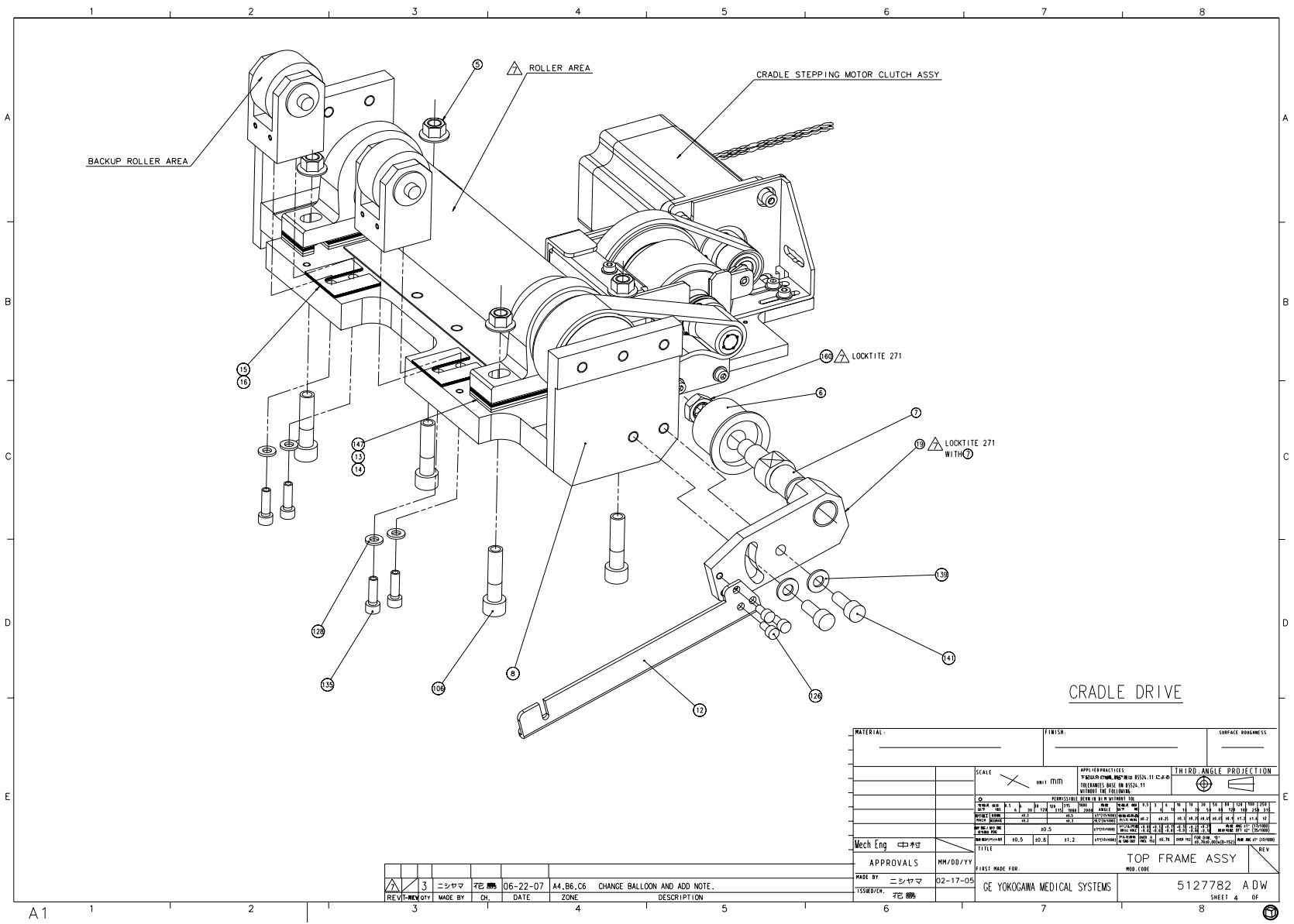


Figure 2-20 Top Frame Assembly (5127782ADW, sheet 4, rev 7)



Section 1.0 — GT 2000 Table (5121647)



2 - GT Table
GT 2000 Table (5121647)

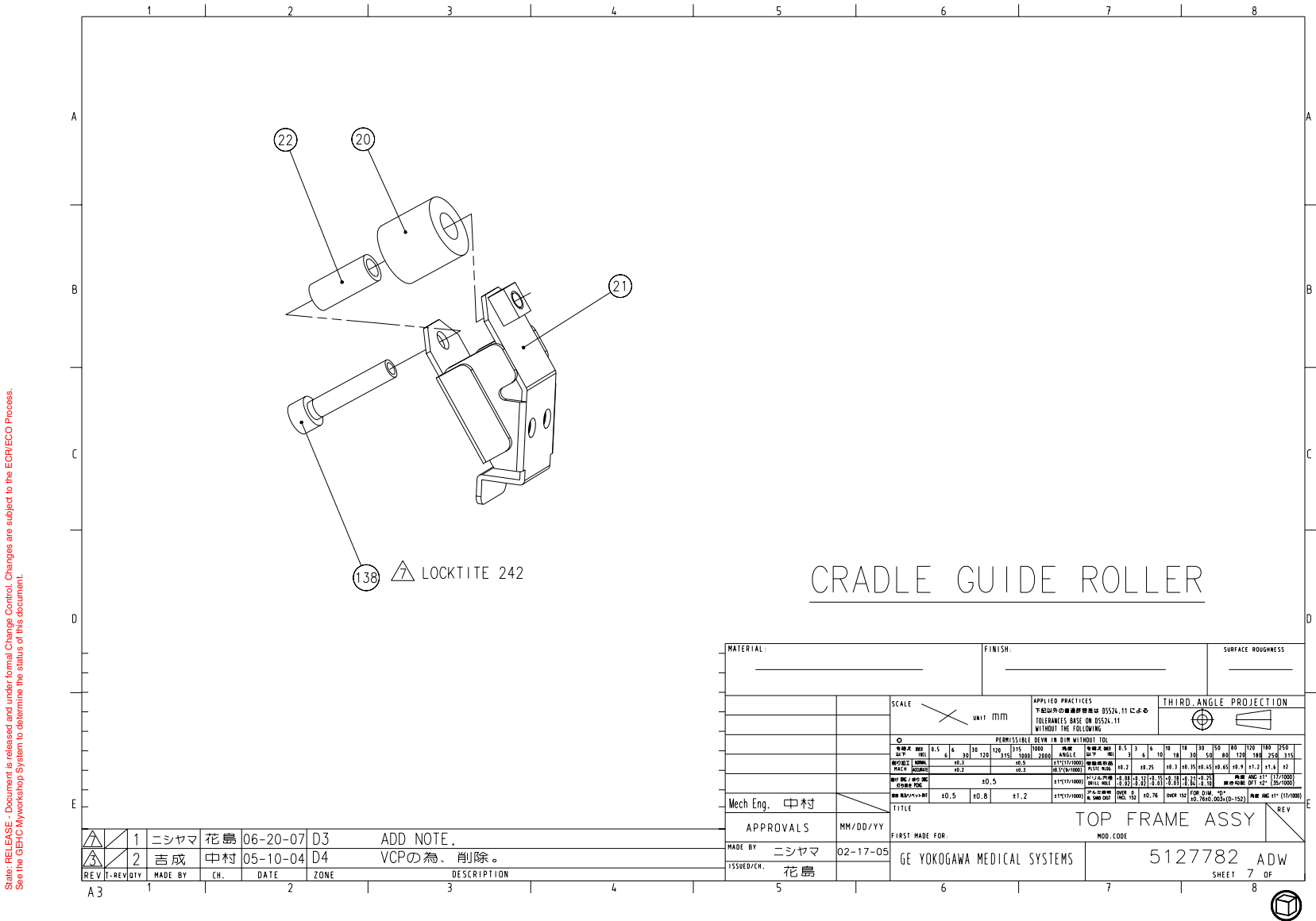


Figure 2-23 Top Frame Assembly (5127782ADW, sheet 7, rev 7)

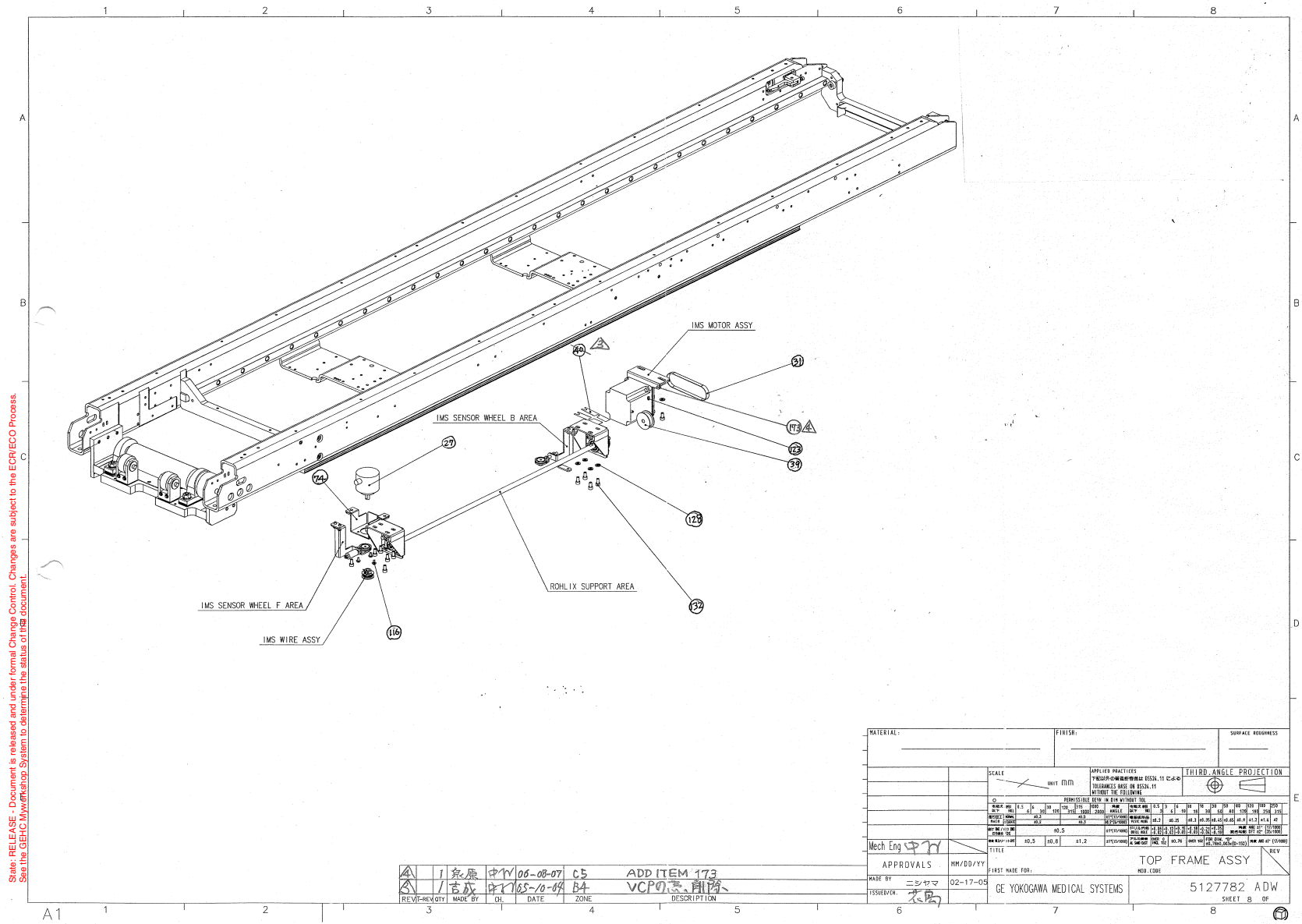


Figure 2-24 Top Frame Assembly (5127782ADW, sheet 8, rev 7)



Section 1.0 — GT 2000 Table (5121647)



2 - GT Table
GT 2000 Table (5121647)

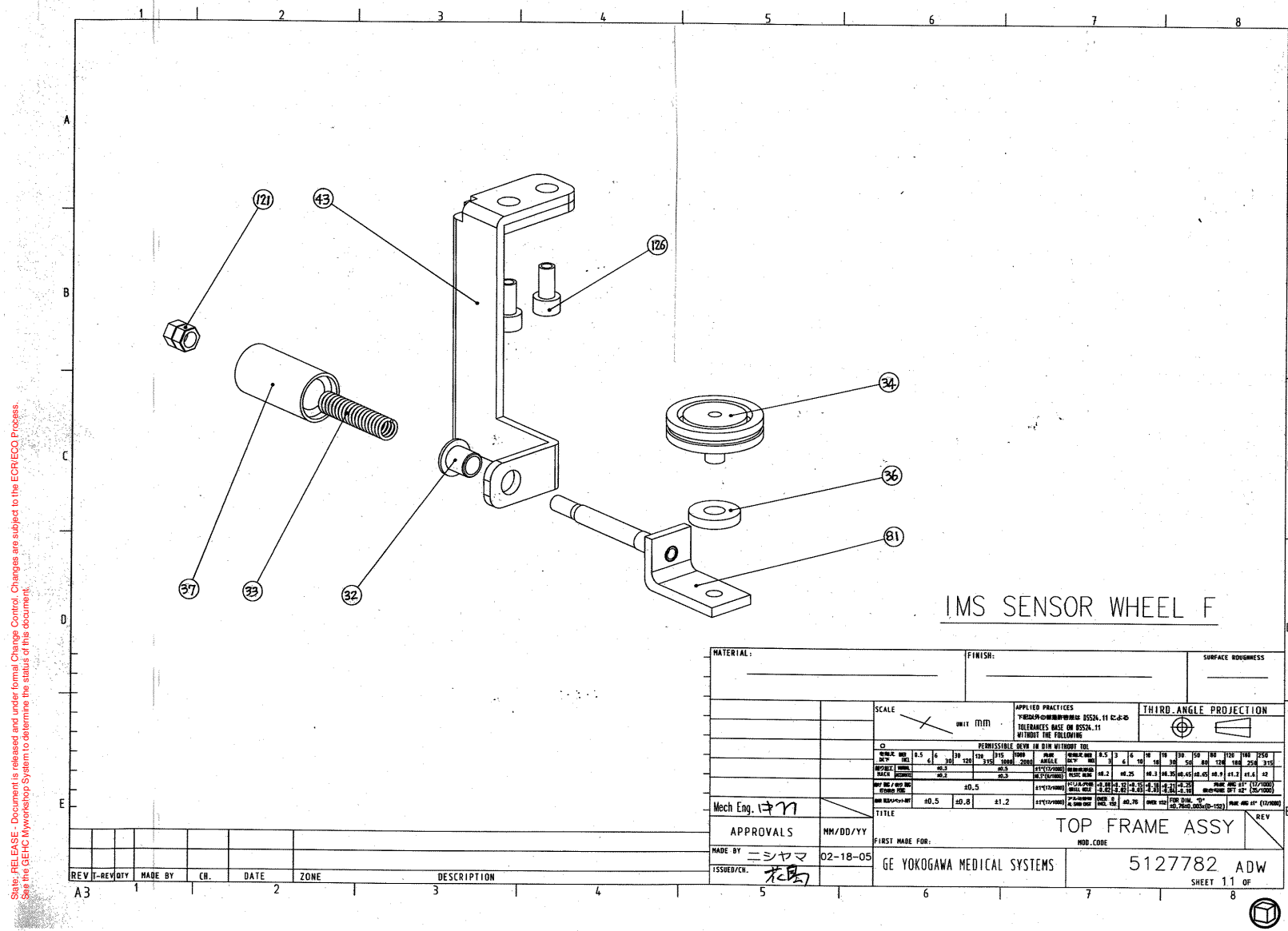


Figure 2-27 Top Frame Assembly (5127782ADW, sheet 11, rev 7)



2 - GT Table
GT 2000 Table (5121647)

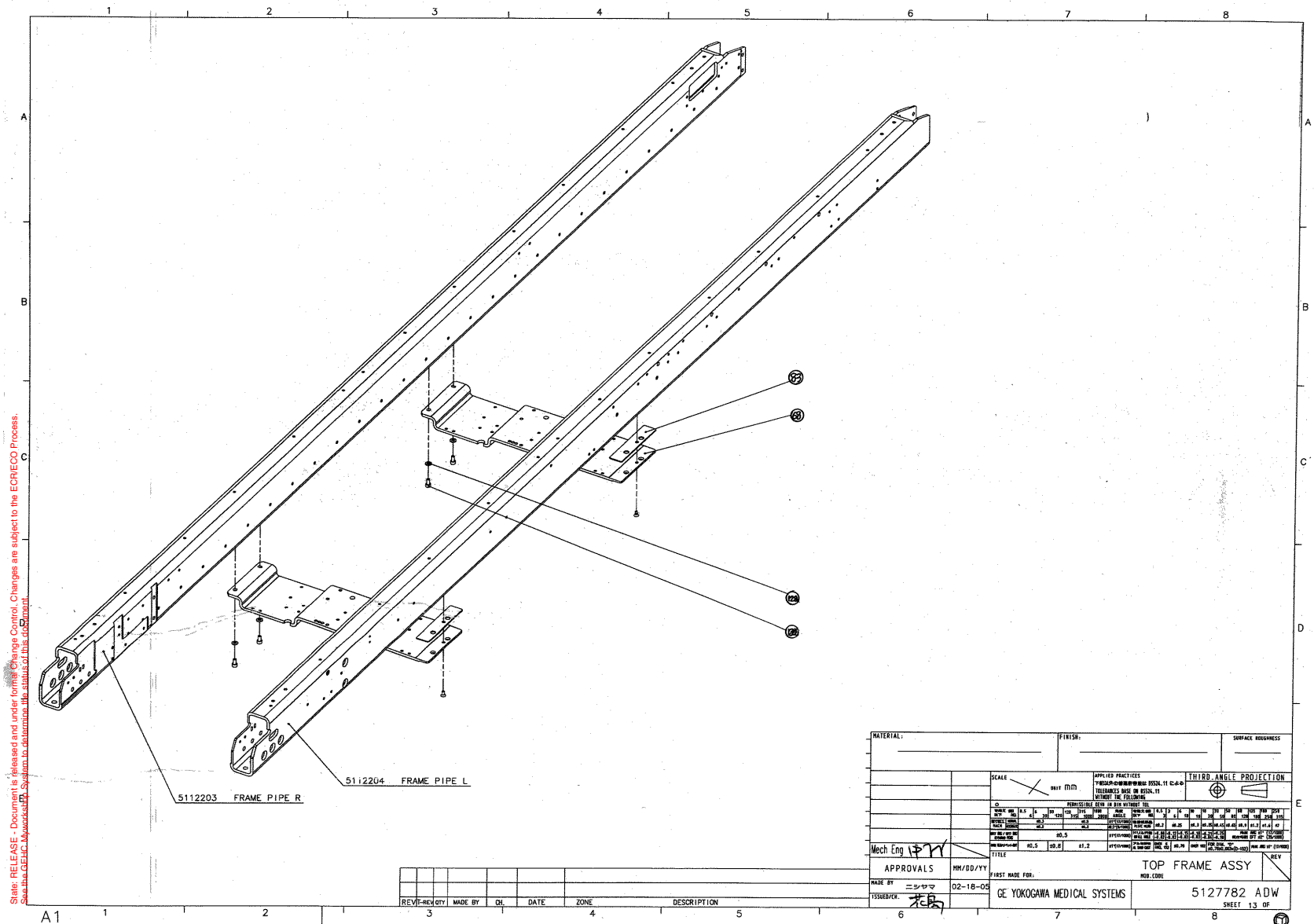


Figure 2-29 Top Frame Assembly (5127782ADW, sheet 13, rev 7)

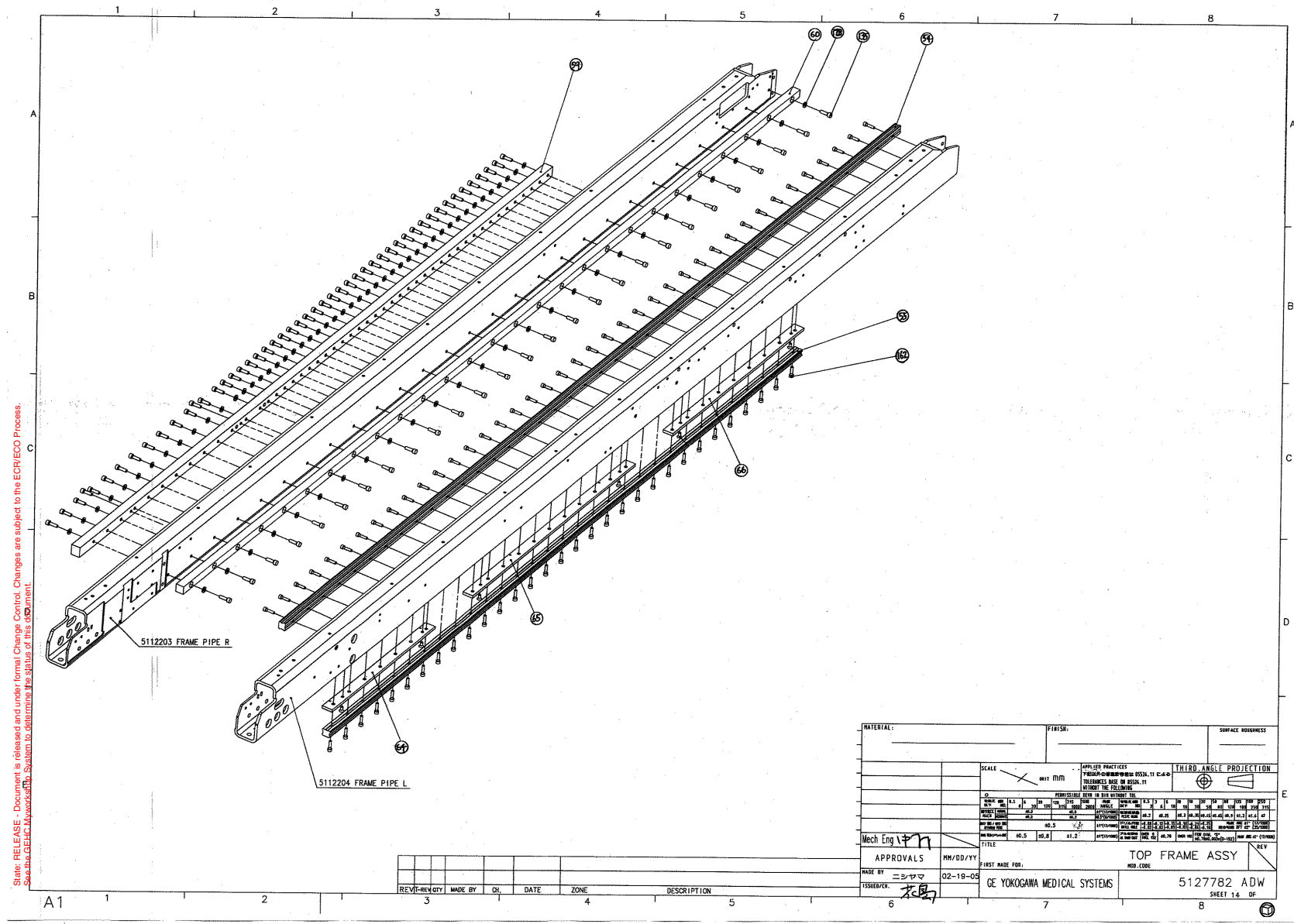


Figure 2-30 Top Frame Assembly (5127782ADW, sheet 14, rev 7)



Section 1.0 — GT 2000 Table (5121647)

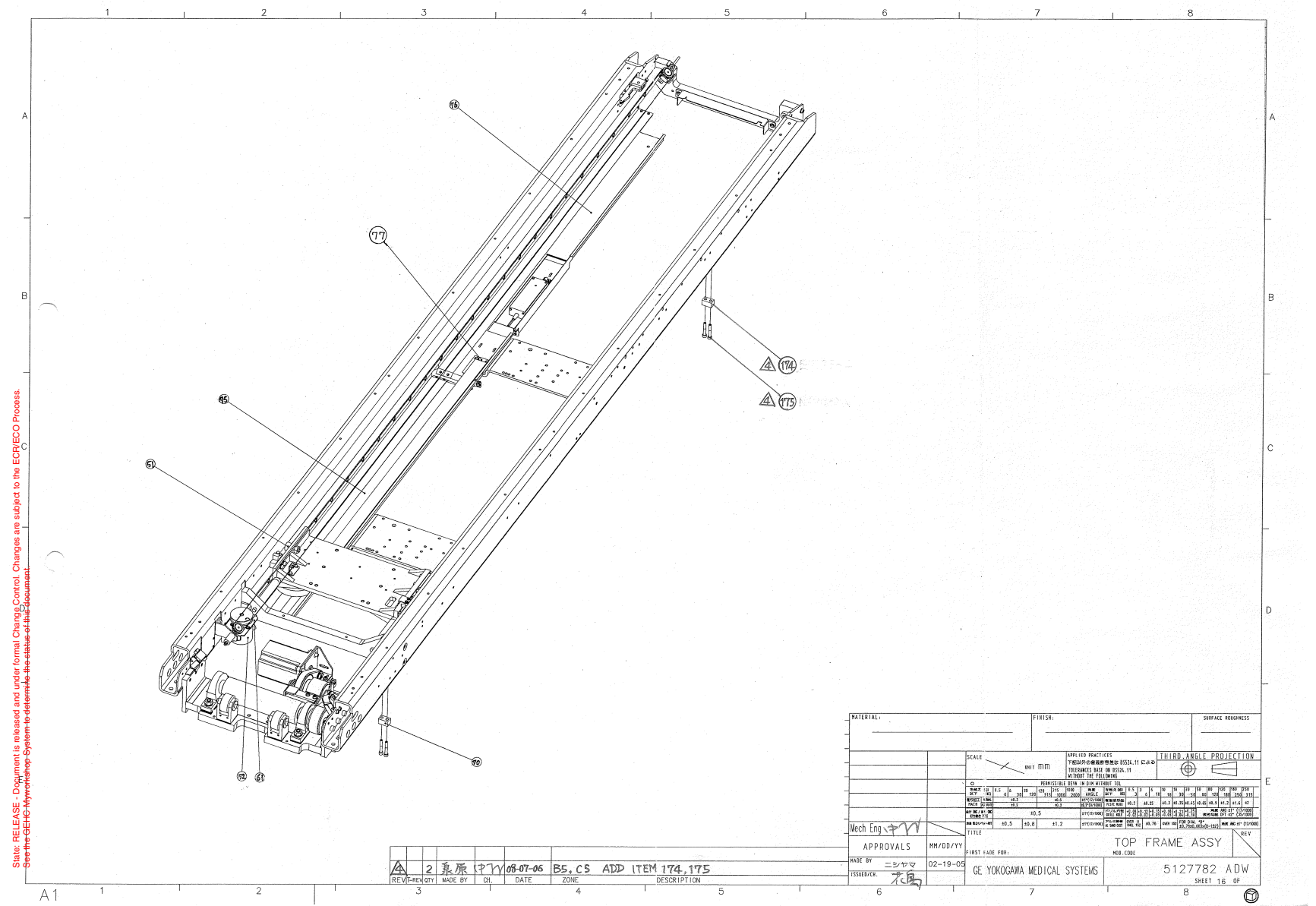


Figure 2-32 Top Frame Assembly (5127782ADW, sheet 16, rev 7)

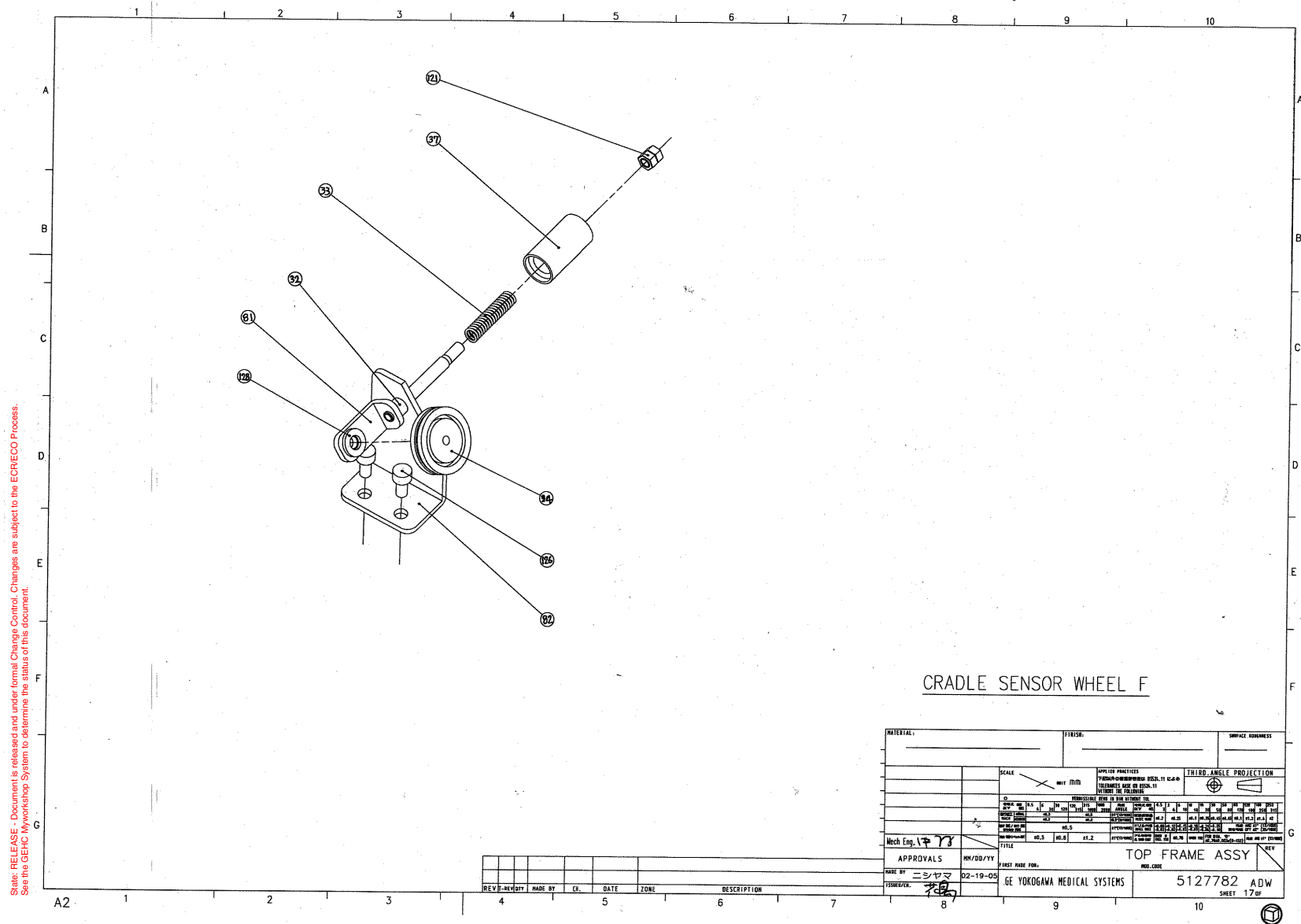


Figure 2-33 Top Frame Assembly (5127782ADW, sheet 17, rev 7)

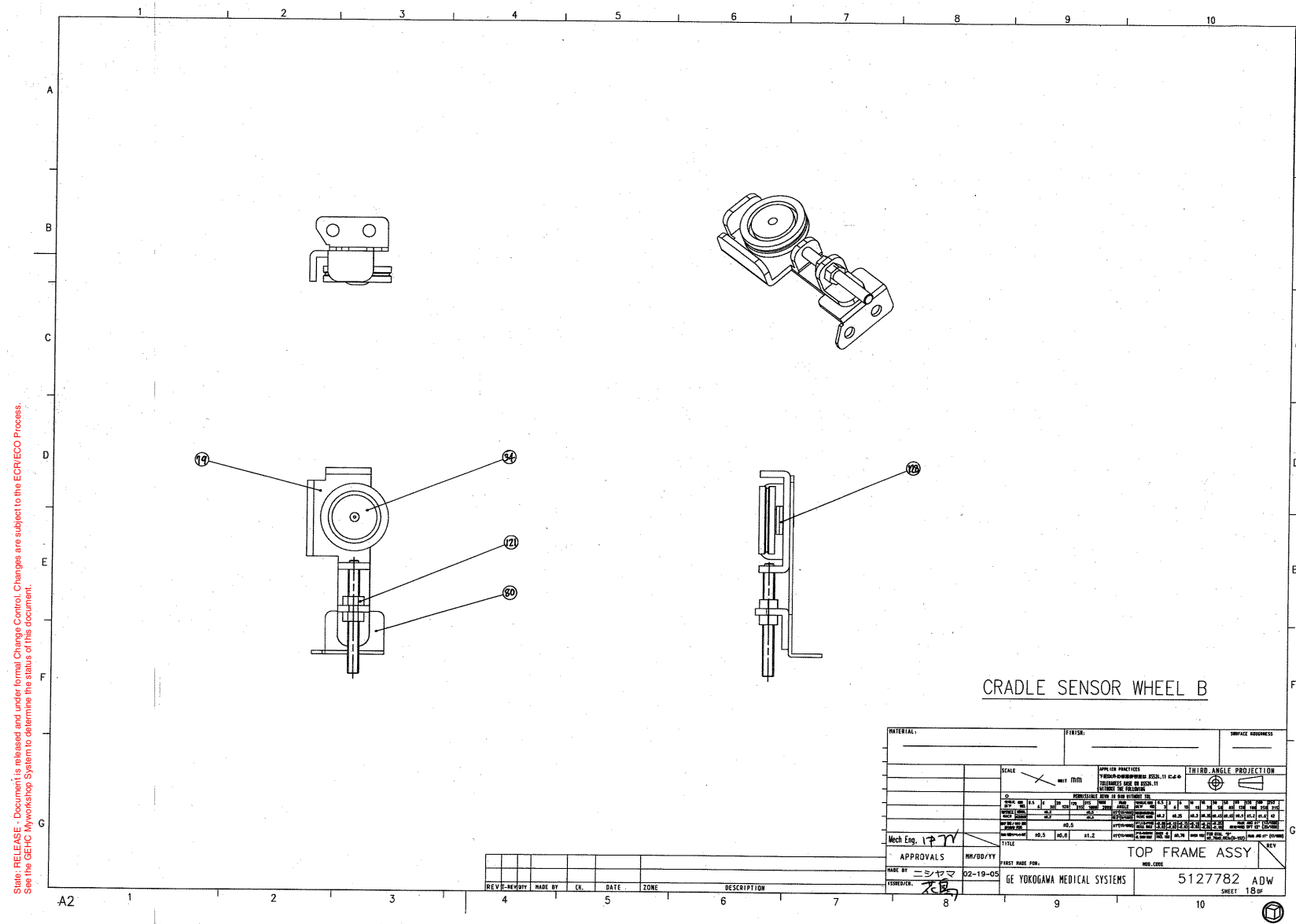


Figure 2-34 Top Frame Assembly (5127782ADW, sheet 18, rev 7)

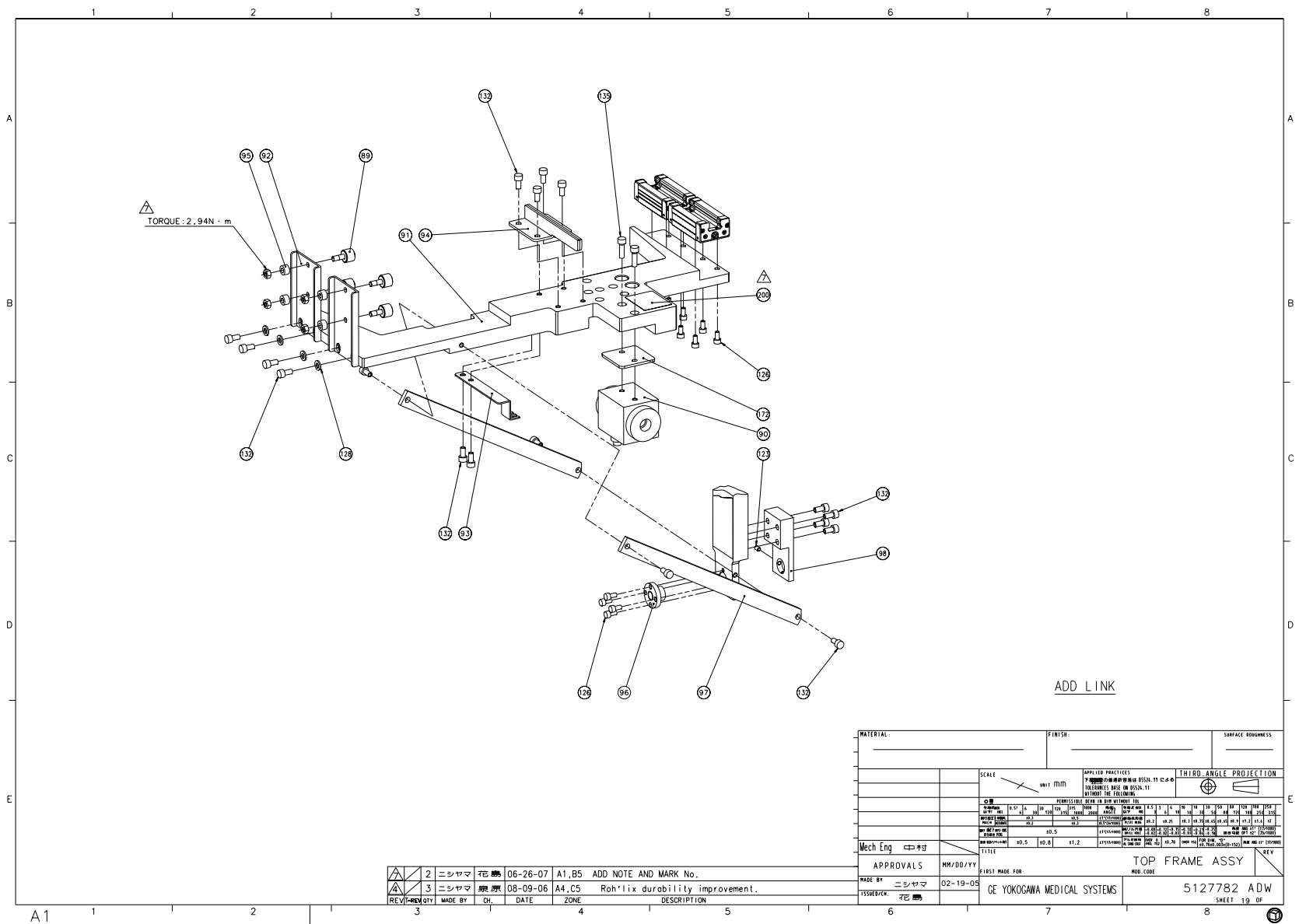


Figure 2-35 Top Frame Assembly (5127782ADW, sheet 19, rev 7)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2201853	2.0	Yes	U0226AH BEARING BALL,DOUBLE RO 20MM
23	P9145CZ	4.0	Yes	WASHER ROLLER
27	5122208-2	1.0	Yes	ABSOLUTE ENCODER FOR IMS GT
31	5123252	1.0	Yes	BELT 318P3M10
32	U0908AH	2.0	Yes	BUSHING 80F-0608
34	U0009AH	4.0	Yes	CAM FOLLOWER WITH SHAFT 6 MM 30 MM DV-30-B0.5
36	P9145DY	4.0	Yes	SPACER
37	P9145DL	2.0	Yes	CASE SPRING
48	U0443BA	2.0	Yes	TUBE .vk X F4X0.25 hY
50	U0443BA	1.0	Yes	TUBE .vk X F4X0.25 hY
52	5122207-2	1.0	Yes	ABSOLUTE ENCODER FOR CRADLE GT
53	5122210	1.0	Yes	MAGNET SCALE GT
56	U0125AE	4.0	Yes	CUSHION CK-17
89	U0402AH	4.0	Yes	CAM FOLLOWER WITH SHAFT 6 MM 16 MM CF 6 VBUU
99	5127495	1.0	Yes	IMS WIRE ASSY Positioing GT
100	5127502	1.0	Yes	CRADLE WIRE ASSY Positioing GT
101	5127674	1.0	Yes	LINEAR SCALE ASSY
102	5127493	1.0	Yes	CRADLE STEEPING MOTOR CLUTCH 2000 ASSY Positioning GT {For further detail, see Section 1.1.2.1, beginning on page 84. }
103	5127778	1.0	Yes	CRADLE LOCK ASSY
104	5127247-2	1.0	Yes	IMS Motor Assy Positioning GT
105	N9020YX	5.0	Yes	SCREW HEXAGON SOCKET 10 MM 20 MM CAP SCREW 10X20 FE
106	N9040YX	4.0	Yes	SCREW HEXAGON SOCKET 10 MM 40 MM CAP SCREW M10X40 FE
111	N9306ZX	4.0	Yes	SCREW HEXAGON SOCKET 3 MM 6 MM CAP SCREW 3X6 FE
116	N9408HR	2.0	Yes	SCREW TP
118	N9408ZX	2.0	Yes	SCREW HEXAGON SOCKET 4 MM 8 MM CAP SCREW 4X8 FE
121	N9501BS	8.0	Yes	NUT HEXAGON 5 MM NUT M5 FE ZN

Table 2-4 TOP FRAME ASSY (5127782 BOM, Rev 24)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
122	N9501WS	2.0	Yes	WASHER PLAIN - LARGE 5.3 MM 10 MM WASHER M5
123	N9508SJ	2.0	Yes	SCREW SET SCREW 5 MM 8 MM SET SCREW M5X0.8X8htXvsg
125	N9510LK	24.0	Yes	SCREW PAN HEAD 5 MM 10 MM
126	N9510ZX	8.0	Yes	SCREW HEXAGON SOCKET 5 MM 10 MM CAP SCREW 5X10 FE
128	N9601WS	11.0	Yes	WASHER PLAIN - LARGE 6.4 MM 12.5 MM WASHER 6.4X12.5X1.6FEZN
132	N9612YX	33.0	Yes	SCREW HEXAGON SOCKET 6 MM 12 MM CAP SCREW 6X12 FE
135	N9620YX	111.0	Yes	SCREW HEXAGON SOCKET 6 MM 20 MM CAP SCREW 6X20 FE
138	N9630YX	2.0	Yes	SCREW HEXAGON SOCKET 6 MM 30 MM CAP SCREW 6X30 FE
139	N9801WS	2.0	Yes	WASHER PLAIN - LARGE 8.4 MM 17 MM
142	N9825YX	4.0	Yes	SCREW CAP SCREW 8X25 Fe
161	N9612YX	6.0	Yes	SCREW HEXAGON SOCKET 6 MM 12 MM CAP SCREW 6X12 FE
162	N9520ZX	69.0	Yes	SCREW HEXAGON SOCKET 5 MM 20 MM CAP SCREW 5X20 FE
170	U0009AP	2.0	Yes	CRIMP STYLE TERMINAL B-4(CU2.5X1.5X3.1)
173	5184277	1.0	Yes	Rotary Dumper U0122AQ
181	N9412ZX	4.0	Yes	SCREW HEXAGON SOCKET 4 MM 12 MM SCREW
210	46-208718P3	1.0	Yes	CLAMP, CUSHIONED, ID .375 LOOP TYPE TA MFG CORP TA713D6
2	5121923	1.0	No	Parallel key U0163AP
3	5121921	2.0	No	BEARING UNIT U0447AH
4	5125774	4.0	No	FLANGE SOCKET U0478AA
5	5125773	4.0	No	FLANGE NUT U0276AA
6	5125253	1.0	No	IDLER NIPS3050 Positioning GT
7	5125254	1.0	No	IDLER PIN PID1522 Positioning GT
8	5112206	1.0	No	CRADLE-DRIVE-BASE Positioning GT
9	5113224	1.0	No	ROLLER positioning GT
10	5113219	2.0	No	BRKT ROLLER positioning GT
11	5113220	1.0	No	ROLLER PULLY positioning GT
12	5118619	1.0	No	TENSION BAR
13	5118587	10.0	No	SPACER 01 PILLOW

Table 2-4 TOP FRAME ASSY (5127782 BOM, Rev 24) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
14	5118588	6.0	No	SPACER 05 PILLOW
15	5118592	10.0	No	SPACER 01 BASE
16	5118591	6.0	No	SPACER 05 BASE
17	2199990	2.0	No	P9235JB PIN
18	2199991	4.0	No	P9235LA SPACER INNER
19	5117121	1.0	No	IDLER BASE positioning GT
20	5118181	2.0	No	CRADLE GUIDE ROLLER Positioning GT
21	5118183	2.0	No	CRADLE GUIDE ROLLER SPRT Positioning GT
22	5118182	2.0	No	SPACER BRASS Positioning GT
28	5126434	1.0	No	BEARING UNIT Positioning GT
29	2201842	1.0	No	U0428AH BEARING ATTACHMENT BEARING UNIT
30	5125163	1.0	No	IMS WHEEL SCREW BRKT Positioning GT
33	U0004AQ	2.0	No	SPRING COMPRESSION 47 MM SC-085-120-0470 4017
35	5115688	1.0	No	ROHLIX-SHAFT positioning GT
39	5115675	1.0	No	40P3M10BF-ROHLIX positioning GT
41	5115679-2	1.0	No	IMS SUP SHAFT FRONT
42	5115678-2	1.0	No	IMS SUP SHAFT BACK
43	5115682	1.0	No	IMS-SENSOR-BRKT positioning GT
44	5115683	1.0	No	IMS-MOTOR-PLATE-FF10 positioning GT
45	5115684	1.0	No	IMS-MOTOR-PLATE positioning GT
46	5115687	1.0	No	IMS-DRIVE-WASHER positioning GT
47	5117312	1.0	No	IMS SENSOR BRKT REAR Positioning GT
54	5112220	1.0	No	CRADLE RAIL Positioning GT
55	5112223	1.0	No	IMS RAIL Positioning GT
57	5112203	1.0	No	FRAME-PIPE-R Positioning GT
58	5112204	1.0	No	FRAME-PIPE-L Positioning GT
59	5167153	1.0	No	FRAME RAIL MACHINING positioning GT
60	5167123	1.0	No	PLATE RAIL MACHINING positioning GT

Table 2-4 TOP FRAME ASSY (5127782 BOM, Rev 24) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
61	5113114	1.0	No	ENCODER SUP CRADLE positioning GT
62	5113118	1.0	No	BRACE MID positioning GT
63	5113119	1.0	No	BRACE REAR positioning GT
64	5116868	1.0	No	IMS RAIL PLATE FRONT
65	5116870	1.0	No	IMS RAIL PLATE MID
66	5116871	1.0	No	IMS RAIL PLATE REAR
67	5113230	1.0	No	GTCN BRKT positioning GT
68	5113241	2.0	No	REIN MID positioning GT
69	5117361	1.0	No	TPB BRKT
70	5113242	1.0	No	STOPPER positioning GT
71	5113243	1.0	No	TOP PLATE BACK positioning GT
72	5117371	1.0	No	TOP PLATE MID
73	5113244	1.0	No	TOP PLATE FRONT positioning GT
74	5116958	1.0	No	ENCODER SUP IMS Positioning GT
75	5117356	1.0	No	CABLE DUCT M
76	5113247	1.0	No	CABLE DUCT B positioning GT
77	5117355	1.0	No	CABLE-GUIDE-BRKT
78	5118631	2.0	No	STOPPER SHAFT
79	5117982	1.0	No	WHEEL-SCREW-BRKT-REAR Positioning GT
80	5117983	1.0	No	REAR-WHEEL-ADJUST-BRKT Positioning GT
81	5117984	2.0	No	WHEEL-SCREW-BRKT-FRONT Positioning GT
82	5117985	1.0	No	WHEEL-SPRING-STAY Positioning GT
83	5117163	2.0	No	REIN MID SPACER
84	5122747	1.0	No	NUT-PLATE positioning GT
88	5126275	1.0	No	TOP-FRONT-BRKT-3L positioning GT
90	5182656	1.0	No	Linear Actuator
91	5112925	1.0	No	IMS BASE GT positioning GT
92	5112921	2.0	No	UPPER SUP ADD positioning GT

Table 2-4 TOP FRAME ASSY (5127782 BOM, Rev 24) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
93	5112926	1.0	No	IMS CABLE HOLDER 1 positioning GT
94	5112927	1.0	No	IMS STOPPER positioning GT
95	5116960	4.0	No	CUMFOLLOWER SPACER
96	5116947	1.0	No	ADDLINK-NUT Positioning GT
97	5124243	2.0	No	IMS BASE SUP Positioning GT
98	5121568	1.0	No	LINK BAR SUPPORT Positioning GT
112	N9406ES	5.0	No	SCREW SARA 4X6 Fe Zn
130	N9610SJ	3.0	No	SET SCREW
141	N9820YX	9.0	No	SCREW HEXAGON SOCKET 8 MM 20 MM CAP SCREW 8X20 FE
147	5118590	4.0	No	SPACER 40 PILLOW
149	5112471-3	1.0	No	GTCN BOARD ASSY Positioning GT
160	N9121CS	1.0	No	NUT M12 3SYU FeZn
171	5132756	1.0	No	SHIM-IMS-MOTOR positioning GT
172	5183123	1.0	No	SPACER ROHLIX ACTUATOR
174	5177564	1.0	No	IMS STOPPER 1700 positioning GT
175	N9616YX	2.0	No	SCREW HEXAGON SOCKET 6 MM 16 MM CAP SCREW 6X16 FE
180	N9310ZX	4.0	No	SCREW HEXAGON SOCKET 3 MM 10 MM CAP SCREW 3X10 FE
182	2299437	9.0	No	U0472AA_FLANGE_SOCKET
200	5135604	1.0	No	SHAFT LABEL Positioning GT
211	2379698	0.001	No	U0064BA_EDGING_TAPE

Table 2-4 TOP FRAME ASSY (5127782 BOM, Rev 24) (Continued)

FRU = YES

Figure 2-36 CRADLE STEEPING MOTAR CLUTCH ASSY Positioning GT (5127493ADW, sheet 1, rev 3)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	5123238	1.0	Yes	BELT 375 UP5M25 Positioning GT
2	5118348	1.0	No	PULLEY SUP BASE
3	2299437	7.0	No	U0472AA_FLANGE_SOCKET
5	5127775	1.0	No	CLUTCH ASSY
6	5127776	1.0	No	CRADLE MOTOR ASSY Positioning GT

Table 2-5 CRADLE STEEPING MOTOR CLUTCH 2000 ASSY Positioning GT (5127493 BOM, Rev 3)

1.1.3 FOOT SW ASSY (5126796, 5126796-2)

FRU = NO

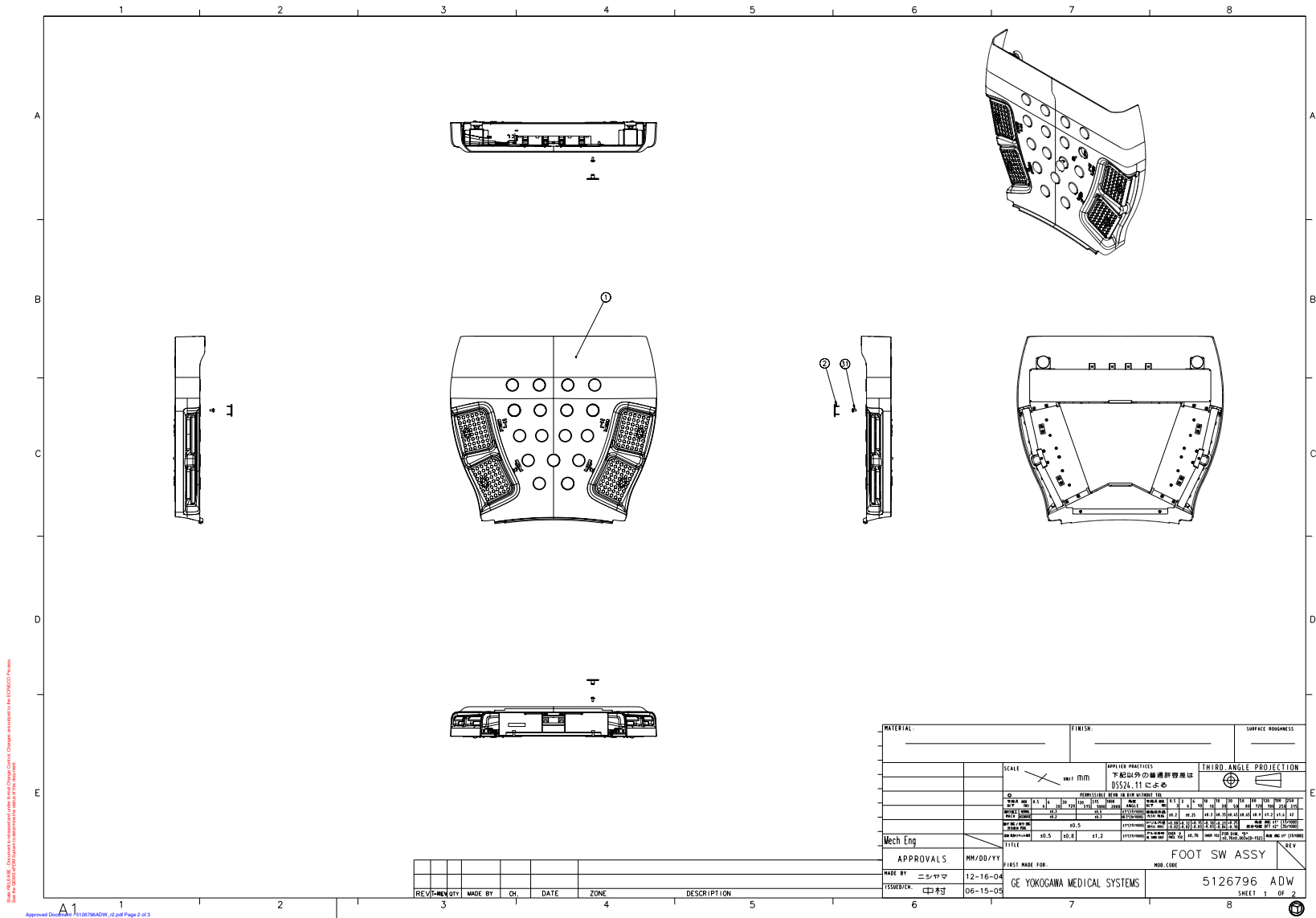


Figure 2-37 Foot Switch Assembly (5126796ADW, sheet 1, rev 2)



MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	5115661-2	1.0	Yes	FS COVER positioning GT
2	5118050-2	3.0	Yes	FS COVER LID positioning GT
3	5120586-2	4.0	Yes	pedal cover positioning GT
4	5120587-2	4.0	Yes	pedal cover lid positioning GT
5	U0023AA	16.0	Yes	SCREW
7	U0233AA	4.0	Yes	NUT HEXAGON 6 MM M3-17-9
9	P9155EQ	16.0	Yes	STOP RUBBER
18	2174553	4.0	Yes	P9235JU PEDAL
22	5121380	4.0	Yes	FOOT SW END ASSY Positioning GT
23	5121381	1.0	Yes	FOOT SW L CABLE ASSY Positioning GT
24	5121382	1.0	Yes	FOOT SW R CABLE ASSY Positioning GT
25	N9408JS	24.0	Yes	SCREW NABE 4X8 Fe Zn
27	N9412ZX	16.0	Yes	SCREW HEXAGON SOCKET 4 MM 12 MM SCREW
28	N9520JS	8.0	Yes	SCREW PAN HEAD 5 MM 20 MM SCREW osX 5X20 FE ZN
31	N9510HR	2.0	Yes	SCREW TP
32	N9612YX	10.0	Yes	SCREW HEXAGON SOCKET 6 MM 12 MM CAP SCREW 6X12 FE
33	N9630YX	2.0	Yes	SCREW HEXAGON SOCKET 6 MM 30 MM CAP SCREW 6X30 FE
34	N9510ZX	16.0	Yes	SCREW HEXAGON SOCKET 5 MM 10 MM CAP SCREW 5X10 FE
37	N9501WS	8.0	Yes	WASHER PLAIN - LARGE 5.3 MM 10 MM WASHER M5
38	N9401WS	16.0	Yes	WASHER PLAIN _ LARGE 4.3 MM 9 MM WASHER 4.3x9x0.8 FE ZN
39	46-297971P1	4.0	Yes	TERMINAL GROUND LUG FOR #4 TO 2/0 CABLE. ALUMINUM BASE, .25~ BOLT HOLE. SCREW COLORED GREEN FOR GROUND IDENT IFICATION. (DWG REV 12/92, JFK - OK) .
40	U0030AC	4.0	Yes	CLAMP OR CLIP CLAMP
41	5127867	1.0	Yes	TABLE GND TO RACEWAY GND BAR CABLE Positioning GT
42	U0507AC	3.12	Yes	PROTECTOR EDGING CE-016
6	U0004AQ	8.0	No	SPRING COMPRESSION 47 MM SC-085-120-0470 4017
8	5123636	2.0	No	ADJUSTER Foot U0165AE

Table 2-6 FOOT SW ASSY (5126796-2 BOM, Rev 4)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
10	P9225EN	8.0	No	SHAFT
11	P9185ES	8.0	No	SW_BRKT2
12	P9225EM	4.0	No	PLATE_WASHER
13	P9225EP	8.0	No	STOPPER
14	P9225ER	4.0	No	SW_BRKT1
15	P9225EW	4.0	No	SW_BRKT3
16	2174552	4.0	No	P9235JT SW PEDAL
17	5118061	1.0	No	FS COVER SUP PLATE positioning GT
19	5112987	1.0	No	BASE FRAME Positioning GT
20	5112989	4.0	No	SW BASE Positioning GT
21	5112992	1.0	No	GROUND Positioning GT
26	N9512LK	3.0	No	SCREW BIND UTIBATUKI 5X12 BsNi
30	N9600YX	4.0	No	SCREW HEXAGON SOCKET 6 MM 100 MM CAP SCREW 6X100 FE
35	N9514ZX	2.0	No	SCREW HEXAGON SOCKET 5 MM 14 MM
36	N9314JS	16.0	No	SCREW PAN HEAD 3 MM 14 MM SCREW osX 3X14 FE ZN

Table 2-6 FOOT SW ASSY (5126796-2 BOM, Rev 4) (Continued)

FRU = YES



Section 1.0 — GT 2000 Table (5121647)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
3	5133731	2.0	Yes	TAPE SWITCH 2 Positioning GT
4	5133732	2.0	Yes	TAPE SWITCH 3 Positioning GT
6	N9510HR	4.0	Yes	SCREW TP
8	N9501BS	2.0	Yes	NUT HEXAGON 5 MM NUT M5 FE ZN
1	5173720	1.0	No	RIM TOP FRONT positioning GT
2	5113608	2.0	No	TOP FRONT BRKT 1 positioning GT
5	5120573	1.0	No	TS_TOP_FRONT Positioning GT
7	N9514ZX	2.0	No	SCREW HEXAGON SOCKET 5 MM 14 MM
9	N9406ES	4.0	No	SCREW SARA 4X6 Fe Zn
10	5123637	2.0	No	Ball Catch U0441AE

Figure 2-40 RIM TOP COVER FRONT ASSY positioning GT (5127465-3 BOM, Rev 3)

1.1.5 TOP COVER REAR ASSY (5154872)

FRU = YES

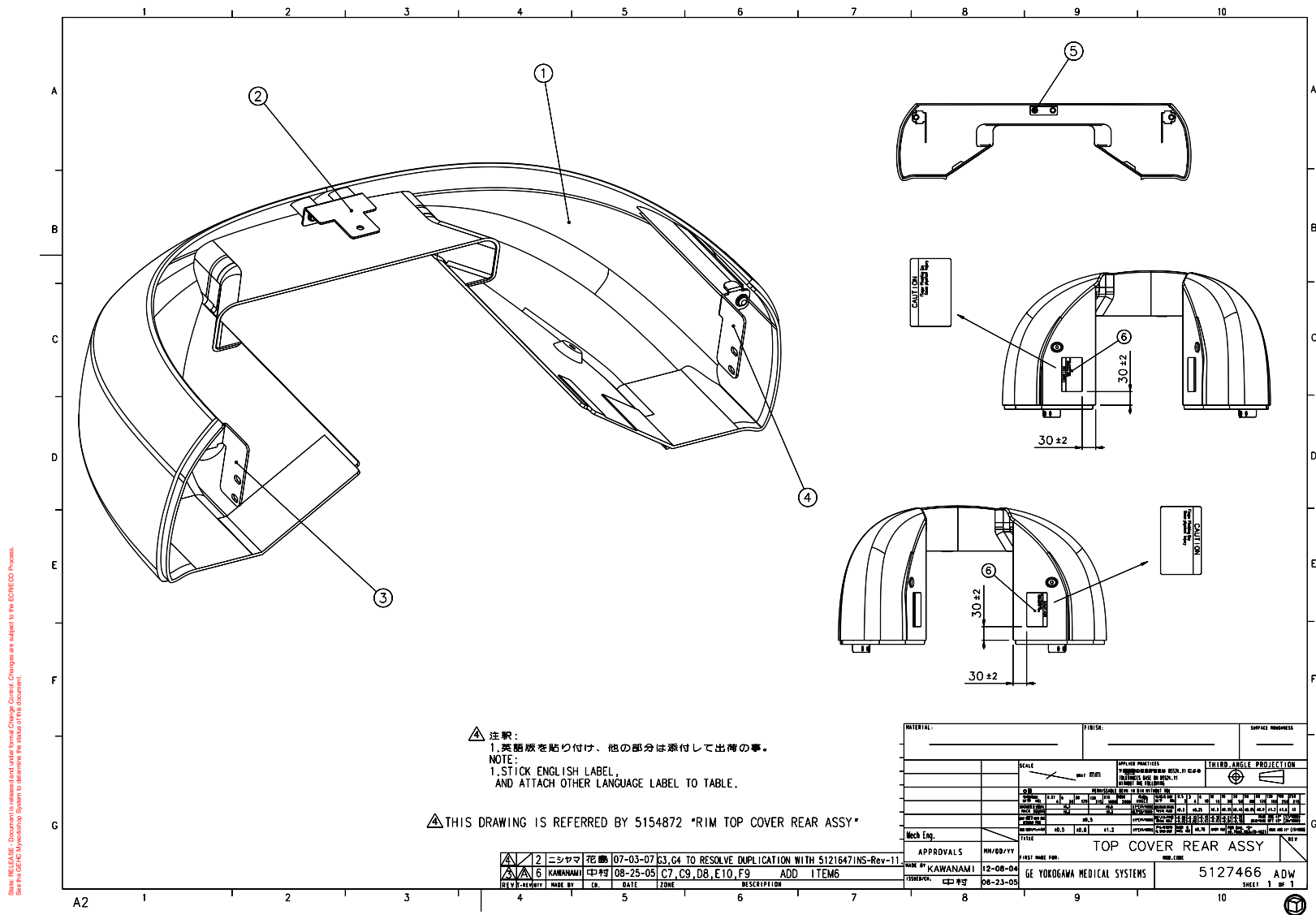


Figure 2-41 TOP COVER REAR ASSY (5127466ADW, sheet 1, rev 4)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
5	N9510HR	4.0	Yes	SCREW TP
6	5129092	1.0	Yes	FINGER PINCHING LABEL Positioning GT
1	5155285	1.0	No	RIM TOP COVER REAR positioning GT
2	5113614	1.0	No	TOP REAR BRKT 2 positioning GT
3	5113615	1.0	No	TOP REAR BRKT 3L positioning GT
4	5113616	1.0	No	TOP REAR BRKT 3R positioning GT

Table 2-7 TOP COVER REAR ASSY (5154872 BOM, Rev 2)

FRU = YES



Section 1.0 — GT 2000 Table (5121647)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
5	N9510HR	18.0	Yes	SCREW TP
6	2199652	1.0	Yes	P9235LZ LABEL
1	5112211-2	1.0	No	TOP SIDE COVER L Positioning GT
2	5113622	2.0	No	TOP COVER SIDE BRKT 9A positioning GT
4	5115665	1.0	No	FINGER PROTCTER L positioning GT

Table 2-8 TOP COVER SIDE L ASSY (5127467 BOM, Rev 3)

1.1.7 TOP COVER SIDE R ASSY (5127468)

FRU = YES

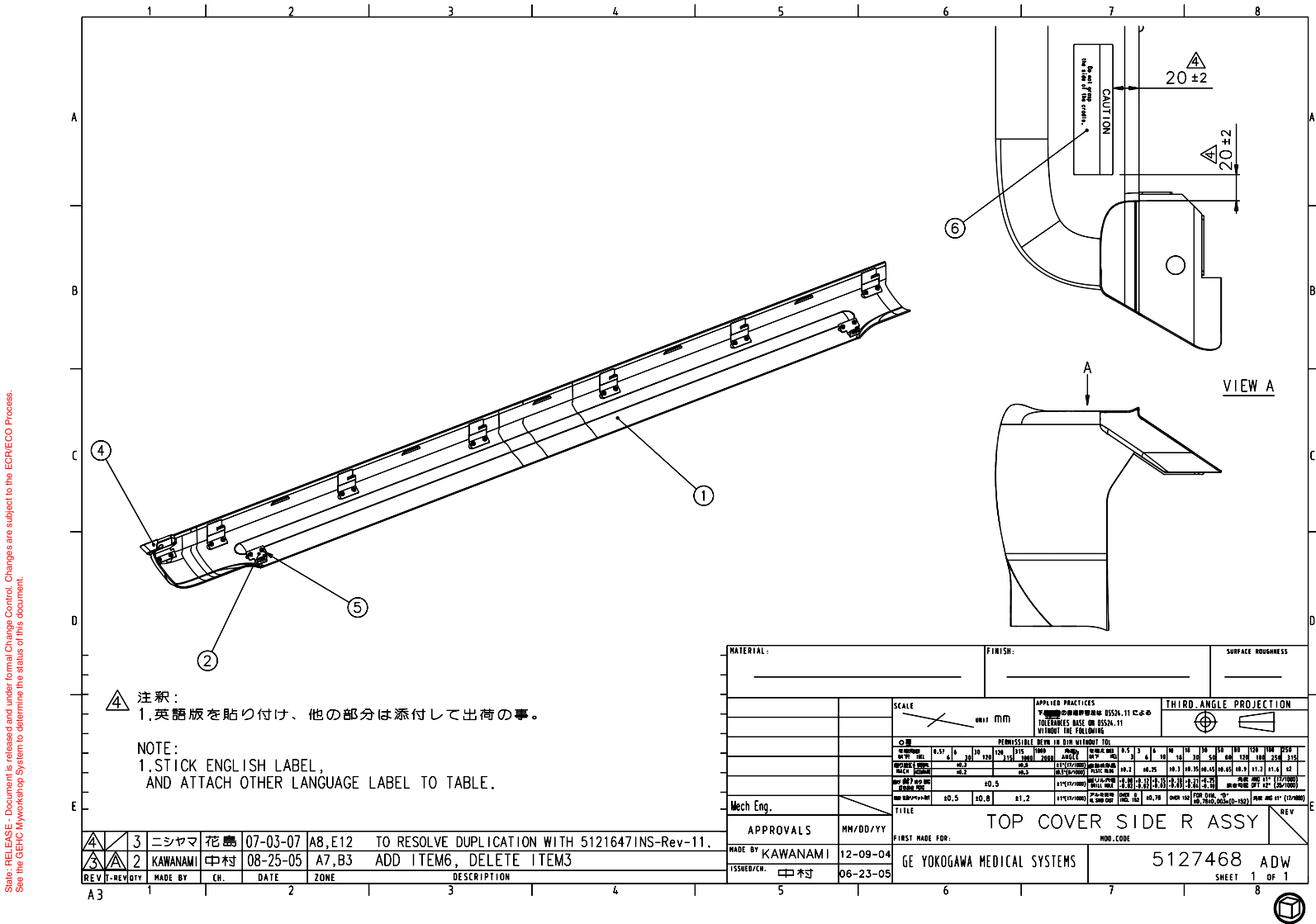


Figure 2-43 TOP COVER SIDE R ASSY (5127468ADW, sheet 1, rev 4)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
5	N9510HR	18.0	Yes	SCREW TP
6	2199652	1.0	Yes	P9235LZ LABEL
1	5112212-2	1.0	No	TOP SIDE COVER R Positioning GT
2	5113622	2.0	No	TOP COVER SIDE BRKT 9A positioning GT
4	5115664	1.0	No	FINGER PROTCTER R positioning GT

Table 2-9 TOP COVER SIDE R ASSY (5127468 BOM, Rev 3)

1.1.8 SIDE-PANEL-FRONT-LEFT-ASSY, Positioning, GT (5130814)

FRU = YES

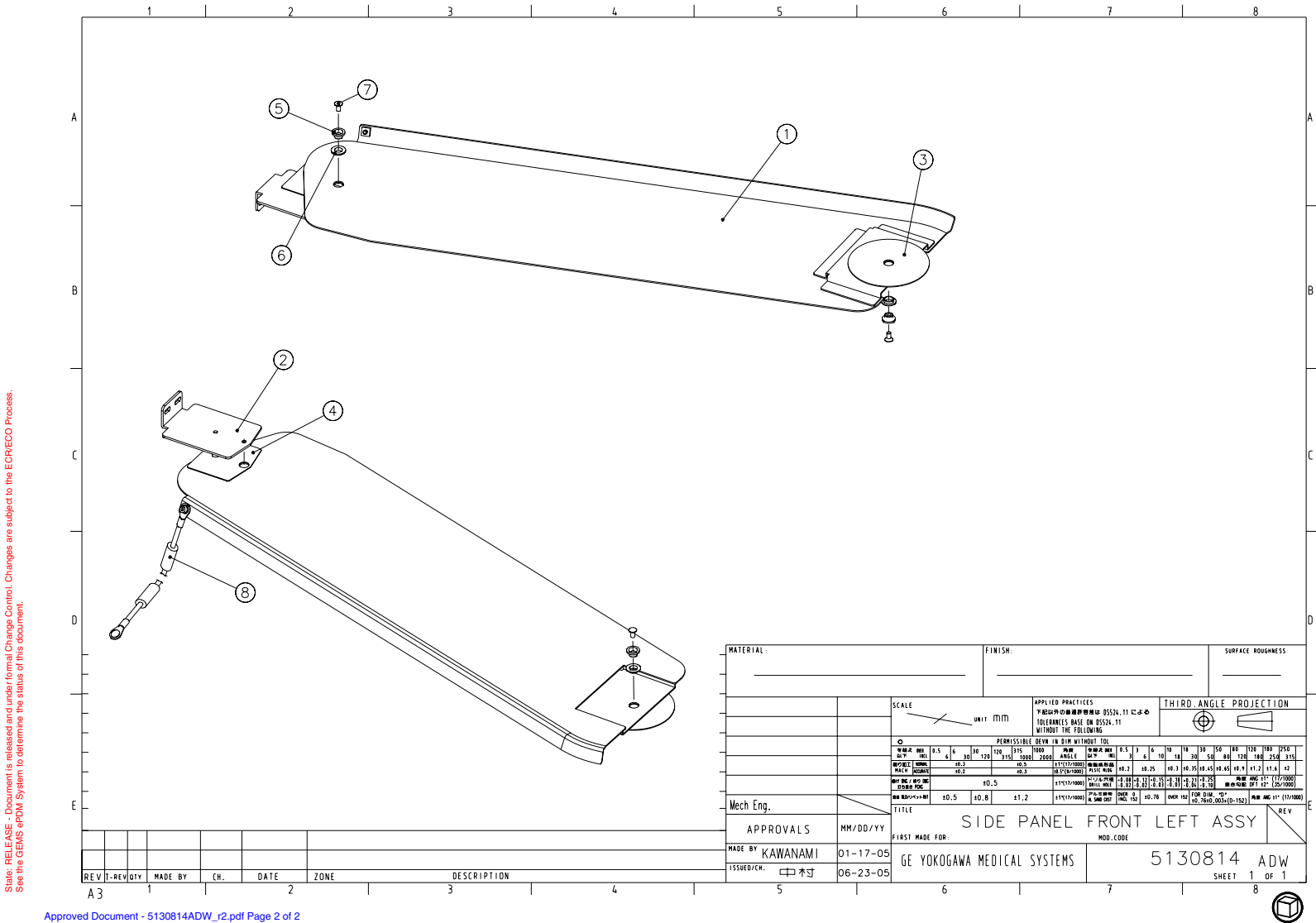


Figure 2-44 Side Panel, Front Left Assembly (5130814ADW, sheet 1, rev 2)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	5113110	1.0	Yes	SIDE PANEL FRONT LEFT positioning GT
8	5126552	1.0	Yes	GND WIRE COVER Positioning GT
10	U0037BA	0.3	Yes	TAPE VHB Y_4950 15M_M 11M
2	5118067	1.0	No	SIDE PANEL BRKT UPR F positioning GT
3	5120589	2.0	No	SIDE-PANNEL-SPCR-LWR positioning GT
4	5120590	2.0	No	SIDE PANEL SPCR UPR positioning GT
5	5112207	2.0	No	SIDE PANEL SHAFT
6	5123639	2.0	No	OILLESS BUSH Bearing Special Parts U0935AH
7	N9510ES	2.0	No	SCREW COUNTERSUNK HEA 5 MM 10 MM SCREW s[5X10 FE ZN
9	N9508ZX	1.0	No	SCREW CAP SCREW 5X8 Fe * CAP SCREW 5X8 FE

Table 2-10 SIDE-PANEL-FRONT-LEFT-ASSY, Positioning, GT (5130814 BOM, Rev 3)

1.1.9 SIDE-PANEL-REAR-LEFT-ASSY, Positioning, GT (5130815)
FRU = YES

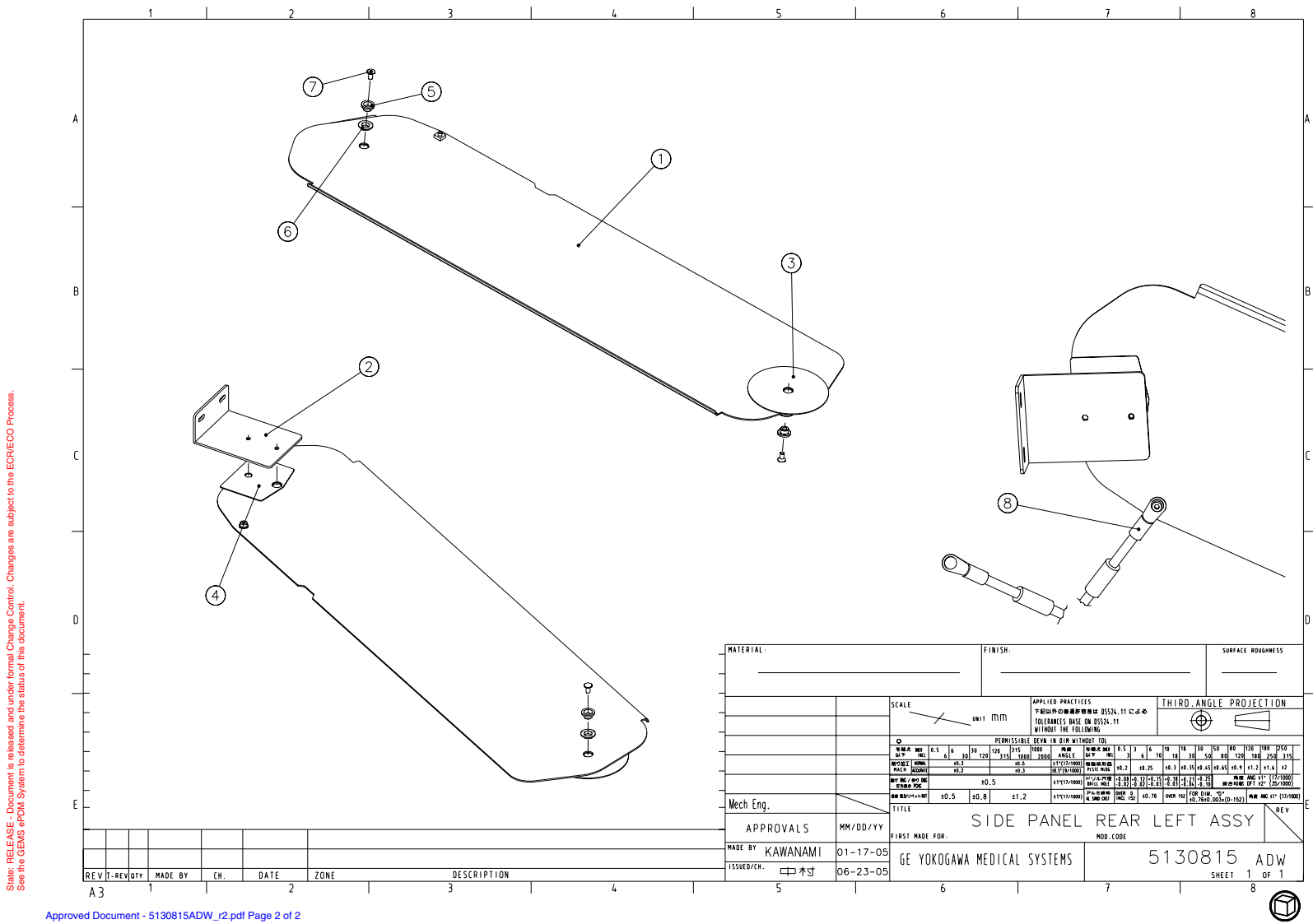


Figure 2-45 Side Panel, Rear Left Assembly (5130815ADW, sheet 1, rev 2)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	5113111	1.0	Yes	SIDE PANEL REAR LEFT positioning GT
8	5126552	1.0	Yes	GND WIRE COVER Positioning GT
10	U0037BA	0.3	Yes	TAPE VHB Y_4950 15M_M 11M
2	5118068	1.0	No	SIDE PANEL BRKT UPR R positioning GT
3	5120589	2.0	No	SIDE-PANNEL-SPCR-LWR positioning GT
4	5120590	2.0	No	SIDE PANEL SPCR UPR positioning GT
5	5112207	2.0	No	SIDE PANEL SHAFT
6	5123639	2.0	No	OILLESS BUSH Bearing Special Parts U0935AH
7	N9510ES	2.0	No	SCREW COUNTERSUNK HEA 5 MM 10 MM SCREW s[5X10 FE ZN
9	N9508ZX	1.0	No	SCREW CAP SCREW 5X8 Fe * CAP SCREW 5X8 FE

Table 2-11 SIDE-PANEL-REAR-LEFT-ASSY, Positioning, GT (5130815 BOM, Rev 3)

1.1.10 SIDE-PANEL-FRONT-RIGHT-ASSY, Positioning, GT (5130816)

FRU = YES

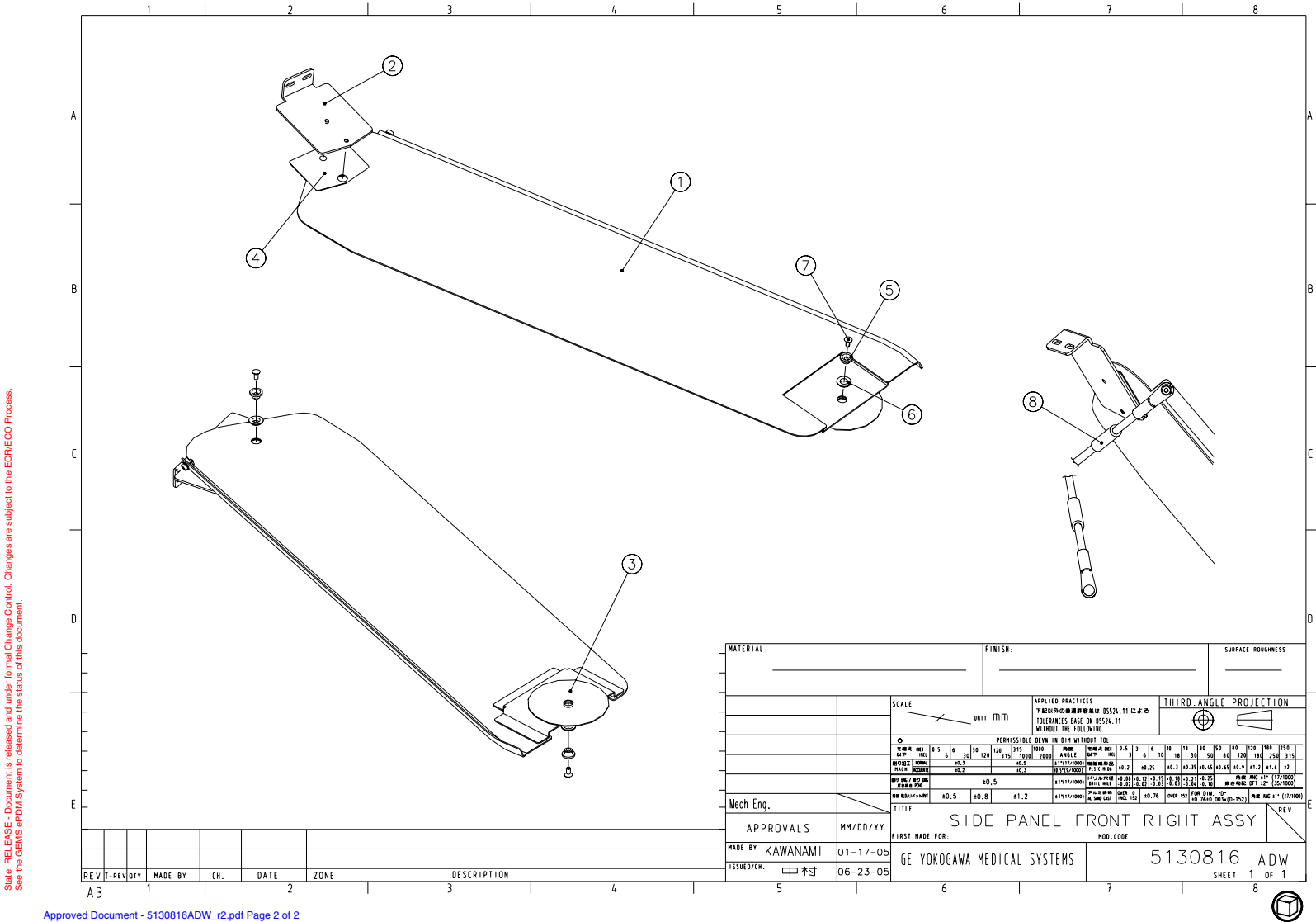


Figure 2-46 Side Panel, Front Right Assembly (5130816ADW, sheet 1, rev 2)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	5113112	1.0	Yes	RIGHT FRONT SIDE PANEL positioning GT
8	5126552	1.0	Yes	GND WIRE COVER Positioning GT
10	U0037BA	0.3	Yes	TAPE VHB Y_4950 15M_M 11M
2	5118067	1.0	No	SIDE PANEL BRKT UPR F positioning GT
3	5120589	2.0	No	SIDE-PANNEL-SPCR-LWR positioning GT
4	5120590	2.0	No	SIDE PANEL SPCR UPR positioning GT
5	5112207	2.0	No	SIDE PANEL SHAFT
6	5123639	2.0	No	OILLESS BUSH Bearing Special Parts U0935AH
7	N9510ES	2.0	No	SCREW COUNTERSUNK HEA 5 MM 10 MM SCREW s[5X10 FE ZN
9	N9508ZX	1.0	No	SCREW CAP SCREW 5X8 Fe * CAP SCREW 5X8 FE

Table 2-12 SIDE-PANEL-FRONT-RIGHT-ASSY, Positioning, GT (5130816 BOM, Rev 3)

1.1.11 SIDE-PANEL-REAR-RIGHT-ASSY, Positioning, GT (5130817)

FRU = YES


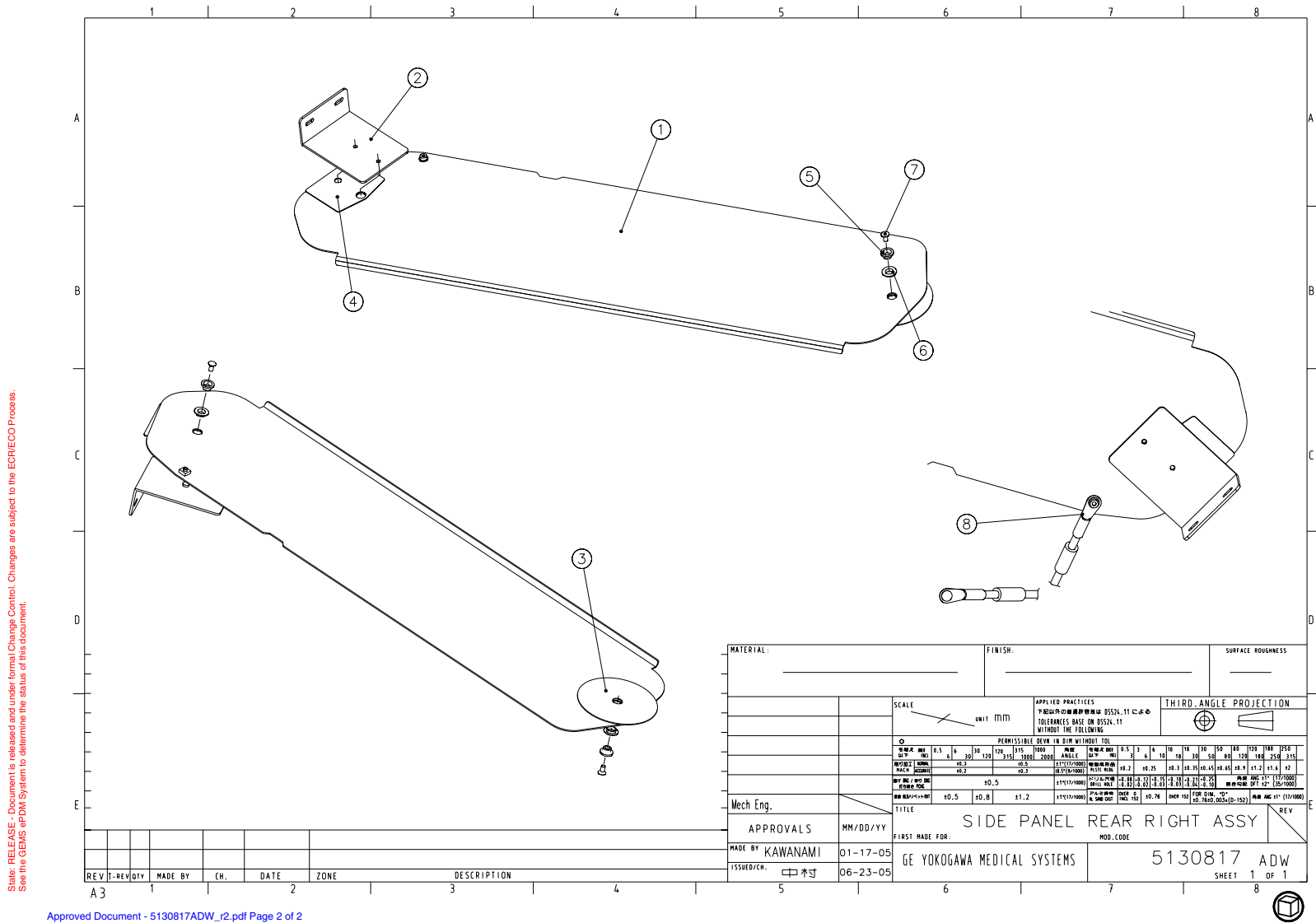


Figure 2-47 Side Panel, Right Rear Assembly (5130817ADW, sheet 1, rev 2)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	5113113	1.0	Yes	SIDE PANEL REAR RIGHT positioning GT
8	5126552	1.0	Yes	GND WIRE COVER Positioning GT
10	U0037BA	0.3	Yes	TAPE VHB Y_4950 15M_M 11M
2	5118068	1.0	No	SIDE PANEL BRKT UPR R positioning GT
3	5120589	2.0	No	SIDE-PANNEL-SPCR-LWR positioning GT
4	5120590	2.0	No	SIDE PANEL SPCR UPR positioning GT
5	5112207	2.0	No	SIDE PANEL SHAFT
6	5123639	2.0	No	OILLESS BUSH Bearing Special Parts U0935AH
7	N9510ES	2.0	No	SCREW COUNTERSUNK HEA 5 MM 10 MM SCREW s[5X10 FE ZN
9	N9508ZX	1.0	No	SCREW CAP SCREW 5X8 Fe * CAP SCREW 5X8 FE

Table 2-13 SIDE-PANEL-REAR-RIGHT-ASSY, Positioning, GT (5130817 BOM, Rev 3)

FRU = YES

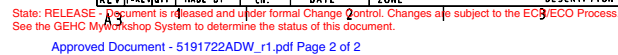

Section 1.0 — GT 2000 Table (5121647)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
6	5133730	1.0	Yes	TAPE SWITCH 1 Positioning GT
7	5133732	1.0	Yes	TAPE SWITCH 3 Positioning GT
8	5181038	1.0	Yes	TAPE SW-5 Positioning GT
10	U0087AA	3.0	Yes	SCREW
12	N9406JS	1.0	Yes	SCREW NABE 4X6 FE ZN
13	N9510HR	8.0	Yes	SCREW TP
1	5173144	1.0	No	RIM IMS COVER SIDE L Positioning GT
2	5175786	1.0	No	SUPPORT BRKT IMS SIDE L Positioning GT
3	5113597	1.0	No	IMS COVER BRKT 8 positioning GT
4	5113600	1.0	No	IMS SIDE BRKT 1L positioning GT
5	5135990	1.0	No	IMS-COVER-BRKT-25L positioning GT
9	5120578	1.0	No	TS_IMS_LEFT Positioning GT
11	U0246AA	3.0	No	NUT a.p NTM4

Table 2-14 RIM IMS COVER SIDE LH ASSY positioning GT (5191722 BOM, Rev 1)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
6	5133730	1.0	Yes	TAPE SWITCH 1 Positioning GT
7	5133732	1.0	Yes	TAPE SWITCH 3 Positioning GT
8	5181038	1.0	Yes	TAPE SW-5 Positioning GT
10	U0087AA	3.0	Yes	SCREW
12	N9406JS	1.0	Yes	SCREW NABE 4X6 FE ZN
13	N9510HR	12.0	Yes	SCREW TP
1	5173123	1.0	No	RIM IMS COVER SIDE R Positioning GT
2	5176056	3.0	No	SUPPORT BRKT IMS SIDE R Positioning GT
3	5113597	1.0	No	IMS COVER BRKT 8 positioning GT
4	5113601	1.0	No	IMS SIDE BRKT 1R positioning GT
5	5135987	1.0	No	IMS-COVER-BRKT-25R positioning GT
9	5120576	1.0	No	TS_IMS_RIGHT Positioning GT
11	U0246AA	3.0	No	NUT a.p NTM4

Table 2-15 RIM IMS COVER SIDE RH ASSY positioning GT (5191736 BOM, Rev 1)

1.2 Table Ship Collector (5129558)

FRU = NO

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2115990-3	1.0	Yes	LCC HSA METALLESS SUPINE CORNAL HEADHOLDER ASM WITH RATING PLATE AND SELF-CATCHING LATCH
2	2115993-2	1.0	Yes	LCC METALLESS CRADLE EXTENDER HSA METALLESS FOOT EXTENDER WITH RATING PLATE AND SELF-CATCHING LATCH
3	2115996-3	1.0	Yes	AXIAL HEADHOLDER ASM WITH RATING PLATE AND SELF-CATCHING LATCH
4	2169679	1.0	Yes	P9150TS BODY STRAP A
5	2169680	1.0	Yes	P9150TT BODY STRAP B
6	2169681	1.0	Yes	P9150TU BODY STRAP C
8	2331933	1.0	Yes	METALLESS COMPATIBLE PHANTOM HOLDER
11	P9140HD	1.0	Yes	BALLOON CLAMP
12	P9140TL	1.0	Yes	ARM SUPPORT ASSY
13	2152502	1.0	Yes	P9150SN BODY STRAP A
14	2152503	1.0	Yes	P9150SP BODY STRAP B
15	2152504	1.0	Yes	P9150SQ BODY STRAP C
20	5193132	1.0	Yes	FIELD EMERGENCY KIT Positioning GT
9	5121757-2	1.0	No	CRADLE PAD GT LONG Positioning GT
10	5122945-2	1.0	No	EXTENDER PAD 400 Positioning GT
16	5131377	1.0	No	MSDS GT TABLE Positioning GT

Table 2-16 Table Ship Collector (5129558 BOM, Rev 8)

Section 2.0 — GT 1700 Table (5122080)

FRU = NO (see [Figure 2-1, beginning on page 30.](#))

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
	5143974	1.0	Yes	GTCN SET
	5164627	1.0	Yes	TOP COVER ASSY 1700 Positioning GT
	5168472-2	1.0	Yes	IMS DRIVE REPLACEMENT SET Positioning GT
	5192989	1.0	Yes	IMS MOTOR ASSY Collector
	5237567	1.0	Yes	Lubricant U0303BD-2
2	46-303268P1	1.0	Yes	PRODUCT LOCATOR CARD FORM FOR RATING PLATES, PRINTED BY DUPLITECH
1	5134845	1.0	No	TABLE ASSY MID positioning GT {For further detail, see Section 2.1, beginning on page 112. }
3	5138509	1.0	No	SHIP COLLECTOR GT TABLE 1700 {For further detail, see Section 2.2, beginning on page 162. }

Table 2-17 Positioning GT MIDDLE TABLE ASSY (5122080 BOM, Rev 18)

2.1 TABLE ASSY MID positioning GT (5134845)

FRU = NO

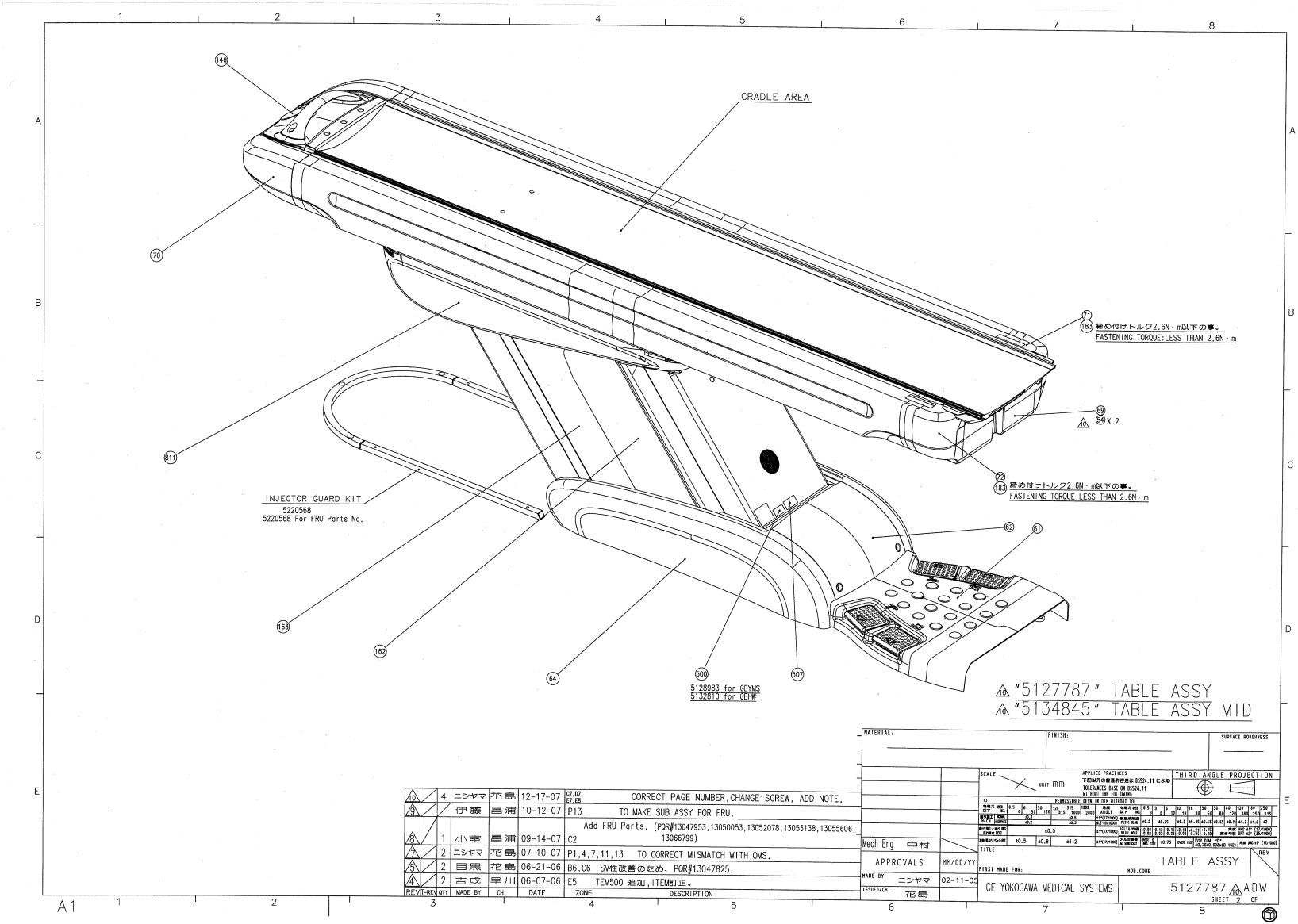


Figure 2-50 Table Assembly (5127787ADW, sheet 2, rev 10)

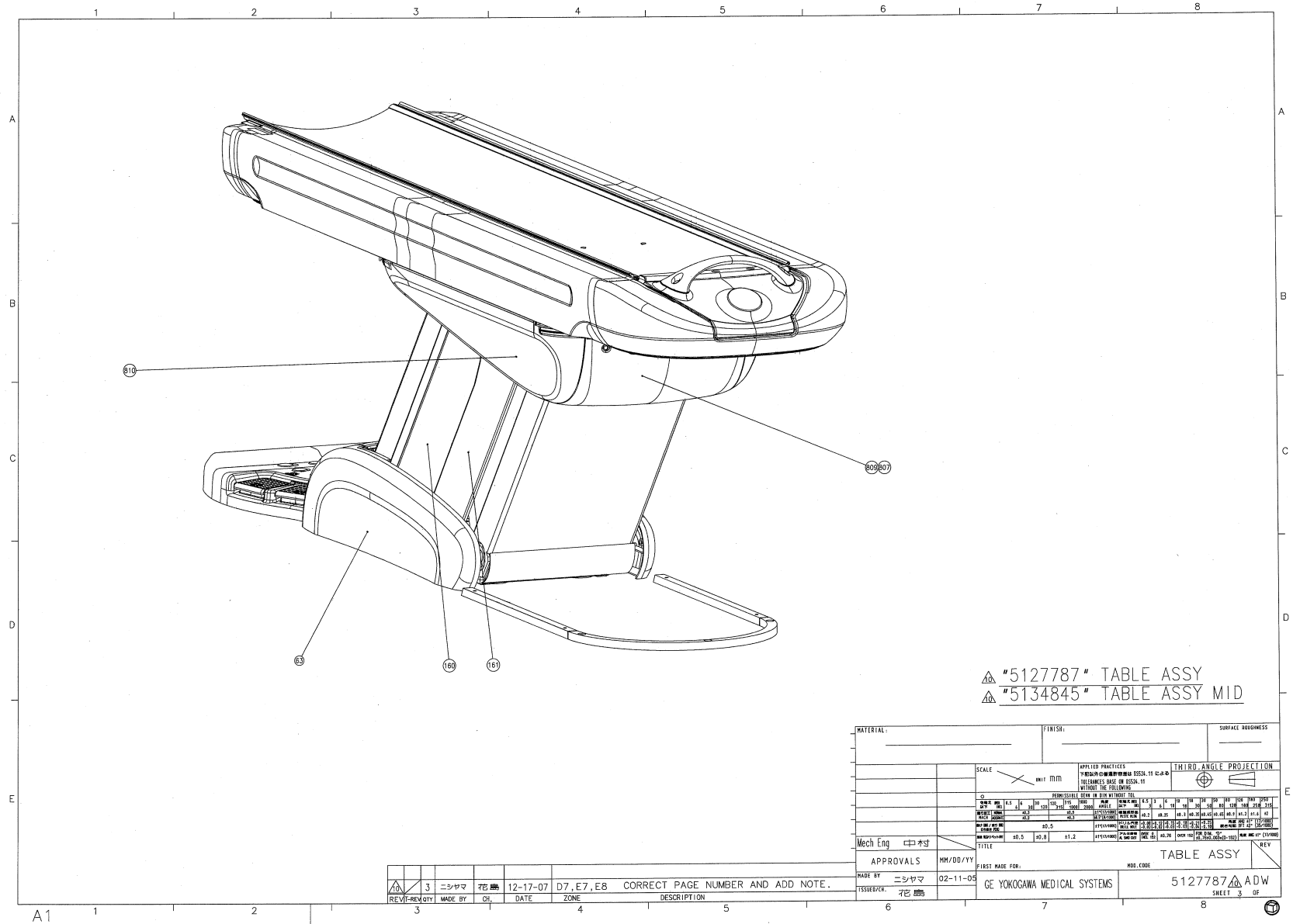
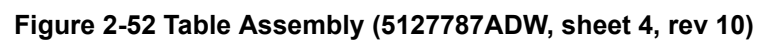
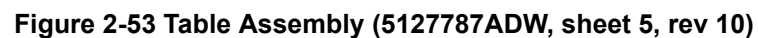
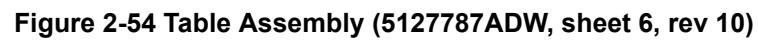


Figure 2-51 Table Assembly (5127787ADW, sheet 3, rev 10)







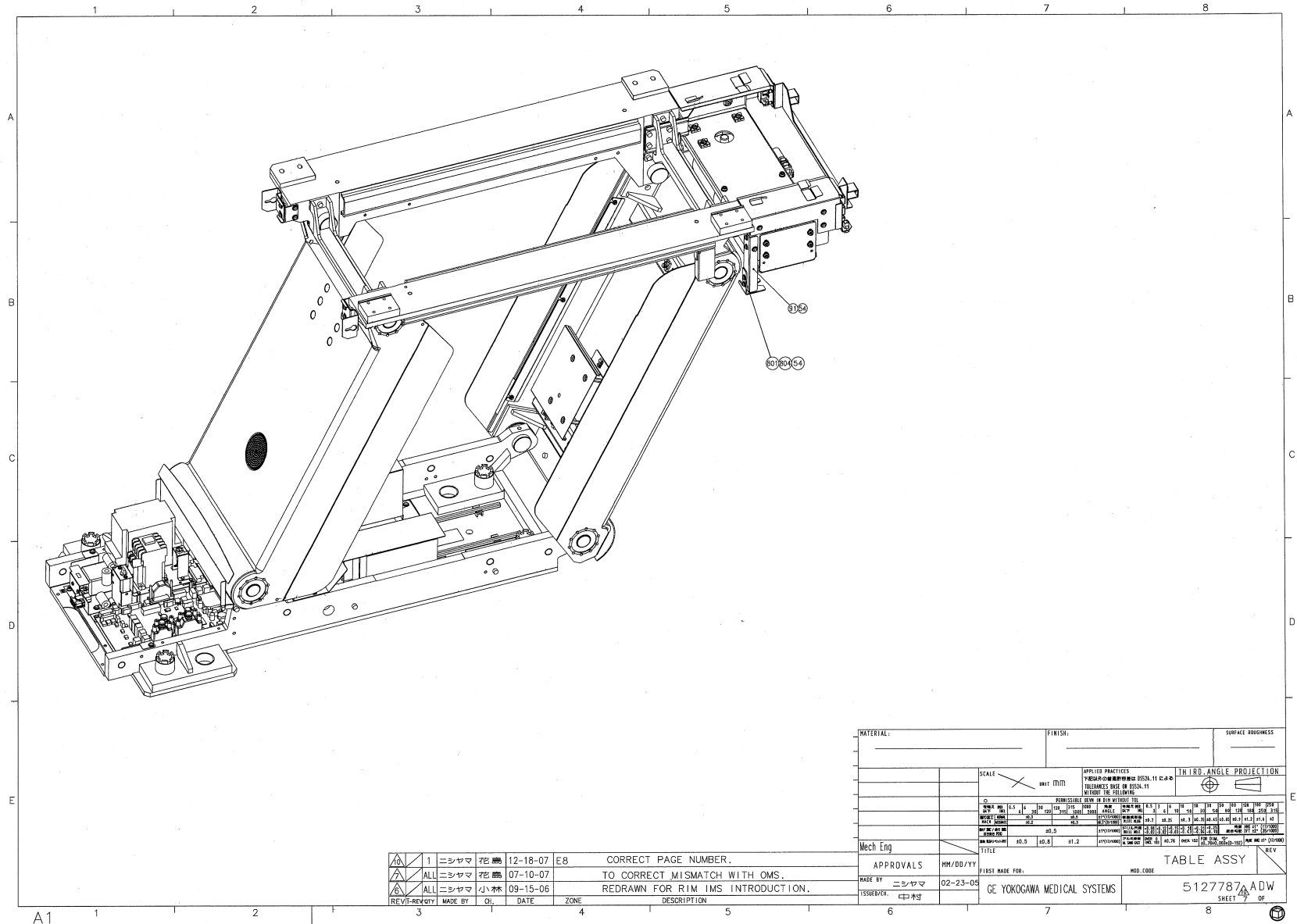


Figure 2-55 Table Assembly (5127787ADW, sheet 7, rev 10)

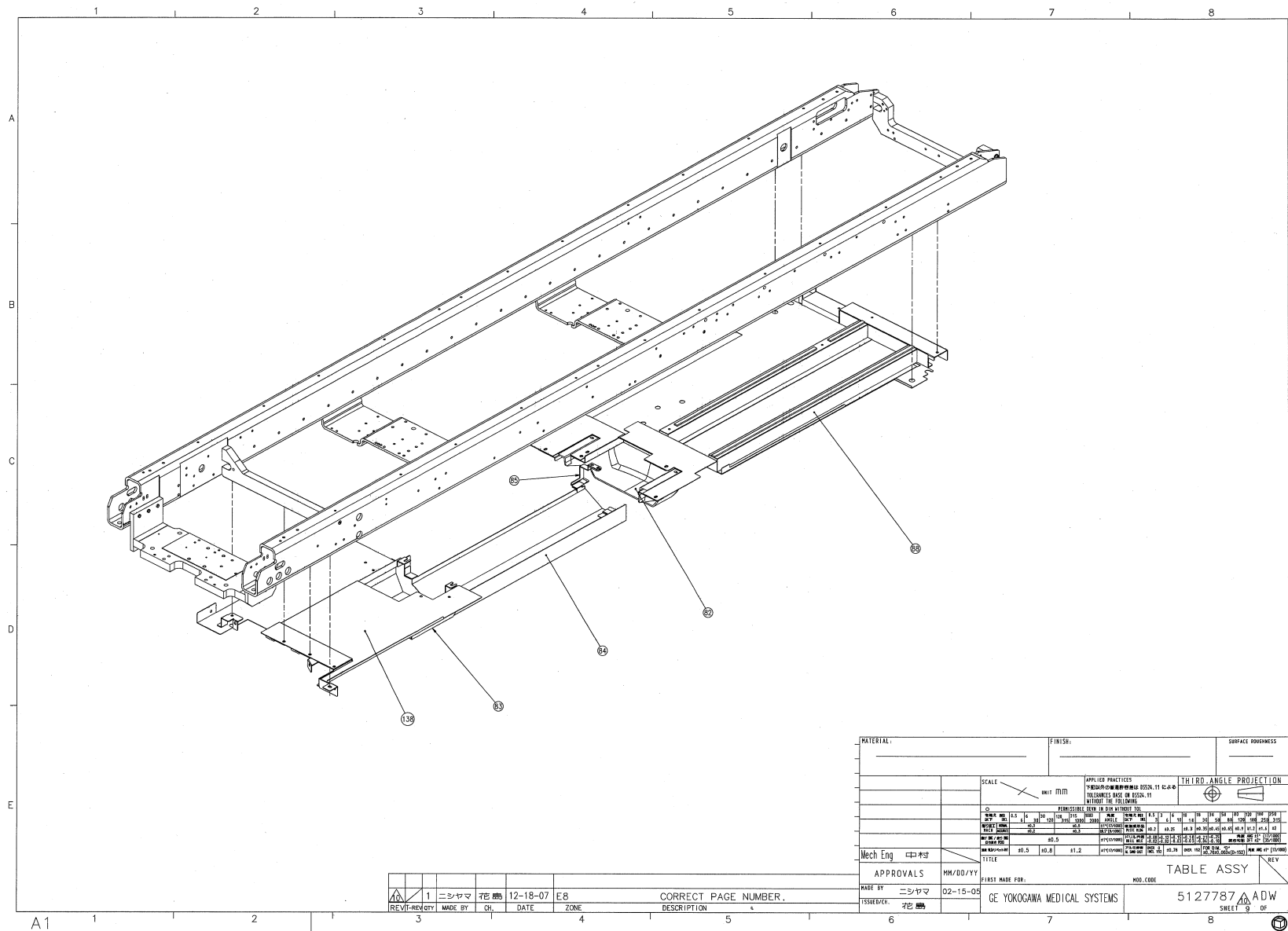


Figure 2-56 Table Assembly (5127787ADW, sheet 9, rev 10)


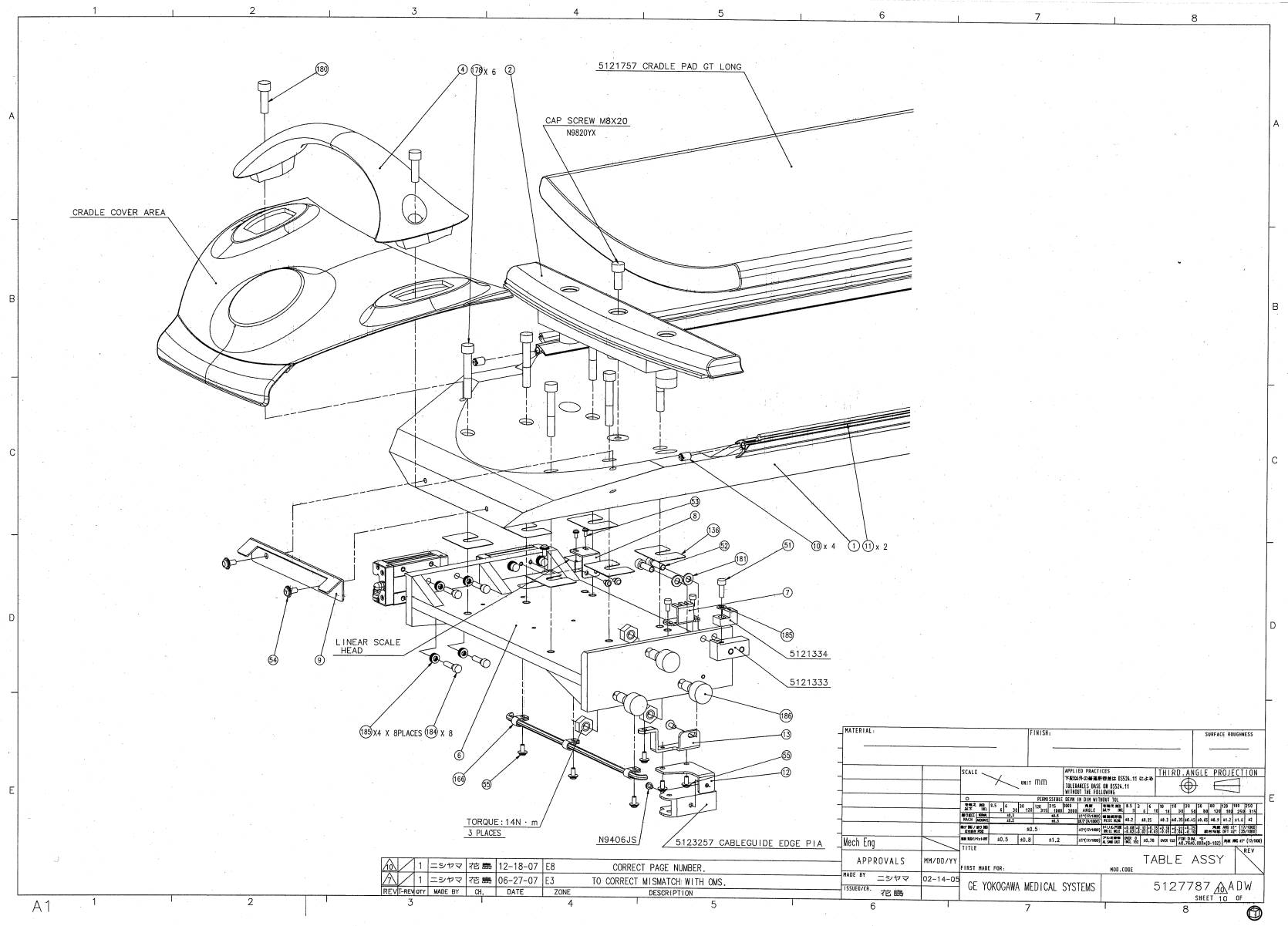


Figure 2-57 Table Assembly (5127787ADW, sheet 10, rev 10)

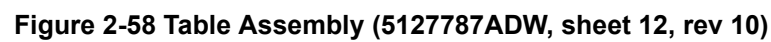

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
	5127282	1.0	Yes	HYDRAULIC CYLINDER ASSY Positoning GT
	5127283	1.0	Yes	HYDRAULIC RETURN HOSE ASSY Positoning GT
	5127284	1.0	Yes	PROPORTIOAL VALVE ASSY Positoning GT
	5127285	1.0	Yes	SHUT OFF VALVE ASSY Positoning GT
	5127286	1.0	Yes	HYDRAULIC HIGH PRESUURE HOSE Positoning GT
	5127287	1.0	Yes	ELEVATION HYDRAULIC PUMP Positoning GT
	5136012	1.0	Yes	IMS FINGER PINCHING LABEL positioning GT
	5138908	1.0	Yes	TABLE SUPPORT JIG ASSY
	5164938	1.0	Yes	TAPE SWITH 4 ASSY Positioning GT
	5166082-2	1.0	Yes	RIM BASE COVER LH REPLACEMENT SET positioning GT
	5166242-2	1.0	Yes	RIM BASE COVER RH REPLACEMENT SET positioning GT
	5181038	1.0	Yes	TAPE SW-5 Positioning GT
	5191991	1.0	Yes	RIM IMS COVER SET positioning GT
	5220568	1.0	Yes	INJECTOR GUARD KIT Positioning GT
1	5173055	1.0	Yes	CRADLE MID ASSY Positioning GT
10	5123613	4.0	Yes	PLASTIC SET SCREW Positioning GT
11	5133822	2.0	Yes	ATTACHMENT RAIL MID positioning GT
25	5117286	1.0	Yes	GTCN J7 TO J221 Positioning GT
26	5114080	1.0	Yes	PE Positioning GT
27	5114118	1.0	Yes	GTCBJ3 TO J104 Positioning GT
28	5114120	1.0	Yes	GTCBJ10 TO SPEAKER Positioning GT
29	5114122	1.0	Yes	GTCBJ11 TO J200 Positioning GT
30	5114124	1.0	Yes	GTCBJ12 TO J201 Positioning GT
31	5114125	1.0	Yes	GTCBJ6 TO IMSDRIVER Positioning GT
32	5114126	1.0	Yes	GTCBJ5 TO CRADLED RIVER Positioning GT
33	5114128	1.0	Yes	CRADLED RIVER TO J204
34	5114129	1.0	Yes	IMSDRIVER TO J203 Positioning GT

Table 2-18 TABLE ASSY MID positioning GT (5134845 BOM, Rev 31)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
35	5114135	1.0	Yes	GTCN TO J207 Positioning GT
36	5114146	1.0	Yes	GTCN TO J218 Positioning GT
37	5114148	1.0	Yes	GTCN TO LINERSCALE Positioning GT
38	5118961	1.0	Yes	LMIT SW CABLE Positioning GT
39	5119130	1.0	Yes	J102 TO J112 Positioning GT
40	5120569	2.0	Yes	4PIN470DUM Positioning GT
41	5120570	1.0	Yes	6PIN470DUM Positioning GT
42	5120572	1.0	Yes	9PIN470DUM Positioning GT
43	5120574	1.0	Yes	TS_TOP_FRONT Positioning GT
44	5120577	1.0	Yes	TS_IMS_RIGHT Positioning GT
45	5120579	1.0	Yes	TS_IMS_LEFT Positioning GT
46	5120581	1.0	Yes	TS_IMS_REAR Positioning GT
47	5120773-2	1.0	Yes	GTCB TO GTTS Cable For EMC Positioning GT
48	5114146-2	1.0	Yes	LIMIT SW Positioning GT
50	N9408ZX	4.0	Yes	SCREW HEXAGON SOCKET 4 MM 8 MM CAP SCREW 4X8 FE
52	N9620YX	2.0	Yes	SCREW HEXAGON SOCKET 6 MM 20 MM CAP SCREW 6X20 FE
54	N9510HR	162.0	Yes	SCREW TP
55	N9408HR	5.0	Yes	SCREW TP
57	5127780-2	1.0	Yes	IMS CABLE BEAR ASSY
62	5127458-2	1.0	Yes	BASE COVER CTR ASSY positioning GT
69	5127465-3	1.0	Yes	RIM TOP COVER FRONT ASSY positioning GT {For further detail, see Section 1.1.4, beginning on page 90. }
70	5133505	1.0	Yes	TOP COVER REAR-1700 ASSY positioning GT {For further detail, see Section 2.1.11, beginning on page 156. }
71	5133507	1.0	Yes	TOP SIDE COVER L-1700 ASSY positioning GT {For further detail, see Section 2.1.12, beginning on page 158. }
72	5133506	1.0	Yes	TOP SIDE COVER R-1700 ASSY positioning GT {For further detail, see Section 2.1.13, beginning on page 160. }
86	5113088	1.0	Yes	COVER-PLATE-L positioning GT

Table 2-18 TABLE ASSY MID positioning GT (5134845 BOM, Rev 31) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
87	5113091	1.0	Yes	COVER-PLATE-R positioning GT
102	5130897	1.0	Yes	TOP SIDE BRKT RF-2 positioning GT
103	5129667	1.0	Yes	TOP SIDE BRKT RR MID positioning GT
104	5130899	1.0	Yes	TOP SIDE BRKT LF-2 positioning GT
105	5129668	1.0	Yes	TOP SIDE BRKT LR MID positioning GT
109	5113617	2.0	Yes	TOP COVER SIDE BRKT 11 positioning GT
119	N9612YX	16.0	Yes	SCREW HEXAGON SOCKET 6 MM 12 MM CAP SCREW 6X12 FE
121	5134666	1.0	Yes	TS MID SIDE-R ASSY positioning GT
122	5134665	1.0	Yes	TS MID SIDE-L ASSY
131	5125112	1.0	Yes	GT REAR TOUCH SWITCH CABLE Positioning GT
146	5219979	1.0	Yes	CRADLE COVER ASSY WITH LOGO LABEL Positioning GT of FRU
148	U0442BA	10.0	Yes	TUBE 3X0.25 BLACK
149	U0030AC	4.0	Yes	CLAMP OR CLIP CLAMP
150	U0031AC	12.0	Yes	CLAMP
151	U0032AC	2.0	Yes	CLAMP OLD CLAMP
154	5191755	1.0	Yes	TOUCH PLATE REAR ASSY
160	5130814	1.0	Yes	SIDE-PANEL-FRONT-LEFT-ASSY, Positioning, GT {For further detail, see Section 1.1.8, beginning on page 98. }
161	5130815	1.0	Yes	SIDE-PANEL-REAR-LEFT-ASSY, Positioning, GT {For further detail, see Section 1.1.9, beginning on page 100. }
162	5130816	1.0	Yes	SIDE-PANEL-FRONT-RIGHT-ASSY, Positioning, GT {For further detail, see Section 1.1.10, beginning on page 102. }
163	5130817	1.0	Yes	SIDE-PANEL-REAR-RIGHT-ASSY, Positioning, GT {For further detail, see Section 1.1.11, beginning on page 104. }
164	U0034BA	150.0	Yes	TAPE #5421 16Mug
165	U0047BA	80.0	Yes	TAPE
166	U0001AC	4.0	Yes	CLAMP NK-4N CLAMP
176	N9508JS	21.0	Yes	SCREW NABE 5X8 FE ZN
180	N9825YX	2.0	Yes	SCREW CAP SCREW 8X25 Fe

Table 2-18 TABLE ASSY MID positioning GT (5134845 BOM, Rev 31) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
181	N9601WS	4.0	Yes	WASHER PLAIN - LARGE 6.4 MM 12.5 MM WASHER 6.4X12.5X1.6FEZN
182	5126552	1.0	Yes	GND WIRE COVER Positioning GT
184	N9520ZX	8.0	Yes	SCREW HEXAGON SOCKET 5 MM 20 MM CAP SCREW 5X20 FE
185	N9501WS	32.0	Yes	WASHER PLAIN - LARGE 5.3 MM 10 MM WASHER M5
302	U0405BD	0.313	Yes	ADHESIVE YhjbN 271(af)
304	2221887	0.00625	Yes	U0427BD ADHESIVE
307	N9306ES	21.0	Yes	SCREW COUNTERSUNK HEA 3 MM 6 MM
308	N9306ZX	3.0	Yes	SCREW HEXAGON SOCKET 3 MM 6 MM CAP SCREW 3X6 FE
309	N9408JS	35.0	Yes	SCREW NABE 4X8 Fe Zn
310	N9414JS	2.0	Yes	SCREW PAN HEAD 4 MM 14 MM SCREW osX 4X14 FE ZN
311	N9510JS	8.0	Yes	SCREW NABE 5X10 Fe Zn
313	U0201AC	6.0	Yes	CLAMP CV_100
314	U0202AC	28.0	Yes	CLAMP CV_150
315	U0203AC	29.0	Yes	CLAMP CV_200
316	U0208AC	3.0	Yes	CLAMP TC-142 TAIRATSUPU KOTEI
317	U0212AC	23.0	Yes	CLAMP TMS_20
324	U0506AC	45.0	Yes	PROTECTOR EDGING CE-012
325	U0509AC	0.35	Yes	PROTECTOR CE_032
809	5142244	1.0	Yes	IMS-COVER-CTR-RIM positioning GT
810	5191722	1.0	Yes	RIM IMS COVER SIDE LH ASSY positioning GT {For further detail, see Section 1.1.12, beginning on page 106. }
811	5191736	1.0	Yes	RIM IMS COVER SIDE RH ASSY positioning GT {For further detail, see Section 1.1.13, beginning on page 108. }
2	5112083	1.0	No	CRADLE ATT BASE Positioning GT
4	5120588	1.0	No	cradle handle positioning GT
6	5129836	1.0	No	CRADLE-CARRIAGE-1700 positioning GT
7	5112420	1.0	No	CARRIAGE WIRE HOLDER Positioning GT
8	5112424	1.0	No	CARRIAGE SENSOR SUP Positioning GT

Table 2-18 TABLE ASSY MID positioning GT (5134845 BOM, Rev 31) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
9	5121378	1.0	No	CRADLE COVER BRKT positioning GT
12	5129840	1.0	No	CARRIAGE-CABLE-DUCT-SUP-MID positioning GT
51	N9514ZX	9.0	No	SCREW HEXAGON SOCKET 5 MM 14 MM
56	N9510ES	4.0	No	SCREW COUNTERSUNK HEA 5 MM 10 MM SCREW s[5X10 FE ZN
58	5127773	1.0	No	ELEVATION ASSY Positioning GT {For further detail, see Section 1.1.1, beginning on page 50. }
59	5117455	1.0	No	LINK BAR ASSEMBLY Positioning GT
60	5133502	1.0	No	TOP FRAME-1700 ASSY Positioning GT {For further detail, see Section 2.1.10, beginning on page 130. }
61	5126796-2	1.0	No	FOOT SW ASSY Positioning GT {For further detail, see Section 1.1.3, beginning on page 86. }
63	5154867-2	1.0	No	RIM BASE COVER LH ASSY positioning GT
64	5155337-2	1.0	No	RIM BASE COVER RH ASSY positioning GT
76	5113566	1.0	No	BASE BRKT 2R positioning GT
77	5113567	1.0	No	BASE BRKT 2L positioning GT
78	5113573	1.0	No	BRKT BASE FRAME L positioning GT
79	5113574	1.0	No	BRKT BASE FRAME R positioning GT
80	5159250	1.0	No	SIDE PANEL LWR BRKT L COMMON positioning GT
81	5159249	1.0	No	SIDE PANEL LWR BRKT R COMMON positioning GT
82	5112108-3	1.0	No	IMS Driving Cover Positioning GT
83	5112217	1.0	No	IMS D COVER FRONT Positioning GT
84	5112218	1.0	No	IMS D COVER MID Positioning GT
85	5118042	1.0	No	IMS D CV M BRKT 1 positioning GT
88	5132406	1.0	No	IMS D COVER MID REAR Postioning GT
90	5136785	2.0	No	IMS COVER BRKT 18 positioning GT
91	5113591	1.0	No	IMS COVER BRKT 13L positioning GT
92	5113592	1.0	No	IMS COVER BRKT 13R positioning GT
95	5113596	2.0	No	IMS COVER BRKT 3 positioning GT
98	5136786	1.0	No	IMS COVER BRKT 19L positioning GT

Table 2-18 TABLE ASSY MID positioning GT (5134845 BOM, Rev 31) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
99	5136787	1.0	No	IMS COVER BRKT 19R positioning GT
106	5113610	1.0	No	TOP FRONT BRKT 2L positioning GT
107	5113611	1.0	No	TOP FRONT BRKT 2R positioning GT
108	5113612	1.0	No	TOP FRONT BRKT 3 positioning GT
111	5119714	1.0	No	CABLE BEAR GUIDE positioning GT
114	N9408ES	4.0	No	SCREW COUNTERSUNK HEA 4 MM 8 MM
117	N9512LK	8.0	No	SCREW BIND UTIBATUKI 5X12 BsNi
129	5119697	2.0	No	NUT-PLATE LS positioning GT
136	5127957	12.0	No	CARRIAGE SHIM02 Positioning GT
138	5132405	1.0	No	IMS D COVER 1700 FRONT Positioning GT
139	5129853	1.0	No	TRAY-F-MID positioning GT
140	5129854	1.0	No	TRAY-M positioning GT
141	5127761	1.0	No	AMP-COVER positioning GT
145	U0200AK	0.215	No	FMI GASKET
152	5130902	3.0	No	WASHER positioning GT
153	5134650	2.0	No	STUD M8x36
155	5129674	1.0	No	REAR SW BRKT-2 positioning GT
156	5176626	1.0	No	REAR SW BRKT 104 positioning GT
157	5180589	1.0	No	REAR SW BRKT 1700 5L
158	5130900	1.0	No	REAR SW BRKT-4R positioning GT
159	5132685	1.0	No	REAR SW BRKT-101 positioning GT
171	2262896	10.0	No	U0465AA FLANGE SOCKET
175	N9430JS	2.0	No	SCREW PAN HEAD 4 MM 30 MM SCREW osX 4X30 FE ZN
177	N9801BS	4.0	No	NUT HEXAGON 8 MM FE ZN
178	N9850YX	6.0	No	SCREW HEXAGON SOCKET 8 MM 50 MM CAP SCREW 8X50 FE
179	U0241AA	10.0	No	NUT HEXAGON 4 MM NUT
183	2299437	36.0	No	U0472AA_FLANGE_SOCKET
186	5121904	3.0	No	CAM FOLLOWER BEARING U0443AH

Table 2-18 TABLE ASSY MID positioning GT (5134845 BOM, Rev 31) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
187	U0242AA	3.0	No	NUT HEXAGON 5 MM
189	N9308ES	1.0	No	SCREW 3X8 SARA FE ZN
190	N9406ES	5.0	No	SCREW SARA 4X6 Fe Zn
191	5125773	3.0	No	FLANGE NUT U0276AA
210	5139862	1.0	No	IMS-COVER-BRKT-24L positioning GT
211	5139863	1.0	No	IMS-COVER-BRKT-24R positioning GT
212	5151800	1.0	No	TOP REAR BRKT COMMON 6L Positioning GT
213	5151801	1.0	No	TOP REAR BRKT COMMON 6R Positioning GT
214	5151802	1.0	No	TOP REAR BRKT COMMON 7L Positioning GT
215	5151803	1.0	No	TOP REAR BRKT COMMON 7R Positioning GT
300	5160028	0.0025	No	LUBRI PLATE,U0304BD
301	U0303BD	6.25	No	MORICORT M-77 PASTE
303	U0409BD	0.313	No	ADHESIVE Yhjb #601
305	U0046BA	180.0	No	TAPE
306	5160029	0.00625	No	KURE CRC,U0305BD
500	5128983	1.0	No	ETL LABEL Positioning GT
502	5129157	1.0	No	PLACE LABEL Potisioning GT
503	5129016	2.0	No	SHIPPING LABEL Positioning GT
506	U0036BN	12.0	No	NAME PLATE
507	5139485	1.0	No	WEEE LABEL WITH BIN ONLY
801	5142259	1.0	No	IMS COVER BRKT 31L positioning GT
802	5142261	1.0	No	IMS COVER BRKT 31R positioning GT
803	5142269	1.0	No	IMS COVER BRKT 33 positioning GT
804	5175971	2.0	No	PLATE NUT Positioning GT
805	5175384	1.0	No	IMS COVER BRKT 33R Positioning GT
806	5175137	1.0	No	IMS COVER BRKT 33L Positioning GT
807	5120580-2	1.0	No	TS_IMS_REAR Positioning GT
808	5189420	1.0	No	IMS DUMMY PLUG LABEL positioning GT

Table 2-18 TABLE ASSY MID positioning GT (5134845 BOM, Rev 31) (Continued)

2.1.1 Elavation Assembly (5127773)

FRU = NO For further detail, see [Section 1.1.1, beginning on page 50.](#)

2.1.2 Foot Switch Assembly (5126796-2)

FRU = NO For further detail, see [Section 1.1.3, beginning on page 86.](#)

2.1.3 RIM Top Cover Assembly Positioning GT (5127465-3)

FRU = YES) For further detail, see [Section 1.1.4, beginning on page 90.](#)

2.1.4 Side-Panel-Front-Left-Assembly Positioning GT (5130814)

FRU = YES For further detail, see [Section 1.1.8, beginning on page 98.](#)

2.1.5 Side-Panel-Rear-Left-Assembly Positioning GT (5130815)

FRU = YES For further detail, see [Section 1.1.9, beginning on page 100.](#)

2.1.6 Side-Panel-Front-Right-Assembly Positioning GT (5130816)

FRU = YES For further detail, see [Section 1.1.10, beginning on page 102.](#)

2.1.7 Side-Panel-Rear-Right-Assembly Positioning GT (5130817)

FRU = YES For further detail, see [Section 1.1.11, beginning on page 104.](#)

2.1.8 RIM IMS Cover Side LH Assembly Positioning GT (5191722)

FRU = YES For further detail, see [Section 1.1.12, beginning on page 106.](#)

2.1.9 RIM IMS Cover Side RH Assembly Positioning GT (5191736)

FRU = YES For further detail, see [Section 1.1.13, beginning on page 108.](#)

2.1.10 TOP FRAME-1700 ASSY Positioning GT (5133502)

FRU = NO

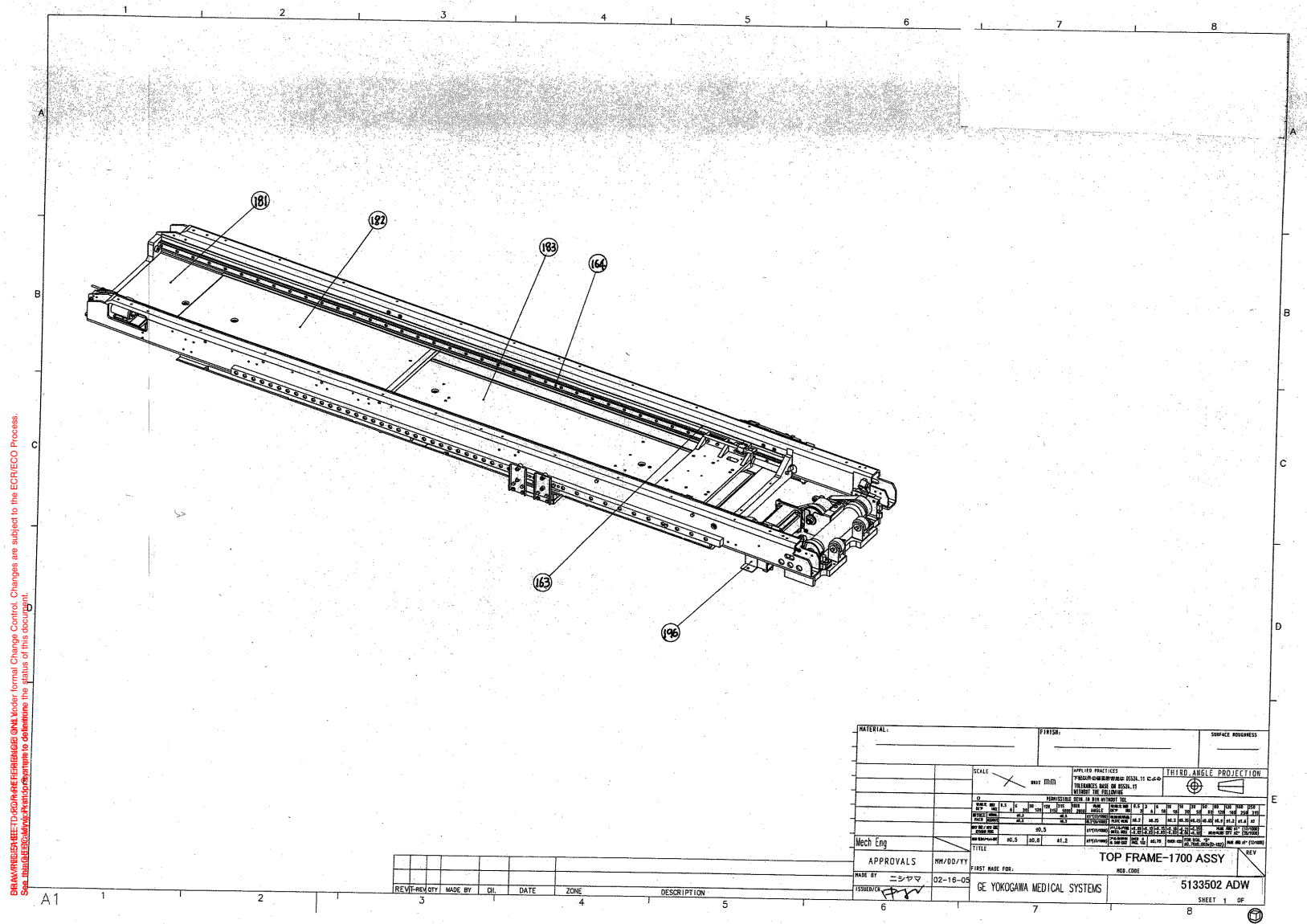


Figure 2-59 Top Frame 1700 Assembly (5133502ADW, sheet 1, rev 5)


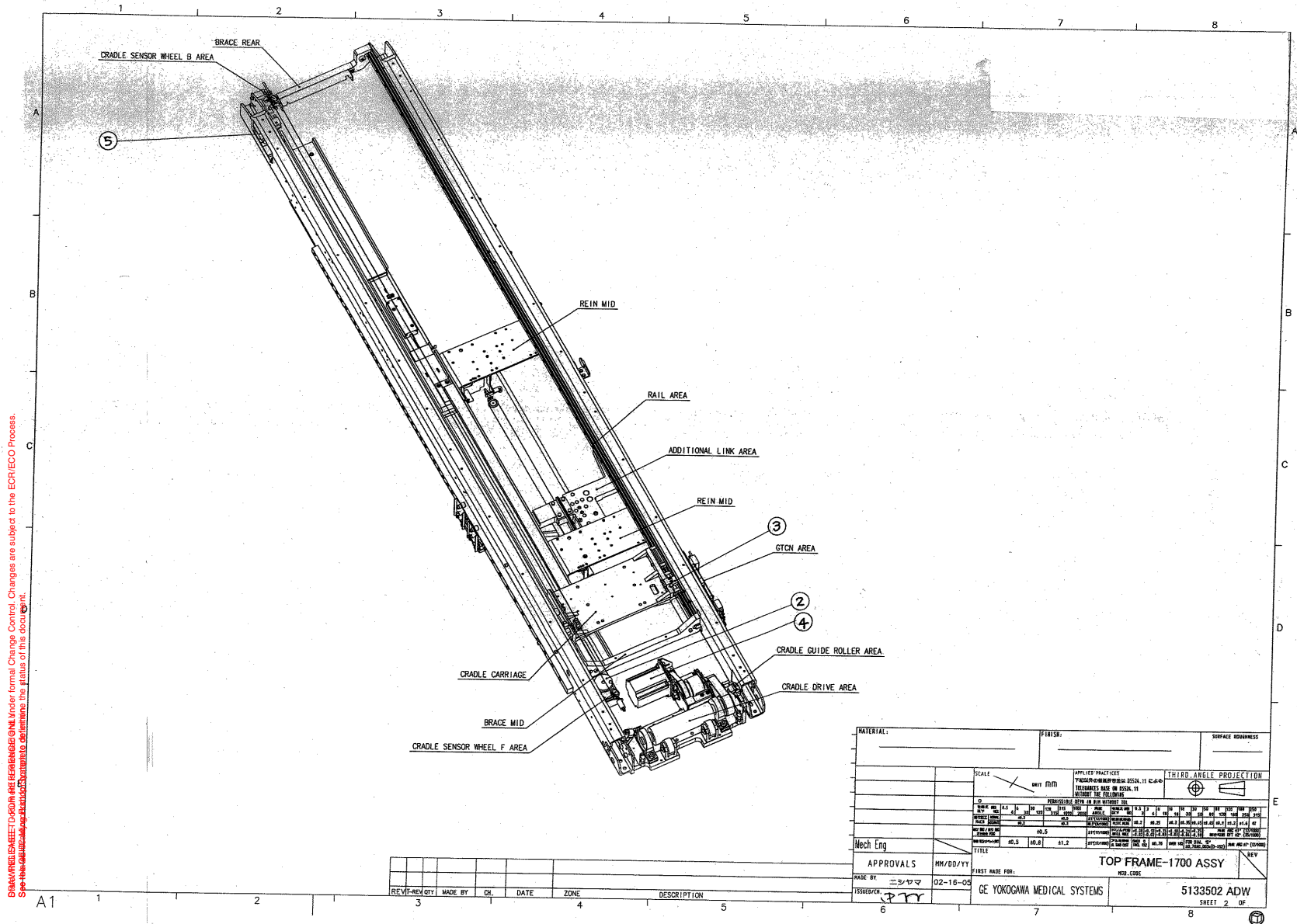


Figure 2-60 Top Frame 1700 Assembly (5133502ADW, sheet 2, rev 5)



Section 2.0 — GT 1700 Table (5122080)



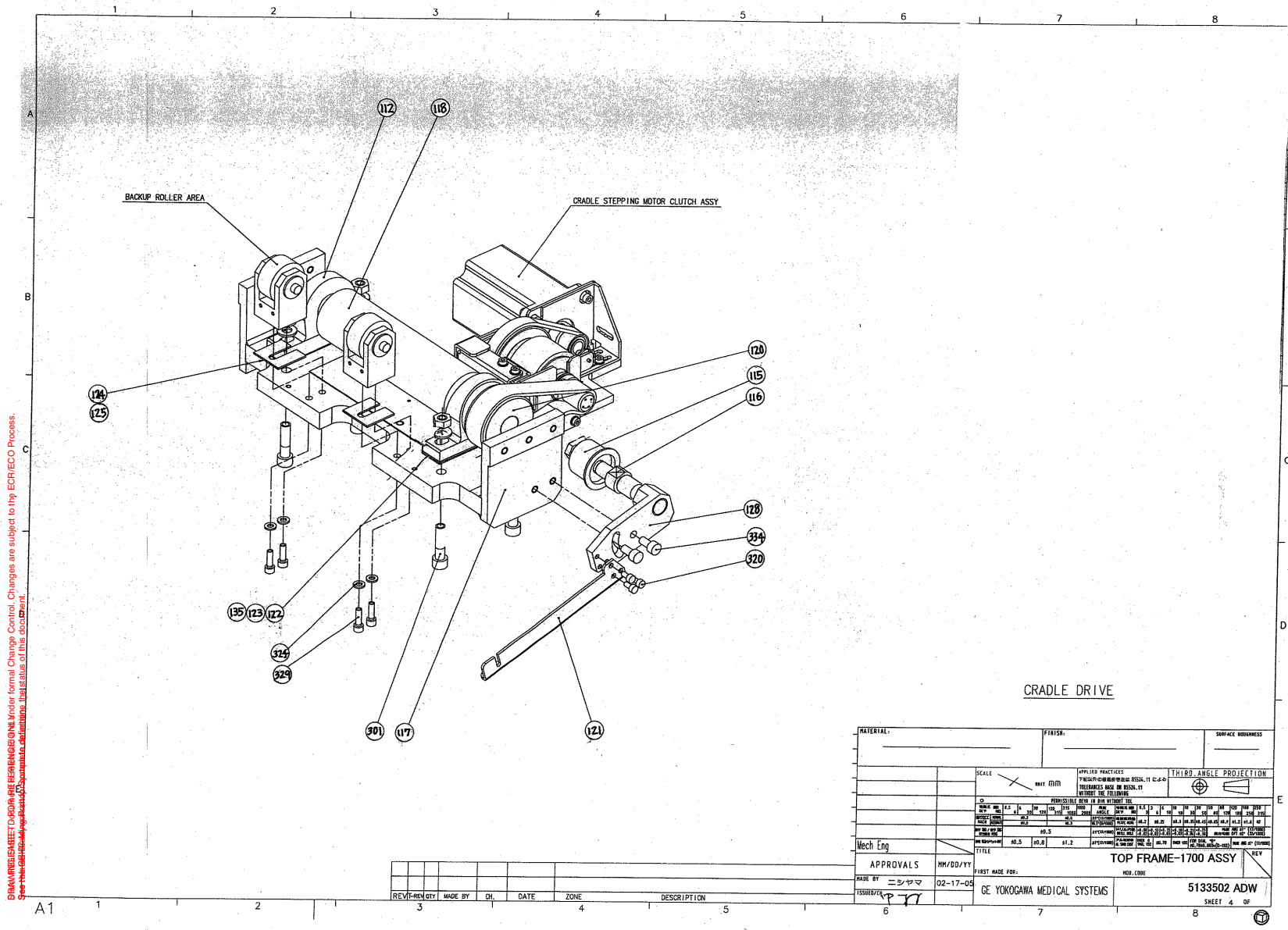

Section 2.0 — GT 1700 Table (5122080)

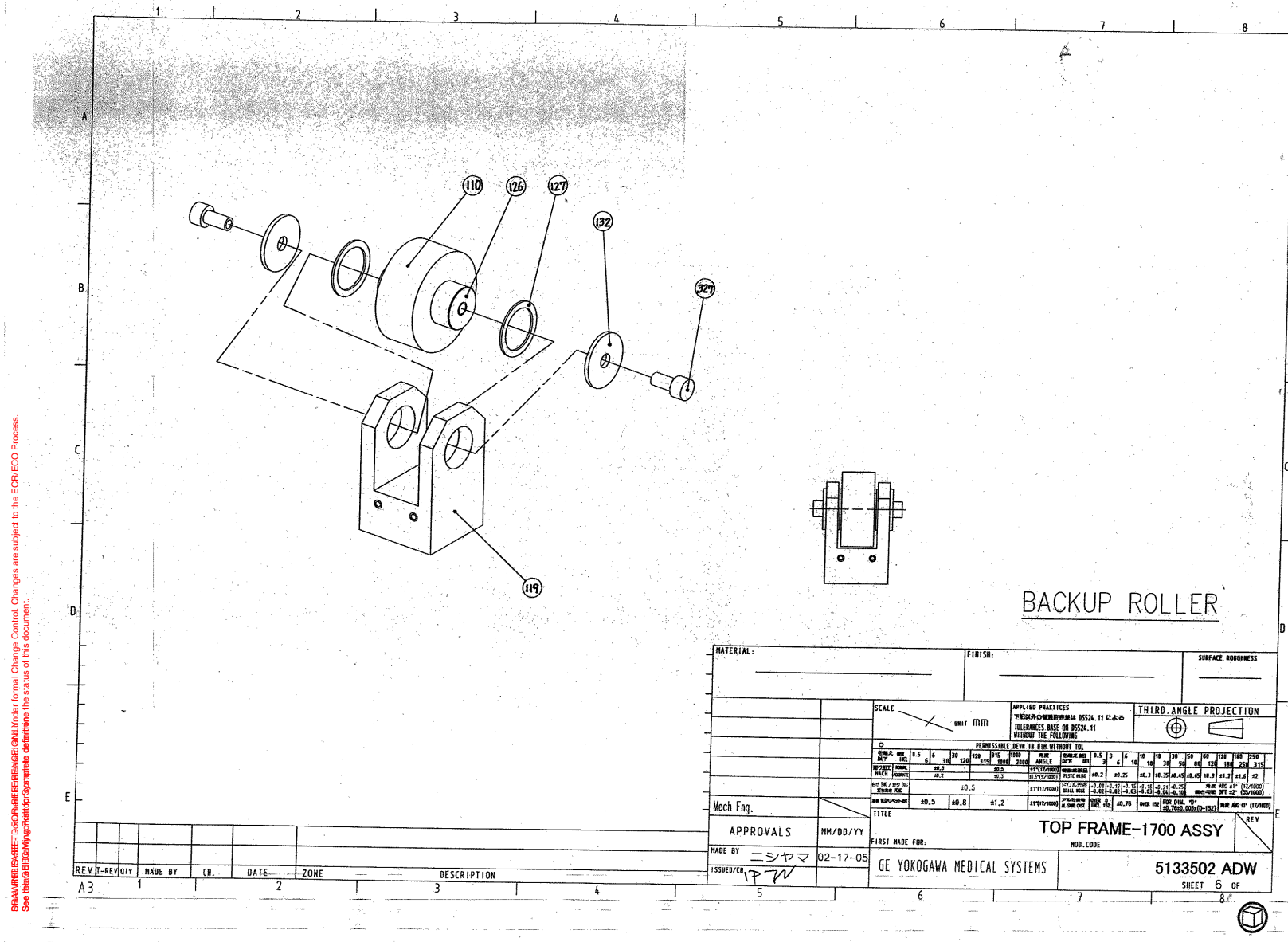
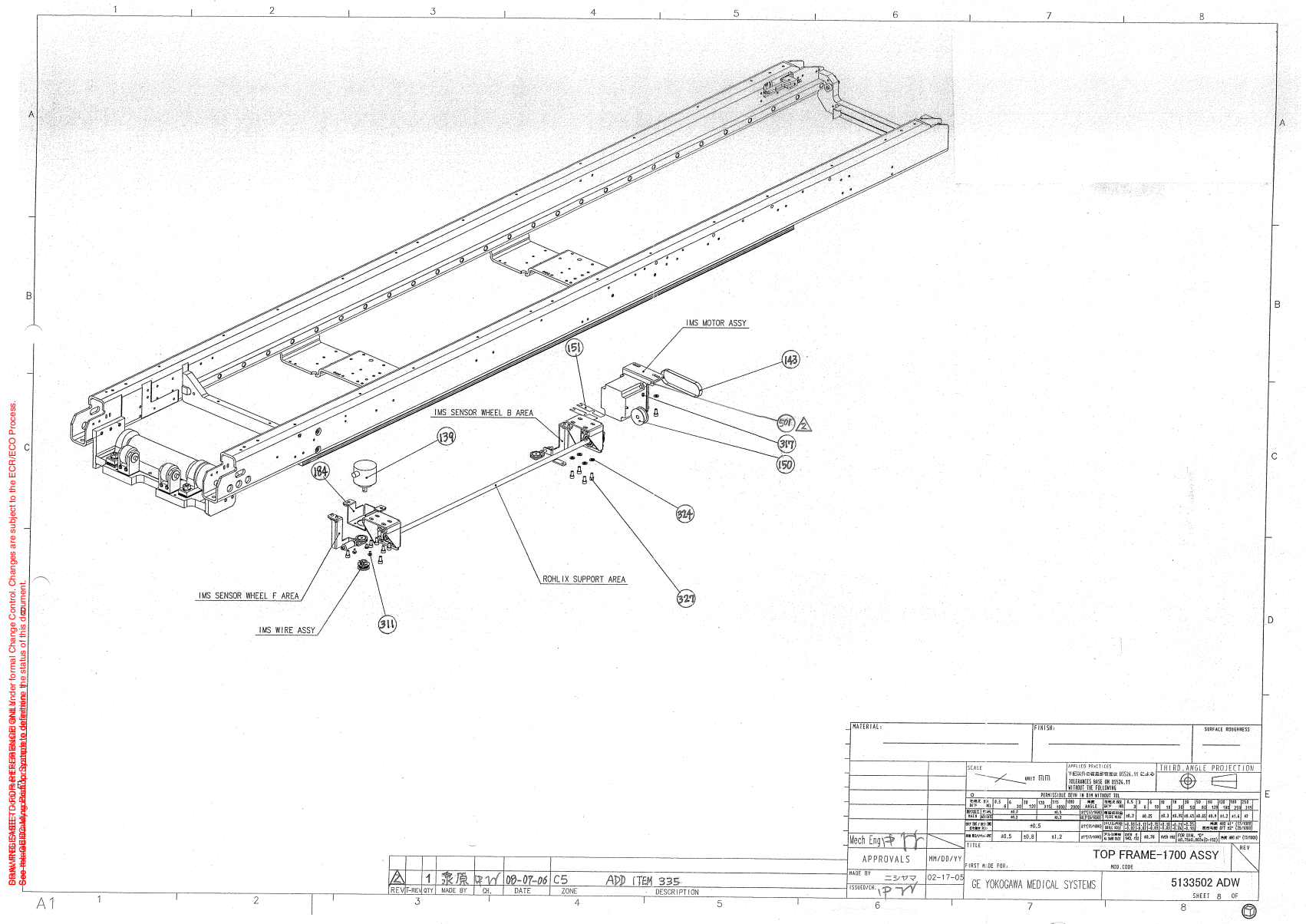


Figure 2-64 Top Frame 1700 Assembly (5133502ADW, sheet 6, rev 5)





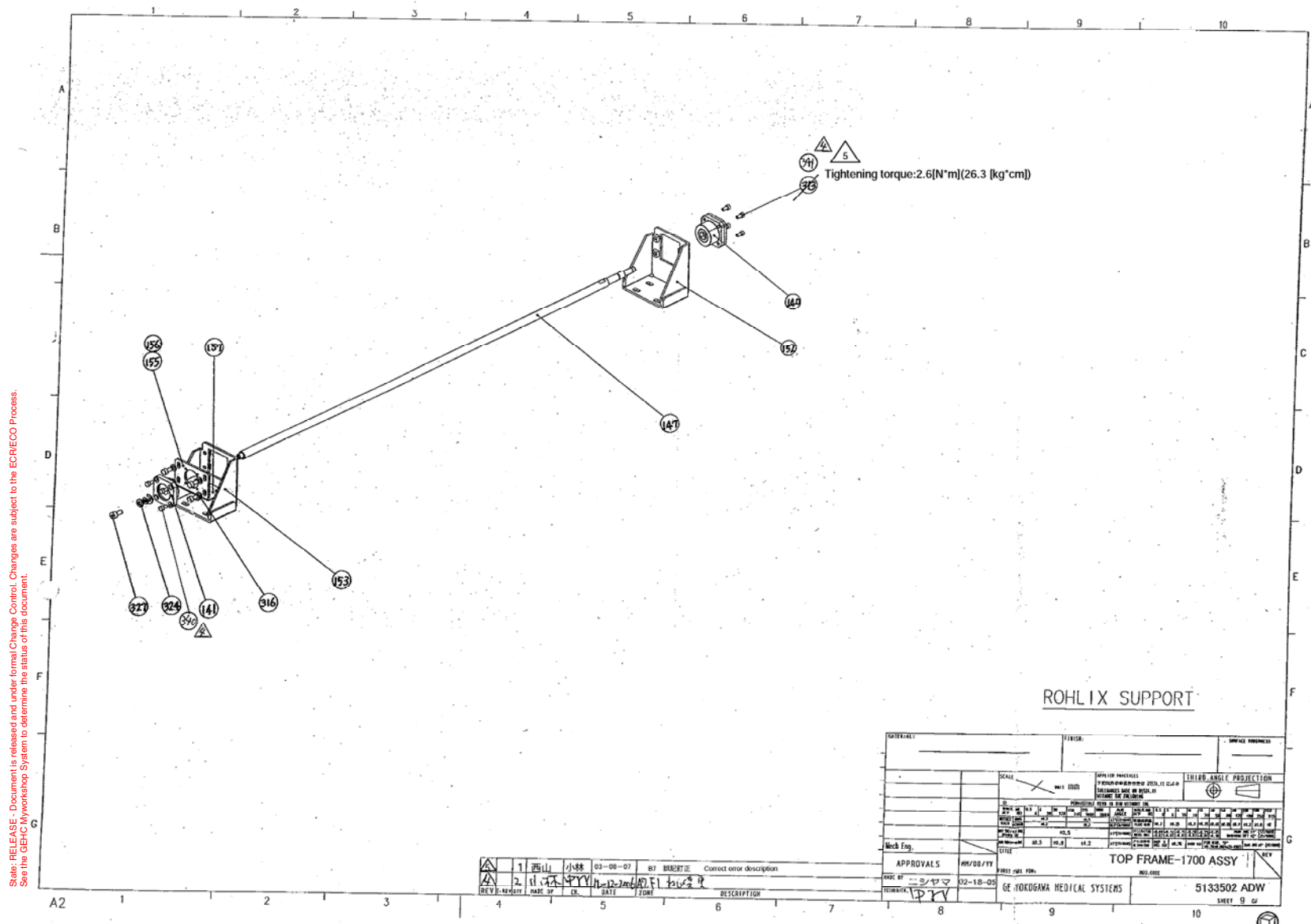


Figure 2-67 Top Frame 1700 Assembly (5133502ADW, sheet 9, rev 5)

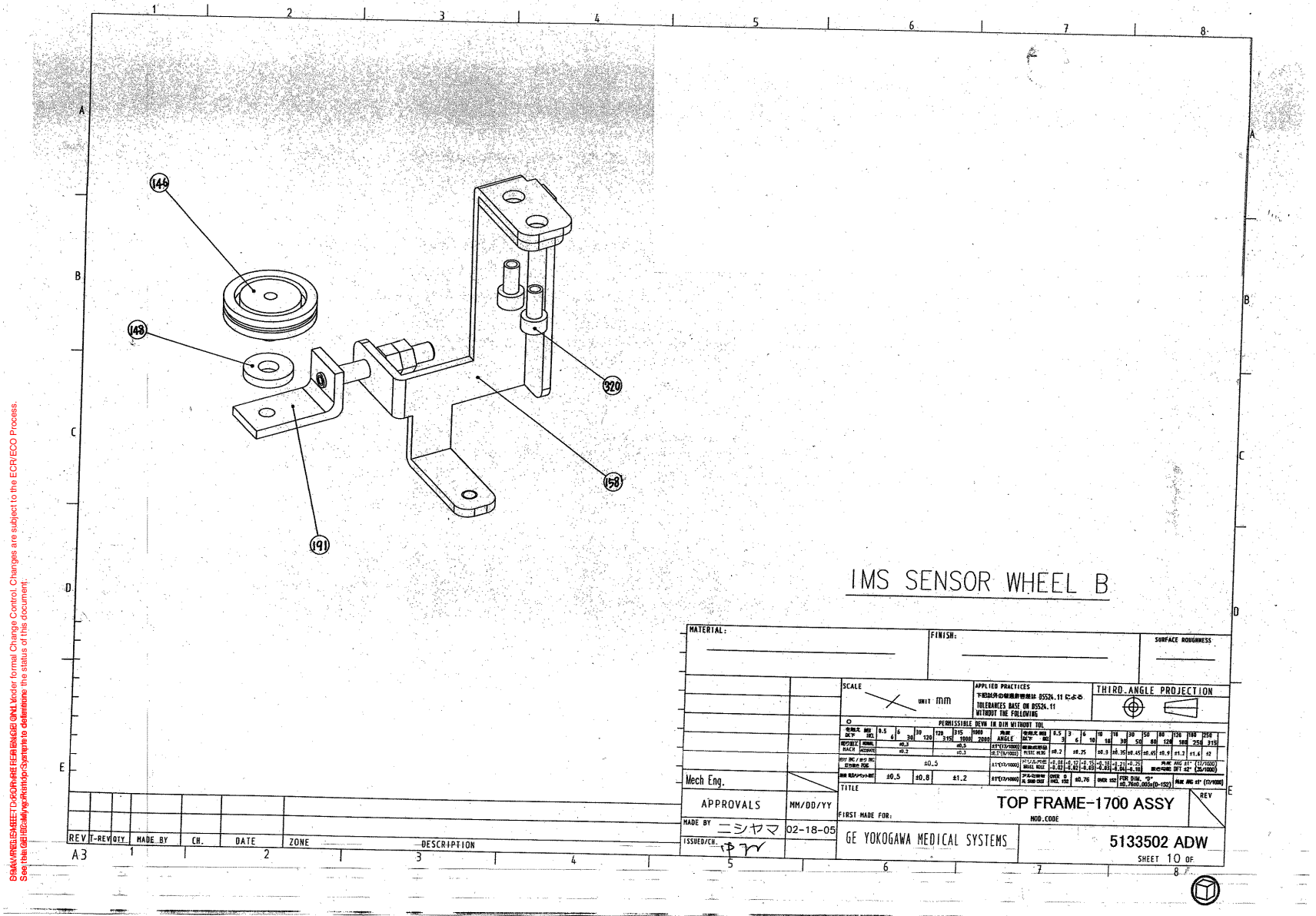


Figure 2-68 Top Frame 1700 Assembly (5133502ADW, sheet 10, rev 5)



Section 2.0 — GT 1700 Table (5122080)



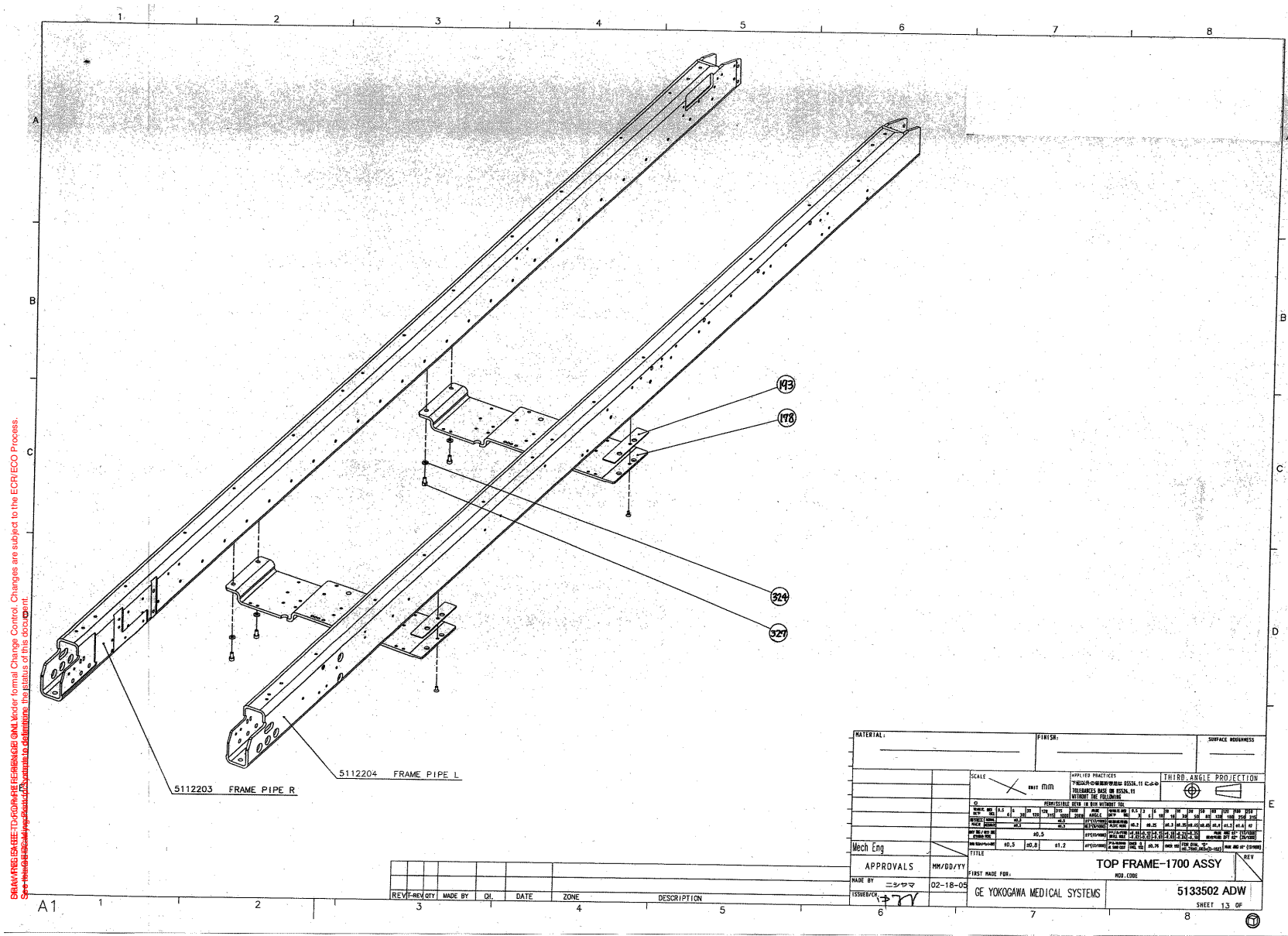


Figure 2-71 Top Frame 1700 Assembly (5133502ADW, sheet 13, rev 5)





Section 2.0 — GT 1700 Table (5122080)



2 - GT Table
GT 1700 Table (5122080)

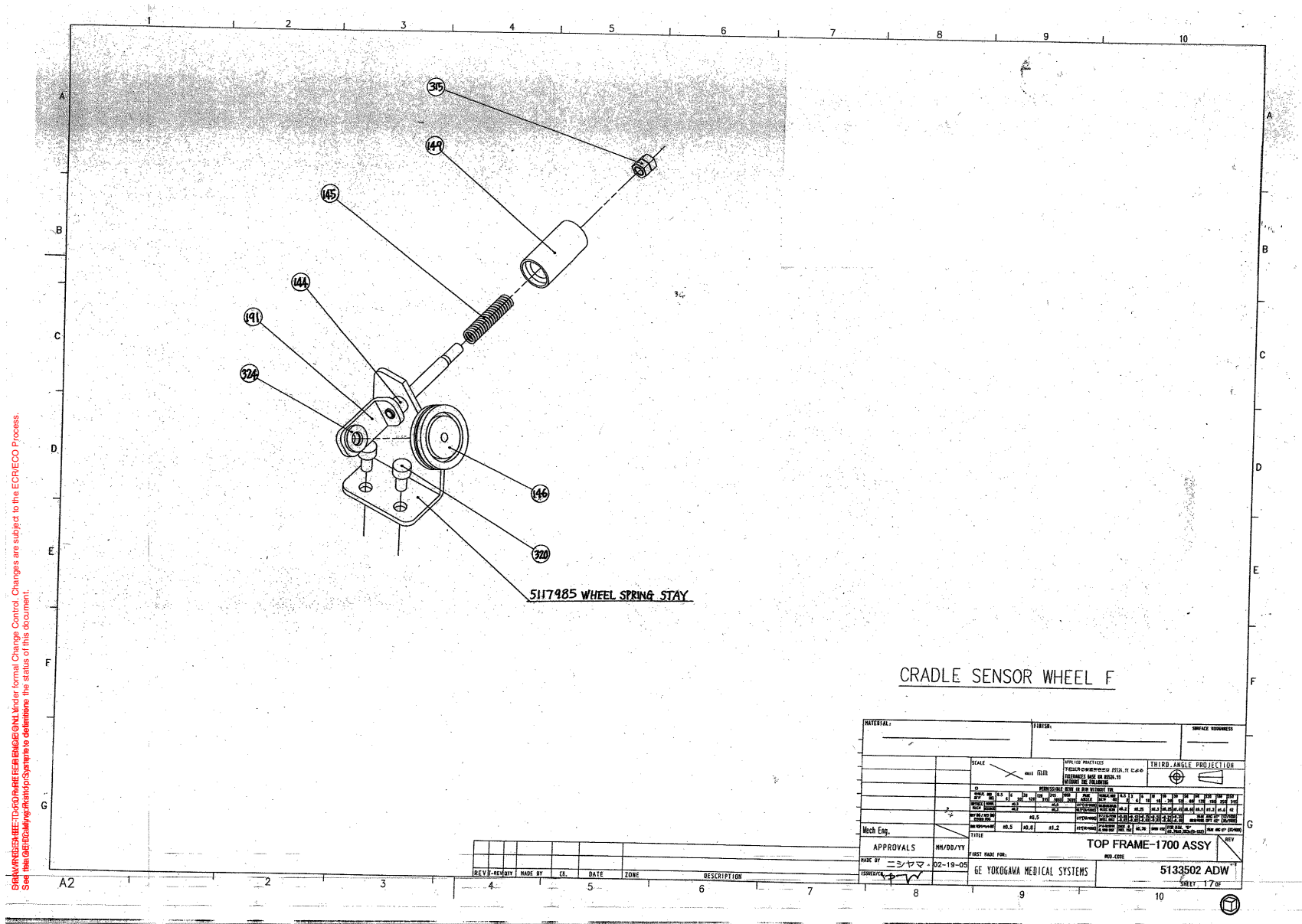


Figure 2-75 Top Frame 1700 Assembly (5133502ADW, sheet 17, rev 5)

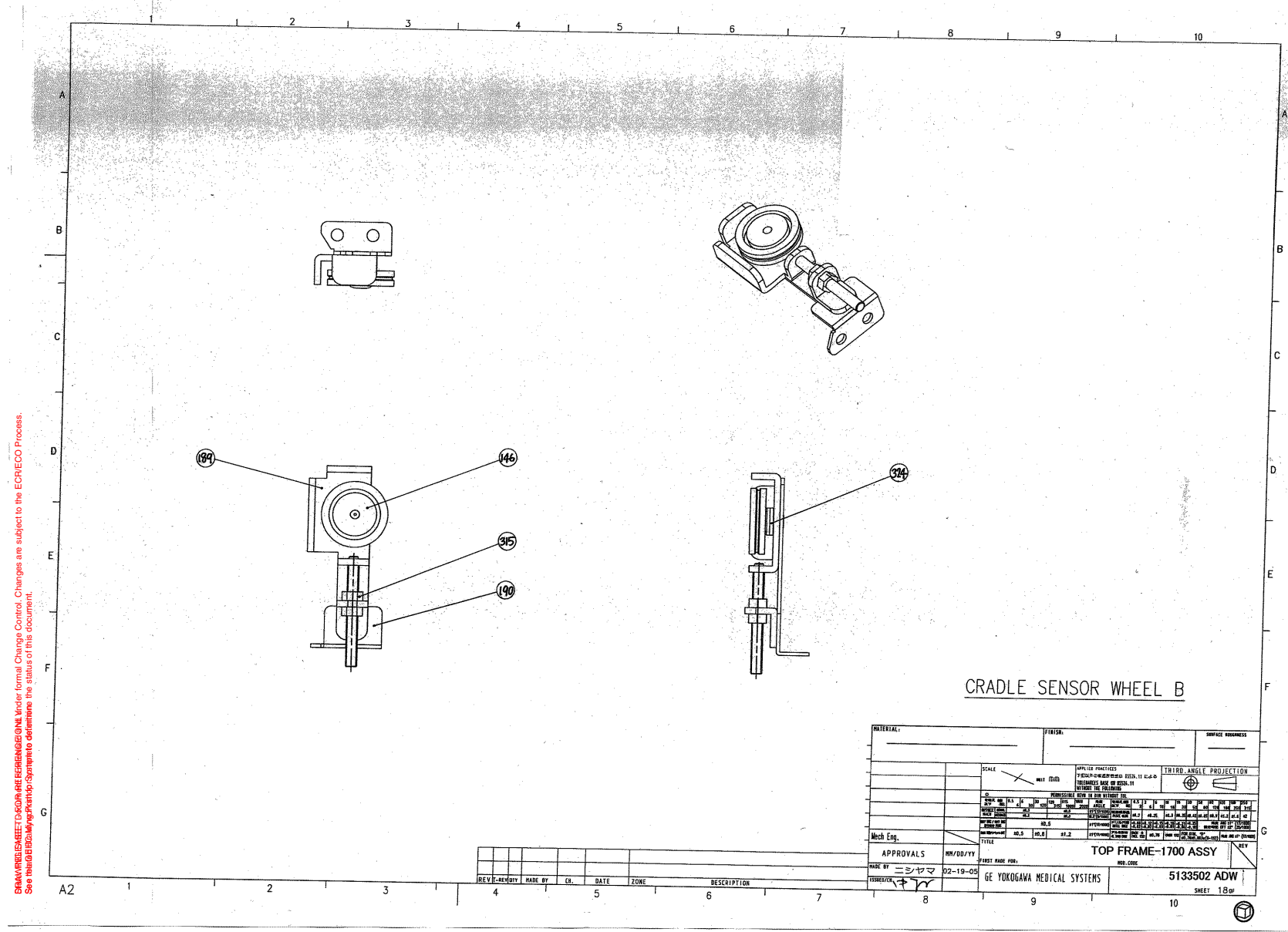


Figure 2-76 Top Frame 1700 Assembly (5133502ADW, sheet 18, rev 5)

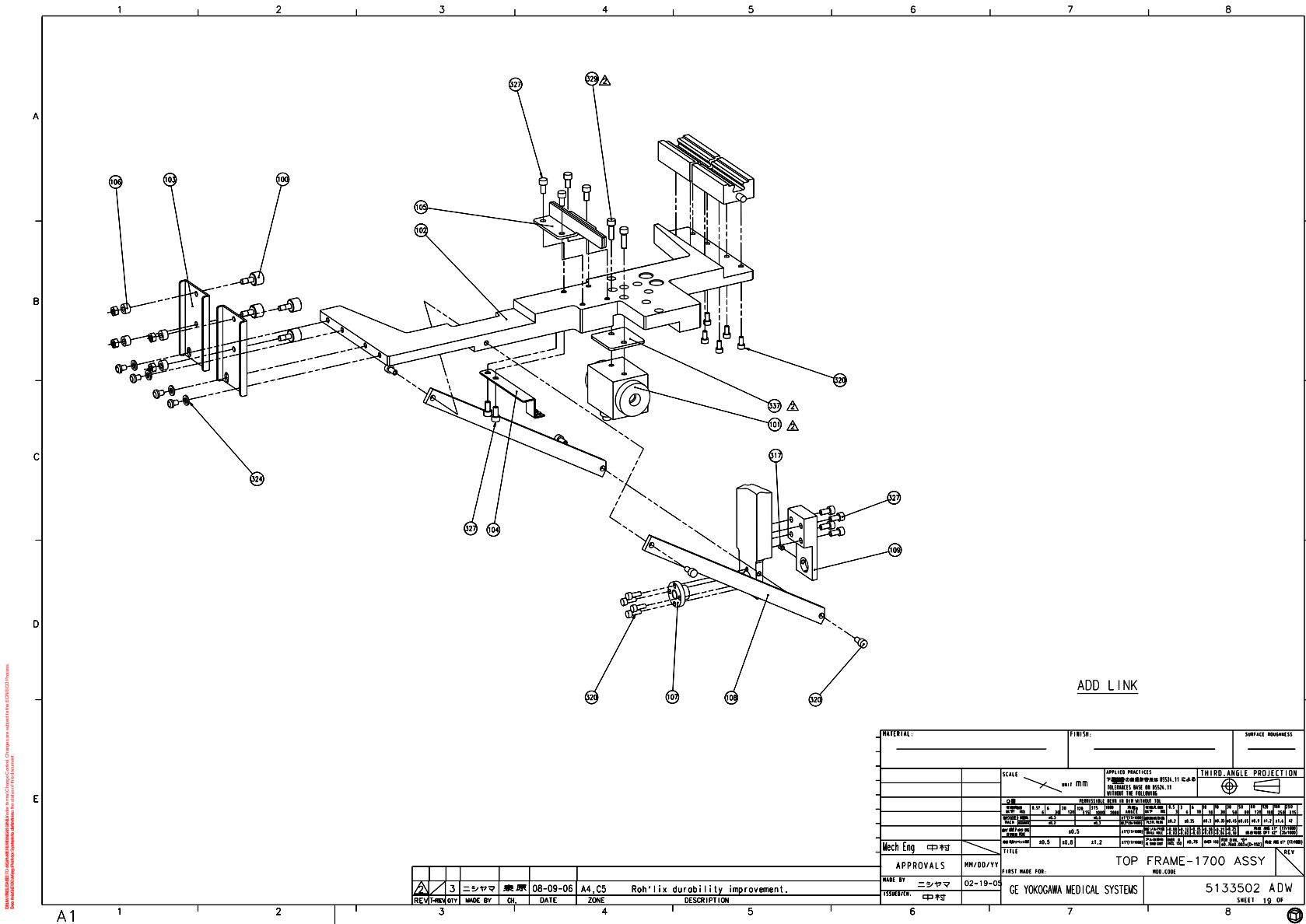


Figure 2-77 Top Frame 1700 Assembly (5133502ADW, sheet 19, rev 5)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2201853	2.0	Yes	U0226AH BEARING BALL,DOUBLE RO 20MM
23	P9145CZ	4.0	Yes	WASHER ROLLER
27	5122208-2	1.0	Yes	ABSOLUTE ENCODER FOR IMS GT
31	5123252	1.0	Yes	BELT 318P3M10
32	U0908AH	2.0	Yes	BUSHING 80F-0608
34	U0009AH	4.0	Yes	CAM FOLLOWER WITH SHAFT 6 MM 30 MM DV-30-B0.5
36	P9145DY	4.0	Yes	SPACER
37	P9145DL	2.0	Yes	CASE SPRING
48	U0443BA	5.0	Yes	TUBE .vk X F4X0.25 hY
52	5122207-2	1.0	Yes	ABSOLUTE ENCODER FOR CRADLE GT
53	5122210	1.0	Yes	MAGNET SCALE GT
56	U0125AE	3.0	Yes	CUSHION CK-17
89	U0402AH	4.0	Yes	CAM FOLLOWER WITH SHAFT 6 MM 16 MM CF 6 VBUU
99	5127495	1.0	Yes	IMS WIRE ASSY Positioing GT
100	5138091	1.0	Yes	CRADLE WIRE ASSY 1700 Positioning GT
101	5133508	1.0	Yes	LINEAR SCALE-1700 ASSY positioning GT
102	5135265	1.0	Yes	CRADLE STEEPING MOTOR CLUTCH 1700 ASSY Positioning GT {For further detail, see Section 2.1.10.1, beginning on page 154. }
103	5133504	1.0	Yes	CRADLE LOCK-1700 ASSY positioning GT
104	5127247-2	1.0	Yes	IMS Motor Assy Positioning GT
105	N9020YX	5.0	Yes	SCREW HEXAGON SOCKET 10 MM 20 MM CAP SCREW 10X20 FE
106	N9040YX	8.0	Yes	SCREW HEXAGON SOCKET 10 MM 40 MM CAP SCREW M10X40 FE
111	N9306ZX	4.0	Yes	SCREW HEXAGON SOCKET 3 MM 6 MM CAP SCREW 3X6 FE
118	N9408ZX	5.0	Yes	SCREW HEXAGON SOCKET 4 MM 8 MM CAP SCREW 4X8 FE
121	N9501BS	9.0	Yes	NUT HEXAGON 5 MM NUT M5 FE ZN
122	N9501WS	2.0	Yes	WASHER PLAIN - LARGE 5.3 MM 10 MM WASHER M5
123	N9508SJ	2.0	Yes	SCREW SET SCREW 5 MM 8 MM SET SCREW M5X0.8X8htXvsg

Table 2-19 TOP FRAME-1700 ASSY Positioning GT (5133502 BOM, Rev 16)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
125	N9510LK	14.0	Yes	SCREW PAN HEAD 5 MM 10 MM
126	N9510ZX	9.0	Yes	SCREW HEXAGON SOCKET 5 MM 10 MM CAP SCREW 5X10 FE
128	N9601WS	55.0	Yes	WASHER PLAIN - LARGE 6.4 MM 12.5 MM WASHER 6.4X12.5X1.6FEZN
132	N9612YX	39.0	Yes	SCREW HEXAGON SOCKET 6 MM 12 MM CAP SCREW 6X12 FE
135	N9620YX	41.0	Yes	SCREW HEXAGON SOCKET 6 MM 20 MM CAP SCREW 6X20 FE
138	N9630YX	2.0	Yes	SCREW HEXAGON SOCKET 6 MM 30 MM CAP SCREW 6X30 FE
139	N9801WS	4.0	Yes	WASHER PLAIN - LARGE 8.4 MM 17 MM
142	N9825YX	2.0	Yes	SCREW CAP SCREW 8X25 Fe
162	N9520ZX	69.0	Yes	SCREW HEXAGON SOCKET 5 MM 20 MM CAP SCREW 5X20 FE
170	U0009AP	2.0	Yes	CRIMP STYLE TERMINAL B-4(CU2.5X1.5X3.1)
173	5184277	1.0	Yes	Rotary Dumper U0122AQ
181	N9412ZX	4.0	Yes	SCREW HEXAGON SOCKET 4 MM 12 MM SCREW
184	N9625YX	4.0	Yes	SCREW CAP 6X25
210	46-208718P3	1.0	Yes	CLAMP, CUSHIONED, ID .375 LOOP TYPE TA MFG CORP TA713D6
2	5121923	1.0	No	Parallel key U0163AP
3	5121921	2.0	No	BEARING UNIT U0447AH
4	5125774	8.0	No	FLANGE SOCKET U0478AA
5	5125773	4.0	No	FLANGE NUT U0276AA
6	5125253	1.0	No	IDLER NIPS3050 Positioning GT
7	5125254	1.0	No	IDLER PIN PID1522 Positioning GT
8	5112206	1.0	No	CRADLE-DRIVE-BASE Positioning GT
9	5113224	1.0	No	ROLLER positioning GT
10	5113219	2.0	No	BRKT ROLLER positioning GT
11	5113220	1.0	No	ROLLER PULLY positioning GT
12	5118619	1.0	No	TENSION BAR
13	5118587	10.0	No	SPACER 01 PILLOW
14	5118588	6.0	No	SPACER 05 PILLOW
15	5118592	10.0	No	SPACER 01 BASE

Table 2-19 TOP FRAME-1700 ASSY Positioning GT (5133502 BOM, Rev 16) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
16	5118591	6.0	No	SPACER 05 BASE
17	2199990	2.0	No	P9235JB PIN
18	2199991	4.0	No	P9235LA SPACER INNER
19	5117121	1.0	No	IDLER BASE positioning GT
20	5118181	2.0	No	CRADLE GUIDE ROLLER Positioning GT
21	5118183	2.0	No	CRADLE GUIDE ROLLER SPRT Positioning GT
22	5118182	2.0	No	SPACER BRASS Positioning GT
28	5126434	1.0	No	BEARING UNIT Positioning GT
29	2201842	1.0	No	U0428AH BEARING ATTACHMENT BEARING UNIT
30	5125163	1.0	No	IMS WHEEL SCREW BRKT Positioning GT
33	U0004AQ	2.0	No	SPRING COMPRESSION 47 MM SC-085-120-0470 4017
35	5115688	1.0	No	ROHLIX-SHAFT positioning GT
39	5115675	1.0	No	40P3M10BF-ROHLIX positioning GT
41	5115679-2	1.0	No	IMS SUP SHAFT FRONT
42	5115678-2	1.0	No	IMS SUP SHAFT BACK
43	5115682	1.0	No	IMS-SENSOR-BRKT positioning GT
44	5115683	1.0	No	IMS-MOTOR-PLATE-FF10 positioning GT
45	5115684	1.0	No	IMS-MOTOR-PLATE positioning GT
46	5115687	1.0	No	IMS-DRIVE-WASHER positioning GT
47	5117312	1.0	No	IMS SENSOR BRKT REAR Positioning GT
54	5133087	1.0	No	CRADLE LINEAR GUIDE positioning GT
55	5133088	1.0	No	IMS LINEAR GUIDE positioning GT
57	5129820	1.0	No	PIPE-R-MID positioning GT
58	5129818	1.0	No	PIPE-L-MID positioning GT
59	5128682	1.0	No	FRAME RAIL MID positioning GT
60	5128683	1.0	No	PLATE RAIL MID Positioning GT
61	5129848	1.0	No	CRADLE-ENCODER-BRKT-MID positioning GT
63	5129677	1.0	No	BRACE REAR 1700 positioning GT

Table 2-19 TOP FRAME-1700 ASSY Positioning GT (5133502 BOM, Rev 16) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
64	5129838	1.0	No	IMS-RAIL-PLATE-MID-FRONT positioning GT
65	5116870	1.0	No	IMS RAIL PLATE MID
66	5129837	1.0	No	IMS-RAIL-PLATE-MID-REAR positioning GT
67	5113230	1.0	No	GTCN BRKT positioning GT
68	5113241	2.0	No	REIN MID positioning GT
69	5129852	1.0	No	SUP-TOP-PLATE-REAR-1700 positioning GT
70	5113242	1.0	No	STOPPER positioning GT
71	5129846	1.0	No	TOP-PLATE-REAR-1700 positioning GT
73	5129845	1.0	No	TOP-PLATE-FRONT-MID positioning GT
74	5116958	1.0	No	ENCODER SUP IMS Positioning GT
75	5129841	1.0	No	CABLE-DUCT-FRONT-MID positioning GT
76	5129842	1.0	No	CABLE-DUCT-REAR-MID positioning GT
77	5117355	1.0	No	CABLE-GUIDE-BRKT
79	5117982	1.0	No	WHEEL-SCREW-BRKT-REAR Positioning GT
80	5117983	1.0	No	REAR-WHEEL-ADJUST-BRKT Positioning GT
81	5117984	2.0	No	WHEEL-SCREW-BRKT-FRONT Positioning GT
83	5117163	2.0	No	REIN MID SPACER
84	5134055	1.0	No	NUT PLATE 1700 positioning GT
88	5126275	1.0	No	TOP-FRONT-BRKT-3L positioning GT
90	5182656	1.0	No	Linear Actuator
91	5112925	1.0	No	IMS BASE GT positioning GT
92	5112921	2.0	No	UPPER SUP ADD positioning GT
93	5112926	1.0	No	IMS CABLE HOLDER 1 positioning GT
94	5112927	1.0	No	IMS STOPPER positioning GT
95	5116960	4.0	No	CUMFOLLOWER SPACER
96	5116947	1.0	No	ADDLINK-NUT Positioning GT
97	5124243	2.0	No	IMS BASE SUP Positioning GT
98	5121568	1.0	No	LINK BAR SUPPORT Positioning GT

Table 2-19 TOP FRAME-1700 ASSY Positioning GT (5133502 BOM, Rev 16) (Continued)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
130	N9610SJ	1.0	No	SET SCREW
141	N9820YX	11.0	No	SCREW HEXAGON SOCKET 8 MM 20 MM CAP SCREW 8X20 FE
147	5118590	4.0	No	SPACER 40 PILLOW
149	5112471-3	1.0	No	GTCN BOARD ASSY Positioning GT
160	N9121CS	1.0	No	NUT M12 3SYU FeZn
172	5183123	1.0	No	SPACER ROHLIX ACTUATOR
174	5177564	1.0	No	IMS STOPPER 1700 positioning GT
175	N9616YX	2.0	No	SCREW HEXAGON SOCKET 6 MM 16 MM CAP SCREW 6X16 FE
180	N9310ZX	4.0	No	SCREW HEXAGON SOCKET 3 MM 10 MM CAP SCREW 3X10 FE
182	2299437	9.0	No	U0472AA_FLANGE_SOCKET
183	N9514ZX	3.0	No	SCREW HEXAGON SOCKET 5 MM 14 MM
200	5135604	1.0	No	SHAFT LABEL Positioning GT
201	5129843	1.0	No	SUP-CABLE-DUCT-FRONT-MID positioning GT
202	5129839	1.0	No	STOPPER-CARRIAGE-FRONT positioning GT
211	2379698	0.001	No	U0064BA_EDGING_TAPE

Table 2-19 TOP FRAME-1700 ASSY Positioning GT (5133502 BOM, Rev 16) (Continued)

2.1.10.1 CRADLE STEEPING MOTOR CLUTCH ASSY (5135265)

FRU = YES

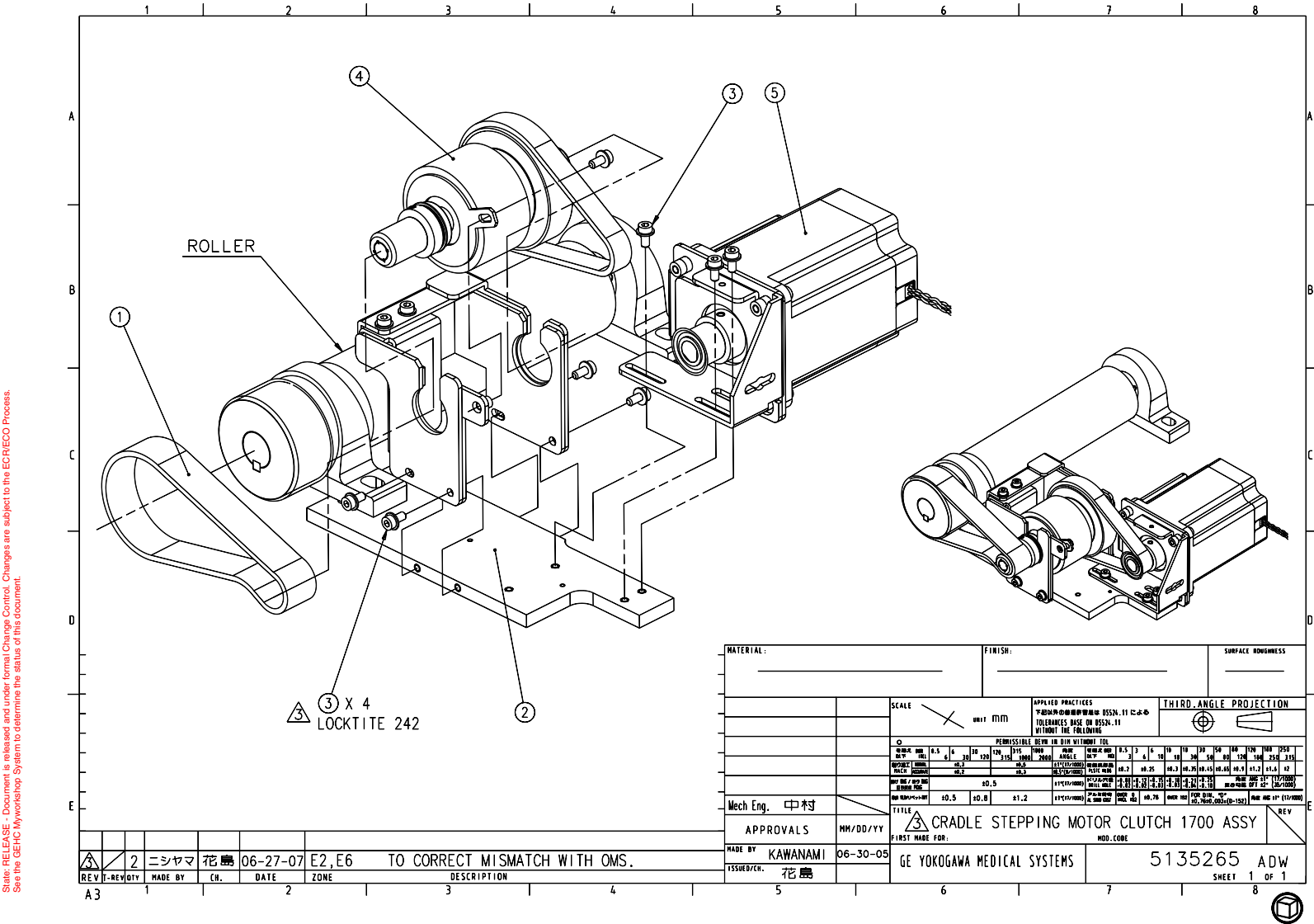


Figure 2-78 CRADLE STEEPING MOTOR CLUTCH ASSY (5135265ADW, sheet 1, rev 3)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	5123238	1.0	Yes	BELT 375 UP5M25 Positioning GT
2	5118348	1.0	No	PULLEY SUP BASE
3	2299437	7.0	No	U0472AA_FLANGE_SOCKET
4	5127775	1.0	No	CLUTCH ASSY
5	5135249	1.0	No	CRADLE MOTOR 1700 ASSY

Table 2-20 CRADLE STEEPING MOTOR CLUTCH 1700 ASSY Positioning GT (5135265 BOM, Rev 2)

2.1.11 TOP COVER REAR-1700 ASSY positioning GT (5133505)

FRU = YES

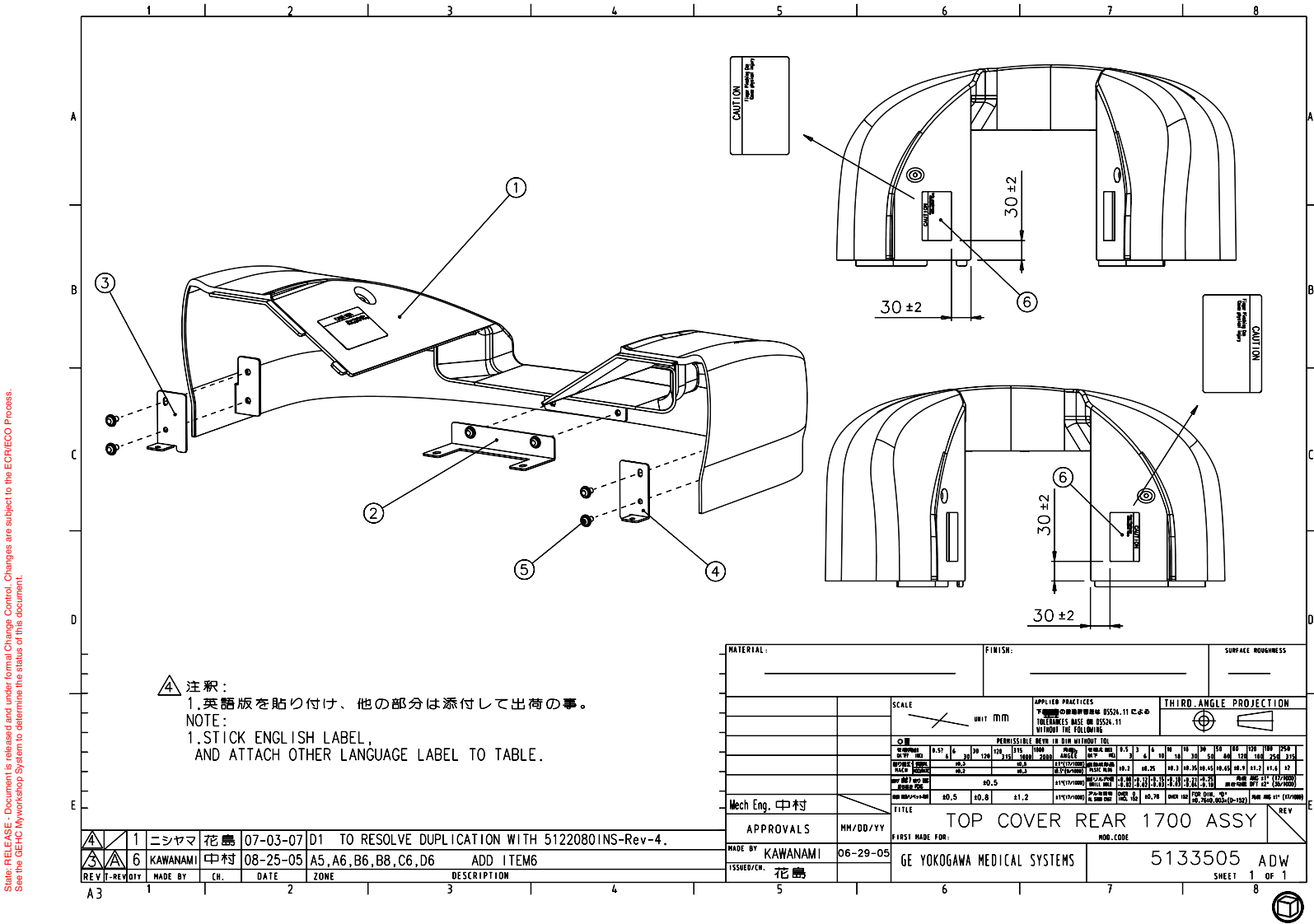


Figure 2-79 TOP COVER REAR-1700 ASSY positioning GT (5133505ADW, sheet 1, rev 4)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
5	N9510HR	6.0	Yes	SCREW TP
6	5129092	1.0	Yes	FINGER PINCHING LABEL Positioning GT
1	5138524	1.0	No	TOP COVER REAR 1700
2	5129664	1.0	No	TOP REAR BRKT-MID-3
3	5132688	1.0	No	TOP REAR BRKT MID-101R positioning GT
4	5132690	1.0	No	TOP REAR BRKT MID-101L positioning GT

Table 2-21 TOP COVER REAR-1700 ASSY positioning GT (5133505 BOM, Rev 4)

FRU = YES



Section 2.0 — GT 1700 Table (5122080)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
5	N9510HR	14.0	Yes	SCREW TP
6	2199652	1.0	Yes	P9235LZ LABEL
7	5136012	1.0	Yes	IMS FINGER PINCHING LABEL positioning GT
1	5129860-2	1.0	No	TOP COVER SIDE L 1700 Positioning GT
2	5113622	1.0	No	TOP COVER SIDE BRKT 9A positioning GT
4	5115665	1.0	No	FINGER PROTCTER L positioning GT

Table 2-22 TOP SIDE COVER L-1700 ASSY positioning GT (5133507 BOM, Rev 4)

2.1.13 TOP SIDE COVER R-1700 ASSY positioning GT (5133506)

FRU = YES

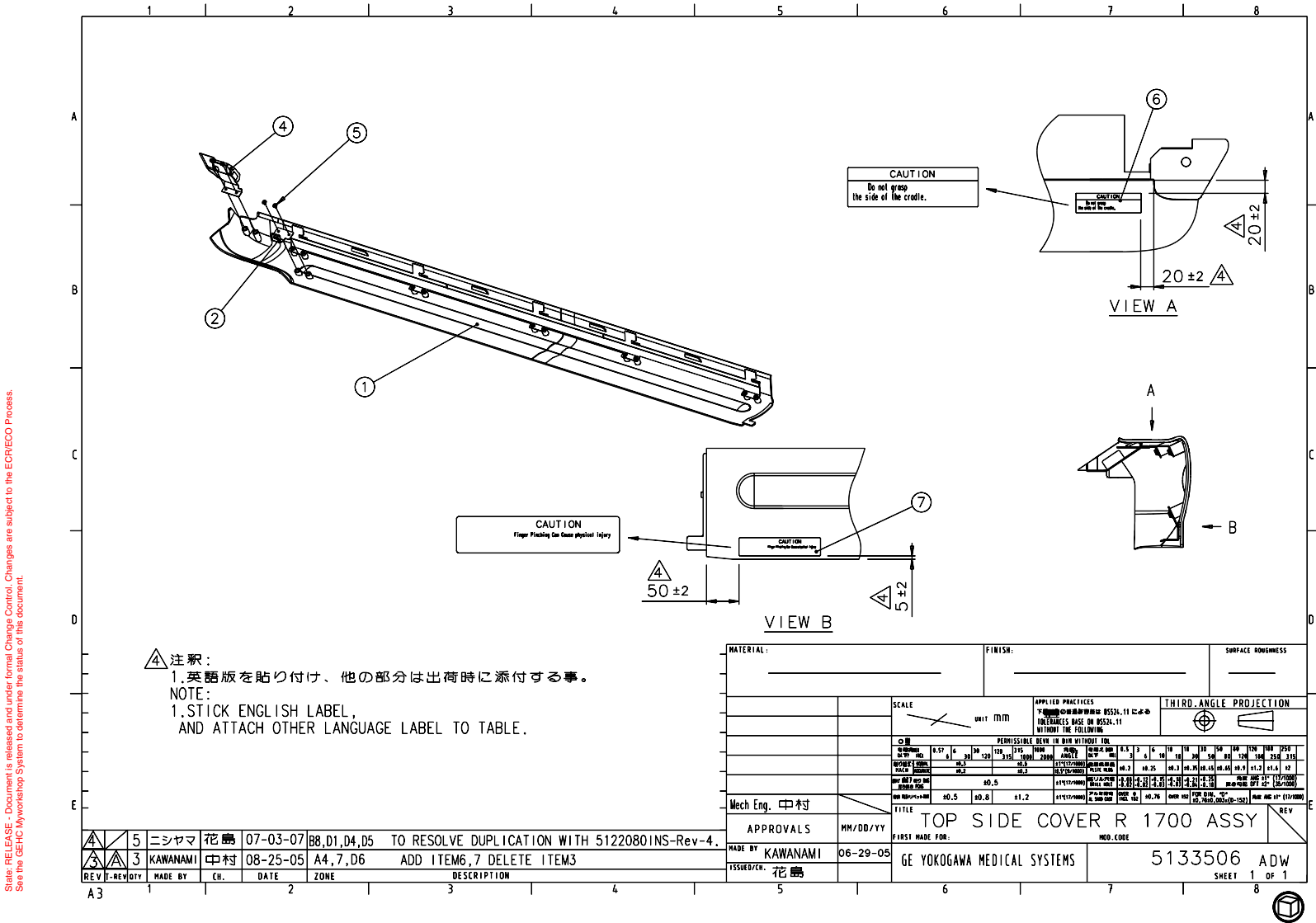


Figure 2-81 TOP SIDE COVER R-1700 ASSY positioning GT (5133506ADW, sheet 1, rev 4)

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
5	N9510HR	14.0	Yes	SCREW TP
6	2199652	1.0	Yes	P9235LZ LABEL
7	5136012	1.0	Yes	IMS FINGER PINCHING LABEL positioning GT
1	5138523-2	1.0	No	TOP COVER SIDE R 1700 positioning GT
2	5113622	1.0	No	TOP COVER SIDE BRKT 9A positioning GT
4	5115664	1.0	No	FINGER PROTCTER R positioning GT

Table 2-23 TOP SIDE COVER R-1700 ASSY positioning GT (5133506 BOM, Rev 5)

2.2 SHIP COLLECTOR GT TABLE 1700 (5138509)

FRU = NO

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2115996-3	1.0	Yes	AXIAL HEADHOLDER ASM WITH RATING PLATE AND SELF-CATCHING LATCH
2	2115990-3	1.0	Yes	LCC HSA METALLESS SUPINE CORNAL HEADHOLDER ASM WITH RATING PLATE AND SELF-CATCHING LATCH
3	2115993-2	1.0	Yes	LCC METALLESS CRADLE EXTENDER HSA METALLESS FOOT EXTENDER WITH RATING PLATE AND SELF-CATCHING LATCH
4	2169681	1.0	Yes	P9150TU BODY STRAP C
5	2169680	1.0	Yes	P9150TT BODY STRAP B
6	2169679	1.0	Yes	P9150TS BODY STRAP A
8	2331933	1.0	Yes	METALLESS COMPATIBLE PHANTOM HOLDER
11	P9140HD	1.0	Yes	BALLOON CLAMP
12	P9140TL	1.0	Yes	ARM SUPPORT ASSY
13	2152504	1.0	Yes	P9150SQ BODY STRAP C
14	2152503	1.0	Yes	P9150SP BODY STRAP B
15	2152502	1.0	Yes	P9150SN BODY STRAP A
20	5193132	1.0	Yes	FIELD EMERGENCY KIT Positioning GT
9	5133823-2	1.0	No	CRADLE PAD GT MID positioning GT
10	5122945-2	1.0	No	EXTENDER PAD 400 Positioning GT
16	5131377	1.0	No	MSDS GT TABLE Positioning GT

Figure 2-82 SHIP COLLECTOR GT TABLE 1700 (5138509 BOM, rev 5)