


System Cable Interconnect

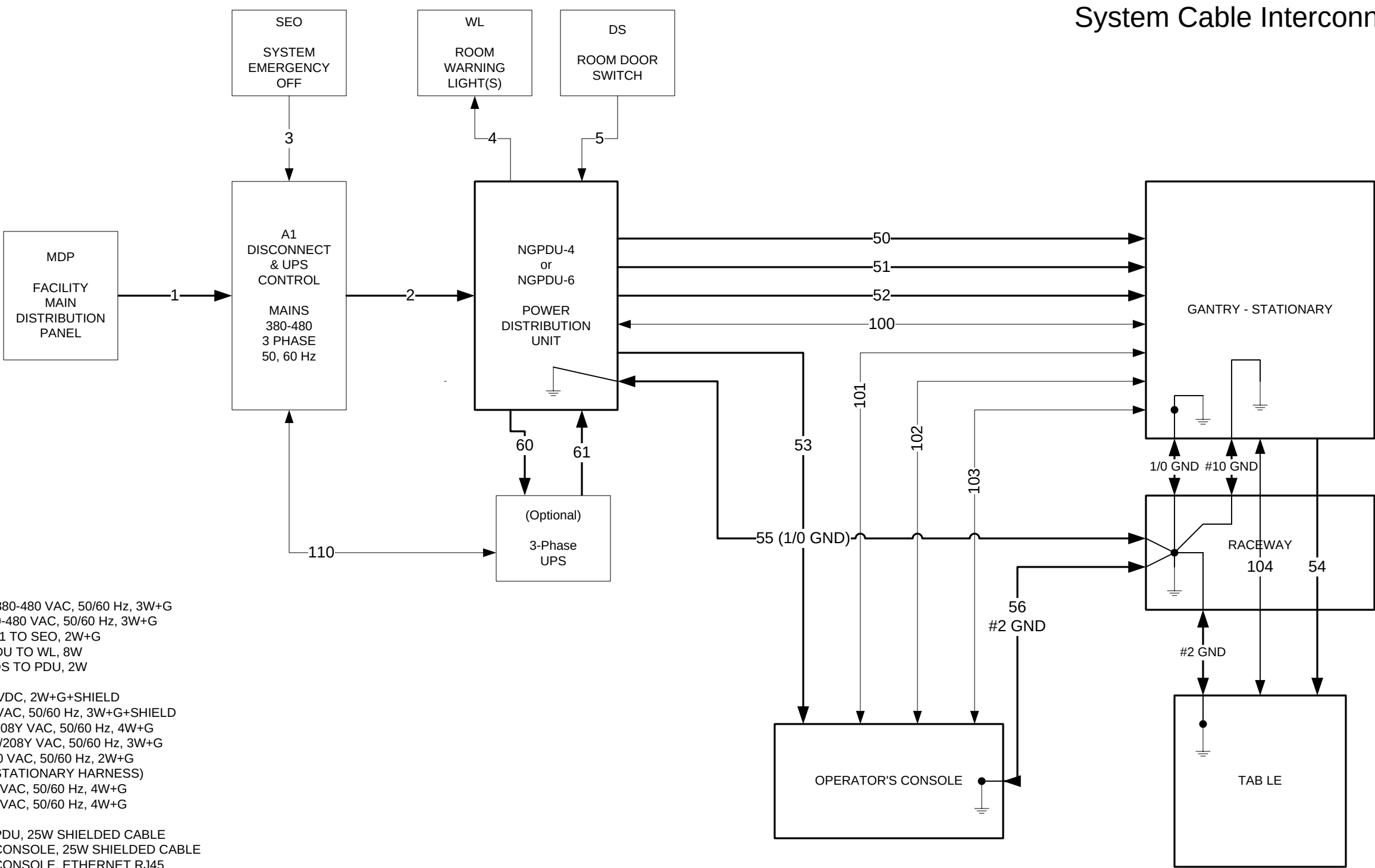

- 1* MAINS FEEDER, MDP TO A1, 380-480 VAC, 50/60 Hz, 3W+G
2* SUB-FEEDER, A1 TO PDU, 380-480 VAC, 50/60 Hz, 3W+G
3* SYSTEM EMERGENCY OFF, A1 TO SEO, 2W+G
4* X-RAY WARNING LIGHT(S), PDU TO WL, 8W
5* DOOR SWITCH INTERLOCK, DS TO PDU, 2W
- 50 HVDC, PDU TO GANTRY, 670 VDC, 2W+G+SHIELD
51 HVAC, PDU TO GANTRY, 440 VAC, 50/60 Hz, 3W+G+SHIELD
52 LVAC, PDU TO GANTRY, 120/208Y VAC, 50/60 Hz, 4W+G
53 LVAC, PDU TO CONSOLE, 120/208Y VAC, 50/60 Hz, 3W+G
54 LVAC, GANTRY TO TABLE, 120 VAC, 50/60 Hz, 2W+G
(PART OF GANTRY STATIONARY HARNESS)
60 LVAC, PDU TO UPS, 120/208Y VAC, 50/60 Hz, 4W+G
61 LVAC, UPS TO PDU, 120/208Y VAC, 50/60 Hz, 4W+G
- 100 SIGNAL CABLE, GANTRY TO PDU, 25W SHIELDED CABLE
101 SIGNAL CABLE, GANTRY TO CONSOLE, 25W SHIELDED CABLE
102 SIGNAL CABLE, GANTRY TO CONSOLE, ETHERNET RJ45
103 DATA CABLE, GANTRY TO CONSOLE, FIBER OPTIC
104 SIGNAL CABLE, GANTRY TO TABLE, 25W SHIELDED CABLE
(PART OF GANTRY STATIONARY HARNESS)
- 110 SIGNAL CABLE, UPS TO A1 , 4W+G+SHIELD

* NOTE: Runs 1, 2, 3, 4 & 5 are customer supplied. All others are supplied by GEHC.