

Effectivity

The following CT Equipment is affected:

2260619 LightSpeed Plus Gantry Assembly

Purpose

The purpose of this FMI is to replace the LightSpeed Plus Gantry base covers (left and right only) with new covers that are Cardiac Option ready.

Furnished Materials

FMI Kit Part Number 2302681 includes the items listed below.

Qty	Description	Part #
1	FMI Document	2302178-100
2	Base Cover Mounting Bracket	2285469
2	Blank Bulkheads	2288882
4	4 mm Hex Screws	46-328417P4
4	Lock Washers	46-328432P2
4	Plain Washers	46-328430P2
4	4 mm Button Head Screws	46-328420P3
4	4mm Hex Nuts	46-328425P1
1	Left Base Cover, with Cardiac Option	2285829-2, or 2285829
1	Right Base Cover, with Cardiac Option	2285830-2, or 2285830

Disposal

U.S. ONLY

- 1.) Return all obsolete parts to:

GEMS Recycling Center
c/o ACE Warehouse, Bldg. 11
2200 College Avenue
Cudahy, WI 53110

- 2.) Destroy and/or dispose of all obsolete software and documentation according to GE policy.

INTERNATIONAL SITES (OUTSIDE THE U.S.)

Destroy and/or dispose of all obsolete hardware, software and documentation according to GE policy.

People Power

Allow 1.0 hour of labor for the completion of this FMI.

Procedure

1 Gantry Preparation - Side Cover Removal

CAUTION
Potential for
injury if covers
removed and
power is left
"ON".

- 1.) Lower table to home (lowest) position.
Always remove the right side cover first, and turn "OFF" power at the STC.
- 2.) Use an 8mm Hex wrench to unlatch the side cover from the front cover. See Figure 1.

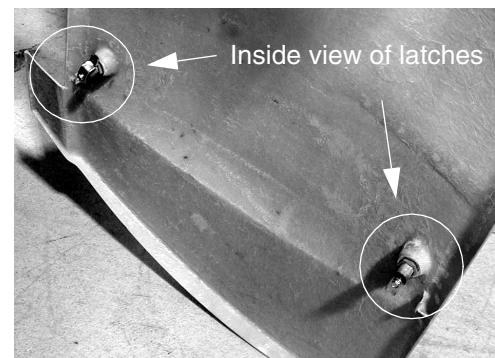
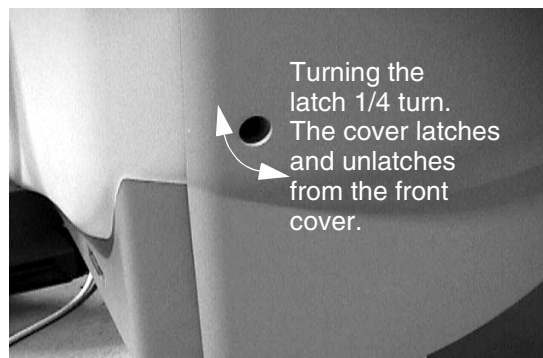


Figure 1 Side Cover Latches

- 3.) Remove the right side cover by lifting it upward to release the two (2) latches, located on the top edge of the cover. Once removed, the STC backplane should be exposed.

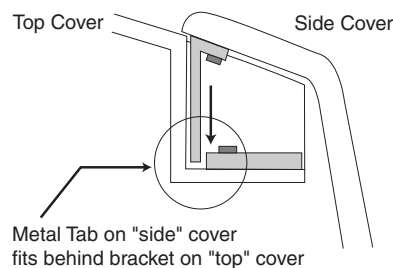


Figure 2 Side and Top Cover Clasp

- 4.) Turn off all three (3) power switches on the STC backplane.

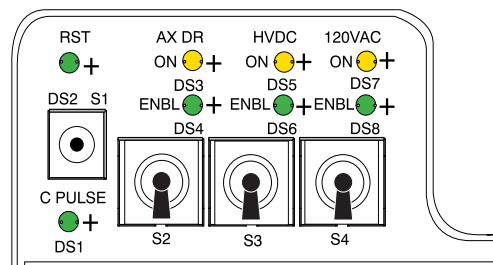


Figure 3 STC Power Switches

2 Base Cover Replacement

- 1.) Remove the existing left and right base covers. Retain the mounting screws for reuse.
- 2.) Remove the existing rear mounting brackets for the left and right base covers. Note the bracket mounting locations—the new brackets will use the same holes.
- 3.) Attach the blank bulkheads (2288882) to the new base cover mounting brackets (2285469), using flat washers and nuts.
- 4.) Install the new base cover mounting brackets (2285469).
- 5.) Install the new left and right base covers.

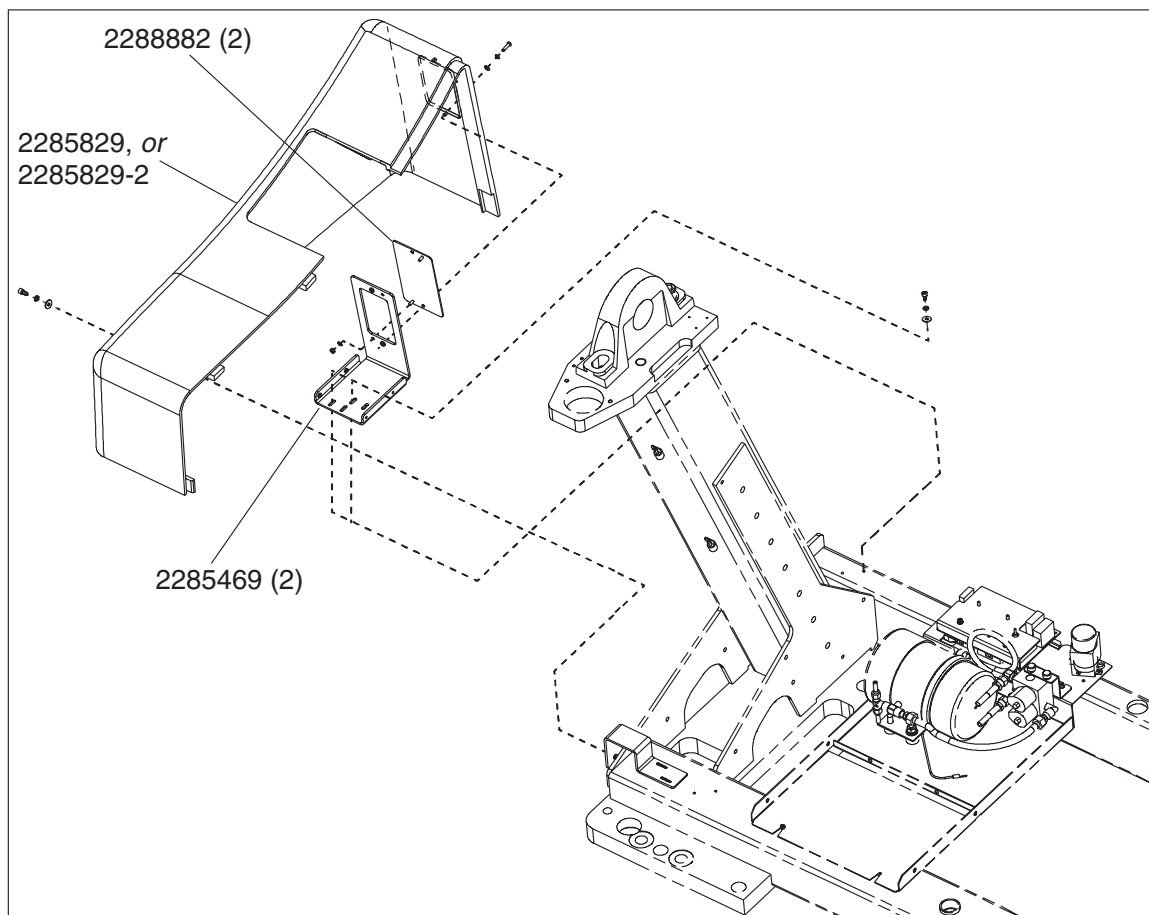


Figure 4 Installation of new left base cover

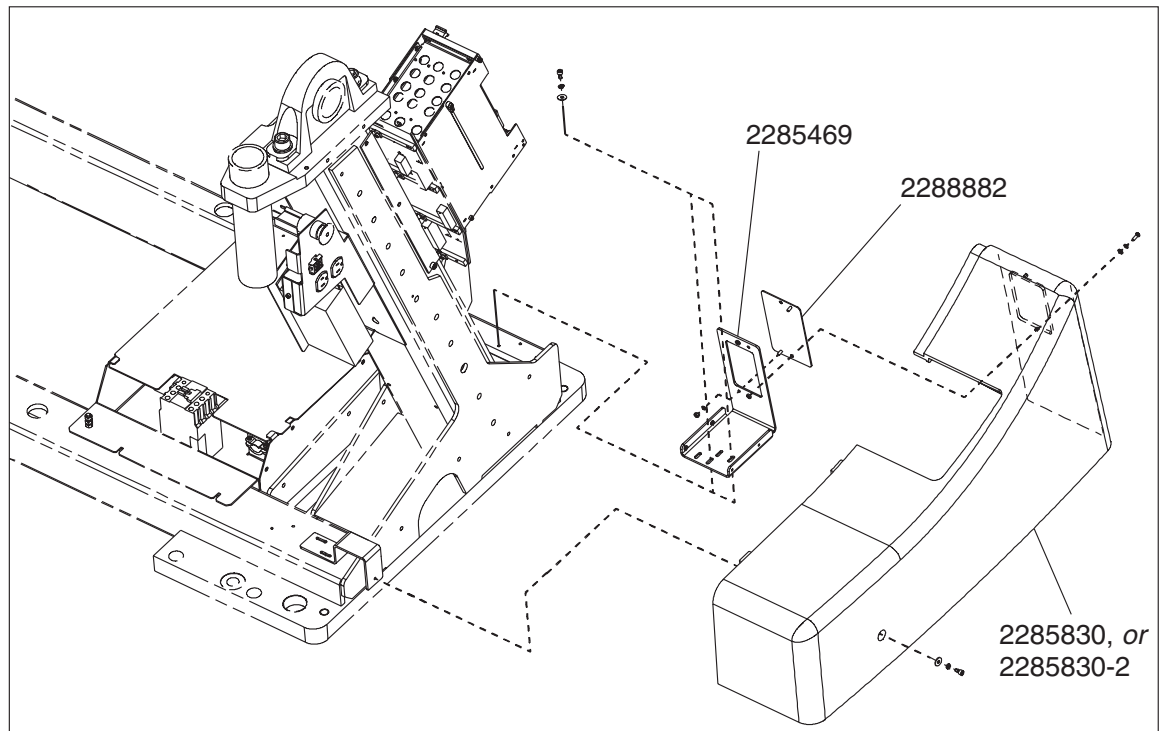


Figure 5 Installation of new right base cover

3 Gantry Reassembly

Reverse the order of the steps in *1 Gantry Preparation - Side Cover Removal*, to return the gantry to its operational-ready state.

End of Document