

**ADMINISTRATIVE INFORMATION****FMI No: 25407****FMI Title: Service Pack 2.0 for Software Versions 07MW18.4,
07MW11.10, and 07BW08.4**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☒ Mandatory Safety	GEMS - AM:	Issue Date:
☐ Mandatory	GEMS - E:	
☐ Mandatory on Request	GEMS - A:	Close Date:

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, Contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center 800-321-7937 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Please contact a Support Engineer in your Region!

Special Instructions:

- Can be installed remotely for the following system types: HiSpeed QX/i, LightSpeed VCT, and LightSpeed 1X – 5X (excluding HP60, RT¹⁶, and Xtra)
- Does not apply to CT/PET systems!
- Does not apply to RT16 systems!
- Does not apply to ATOM systems!

Unit Tracked: Multiple Host Console Model Numbers**FMI Instructions Part No.:** 5264685-1EN, Rev 2**Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center 800-321-7937 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Send to your Country FMI or Service Administrator

FMI No.: 25407

FMI Title: Service Pack 2.0 for Software Versions 07MW18.4,
07MW11.10, 07BW08.4

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not

Customer Refusal or Equipment Altered - Modification Pending.

Installer
Comments

Table of Contents

Table of Contents	1
Revision History.....	2
Effectivity	3
Purpose	4
Prerequisite FMIs	4
Furnished Materials.....	5
Tools Required	5
Personnel Requirements.....	5
Installation Notes.....	7
Remote Installation Procedure	8
1. BEFORE YOU BEGIN.....	8
2. REMOTE INSTALLATION CHECKS [HiSPEED QX/i, VCT, AND LIGHTSPEED 1X – 5X (EXCLUDING HP60, RT ¹⁶ , AND XTRA) ONLY]	8
3. SAVE SYSTEM STATE – BEFORE SP2	9
4. UNINSTALLING PREVIOUS SERVICE PACKS AND PATCHES	10
5. INSTALLING SERVICE PACK 2.0	11
6. SANITY SCANS (ON-SITE ASSISTANCE REQUIRED)	16
7. SAVE SYSTEM STATE – AFTER SP2.0.....	16
Finishing Up.....	17
1. CHECK AUTOVOICE VOLUME SETTINGS	17
2. DOCUMENTATION	17
3. FMI CLOSURE	17
4. VERIFICATION OF FUNCTIONALITY.....	17
Appendix A: Download Troubleshooting	18

Revision History


Revision	Date	Reason for Change
1	6/13/2008	Initial Release
2	8/20/2008	Updated with F3 pilots feedback

Effectivity


This FMI is to be installed on the following **CT-only** systems: HiSpeed QX/I, BrightSpeed (4, 8, and 16-slice), BrightSpeed Select (4, 8, and 16-slice), LightSpeed 1X, 2X, 3X, 4X, 5X (excluding LightSpeed RT¹⁶ and LightSpeed Xtra), 7X, RT4.

The CT system must be configured with a **GRE (GOC3, GOC4, or GOC5)** host computer console or gantry with one of the following tracked model numbers:

- 2377708
- 2377708-4
- 2377708-7
- 2377708-10
- 2377708-13
- 2377708-14
- 2377708-20
- 2377708-21
- 2377708-22
- 2377708-30
- 2377708-31
- 2377708-32
- 2377708-40
- 2377708-41
- 2377708-42
- 2377708-43
- 2377708-50
- 2377708-51
- 2377708-52
- 2377708-53
- 2377708-56
- 2377708-120
- 2377708-520
- 2377708-550
- 2377708-551
- 2377708-552
- 2377708-553
- 5111750
- 5377708-120
- 5377708-520
- 5377708-530
- 5124069
- 5124069-2
- 5124069-3
- 5124069-4
- 5124069-5
- 5124069-6
- 5126093

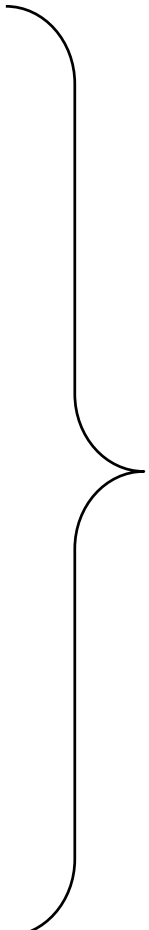

HOST CONSOLE TRACKED MODEL NUMBERS

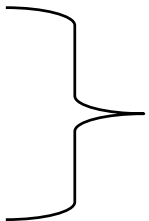

GANTRY TRACKED MODEL NUMBERS

DO NOT INSTALL this FMI on any of the following systems:

- PET/CT
- LightSpeed RT¹⁶
- LightSpeed Xtra
- ATOM

Forward Production cut in date: FW25, 2008.

Purpose

This FMI delivers Service Pack 2.0 to CT-only systems except for RT¹⁶. Service Pack 2.0 provides 8 safety and 108 quality fixes. Additionally, Service Pack 2.0 also provides support for XW8400 host computer, and the new ICOM board. This FMI can be completed remotely for the following system types: HiSpeed QX/i, LightSpeed VCT, and LightSpeed 1X – 5X (excluding HP60, RT¹⁶, and Xtra)

For a more complete list of fixes, please see TWA (Tips and Workarounds document) included in the FMI kit and available online in Common Documentation Library.

Prerequisite FMIs

The CT-only system must be loaded with one of the software releases listed below before installing this FMI:

- 07MW11.10
- 07MW18.4
- 07BW08.4

STOP if your system does not have one of these software releases already loaded. One of the software releases listed above must first be installed before proceeding with this FMI. Please refer to table below to identify the prerequisite FMI based on your system.

Note: If your site has closed FMI 25389 or 25394 as a Code 4 (Customer Refusal), and you are now installing one of those FMI's (25389 or 25394), a new dispatch will need to be opened for those FMI's, then closed as Code 1.

FMI Number	Title
FMI 25389	07MW11.10 for LS 1X2X3X4X5X, HS QX/i, & BSD 4/8/16
FMI 25394	07MW18.4 VCT Quality Release

STOP if your system is a PET/CT system. This FMI is not to be installed on PET/CT systems. Close this FMI dispatch as Code 2 and return the FMI kit.

Furnished Materials

The following items are included in FMI 25407, kit part number 5264685, Rev 2:

Item	Name	Rev	Qty	Description
1	5264685-1EN	2	1	FMI 25407 Instructions - English
2	5234959-20	2	1	Service Pack 2.0 Collector for SW releases 07MW11.10, 07MW18.4, 07BW08.4
3	5309078-1EN	5	1	Service Pack 2.0 Install Instructions
4	5265325-299	1	1	CT Service Pack 2 Common Software Addendum TWA - MULTI-LANGUAGE
5	5265325-1EN	1	1	CT Service Pack 2 Common Software Addendum TWA - ENGLISH
6	5265325-1JA	1	1	CT Service Pack 2 Common Software Addendum TWA - JAPANESE
7	5265325-1ZH-CN	1	1	CT Service Pack 2 Common Software Addendum TWA - CHINESE FOR CHINA
8	2337414	4	1	DVD RAM MEDIA, 9.4 GB

Tools Required

No tools are required for this FMI.

Personnel Requirements

- **REMOTELY OR ONSITE** - VCT, HiSpeed QX/i, and LightSpeed 1X – 5X (excluding HP60, RT¹⁶, Xtra)
 - Travel: 1 Hour
 - Labor: 2 Hour
 - Field Engineer (RSE): 1
- **ONSITE** - BrightSpeed Power, BrightSpeed Select, and LightSpeed 5X HP60
 - Travel: 1 Hour
 - Labor: 3 Hours
 - Field Engineer: 1

Installation Notes

This FMI has two methods of installation:

A. On-site via CD.

For **ON-SITE installation**, use Service Pack 2.0 Install Instructions (part # 5309078-1EN, included in this kit) to install the Service Pack. Skip all remaining steps in Section 1 of this document and proceed directly to section 2 (Finishing Up), to complete FMI.

B. Remotely via ISD (Integrated Service Desktop)

For **REMOTE installation** instructions, proceed to Section 1. Remote installation is applicable to the following systems with Broadband connectivity:

- HiSpeed QX/i
- LightSpeed VCT
- LightSpeed 1X – 5X (excluding HP60, RT¹⁶, and Xtra)

NOTE: Remote installation **IS NOT** applicable to any BrightSpeed, BrightSpeed Select, or HP60's systems, as they require additional calibrations.

NOTE: Service Pack 2 will be automatically downloaded, (but NOT auto-installed) to sites that meet the following criteria:

1. Are listed on FMI 25407 Effectivity List at the time of F4 (full) release.
2. Have functioning broadband connectivity.
3. Have ProDiags task sw_download enabled (can be enabled remotely; see appendix A).
4. Have Auto-Download enabled (can be enabled remotely; see appendix A).

* For those sites that should have Service Pack auto-downloaded, but do not, see **Appendix A** for troubleshooting.

Remote Installation Procedure

1. Before you Begin

- Please note that this FMI can be installed remotely for the following systems **ONLY**:

- HiSpeed QX/i
- LightSpeed VCT
- LightSpeed 1X – 5X (excluding HP60, RT¹⁶, and Xtra)

NOTE: Remote installation **IS NOT** applicable to any BrightSpeed, BrightSpeed Select, or HP60's systems, as they require additional calibrations.

- If this FMI is remotely installed, please make sure to **schedule time with customer in advance**. Also, please let them know that you will need their assistance during the installation to perform the final sanity scans.

2. Remote Installation Checks [HiSpeed QX/i, VCT, and LightSpeed 1X – 5X (excluding HP60, RT¹⁶, and Xtra) only]

Follow the steps listed below to remotely install this FMI:

1. Confirm that the site is InsiteTM connected (remote installation is recommended for sites with Broadband connectivity only).
2. Become **ctuser**
3. Verify that one of the following software versions is installed: 07MW11.10 or 07MW18.4

```
{ctuser@hostname} swwhinfo -all
```

If **not** 07MW11.10 or 07MW18.4, CANNOT do remote install.

If the system is an HP60 (noted by the *HP_P* string), CANNOT do remote install.

4. Was SP2.0 previously installed?

As root, type: [root@hostname ~]# **showprods | grep ServicePack**

Check output to see if SP2.0 is installed.

- If Service Pack 2.0 is already installed on the system, close this FMI dispatch as Code 1 and return the FMI kit.
- If Service Pack 2.0 has not been installed, continue with step 5.

5. Verify that Service Pack 2.0 was successfully downloaded and is located in the following directory. Additionally, verify that the size is 148MB. If the size is not correct, re-download. Please see Appendix A to re-download.

```
ls -l /usr/g/ctuser/sw_download/patches
```

Output Displayed:

```
-rwxrwxrwx 1 ctuser users 148672212 Jun 10 14:07 07MW11.10_H1_H2_GEHCSERVICEPACK2.0-2-26.i386.rpm*
```

If you do NOT see a listing for Service Pack 2.0, refer to **Appendix A** for Download Instructions. Note that you may also see a listing for Service Pack 1.0 and/or Service Pack 1.5. Disregard these two.

If you know that your system **IS NOT** an H-Power system, **OR** if you know that your system has a **TGPU** board instead of a TGP board, skip step 6 and proceed to step 7.

6. Please verify the TGP board part number:

WARNING: Customer CANNOT be scanning or using the system during step 6.

H-Power systems with TGP board part number 2290138-0500 will need a manual gantry reset during the installation. Remote install is **NOT** applicable for systems with this TGP board.

```
{ctuser@hostname} rlogin tgp
type: print_cpu_board_revision
Output displayed:
PartNo: 2290138-0500 (Variation:5 Revision:0)
value = 46 = 0x2e = '.'
type: logout to exit
```

If the output displayed is **PartNo# 2290138-0500**, then Remote Installation of SP2.0 **CANNOT** be completed, load SP2.0 via CD onsite.

7. If the previous steps checked out successfully, then continue on with Save System State in the next section.

3. Save System State – Before SP2

NOTE: Make sure that applications are running before proceeding with Save System State

A current System State, labeled “**Before SP2**”, must be available before attempting any software or hardware modification. A second System State, labeled “**After SP2**”, will be saved following Service Pack 2 load. The following steps will save a copy of System State to the hard drive. On the next visit to the site, the saved version of System State should be put onto a DVD. Those commands are also listed below.

Complete Save State process now using the following steps:

1. Create disk directory. Must be **ctuser**

```
{ctuser@hostname} mkdir /usr/g/systemstate (dir may already exist. That's ok)
{ctuser@hostname} rm -rf /tmp/systemstate
{ctuser@hostname} ln -sf /usr/g/systemstate /tmp/systemstate
```

2. Start the command line version of SaveAll SystemState to /tmp directory:

```
{ctuser@hostname} sysStateABB -saveAll 3
{ctuser@hostname} cd /usr/g/systemstate
```

3. Now tar and compress all files into a single one:

```
{ctuser@hostname} tar zPpcvf /usr/g/ctuser/stateBeforeSP2.tgz . (Include period)
{ctuser@hostname} cd /usr/g/ctuser
```

```
{ctuser@hostname} ls -l state*
```

The **.tgz** file displayed by the command above is the compressed system state file, which can be manually copied to a DVD.

Also a 2nd backup to the DARC (less than 10 seconds)

Type: **rcp -p stateBeforeSP2.tgz ctuser@darc:~;**

4. On the next visit to the site, make sure to save the system state files in the directory above to a system state DVD using the commands below:

➤ Mount the DVD: {ctuser@hostname} **su -**

Password: **#bigguy**

[root@hostname ~]# **mountDVD**

(Note: the following commands assume the DVD has been previously formatted as a systemstate)

- Untar the file: [root@hostname ~]# **tar xpvfz /usr/g/ctuser/stateBeforeSP2.tgz -C /DVD**
- UnMount the DVD: **unmountDVD**

4. Uninstalling Previous Service Packs and Patches

Service Pack 2 will not be able to install successfully if **Service Pack 1.0 / Service Pack 1.5** is installed on the system. Additionally, Application software must be down to install or uninstall the Service Pack.

1. Verify whether any previous Service Packs/Patches are already installed on the system:

```
[root@hostname ~]# showprods | grep ServicePack
```

```
[root@hostname ~]# showprods | grep Patch
```

2. If no output is displayed by the two commands above, proceed to "Installing Service Pack 2.0" section; otherwise proceed to step 3.

3. Via **ISD** pull down menu, shutdown applications. Wait approximately 1-2 minutes for Application software to go down.

Applications are down when the following output is displayed in the telnet window.

cleanMon.oc[12133]: OC applications down, cleanMon.oc exiting

4. In a telnet window.

5. Become root:

```
{ctuser@hostname} su -
```

6. Password: **#bigguy**

7. As root on the Host type:

8. [root@hostname ~]# **patch_install -u**

9. To uninstall all existing patches type: **y** (Answer yes to each existing patch)

10. Verify that all the previous Service Packs and Patches uninstalled correctly:

```
[root@hostname ~]# showprods | grep ServicePack
```

```
[root@hostname ~]# showprods | grep Patch
```

No output should be displayed by these commands.

11. Perform reboot command to restart the console:

```
[root@hostname ~]# reboot
```

NOTE: The **reboot** command will disconnect the telnet window. Please wait 4-5 minutes before re-connecting. After reboot is complete, open a new telnet window and continue.

5. Installing Service Pack 2.0

This section outlines the steps to load Service Pack 2. A reboot is **ALWAYS** necessary after the Service Pack load.

1. In a telnet window.
2. Verify one of the following base software versions is installed:

```
{ctuser@hostname} swwhwinfo -v
```

```
07MW11.10
```

```
07MW18.4
```

3. Verify the current keyboard language. As *ctuser* type:

```
{ctuser@hostname} more /etc/sysconfig/keyboard
```

Output Displayed: KEYBOARDTYPE="pc"

```
KEYTABLE="us"
```

(Sample output only. Output will vary based on keyboard language). **PLEASE NOTE THE KEYBOARD LANGUAGE DENOTED BY THE "KEYTABLE" ENTRY IN THE OUTPUT DISPLAYED. THE KEYBOARD LANGUAGE SELECTION WILL BE NEEDED LATER.**

4. Shutdown applications using the **ISD** pull down menu. Wait approximately 1-2 minutes for Application software to go down.

Applications are down when the following output is displayed in the telnet window.

cleanMon.oc[12133]: OC applications down, cleanMon.oc exiting

5. Become root:

```
{ctuser@hostname} su -
```

Password: **#bigguy**

6. It is critical that the DARC is up prior to loading Service Pack, verify that the DARC is up.

```
[root@hostname ~]# rsh darc
```

```
[root@darc ~]# exit
```

7. Type the Service Pack installation command:

STOP and READ THIS: Read screen install questions carefully; system may ask if you want to install SP1.0 or SP1.5. **DO NOT install SP1.0 or SP1.5!**

Only install SP2.0 at this time!

```
[root@hostname ~]# patch_install -i
```

ATTENTION: System may ask if you want to install SP1.0 or SP1.5. **DO NOT install SP1.0 or SP1.5!**

Only install SP2.0 at this time!

Wait for the following prompt to display before moving on: "I will install 07MXX.XX_GEHCSERVICEPACK2.0 from CDROM, is this ok? [y/n]"

8. Type: y

9. The load will take approximately **5-10 minutes**. The following output is displayed:

NOTE: Output displayed will vary based on software version.

07MW11.10_Gantry Info_GEHCSERVICEPACK2 update has been successfully installed.

10. Verify that the Service Pack loaded correctly using the following command:

```
[root@hostname ~]# showprods | grep ServicePack
```

NOTE: Output displayed will vary based on software version.

07MW11.10_Gantry Info_GEHCSERVICEPACK2 DATE 07MW11.10_Gantry Info Service Pack 2 ver x rel x.xx

11. Verify the Application software package:

```
{ctuser@hostname} swhwinfo
```

NOTE: Service Pack version is now included in the *swhwinfo* string. Output displayed will vary based on software version.

Example: *07MW11.10_SP2-0-26.Hardware info*

12. Verify that Application software is DOWN.

Applications are down when the following appears in the telnet window.

cleanMon.oc[12133]: OC applications down, cleanMon.oc exiting

NOTE: Entering reconfig, network menu and selecting accept is required for changes made to correct station name issue.

13. Continue with Remote Reconfig in the following steps.

14. Become root:

```
{ctuser@hostname} su -
```

Password: **#bigguy**

14. Set Display so the GUI appears on your laptop.

- Remove previous uses of CSD

```
[root@hostname ~]# rm /tmp/sdm_temp*
```

Answer yes to remove each temp file listed.

- Start CSD by selecting the Common Service Desktop (CSD) button in telnet window of ISD. Login as "Insite" and enter CSD password (CSD Redirect was avoided because not all users have access to Exceed and CSD Redirect). Leave CSD window open, as it will be needed later.
- In the Telnet window that was previously open, view contents of newly created files:

```
[root@hostname ~]# cat /tmp/sdm_temp*
```

Example output:

150.2.12.158

125.0

- Create the command line with the above info. See example below:

```
[root@hostname ~]# setenv DISPLAY 150.2.12.158:125.0
```

15. Start the reconfig GUI (GUI will open in the CSD window opened previously)

```
[root@hostname ~]# reconfig
```

NOTE: An Error in **tkerror** will occur. That is normal.

16. From the open GUI, select **[Config]**

17. It is now necessary to set the appropriate keyboard language based on the setting before Service Pack 2. The **[Preferences]** tab will be displayed. Set the **Keyboard Selection** field based on the keyboard language recorded in section 5, step 3. Select the appropriate keyboard language from the **[Preferences]** tab GUI by referring to the table below:

English	"us"
French	"fr"
German	"de"
Italian	"it"
Spanish	"es"
Swedish	"sv"
Portugues e	"br"

(Brazilian)	
-------------	--

WARNING: DO NOT SELECT ONE OF THE FOLLOWING KEYBOARD SELECTIONS UNLESS YOU HAVE THE CORRESPONDING KEYBOARD: DANISH, EUROPEAN PORTGUESE, DUTCH, AND NORWEGIAN

18. Select the **[Network]** tab

At this time, you can set the **Station Name**. **Station Name** is now an editable field. If your customer requires a specific *Station Name* for network purposes, it should be set now. Consult with your customer and set accordingly using the guidelines below. For most cases the *Station Name* default value (Host Name) is used. If no special *Station Name* setting is required, select **[Accept]** to set to default, and proceed to step 20.

Station Name guidelines:

- The *Station Name* identifies the text that will be stored in the DICOM attribute "Station Name" (DICOM TAG 0008,1010) in all CT images created. Some pacs systems may use the Station Name attribute or system tracking for identification of the scanner where the CT images were generated. Typically, the Host Name and Station Name will be the same.
- *Station Name* MUST NOT exceed 16 Characters
- *Station Name* MUST only contain the following characters: a through z, A through Z, 0 through 9, - and _
- *Station Name* will default to the Host Name if left blank.

19. Select the **[Accept]** tab (Long pause. Be patient, wait until reconfig completes)

20. When prompted to reboot, select **[Yes]**

NOTE: Please wait until reboot is completely finished before moving on to next step. The **reboot** will disconnect the telnet window. Please wait 4-5 minutes before re-connecting. After reboot is complete, open a new telnet window and continue.

21. Become **root**:

```
{ctuser@hostname} su -
```

Password: **#bigguy**

The following command will restore the configuration file that holds the users HIPPA account information.

NOTE: This file is different for different systems. If not sure of system, type `swhwinfo` to get software and gantry info.

for 07MW18.4 VCT :

```
[root@hostname ~]# cd /etc
```

```
[root@hostname ~]# cp _ea3.config.pre_07MW18.4_GEHCSERVICEPACK2.0-2-26 ea3.config
```

Output Displayed: **cp: overwrite 'ea3.config'?**

Type: **y**

for 07MW11.10 (H1, H2):

```
[root@hostname ~]# cd /etc
```



```
[root@hostname ~]#
```

```
cp _ea3.config.pre_07MW11.10_H1_H2_GEHCSERVICEPACK2.0-2-26 ea3.config
```

Output Displayed: **cp: overwrite 'ea3.config'?**

Type: **y**

for 07MW11.10 & 07BW08.4:

```
[root@hostname ~]# cd /etc
```

```
[root@hostname ~]#
```

```
cp _ea3.config.pre_07MW11.10_07BW08.4_GEHCSERVICEPACK2.0-2-26 ea3.config
```

Output Displayed: **cp: overwrite 'ea3.config'?**

Type: **y**

22. Check if **DENTASCAN** option is installed on the system:

```
{ctuser@hostname} swokinstall -p
```

The list of all options installed on the system will be displayed. Please look for **DENTASCAN**. If **DENTASCAN** is NOT installed on the system, skip the rest of this step and move on to step 23; if **DENTASCAN** is installed, perform the following command:

```
{ctuser@hostname} installcamera.denta
```

23. As ctuser type: **startup**

This brings applications back up after a reconfig. (Wait at least 5 minutes for Apps to come up)

24. Close pop-up message windows

This command will close any pop-up windows that appear on Console Monitor after reboot.

```
{ctuser@hostname} slay -f attention
```

25. From the Common Service Desktop icon, select **[Utilities]**

Select **[Install]**

Select **[Flash Download Tool]**

26. Select **[Query]** (Wait for completion, Query command will indicate which controllers need to be flashed. Continue with update even if a controller fails).

27. Select **[Update]** (Please be patient; the update may take several minutes).

28. If Flash Download fails on the first pass, try running the Update again before attempting a hardware reset via the Common Service Desktop.

29. Once Flash Download is completed successfully, select **[Dismiss]**

6. Sanity Scans (ON-SITE ASSISTANCE REQUIRED)

ATTENTION: At this point, you will need the customer's assistance to perform the sanity scans. **DO NOT** perform the following steps without being in DIRECT contact (i.e. via phone) with your customer.

1. Ask the customer to place a phantom on the CT table
2. Confirm that the customer has left the scanning room and the door is closed.
3. Ask the customer to perform at a minimum one Scout and a series of 20 Helical or Axial scans
4. Monitor the error log for any unexpected messages during or after the scans
5. Ask the customer to confirm that all scans are completed successfully.
6. Please note the exam number of the sanity scans. This will be recorded during debrief.

NOTE: It is **NOT** required to check the Image Quality of the scans taken. Only confirm that all scans were completed successfully with no issues.

Customer assistance is not required for the rest of this install. Please remember to thank them for their help.

7. Save System State – After SP2.0

Perform the following steps to save a second System State on the disk, labeled "**After SP2**". On the next visit to the site, the saved version of System State should be put onto a DVD. Those commands are also listed below.

Complete Save State process now using the following steps:

1. Create disk directory:

```
{ctuser@hostname} mkdir /usr/g/systemstate  
{ctuser@hostname} rm -rf /tmp/systemstate  
{ctuser@hostname} ln -sf /usr/g/systemstate /tmp/systemstate
```

2. Start the command line version of SaveAll SystemState to /tmp directory:

```
{ctuser@hostname} sysStateABB -saveAll 3
```

Now, *cd* into /usr/g/ctuser directory

```
{ctuser@hostname} cd /usr/g/systemstate
```

3. Now tar and compress all files into a single one:

```
{ctuser@hostname} tar zPpcvf /usr/g/ctuser/stateAfterSP2.tgz . (include period)
```

```
{ctuser@hostname} cd /usr/g/ctuser
```

```
{ctuser@hostname} ls -l state*
```

The **.tgz** file displayed by the command above is the compressed system state file, which can be manually copied to a DVD.

Also a 2nd backup to the DARC (less than 10 seconds)

Type: **rcp -p stateAfterSP2.tgz ctuser@darc:~;**

- On the next visit to the site, inform primary FE to perform a normal **SystemStateDVD SaveALL** per system documentation. The above processes can be retrieved if necessary. There are two copies of the **stateAfterSP2.tgz** file, one on the HP and one on the DARC. Please see FMI need to know document for retrieval process.

Finishing Up

1. Check Autovoice Volume Settings

Service Pack 2 will increase the *Autovoice Volume* setting. If the autovoice volume is too high, try the following:

1. Adjust the volume using the knob at the operator console. If this FMI is being done remotely, have the customer adjust the knob.
2. If the above step does NOT fix the issue, select the “*Autovoice Volume*” (“audioSetting”) tool from the *Tool Chest* menu to decrease the volume as needed.

2. Documentation

The following set of documents should be left with the customer or delivered to the customer on your next visit:

Item	Name	Rev	Qty	Description
1	5234959-20	2	1	Service Pack 2.0 Collector for SW releases 07MW11.10, 07MW18.4, 07BW08.4
2	5309078-1EN	5	1	Service Pack 2.0 Install Instructions
3	5265325-299	1	1	CT Service Pack 2 Common Software Addendum TWA - MULTI-LANGUAGE
4	5265325-1EN	1	1	CT Service Pack 2 Common Software Addendum TWA - ENGLISH
5	5265325-1JA	1	1	CT Service Pack 2 Common Software Addendum TWA - JAPANESE
6	5265325-1ZH-CN	1	1	CT Service Pack 2 Common Software Addendum TWA - CHINESE FOR CHINA
7	2337414	4	1	DVD RAM MEDIA, 9.4 GB

3. FMI Closure

Note the Sanity Scans Exam number. Close FMI per regional procedure.

4. Verification of Functionality

Add the Sanity Scans Exam number to the debrief comments as follows:

Exam #: xxxxx completed as verification of functionality

NOTE: Due to changes with the Voxel SW, customer may need to reset film save option.

Appendix A: Download Troubleshooting

This section is for those sites that should have Service Pack 2.0 auto-downloaded on the system, but **DO NOT**. The steps below describe the process of enabling the ProDiags and Auto-Download features:

1. Confirm that site is on the FMI 25407 Effectivity List
2. Verify that site is Broadband connected with Insite™ connectivity
3. Become root
4. Check if ProDiags is already enabled:

```
[root@hostname ~]# crontab -l | grep -i prodiags
```

If an output similar to the following is displayed, then ProDiags is already enabled. Skip steps 5 and 6, and proceed with step 7. If you **DO NOT** see a similar output, then ProDiags is **NOT** enabled; proceed to step 5:

```
1,16,31,46 * * * * /usr/g/insite/bin/cronWrapper /usr/g/insite/ProDiags/bin/proDiagsCheck >/dev/null 2>&1
```

5. Enable ProDiags to use the default schedule and task settings:

```
[root@hostname ~]# cd ~insite/ProDiags/bin
```

```
[root@hostname ~]# pd_install -auto
```

6. Perform step 4 again to ensure that ProDiags was enabled successfully.
7. Open the Common Service Desktop
8. Select the **[Utilities]** tab
9. Select **[Update Management]**


NOTE: Service Key is required for onsite users to access Update Management Menu. Remote users coming in as insite do not need a Service Key.

10. Toggle the **[Enable/Disable Auto Download]** button until Auto-Download is enabled. The following output will be displayed when Auto-Download is enabled:

“Auto Download of updates is enabled in this scanner”

11. Select **[Query Updates]**
12. If the system displays the Service Pack 2.0 update, select **[Download Updates]** to download Service Pack 2. After the download is completed, close the Common Service Desktop and return to section 1.

13. If after selecting **[Query Updates]**, you get an error message similar to the one below, proceed to the next steps to resolve it. NOTE: this error will be fixed with Service Pack 2.

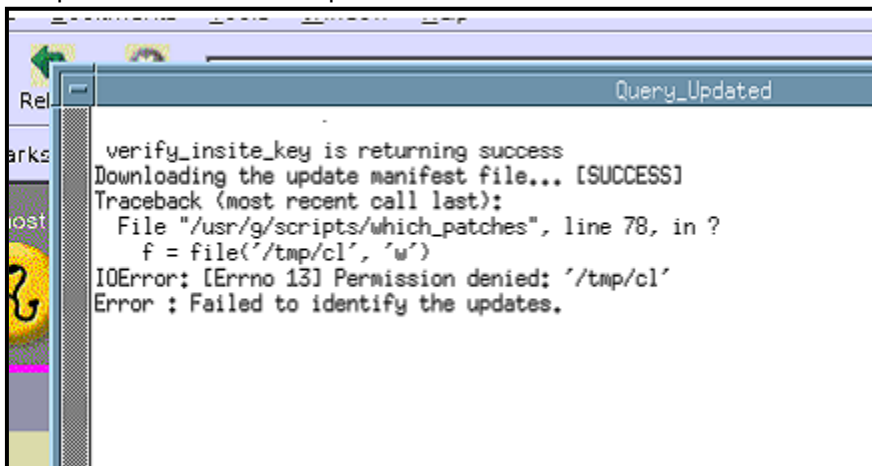

14. You will need to change ownership of the following files:

/tmp/cl

/tmp/cl1

/tmp/cl2

Become root:

```
[insite@ncict1 patches]$ su - root
```

Password: **enter root password**

Change directory to /tmp and list cl files. Note that ownership is likely root.*

```
[root@ncict1 ~]# cd /tmp
```

```
[root@ncict1 /tmp]# ls -al cl*
```

```
-rw-r--r-- 1 root root 9 Jun 30 11:16 cl
```

```
-rw-r--r-- 1 root root 2 Jun 30 11:16 cl1
```

```
-rw-r--r-- 1 root root 27 Jun 30 11:16 cl2
```

```
-rw-rw-rw- 1 insite root 733 Jul 1 10:21 clientdrv.trace
```

Remove files cl cl1 cl2, answer yes to remove each one.

```
[root@ncict1 /tmp]# rm -r cl cl1 cl2
```

```
rm: remove regular file `cl'? y
```

```
rm: remove regular file `cl1'? y
```

```
rm: remove regular file `cl2'? y
```

Recreate cl, cl1, cl2 with ctuser ownership

```
[root@ncict1 /tmp]# su - ctuser
```

```
{ctuser@ncict1}[1] touch /tmp/cl
```

```
{ctuser@ncict1}[2] touch /tmp/cl1
```

```
{ctuser@ncict1}[3] touch /tmp/cl2
```

Now note that ownership has changed to ctuser

```
{ctuser@ncict1}[4] cd /tmp
```

```
{ctuser@ncict1}[5] ls -al cl*
```

```
-rw-rw-r-- 1 ctuser users 0 Jul 1 15:51 cl
```

```
-rw-rw-r-- 1 ctuser users 0 Jul 1 15:51 cl1
```

```
-rw-rw-r-- 1 ctuser users 0 Jul 1 15:51 cl2
```

```
-rw-rw-rw- 1 insite root 733 Jul 1 10:21 clientdrv.trace
```

15. Now retry to **[Query Updates]**

16. You may get a message stating that “**there are no relevant updates available**”. This may mean that:

- The software revision of the scanner is not compatible with the Service Pack, not current revision or not applicable.
- The Service Pack has already been loaded onto the scanner... check **/user/g/ctuser/sw_download/patches** directory to see if it is already there.
- The System ID was not loaded into the download server. This may be because the System ID was not on the effectivity list. It may also be that the site is not truly Broadband connected. Finally, the System ID may have been accidentally missed during download.
- If you think the System ID should have been downloaded to the server but it isn't, contact Gerry Bass or Abdallah Elsayed.