

GE Healthcare

LightSpeed 4.X

Illustrated Parts List

OPERATING DOCUMENTATION



2341961-100
Rev 8

The information contained in this manual applies to the following LightSpeed 4.X (LightSpeed16) catalog numbers:

LightSpeed 16 Cat Numbers

- B7868JA
- B7868JB
- B7868JC
- B7868JD
- B7868JE
- B7868JF

IMPORTANT PRECAUTIONS

LANGUAGE

WARNING

- THIS SERVICE MANUAL IS AVAILABLE IN ENGLISH ONLY.
- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR OR PATIENT FROM ELECTRIC SHOCK, MECHANICAL OR OTHER HAZARDS.

AVERTISSEMENT

- CE MANUEL DE MAINTENANCE N'EST DISPONIBLE QU'EN ANGLAIS.
- SI LE TECHNICIEN DU CLIENT A BESOIN DE CE MANUEL DANS UNE AUTRE LANGUE QUE L'ANGLAIS, C'EST AU CLIENT QU'IL INCOMBE DE LE FAIRE TRADUIRE.
- NE PAS TENTER D'INTERVENTION SUR LES ÉQUIPEMENTS TANT QUE LE MANUEL SERVICE N'A PAS ÉTÉ CONSULTÉ ET COMPRIS.
- LE NON-RESPECT DE CET AVERTISSEMENT PEUT ENTRAÎNER CHEZ LE TECHNICIEN, L'OPÉRATEUR OU LE PATIENT DES BLESSURES DUES À DES DANGERS ÉLECTRIQUES, MÉCANIQUES OU AUTRES.

WARNUNG

- DIESES KUNDENDIENST-HANDBUCH EXISTIERT NUR IN ENGLISCHER SPRACHE.
- FALLS EIN FREMDER KUNDENDIENST EINE ANDERE SPRACHE BENÖTIGT, IST ES AUFGABE DES KUNDEN FÜR EINE ENTSPRECHENDE ÜBERSETZUNG ZU SORGEN.
- VERSUCHEN SIE NICHT, DAS GERÄT ZU REPARIEREN, BEVOR DIESES KUNDENDIENST-HANDBUCH ZU RATE GEZOGEN UND VERSTANDEN WURDE.
- WIRD DIESE WARNUNG NICHT BEACHTET, SO KANN ES ZU VERLETZUNGEN DES KUNDENDIENSTTECHNIKERS, DES BEDIENERS ODER DES PATIENTEN DURCH ELEKTRISCHE SCHLÄGE, MECHANISCHE ODER SONSTIGE GEFAHREN KOMMEN.

AVISO

- ESTE MANUAL DE SERVICIO SÓLO EXISTE EN INGLÉS.
- SI ALGÚN PROVEEDOR DE SERVICIOS AJENO A GEMS SOLICITA UN IDIOMA QUE NO SEA EL INGLÉS, ES RESPONSABILIDAD DEL CLIENTE OFRECER UN SERVICIO DE TRADUCCIÓN.
- NO SE DEBERÁ DAR SERVICIO TÉCNICO AL EQUIPO, SIN HABER CONSULTADO Y COMPRENDIDO ESTE MANUAL DE SERVICIO.
- LA NO OBSERVANCIA DEL PRESENTE AVISO PUEDE DAR LUGAR A QUE EL PROVEEDOR DE SERVICIOS, EL OPERADOR O EL PACIENTE SUFRAN LESIONES PROVOCADAS POR CAUSAS ELÉCTRICAS, MECÁNICAS O DE OTRA NATURALEZA.

ATENÇÃO

- ESTE MANUAL DE ASSISTÊNCIA TÉCNICA SÓ SE ENCONTRA DISPONÍVEL EM INGLÊS.
- SE QUALQUER OUTRO SERVIÇO DE ASSISTÊNCIA TÉCNICA, QUE NÃO A GEMS, SOLICITAR ESTES MANUAIS NOUTRO IDIOMA, É DA RESPONSABILIDADE DO CLIENTE FORNECER OS SERVIÇOS DE TRADUÇÃO.
- NÃO TENTE REPARAR O EQUIPAMENTO SEM TER CONSULTADO E COMPREENDIDO ESTE MANUAL DE ASSISTÊNCIA TÉCNICA.
- O NÃO CUMPRIMENTO DESTE AVISO PODE POR EM PERIGO A SEGURANÇA DO TÉCNICO, OPERADOR OU PACIENTE DEVIDO A' CHOQUES ELÉTRICOS, MECÂNICOS OU OUTROS.

AVVERTENZA

- IL PRESENTE MANUALE DI MANUTENZIONE È DISPONIBILE SOLTANTO IN INGLESE.
- SE UN ADDETTO ALLA MANUTENZIONE ESTERNO ALLA GEMS RICHIEDE IL MANUALE IN UNA LINGUA DIVERSA, IL CLIENTE È TENUTO A PROVVEDERE DIRETTAMENTE ALLA TRADUZIONE.
- SI PROCEDA ALLA MANUTENZIONE DELL'APPARECCHIATURA SOLO DOPO AVER CONSULTATO IL PRESENTE MANUALE ED AVERNE COMPRESO IL CONTENUTO.
- NON TENERE CONTO DELLA PRESENTE AVVERTENZA POTREBBE FAR COMPIERE OPERAZIONI DA CUI DERIVINO LESIONI ALL'ADDETTO ALLA MANUTENZIONE, ALL'UTILIZZATORE ED AL PAZIENTE PER FOLGORAZIONE ELETTRICA, PER URTI MECCANICI OD ALTRI RISCHI.

警告

このサービスマニュアルには英語版しかありません。
GEMS以外でサービスを担当される業者が英語以外の言語を要求される場合、翻訳作業はその業者の責任で行うものとさせていただきます。
このサービスマニュアルを熟読し理解せずに、装置のサービスを行わないで下さい。
この警告に従わない場合、サービスを担当される方、操作員あるいは患者さんが、感電や機械的又はその他の危険により負傷する可能性があります。

注意:

本维修手册仅存有英文本。
非 GEMS 公司的维修员要求非英文本的维修手册时，
客户需自行负责翻译。
未详细阅读和完全了解本手册之前，不得进行维修。
忽略本注意事项会对维修员，操作员或病人造成触电，机械伤害或其他伤害。

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation "damage in shipment" written on all copies of the freight or express bill before delivery is accepted or "signed for" by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage MUST be reported to the carrier immediately upon discovery, or in any event, within 14 days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this 14 day period.

To file a report:

- Call 1-800-548-3366 and use option 8.
- Fill out a report on <http://us44hdd21/sctq/InstallFulfill/InstalFulfillment.htm>
- Contact your local service coordinator for more information on this process.

Rev. Jan. 5, 2005

CERTIFIED ELECTRICAL CONTRACTOR STATEMENT

All electrical Installations that are preliminary to positioning of the equipment at the site prepared for the equipment shall be performed by licensed electrical contractors. In addition, electrical feeds into the Power Distribution Unit shall be performed by licensed electrical contractors. Other connections between pieces of electrical equipment, calibrations and testing shall be performed by qualified GE Healthcare personnel. The products involved (and the accompanying electrical installations) are highly sophisticated, and special engineering competence is required. In performing all electrical work on these products, GE will use its own specially trained field engineers. All of GE's electrical work on these products will comply with the requirements of the applicable electrical codes.

The purchaser of GE equipment shall only utilize qualified personnel (i.e., GE's field engineers, personnel of third-party service companies with equivalent training, or licensed electricians) to perform electrical servicing on the equipment.

IMPORTANT...X-RAY PROTECTION

X-ray equipment if not properly used may cause injury. Accordingly, the instructions herein contained should be thoroughly read and understood by everyone who will use the equipment before you attempt to place this equipment in operation. The General Electric Company, Medical Systems Group, will be glad to assist and cooperate in placing this equipment in use.

Although this apparatus incorporates a high degree of protection against x-radiation other than the useful beam, no practical design of equipment can provide complete protection. Nor can any practical design compel the operator to take adequate precautions to prevent the possibility of any persons carelessly exposing themselves or others to radiation.

It is important that anyone having anything to do with x-radiation be properly trained and fully acquainted with the recommendations of the National Council on Radiation Protection and Measurements as published in NCRP Reports available from NCRP Publications, 7910 Woodmont Avenue, Room 1016, Bethesda, Maryland 20814, and of the International Commission on Radiation Protection, and take adequate steps to protect against injury.

The equipment is sold with the understanding that the General Electric Company, Medical Systems Group, its agents, and representatives have no responsibility for injury or damage which may result from improper use of the equipment.

Various protective materials and devices are available. It is urged that such materials or devices be used.

LITHIUM BATTERY CAUTIONARY STATEMENT

CAUTION



Risk of Explosion.

Danger of explosion if battery is incorrectly replaced.

Replace only with the same or equivalent type recommended by the manufacturer. Discard used batteries according to the manufacturer's instructions.

OMISSIONS & ERRORS

Customers, please contact your GE Sales or Service representatives.

GE personnel, please use the GE Healthcare PQR Process to report all omissions, errors, and defects in this publication.

Revision History

Revision	Date	Reason for change
8	2/21/07	Added NGPDU information (applies to some CT and PET/CT systems shipped in 2006 or later) Added Cross-References to Dir 5134893-100 (GOC4 IPL) Updated formatting to comply with new GE Healthcare standards
7	12/22/04	Chapter 4: Removed n/a Tube Plate Assembly Added Cross-References to Dir 2394950-100 (GOC3 IPL)
6	05/25/04	Chapter 4: - New P/N for Performix Tube. - Removed non-applicable Detector
5	04/06/04	Chapter 4: Corrected FRU Code for p/n 46-328417P16 in Table 4-47 (Slipring BOM 2333608)
4	01/19/04	CQA 13011416: Updated Table 4-47 (Slipring BOM 2333608) Chapter 4: Added Section 16.0 — DAS Fiber Optic Cables Chapter 6: Updated drawing 2341924
3	04/25/03	CQAs 1030542 & 1030544: Added AILs to Table 4-29 DAS Assembly, H16 (2332633 BOM, rev 1) Updated BOMs and drawings throughout.
2	03/13/03	CQA 1030541: Added LCD Display to Chapter 5. Added Chapter 6 Global Console - Linux (GOC2)
1	11/21/02	CQA 10210818: FRU Codes for Console Cables Chapter 4: - Updated drawing 2217930 - Added Section 15.0 — Gantry Service Items Chapter 5: - Updated various BOMs and drawings - Added Section 1.3 Octane2 FRUs
0	08/21/02	New Product Introduction (LightSpeed 16)

End of Revision History.

Table of Contents

Chapter 1

System Matrices & Manual Conventions	27
Section 1.0 — Page Layout, Including Finger Tabs	27
Section 2.0 — Manual Organization.....	28
Section 3.0 — System Matrices.....	28
3.1 LightSpeed 4.X System Matrix	29
LightSpeed 4.X System Catalog Tree	29
3.2 Table Sub-System Matrix.....	30
Table BOM Tree	30
3.3 Gantry Sub-System Matrices	31
Gantry BOM Tree (sheet 1 of 2)	31
Gantry BOM Tree (sheet 2 of 2)	32
3.4 Console Sub-System Matrix - GOC1 (Global Operator Console - Octane2)	33
GOC1 Console BOM Tree	33
3.5 Console Sub-System Matrix - GOC2 (Global Operator Console - Linux).....	34
GOC2 Console BOM Tree	34
3.6 GRE Xstream Console (GOC3 or GOC4).....	35

Chapter 2

Table	37
Table BOM Tree.....	38
Table w/ Footswitch (2276776-2 BOM, rev 2)	39
Section 1.0 — Table Assembly - LCC (2269165).....	40
Table Assembly - LCC (2269165, sheet 1, rev 3)	40
Table Assembly - LCC (2269165, sheet 2, rev 3)	41
Table Assembly - LCC (2269165, sheet 3, rev 0)	42
Table Assembly - LCC (2269165 BOM, rev 9)	43
1.1 Side Panel Kit (2269167-2).....	48
Side Panel Kit (2269167-2 BOM, rev 0)	48

1.2	Right Base Cover Assembly (2271244)	50
	Right Base Cover Assembly (2271244, sheet 1, rev 1)	50
	Right Base Cover Assembly (2271244 BOM, rev 3)	51
1.3	Left Base Cover Assembly (2271247)	52
	Left Base Cover Assembly (2271247, sheet 1, rev 1)	52
	Left Base Cover Assembly (2271247 BOM, rev 3)	53
1.4	Cover Table Top Left (2372371)	54
	Cover Table Top Left (2372371 sheet 1, rev 1)	54
	Cover Table Top Left (2372371 BOM, rev 1)	55
1.5	Cover Table Top Right (2372372)	56
	Cover Table Top Right (2372372 sheet 1, rev 1)	56
	Cover Table Top Left (2372372 BOM, rev 1)	57
1.6	Metalless Cradle (2115992-2)	58
	Metalless Cradle (2115992-2 BOM, rev 0)	58
1.7	Table Subassembly (2269166)	60
	Table Subassembly (2269166, sheet 1, rev 0)	60
	Table Subassembly (2269166, sheet 2, rev 0)	61
	Table Subassembly (2269166, sheet 3, rev 0)	62
	Table Subassembly (2269166, sheet 4, rev 0)	63
	Table Subassembly (2269166 BOM, rev 2)	64
	1.7.1 – Carriage Assembly (2269168)	66
	Carriage Assembly (2269168, sheet 1, rev 1)	66
	Carriage Assembly (2269168 BOM, rev 1)	67
	1.7.2 – Latch Assembly (2269177)	68
	Latch Assembly (2269177, sheet 1, rev 0)	68
	Latch Assembly (2269177 BOM, rev 0)	69
	1.7.3 – Cradle Drive Assembly (2269169)	70
	Cradle Drive Assembly (2269169 sheet 1, rev 0)	70
	Cradle Drive Assembly (2269169 sheet 1, rev 0)	71
	Cradle Drive Assembly (2269169 BOM, rev 2)	72
	1.7.4 – Longitudinal Encoder Asm (2269175)	76
	Longitudinal Encoder Asm (2269175, sheet 1, rev 0)	76
	Longitudinal Encoder Asm (2269175, sheet 2, rev 0)	77
	Longitudinal Encoder Asm (2269175 BOM, rev 1)	78

1.7.5 – Servo Assembly (46-297059G1)	80
Servo Assembly (46-297059, sheet 1, rev 2)	80
Servo Assembly (46-297059G1 BOM, rev 3)	81
1.7.6 – Power Supply Assembly (2269172)	82
Power Supply Assembly (2269172, sheet 1, rev 0)	82
Power Supply Assembly (2269172, sheet 2, rev 0)	83
Power Supply Assembly (2269172 BOM, rev 0)	84
1.7.7 – Elevation Encoder Assembly (2269173)	86
Elevation Encoder Assembly (2269173, sheet 1, rev 1)	86
Elevation Encoder Assembly (2269173, sheet 2, rev 1)	87
Elevation Encoder Assembly (2269173 BOM, rev 1)	88
1.7.8 – Frame Assembly (2269174)	90
Frame Assembly (2269174, sheet 1, rev 0)	90
Frame Assembly (2269174, sheet 2, rev 0)	91
Frame Assembly (2269174, sheet 3, rev 0)	92
Frame Assembly (2269174, sheet 4, rev 0)	93
Frame Assembly (2269174, sheet 5, rev 0)	94
Frame Assembly (2269174 BOM, rev 4)	95
Section 2.0 — Foot Switch Assembly (2330508)	100
Foot Switch Assembly (2330508, sheet 1, rev 2)	100
Footswitch Assembly (2330508 BOM, rev 5)	101
 Chapter 3	
Power Distribution Unit (PDU)	103
Section 1.0 — CPDU (2269902)	104
Compact PDU (2269902, front view)	104
Compact PDU (2269902, rear view)	105
Compact PDU (2269902)	106
Section 2.0 — NGPDU (2326492-2)	108
NGPDU (2326492-2, sheet 1, rev 0)	108
NGPDU - 90kVA (2326492-2 BOM, Rev 3)	109
NGPDU Assembly2 (2348075, sheet 1, rev 0)	110
NGPDU Assembly2 (2348075, sheet 2, rev 0)	111
NGPDU Assembly2 (2348075, sheet 3, rev 0)	112

NGPDU ASSEMBLY2 (2348075 BOM, Rev 18)	113
---	-----

Chapter 4

Gantry..... 125

Gantry BOM Tree (sheet 1 of 2)	125
Gantry BOM Tree (sheet 2 of 2)	126
Gantry (2339985, sheet 1, rev 1)	127
Gantry (2339985, sheet 2, rev 0)	128
Gantry (2339985, sheet 3, rev 0)	129
Gantry (2339985, sheet 4, rev 1)	130
Gantry (2339985, sheet 5, rev 0)	131
Gantry (2339985, sheet 6, rev 0)	132
Gantry (2339985 BOM, rev 1)	133

Section 1.0 — Stationary Base Assembly (2260616)..... 136

Stationary Base Assembly (2260616, sheet 1, rev 6)	136
Stationary Base Assembly (2260616, sheet 2, rev 6)	137
Stationary Base Assembly (2260616, sheet 3, rev 5)	138
Stationary Base Assembly (2260616 BOM, rev 10)	139

1.1 STC Assembly (2266927)..... 142

STC Assembly (2266927, sheet 1, rev 3)	142
STC Assembly (2266927 BOM, rev 6)	143

1.2 Hydraulic Tilt Motor Assembly (2219018) 146

Hydraulic Tilt Motor Assembly (2219018, sheet 1, rev 2).....	146
Hydraulic Tilt Motor Assembly (2219018, sheet 2, rev 2).....	147
Hydraulic Tilt Motor Assembly (2219018, sheet 3, rev 2).....	148
Hydraulic Tilt Motor Assembly (2219018 BOM, rev 4)	149

1.2.1 – Hydraulic Tilt Parts (2225778) 150

Collector Hydraulic Tilt Parts (2225778 BOM, rev 0).....	150
1.2.1.1 – Hydraulic Pump (2262276).....	152
Hydraulic Pump (2262276, sheet 1, rev 0)	152
Tilt (2262276 BOM, rev 1)	153

1.3 E-Stop Assembly (2266187) 154

E-Stop Assembly (2266187, sheet 1, rev 0).....	154
E-Stop Assembly (2266187 BOM, rev 2)	155

1.4	Power Pan (2275608-3)	156
	Power Pan (2275608, sheet 1, rev 4).....	156
	Power Pan (2275608, sheet 2, rev 4).....	157
	Power Pan (2275608-3 BOM, rev 7)	158
	Section 2.0 — Tilting Frame (2339984)	160
	Tilting Frame (2339984, sheet 1, rev 1)	160
	Tilting Frame (2339984, sheet 2, rev 0)	161
	Tilting Frame (2339984 BOM, rev 2)	162
2.1	Indexer/Rotating Lock Assembly (2181221)	166
	Indexer/Rotating Lock Assembly (2181221, sheet 1, rev 0).....	166
	Indexer/Rotating Lock Assembly (2181221 BOM, rev 0)	167
2.2	Home Flag Assembly (2227378)	168
	Home Flag Assembly (2227378, sheet 1, rev 0)	168
	Home Flag (2227378 BOM, rev 1)	169
2.3	Gantry Fan Assembly - RH (2315618)	170
	Gantry Fan Assembly - RH (2315618, sheet 1, rev 1)	170
	Fan Assembly - Gantry, H3 RH (2315618 BOM, rev 1)	171
2.4	Motor Mount/Belt Tensioner (2217928-2)	172
	Motor Mount/Belt Tensioner (2217928, sheet 1, rev 1).....	172
	Motor Mount/Belt Tensioner (2217928-2 BOM, rev 3)	173
2.5	Axial Dynamic Braking Module (2266934-2)	174
	Axial Dynamic Braking Module (2266934, sheet 1, rev 3)	174
	Axial Dynamic Braking Module (2266934, sheet 2, rev 1)	175
	Axial Dynamic Brake Module, w/cvr mnt hls (2266934-2 BOM, rev 3).....	176
2.6	Gantry Fan Assembly - LH (2315619)	178
	Gantry Fan Assembly - LH (2315619, sheet 1, rev 2).....	178
	Fan Assembly - Gantry, H3 LH (2315619 BOM, rev 1).....	179
2.7	Drive Assembly, Axial (2217930)	180
	Drive Assembly, Axial (2217930, sheet 1, rev 3)	180
	Drive Assembly, Axial (2217930, sheet 2, rev 2)	181
	Drive Assembly, Axial (2217930 BOM, rev 6)	182
	2.7.1 – Axial Drive Motor Brake Assembly (2277225)	184
	Axial Drive Motor Brake Assembly (2277225, sheet 1, rev 0).....	184
	Brake Assembly, Axial Drive Motor (2277225 BOM, rev 2).....	185

2.7.2 – Axial Drive Motor Controller w/Corers Assembly (2233402)	186
Axial Drive Motor Controller w/Corers Assembly (2233402, sheet 1, rev 1)	186
Axial Drive Motor Controller w/Corers Assembly (2233402, sheet 2, rev 1)	187
Axial Drive Motor Controller w/Corers Assembly (2233402 BOM, rev 2)	188
2.8 Encoder Mount Assembly (2175817)	190
Encoder Mount Assembly (2175817, sheet 1, rev 1)	190
Encoder Mount Assembly (2175817 BOM, rev 2)	191
2.9 Brace Assembly, Upper (2273029)	192
Brace Assembly, Upper (2273029, sheet 1, rev 1)	192
Brace, Gantry Cover Upper (2273029 BOM, rev 0)	193
2.10 Brace Assembly, Lower (2273030)	194
Brace Assembly, Lower (2273030, sheet 1, rev 1)	194
Brace, Gantry Cover Lower (2273030 BOM, rev 0)	195
Section 3.0 — Rotating Base Assembly (2334114)	196
Rotating Base Assembly (2334114, sheet 1, rev 1)	196
Rotating Base Assembly (2334114, sheet 2, rev 0)	197
Rotating Base Assembly (2334114, sheet 3, rev 0)	198
Rotating Base Assembly (2334114, sheet 4, rev 0)	199
Rotating Base Assembly (2334114, sheet 5, rev 0)	200
Rotating Base Assembly (2334114, sheet 6, rev 1)	201
Rotating Base Assembly (2334114 BOM, rev 3)	202
3.1 Collimator w/Plate (2214768-3)	206
Collimator w/Plate (2214768, sheet 3, rev 3)	206
Collimator w/Plate (2214768-3 BOM, rev 0)	207
3.1.1 – Collimator (2120100-3)	208
Collimator (2120100, sheet 5, rev 8)	208
Collimator (2120100, sheet 6, rev 7)	209
Collimator (2120100-3 BOM, rev 5)	210
3.2 OBC Assembly (2231831)	212
OBC Assembly (2231831, sheet 1, rev 1)	212
OBC Assembly (2231831 BOM, rev 8)	213
3.2.1 – OBC Chassis Assembly (2232299)	216
OBC Chassis Assembly (2232299, sheet 1, rev 3)	216
OBC Chassis Assembly (2232299 BOM, rev 4)	217

3.3	HEMRC Assembly (2208043)	220
	HEMRC Assembly (2208043, sheet 1, rev 4)	220
	HEMRC Assembly (2208043, sheet 2, rev 1)	221
	HEMRC Assembly (2208043, sheet 3, rev 1)	222
	HEMRC Assembly (2208043 BOM, rev 5)	223
3.4	Detector Plate Assembly (2332635)	226
	Detector Plate Assembly (2332635, sheet 1, rev 2)	226
	Detector Plate Assembly (2332635 BOM, rev 2)	227
	3.4.1 – DAS Assembly (2332633)	228
	DAS Assembly (2332633, sheet 1, rev 2)	228
	DAS Assembly, H16 (2332633 BOM, rev 1)	229
	<i>3.4.1.1 – Backplane/Chassis, Left Assembly (2332632)</i>	<i>230</i>
	Backplane/ Chassis Assembly, Left (2332632, sheet 1, rev 0)	230
	Backplane/Chassis, Left Assembly (2332632 BOM, rev 0).....	231
	<i>3.4.1.2 – Backplane/Chassis, Center Assembly (2332631)</i>	<i>232</i>
	Backplane/Chassis, Center Assembly (2332631, sheet 1, rev 0)	232
	Backplane/Chassis, Center Assembly (2332631 BOM, rev 0)	233
	<i>3.4.1.3 – Backplane/Chassis, Right Assembly (2332630)</i>	<i>234</i>
	Backplane/Chassis, Right Assembly (2332630, sheet 1, rev 0).....	234
	Backplane/Chassis, Right Assembly (2332630 BOM, rev 0)	235
	3.4.2 – Detector Plate Sub-Assembly (2333652)	236
	Detector Plate Sub-Assembly (2333652, sheet 1, rev 1)	236
	Detector Plate Sub-Assembly (2333652 BOM, rev 1)	237
	3.4.3 – Detector Assembly (2323504)	238
	Detector Assembly (2323504, sheet 1, rev 0)	238
	Watson16 Detector Assembly (2323504 BOM, rev 1)	239
	<i>3.4.3.1 – Detector Memory Board Assembly (2296076)</i>	<i>242</i>
	Detector Memory Board Assembly (2296076, sheet 1, rev 1)	242
	Detector Memory Board Assembly (2296076 BOM, rev 3)	243
	3.4.4 – DAS Duct Fan Assembly (2333651)	244
	DAS Duct Fan Assembly (2333651, sheet 1, rev 1).....	244
	DAS Duct Fan Assembly (2333651 BOM, rev 1)	245
3.5	Power Supply, 24VDC (2334457)	246
	Power Supply Assembly, 24VDC (2334457, sheet 1, rev 0).....	246
	Power Supply Assembly, 24VDC (2334457 BOM, rev 0)	247

3.6	Fuse Box Assembly - H16 Upgraded (2331963-3)	248
	Fuse Box Assembly - H16 Upgraded (2331963-2, sheet 1, rev 1)	248
	Fuse Box Assembly - H16 Upgraded (2331963-2, sheet 2, rev 0)	249
	Fuse Box Assembly - H16 Upgraded (2331963-3 BOM, rev 0)	250
3.6.1	Harness (2347052)	252
	Harness (2347052, sheet 1, rev 0)	252
	Harness (2347052 BOM, rev 0)	253
3.7	Power Supply, 12VDC (2334459)	254
	Power Supply, 12VDC (2334459, sheet 1, rev 0)	254
	Power Supply, 12VDC (2334459 BOM, rev 0)	255
3.7.1	Power Supply Assembly, 12V DC (2334456)	256
	Power Supply Assembly, 12V DC (2334456, sheet 1, rev 0)	256
	Power Supply Assembly (2334456 BOM, rev 0)	257
3.8	DHCB Assembly (2294856)	258
	DHCB Assembly (2294856, sheet 1, rev 1)	258
	DHCB Assembly (2294856 BOM, rev 1)	259
3.9	Anode Inverter Harness (2240779)	260
	Anode Inverter Harness (2240779, sheet 1, rev 1)	260
	Anode Inverter Harness (2240779 BOM, rev 1)	261
3.10	Cathode Inv/OBC Harness (2253952)	262
	Cathode Inv/OBC Harness (2253952, sheet 1, rev 1)	262
	Cathode Inv/OBC Harness (2253952 BOM, rev 3)	263
3.11	Power Supply Assembly, 5VDC (2334455)	264
	Power Supply Assembly, 5VDC (2334455, sheet 1, rev 0)	264
	Power Supply Assembly, 5VDC (2334455 BOM, rev 0)	265
3.12	Power Supply, 5V Logic (2334458)	266
	Power Supply, 5V Logic (2334458, sheet 1, rev 0)	266
	Power Supply, 5V Logic (2334458 BOM, rev 0)	267
Section 4.0	Slipring Assembly (2333608)	268
	Slipring Assembly (2333608, sheet 1, rev 1)	268
	Slipring Assembly (2333608, sheet 2, rev 1)	269
	Slipring Assembly (2333608 BOM, rev 2)	270
4.1	850MB Slipring Assembly (2304038)	272
	850MB Slipring Assembly (2304038, sheet 1, rev 0)	272

Slipring Assembly (2304038 BOM, rev 0)	273
4.1.1 – 850MB Receiver/Antenna Assembly (2333609)	274
850MB Receiver/Antenna Assembly (2333609, sheet 1, rev 0).....	274
Receiver-Antenna Assembly (2333609 BOM, rev 0)	275
4.1.2 – Brushblock Assembly (2240267)	276
Brushblock Assembly (2240267, sheet 1, rev 1).....	276
Brushblock Assembly (2240267 BOM, rev 1)	277
4.1.3 – Slipring Assembly (2333614).....	278
Slipring Assembly (2333614, sheet 1, rev 0).....	278
Slipring Assembly (2333614 BOM, rev 0)	279
Section 5.0 — S/R Cover Collector (2282516)	280
S/R Cover Collector (2282516 BOM, rev 0)	280
Section 6.0 — Scan Window Assembly (2101967-4)	281
Scan Window Assembly - Gantry (2101967-4 BOM, rev 0)	281
Section 7.0 — Cover Assembly - Front (2272939)	282
Cover Assembly - Front Gantry (2272939, sheet 1, rev 2).....	282
Cover Assembly - Front Gantry (2272939 BOM, rev 4)	283
7.1 Cover Sub Assembly, Front (2273017).....	284
Cover Sub Assembly, Gantry Front (2273017, sheet 1, rev 6)	284
Cover Sub Assembly, Gantry Front (2273017, sheet 2, rev 5)	285
Cover Sub Assembly, Gantry Front (2273017 BOM, rev 7)	286
7.2 Mic Assembly (2111955)	288
Mic Assembly (2111955 BOM, rev 0).....	288
7.3 Patient Sensor (2277405).....	290
Front Sensor Tilt Assembly (2277405, sheet 1, rev 0)	290
Patient Sensor, Small (2277405 BOM, rev 0)	291
Section 8.0 — Cover Assembly - Rear (2273016)	292
Cover Assembly, Rear Gantry (2273016 sheet 1, rev 3)	292
Cover Assembly, Rear Gantry (2273016 BOM, rev 3)	293
8.1 Cover Sub Assembly, Rear (2273018)	294
Cover Sub Assembly, Rear Gantry (2273018, sheet 1, rev 5)	294
Cover Sub Assembly, Rear Gantry (2273018, sheet 2, rev 3)	295
Cover Sub Assembly, Rear Gantry (2273018 BOM, rev 6).....	296

8.2	Cover Sub-Assy, Rear Gantry - FRP (2291369)	298
	Cover Sub-Assy, Rear Gantry - FRP (2291369 sheet 1, rev 0)	298
	Cover Sub-Assy, Rear Gantry - FRP (2291369 sheet 2, rev 0)	299
	Cover Sub-Assy, Rear Gantry - FRP (2291369 BOM, rev 3)	300
8.3	Patient Sensor, Large (2277404)	302
	Patient Sensor, Large (2277404 sheet 1, rev 0)	302
	Patient Sensor, Large (2277404 BOM, rev 0)	303
	Section 9.0 — Gantry Top Covers Kit - RIM (2292881)	304
	Gantry Top Covers Kit - RIM (2292881, sheet 1, rev 1)	304
	Gantry Top Covers Kit - RIM (2292881 BOM, rev 2)	305
9.1	Cover Assembly, Top RH (2273013)	306
	Cover Assembly, Top RH (2273013, sheet 1, rev 3)	306
	Cover Assembly, Top RH (2273013, sheet 2, rev 3)	307
	Cover Assembly, Gantry Top RH (2273013 BOM, rev 5)	308
9.2	Cover Assembly, Top LH (2273014)	310
	Cover Assembly, Gantry Top LH (2273014, sheet 1, rev 3)	310
	Cover Assembly, Gantry Top LH (2273014, sheet 2, rev 3)	311
	Cover Assembly, Gantry Top LH (2273014 BOM, rev 4)	312
9.3	Cover Assembly, Right Side (2276690)	314
	Cover Assembly, Right Side (2276690, sheet 1, rev 2)	314
	Cover Assembly, Side Gantry Right (2276690 BOM, rev 3)	315
9.4	Cover Assembly, Left Side (2276691)	316
	Cover Assembly, Left Side (2276691, sheet 1, rev 2)	316
	Cover Assembly, Side Gantry Left (2276691 BOM, rev 2)	317
	Section 10.0 — Gantry Top Covers Kit - FRP (2292882)	318
	Gantry Top Covers Kit - FRP (2292882, sheet 1, rev 0)	318
	Gantry Top Covers Kit - FRP (2292882 BOM, rev 1)	319
10.1	Cover Assembly, Top RH - FRP (2290538)	320
	Top Cover Assembly, RH - FRP (2290538, sheet 1, rev 0)	320
	Cover Assembly, Top RH - FRP (2290538 BOM, rev 1)	321
10.2	Cover Assembly, Top LH - FRP (2290539)	322
	Top Cover Assembly, LH - FRP (2290539, sheet 1, rev 0)	322
	Cover Assembly, Top LH - FRP (2290539 BOM, rev 0)	323

10.3	Cover Assembly, Side RH - FRP (2290540)	324
	Side Cover Assembly, RH - FRP (2290540, sheet 1, rev 0)	324
	Cover Assembly, Side RH - FRP (2290540 BOM, rev 0)	325
10.4	Cover Assembly, Side LH - FRP (2290541)	326
	Side Cover Assembly, LH - FRP (2290541, sheet 1, rev 0)	326
	Cover Assembly, Side LH - FRP (2290541 BOM, rev 1).....	327
	Section 11.0 — Tilt Pot Assembly (2178677)	328
	Tilt Pot Assembly (2178677, sheet 1, rev 1)	328
	Tilt Pot Assembly (2178677 BOM, rev 4)	329
11.1	Tilt Pot Harness Assembly (46-297036G1)	330
	Tilt Pot Harness Assembly (46-297036, sheet 1, rev 1)	330
	Tilt Pot Harness Assembly (46-297036G1 BOM, rev 0).....	331
11.2	Tilt Switch Assembly Harness (2277131)	332
	Switch Harness Assembly (2277131, sheet 1, rev 1).....	332
	Switch Harness Assembly (2277131, sheet 2, rev 1).....	333
	Switch Harness Assembly (2277131 BOM, rev 1)	334
	Section 12.0 — Base Cover Collector (2290278, 2290278-2)	336
	Gantry Base Covers (2290278, sheet 1, rev 2).....	336
	Base Covers, Gantry (2290278, -2 BOM, rev 1)	337
12.1	Cover Assembly, Base Rear - LH (2275686, 2275686-2)	338
	Cover Assembly, Base Rear - LH (2275686, sheet 1, rev 2)	338
	Cover Assembly, Base Rear - LH (2275686 BOM, rev 2).....	339
	Cover Assembly, Base Rear - LH (FRP) (2275686-2 BOM, rev 0)	339
12.2	Cover Assembly, Base Rear - RH (2275685, 2275685-2)	340
	Cover Assembly, Base Rear - RH (2275685, sheet 1, rev 2).....	340
	Cover Assembly, Base Rear - RH (2275685 BOM, rev 2)	341
	Cover Assembly, Base Rear - RH (2275685-2 BOM, rev 0)	341
	Section 13.0 — Cover Dollies	342
	Cover Dollies (rev 03/19/02).....	342
	Section 14.0 — Phantoms (2180544)	342
	Phantoms (2180544 BOM, rev 0).....	342
	Section 15.0 — Gantry Service Items	343
	Gantry Service Items	343

Section 16.0 — DAS Fiber Optic Cables.....	343
DAS Fiber Optic Cables	343

Chapter 5

Global Console - Octane (GOC1).....	345
H16 Console (GOC1), BOM Tree	345
Section 1.0 — Console (2341923, 2341923-2)	346
Console (2341923, sheet 1, rev 1)	346
Console (2341923, sheet 2, rev 0)	347
Console (2341923, sheet 3, rev 1)	348
Console (2341923, sheet 4, rev 0)	349
Console Assembly (2341923, 2341923-2 BOMs, rev 6)	350
1.1 Console Table Top Collector (2330564)	352
Console Table Top Collector (2330564 BOM, rev 0)	352
1.1.1 – Table Top Assembly, Monitor (2315657).....	354
Table Top Assembly, Monitor (2315657, sheet 1, rev 0)	354
Table Top Assembly, Monitor (2315657 BOM, rev 0)	355
1.1.2 – Table Top Assembly, Keyboard (2330566)	356
Table Top Assembly, Keyboard (2330566, sheet 1, rev 0)	356
Table Top Assembly, Keyboard (2330566 BOM, rev 0).....	357
1.2 Host Computer - Octane2 (2352078, 2352078-2)	358
Octane2 Computer Assembly (2295618, sheet 1, rev 0)	358
Octane2 18/73 Computer Assembly (2352078 BOM, rev 0).....	359
Octane2 18/73 Computer Assembly, w/D3D (2352078-2 BOM, rev 0)	359
1.3 Octane2 FRUs.....	360
Octane2 FRUs, Rev 2002-11-20	360
1.4 Recon Box Assembly (2341924)	362
Recon Box Assembly (2341924, sheet 1, rev 2)	362
Recon Box Assembly (2341924, sheet 2, rev 0)	363
Recon Box Assembly (2341924, sheet 3, rev 0)	364
Recon Box Assembly (2341924, sheet 4, rev 1)	365
Recon Box Assembly, H16 (2341924 BOM, rev 4)	366
1.4.1 – Chassis, VME Recon Box (2333179).....	368
VME Recon Box Chassis (2333179 sheet 1, rev 1)	368

	VME Recon Box Chassis (2333179 sheet 2, rev 1)	369
	VME Recon Box Chassis (2333179 sheet 3, rev 1)	370
	Chassis, VME Recon Box (2333179 BOM, rev 4)	371
	1.4.1.1 – Fan Assembly, VME Chassis (2324186).....	374
	Fan Assembly, VME Chassis (2324186, sheet 1, rev 0)	374
	Fan Assembly, VME Chassis (2324186 BOM, Rev 0)	375
	1.4.1.2 – VME Console Power Supply (2260421-2).....	376
	VME Console Power Supply (2260421, sheet 1, rev 1)	376
	VME Console Power Supply (2260421-2 BOM, rev 1)	377
	1.4.2 – VME Chassis Assembly (2307794)	378
	VME Chassis Assembly (2307794, sheet 1, rev 2)	378
	VME Chassis Assembly (2307794, sheet 2, rev 2)	379
	VME Chassis Assembly (2307794, sheet 3, rev 0)	380
	VME Chassis Assembly (2307794 BOM, rev 0)	381
1.5	Console Chassis Asm, Octane2 (2341851)	384
	Console Chassis Assembly, Octane2 (2341851 sheet 1, rev 0)	384
	Console Chassis Assembly, Octane2 (2341851 sheet 2, rev 0)	385
	Console Chassis Assembly, Octane2 (2341851 BOM, rev 0)	386
	1.5.1 – Outlet Box Assembly (2341413)	388
	Outlet Box Assembly (2341413, sheet 1, rev 1)	388
	Outlet Box Assembly (2341413 BOM, rev 0)	389
	<i>Breaker Assembly (2341432)</i>	<i>390</i>
	Breaker Assembly (2341432, sheet 1, rev 0)	390
	Breaker Assembly (2341432 BOM, rev 1)	391
	1.5.2 – Fan Assembly, Global Console Chassis (2323943)	392
	Fan Assembly, Console Chassis (2323943, sheet 1, rev 0)	392
	Fan Assembly, Global Console Chassis (2323943 BOM, rev 0)	393
1.6	Console Chassis Assembly (2341103)	394
	Console Chassis Assembly (2341103, sheet 1, rev 1)	394
	Console Chassis Assembly (2341103, sheet 2, rev 1)	395
	Chassis Assembly (2341103 BOM, rev 1)	396
1.7	SCSI Tower Assembly (2351055)	398
	SCSI Tower Assembly (2351055, sheet 1, rev 0)	398
	SCSI Tower Assembly (2351055 BOM, rev 2)	399

1.8	Global Modem Kit (2245794)	400
	Global Modem Kit (2245794 BOM, rev 3)	400
1.9	Side Covers (2337534)	400
	Side Covers (2337534 BOM, rev 0)	400
1.10	Seismic Kit (2336557)	400
	Seismic Kit (2336557 BOM, rev 1)	400
Section 2.0 — Console Display Monitors		402
2.1	21" Color CRT Monitors (B7850MM)	402
	Console Display Monitors (B7850MM BOM, rev 1).....	402
2.2	LCD Monitor & Cable (2309592)	402
	LCD Monitor & Cable (2309592 BOM, rev 1)	402
Section 3.0 — DASMs		403
	DASMs	403
Section 4.0 — Service Kits		404
	GOC Service Kits	404
4.1	SCSI Cable Kit (2333829)	405
	SCSI Cable Kit (2333829 BOM, rev 2)	405
4.2	Communication Cable Kit (2333831)	405
	Communication Cable Kit (2333831 BOM, rev 0)	405
4.3	Video/Power Cables (2333832)	406
	Video/Power Cables (2333832 BOM, rev 0)	406
4.4	Electrical Components (2333833)	406
	Electrical Components (2333833 BOM, rev 1)	406
4.5	Brackets & Straps (2333834)	406
	Brackets & Straps (2333834 BOM, rev 0)	406
4.6	Miscellaneous PM Hardware (2333835)	407
	Miscellaneous PM Hardware (2333835 BOM, rev 0)	407
4.7	Hardware Kit (2333836)	407
	Hardware Kit (2333836 BOM, rev 0)	407
4.8	Power Switch & Cables (2333837)	409
	Power Switch & Cables (2333837 BOM, rev 0)	409

4.9	Slide Service Kit (2336912)	409
	Slide Service Kit (2336912 BOM, rev 0)	409
	Section 5.0 — Keyboard w/SCIM, Collectors	410
	Keyboard & SCIM Assembly (Rev 11/15/00 - English shown)	410
5.1	English (2275758)	411
	Keyboard w/SCIM Collector - English (2275758 BOM, rev 1)	411
	<i>SCIM Collector – English (2275752)</i>	411
	SCIM Collector - English (2275752 BOM, rev 0)	411
5.2	French (2275758-2)	412
	Keyboard w/SCIM Collector - French (2275758-2 BOM, rev 1)	412
	<i>SCIM Collector - French (2275752-2)</i>	412
	SCIM Collector - French (2275752-2 BOM, rev 0)	412
5.3	German (2275758-3)	413
	Keyboard w/SCIM Collector - German (2275758-3 BOM, rev 1)	413
	<i>SCIM Collector - German (2275752-3)</i>	413
	SCIM Collector - German (2275752-3 BOM, rev 0)	413
5.4	Swedish (2275758-4)	414
	Keyboard w/SCIM Collector - Swedish (2275758-4 BOM, rev 1)	414
5.5	International (2275758-5)	415
	Keyboard w/SCIM Collector - Swedish (2275758-5 BOM, rev 1)	415
	<i>SCIM Collector - International (2275752-4)</i>	415
	SCIM Collector - International (2275752-4 BOM, rev 0)	415

Chapter 6

Global Console - Linux (GOC2)	417
Global Console - Linux (GOC2) BOM Tree	417
Global Console Assembly (2366438, sheet 1, rev 0)	418
Global Console Collector (2366438 BOM, rev 0)	419
Section 1.0 — Console Assembly (2350046)	420
Console Assembly (2350046, sheet 1, rev 1)	420
Console Assembly (2350046, sheet 2, rev 0)	421
Console Assembly (2350046, sheet 3, rev 1)	422
Console Assembly (2350046, sheet 4, rev 1)	423

	Console Assembly (2350046 BOM, rev 3)	424
1.1	Console Table Top Collector (2330564)	426
	Console Table Top Collector (2330564 BOM, rev 0)	426
	1.1.1 – Table Top Assembly, Monitor (2315657).....	428
	Table Top Assembly, Monitor (2315657, sheet 1, rev 0)	428
	Table Top Assembly, Monitor (2315657 BOM, rev 0)	429
	1.1.2 – Table Top Assembly, Keyboard (2330566)	430
	Table Top Assembly, Keyboard (2330566, sheet 1, rev 0)	430
	Table Top Assembly, Keyboard (2330566 BOM, rev 0).....	431
1.2	Recon Box Assembly (2341924)	432
	Recon Box Assembly (2341924, sheet 1, rev 3)	432
	Recon Box Assembly (2341924, sheet 2, rev 0)	433
	Recon Box Assembly (2341924, sheet 3, rev 3)	434
	Recon Box Assembly (2341924, sheet 4, rev 1)	435
	Recon Box Assembly (2341924 BOM, rev 4).....	436
	1.2.1 – Chassis, VME Recon Box (2333179).....	438
	Chassis, VME Recon Box (2333179 sheet 1, rev 1)	438
	Chassis, VME Recon Box (2333179 sheet 2, rev 1)	439
	Chassis, VME Recon Box (2333179 sheet 3, rev 1)	440
	Chassis, VME Recon Box (2333179 BOM, rev 4).....	441
	1.2.1.1 – Fan Assembly, VME Chassis (2324186).....	444
	Fan Assembly, VME Chassis (2324186, sheet 1, rev 0)	444
	Fan Assembly, VME Chassis (2324186 BOM, Rev 0)	445
	1.2.1.2 – VME Console Power Supply (2260421-2).....	446
	VME Console Power Supply (2260421, sheet 1, rev 1)	446
	VME Console Power Supply (2260421-2 BOM, rev 1)	447
	1.2.2 – VME Chassis Assembly (2307794)	448
	VME Chassis Assembly (2307794, sheet 1, rev 2)	448
	VME Chassis Assembly (2307794, sheet 2, rev 2)	449
	VME Chassis Assembly (2307794, sheet 3, rev 0)	450
	VME Chassis Assembly (2307794 BOM, rev 0).....	451
1.3	Console Chassis Assembly (2333194)	454
	Console Chassis Assembly (2333194 sheet 1, rev 1).....	454
	Console Chassis Assembly (2333194 sheet 2, rev 1).....	455
	Chassis Assembly (2333194 BOM, rev 2)	456

1.3.1 – Outlet Box Assembly (2174159-2)	458
Outlet Box Assembly (2174159-2, sheet 1, rev 1).....	458
Outlet Box Assembly (2174159-2 BOM, rev 0)	459
Power Switch Box (2190750-2)	460
Power Switch Box (2190750-2, sheet 1, rev 1)	460
Power Switch Box (2190750-2 BOM, rev 2).....	461
1.3.2 – Fan Assembly, Global Console Chassis (2323943)	462
Fan Assembly, Console Chassis (2323943, sheet 1, rev 0)	462
Fan Assembly, Global Console Chassis (2323943 BOM, rev 0).....	463
1.4 Console Chassis Assembly (2317687)	464
Console Chassis Assembly (2317687, sheet 1, rev 1).....	464
Console Chassis Assembly (2317687, sheet 2, rev 1).....	465
Console Chassis Assembly (2317687 BOM, rev 1)	466
1.5 SCSI Tower Assembly (2368838)	468
SCSI Tower Assembly (2368838, sheet 1, rev 0)	468
SCSI Tower, DVD-RAM, MOD (2368838 BOM, rev 0)	469
1.6 Global Modem Kit (2245794)	470
Global Modem Kit (2245794 BOM, rev 3)	470
Section 2.0 — Side Covers (2337534)	471
Side Covers (2337534 BOM, rev 0)	471
Section 3.0 — Miscellaneous (Associated Items)	472
Miscellaneous Items (2366438 AILs, rev 2003-04-25)	472
3.1 Seismic Kit (2336557)	473
Seismic Kit (2336557 BOM, rev 1)	473
Section 4.0 — Service Kits (Associated Items)	474
Service Kits (2350046 AILs, rev. 2003-04-25)	474
4.1 SCSI Cable Kit (2333829)	475
SCSI Cable Kit (2333829 BOM, rev 2).....	475
4.2 Ethernet / Communication Cable Kit (2333831)	475
Communication Cable Kit (2333831 BOM, rev 0)	475
4.3 Video / Power Cables (2333832)	476
Video/Power Cables (2333832 BOM, rev 0)	476

4.4	Electrical Components Kit (2333833)	476
	Electrical Components (2333833 BOM, rev 1)	476
4.5	Brackets and Straps (2333834)	476
	Brackets & Straps (2333834 BOM, rev 0)	476
4.6	Miscellaneous PM Hardware (2333835)	477
	Miscellaneous PM Hardware (2333835 BOM, rev 0)	477
4.7	Console Hardware Kit (2333836)	477
	Hardware Kit (2333836 BOM, rev 0)	477
4.8	Power Switch & Cables (2333837)	479
	Power Switch & Cables (2333837 BOM, rev 0)	479
4.9	Slide Service Kit (2336912)	479
	Slide Service Kit (2336912 BOM, rev 0).....	479
	Section 5.0 — Keyboard w/SCIM, Collectors	480
	Keyboard & SCIM Assembly (Rev 11/15/00)	480
5.1	English (2275758)	481
	Keyboard w/SCIM Collector - English (2275758 BOM, rev 1).....	481
	<i>SCIM Collector – English (2275752)</i>	481
	SCIM Collector - English (2275752 BOM, rev 0).....	481
5.2	French (2275758-2)	482
	Keyboard w/SCIM Collector - French (2275758-2 BOM, rev 1)	482
	<i>SCIM Collector - French (2275752-2)</i>	482
	SCIM Collector - French (2275752-2 BOM, rev 0)	482
5.3	German (2275758-3)	483
	Keyboard w/SCIM Collector - German (2275758-3 BOM, rev 1)	483
	<i>SCIM Collector - German (2275752-3)</i>	483
	SCIM Collector - German (2275752-3 BOM, rev 0)	483
5.4	Swedish (2275758-4)	484
	Keyboard w/SCIM Collector - Swedish (2275758-4 BOM, rev 1).....	484
5.5	International (2275758-5)	485
	Keyboard w/SCIM Collector - Swedish (2275758-5 BOM, rev 1).....	485
	<i>SCIM Collector - International (2275752-4)</i>	485
	SCIM Collector - International (2275752-4 BOM, rev 0)	485

Chapter 1

System Matrices & Manual Conventions

Section 1.0 — Page Layout, Including Finger Tabs

Headers and footers in this publication are designed to allow you to quickly identify your location. The document's part number, revision number and title appear in the header on every page. All page footers display the current page number, at left (when viewed landscape). Odd numbered page footers include a reversed "finger tab" to help the user quickly flip through this manual. These tabs include chapter number and name, plus the name of the current major (level one) section heading. Where applicable, even page footers will show the current minor (level two) section number and name (e.g., 3.3 Gantry Sub-System Matrices) - otherwise, they repeat the current major (level one) section heading and name (e.g., Section 1.0 — Page Layout, Including Finger Tabs).

Section 2.0 — Manual Organization

This document is organized according to product structure. That is, the order in which the contents appear is based on the BOMs (Bills of Material) used herein. This means that the user may be required to “dig down” through the product structure in order to find the part they are looking for. [For pdf users, hyperlinks are provided for quick navigation downward (the “back” arrow may be used to move back up)].

For example:

Knowing that the DCB is on the Rotating Base, a user might look up “Rotating Base” in the table of contents. According to the TOC, information on the Rotating Base begins on [Page 196](#). But, because this manual follows the product structure defined by the BOMs, the DCB will not be found anywhere near [Page 196](#). The user will find, instead, a reference in the Rotating Base BOM to further detail on the Detector Plate Assembly. Turning to the directed page ([Page 226](#)), the user will find a further reference to detail on the DAS Collector. Digging down this one more level brings the user to the DAS Assembly drawing and BOM (beginning on [Page 228](#)), and there the user will find the DCB (mark 6).

WHAT IS A BOM?

A BOM is a Bill (list) Of Materials (parts) used to produce an item. The Parent Item is the thing that is made from the parts listed on the BOM. Because BOMs describe the materials required for a specific manufacturing process, they may not necessarily list every single item used to produce the parent item. The parts on a BOM themselves may be the result of a separate process, and thus have their own BOMs, listing the secondary items required to produce these primary parts.

BOM TREES

As indicated above, this manual is organized according to the overall Product Structure of the LightSpeed¹⁶ CT System's major components: Table, PDU, Gantry and Console. The many levels of BOMs contained in the Product Structures for the Table, Gantry and Console are represented on the next few pages in Tree form, revealing which BOM parts are themselves parent items to secondary BOMs. The trees list only key items, providing the Mark number (in parentheses) from the BOM of its parent item.


Section 3.0 — System Matrices

The BOM Trees shown on the following pages act as **SYSTEM MATRICES**, providing visual cues to this manual's organization. No BOM tree is shown for the Compact PDU or NGPDU.

Note: As you use these trees, keep in mind that BOMs are not necessarily organized by function, but rather by manufacturing process.

3.1 LightSpeed 4.X System Matrix


The information in this manual applies to LightSpeed 4.X (LightSpeed¹⁶) systems with the following catalog numbers: B7868JA, B7868JB, B7868JC, B7868JD, B7868JE, B7868JF and B7868MB. The matrix shown below shows the catalog structure of these systems, and includes page numbers to point the user to the relevant portions of this manual. Please note that not all catalog items are shown in this tree, only those that are relevant to this document. Furthermore, not every catalog item applies to every system — applicability notes indicate inclusions and exclusions (e.g., “Cat: ALL” indicates that this item applies to all of the listed LightSpeed 4.X system catalogs).

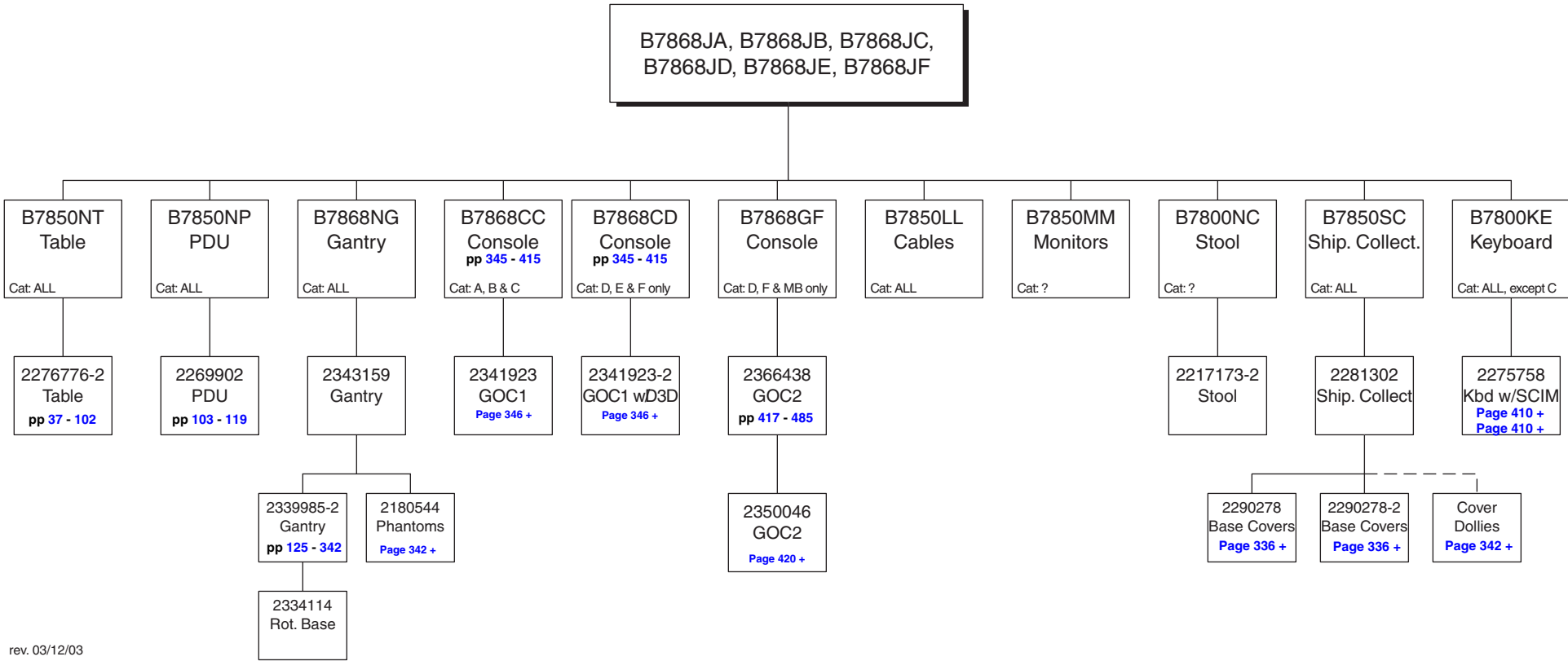

rev. 03/12/03

Figure 1-1 LightSpeed 4.X System Catalog Tree

3.2 Table Sub-System Matrix

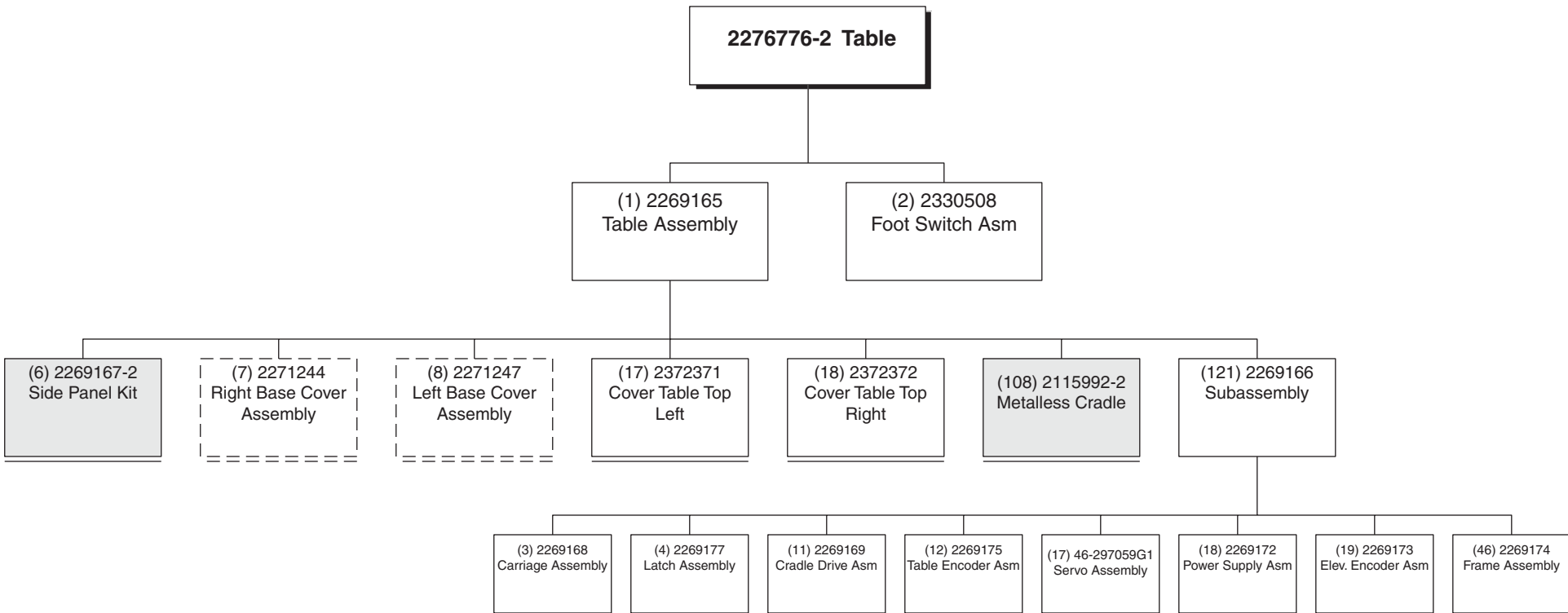

KEY:

= No FRUs at this level

= No drawing available

= No structure below this level

(#) = Mark number (from parent BOM)

Page ## + = Starting page number of info

rev. 04/23/03

Figure 1-2 Table BOM Tree

3.3 Gantry Sub-System Matrices

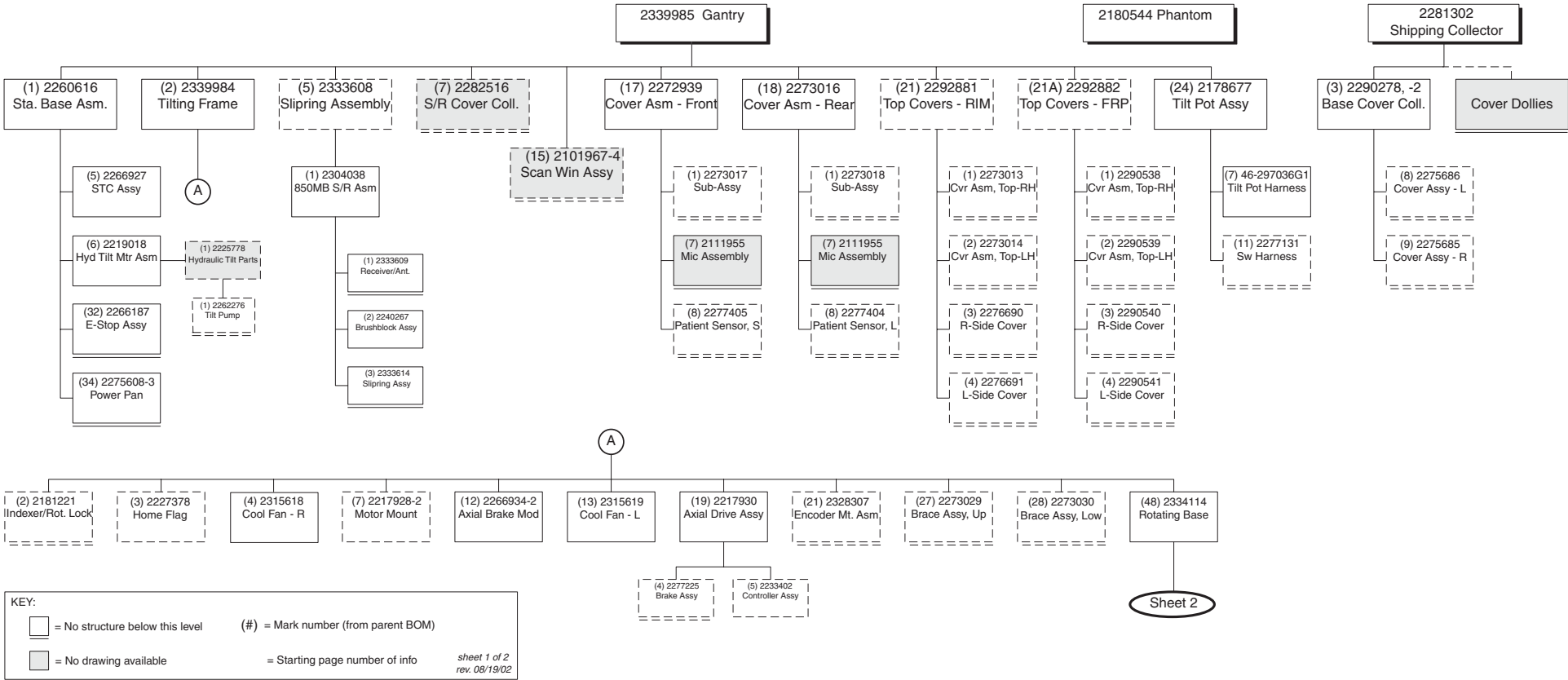


Figure 1-3 Gantry BOM Tree (sheet 1 of 2)

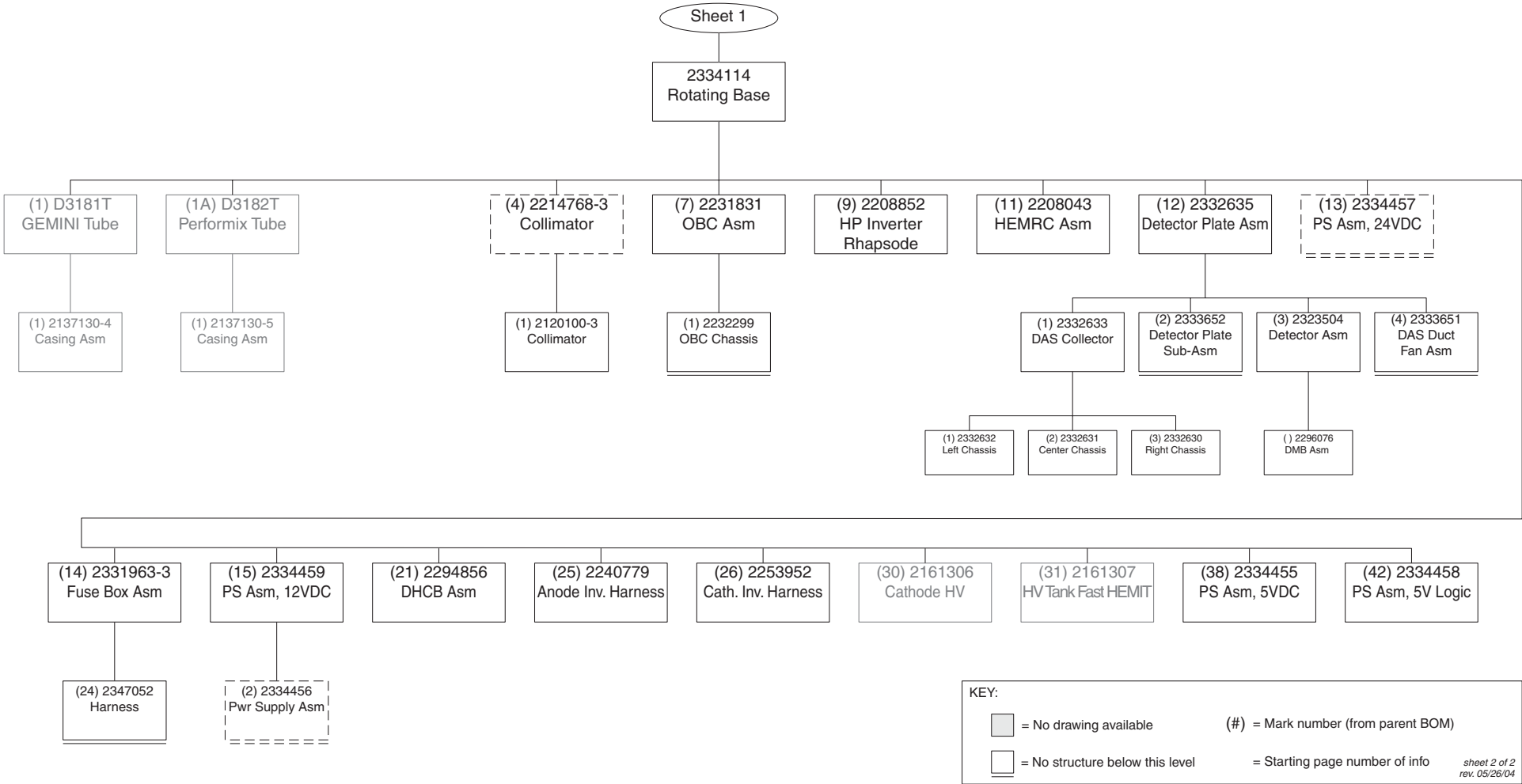


Figure 1-4 Gantry BOM Tree (sheet 2 of 2)

3.4 Console Sub-System Matrix - GOC1 (Global Operator Console - Octane2)

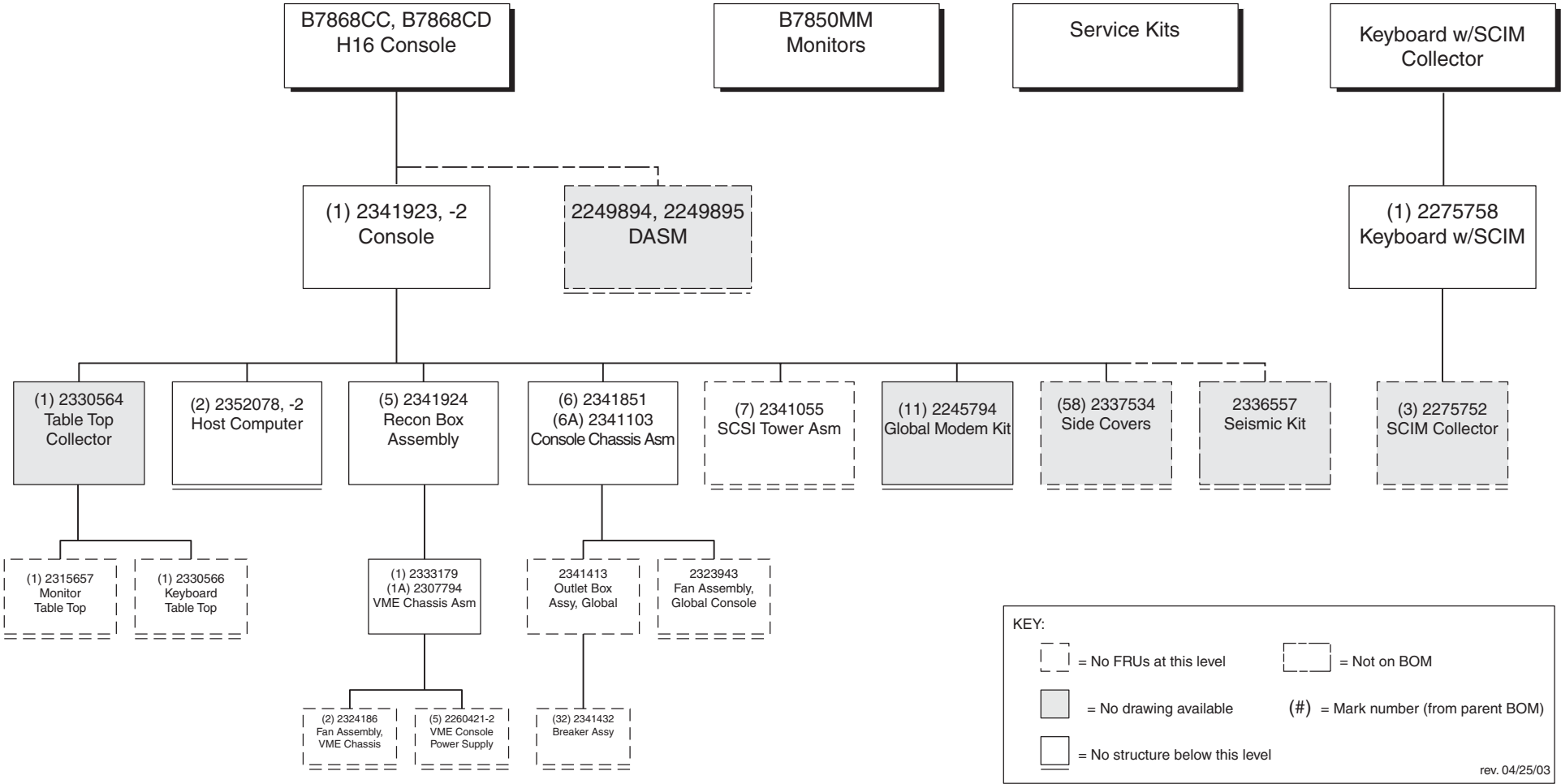


Figure 1-5 GOC1 Console BOM Tree

3.5 Console Sub-System Matrix - GOC2 (Global Operator Console - Linux)

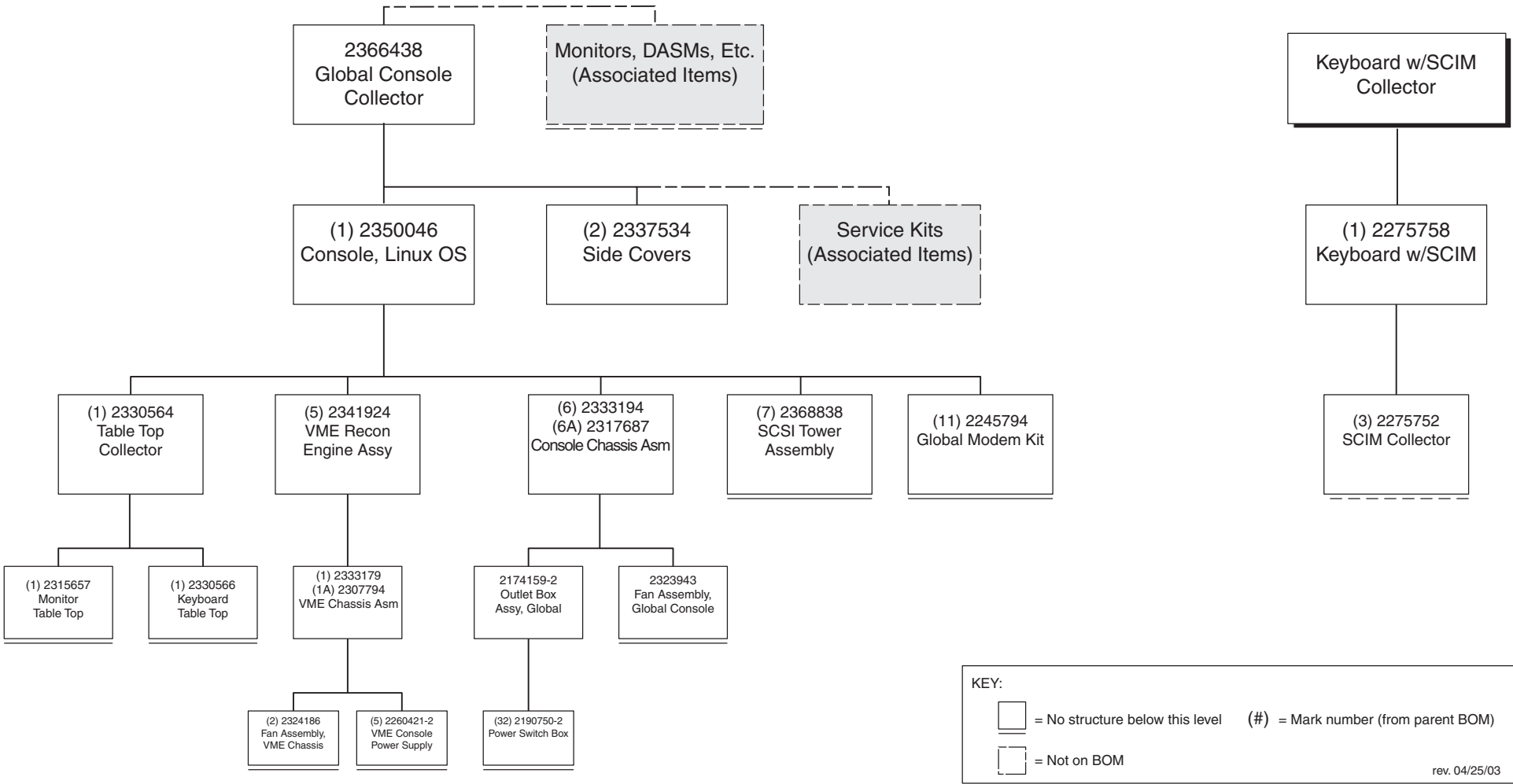


Figure 1-6 GOC2 Console BOM Tree

3.6 GRE Xstream Console (GOC3 or GOC4)

For GOC3 (Global Operator Console, version 3) Parts Information: Please refer to the GOC3 Illustrated Parts List, Direction 2394950-100. This IPL is available for order in paper, or may be found in PDF form on the GOC3 Installation and Service Information CD-ROM (2384451-200) that came with the GOC3 Console. The PDF may also be found on Revision 6 or later of the *LightSpeed 4.X Service Information* CD-ROM (2353846-200) and on GE's Common Documentation Library.

For GOC4 (Global Operator Console, version 4) Parts Information: Please refer to the GOC4 Illustrated Parts List, Direction 5134893-100. This IPL is available for order in paper, or may be found in PDF form on the GOC4 Installation and Service Information CD-ROM (5134894-200) that came with the GOC4 Console. The PDF may also be found on Revision x or later of the *LightSpeed 4.X Service Information* CD-ROM (2353846-200) and on GE's Common Documentation Library.

End of Chapter.