

Installing the License Key for Context Vision Filter on Leonardo Workstation


It is absolutely essential to use the right accounts (Administrator, Meduser) otherwise the created license file will not have the right permission!

1. Login as Administrator

2. Open Windows Explorer → Navigate to C:\Leonardo\MriSiteData\ContextVision

Select Tools → Folder Options → View → select Show hidden files and folders

Select OK to save setting.

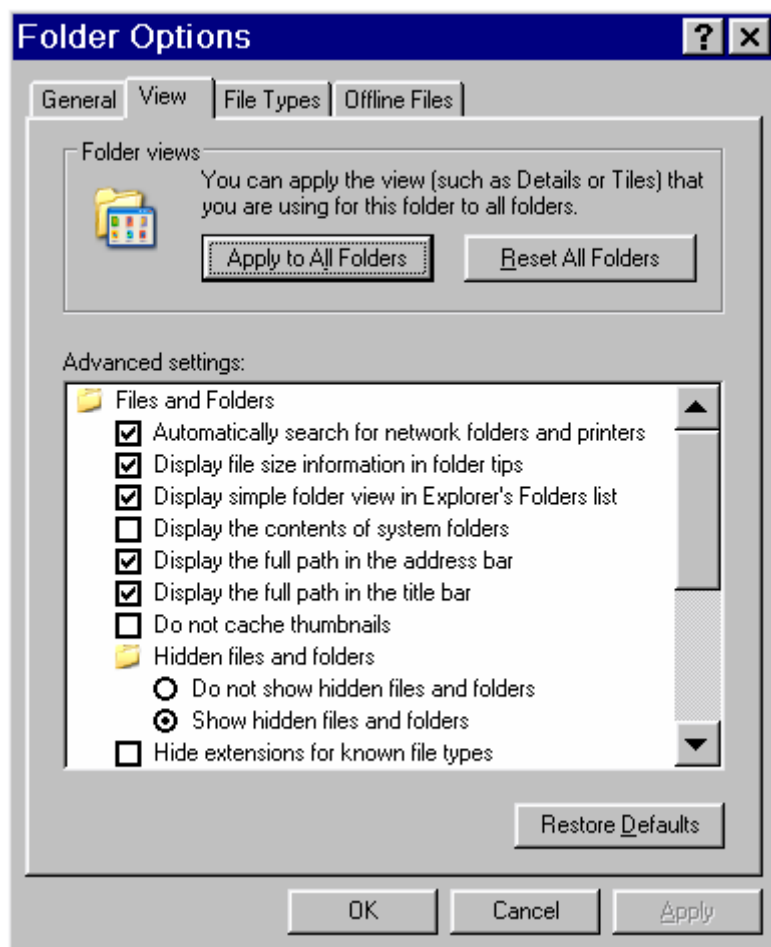


3. Delete C:\Leonardo\MriSiteData\ContextVision\.\siemensIRU-PCNT-init (the file name has a "dot" in front)

4. Log off Administrator

5. Login as Meduser


6. Switch to advanced user → Ctrl/esc → Run → cmd (open comdtol)

7. Enter cvlicense+your key + enter (key is just an example)

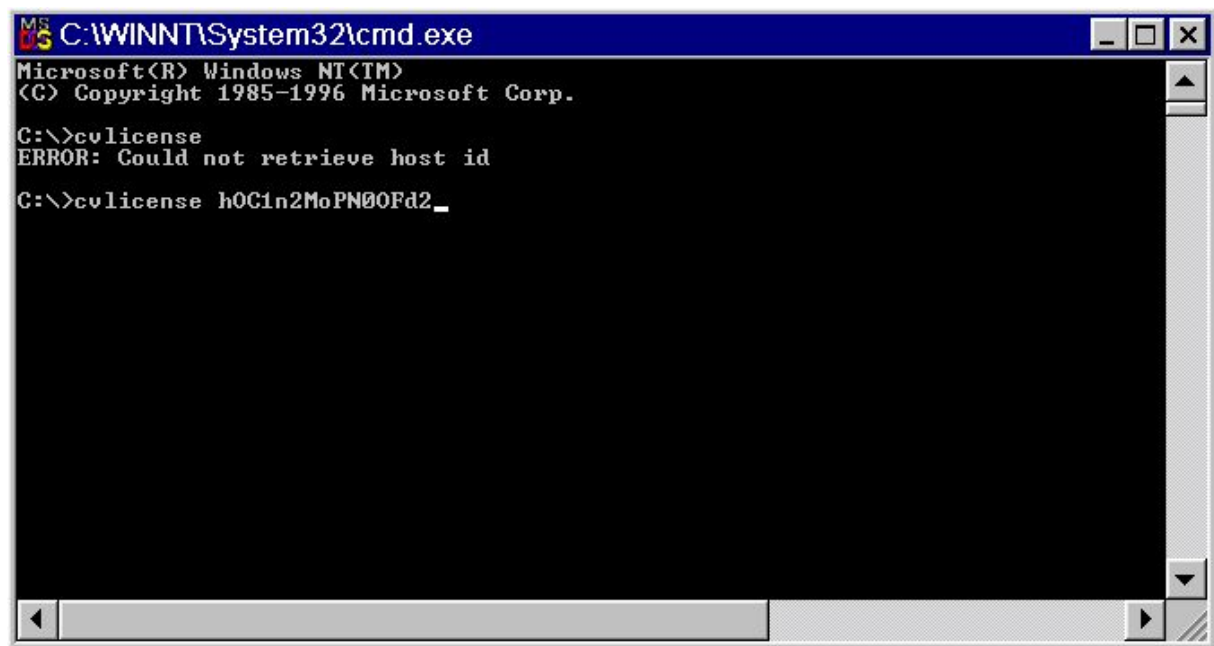


8. Please check the file properties for file (this step is not necessary, just a check)

C:\Leonardo\MriSiteData\ContextVision\.siemensIRU-PCNT-init Read only - has to be unchecked (the file name has a "dot" in front)

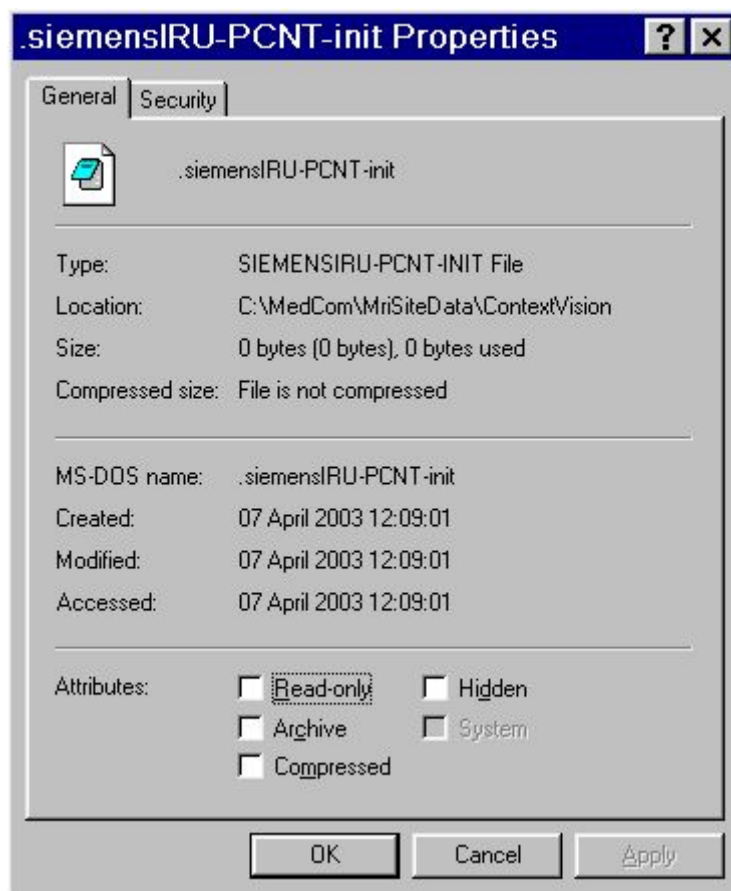

Check also in Security Permissions, that owner SiemensMed has Change(RWXD) permission.

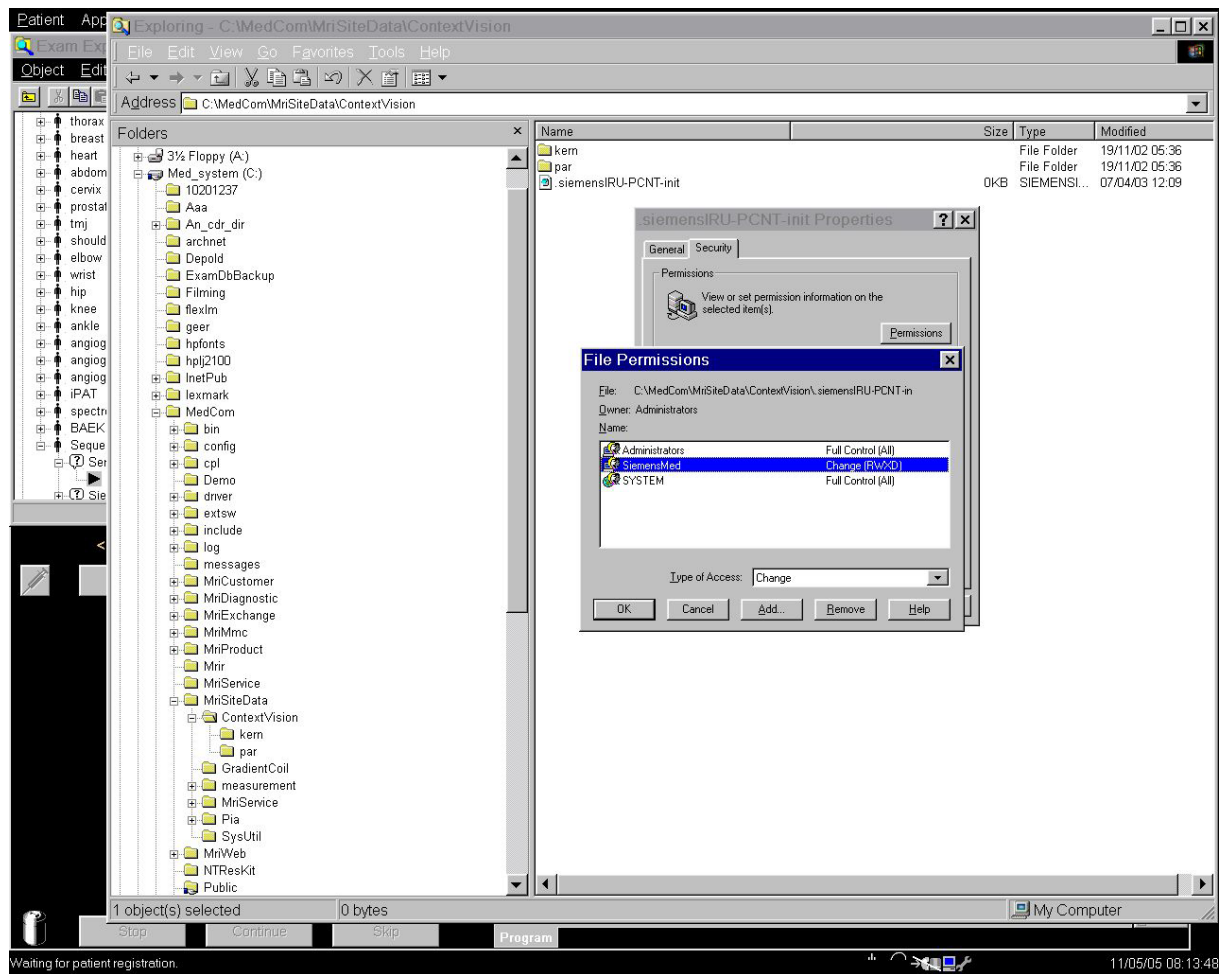

9. Reboot System