




1-1 Cable map (Type A Detachable Table)

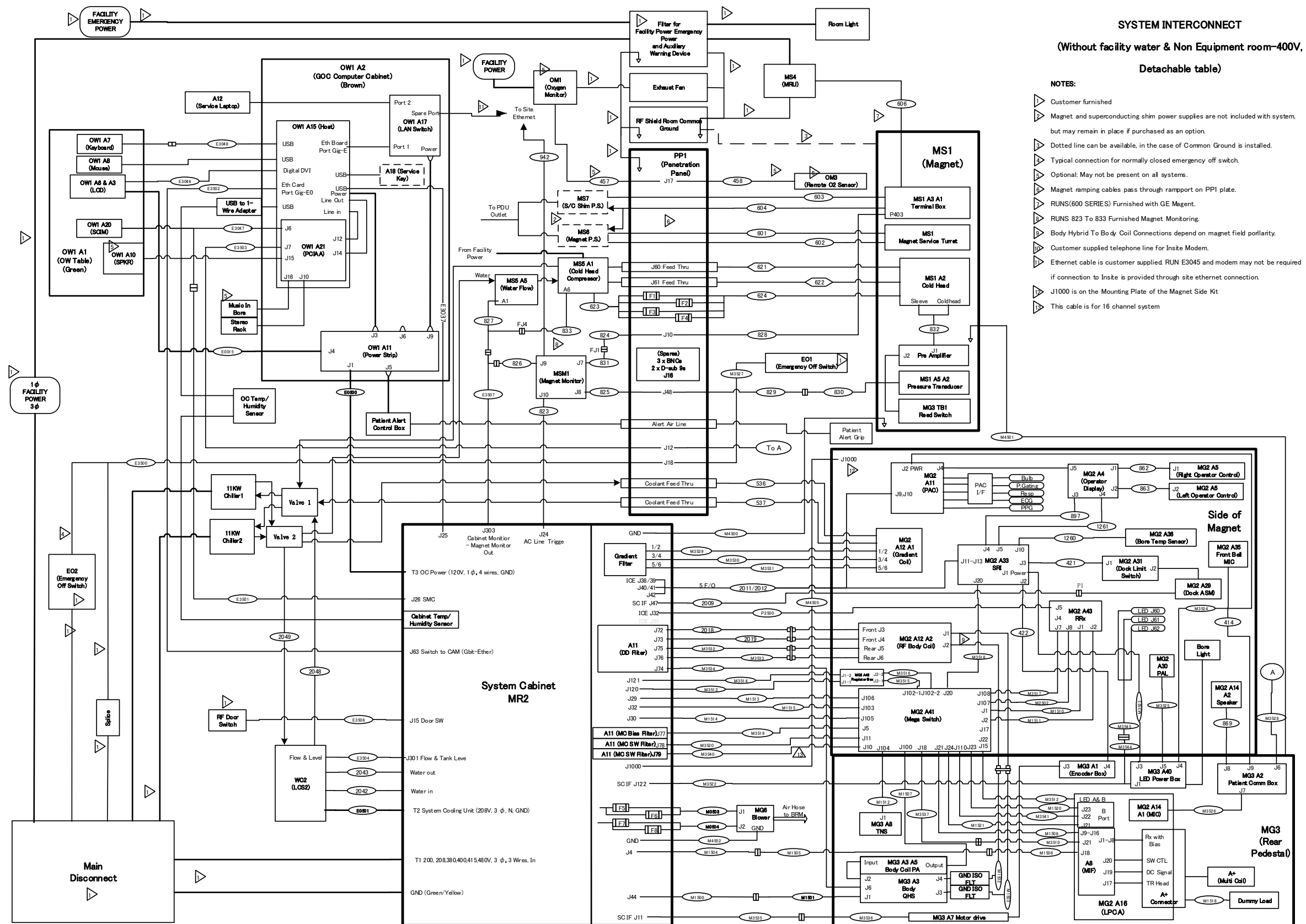


1-2 Cable map (Type B Detachable Table)

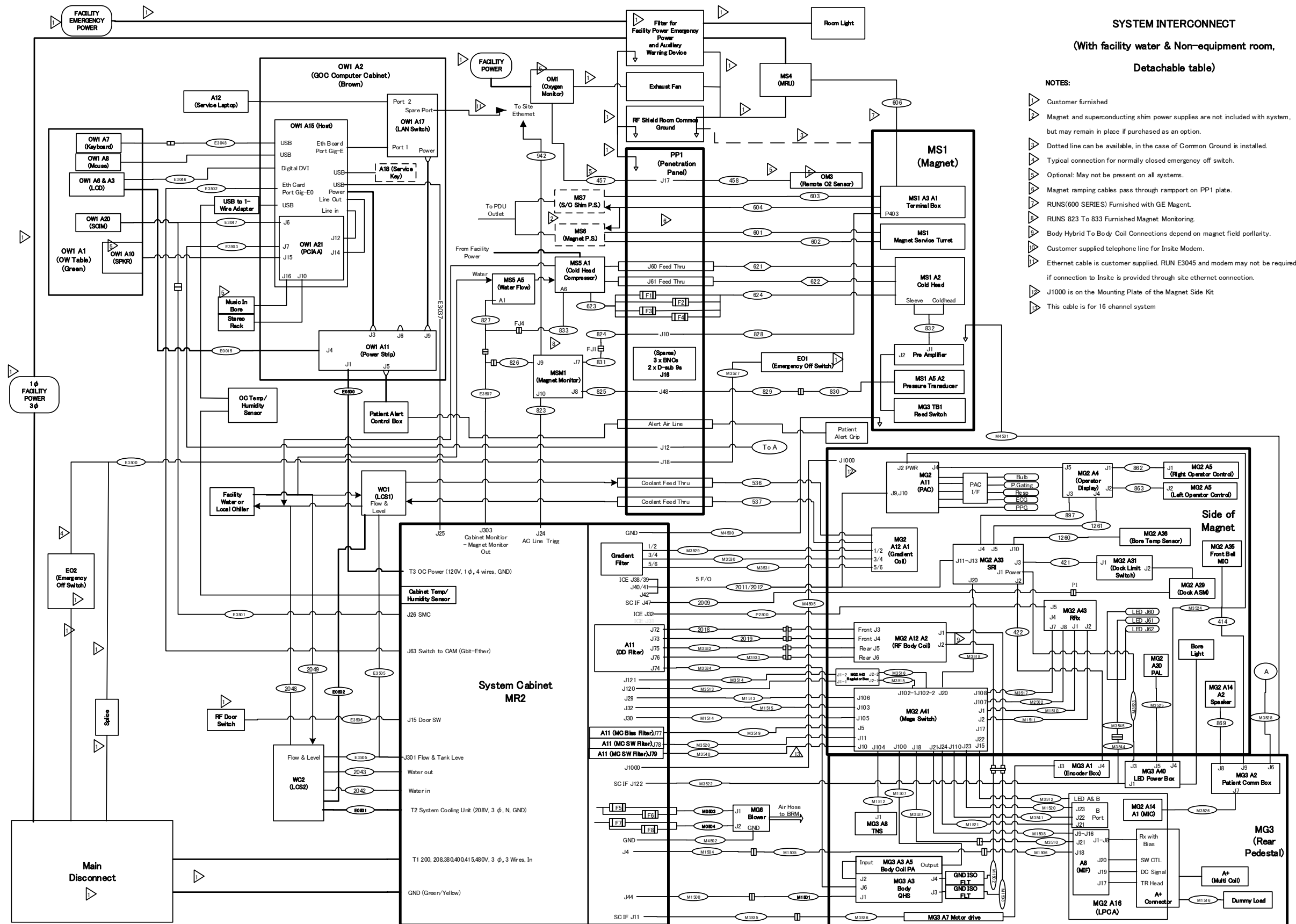


1-3 Cable map (Type C Detachable Table)

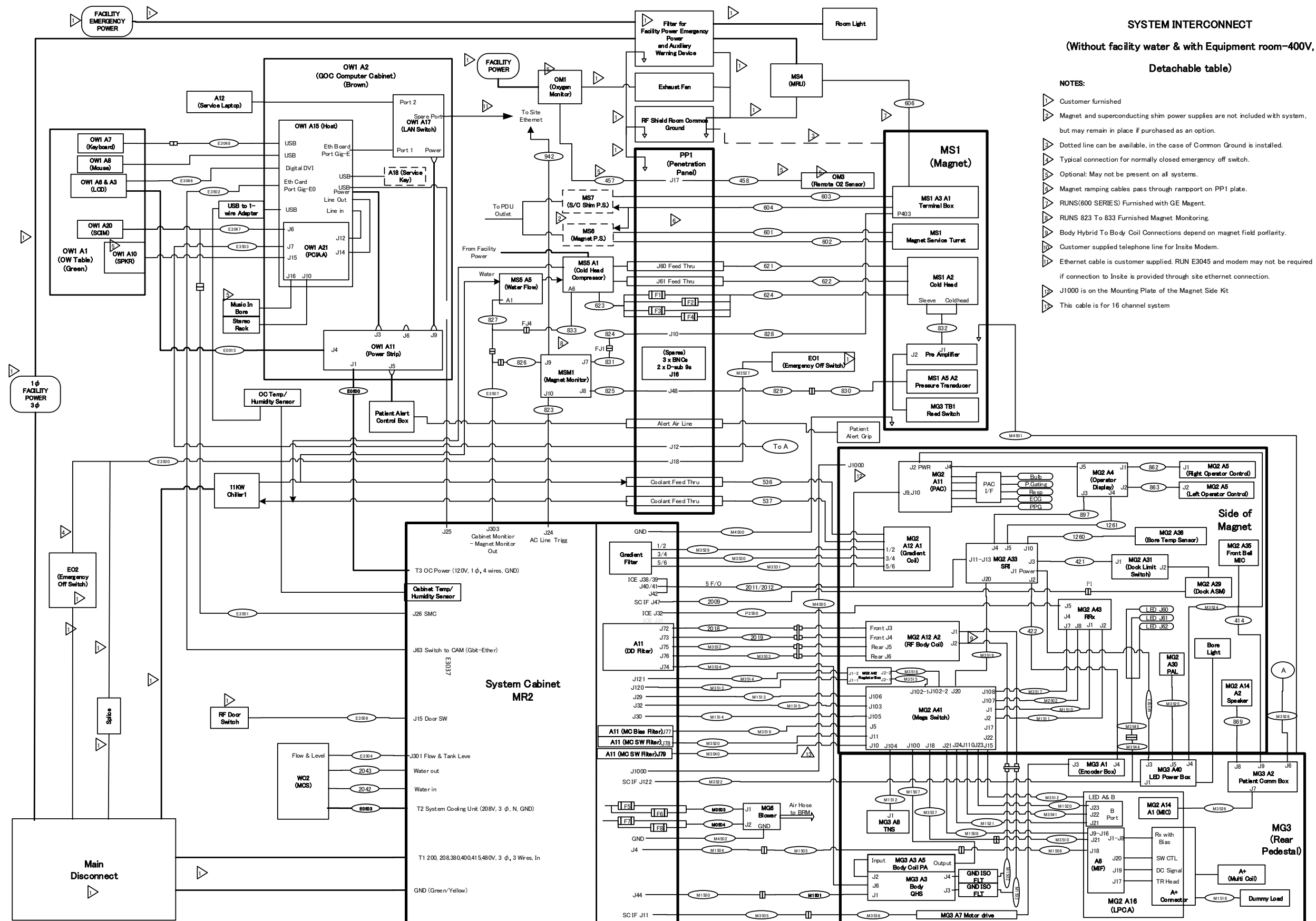

1-4 Cable map (Type D Detachable Table)

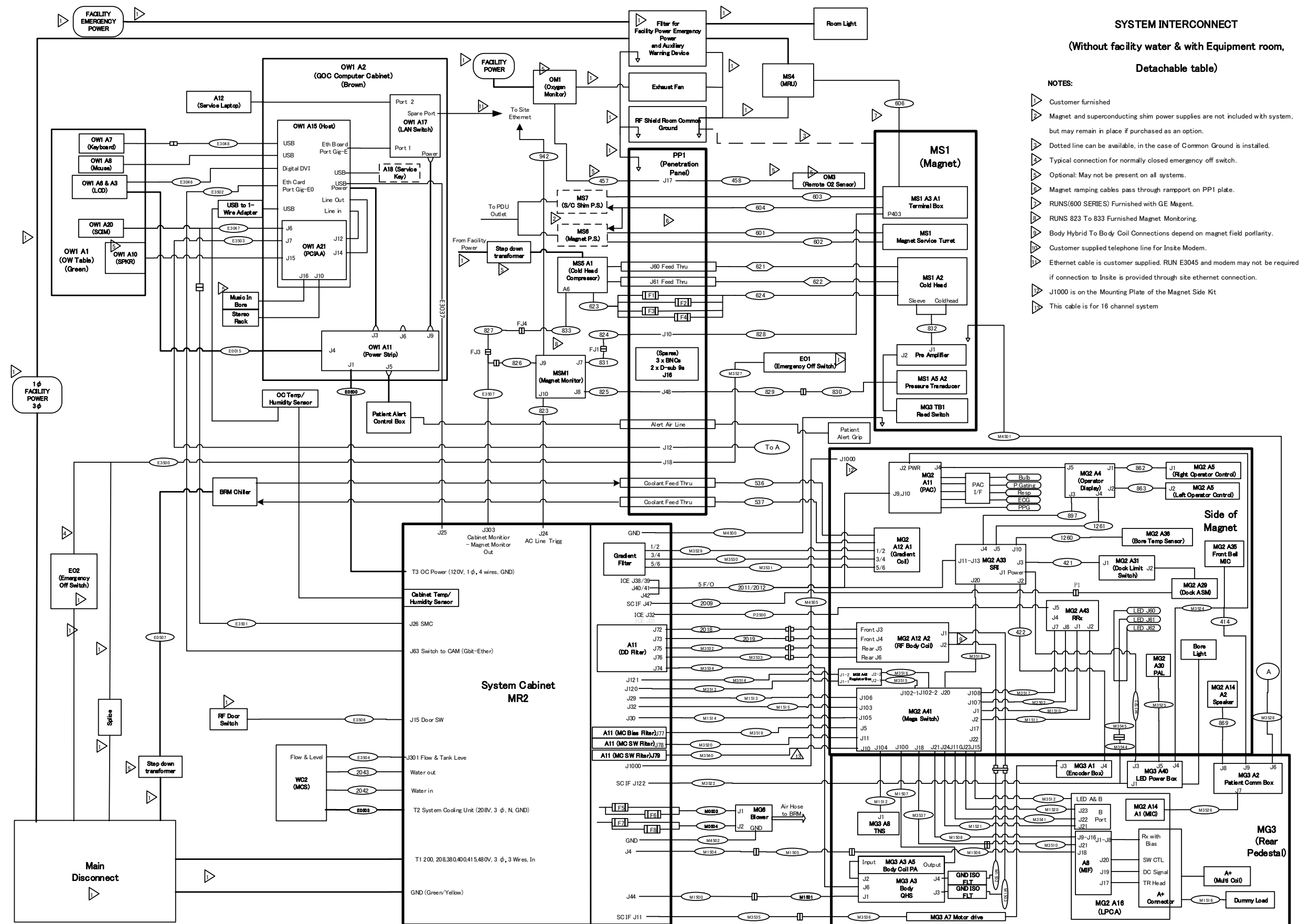

1-5 Cable map (Type E Detachable Table)

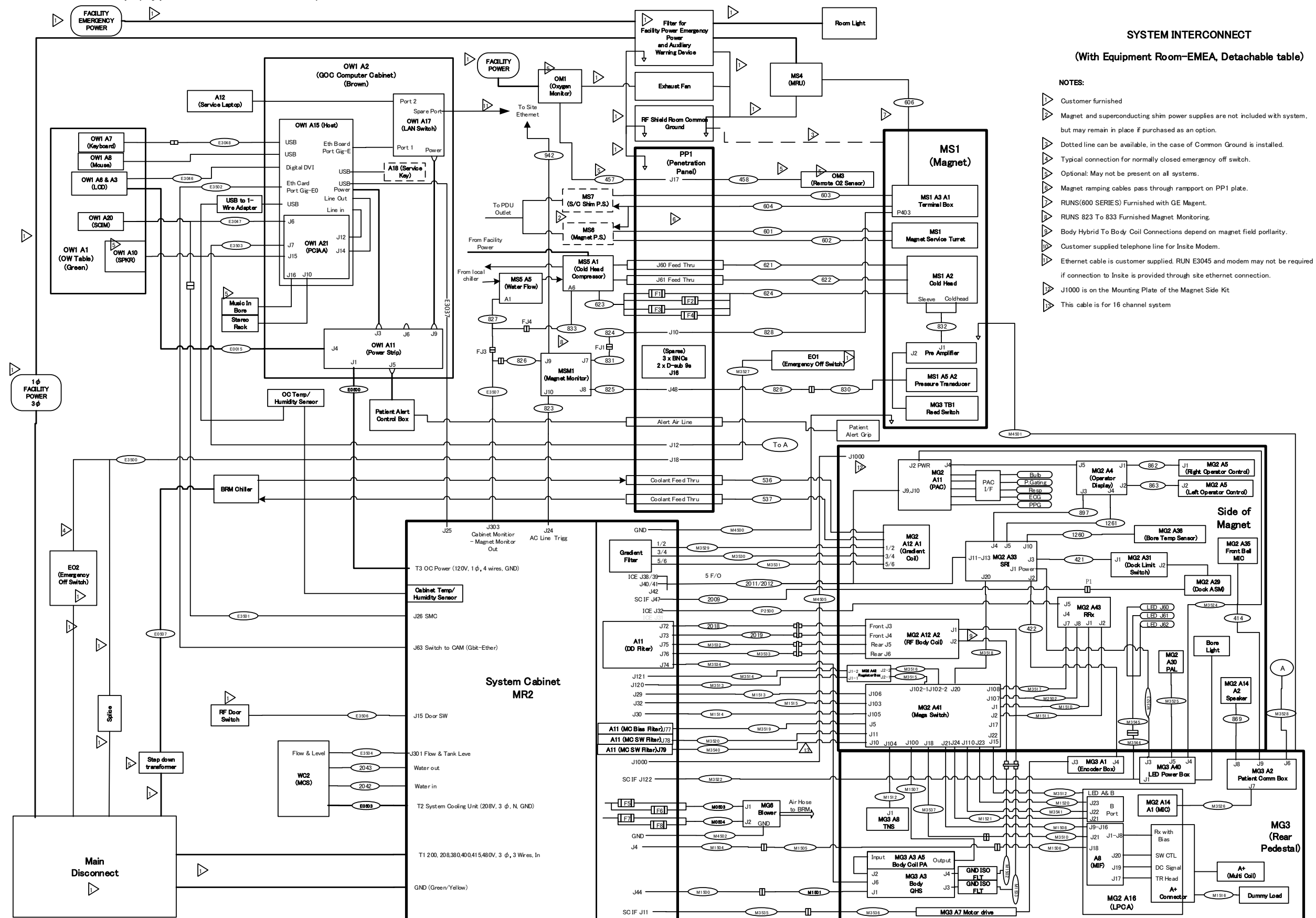

SYSTEM INTERCONNECT

(Without facility water & Non Equipment room-400V,

Fixed table)


NOTES:

- 1 Customer furnished
- 2 Magnet and superconducting shim power supplies are not included with system, but may remain in place if purchased as an option.
- 3 Dotted line can be available, in the case of Common Ground is installed.
- 4 Typical connection for normally closed emergency off switch.
- Optional: May not be present on all systems.
- 6 Magnet ramping cables pass through rampout on PP1 plate.
- 7 RUNS(600 SERIES) Furnished with GE Magent.
- 8 RUNS 823 To 833 Furnished Magnet Monitoring.
- 9 Body Hybrid To Body Coil Connections depend on magnet field polarity.
- 10 Customer supplied telephone line for Insite Modem.
- 11 Ethernet cable is customer supplied. RUN E3045 and modem may not be required if connection to Insite is provided through site ethernet connection.
- 12 J1000 is on the Mounting Plate of the Magnet Side Kit
- 13 This cable is for 16 channel system


Fixed table)

- Customer furnished
- Magnet and superconducting shim power supplies are not included with system, but may remain in place if purchased as an option.
- Dotted line can be available, in the case of Common Ground is installed.
- Typical connection for normally closed emergency off switch.
- Optional: May not be present on all systems.
- Magnet ramping cables pass through rampport on PP1 plate.
- RUNS(600 SERIES) Furnished with GE Magnet.
- RUNS 823 To 833 Furnished Magnet Monitoring.
- Body Hybrid To Body Coil Connections depend on magnet field polarity.
- Customer supplied telephone line for Insite Modem.
- Ethernet cable is customer supplied. RUN E3045 and modem may not be required if connection to Insite is provided through site ethernet connection.
- J1000 is on the Mounting Plate of the Magnet Side Kit
- This cable is for 16 channel system


1-7 Cable map (Type B Fixed Table)

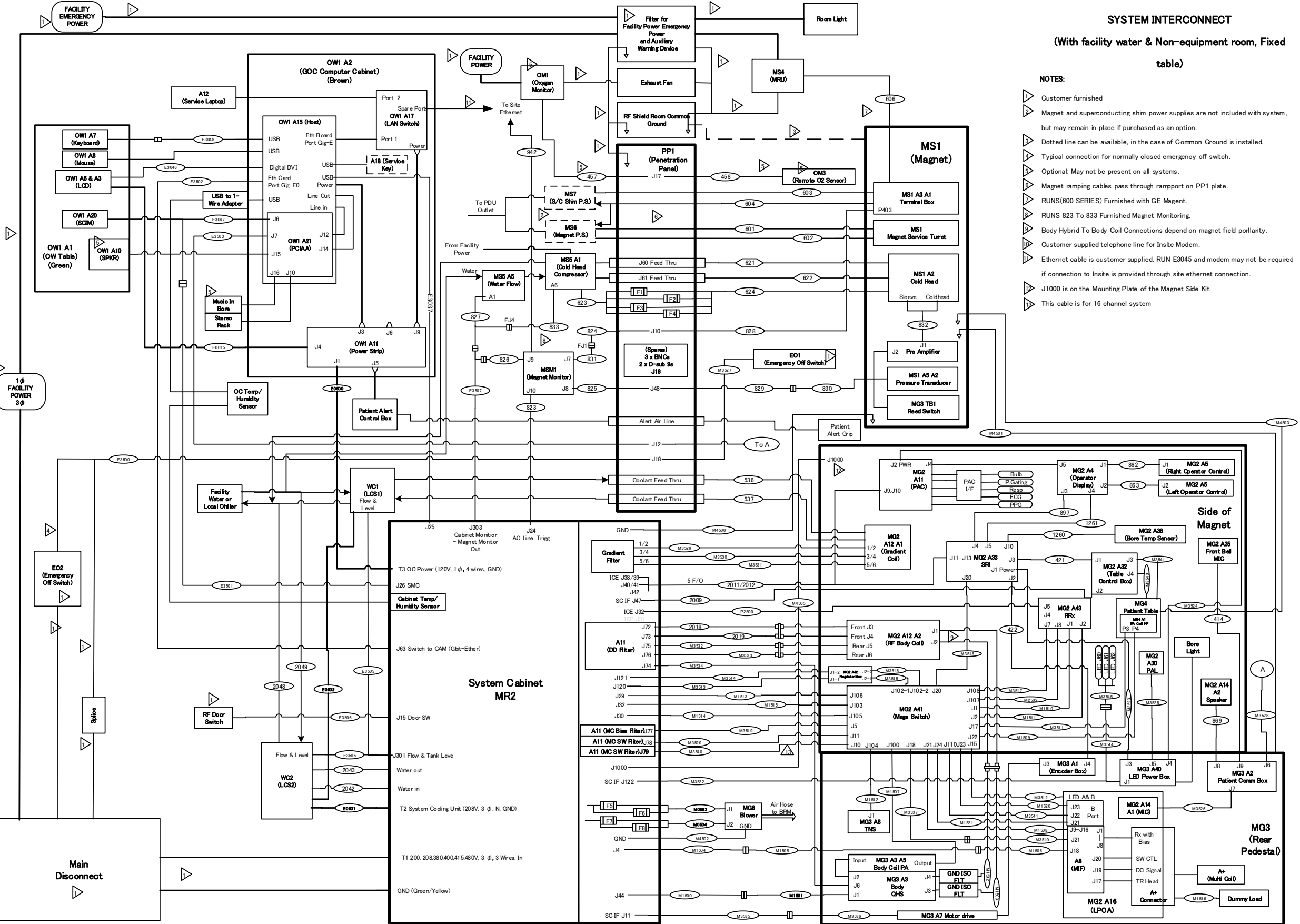

1-8 Cable map (Type C Fixed Table)

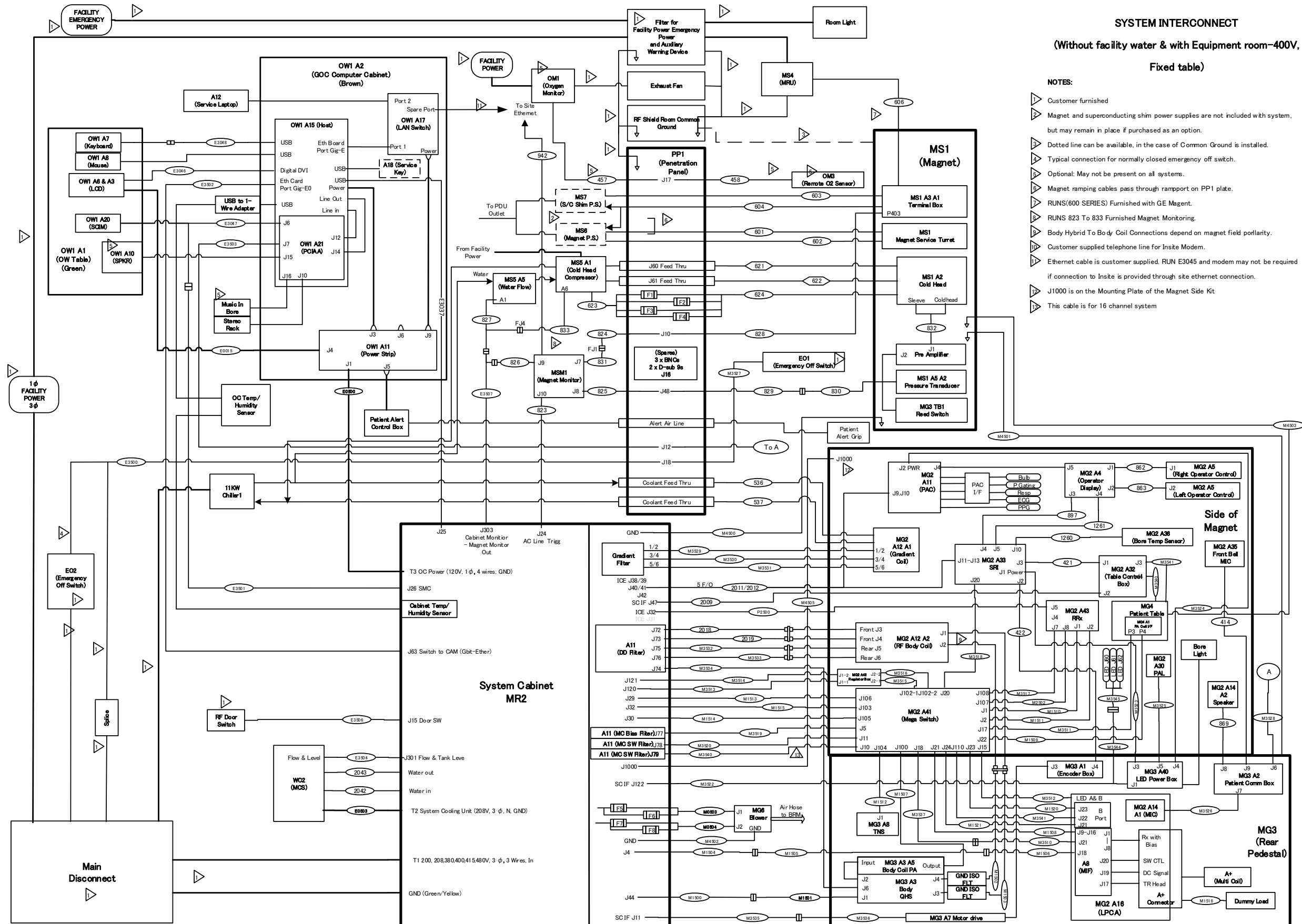

1-9 Cable map (Type D Fixed Table)

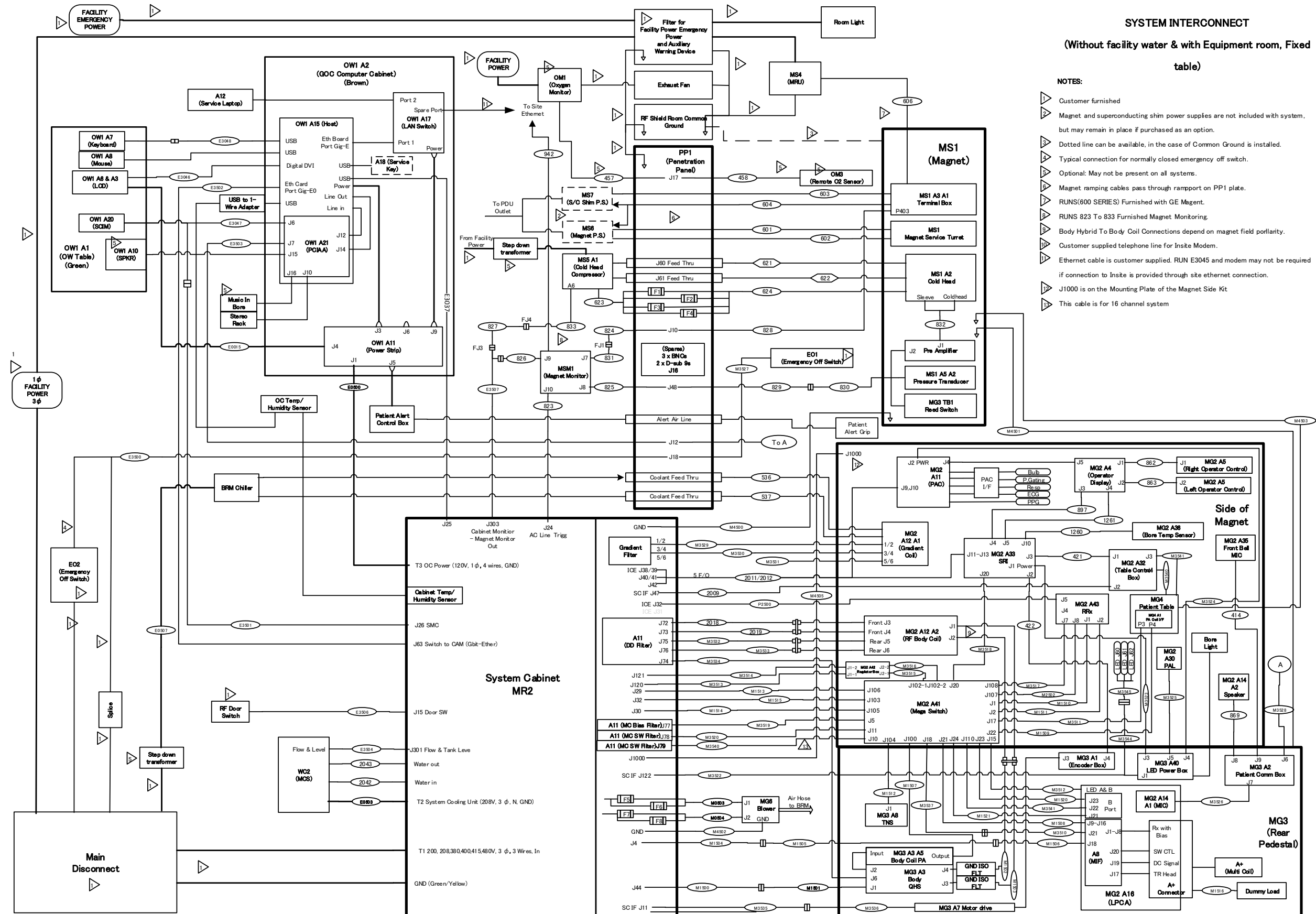

1-10 Cable map (Type E Fixed Table)

