

CREATE A SITE TAPE

You should create a site tape at this time. This is done via the Main Startup Menu. Click in the lower left corner and type the word “site” (lower case) and press <CR>. After the Save Site / Restore Site menu appears, push the tape in, wait for the green light to stop flashing and select “1” from the menu to Save Site. The system will eject the tape and reboot the system when finished. (Restore Site is done in a similar fashion.)

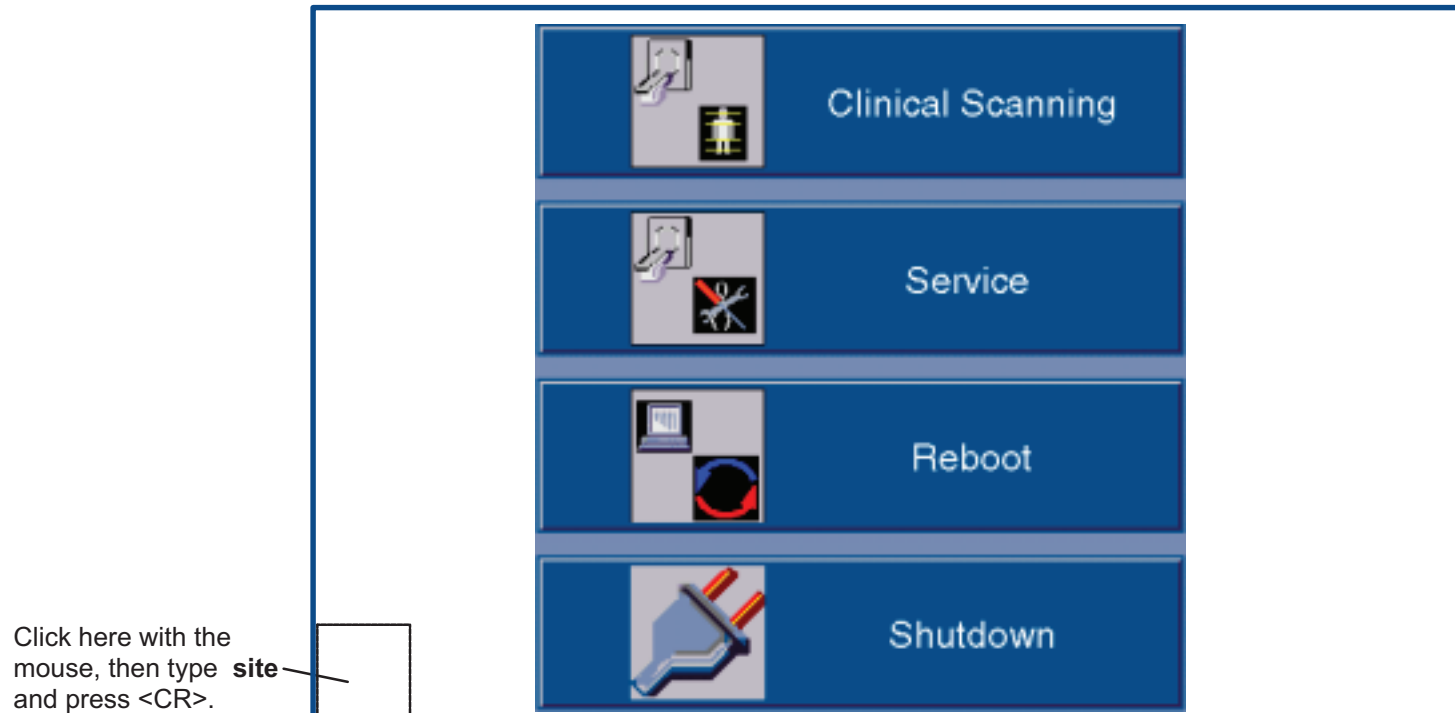


FIGURE 61 Starting the Save Site / Restore Site Menu

ALIGN LASERS

NOTE

[Laser Alignment](#) is performed in the Mechanical Calibration manual.