

GE Healthcare

Signa HDx 3.0T Series
Field Replaceable Units (FRU)

Technical Publication

© 2006 – 2011 General Electric Company
GE Medical Systems, a General Electric Company, going to market as GE Healthcare
3200 N. Grandview Boulevard
Waukesha, Wisconsin 53188
USA

www.gehealthcare.com



5166013
Rev 14

**COVER
BOTTOM
CROP LINE**

**COVER
BOTTOM
CROP LINE**

**5/16" RADIUS
CORNER TRIM
TYPICAL**



**COVER
LEFT
CROP LINE**

**COVER
FOLD
LINE** | **(2) CENTER STAPLES
FOR BINDING**

**COVER
RIGHT
CROP LINE**



**5/16" RADIUS
CORNER TRIM
TYPICAL**

SEND COMMENTS ABOUT THIS MANUAL TO:

GE MEDICAL SYSTEMS
MANAGER, SERVICE ENGINEERING
W-808
3200 N. GRANDVIEW BLVD.
WAUKESHA, WISCONSIN 53188
USA

This manual belongs to:

INTRODUCTION TO FIELD REPLACEABLE UNIT (FRU)

This direction is a guide to the FRU's for HDxt and HDx.

Parts designated as a FRU are critical to maintain the safety and/or functionality of the system. FRU's are available on a Priority 1 basis from Service Parts Operation.

As a reference aid, the FRU's are listed under an upper-level part number where they are used. In most cases the upper-level part number is not a FRU. In some cases the parts listed under an upper-level assembly may have been nested in intermediate sub-assemblies which are not shown.

Where possible, FRU's which are part of Option Catalogs are shown under the area where they would be installed, rather than being listed in the Coils or Options Sections (SEE TABLE OF CONTENTS).

LIST OF EFFECTIVE PAGES

<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>
Cover Page .	14	22	14	31	11
1	14	23	11	32 to 39	13
2	13	24	14	40 to 49	11
3	14	25	12	50 to 61	13
4 to 15	11	26	14	62 to 63	14
16	14	27 to 28	11	64	13
17 to 21	11	29 to 30	14		

TITLE	PAGE
INTRODUCTION TO FIELD REPLACEABLE UNIT(FRU)	1
LIST OF EFFECTIVE PAGES	1
TABLE OF CONTENTS	2
PICTORIAL TABLE OF CONTENTS	3
MAGNET AND ENCLOSURE (ASSEMBLIES)	
FRONT COVER, MISC. PARTS, & TRM	5
BRM & TRM LIST (OF PARTS NOT SHOWN)	5
ENCLOSURE & VACUUM RELATED	6
FRONT & REAR ENDBELL, BRIDGE SUPPORT	6
3T94 ENCLOSURE COMPONENTS	7 & 8
REAR PEDESTAL PARTS	9
SPLIT BRIDGE & DOCKING PARTS	10
LOW PROFILE CARRIAGE ASSEMBLY	11
SRI/PAC & TOUCH-UP PAINTS	12
MRU, MONITORS, & TRANSFORMER	13
HDx MAGNET ENCLOSURE	14 & 15
HDx BRIDGE PARTS	16
HDx LPCA	17
PATIENT SUPPORT TABLE	18 & 19
CHILLER	20 & 21
PENETRATION PANEL	22
OPERATOR WORKSPACE	24
GOC COMPUTER CABINET, SCIM AND KEYBOARD	25
MAIN DISCONNECT PANEL	28
HFD-S/PDU & HFD/PDU CABINET	29, 30
RFS CABINET	32, 33, & 34
SYSTEM CABINET	35
CAM, MGD, RRF, & ASC MODULES	36, 37, 38, 39, 40
HDx 32 CHANNEL CABINET (OPTION)	41
RF ACCESSORY CABINET	43
NB RF AMP CABINET	44 & 45
8KW MNS BB RF AMP CABINET (OPTIONAL)	46
TWIN ACCESSORY CABINET (TAC)	47, 48, 49
REPLACEMENT CABLES	51 & 52
UNIFIED SURFACE COIL PHANTOM SET	53
SERVICE TOOLS	54
TOOLS, REPAIR KITS, AND HEAD COILS	55
COILS	56, 57 & 58
BrainWave OPTIONS	58
MNS OPTIONS	60
SOFTWARE OPTION KEYS	61
REVISION HISTORY	62

INTRO TO FRU'S	PAGE 1
MRU & MONITORS	PAGE 13
CHILLER PARTS	PAGE 20
CABLES	PAGE 51
COILS	PAGE 57
MNS OPTIONS	PAGE 60
S/W OPTIONS	PAGE 58
BrainWave	PAGE 58

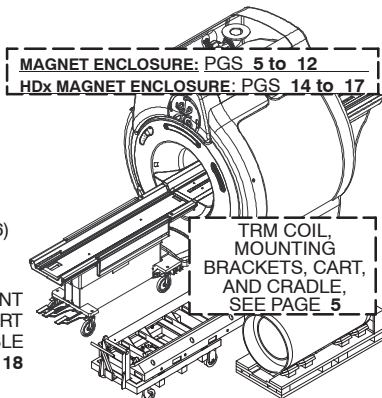


PATIENT AIR SUPPLY
BLOWER CABINET (MG6)
3T94 – PAGE 8
G3 – PAGE 16



PENETRATION
PANEL
PAGE 22

PATIENT
SUPPORT
TABLE
PAGE 18



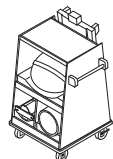
MAGNET ENCLOSURE: PGS 5 to 12

HDx MAGNET ENCLOSURE: PGS 14 to 17

TRM COIL,
MOUNTING
BRACKETS, CART,
AND CRADLE,
SEE PAGE 5

MAGNET ROOM

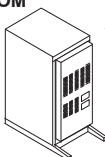
OPERATOR/EQUIPMENT ROOM



SERVICE TOOLS
PAGE 54



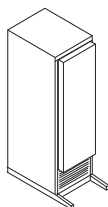
OPERATOR
WORKSPACE
PAGES 24 & 25



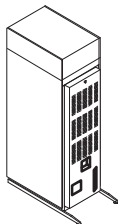
ACCESSORY
CABINET
PAGES
47 THRU 49



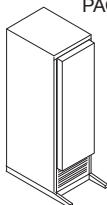
MAIN
DISCONNECT
PANEL
PAGE 28



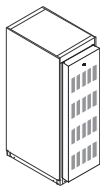
32 CHANNEL
CABINET
PAGE 41



HFD-S/PDU
PAGE 29
HFD/PDU
PAGE 30



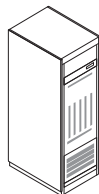
RFS &
SYSTEM
CABINET
PAGES
33, THRU 40



RF ACCESSORY
CABINET
PAGE 43



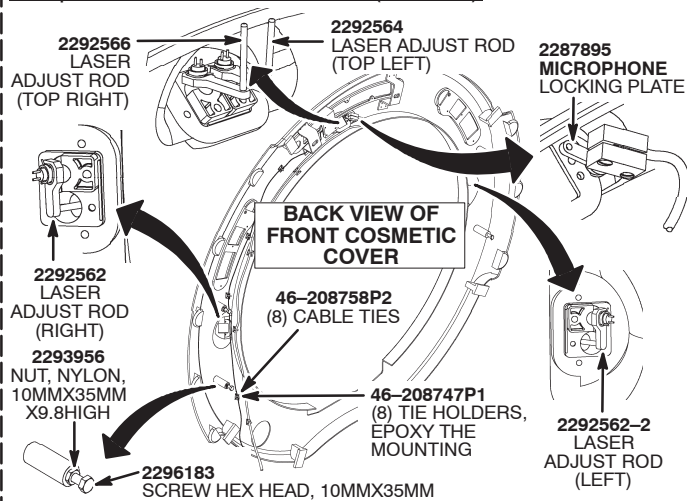
NB RF AMP
CABINET
PAGES
44 & 45



8KW BB RF
AMP CABINET
PAGE 46

PICTORIAL TABLE OF CONTENTS

[illegible]

WideOpen FRONT COSMETIC COVER & PARTS (BACK VIEW):

NOTE: The following LCC 3.0T FIXED W/TRM & Enclosure related parts **Not Shown**;
MISCELLANEOUS:

2203803	SPLIT O-RING
2242756	RELIEF VALVE, 3/8 NPT X 7 PSIG, GEN
2299843	MAGNET PRESSURE TRANSDUCER, SUNX DP2-41N
46-294231G1	MAGNET RUNDOWN UNIT
46-294231P1	BUSSMAN, F1, 110VAC .400 GDC

TRM:

5113878*	RF SHIELD FOR TRM COIL
2144093-9	TRM WITH MOUNTING BRACKETS, CART, AND CRADLE
2208999-7	3.0T RF BODY COIL AND CASE
2208999-12	HDx new color 3.0T RF BODY COIL & Purple/Black CASE
2218148	COLD SHRINK (FOR REAR ENDBELL AND TRM COIL)
2299498	COOLING CIRCUIT REPAIR KIT
2383070 or 5250104	DIODE SLOT BOARD
5123812	(1260) BORE TEMPERATURE SENSOR CABLE
5123813	(1261) OPERATOR DISPLAY TO SRI CABLE
5127441	GRADIENT CABLE CLAMP BLOCK KIT
2293612	MALE CONNECTOR 3/4" HOSE TO 3/4" NPT
2324086-2	5 GALLON PAIL w/ETHYLENE GLYCOL WARNING LABEL

***2389049** RF Shield Change Out Kit – MUST be ordered for the RF Shield replacement

MAGNET AND ENCLOSURE (ASSEMBLIES)

ENCLOSURE & VACUUM RELATED PARTS:

2293771
RADIAL SHIM SPACER, 1.5MM
2293771-2
RADIAL SHIM SPACER, 0.8MM
2241176
RUBBER ISOLATION PAD

**FRONT ENDBELL PARTS:**
(TRM ONLY)

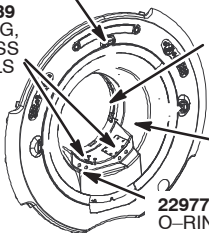
2162554-4
FRONT COVER LASER

2295639
O-RING,
ACCESS
PANELS

GASKET,
ENDBELL
TO RF TUBE
2294358, 4MM
2294358-2, 6MM

FRONT
ENDBELL ASM
2296177

2297761
O-RING, WATERLINE
FEEDTHROUGH

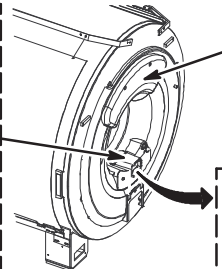
**BRIDGE SUPPORT PARTS:**
(TRM ONLY)

2293539
BRIDGE
SUPPORT
BUSHING

2293956
NYLON HEX NUT,
M10 X 1.5MM
X 9.5MM HIGH

**REAR ENDBELL ASSEMBLY:**
(TRM ONLY)

2284170
REAR ENDBELL
VACUUM ASSEMBLY;
& 2218148
COLD SHRINK
(FOR REAR ENDBELL &
TRM COIL)
& 2218149
COMPRESSION
SPLICER
(FOR REAR ENDBELL)



2358470
COMFORT MODULE
ASSEMBLY

2206437
M4 X 12mm NYLON
SCREW FOR
SPEAKER GRILL

GRADIENT GASKETS:
(TRM ONLY)

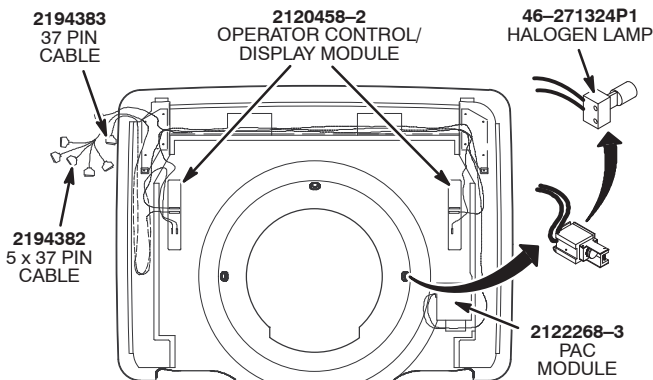
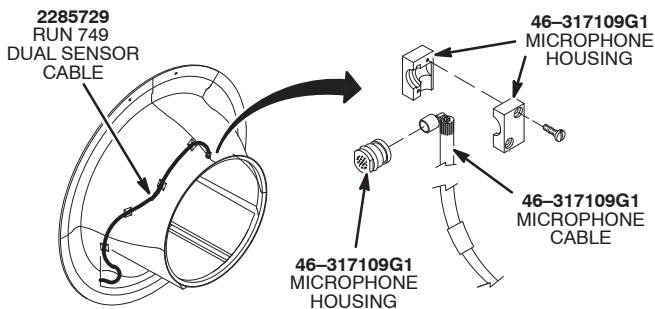
2297540
(4MM)
2297540-2
(6MM)



NOTE: The following parts are Not Shown.

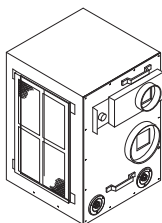
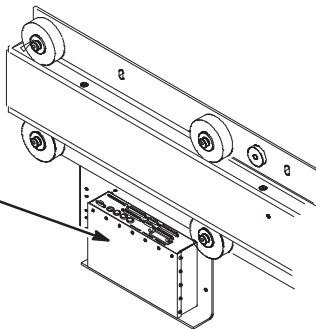
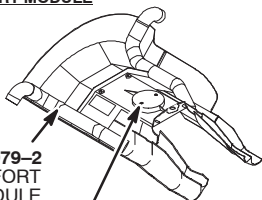
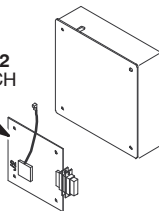
2354050	SPECTRO T/R ASM (OPTION)
2381249	MNS LPCA CABLE HARNESS (OPTION)
2297760	VACUUM GREASE, MSDS# 8366442
2297672	GLYCOL, 1 GALLON/3.8 LITERS

6 MAGNET AND ENCLOSURE (ASSEMBLIES)

3T94 FRONT COVER & PARTS (BACK VIEW):**FRONT END BELL (BACK VIEW):****NOTE: The following parts Not Shown:**

5113881*	3T CRM RF Shield and FRU Container
5109260	3T Eclipse Body Coil FRU (in case)
2319599	3T/94 RT Shim Kit
2370499	3T Eclipse CRM B0 Coil FRU

***2389049** RF Shield Change Out Kit – MUST be ordered for the RF Shield replacement

3T94 COMPONENTS:**BLOWER CABINET****2121199-2**
COMFORT MODULE
BLOWER**5106787**
BODY COIL
BLOWER**SRI MODULE****5121793**
SRI-3 MODULE
5123812
RUN 1260
5123813
RUN 1261
5115910
SRI-2 MODULE**COMFORT MODULE****2220079-2**
COMFORT
MODULE**2206437**
M4 X 12mm NYLON
SCREW FOR
SPEAKER GRILL**MAGNET INTERFACE BOARD****46-226552G2**
LIMIT SWITCH
BOARD

REVISION 11

DIRECTION 5166013

REAR PEDESTAL**5224109**

LONGITUDINAL ENCODER

5250136 or 6250012
REMOTE INTERCOM
BOARD, MG3 A2 A1

46-226552G2
ENCODER LIMIT
SWITCH BD

46-282363P1
FLEXIBLE COUPLING



46-282359P2
GEAR 16 MM

46-282360P1
MOTOR, MG3-A7
(NO LABELS)

5113942
3.0T 8W8
CABLE ASM

5109645-2
3.0T 16 CHANNEL
SWITCH ASM
or
5400021
3.0T DEPOPULATED
16 CH SWITCH ASM
or
5109645-4
3.0T 16 CHANNEL
SWITCH FOR 32
CHANNEL (HDxt)
(NOT SHOWN)

5110428
SWITCH MOUNT PL.

5137348-2
HDxt 3.0T IMAGE
REJECT FILTER

3T94 3.0T HYBRID SPLITTER:

2364795
RUN 893
COAX CABLE

2214448
3T PREAMP

2302305
3T BODY
HYBRID ASM

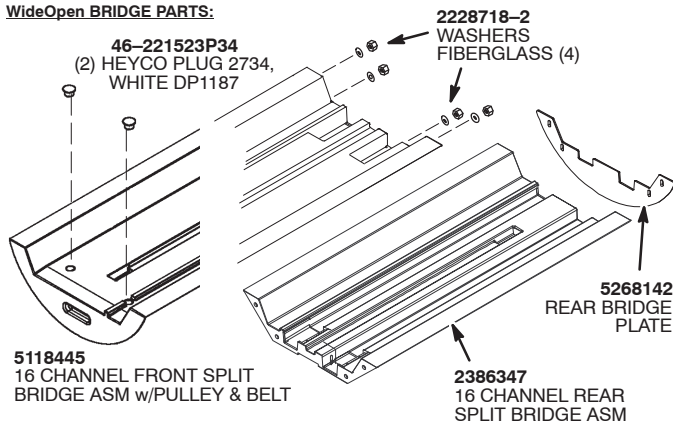
G3 3.0T HYBRID SPLITTER:

2370221
3T DUAL BODY
HYBRID

2369611-2
IMAGE REJECT
FILTER

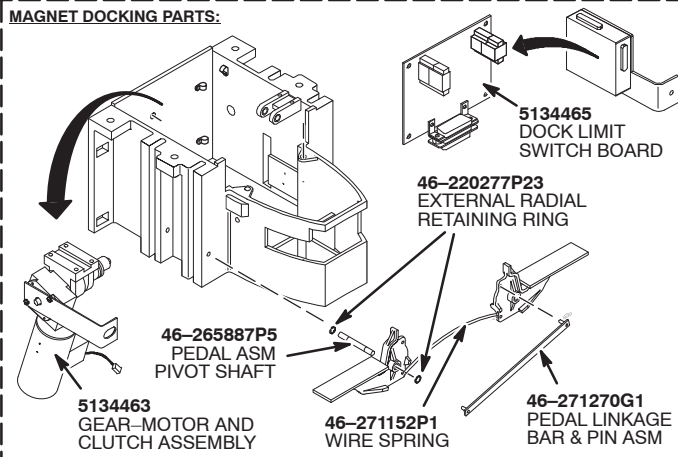
2377703
3T BODY PREAMP
& COMBINER
MODULE

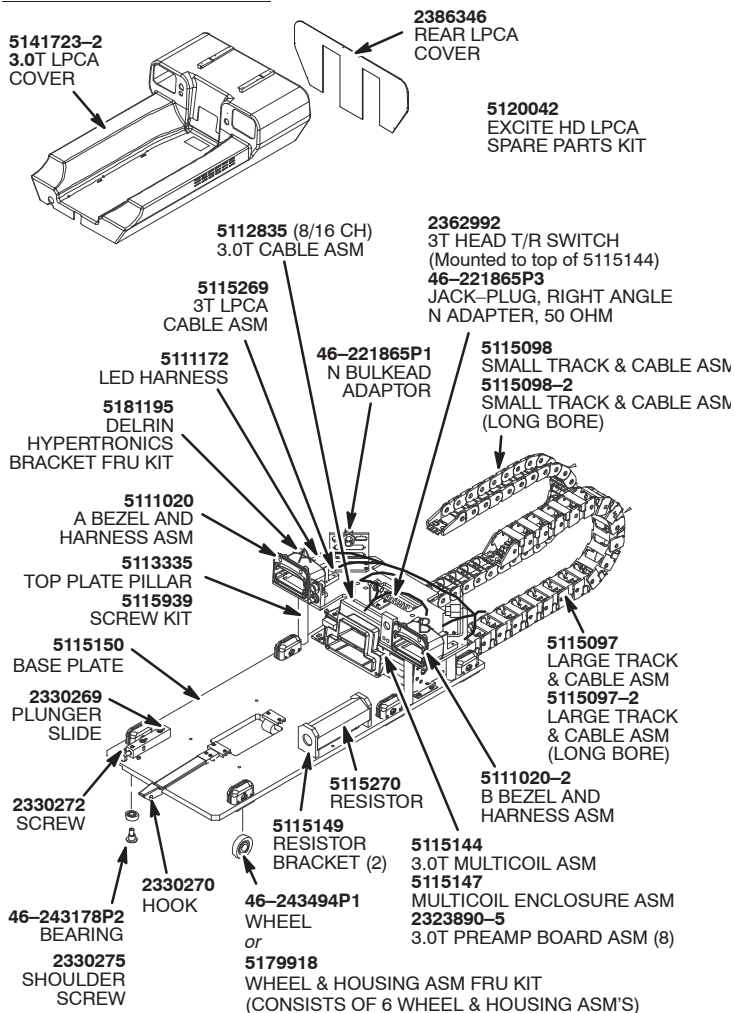
MAGNET AND ENCLOSURE (ASSEMBLIES)

WideOpen BRIDGE PARTS:

NOTE: The following part is Not Shown.

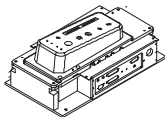
2386348	FRONT BRIDGE FIBERGLASS
5117304	FRONT BRIDGE PULLEY ASM
2293538	BRIDGE EDGE FOAM
46-282436P1	CARRIAGE TRAVEL SENSOR FLAG

MAGNET DOCKING PARTS:

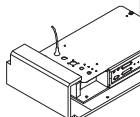
HD LOW PROFILE CARRIAGE ASM**MAGNET AND ENCLOSURE (ASSEMBLIES)**

SRI/PAC PLATFORM:

5112668
NON-LCC
PAC-2B MODULE



5122500
ECG LEAD KIT



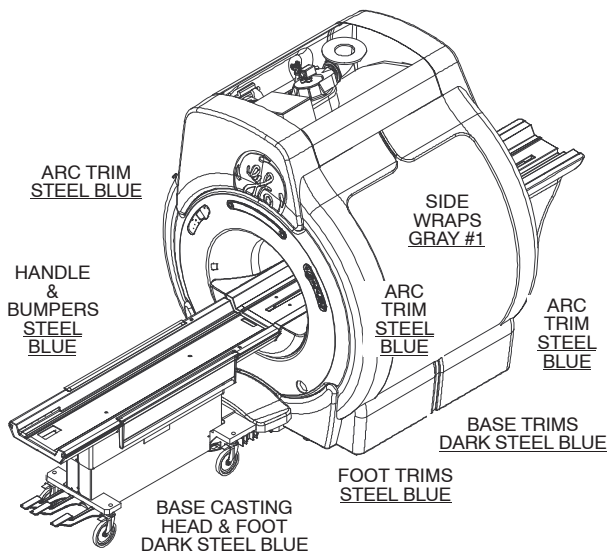
5105270
SRI-3
MODULE

5121793
SRI-3 U/G KIT
w/CABLES

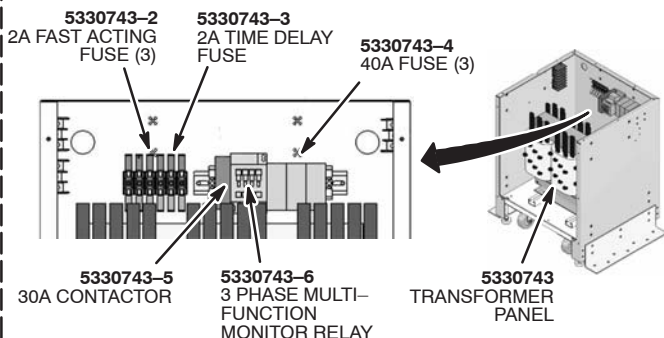
5112668
PAC-2B MODULE

TOUCH-UP PAINTS:

2183491-3	GRAY #1	0.6 OZ. NON-AEROSOL
46-303460P3	STEEL BLUE	0.6 OZ. NON-AEROSOL
46-303460P4	DARK STEEL BLUE	0.6 OZ. NON-AEROSOL



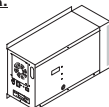
NOTE: To obtain enclosure part numbers not found in this manual, refer to Direction 2138247, *Signa Illustrated Parts List*, for enclosure covers, endbells, and other parts.

CRYOGEN PANEL STEP-DOWN TRANSFORMER FRU'S:**MAGNET RUNDOWN UNIT (MRU) FRU'S:**

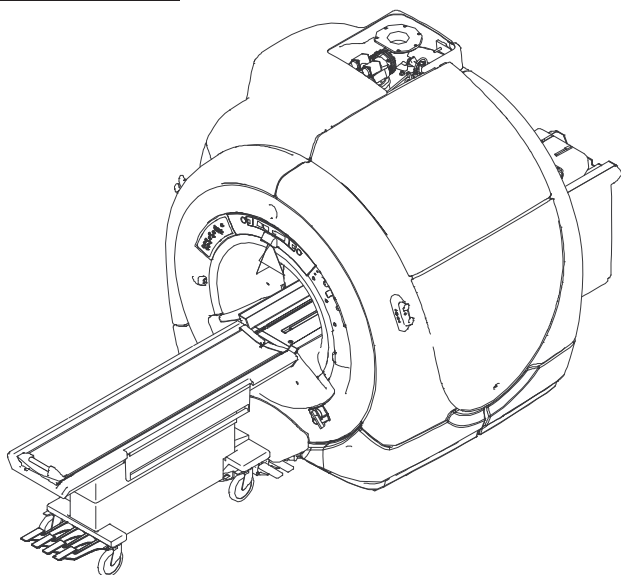
5196918	MAGNET RUNDOWN UNIT (Primary)
5196919	MAGNET RUNDOWN UNIT (Remote)
5363631	SPARE FUSE KIT includes: 1.25A, 250V Fuse; 2.5A, 250V Fuse; 400ma, 120V Fuse; 200ma, 250V fuse
5199490	BATTERY PACK ASM
5180324	ADAPTER PLATE
5196921	REMOTE TO PEN PANEL CABLE
5196920	PRIMARY TO PEN PANEL CABLE
2128699	MAGNET TO MRU CABLE, RUN 606

OXYGEN MONITOR (O2) FRU'S:

5390676	OXYGEN MONITOR FRU – Contains all of the parts included with the O2 monitor, box, sensor and cables
5399042	OXYGEN MONITOR CALIBRATION KIT – Contains all of the parts needed to calibrate either the nw or old O2 monitor
5390678-3	REMOTE 5-YEAR SENSOR – Sensor for the ISA-60M O2 monitor

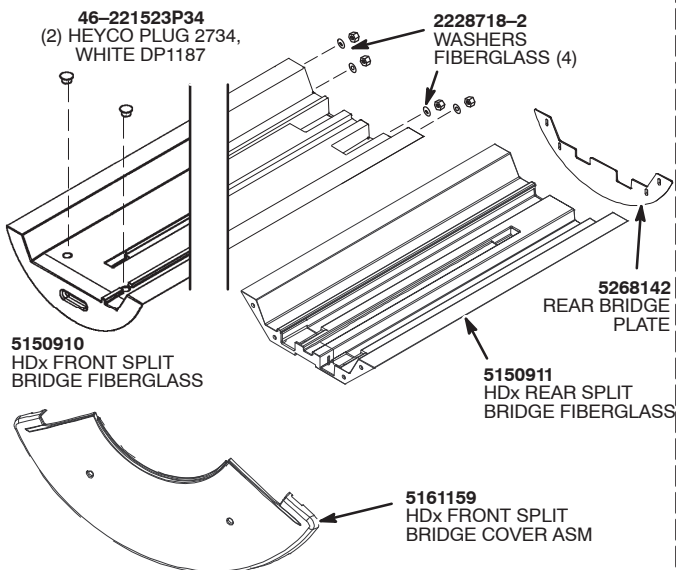
MAGNET MONITOR:

2219400
MAGNET MONITOR (MSM)

HDx MAGNET ENCLOSURE**3.0T HDx ENCLOSURE PARTS**

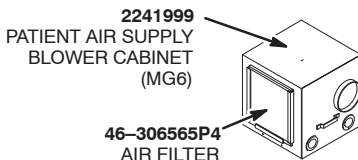
5164786	HDx Enclosure Hardware Spares Kit
5161262	Turtle Bracket Gasket TRM
5150552	Rear Turtle Arm, TRM Single Boss
5150552-2	Rear Turtle Arm, TRM Double Boss
5150552-3	Rear Turtle Arm, TRM Single Boss with Relief Cut
5160997	Front End Bell TRM Assembly
5160999	Filler Cover Left, Front End Bell TRM
5160999-2	Filler Cover Right, Front End Bell TRM
5161000	Rear End Bell TRM Assembly
5161003	Filler Cover Left, Rear End Bell TRM
5161003-2	Filler Cover Right, Rear End Bell TRM
5143732-2	Rear Cover Assembly, TRM
5146696	TRM Comfort Module Assembly
2222549-4	Twin Lumitex Light Panel, 1085mm

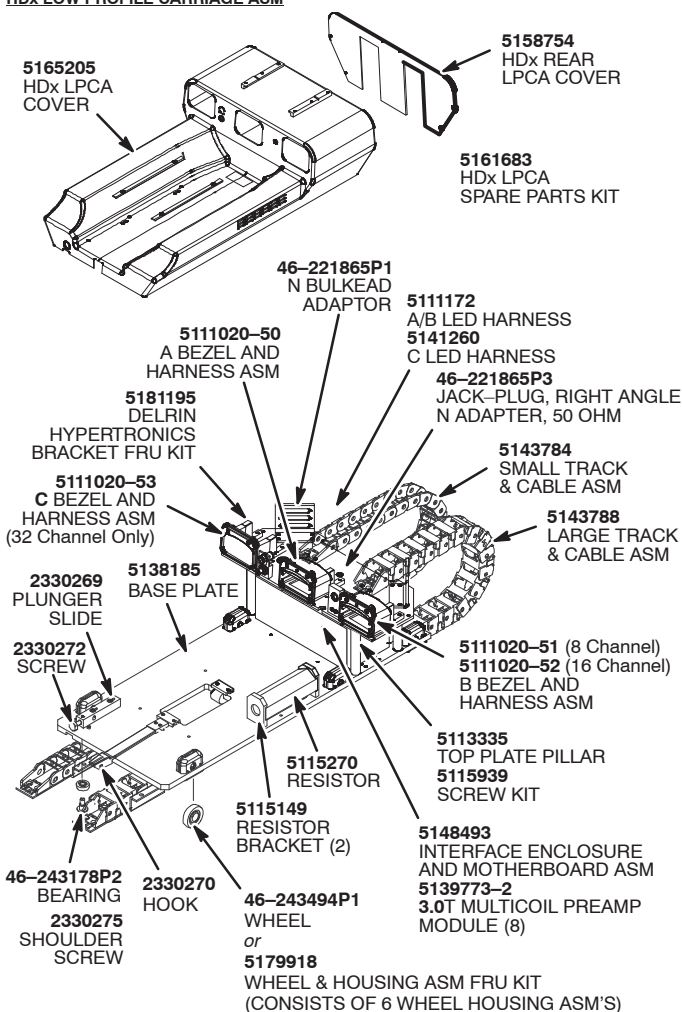
3.0T HDx ENCLOSURE PARTS	
5151721-2	Front Right Arc Cover Assembly
5151722-2	Front Left Arc Cover Assembly
5151723-2	Rear Right Arc Cover Assembly
5151725-2	Rear Left Arc Cover Assembly
5151718-2	Right Side Top Cover, White
5151771-2	Right Side Bottom Cover, White
5161378-2	Right Latch Hook
5151720-2	Left Side Bottom Cover, White
5161378	Left Latch Hook
5145172	Left Side Top Cover, White
5146812	Left Side Top Cover, Rear Corner 3.0T
5151729	Front Left Skirt Cover
5151727	Front Right Skirt Cover
5151731	Rear Right Skirt Cover
5151732	Rear Left Skirt Cover
5143554	Left Rear Pedestal Cover
5143555	Right Rear Pedestal Cover
2215716-51	Center Display Panel With Overlay
2215716-10	Center Display Panel
2215716-60	Left Control Panel with all Overlays
2215716-11	Left Control Panel
2215716-70	Right Control Panel with all Overlays
2215716-12	Right Control Panel
5160276	Bracket I and Q Lines
5150097	Turtle Bracket Front Arm
5160991	Dock Bracket Assembly
5150410	Front Trim Ring Assembly
5150411	Rear Trim Ring Assembly
5143642	Front Cover Assembly
5141500	Microphone Block and Cable Assembly
5142985	Vacuum Port Cover Assembly
5149066	Remote PAC Base
5149068	Remote PAC Inner Mounting Tray
5149065	Remote PAC Accessory Hook
2344509-2	Bore Light End Cap
5160206	Side Cover Mounting Bracket Assembly RH
5160206-2	Side Cover Mounting Bracket Assembly LH
5160992	Rear Bracket Assembly, Left
5160992-2	Rear Bracket Assembly, Right
2206437	M4 x 12mm Nylon Screw for Speaker Grill

HDx BRIDGE PARTS:

NOTE: The following part is **Not Shown**.

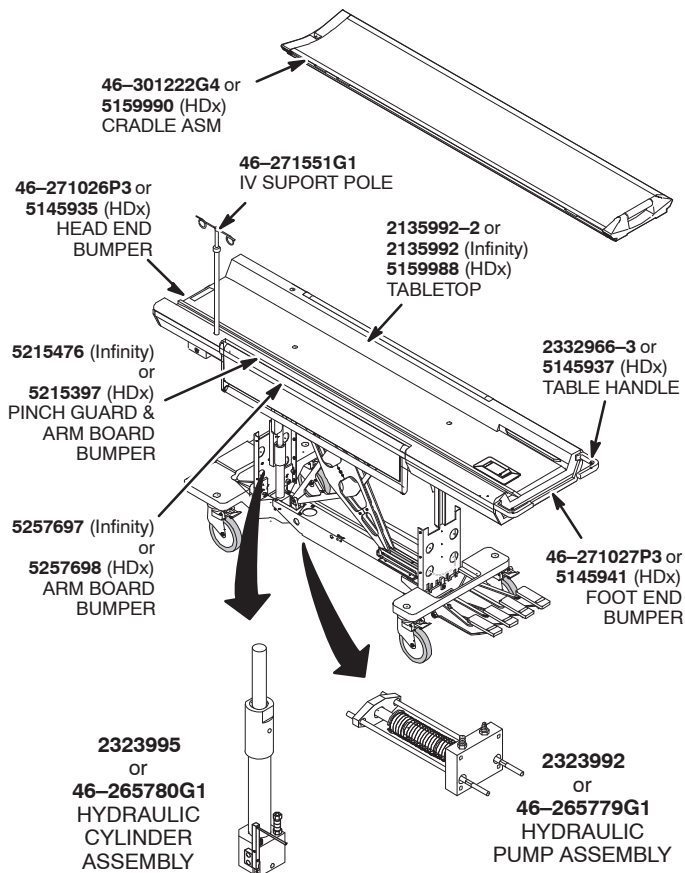
5117304	FRONT BRIDGE PULLEY ASM
2293538-5	BRIDGE EDGE FOAM
46-282436P1	CARRIAGE TRAVEL SENSOR FLAG

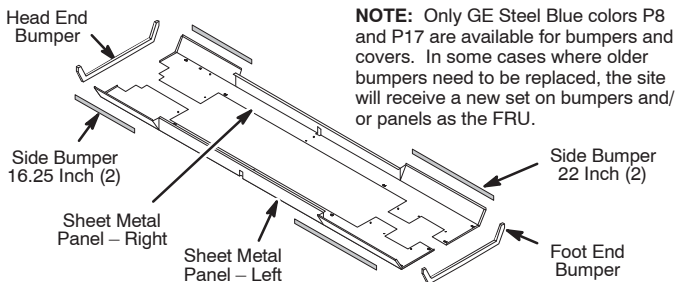
BLOWER CABINET

HDx LOW PROFILE CARRIAGE ASM**MAGNET AND ENCLOSURE (ASSEMBLIES)**

PATIENT TABLE – INFINITY and HDx

NOTE: Refer to page 19 for information on bumper FRU's for HDx, Infinity, and Pre-Infinity patient tables.



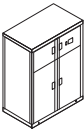
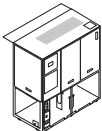
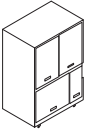
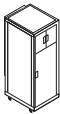
PATIENT TABLE – BUMPERS AND PANELS

HDx EXISTING PART	HDx REPLACEMENT FRU
5145946, Panel – Left 5145947, Panel – Right	5179590 , New panels and side bumpers FRU kit for HDx table
5145939, 22 Side Bumper	5161172 , Side Bumper 22, Color P17
5145939–2, 16.25 Side Bumper	5161172–2 , Side Bumper 16.25, P17
5145941, Foot End Bumper	5145941 , Foot End Bumper
5145935, Head End Bumper	5145935 , Head End Bumper

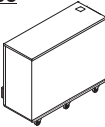
INFINITY EXISTING PART	INFINITY REPLACEMENT FRU
46–271022G1, 2393106, Panel – Left 46–271021G1, 2393106–2, Panel – Right	5179706 , New panels and side bumpers FRU kit for Infinity table
5154908–4, Panel – Left 5154908–3, Panel – Right	5179706 , New panels and side bumpers FRU kit for Infinity table
46–271025P5, 22 Side Bumper	5161172–3 , Side Bumper 22, Color P17
46–271025P6, 16.25 Side Bumper	5161172–4 , Side Bumper 22, Color P17
46–271027P3, Foot End Bumper	46–271027P3 , Foot End Bumper
46–271026P3, Head End Bumper	46–271026P3 , Head End Bumper

Pre-INFINITY EXISTING PART	Pre-INFINITY REPLACEMENT FRU
5154908–4, Panel – Left (Grey) 5154908–3, Panel – Right (Grey)	5179706 , New panels and side bumpers FRU kit for Infinity table
22 Side Bumper (Accent Grey or Blue)	5161172–3 , Side Bumper 22, Color P8
16.25 Side Bumper (Accent Grey or Blue)	5161172–4 , Side Bumper 16.25, Color P8
Foot End Bumper (Accent Grey or Blue)	46–271027P3 , Foot End Bumper
Head End Bumper (Accent Grey or Blue)	46–271026P3 , Head End Bumper

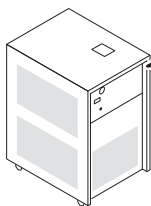
G3 SYSTEM COOLING CABINETS:REFER TO **ELLIS & WATTS** SERVICE MANUALS FOR FRU INFORMATION

OUTDOOR TSCC   REMOTE CONTROL PANEL	INDOOR WATER COOLED TSCC 	INDOOR AIR COOLED TSCC 
OUTDOOR TGWC   REMOTE CONTROL PANEL	INDOOR WATER COOLED TGWC   REMOTE CONTROL PANEL	INDOOR AIR COOLED TGWC 

REFER TO **AIRSYS** SERVICE MANUALS FOR FRU INFORMATION

GWHX 	MRCC   REMOTE CONTROL PANEL
--	---

2343325 AIRSYS HEAT EXCHANGER (GWHX)**2341077-2** AIRSYS QUICK DISCONNECT KIT**2341077-3** AIRSYS REMOTE CONTROL PANEL**2298100-4** E & W MODEL LTL10 HEAT EXCHANGER CABINET**2298100-5** CONTROL PANEL FOR LTL10 ONLY**2297672** GLYCOL USED WITH E & W OR AIRSYS

3T94 GRADIENT WATER CHILLER

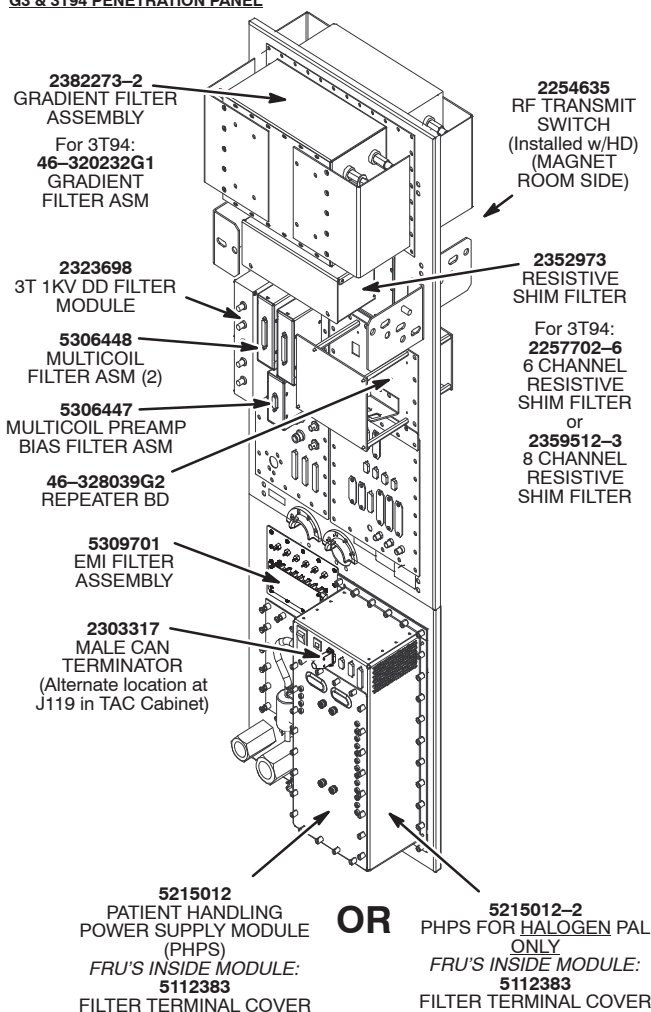
2222564-2
60 Hz GRADIENT CHILLER

OR

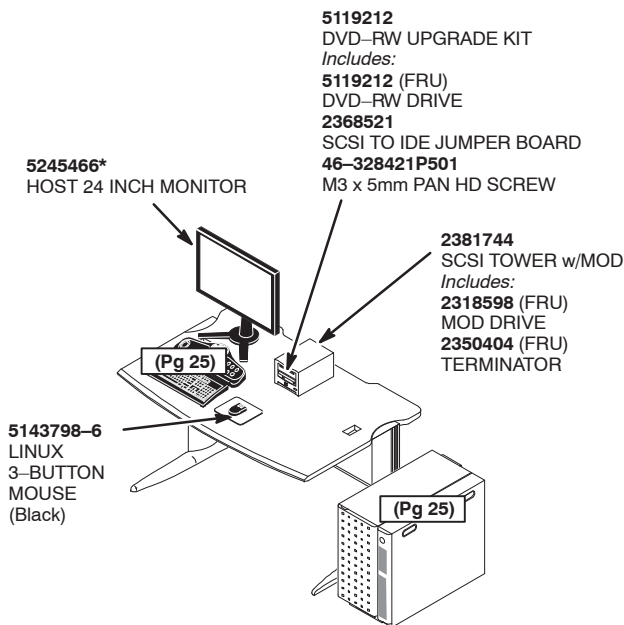
2222564-21
50 Hz GRADIENT CHILLER

NOTE: The following WATER CHILLER parts are Not Shown:

2222564-6	LYTRON 50/60 Hz PUMP MOTOR
2222564-8	LYTRON 0.25A/250V FUSE
2222564-18	LYTRON 60Hz FAN MOTOR
2222564-19	LYTRON 50Hz FAN MOTOR
2222564-31	LYTRON NYLON MOTOR COUPLING
2138791	BODY COIL COOLANT, 1 GALLON

G3 & 3T94 PENETRATION PANEL

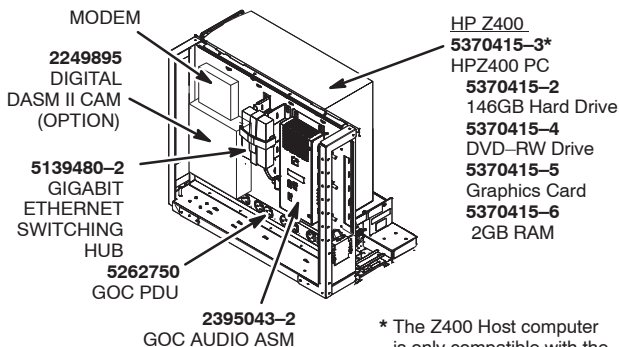
[illegible]

OPERATOR WORKSPACE – HDx

* When ordering a monitor replacement, be aware that models may be updated without notice.

NOTE: The following parts are Not Shown.

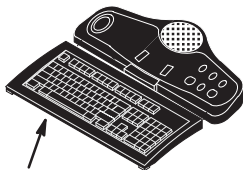
2156158–73	S/W COLLECTOR FOR 15.0 M4
2156158–82	S/W COLLECTOR FOR 15.0 M4A
2156158–84	S/W COLLECTOR FOR 15.0 M4B
2156158–86	S/W COLLECTOR FOR 16.0 M3
2156158–89	S/W COLLECTOR FOR 16.0v2
2156158–91	S/W COLLECTOR FOR HD23.0
2389540	SCAN TOOLS MOD

GOC COMPUTER CABINET WITH HP Z400:

* The Z400 Host computer is only compatible with the 15M4B software or higher.

NOTE: The following parts are Not Shown.

- 5373029** Cable FRU kit for GOC contains:
Run E3016, Run E3046, Run E3047, Run E3048, Host to IRDC DVI-D cable, Host to IRDC USB cable, IRDC Power cable
- 5376172** Dust Filter FRU Kit for Z400
- 5364704** One Wire Collector for Z400

SCIM & KEYBOARD:**NOTE:**

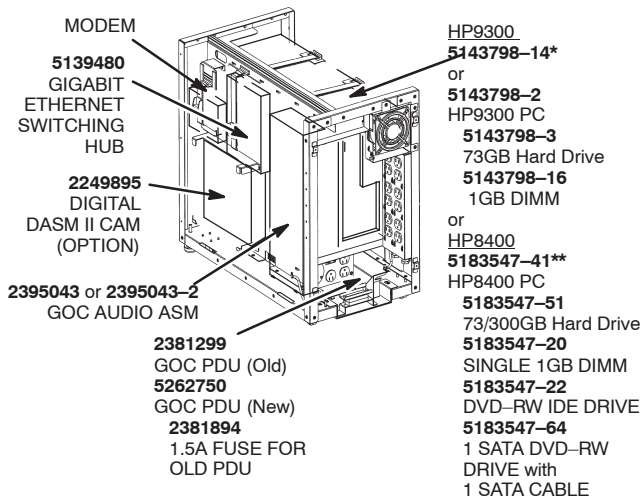
If Keyboard is to be replaced, also order **46-317475P2**, 2 inch x 4 inch adhesive backed velcro to attach to bottom of new keyboard.

KEYBOARDS ONLY:

- 5143798-7** (English)
5143798-8 (French)
5143798-9 (German)
5143798-10 (Spanish)
5143798-11 (Italian)
5143798-12 (Portuguese)
5143798-13 (Scandinavian)

NOTE: The following parts are Not Shown.

- 2262158-9** CABLE, 50 PIN SCIM KEYBOARD
- 5148389-2** SCIM ^W/ENGLISH OVERLAY (Black)
- 5148389-9** MR SCIM GLOBAL OVERLAY SET (Black)

GOC COMPUTER CABINET WITH HP9300 OR HP8400:

* **5143798-14** is only compatible with 14M5 software or later

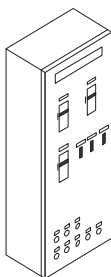
** Sites with SCSI Tower MUST order 5183547-41 as a replacement.
 Sites without a SCSI Tower can order 5183547-41 or 5183547-60 as
 a replacement.

NOTE: The following parts are Not Shown.

5151612 USB CABLE FOR OPERATOR WORKSPACE
5141564 AUTOVOICE CABLE

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

NOTE: Refer to the Directions on applicable MR Service Methods DVD for additional product information.



NOTE: The following items are Not Shown.

5338648 Upgrade kit for MDP Battery and Charger

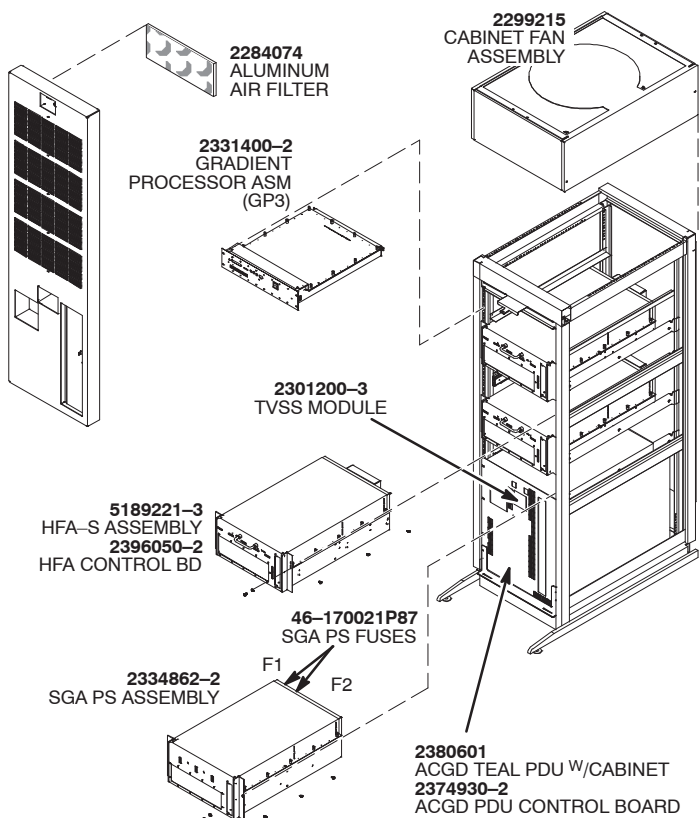
5375947 Power Sonic Battery FRU

5316807 Power Supply FRU

LOCATION	AMPS	VOLTS	GEMS PART NUMBER	MANUFACTURER	
				BUSS	G-S
FUSE F1	1.25	480	2276547-2	SNQ-R-1-1/4	ATQR1-1/4
FUSE F2	3	120	2276547-3	FNQ-R-3	ATQR3
FUSE F3	1	12VDC	2276547-8	LP-CC-1	ATDR1
FUSE F4	0.12	480	2276547-9	SNQ-R-8	ATQR1/8
FUSE F5	125	600	2276547-19	LPJ125SP	AJT125
FUSE F6	50	600	2276547-20	LPJ-50SP	AJT50
FUSE F7	20	600	2276547-22	LPJ-20SP	AJT20
FUSE F8	5	600	2276547-21	FRS-R-5	TRS5R
FUSE F9	8	400	2276547-6	FRS-R-8	TRS8R
FUSE F10	10	120	2276547-11	FNQ-R-10	ATQR10
FUSE F11	10	120	2276547-11	FNQ-R-10	ATQR10

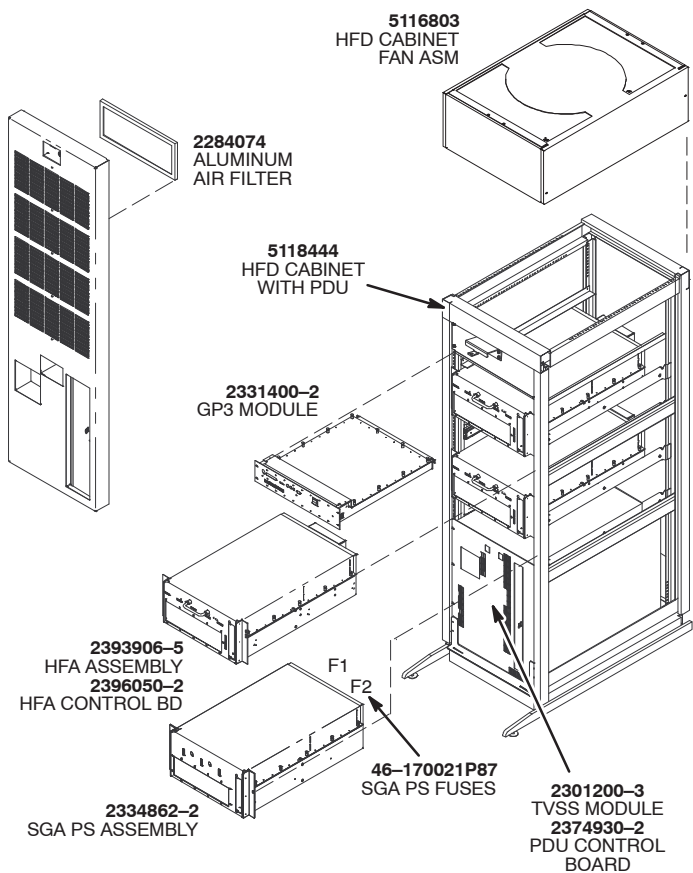
2276547-14 12 VDC BATTERY

2276547-16 12 VDC SUPPLY

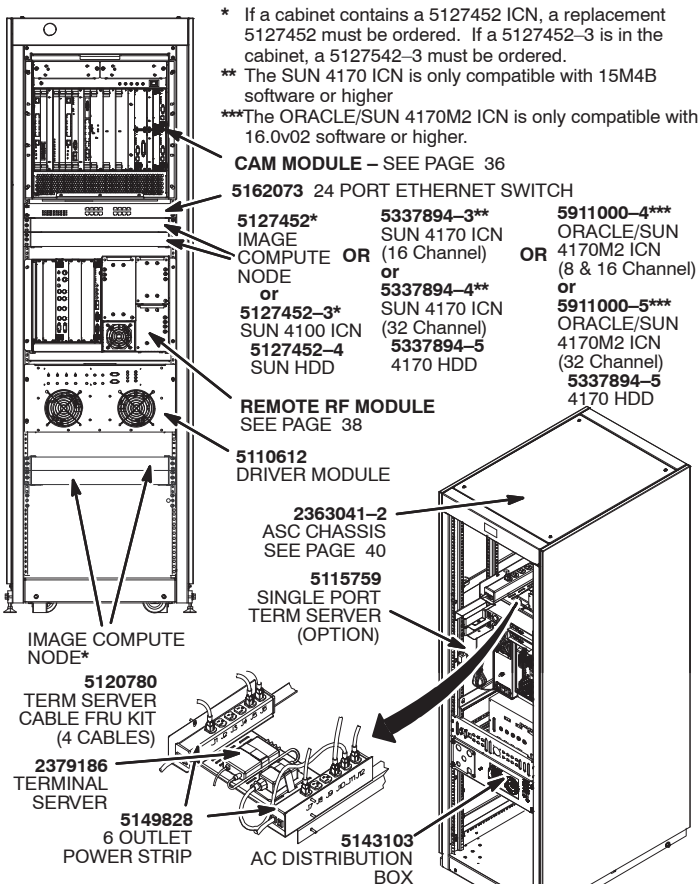
HFD-S/PDU CABINET ASSEMBLY W/SGA/SGAPS/GP

NOTE: The following Part is Not Shown.

2303317 CAN Terminator

HFD/PDU CABINET

[illegible]

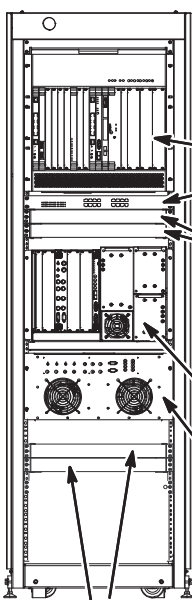
HDx RFS CABINET with CAM Module

NOTE: The following items are Not Shown.

5272613 NMR SERVICE CABLE KIT

46-296817P150 OHM BNC TERMINATOR

5306523 DVMR SERVICE CABLE KIT (Supports IRF2b, RRF-DIF3 & later
(Order if DVMR Loopback adaptor is needed; Supports HD 16.0v02 or later)

HDx RFS CABINET

* If a cabinet contains a 5127452 ICN, a replacement 5127452 must be ordered. If a 5127452-3 is in the cabinet, a 5127542-3 must be ordered.

** The SUN 4170 ICN is only compatible with 15M4B software or higher

***The ORACLE/SUN 4170M2 ICN is only compatible with 16.0v02 software or higher.

MGD MODULE – SEE PAGE 37

5162073 24 PORT ETHERNET SWITCH

5127452*

IMAGE
COMPUTE
NODE
OR

5127452-3*
SUN 4100 ICN
5127452-4
SUN HDD

5337894-3**

SUN 4170 ICN
(16 Channel)

OR

5337894-4**
SUN 4170 ICN
(32 Channel)
5337894-5
4170 HDD

5911000-4***

ORACLE/SUN
4170M2 ICN
(8 & 16 Channel)

OR
5911000-5***
ORACLE/SUN
4170M2 ICN
(32 Channel)
5337894-5
4170 HDD

REMOTE RF MODULE
SEE PAGE 38

5110612
DRIVER MODULE

2363041-2
ASC CHASSIS
SEE PAGE 40

5125176
UTNS
SPLITTER
BOX

IMAGE COMPUTE
NODE *

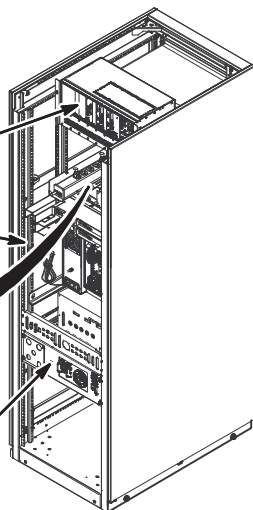
5120780
TERM SERVER
CABLE FRU KIT
(4 CABLES)

2379186
TERMINAL
SERVER

5115759
SINGLE PORT
TERM SERVER
(OPTION)

5149828
6 OUTLET
POWER STRIP

5143103
AC DISTRIBUTION
BOX



NOTE: The following items are Not Shown.

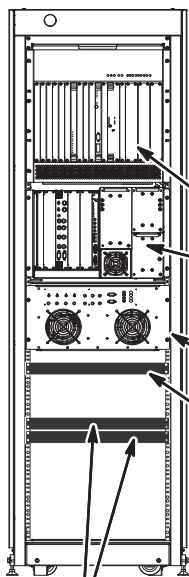
5272613 NMR SERVICE CABLE KIT

46-296817P150 OHM BNC TERMINATOR

5306523 DVMR SERVICE CABLE KIT (Supports IRF2b, RRF-DIF3 & later
(Order if DVMR Loopback adaptor is needed; Supports HD 16.0v02 or later)

RFS CABINET

33

HD RFS CABINET UPGRADED TO HDx

- * If a cabinet contains a 5127452 ICN, a replacement 5127452 must be ordered. If a 5127452-3 is in the cabinet, a 5127542-3 must be ordered.
- ** The SUN 4170 ICN is only compatible with 15M4B software or higher
- ***The ORACLE/SUN 4170M2 ICN is only compatible with 16.0v02 software or higher.

MGD MODULE

SEE PAGE 37

REMOTE RF MODULE

SEE PAGE 38

**5110612
DRIVER
MODULE****5162073
24 PORT
ETHERNET
SWITCH****5111201
AC 120V
DISTRIBUTION
MODULE****5115474
UTNS SPLITTER BOX****2363041-2
ASC CHASSIS
SEE PAGE 40****2379186
TERMINAL
SERVER****2272740
ETHERNET HUB
SERVICE KIT****5120780
TERM SERVER
CABLE FRU KIT
(4 CABLES)****5115759
SINGLE PORT
TERM SERVER
(OPTION)****5139347
AC
DISTRIBUTION
BOX****5127452***
IMAGE
COMPUTE
NODE

or

5127452-3*
SUN 4100 ICN
5127452-4
SUN HDD

OR

5337894-3**
SUN 4170 ICN
(16 Channel)
5337894-5
4170 HDD

or

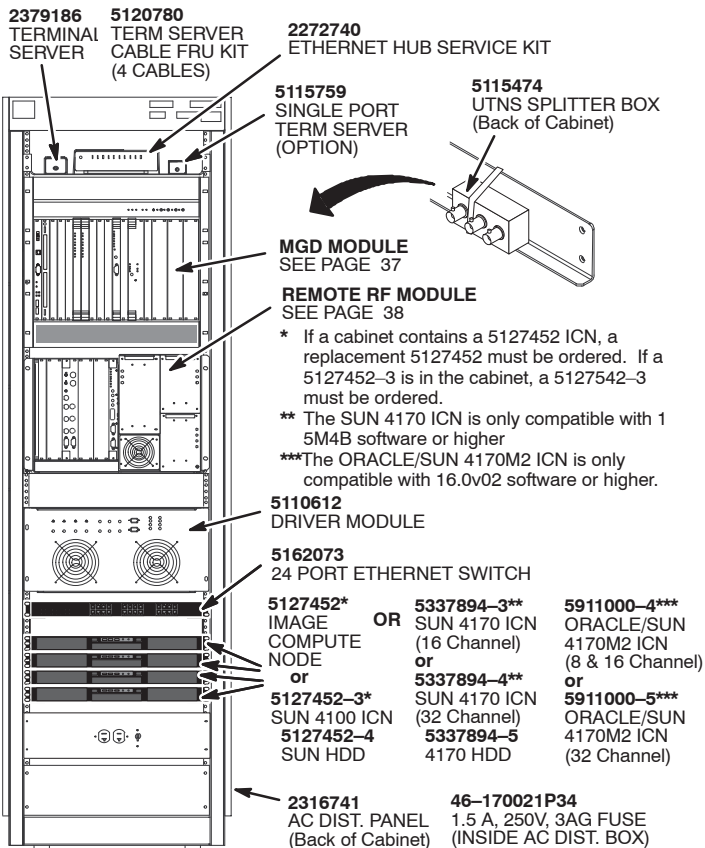
5337894-4**
SUN 4170 ICN
(32 Channel)
5337894-5
4170 HDD

OR

5911000-4***
ORACLE/SUN
4170M2 ICN
(8 & 16 Channel)
5337894-5
4170 HDD

or

5911000-5***
ORACLE/SUN
4170M2 ICN
(32 Channel)
5337894-5
4170 HDD**NOTE: The following items are Not Shown.****5272613** NMR SERVICE CABLE KIT**46-296817P150** OHM BNC TERMINATOR**5306523** DVMR SERVICE CABLE KIT (Supports IRF2b, RRF-DIF3 & later
(Order if DVMR Loopback adaptor is needed; Supports HD 16.0v02 or later)

EXISTING SYSTEM CABINET UPGRADED TO HDx

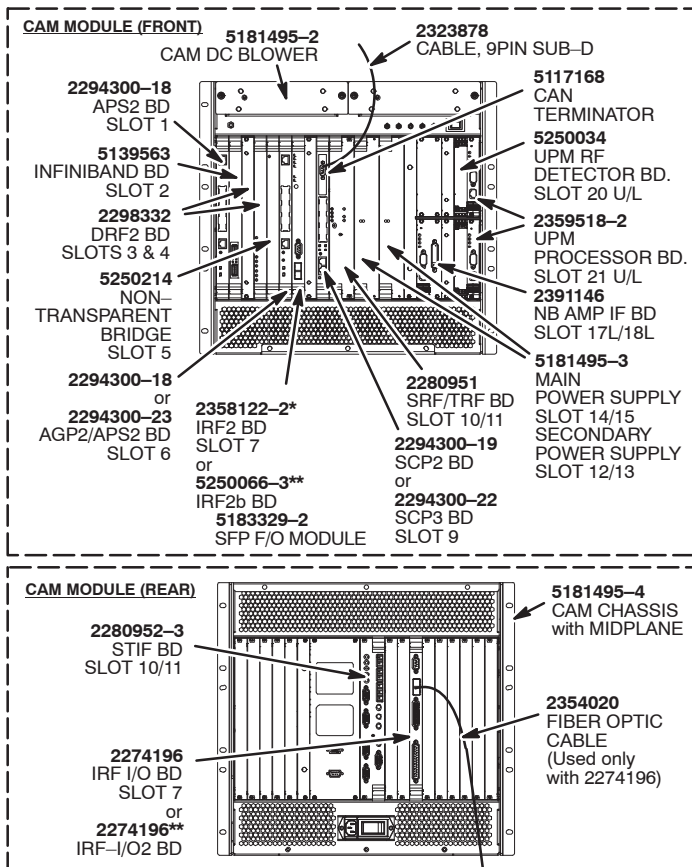
NOTE: The following items are **Not Shown**.

5272613 NMR SERVICE CABLE KIT

46-296817P150 OHM BNC TERMINATOR

5306523 DVMR SERVICE CABLE KIT (Supports IRF2b, RRF-DIF3 & later
(Order if DVMR Loopback adaptor is needed; Supports HD 16.0v02 or later)

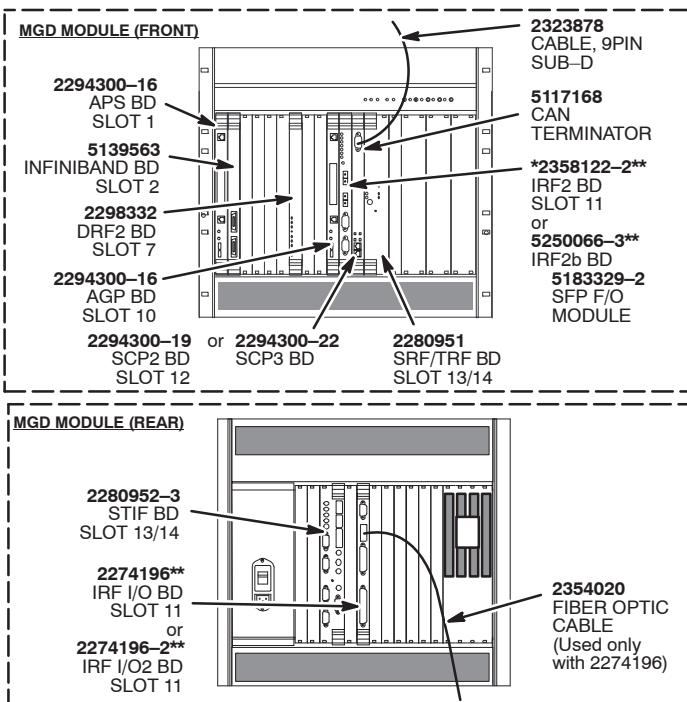
SYSTEM CABINET



* For 32 Channel systems, use **2358122** IRF2 Bd.

** Replace the following boards only with the same P/N board:
2358122-2, 5250066-3, 2274196, and 2274196-2.

NOTE: In the MGD, CAM, and RRF Chassis, the RRF-DIF3 (6250184), IRF2b (5250066-3), and IRF-I/O2 (2274196-2) Boards only work with each other in the system. Do not intermingle other P/N boards with these. Also, they are only compatible with 16.0v02 software and higher.



* For 32 Channel systems, use **2358122** IRF2 Bd.

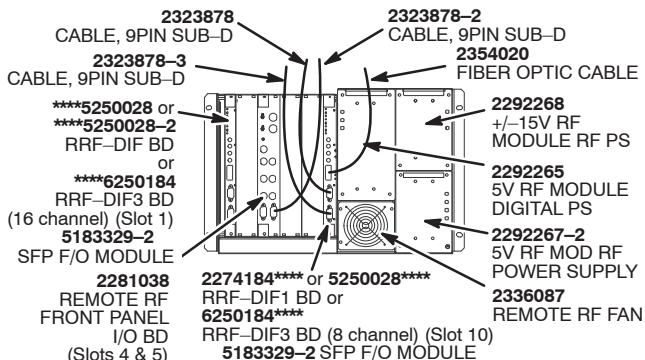
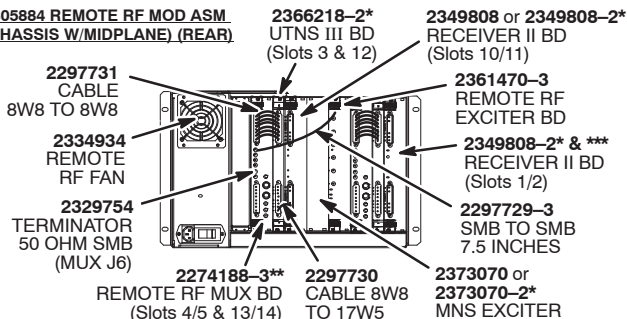
** Replace the following boards only with the same P/N board:
2358122-2, 5250066-3, 2274196, and 2274196-2.

NOTE: In the MGD, CAM, and RRF Chassis, the RRF-DIF3 (6250184), IRF2b (5250066-3), and IRF-I/O2 (2274196-2) Boards only work with each other in the system. Do not intermingle other P/N boards with these. Also, they are only compatible with 16.0v02 software and higher.

NOTE: The following MGD Module parts are **Not Shown**.

2294300-5	BRIDGE BOARD
2294300-6	MGD CHASSIS W/BACKPLANE & AC ENTRY (No Boards)
2294300-7	MGD FAN
2294300-8 or	
5181495-3	POWER SUPPLIES
2294300-9	INDICATOR BOARD

MGD MODULE

2305884 REMOTE RF MOD ASM (CHASSIS W/MIDPLANE) (FRONT)**2305884 REMOTE RF MOD ASM (CHASSIS W/MIDPLANE) (REAR)**

NOTE: In the MGD, CAM, and RRF Chassis, the RRF-DIF3 (6250184), IRF2b (5250066-3), and IRF-I/O2 (2274196-2) Boards only work with each other in the system. Do not intermingle other P/N boards with these. Also, they are only compatible with 16.0v02 software and higher.

NOTE: Additional REMOTE RF MOD ASM Parts/Cables are located on page 39.

* A part that does not have a "-2" can be replaced with either the same part number or the same number ending with a "-2". If replacing the "-2" version, then only the part number with a "-2" can be used.

** 2274188 and 2274188-2 are replaced by 2274188-3. Also, all 2274188 series MUX boards can be replaced by backwards compatible MUX2 board, 5250106.

*** For 32 Channel systems, use 5250144 Exciter 2.1, in main RFS Cabinet only (not in 32 Channel cabinet)

**** Replace the following boards only with the same P/N board:
2274184, 5250028, 5250028-2, and 6250184.

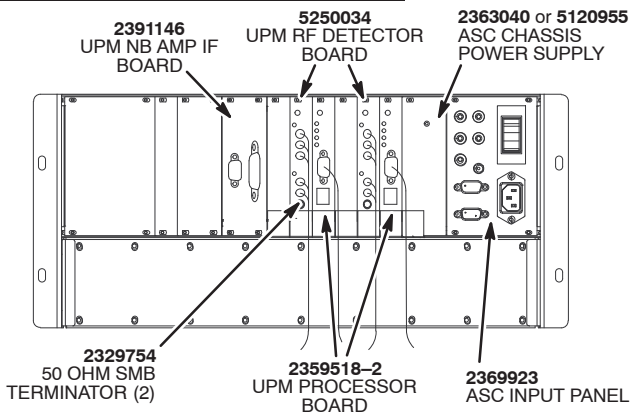
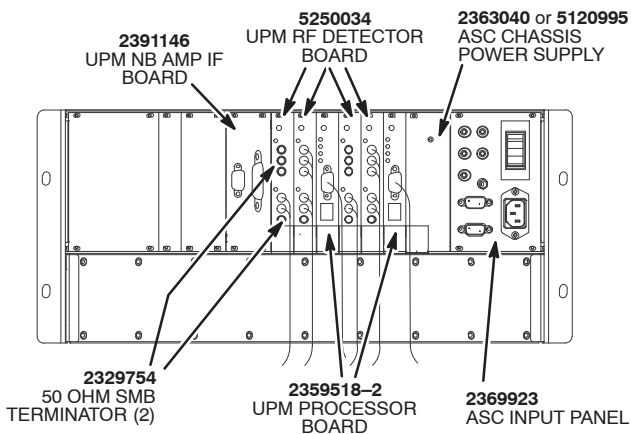
NOTE: The following RRF MODULE (page 38) Parts/Cables are Not Shown.

46-306565P1	CABINET FILTER
2294728	COMPACT PCI FRONT PANEL BLANK (8HP) (2 SLOT)
2319152	COMPACT PCI FRONT PANEL BLANK (4HP) (1 SLOT)
2308729	BNC TO SMB CABLE (MUX J8 TO CAB I/F J1)
2308728	8W8 CABLE ASSEMBLY (MUX J5)
2322578	J15/23/25/26 DATA INPUT CABLE (CAB I/F)
2322580	CABLE J18, UNBLANK (CAB I/F)
2322581	CABLE J24, MAGNET MONITOR (CAB I/F)
2322582	CABLE J30/31, SERIAL EXP (CAB I/F)
2322583	CABLE J62 (15 PIN SUB-D) (CAB I/F)
2322586	CABLE J61 (37 PIN SUB-D) (CAB I/F)
2308729-2	BNC TO SMB CABLE (MUXJ13 TO CAB I/F J2)
2308729-3	BNC TO SMB CABLE (MUXJ14 TO CAB I/F J3)
2394870-3	RRF slot14R-J5 I/O J200
5110796	RRF slot14R-J13 I/O J202
5110795	RRF slot 14R-J13 I/O J202
5115417-3	RRF slot14R-J8 RF Amp J14
5115417	RRF slot12R-J31 Splitter-1
5113193	RRF slot12R-J7 RRF slot11R-J6/7
2394870-2	RRF slot5R-J5 I/O J201
5113193-2	RRF slot3R-J7 RRF slot2R-J6/7
5115417-2	RRF slot3F-J31 Splitter-2
5115630	J7 MUX (2) TO J6 MUX (1) CABLE
5115630-2	RRF slot14R-J11 RRF slot9R-J6
5116213-2	RF Extension cable to Run 1238
5121269	RF Extension cable to Run 1259
5145157	Fiber Optic cable, MGD Slot 11R-J20 to RF Slot 10F-J21
5145157-2	Fiber Optic cable, MGD Slot 11F-J7 to RF Slot 1F-J21
5150527	Y Power cable, 10 feet
5160775	HDx, FRU ICN Ethernet Cable Collector

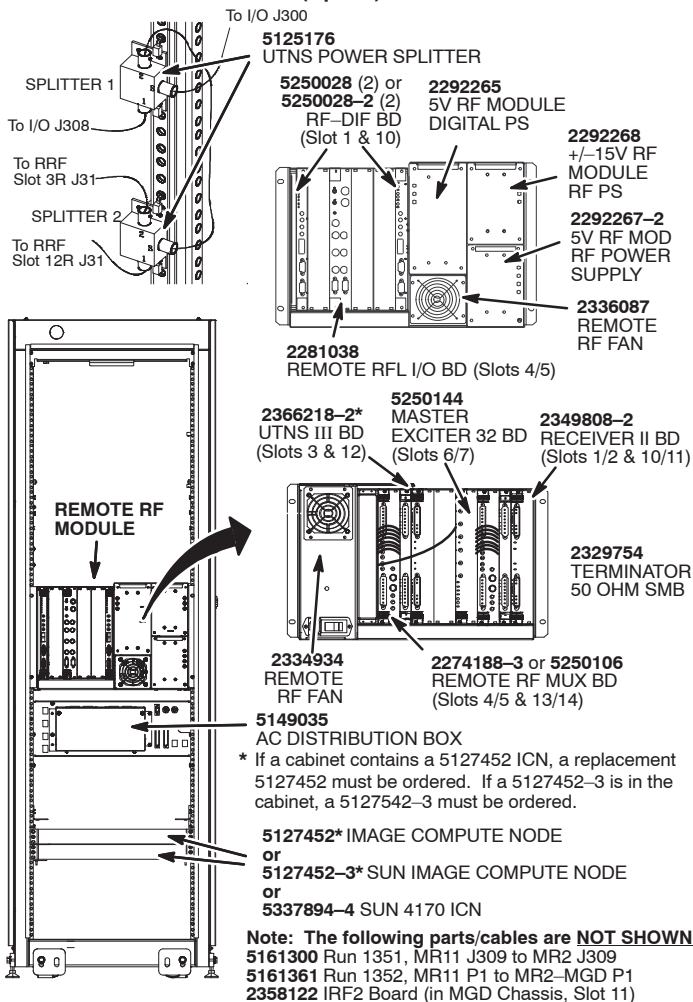
5429284	Fiber Optic Cable FRU Kit for RRF-DIF3, IRF2b & IRF-I/O2
	contains:
	5257686-2, MGD Slot 11F J10 to RF Slot 10F J21
	5416142-2, MGD Slot 11F J8 to RF Slot 10F CLK
	5257686-3, MGD Slot 11F J11 to RF Slot 1F J21

NOTE: The following ASC MODULE (page 40) Parts/Cables are Not Shown.

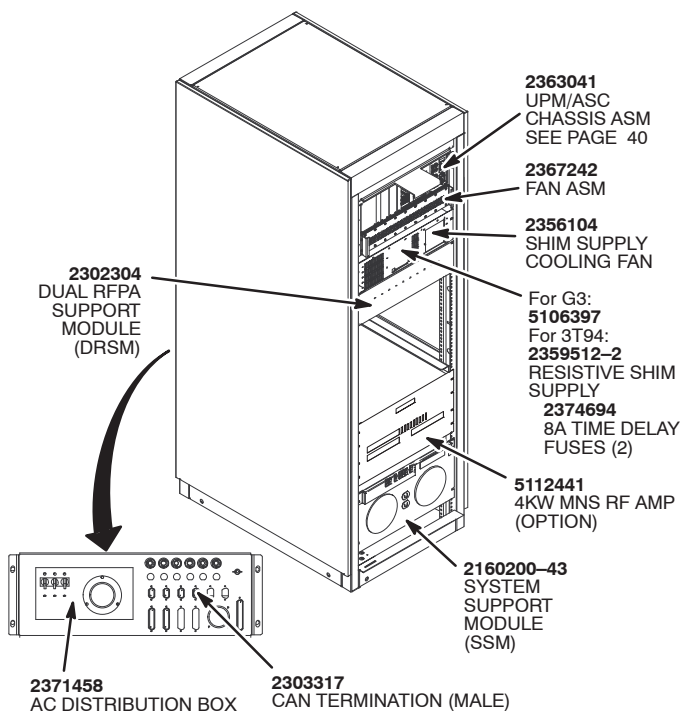
5115417-4	ASC slot14-J1 RF Amp J8
5115417-5	ASC slot14-J4 RF Amp J7
5115417-6	ASC slot11-J1 RF Amp J6
5120955	ASC POWER SUPPLY
2391148	BB2 AMP I/F ASM

2363041 or 2363041-2* ASC CHASSIS ASSEMBLY**2363041 or 2363041-2* ASC CHASSIS ASSEMBLY WITH MNS**

NOTE: Additional ASC MODULE Parts/Cables are located on page 39.

HDx 32 CHANNEL CABINET (Option)

[illegible]



NOTE: The following RF ACCESSORY CABINET parts Not Shown:

46-306565P1 10" FRONT COVER FILTER

2365797 CABLE – I/O J127 TO ASC SLOT 11/14 J1/2/4/5

2368552 CABLE – ASC SLOT 14–J3 TO I/O J100

2369552-8 CABLE – ASC SLOT 13–J4 TO I/O J102

2369552-9 CABLE – ASC SLOT 13–J5 TO I/O J103

2369552-10 CABLE – ASC SLOT 10–J4 TO I/O J104

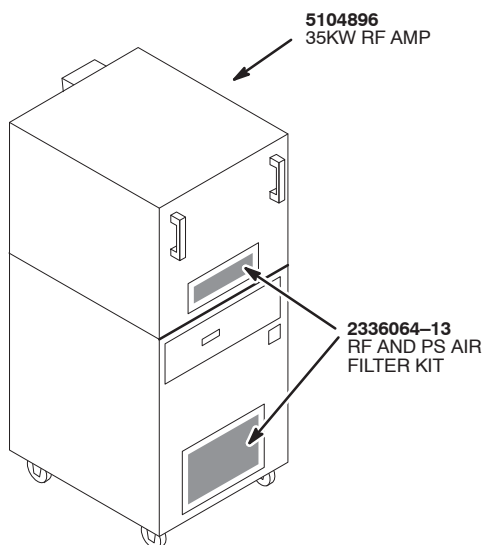
2369552-11 CABLE – ASC SLOT 10–J5 TO I/O J105

2379507 DRSM JUMPER PLUG

2383066 DRSM JUMPER CABLE (ASC–DRSM–I/O)

This page for G3 MAGNET (Short Bore)

35KW NARROW BAND RF AMP



NOTE: The following NB RF AMP CABINET parts are **Not Shown**:

2336064-8 35KW RF AMP FUSE KIT

2336064-15 35KW RF AMP TEMPERATURE SENSOR

NOTE: The following RF TRANSMIT FLEX CABLES are **Not Shown**:

2384850 RUN 1156 – MR8-J4 TO CONNECT TO RUN 1125

2384850-2 RUN 1160 – MG3-A3-J1 TO CONNECT TO RUN 1127

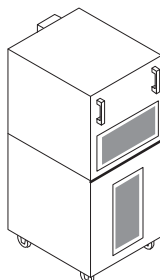
2384851 RUN 1157 – MR8-J3 TO CONNECT TO RUN 1126

2384851-2 RUN 1158 – MR1-CPC-J2 TO CONNECT TO RUN 1137

2384851-3 RUN 1159 – MR9-J4 TO CONNECT TO RUN 1146

This page for 3T94 MAGNET (Long Bore)*

25kW NB RF AMPLIFIER CABINET (for sites WITHOUT LCD Display)



2278810-4
NB RF AMP
CABINET

NOTE: The following 25kW NB RF AMP CABINET (without LCD) parts are Not Shown:

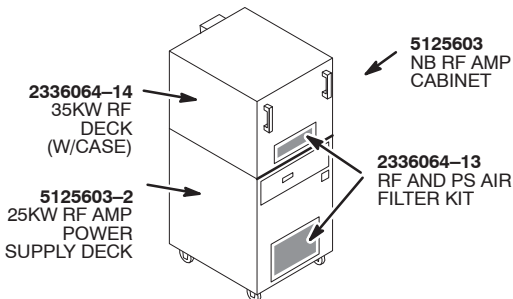
2336064-4 35KW RF AMP IPA TUBE & TUBE PULLER KIT

2336064-5 35KW RF AMP PA TUBE & HEX WRENCH KIT

5111003 55-50dB DUAL DIRECTIONAL BODY RF COUPLER FOR E2

5111004 45-40dB DUAL DIRECTIONAL HEAD RF COUPLER FOR E2

25/35kW NB RF AMPLIFIER CABINET (for sites WITH LCD Display)



2336064-14
35KW RF
DECK
(W/CASE)

5125603-2
25KW RF AMP
POWER
SUPPLY DECK

5125603
NB RF AMP
CABINET

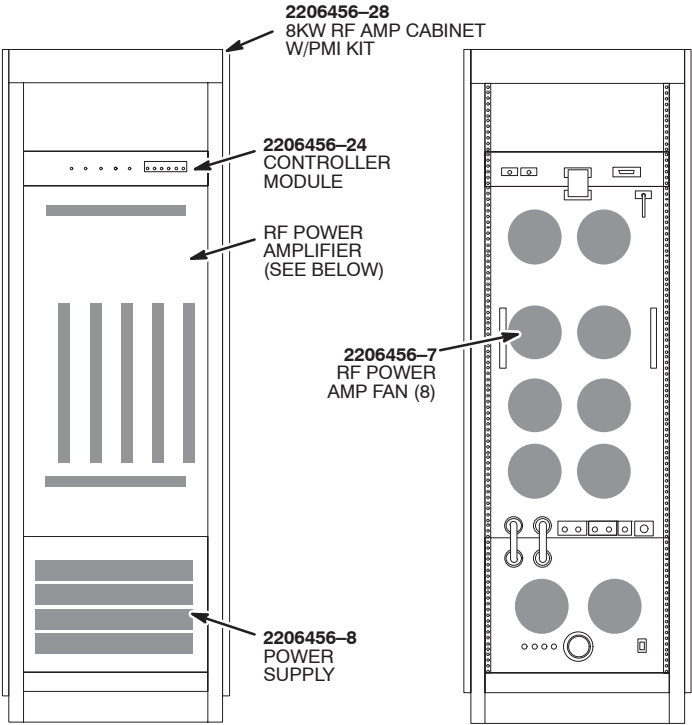
2336064-13
RF AND PS AIR
FILTER KIT

NOTE: The following 25/35kW NB RF AMP CABINET parts are Not Shown:

2336064-8 35KW RF AMP FUSE KIT

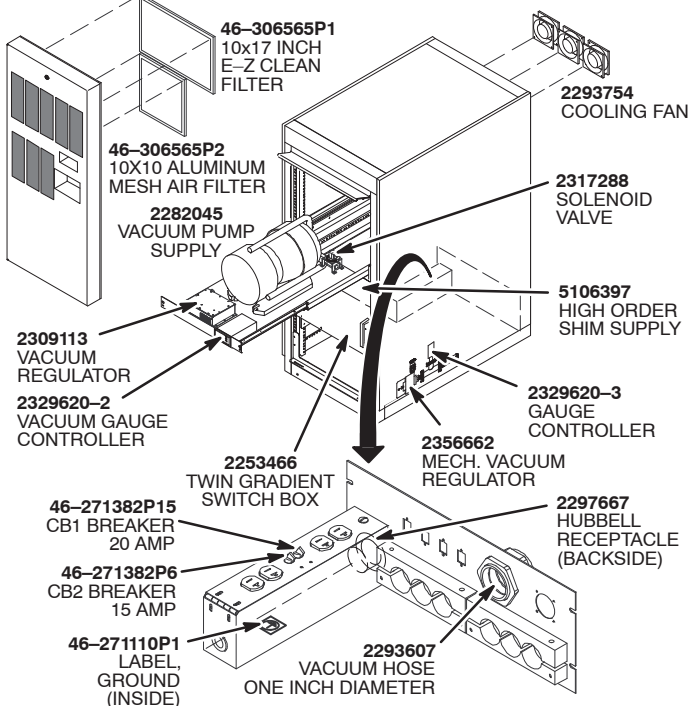
2336064-15 35KW RF AMP TEMPERATURE SENSOR

* **Note:** If the site has a 3T94 Long Bore magnet, it MUST use a 25KW NB RF Amplifier on this page.



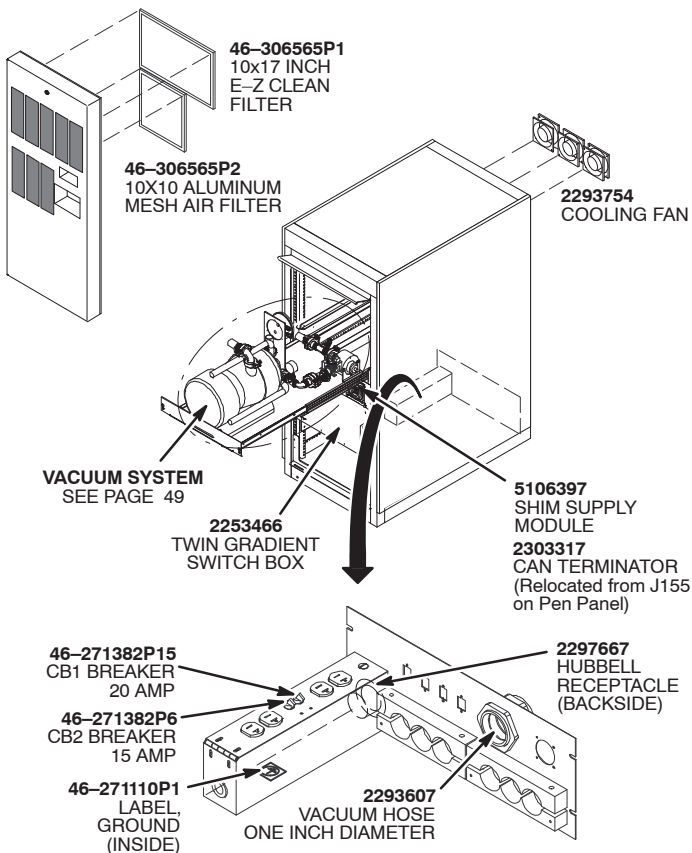
NOTE: The following RF POWER AMP internal parts are **Not Shown**:

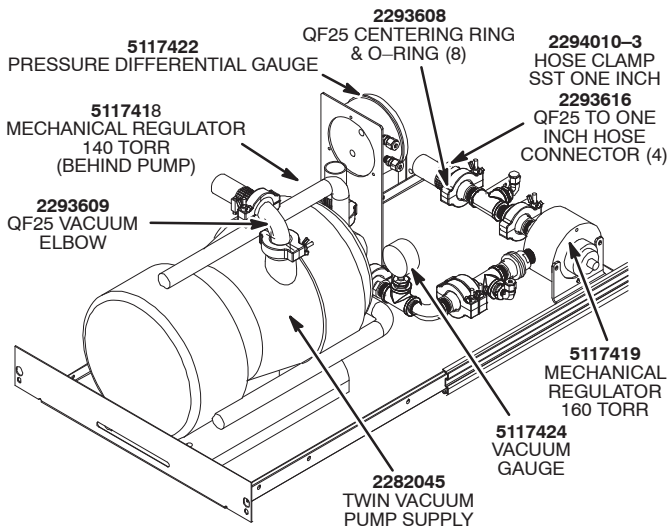
2206456-2	RF POWER AMP PREAMP/SPLITTER
2206456-3	RF AMPLIFIER (5)
2206456-4	RF POWER AMP COMBINER/COUPLER
2206456-5	RF POWER AMP COUPLER (4)
2206456-6	RF POWER AMP BIAS TEE
2206456-13	CABLE – SPLITTER TO AMPLIFIER (5)
2206456-20	COUPLER-BIAS TEE
2206456-21	POWER BIAS TEE
2206456-22	PMI INTERFACE KIT
46-170021P29	5 AMP FUSES (3)

2239500 TAC

NOTE: The following parts are Not Shown.

2282045-2	VACUUM PUMP TIP SEAL KIT
2293608	QF25 CENTERING RING & O-RING
2293608-2	QF16 CENTERING RING & O-RING
2293609	QF25 VACUUM ELBOW
2293616	QF25 ADAPTER TO ONE INCH HOSE
2293617	VACUUM REDUCING 'T'
2293618	QF16 CLAMP
2293618-2	QF25 CHAIN CLAMP
2294010-3	HOSE CLAMP, SST ONE INCH
2297760	VACUUM GREASE, MSDS# 8366442
2329620-4	VACUUM GAUGE CONTROLLER CABLE
2284301	VACUUM GAUGE TO TAC J2 CABLE

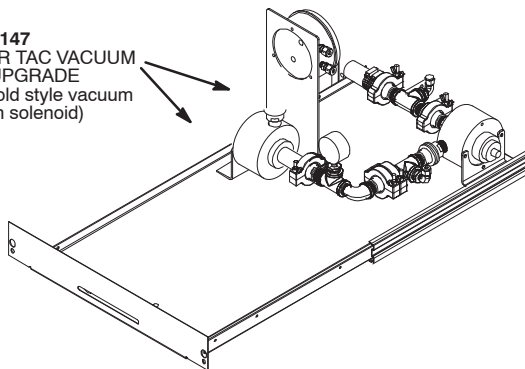
5118162 TAC



NOTE: The following parts are Not Shown.

2282045-2	VACUUM PUMP TIP SEAL KIT
5118818	CLAMP NW 20-25
2297760	VACUUM GREASE, MSDS# 8366442

5127147
SERVICE KIT FOR TAC VACUUM
SYSTEM UPGRADE
(Used to replace old style vacuum
system, with solenoid)



TAC VACUUM SYSTEM

[illegible]

The table below provides information for generic replacement cables. If any of the cables are ordered, make sure to order write-on labels to be used for identifying cable part number and connection designations at the end of each cable.

GENERIC REPLACEMENT CABLES

Description	Length	
	65 Feet	100 Feet
9 Pin Sub-D Male/Female	46-244714G30	46-244714G37
9 Pin Sub-D Female/Female	46-244713G30	46-244713G37
15 Pin Sub-D Male/Female	46-244712G30	46-244712G37
15 Pin Sub-D Female/Female	46-244711G30	46-244711G37
25 Pin Sub-D Male/Female	46-244710G30	46-244710G37
25 Pin Sub-D Female/Female	46-244709G30	46-244709G37
37 Pin Sub-D Male/Female	46-258877G30	46-258877G37
	50 Feet	65 Feet
BNC RG 223	46-251710G1	46-251710G18
Cable Write-on Labels (60 labels per booklet)	46-208781P3	

DESIGNATED CABLE REPLACEMENTS

CABLE PART #	CABLE RUN
46-244710G30	Run 1227 Run 1236
46-244712G30	Run 1188
46-244713G30	Run 1181
46-244714G30	Run 1180 Run 1181 Run 1233
46-258877G30	Run 1189 Run 1190

CABLE KITS

FRU Kit Description	Part #
Receive Cables for Equipment Room, 34.5 ft (10.56m)	5119953
Receive Cables for Equipment Room, 65.6 ft (20m)	5119953-2
Receive Cables for Magnet Room, 28.2 ft (8.6m)	5119954
Receive Cables for Magnet Room, 49.2 ft (15m)	5119954-2
Transportable/Relocatable MC Rcv Cable Kit (Eqmt. Rm)	5119953-3
Body Receive Cables (both rooms), 34.5 ft (10.52m)	5119956
Body Receive Cables (both rooms), 65.6 ft (20m)	5119956-2
Term Server Cables (4 cables)	5120780

FIBER OPTIC CABLE KIT – FIXED SITE

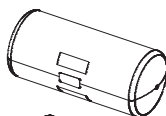
Description	Part #
Includes clamps, cable ties, mounts, and Runs 1045, 1083, & 1179.	5115713

TIMES MICROWAVE CABLE PARTS

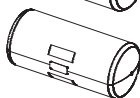
Description	Part #
Straight Male N Connector for LMR 600	2352195
7/8 EIA Right Angle Plug for LMR 1200 & Center Pin	2384859

VRE CABLE KITS

Description	Part #
Infiniband cables – Infiniband to MGD and ICN	5160774
Ethernet cables – ICN to Ethernet, Runs 1371 & 1345	5160775

UNIFIED SURFACE COIL PHANTOM SET

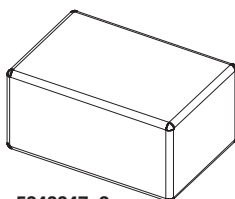
5342679-2
Large Cylindrical
Unified Phantom



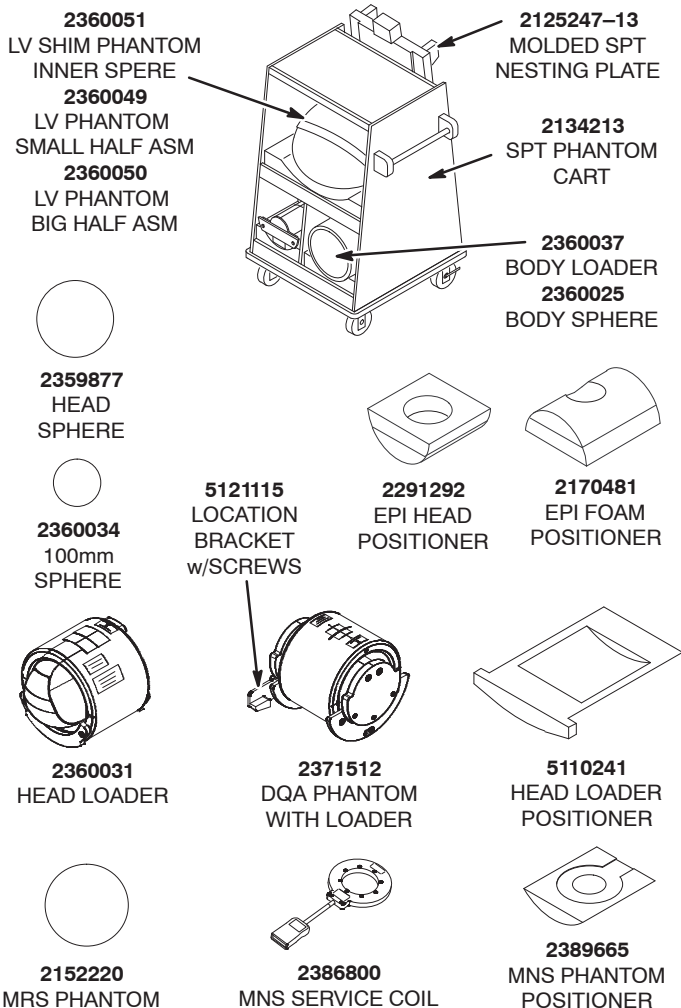
5342680
Small Cylindrical
Unified Phantom



5342681
Cubical
Unified Phantom



5343347-2
TL Unified Phantom

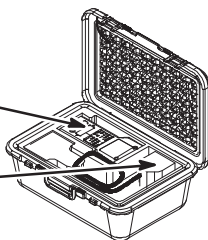


MCR3 TOOL:

For use with systems that have Hypertronics connectors only in LPCA.

5146813
3.0T MCR3
ENCLOSURE TOOL

5146208
MCR3 CABLE ASM
w/CONNECTOR

**RF TRANSMIT TERMINATION KIT:**

5131924 1.5T AND 3.0T RF TRANSMIT CABLE SERVICE KIT

GRAFIDY KIT:

2270695
CuSO4 SAMPLE VIAL

HDx TEMPERATURE AND SENSOR KIT:

5147136-2 ENVIRONMENTAL REPAIR KIT

HEAD COILS:**HDx SINGLE CHANNEL SURFACE COIL ADAPTOR**

5145746 SURFACE COIL SINGLE CHANNEL ASSEMBLY

3T SPLIT TOP HEAD COIL

2376114 SPLIT HEAD COIL ASM

5191687 HEAD COIL LATCHING FRU KIT

2145391 MIRROR ASM

2376104 QD BOX ASM

2355094 REPLACEMENT SPRING CONTACTS KIT

5140887 HDx HEAD COIL QHTR ASSEMBLY

3T94 HEAD COIL

46-282118G8 3T/94 HEAD COIL

2376104-2 HEAD QUICK DISCONNECT CABLE FOR ECLIPSE 2

G3 & 3T94 HDx COILS:**3T 8 CHANNEL BODY BREAST BIOPSY COIL M3335LA**

2417367	BREAST COIL w/o CABLE
2417366	OUTPUT CABLE ASM
2417234	PHANTOM ASM
2414541	BIOPSY MOUNT
2412285	BIOPSY ARM
2412281	BREAST COIL HEADREST

3T 8 CHANNEL TORSO COIL (HD CONNECTOR) M3335LN

2418063	ANTERIOR COIL
2418064	POSTERIOR COIL
2416958	CABLE ASM

3T 8 CHANNEL CTL COIL (HD CONNECTOR) M3335LM

2417540	COIL ASM
2417542	CABLE ASM
2414538	EAR PAD SET
2381682-9	CTL CT PHANTOM
2381682-10	TL PHANTOM

3T 3 CHANNEL SHOULDER COIL M3335LR

2416878	COIL ASM
2416877	CABLE ASM
2375136-5	PATIENT PAD
2375136-3	PHANTOM POSITIONER

3T NV ARRAY COIL

2414383	E3 NV ARRAY W/O CABLE
2414385	E3 NV SYSTEM COIL
2414386	ANTERIOR CHEST WITH CHEST CABLE
2414389	NV PHANTOM SET
2414391	E2 NV SYSTEM CABLE

3T BRAIN/SPINE COIL M3335LC (FP) & M3335LD (U/G)

2417351	HEAD NECK UNIT (A & P)
2418075	HNU ADAPTOR BLOCK
2417350	NECK CHEST UNIT
2417349	TL UNIT
2417352	OUTPUT CABLE – HYPERTRONICS (UPGRADE)
2417353	OUTPUT CABLE – BENDIX (UPGRADE)
2417355	OUTPUT CABLE (FORWARD PRODUCTION)
2417354	TL UNIT CABLE
2418080	MIRROR

1.5T BRAIN/SPINE ACCESSORIES

2418077	TL PHANTOM KIT
2418079	HN PHANTOM KIT
2418081	PAD KIT
5344844	PHANTOM POSITIONER, HNS-UNIFIED

COILS CONTINUED ON NEXT PAGE.

G3 & 3T94 HD COILS:**3T 8 CHANNEL BRAIN COIL**

2380637-2	BRAIN COIL/MIRROR/BASEPLATE ASM
------------------	---------------------------------

2317112-3	PHANTOM POSITIONER
------------------	--------------------

2380637-5	CABLE ASM
------------------	-----------

5146634-6	MIRROR ASM
------------------	------------

3T QUAD EXTREMITY COIL

2380635-2	KNEE/FOOT ASM W/CABLE
------------------	-----------------------

2380635-5	KNEE/FOOT CABLE ASM
------------------	---------------------

2357819-10	EXTREMITY REPAIR KIT
-------------------	----------------------

2357819-11	PHANTOM POSITIONER
-------------------	--------------------

2357819-12	PHANTOM
-------------------	---------

3T GP FLEX COIL

2382817	GP FLEX COIL
----------------	--------------

2389979	GP FLEX COIL CABLE KIT
----------------	------------------------

G3 HD COILS:**3T 8 CHANNEL TORSO ARRAY COIL**

2381683-2	COIL WITH CABLE
------------------	-----------------

2381683-3	PHANTOM
------------------	---------

2382683-4	PHANTOM POSITIONER
------------------	--------------------

3T CTL ARRAY COIL

2413107-1	CTL ARRAY W/CABLE
------------------	-------------------

2413107-4	CTL ARRAY CABLE
------------------	-----------------

2381682-9	CTL PHANTOM CT SECTION
------------------	------------------------

2381682-10	CTL PHANTOM TL SECTION
-------------------	------------------------

3T94 HD COILS:**3.0T 4 CHANNEL TORSO ARRAY COIL**

2366942	3T TORSO ARRAY ASSEMBLY
----------------	-------------------------

5110897 or -2	3T/94 EXCITE 4 CHANNEL TORSO COIL CABLE
----------------------	---

3.0T NV ARRAY COIL

2383478-2	NV ARRAY W/COIL
------------------	-----------------

2383477-3	NV COIL PHANTOM
------------------	-----------------

2383477-4	NV COIL PHANTOM POSITIONER
------------------	----------------------------

2383478-5	NV CABLE ASM
------------------	--------------

3.0T 4 CHANNEL CTL COIL

2380634-2	3T/94 CTL ARRAY
------------------	-----------------

2377326-3	USA1 8 CHANNEL COIL PHANTOM
------------------	-----------------------------

5110896 or -2	3T/94 EXCITE 4 CHANNEL CTL ARRAY CABLE
----------------------	--

3.0T 8 CHANNEL CTL ARRAY COIL

2413107-1	CTL ARRAY W/CABLE
------------------	-------------------

2413107-4	CTL OUTPUT CABLE
------------------	------------------

COILS CONTINUED ON NEXT PAGE.

3.0T 8 CHANNEL WRIST COIL		M3340CA
5268024-2	3T HD 8 CH WRIST COIL W/ BASEPLATE & CABLE	
5268024-3	3T HD 8 CH WRIST ARRAY COIL CABLE	
3.0T 8 CHANNEL FOOT ANKLE COIL		M3340CB
5273235-2	3T HD 8CH FOOT ANKLE COIL W/ CABLE. ELEC ASSY	
5273235-3	3T HD 8CH FOOT ANKLE COIL CABLE	
5273235-5	HD 8 CH FAC PHANTOM ASSEMBLY	
5273235-6	3T HD 8 CH FOOT ANKLE COIL BASEPLATE ASSY	
3.0T 8 CHANNEL SHOULDER COIL		M3340CC
5271265-2	3T HD 8 CH SHOULDER COIL ASSY WITH CABLE	
5271265-3	HD 8 CH SHOULDER COIL CABLE ASSEMBLY	
5271265-5	8 CH SHOULDER COIL PHANTOM HOLDER ASSY	

BrainWave OPTIONS:

2363727	15" LCD Display with 5 foot VGA Cable for Computer connection
2394034-5	(Cedrus) Response Pad Asm. Left & right hand units with two buttons for each hand. Includes fiber optic cable (20 ft) and optical converter box
2394034-6	(Cedrus) "Lumina" response controller box
2394034-7	(Cedrus) AC Power module for Lumina Controller Box
2394034-9	FRU usage only Kit containing all cables [Cedrus-4, -5, -6, -7, 2xxxxxx-6, 46-265505P1 (qty 2)]
2394034-10	Penetration Sub-Panel Asm
2394034-16	PC compatible Keyboard (narrow)
2394034-17	PC compatible Mouse
5160492	BrainWave 2.0 Application Software CD-ROM
5173371-2	PC and Platform Software Asm (For HDx Only)
5110968	Patch CD for BrainWave Stimulus PC

MNS OPTIONS:

Description	Part #
3.0T 13C TR Switch	2354050-6
3.0T 13C MNS Preampl	5114799-2
MNS QA Probe 3.0T 13C	5114403-3
MNS QA Probe 3.0T 31P	5114978-4
3.0T 13C Sample Vial	5114978-3
3.0T 31P Sample Vial	5114978-4

OPTIONS:

CADSTREAM 4.1 w/SURELOC **M3335KL**
5164701 CONFIRMA 4.1 SERVER PC & S/W w/MONITOR & MOUSE
5160894 CADSTREAM S/W & DOCS CD-ROM
5154402 CADSTREAM 3.X TO 4.1 UPGRADE

CONFIRMA MAMMO CAD **M3033NF**
5117875 COMPLETE PC WITH SOFTWARE

BAR CODE READER KIT **M1090PM**
2180690-2 READER ONLY
2180690-3 Y CABLE ONLY

ReportCARD HARDWARE
5117866-26 ReportCARD COMPUTER
5160890 ReportCARD SOFTWARE COLLECTOR
2390369 ReportCARD OPTON KEY

FUNCTOOL
5160880 FUNCTOOL 4.0 APPLICATIONS CD

EXCITE SOFTWARE OPTION KEYS:

2156158-64	Release 15.x Software Collector HDxt (M3)
2156158-65	Release 14.x Software Collector HDx (M5)
2383494	Fiesta C
2390037	Asset
2390040	ConnectPro
2390075	3D FRFSE
2390078	EchoPlus
2390079	Diffusion Tensor
2390084	Probe 2000
2390085	Probe 200 Upgrade
2390086	Sage 7
2390091	2D Delayed Enhancement
2390095	iDrive ProPlus
2390097	SmartPrep 2000
2390098	Flouro Trigger MRA
2390103	Performed Proc Step
5110228	Single Voxel Probe-P
5110229	Single Voxel Probe-S and Probe-P
5110230	Probe 2D – CSI
5110231	Probe 3D – CSI
5110232	MNS, SAGE, Probe-P
5110319	Research Proton Spectro
5120544	Release 12.0 M4 Apps DVD
5121548	ReportCard on AW
5121550	Tractography on AW
5122363	LAVA
5122364	PFO Evaluation
5122365	PROPELLOR Brain Suite
5122366	Tractography
5122367	2D FatSat FIESTA
5122368	3D Flair
5122370	RealTime FIESTA
5122371	VIBRANT-II
5122372	8 Channel HDx
5122373	16 Channel HDx
5122375	3.0T HDx ScanTools 14.0

SOFTWARE OPTIONS

REVISION HISTORY

REV/DATE REASON FOR CHANGE

- 1 Mar 03, 06 Initial release
- 2 Jun 26, 06 Added HDx new color RF Coil. Updated ethernet switch part number to 5162073 for 24 port.
- 3 Jan 26, 07 Added RF Shield Change Out Kit, Speaker Grill nylon screw, LPCA Wheel Housing FRU kit (MRIhc21558), & Dock pedal FRU parts (MRIhc17629). Corrected Linux 3-button mouse part number. For RFS cabinet, added Term Server cable FRU kit for and RRF Module MUX2 board. Removed Infiniband Switch and MCR2 Tool.
- 4 Mar 13, 07 Added front cover laser, 5 gallon coolant pail, and LPCA LED cable. Added note for referring to 2138247 for parts not in this manual. Updated Ethernet cable info. Added center pin to description of Times Microwave right angle connector (CQA 13099582). Removed ECMT kit (CQA 13097585). Added Head Coil Latch Kit. Updated BrainWave s/w asm part # to 5173371.
- 5 Jun 7, 07 Removed erroneous LED cable from LPCA. Updated APS & APG Module parts numbers in MGD to 2294300-16.
- 6 Sep 17, 07 Moved Revision History to last page. Updated HFA-S Asm part # in HFD-S cabinet (ECO 20428110). Added RFS Cabinet with CAM Module and CAM Module board details. Revised MGD & ASC Module part #'s. Added LPCA bracket, pillar, and screws FRU kit (PQR 13121160), Patient Table bumper and panel FRU kit (PQR 13127363), Grafidy CuSO4 vial (PQR 13098795), and Brain Spine accessories. Added NV coil FRU parts (PQR 13113173)
- 7 May 28, 08 Changed name of manual to *Signa HDx 3.0T Series FRU's*. Rear Pedestal – updated Longitudinal Encoder part number (PQR 13175398) and added Encoder board part # 46-226552G2. Added Lytron pump coupler part # (PQR 13136000). Corrected 25kw NB RF Amp part # (PQR13142203). HFD: Updated FHA Control Bd part #. RFS Cab: Added Sun ICN part # and replacement note. Updated SCP2 Bd part #. Added 5181495-3 MGD P.S. NB RF Amp: Corrected MKS 25kw Amp part # (PQR 13142203). Tools: Added RF Transmit Termination kit. Coils: Added Mirror assembly to Brain Coil (ECR 203882). Options: Updated BrainWave software part #.
- 8 Sep 23, 08 Added HDxt parts to OW, coils, and software sections. Updated Remote Intercom Board part number. Added EMI Filter Asm to Pen Panel (ECO 2060894). Updated PHPS FRU's on Pen Panel (ECO 2040374). Added Linux PC module for 14M5 s/w. Added Sun hard drive FRU. Updated part # for Environmental Repair Kit.
- 9 Jun 15, 09 Added Flexible Coupling (PQR 13206931) and Depopulated 16-CH Switch to Rear Pedestal (ECO 2063089). Updated Filter FRU's on PEN Panel (ECO 2084029). Updated Monitor offerings on OW (ECO 2047560). Also added HP9300 DIMM (ECR 2078642), HP8400 DIMM (ECR 2038518) & DVD-RW (ECR 2059074). In RRF Chassis, added 5250028-2 RF-DIF (ECO 2080729).
- 10 Nov 3, 09 Updated OW LCD monitor to 5343103 (ECO 2083909) and added new GOC AA 2395043-2 (ECO 2084000). Added 2294300-22 SCP3 Bd. to CAM and MGD modules (ECO 2064203, MRIhc42663). Added 2294300-23 APS2 bd to CAM (ECO 2064898).
- 11 Jun 28, 10 Added Rear Bridge Plate Cover (ECO 20635630. Added O2 Monitor and MRU (ECO 2076814) and Step-Down Transformer (ECO 2091712). Updated Host Monitor to 5245466. Added HP Z400 computer (ECO 2056587). Added HP8400 alternate replacement (ECO 2101300). Added MDP battery and charger FRU's (ECO 2097255). Added new 16 & 32 channel ICN's (ECO 2101939). Added Unified Surface Coil Phantom set (ECO 2081133).

CONTINUED ON NEXT PAGE.

REVISION HISTORY (Continued)

REV/DATE REASON FOR CHANGE

- | | | |
|----|------------|--|
| 12 | Jan 20, 11 | Removed BrainWave option key, 5110736 (ECO 2109423). Added 2156158–86 S/W collector for 16.0 M3 (DOC0853690). Added S/W collectors for 15.x. Added notes for Z400 computer and 4170 ICN are only compatible for use with 15M4B s/w or higher (MRIhc54664). |
| 13 | Aug 2, 11 | Changes made to reflect HD 16.0v2 update, ICN end-of-life, and IRF2b & RRF–DIF3 hardware upgrade configurations. |
| 14 | Dec 5 11 | Changes made to reflect HD23.x update (SPR HCSMD00101291). Revised Blower Box FRUs (PCN 207370, ECO 2123843). Updated PEN panel FRUs (ECR 2113479). Updated HFA Assembly FRU (ECR 2115417) |

[illegible]